

matrx[®]
SEATING SERIES

Flo-tech[™] **Solution Xtra**[™]



Cushion & Modular System



Combines excellent pressure-reducing properties with the promotion of active correction and management of posture.

The *Invacare* **Matrix Flo-tech Solution Xtra Cushion and Modular System** has been developed alongside clinicians to address complex seating needs.

Through using this **modular system**, therapeutic interventions can be more effectively addressed with regards to postural maintenance by enhancing function, maximizing comfort and preserving skin integrity.

Designed to promote active correction and management of posture

The modular components are manufactured from high-density CMHR foam, to achieve active correction and posture management. The wedges feature witness marks and Hook&Loop fixings, enabling them to be cut to meet the needs of individual users. The application of wedges beneath the cushion combined with the deployment of modular components on the top surface deliver a highly versatile range of options for therapists seeking to achieve specific postural management interventions.



Large dual fluid sac covering the entire top surface of the cushion

The contact area between the user and the foam cushion is protected, as the fluid re-distributes pressure away from bony prominences. The additional layer of fluid is incorporated under the ischial and sacral area to further enhance the pressure reducing qualities of this cushion, making this cushion suitable for Very High Risk users.

Cover Design

The cover is manufactured in an ingress-resistant material, which inhibits fluids from wetting the cover fabric. In addition, an innovative design has been incorporated into the cover, allowing it to expand and thereby increasing a user's immersion into the cushion once *Invacare* **Matrix Flo-tech Solution Xtra** components have been added. The cover features two zips, which have colour coded sliders: Blue to open the cover and green to expand. The cover can be opened either from the top or the bottom to adjust the components as necessary.

Features and options



Pelvic Obliquity pad

The pelvic obliquity pad can be inserted into either side of the ischial well under the gel pad to help correctable and fixed pelvic obliquities.



Hip Aligner

The Hip Aligner can be placed on either side on the rear edge of the cushion, secured by the Hook&Loop, and used either unilaterally or bilaterally to encourage mid-lining and hip alignment.



Thigh Aligner

The Thigh Aligner can be placed on either side at the front of the femoral shelf and used either unilaterally or bilaterally to reduce external rotation and abduction of the thighs.

Features and options



Pommel - Standard or Deep

The pommel can be placed on the front of the cushion to increase the existing pommel height to help femoral alignment, maintains hip joints and may assist with spasticity.



Rigidisers

The Rigidiser is to be placed under the whole cushion (including wedges). This maintains a firm base on slung canvases to maximise effect of component parts.



Fluid Pack

A Fluid Pack can be inserted into the ischial well where required, and can be doubled for specific placing. Specific placement will help prevent forward sliding, assist with sacral pressure and help mild pelvic obliquity.



Front Wedge

Front Wedge is to be placed on the underneath of the cushion base and will need to be trimmed to size. Helps assist with pelvic stability and accommodates posterior tilt whilst reducing anterior tilt.



Side (tapered) Wedges

Side wedges can be placed on any corner of the cushion base. Can help assist lateral hip rotation and thigh alignment. Can be used under one thigh for mild hip flexion discrepancies and for milder postural positioning.



Lateral Wedges

Lateral Wedges can be placed underneath the base of the cushion and must be trimmed to the size of the cushion. Helps facilitate lateral pelvic and hip stability and encourages more contoured positioning.








Ultra-fresh™ added to all cushions for antimicrobial and odour protection

Antimicrobial treatment that has been added to all **Flo-tech** foam to help make it resistant to the growth of odour causing bacteria. Helps provide longer lasting freshness.

Technical data

For more comprehensive pre-sales information about this product, including the product's user manual, please see your local Invacare website.

Size range	Height	Total product weight	Max. user weight ¹	Max. cleaning temperature
 44	 78 mm	 3.08 kg	 140 kg	 80°

**Matrix
Flo-tech Xtra**



Bariatric use


1. Guide maximum user weight limits are set assuming that the cushion selected is of appropriate width and depth for the intended user/s.

Manufactured to European Standards:

BS EN 1021-1 Cig, BS EN 1021-2 Match
and BS 7176 (Crib 5)

Colours

Black smooth

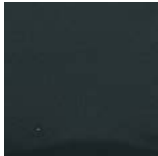

Black spacer
fabric



Please notice colors may vary slightly from those displayed above.

QR
Code

Invacare Ltd
Pencoed Technology Park - Pencoed
Bridgend CF35 5AQ - United Kingdom
Tel +44 1 656 776 222
E-mail: uk@invacare.com
www.invacare.co.uk

Invacare Ireland Ltd
Unit 5 - Seatown Business Campus
Seatown Road - SWORDS, County Dublin - Ireland
Tel. +353 1 8107084
Email: ireland@invacare.com
www.invacare.ie



© 2018 Invacare International GmbH. All rights reserved. All information quoted is believed to be correct at time of print. Invacare reserves the right to alter product specifications without prior consultation. Extra - 1628871 UK - 10/2018