



INVACARE® TOP END®
T-5 7000 SERIES
TENNIS WHEELCHAIR




The **Invacare® Top End® T-5 7000 Series Tennis Wheelchair** is faster, lighter, stronger and better than ever. In 2011 it hit the market with new material, enhanced headtube design and axle securing system, reducing overall weight by two pounds and improving performance. Each Top End T-5 7000 Series Tennis Chair is custom built to each individual athlete's specifications and is overseen by Paul Schulte who is the design engineer for Top End and the Captain of the USA National Men's Wheelchair Basketball Team.

Features

Custom-built fixed frame in 7005 aluminum
Adjustable height fixed swivel anti-tip
58mm fifth wheel
Fixed back with fixed or adjustable height
Adjustable tension seat and back upholstery
Fixed axle position with fixed wheelbase
Fixed camber: 15°, 18° or 20°
High performance Sun® wheels:
24", 25", 26" or 700C
Four bed offensive wing

Anodized pushrims
High pressure clincher tires
Precision fork caster system
72mm casters
Adjustable height tubular footrest
Powder coat paint
Rear frame protector
Axle securing system
250 lb. weight capacity

Model #: T57

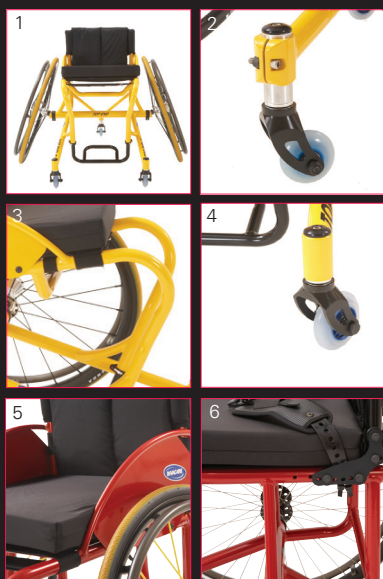

1. Custom front frame angle and tubular footrest
2. Adjustable height fixed swivel anti-tip
3. 7005 aluminum heat treating process
4. Precision fork caster system
5. Optional welded fixed tube sideguards
6. Adjustable center-of-gravity option



Invacare Corporation
www.invacare.com

USA

One Invacare Way
Elyria, Ohio
44035-4190
(800) 333-6900

Canada

570 Matheson Blvd. E.,
Unit 8
Mississauga, Ontario
L4Z 4G4 Canada
(800) 668-5324

Top End

www.topendwheelchair.com
Sports and Recreation
Products
4501 63rd Circle North
Pinellas Park, Florida
33781-5914
(727) 522-8677
(800) 532-8677
Fax: (727) 522-1007

©2014

Invacare Corporation.
All rights reserved.
Trademarks are identified
by the symbols ™ and ®.
All trademarks are owned
by or licensed to Invacare
Corporation unless
otherwise noted. Spinergy
is a registered trademark
of Spinergy Inc. Sun and
Equalizer are registered
trademarks of Sun
Components. Kenda is a
registered trademark of
Kenda Rubber Industrial
Co., Ltd. Specifications are
subject to change without
notification.
Form No. 12-508
Rev. 02/14
140186

Options

Spinergy® SLX wheels
Sun Fusion wheels
Adjustable Center-of-Gravity
Welded seat rigidizer bars with solid seat
Telescoping removable anti-tip bar
Click straps
Bolt-on fixed side guards
Fixed side guards with welded tube
Fixed aluminum or carbon fiber side guards
Titanium axles
Titanium pushrims

Coated pushrims
Footrest covers
Frame protectors
Custom strapping/toe clips
Alignment gauge
Adjustable angle footplate/toe clips
Recessed platform footplate with toe
stops/footrest strap
Kenda® tires
Close mount/Cut pushrim tab
O-ring positioning belt

