

Lift Out Chairs



Ever more rest and comfort with
the new lift out chairs range!

Invacare *Move Up*

NEW

The *Move Up* is the latest addition to the Invacare range of lift out chairs, integrating new features for comfort, mobility and aesthetics.

4 castors for greater comfort


Available as a 1 motor electric version, the *Move Up* has been completely redesigned. Its cosy and inviting style features a soft, padded back and padded backrest. For optimum comfort, side headrests have been integrated into the backrest. The shape of the armrests with rounded corners promotes good upper limb comfort and better blood circulation.

The *Move Up* stands out thanks to its innovative 4-wheel activation mechanism which allows the patient's chair to be moved around the home. Its subtle and elegant look fits easily into most environments.



Unique activation system of the 4 castors

In order to move the chair from one room to another, the *Move Up* is equipped with 4 retractable castors. Activation is effortless by means of a foot pedal on the back of the lift out chair. The user can remain comfortably seated in the chair during the transition. In order to regain its initial stability, the same pedal can be used to retract the castors and ensure safety during use.

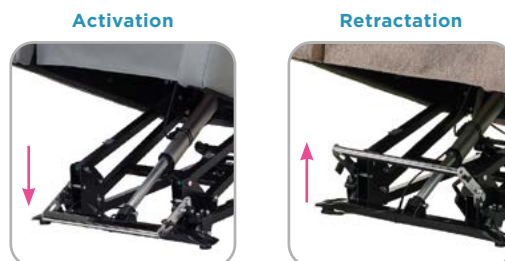

A modern and elegant design

Discover our new Heritage coating, designed with more durable and resistant materials. Its waterproof and water-repellent properties make the fabric more durable. Easy to maintain, stains are easily removed with a damp cloth.

This model is available in two colours:

- ▶ Heritage Camel for a cosy and refined look.
- ▶ PVC/Polyester grey for a chic and trendy design.



Features and Options ▼



Removable coating

The **Move Up** has a removable backrest for machine washing. Two zips at the top and bottom of the backrest are easily accessible to remove the backrest coating.



Side pocket

This chair has a storage space for the remote control so that it can be kept close by. This prevents it from being damaged.



Two button remote control

The **Move Up** comes with an ergonomic, easy-to-use 2 backlit buttons remote control.

Colours ▼



Camel Heritage



Grey PVC/
Polyester



Washing



Chloride



Drying



Ironing



Dry cleaning

Heritage coating maintenance ▼

References and Certifications ▼

1673370 **MOVE UP** 1 motor, Camel Heritage
1673371 **MOVE UP** 1 motor, Grey PVC/Polyester

Complies with fire standard 1021-1
Complies with fire standard 1021-1 and 1021-2

Invacare **Cosy Up**

The lift out chair **Cosy Up** provides maximum comfort thanks to its real duvet seat and the continuity between the seat and the footrest.

Its modern and uncluttered look is suitable for most people and allows it to easily find its place in any type of environment.



Depth-adjustable backrest as standard

The backrest of the **Cosy Up** can be easily adjusted in depth thanks to its two fixing brackets.

Thus, this chair adapts perfectly to all morphologies, from the smallest (with a depth of 450 mm) to the largest (with a depth of 540 mm).



Two seat widths

The **Cosy Up** lift out chair is delivered with a pair of padded armrests. The pair of thin armrests is available as an accessory for two seat widths. Both are attached to the rigid armrests by means of a velcro and zip system. They are therefore easily interchangeable.



To achieve the small seat width of 450 mm, the padded armrest must be installed.



To achieve the large seat width of 530 mm, thin armrest must be installed.

Features and Options ▼



Large side pockets

Lift out chair with two pockets, one on the left and one on the right, to hold magazines or small objects.



Removable coating for machine washing

The armrests, backrest and seat cover are attached to the base of the recliner by means of a zip and velcro system. All these elements are therefore easily removable, allowing machine washing at 30° or 40° depending on the fabric.



The remote control always at hand

Practice, the design of the armrests allows the user to place the chair's remote control (or TV remote control) under the armrest covers.

Options and Configuration ▼



1 motor version

Ergonomic remote control with 2 backlit buttons, easy to use.



2 motors version

Ergonomic remote control with 4 buttons, with integrated USB port as standard, easy to use.



Back-to-the-wall function

Thanks to the back-to-the-wall function, standard on the 1 motor version of the **Cosy Up**, the backrest of the lift out chair can be placed very close to a wall without this hindering the perfect use of all functions. Thanks to a design optimised for small spaces, only 9 cm are needed between the wall and the chair to allow the backrest to recline fully.

Features and Options ▼

EasyClean fabric : Water cleanable!

Antibacterial and water-repellent fabric, specially treated for cleaning with water. EasyClean fabrics allow you to remove everyday stains such as coffee, chocolate and wine with just a clean cloth and water. The EasyClean fabric prevents stains from penetrating the fabric and thus protects the chair permanently.

Practical and aesthetic, this stain-resistant fabric prevents mildew and is fire-resistant in accordance with the fire standards 1021-1 and 1021-2. Machine washable up to 40°C and suitable for tumble drying, the trendy and comfortable EasyClean cushions make the **Cosy Up** the ideal lift out chair for use in institutions, nursing homes and retirement homes!



Ocean EasyClean



Blue EasyClean



Taupe EasyClean



Grey EasyClean



Anthracite EasyClean

Cosy Up: mix, choose and change colours according to your mood!

With its 6 choices of cushion colours and 3 choices of bases, the **Cosy Up** is THE configurable lift out chair. Change the cover as you change your mood and transform your interior in just a few minutes. The **Cosy Up**'s configurability allows the distributor, in a reduced space, to present several colours by changing only the cushion (saving space and inventory!).

You can either order a complete chair or create the chair you want by selecting the base and by choosing the cushions you want: you can now change as many times as you like!

Bases ▾



Remote control base 1 motor PVC/ Polyester Taupe

Réf. 1650169



Remote control base 2 motor PVC/ Polyester Taupe

Réf. 1650205



Remote control base 1 motor PVC/ Polyester Dark Grey

Réf. 1670878

Cushions ▾



Ocean Easyclean

Fire standard 1021-1 and 1021-2

Réf. 1670881



Blue Easyclean

Fire standard 1021-1 and 1021-2

Réf. 1650211



Taupe Easyclean

Fire standard 1021-1 and 1021-2

Réf. 1650212



Anthracite Easyclean

Fire standard 1021-1 and 1021-2

Réf. 1670880



Light Grey Easyclean

Fire standard 1021-1 and 1021-2

Réf. 1670879



Taupe microfiber

Réf. 1650208



Taupe PVC/ Polyester

Réf. 1650209

Invacare® **Fino**, Invacare® **Faro**

The good compromise between comfort and simplicity

Functional and uncluttered look, the **Fino** and **Faro** chairs have a modern design that will easily fit into any environment. The standard headrest offers good support and completes the look of this lift out chair. The new Heritage Camel colour is perfect for modern environments.

The Fino chair has small seat dimensions: width 450 mm x depth 450 mm and has padded armrests.

Back-to-the-wall function: requires only 90 mm of clearance for the 1 motor version.



Faro
Camel Heritage



Fino
Camel Heritage

Several versions available

Both models are available in a 1 motor version that allows the footrest to be raised and the backrest to be tilted simultaneously. The **Faro** chair is also available in a 2 motors version (allows independent adjustment of the backrest and footrest) and in a manual recliner version (progressive tilting backwards is done by simply pressing on the armrests).



Faro
Black PVC/ Polyester

Features and Options ▼



Headrest for better support and maximum comfort. It is removable if needed.



Ergonomic and easy to use remote control with 2 backlit buttons (1 motor version).



Ergonomic and easy to use remote control with 4 buttons (2 motors version).



Side pocket for the remote control or small objects (versions with motor).

Colours ▼



Camel Heritage



Black PVC/ Polyester

References ▼

1644969

Fino 1 motor Camel Heritage

1608090

Faro 1 motor Camel Heritage

1608091

Faro 1 motor Black PVC/ Polyester

1608092

Faro 2 motors Camel Heritage

Invacare® Douro

The meeting of elegance and comfort!

The **Douro** lift out chair offers extra comfort thanks to the continuity between the seat and the footrest. This chair has lateral headrest and a removable fabric headrest. With its uncluttered style and small overall width, it will easily fit into the user's home.

Back-to-the-wall function: requires only 100 mm of setback for the 1 motor version.



Camel Heritage

Several versions available

The **Douro** chair is available in 1 motor version that allows the footrest to be raised and the backrest to be reclined simultaneously.

It is also available in a manual lift out chair version: the progressive tilting backwards is done by simply pressing on the armrests.



Camel Heritage



Chocolate PVC/ Polyester

Features and Options ▼



Lateral headrest for extra comfort.



Ergonomic and easy to use remote control with 2 backlit buttons (1 motor version).



Side pocket for the remote control or small objects (versions with motor).



Delivered with a second fabric headrest cover to allow rotation for washing.

Colours ▼



Camel Heritage



Chocolate PVC/ Polyester

References ▼

1608093

Douro Relax Chocolate PVC/ Polyester

1608094

Douro 1 motor Camel Heritage

1608097

Douro 1 motor Chocolate PVC/ Polyester

Invacare® Porto_{NG}

Sober and classic design

The design of the **Porto_{NG}** lift out chair in Heritage Grey allows it to be integrated in any kind of environment. It offers even more comfort thanks to the continuity between the seat and the footrest and the lateral supports at head level.

Back-to-the-wall function: only 50 mm of backward movement is required for the 1-motor version.

Several versions available

The **Porto_{NG}** is available in a single motor version that allows the footrest to be raised and the backrest to be reclined simultaneously.

The 2 motors version allows independent adjustment of the backrest and footrest.



Grey Heritage



Grey Heritage

Features and Options ▼




Lateral headrest to relieve head and neck.



Ergonomic and easy to use remote control with 2 backlit buttons (1 motor version).



Ergonomic and easy to use remote control with 4 buttons (2 motors version).

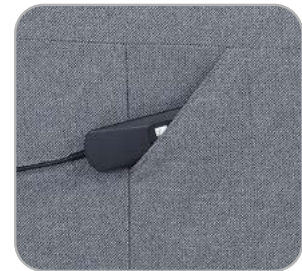

Side pocket for the remote control or small objects (versions with motor).

Colours ▼



Taupé PVC/
Polyester



Chocolate
PVC/
Polyester



Grey Heritage

References ▼

1608134

Porto_{NG} 1 motor Chocolate PVC/ Polyester

1608135

Porto_{NG} 2 motors Grey Heritage

1608136

Porto_{NG} 2 motors Taupé PVC/ Polyester

Characteristics table ▼



Colours ▼

Relax position ▼

Lift Out Chair position ▼

Rest position ▼

Continuity between seat and footrest ▼

Seat depth ▼

Seat height ▼

Seat width ▼

1 MOTOR								
1644969	Fino	Camel Heritage	◆	◆		450 mm	470 mm	450 mm
1608090	Faro	Camel Heritage	◆	◆		540 mm	470 mm	530 mm
1608091		Black PVC/ Polyester	◆	◆				
1608094	Douro	Camel Heritage	◆	◆	◆	540 mm	470 mm	530 mm
1608097		Chocolate PVC/ Polyester	◆	◆	◆			
1608134	Porto NG	Chocolate PVC/ Polyester	◆	◆	◆	540 mm	470 mm	530 mm
1644970		Taupe PVC/ Polyester & Taupe Microfiber	◆	◆	◆			
1650215		Taupe PVC/ Polyester	◆	◆	◆			
1650217	Cosy Up	Taupe PVC/ Polyester & Blue EasyClean	◆	◆	◆	450 or 540 mm	480 mm	450 or 530 mm
1650227		Taupe PVC/ Polyester & Taupe EasyClean	◆	◆	◆			
1660843		Dark Grey PVC/ Polyester & Light Grey EasyClean	◆	◆	◆			
1673370		Move Up	Camel Heritage	◆	◆	◆	480 mm	500 mm
1673371	Grey Heritage		◆	◆	◆			
2 MOTORS								
1608092	Faro	Camel Heritage	◆	◆	◆	540 mm	470 mm	530 mm
1608135	Porto NG	Grey Heritage	◆	◆	◆	540 mm	470 mm	530 mm
1608136		Taupe PVC/ Polyester	◆	◆	◆			
1644971		Taupe PVC/ Polyester & Taupe Microfiber	◆	◆	◆			
1650238	Cosy Up	Taupe PVC/ Polyester	◆	◆	◆	450 or 540 mm	480 mm	450 or 530 mm
1650240		Taupe PVC/ Polyester & Blue EasyClean	◆	◆	◆			
1650241		Taupe PVC/ Polyester & Taupe EasyClean	◆	◆	◆			
MANUAL								
1608093	Douro	Chocolate PVC/ Polyester	◆		◆	540 mm	470 mm	530 mm

The Cosy Up model is also configurable

The different positions ▼



RELAX POSITION:

Comfort position with seat and backrest tilt. This action is combined with the 1 motor chairs and can be desynchronised for a customised position with the 2 motors lift out chairs.



LIFT OUT CHAIR POSITION:

Position to sit and stand up effortlessly.



REST POSITION:

Reclining position for lying down and sleeping.



Max. backrest angle ▼



Total height ▼



Total width ▼



Total depth close ▼



Total depth open ▼



Necessary setback for max. tilt ▼



Floor/armrests height ▼



Max. product weight ▼



Max. user weight ▼



Foam density ▼

1 MOTOR

Fino	120°	1000 mm	730 mm	780 mm	1520 mm	90 mm*	670 mm	53 kg	125 kg	Seat 28kg/m ³ Backrest 18kg/m ³ Armrests 28kg/m ³
Faro	120°	1110 mm	770 mm	880 mm	1580 mm	90 mm*	600 mm	50 kg	125 kg	Seat 28kg/m ³ Backrest 18kg/m ³ Armrests 28kg/m ³
Douro	120°	1080 mm	740 mm	860 mm	1580 mm	100 mm*	640 mm	57 kg	125 kg	Seat 28kg/m ³ Backrest 28kg/m ³ Armrests 28kg/m ³
Porto NG	120°	1020 mm	790 mm	910 mm	1670 mm	50 mm*	600 mm	52 kg	125 kg	Seat 28kg/m ³ Backrest 18kg/m ³ Armrests 28kg/m ³
Cosy Up	120°	1050 mm	770(1)/ 800 mm	910 mm	1520 mm	90 mm*	640(1)/ 700 mm	57 kg	125 kg	Seat 28kg/m ³ Armrests 18kg/m ³
Move Up	120°	1110 mm	770 mm	970 mm	1580 mm	360 mm	500 mm	65 kg	120 kg	Backrest 28kg/m ³ Armrests 23kg/m ³

2 MOTORS

Faro	135°	1080 mm	770 mm	880 mm	1580 mm	480 mm	600 mm	50 kg	125 kg	Seat 28kg/m ³ Backrest 18kg/m ³ Armrests 28kg/m ³
Porto NG	135°	1020 mm	790 mm	860 mm	1670 mm	460 mm	600 mm	52 kg	125 kg	Seat 28kg/m ³ Backrest 18kg/m ³ Armrests 28kg/m ³
Cosy Up	135°	1050 mm	770(1)/ 800 mm	910 mm	1520 mm	90 mm*	640(1)/ 700 mm	57 kg	125 kg	Seat 28kg/m ³ Armrests 18kg/m ³

MANUAL

Douro	120°	1080 mm	740 mm	860 mm	1580 mm	200 mm	640 mm	57 kg	125 kg	Seat 28kg/m ³ Backrest 28kg/m ³ Armrests 28kg/m ³
--------------	------	---------	--------	--------	---------	--------	--------	-------	--------	--

* back-to-the-wall model

Accessories ▼

1663212	2 blue buttons backlit remote control with integrated USB port for lift chairs Cosy Up 1 motor	1677119	Pair of Cosy Up thin armrests, Blue EasyClean
1663213	4 blue button backlit remote control with integrated USB port for Cosy Up 2 motor lift chair	1677123	Pair of Cosy Up thin armrests, Ocean EasyClean
1677095	Pair of Cosy Up thin armrests, Taupe Microfiber	1677127	Pair of Cosy Up thin armrests, Taupe EasyClean
1677099	Pair of Cosy Up thin armrests, Taupe PVC/Polyester	1677131	Pair of Cosy Up thin armrests, Light Grey EasyClean
		1677173	Pair of Cosy Up thin armrests, Anthracite EasyClean

Spare parts ▼

For models with motor(s), common parts :

1577638	Transformer + mains cable
1577640	External battery + cable (allows models with motor(s) to be used without being connected to the mains, occasionally)

For models in 1 motor version :

SP1604803	Motor
SP1604806	Wired remote control, 2 backlit buttons

For models with 2 motors:

SP1604804	Motor for combined lift and footrest functions
SP1604805	Backrest motor
SP1604807	4 button wired remote control

Tel: +44 61 487 70 70
hqueurope@invacare.com
www.invacare.eu.com

Invacare International GmbH
Neuhofweg 51
4147 Aesch BL
Switzerland

