

Dacapo

preventative pressure care



Comfort and choice

Support surface solutions



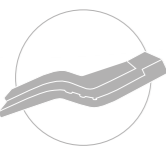
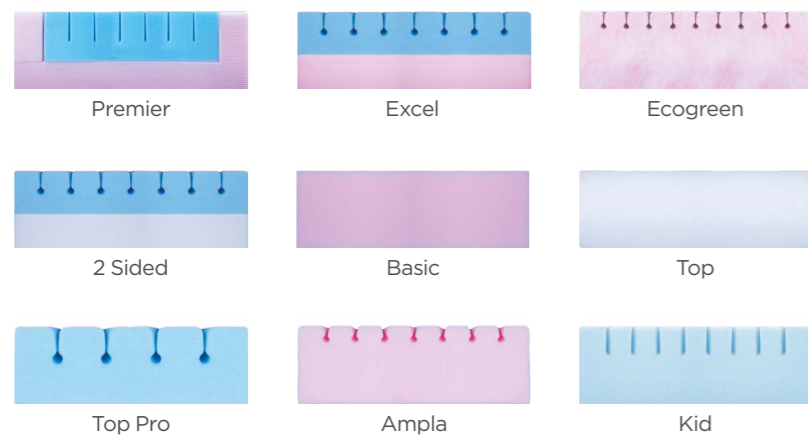
Yes, you can.®

Discover our Dacapo range

A comprehensive range of pressure redistributing foam mattresses


For over 20 years, Dacapo our trusted range of pressure redistributing mattresses have offered comfort and support to patients at risk of pressure injuries.

Expertly designed with clinical insight combined with our design innovation, Dacapo represents a complete range of mattress support solutions.



Pressure Redistribution

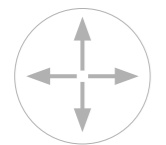All our clinical mattresses offer effective pressure redistribution in every patient position.



Shear & Friction Reduction

Our multi-directional stretch fabric helps to reduce shear and friction.



Microclimate Management

Well-aerated foam structures ensure optimal microclimate management.



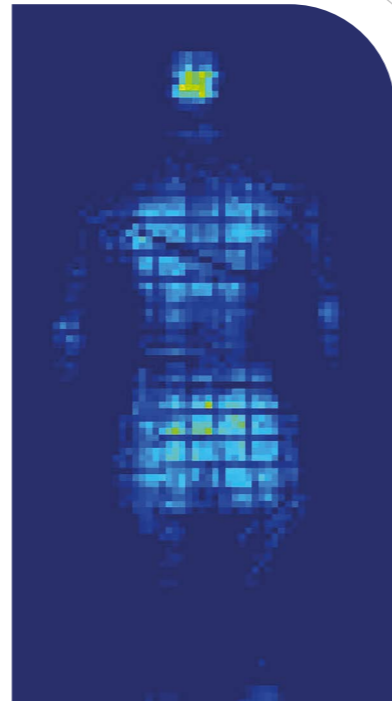
“ Developed to offer optimal cushioning and pressure redistribution ”

Pressure Redistribution

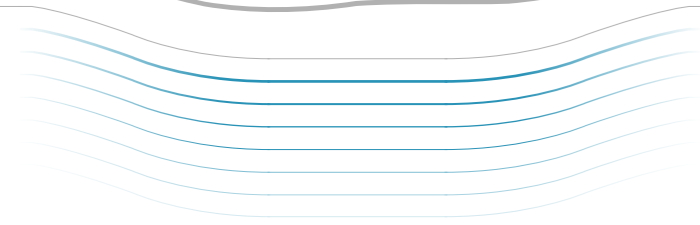
At Invacare, we use standardised pressure mapping technology to inform and better understand the effectiveness of pressure redistribution across our mattresses.

We are confident that each and every mattress design has undergone rigorous analysis and testing. In addition to pressure mapping our testing includes independent durability tests, cytotoxicity testing, fire testing and accelerated aging.




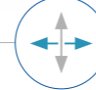


We ensure all our clinical mattresses are quality-assured and provide optimal comfort and support for the patient.



Illustrated above Dacapo Excel with patient in supine position.



Features you can expect from a Dacapo mattress:

- 
Quality-assured high specification foam
- 
Soft and supportive foam
- 
Durable and breathable material
- 
Multi-directional stretch fabric
- 
Machine-washable hygienic cover
- 
Delivered compressed

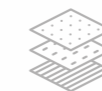
Dacapo Premier

Our flagship, high resilience U-core base mattress, is highly recognised for its positive affect on clinical outcomes.

The U-core design offers stability for patients mobilising in and out of bed and increases their feeling of being safe whilst transferring. A castellated foam insert provides effective pressure redistribution.

Dacapo Premier has outstanding responsive and recovery properties whereby the high resilience foam, recovers extremely well after each use, providing patients continued support and helps to reduce the risk of pressure injuries.


The castellations conform to every curve of the body, in every position and increases contact area with the skin, minimising surface pressure and shear.



> **Multi-layer** high resilience foam.


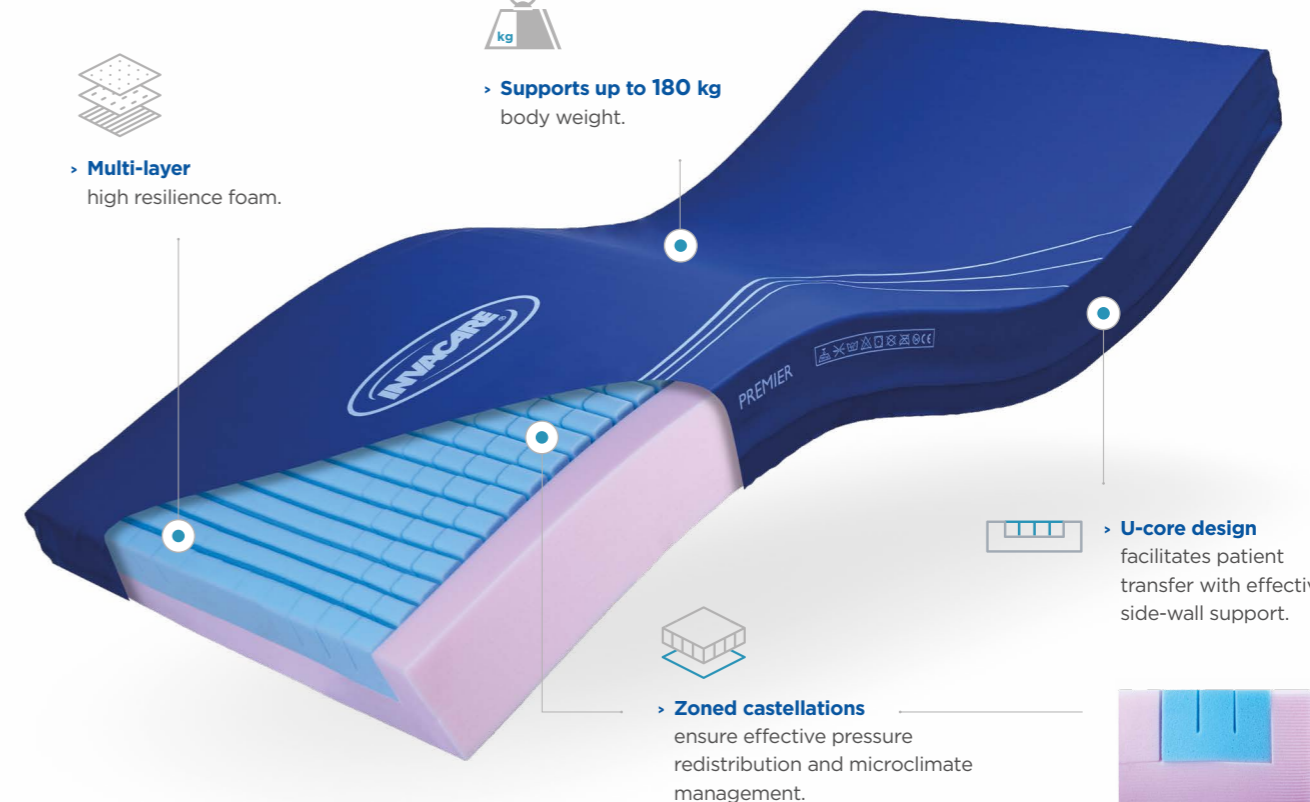
> **Supports up to 180 kg** body weight.



> **Zoned castellations** ensure effective pressure redistribution and microclimate management.



> **U-core design** facilitates patient transfer with effective side-wall support.

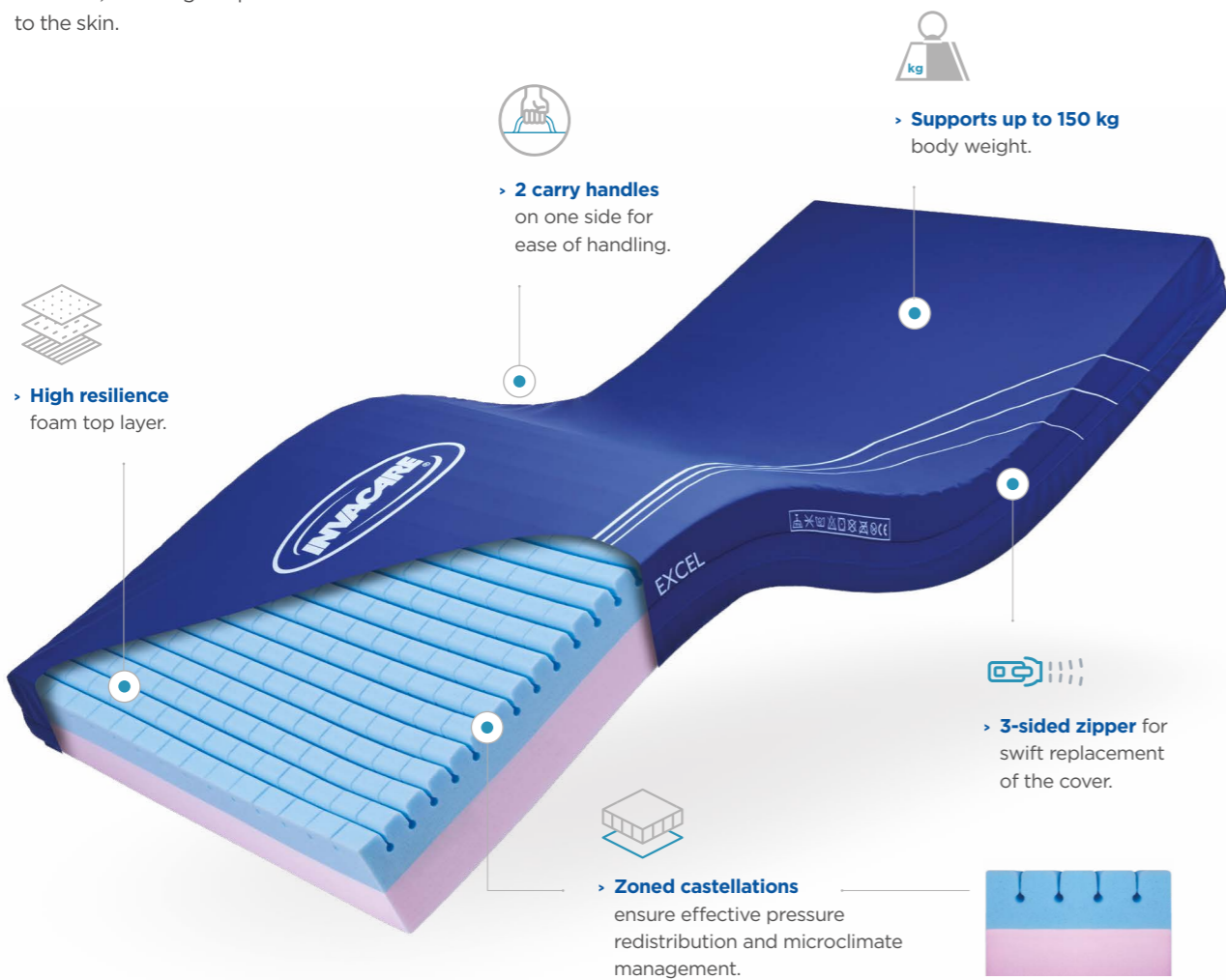


Dacapo Excel

Our latest pressure redistributing mattress has a multi-layer surface to support effective pressure care and microclimate management.

Specially designed to conform to the patient's body whilst ensuring optimal microclimate management. The Dacapo Excel provides outstanding responsive and recovery properties. Allowing the patient's body to immerse and envelop into the mattress, reducing the pressure to the skin.

The high quality foam ensures the mattress recovers extremely well after each use, providing optimal support for long-term care.



Features you can expect from a Dacapo mattress:

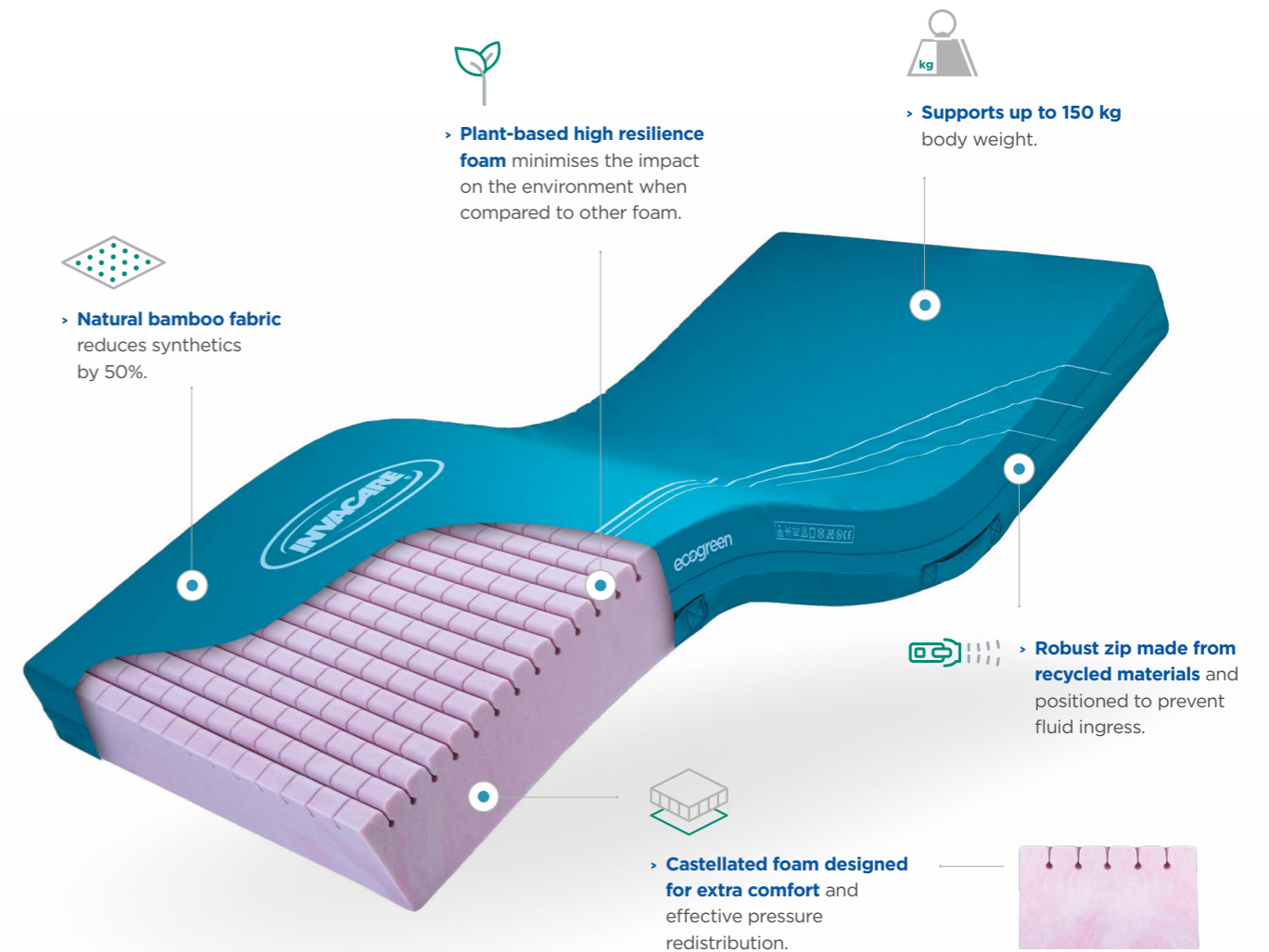
- Quality-assured high specification foam
- Soft and supportive foam
- Durable and breathable material
- Multi-directional stretch fabric
- Machine-washable hygienic cover
- Delivered compressed

Dacapo ecogreen

Responsibly designed, our new pressure redistributing mattress achieves clinical excellence in a sustainable way. Discover the future.

The design focuses on achieving excellent clinical outcomes required by healthcare professionals. Optimal pressure redistribution, reduction of shear and friction, mattress durability as well as patient comfort are all achieved whilst caring for our planet.

Dacapo Ecogreen is also fully compliant to the recommendations of the EPUAP, for high specification foam mattresses.



Dacapo 2Sided

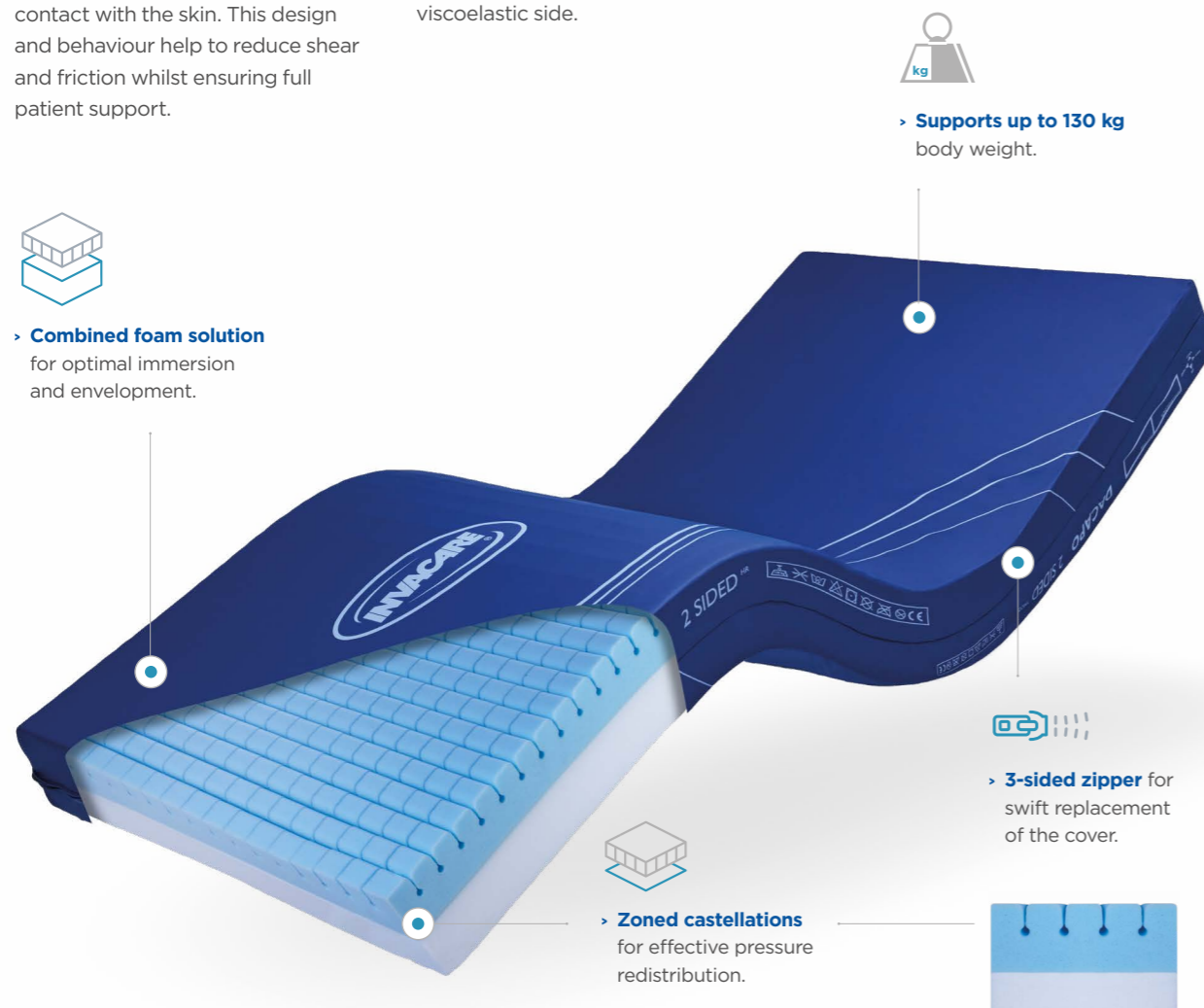
A combination of high resilience and visco elastic foams, that provide excellent immersion and envelopment properties and facilitates pressure redistribution.

Unique castellations enable the mattress to conform to the patient's body in every position with increased contact with the skin. This design and behaviour help to reduce shear and friction whilst ensuring full patient support.

Dacapo 2Sided is versatile, offering a choice of sides, between a high resilience foam side or a viscoelastic side.



> **Combined foam solution** for optimal immersion and envelopment.



Features you can expect from a Dacapo mattress:

- Quality-assured high specification foam
- Soft and supportive foam
- Durable and breathable material
- Multi-directional stretch fabric
- Machine-washable hygienic cover
- Delivered compressed

Dacapo Basic

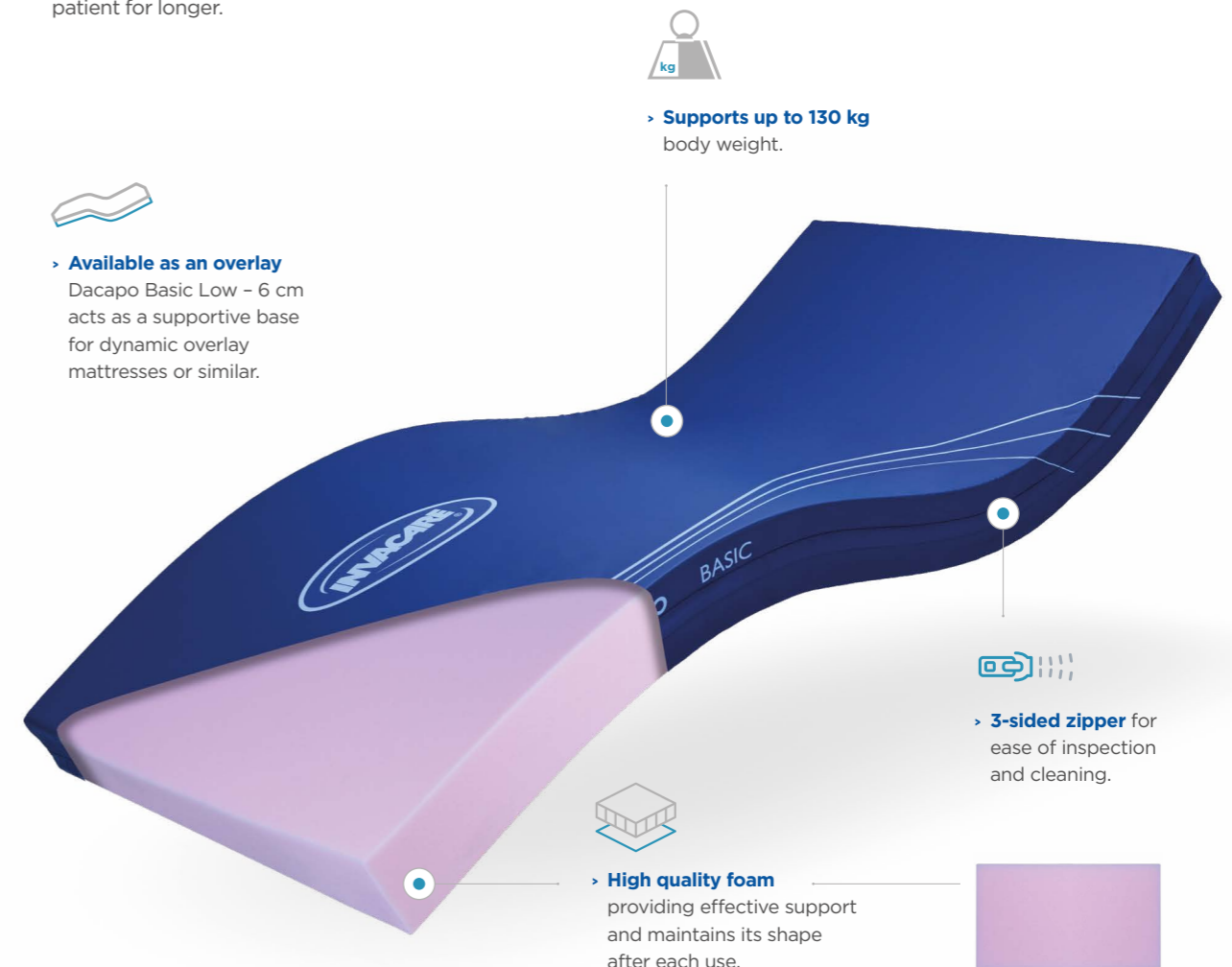
Providing patients, a comfortable and supportive surface whilst helping to reduce shear and friction.

Our high quality foam provides an effective support surface for the patient. Good responsive and recovery properties ensure the mattress restores and maintains its shape after every use, supporting the patient for longer.

Protected by our multi-stretch, vapour-permeable waterproof cover that is designed to be gentle on the skin and reduce shear and friction.



> **Available as an overlay**
Dacapo Basic Low - 6 cm acts as a supportive base for dynamic overlay mattresses or similar.




Dacapo Top

Our high quality, viscoelastic foam top layer offers increased support when used in combination with basic mattresses.

The foam adapts to the patient's temperature and body contour. With an open cell structure, the foam is well ventilated to reduce the build-up of heat and moisture, for effective microclimate management.

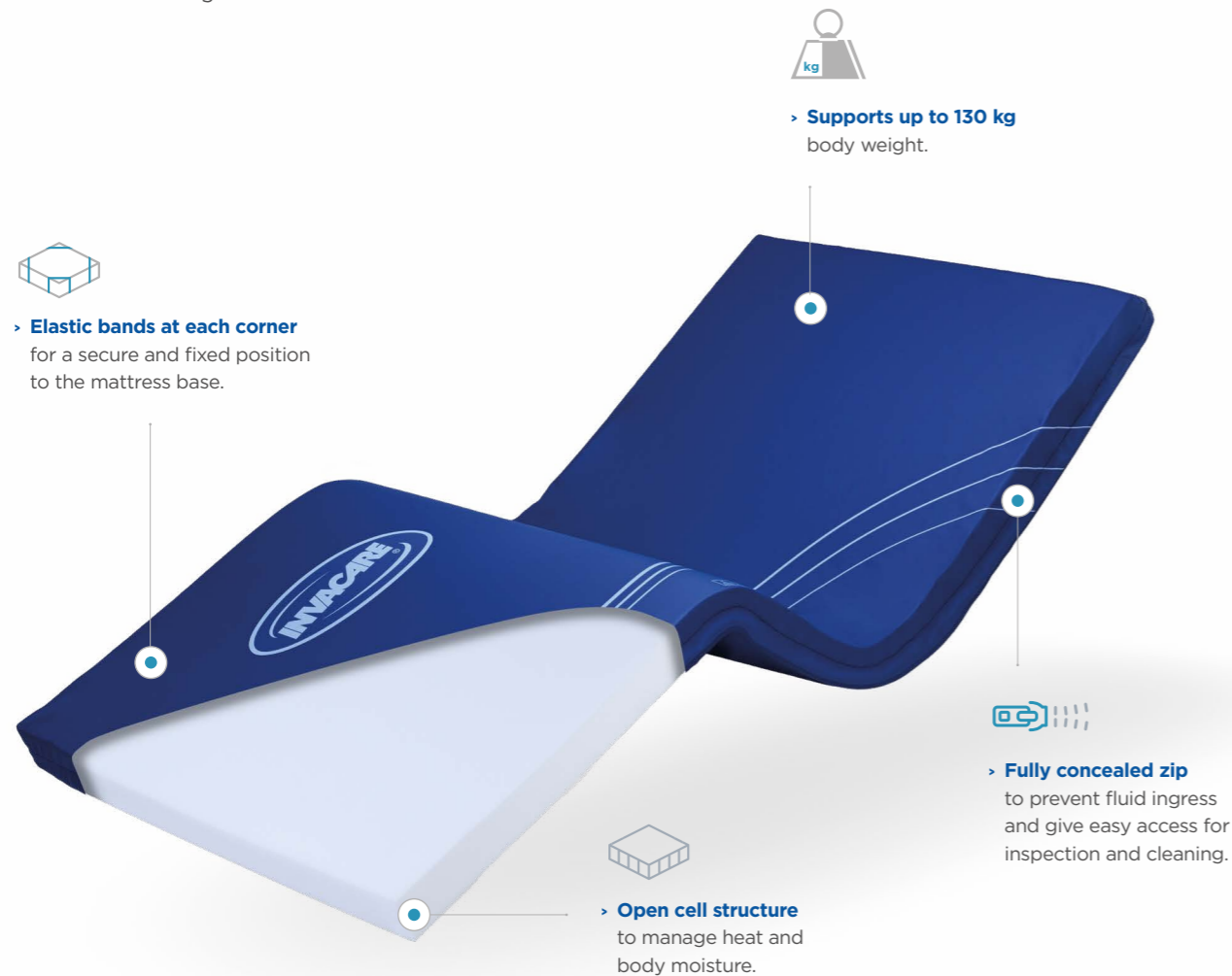Supported by our multi-directional stretch fabric and breathable cover. The Dacapo Top offers an advanced top layer for basic mattresses.




Dacapo Top Pro

A high-density foam mattress overlay, that helps to improve clinical outcomes when used with any other mattress base.

Zoned castellations in three anatomical areas, conform to the body and increase the contact area with the skin. Minimising shear and surface pressure with outstanding pressure redistribution.

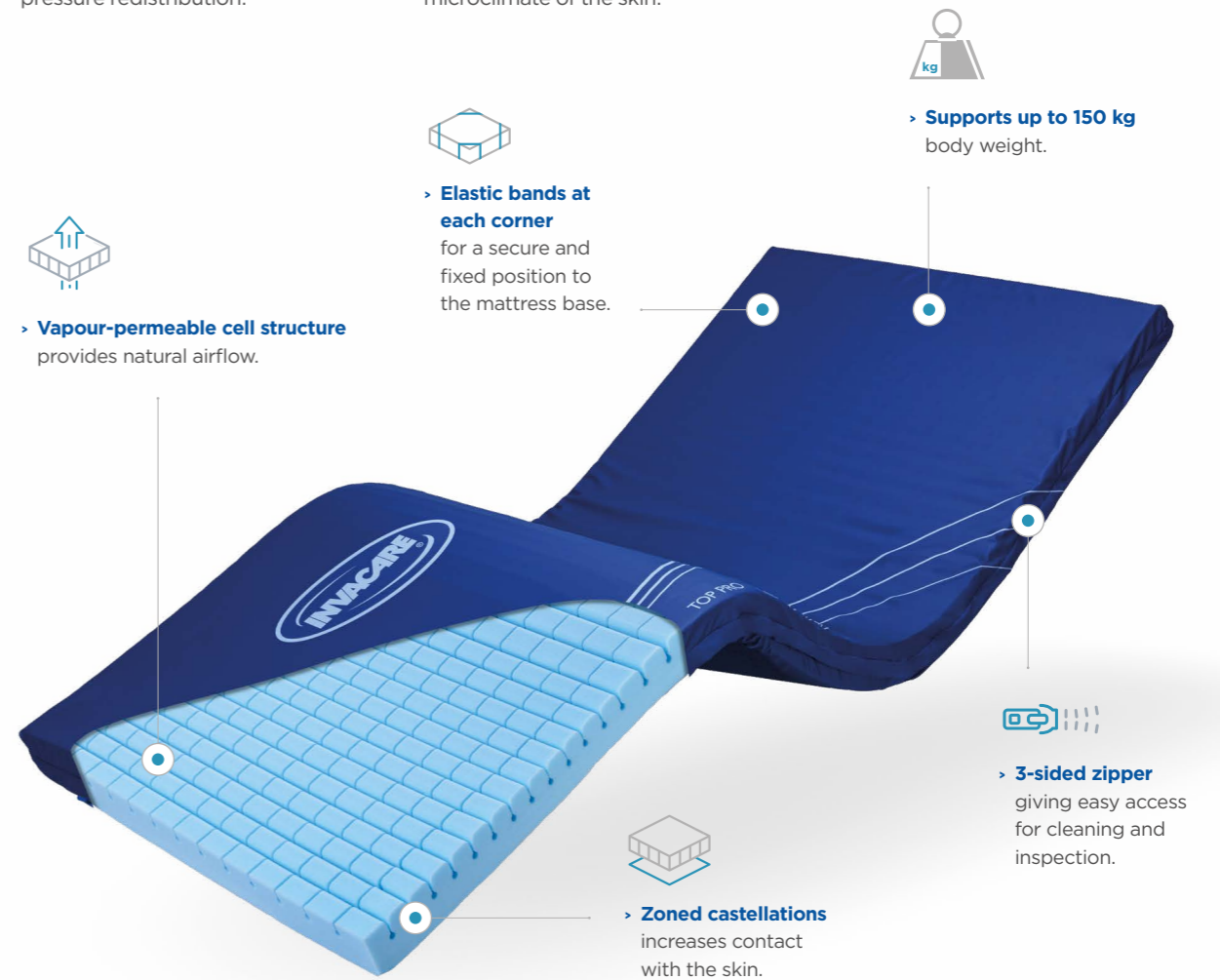Its vapour-permeable cell structure assures a well-aerated interior and in addition our breathable and waterproof cover, the Dacapo Top Pro positively contributes to the microclimate of the skin.



Features you can expect from a Dacapo mattress:

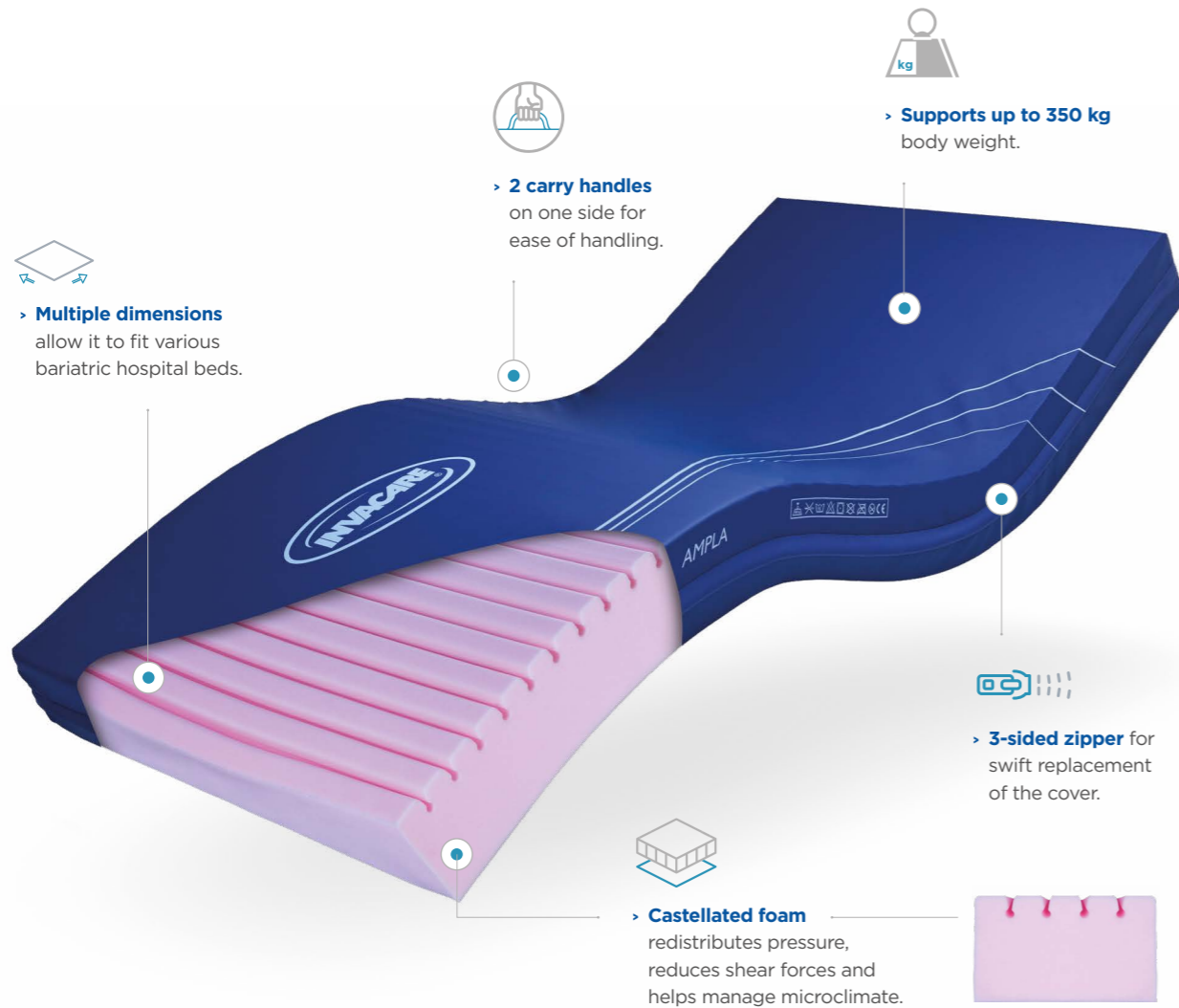
- Quality-assured high specification foam
- Soft and supportive foam
- Durable and breathable material
- Multi-directional stretch fabric
- Machine-washable hygienic cover
- Delivered compressed

Dacapo Ampla

Our new high specification mattress provides effective pressure redistribution and microclimate management for plus-size patients.

The design focuses on achieving excellent clinical outcomes required by healthcare professionals including optimal immersion and envelopment, pressure redistribution and the reduction of shear and friction forces.

The wider dimensions of the mattress offer dignified support and ensure maximum comfort for the patient in every position.



Features you can expect from a Dacapo mattress:

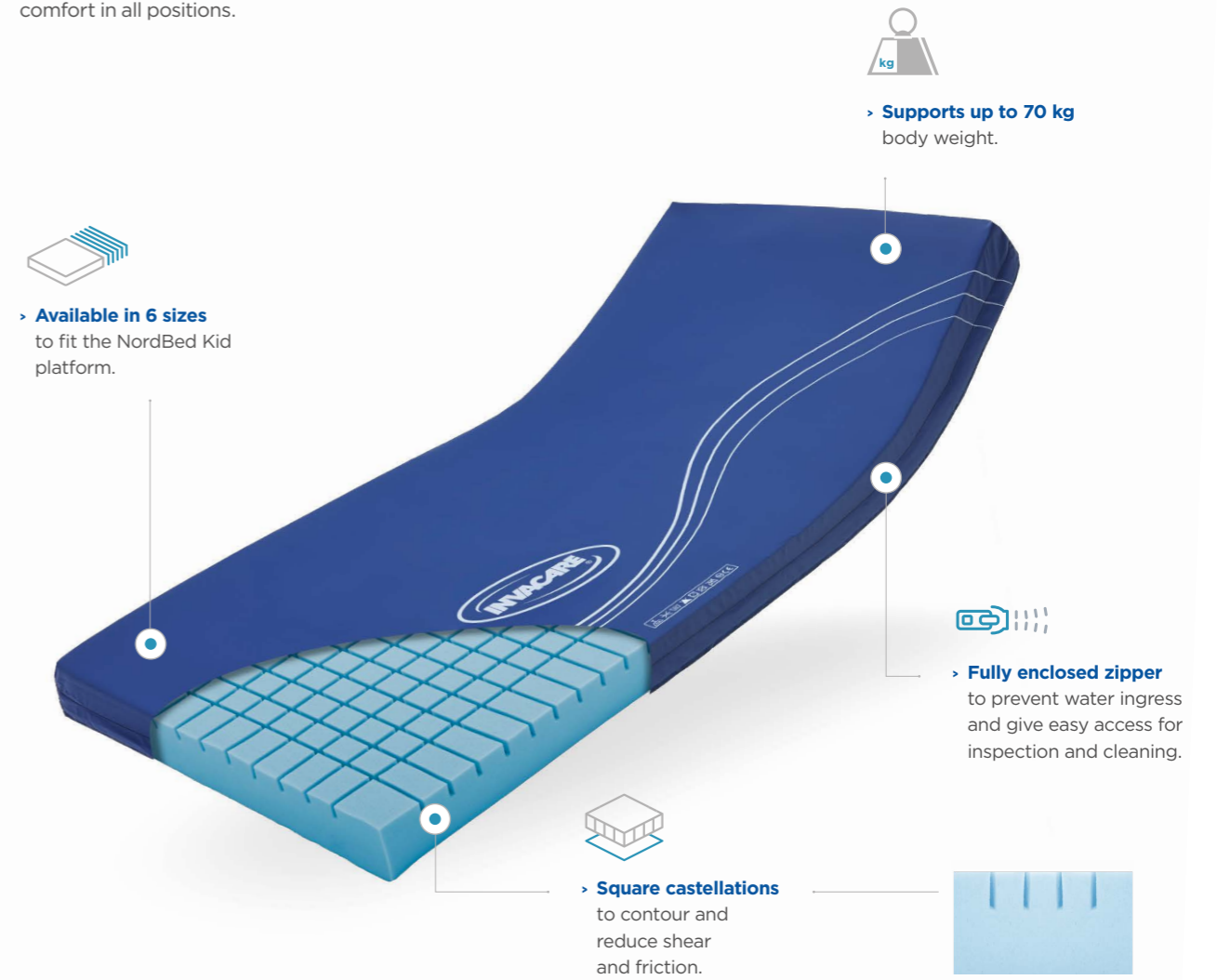
- Quality-assured high specification foam
- Soft and supportive foam
- Durable and breathable material
- Multi-directional stretch fabric
- Machine-washable hygienic cover
- Delivered compressed

Dacapo Kid

Our specially designed mattress for children gives the maximum support and comfort, in combination with our NordBed Kid medical bed.

A single-layer high resilience foam with square castellations ensure increased contact with the skin for effective pressure redistribution. The castellations conform to the contours of the child for increased support and comfort in all positions.











Our ingress-resistant cover is designed with the highest quality and hygiene standards, making the Dacapo Kid a comfortable support for any child.



Clinical Selection Guide

All Dacapo mattresses have been designed to manage varying degrees of pressure injury prevention. The chart below identifies those which are either low, medium or high when considering different areas of pressure injury prevention.

 Low
  Medium
  High

	Comfort & Prevention			Higher prevention			Specific Segments		
	Dacapo Top	Dacapo Top Pro	Dacapo Basic	Dacapo 2 Sided	Dacapo Excel	Dacapo ecogreen	Dacapo Premier	Dacapo Ampla	Dacapo Kid
Pressure									
Shear									
Friction									
Moisture									
Temperature									
Stability									
SWL (kg)	130	150	130	130	130 or 150	150	180	350	70

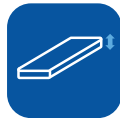
+ From a mattress overlay through to more advanced or segment specific mattresses, there is something for everyone. The Dacapo range can be relied upon for consistent features that actively contribute to pressure injury prevention, such as its durable and breathable material, multi-directional stretch fabric or positive microclimate management.

Quality-assured

We ensure all our clinical mattresses deliver optimal quality. In addition to pressure mapping, our testing includes, for example, independent durability testing, cytotoxicity testing, fire testing and accelerated aging.

Technical data ▼

For more information about this product, including the user manual, please visit our website: www.invacare.eu



Length Range ▼

Width Range ▼

Height Range ▼

Weight ▼

User weight ▼

Clean temperature ▼

Dacapo Premier

2000, 2100 mm

800 – 1200 mm

150 mm

14 kg

180 kg

95 °C

Dacapo Excel

2000, 2100, 2200 mm

800 – 1200 mm

120, 140 mm

10 – 15 kg

130, 150 kg

95 °C

Dacapo Ecogreen

2000 mm

800 - 1200 mm

150 mm

14 kg

150 kg

95 °C

Dacapo 2 Sided

2000, 2050 mm

800, 850 mm

150 mm

11.5 – 17 kg

130 kg

95 °C

Dacapo Basic

2000 mm

800, 850, 900 mm

120 mm

8.5 – 12 kg

130 kg

95 °C

Dacapo Top

1800 - 2200 mm

750 - 1200 mm

70 mm

5 – 7.5 kg

130 kg

95 °C

Dacapo Top Pro

2000, 2150, 2200 mm

800 – 1200 mm

70 mm

4.5 – 5 kg

150 kg

95 °C

Dacapo Ampla

2000 mm

850 – 1200 mm

150 mm

10 – 14 kg

350 kg

95 °C

Dacapo Kid

1358, 1800 mm

600 – 800 mm

100 mm

5 – 6.5 kg

70 kg

95 °C

Tel: +41 61 487 70 70
hqeuropa@invacare.com
www.invacare.eu.com

Invacare International GmbH
 Neuhofweg 51
 4147 Aesch BL
 Switzerland



Yes, you can.®