


Softform[®] Premier

Pressure-redistributing mattress

A high-quality mattress for comfort and pressure redistribution with a contoured insert pad.

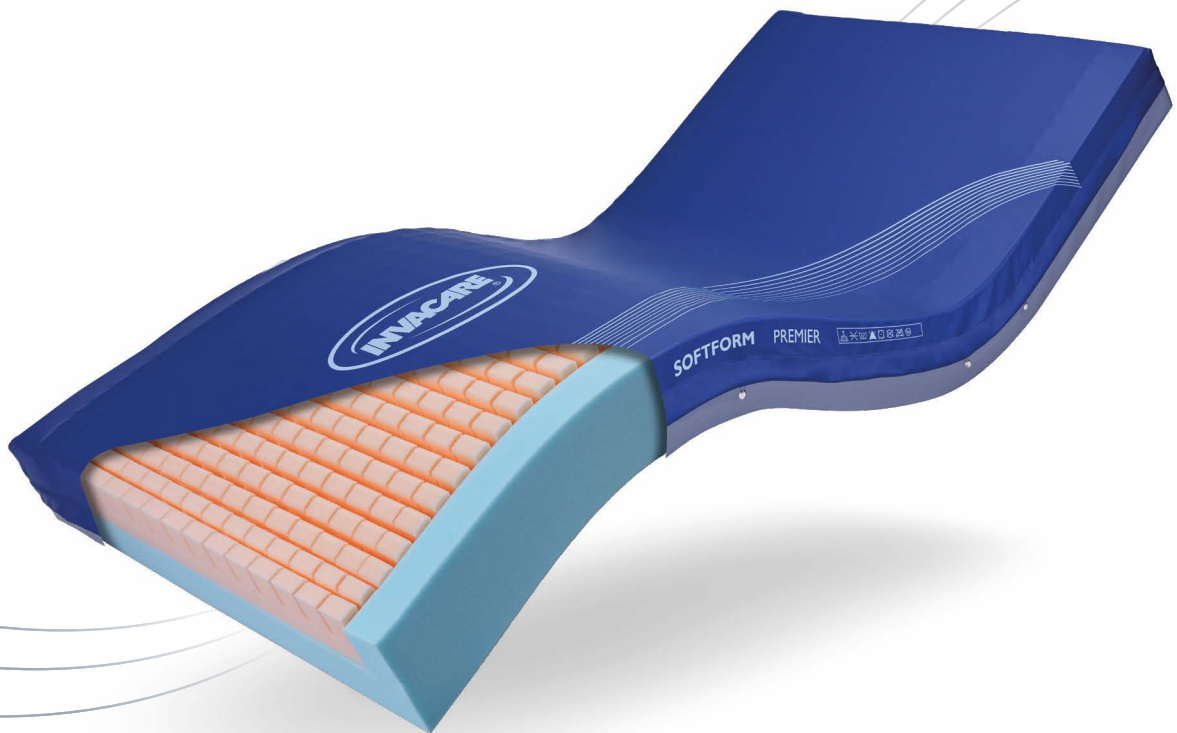

Stability and pressure redistribution

For personalised pressure injury management



Yes, you can.®

Effective pressure injury prevention



Pressure Redistribution

- › High quality foam offers a comfortable and effective pressure-redistributing support surface for patients at high to very high risk of developing pressure injuries.
- › Castellated foam helps redistribute weight over a greater area.
- › Non-glued construction allows good fit with profiled and electric beds.



Shear & Friction Reduction

- › Castellated foam conforms to movement and body position, helping to reduce shear and friction.
- › The top surface of the mattress tapers to relieve pressure on the heels, which are particularly vulnerable to pressure injuries.
- › Two-way stretch, water-resistant, polyurethane cover is durable and helps to reduce shear and friction forces.



Comfort

- › Castellations provide excellent comfort and air channels contributing to microclimate management.
- › Robust side walls and base facilitate stability during transfers and prevents the mattress from bottoming out.
- › Four-sided concealed welded zip reduces potential fluid ingress and allows easy detachment for inspection and laundering.

Technical data ▼

For more information about this product, including the user manual, please visit our website: www.invacare.eu.com

Sizes available ▼

		830 mm	850 mm	880 mm	900 mm	1050 mm	1200 mm	1370 mm
140 mm height	2000 mm	–	•	•	–	–	–	–
	1900 mm	•	•	•	•	–	•	•
150 mm height	1970 mm	•	•	•	•	•	•	•
	2000 mm	•	•	•	•	–	•	•
200 mm height	2000 mm	–	–	•	–	–	–	–



Max. User Weight ▼

247.6 kg



Max. Cleaning Temperature ▼

95 °C

Other sizes may be available upon request. Please ask your local Invacare representative.

Manufactured to BS 7177 (Ignition Source 5); Foam Interior: BS 5852:1982 (Ignition Source 5); Cover BS 7175 (Ignition Source 5).

Essential nursing care is pivotal in pressure injury management. This mattress may positively contribute to the outcome of a pressure injury prevention plan that includes education, clinical judgment and action planning based on vulnerability, which are fundamental factors in the prevention of pressure injuries. A range of assessment scales can be used as a formal method of assessing risk from pressure injury development, and should be used in conjunction with an informal assessment (informed nursing judgement) as part of a comprehensive and individualised plan.

Tel: +41 61 487 70 70
hqeuropa@invacare.com
www.invacare.eu.com

Invacare International GmbH
 Neuhofweg 51
 4147 Aesch BL
 Switzerland

