

Softform[®] Ampla

Supporting plus-size patients



High quality foam mattress

For effective pressure redistribution



Yes, you can.[®]

An overview of challenges for plus-size patients


Worldwide obesity continues to increase in both adults and children. In most EU states, over half of adults according to Eurostat* are declared overweight. Obesity causes a range of everyday challenges from self-care to limited mobility and leaves patients exposed to many adverse health problems.

The importance of access to adequately designed equipment that helps patients live a more active and independent lifestyle, has a direct impact on both mental health and clinical outcomes.

The need for comprehensive plus-size nursing care is increasingly required for global healthcare environments. Maintaining high levels of comfort and support to patients in a dignified way should be inclusive to all.

Our new Softform® Ampla pressure redistributing mattress is the latest addition to our range of products specifically designed for plus-size patients.

* Source: https://ec.europa.eu/eurostat/statistics-explained/index.php?title=Main_Page

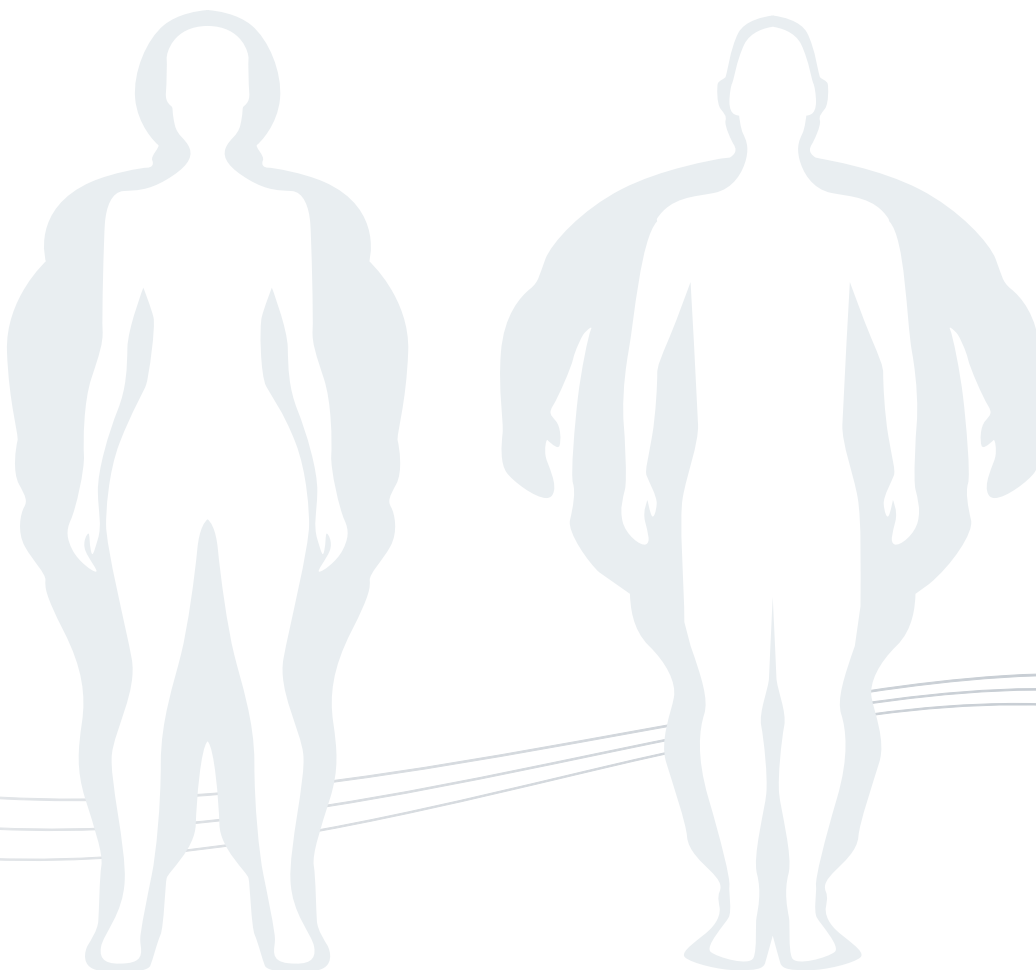


Designed for plus-sized patients

Our new Softform® Ampla high quality mattress provides effective pressure redistribution for plus-size patients.

The design focuses on achieving excellent clinical outcomes required by healthcare professionals including optimal immersion and envelopment, pressure redistribution and the reduction of shear and friction forces.

Maximum comfort in every position to support patient repositioning schedules and care giving routines such as hygiene and personal care.

1

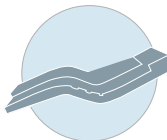

Clinically effective

Presents outstanding pressure redistribution properties for effective pressure injury prevention.

2



Demographically representative

Accommodates the growing need for larger, resilient and clinically effective products for plus-size patients.

3



Quality assured

Delivering optimal quality. Our testing includes pressure mapping, independent durability testing, cytotoxicity testing, fire testing and accelerated aging.



› **High quality medical grade foam** supports plus-size patients up to 340kg safely and effectively.




› **Breathable, water resistant cover** for effective temperature control and microclimate management.



› **Castellated foam** redistributes pressure, reduces shear forces and helps manage the skins' microclimate.




› **Zipper on 4 sides** for easy removal of cover for laundering.

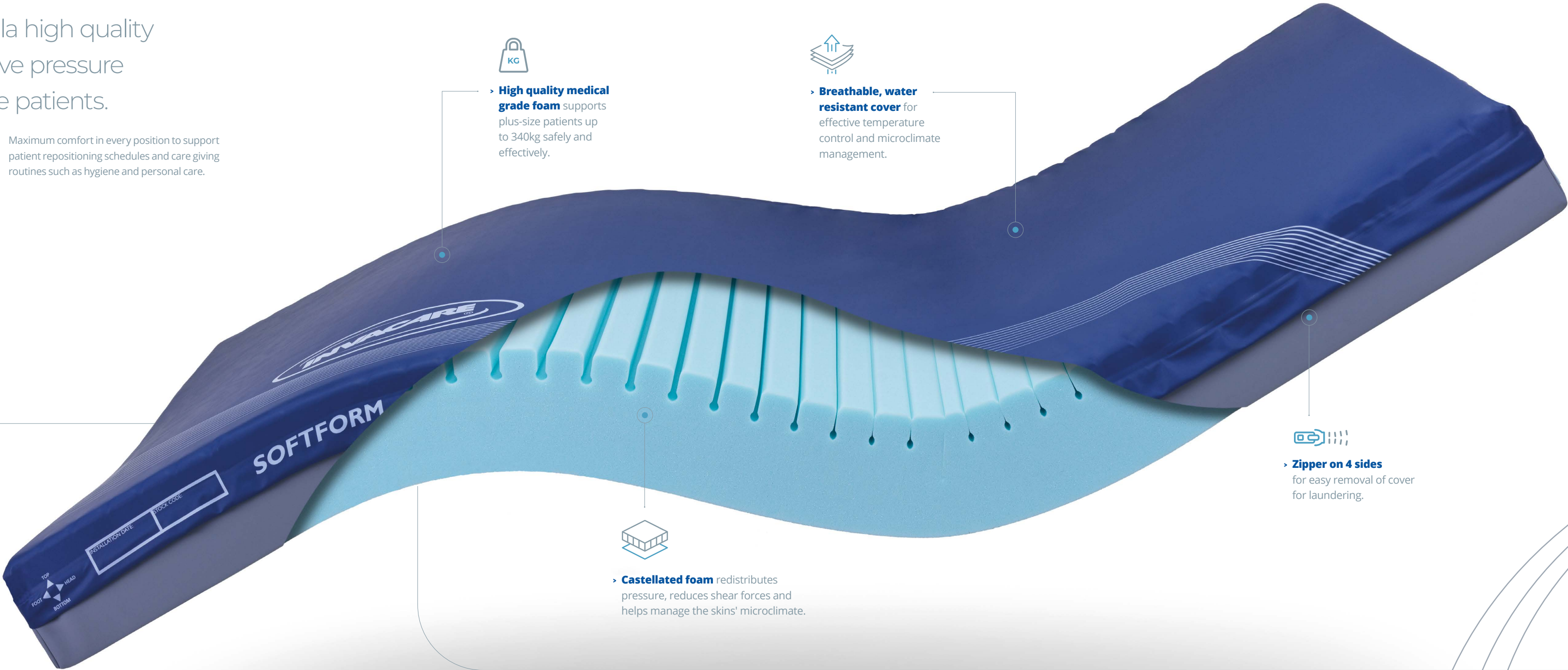

Pressure Redistribution

Softform® Ampla allows for optimal immersion and envelopment of the patient. This translates to a lower tissue interface pressure for vulnerable skin and facilitates a lower risk of developing pressure injuries.

The pressure mapping test results (illustrated right) demonstrates the excellent pressure redistribution of Softform® Ampla.

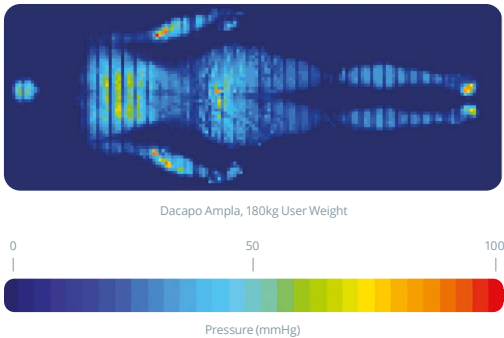

Technical data ▼

For more information about this product, including the user manual, please visit our website: www.invacare.eu.com



Sizes
Available ▼



Weight ▼



Max. User
Weight ▼



Max. Cleaning
Temperature ▼

150 mm height

2000 x 1200 mm

14 kg

340 kg

95 °C

Manufactured to BS 7177 (Ignition Source 5). Foam Interior: BS 5852:1982 (Ignition Source 5); Cover BS 7175 (Ignition Source 5).

Essential nursing care is pivotal in pressure injury management. This mattress may positively contribute to the outcome of a pressure injury prevention plan that includes education, clinical judgment and action planning based on vulnerability, which are fundamental factors in the prevention of pressure injuries. A range of assessment scales can be used as a formal method of assessing risk from pressure injury development, and should be used in conjunction with an informal assessment (informed nursing judgement) as part of a comprehensive and individualised plan.

Tel: +41 61 487 70 70
hqeuropa@invacare.com
www.invacare.eu.com

Invacare International GmbH
Neuhofweg 51
4147 Aesch BL
Switzerland



© 2024 Invacare International GmbH. All rights reserved. All information quoted is believed to be correct at time of print. Invacare reserves the right to alter product specifications without prior consultation. Softform Premier - EU - 11/2024

Yes, you can.®