

Perfecto2 V



Packed with features for today's
demanding lifestyles



Yes, you can.™

The new *Invacare Perfecto2 V* Oxygen Concentrator is an economical addition to the premier Invacare stationary concentrator line. This full 5-litre concentrator is small, quiet and lightweight and incorporates **SensO2** oxygen sensor technology. The **Perfecto2 V** is your cost-effective home oxygen concentrator, and comes with no scheduled maintenance for 3 years.

Quiet and Lightweight

Innovative techniques reduce the compressor and purge sounds to an almost undetectable level. By utilising the latest advancement in sieve material **Perfecto2 V** has reduced weight whilst maintaining the highest oxygen purity at all flow rates.

High-performance, safe and reliable

All of Invacare's rigorously tested safety features are incorporated into this latest evolution of our proven line of home oxygen concentrators. The durable product design and construction ensures years of reliable service with a minimum of fuss.



HomeFill compatible

Simplify your home oxygen business with **HomeFill** – the cylinder refilling system designed with patients in mind.

Features and Options



No scheduled maintenance for 3 years

Provider maintenance is effortless and fast.

Simple for the user

All device controls are on the front panel and within easy reach of the user.












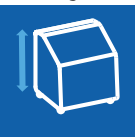


Paediatric flowmeter

50 to 750cc stand alone device for accurate low flow delivery.




Technical data

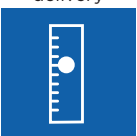For more comprehensive pre-sales information about this product, including the product's user manual, please refer to your local Invacare website.

	Length	Width	Height	Weight	Charger	Power consumption	Sound level	O2 purity
<i>Perfecto₂ V</i>	 305 mm	 380 mm	 585 mm	 18.2 kg	 230 V, 50 Hz (±10%)	 320 W	 < 40dB(A)	 87 - 96%

	Outlet pressure	Pulse dose delivery	Altitude	ISO Compliant
<i>Perfecto₂ V</i>	 34.5 kPa ± 10%	 0.5 - 5.0 LPM	 6000ft (1828 m)	 ISO80601-2-69:2014

SensO₂-System

O ₂ -concentration	
	higher than 85% - green light - Normal
	between 73% and 85% - yellow light - Check filter and flow rate
	lower than 73% - red light & audible alarm - Call service technician

Service interval

No provider routine preventative maintenance for 3 years, replace filters as required.



Scan here to access the website oxygen-therapy.net

Alarms & safety

- Battery-free power-loss alarm
- Resettable circuit breaker
- SensO2 oxygen purity monitoring
- Compressor High Pressure relief valve
- Compressor High Temperature shut down
- High & Low internal pressure alarms
- Self-diagnostic codes with memory for technicians
- Low flow alarm

Tel: +41 61 487 70 70
hqeurope@invacare.com
www.invacare.eu.com

Invacare International GmbH
Benkenstrasse 260
4108 Witterswil
Switzerland



Yes, you can.