

Check O₂ Plus



Oxygen Analyzer



The *Invacare Check O₂ Plus Oxygen Analyzer* is a single economical instrument that reports all the operating parameters necessary for performing oxygen concentrator checks. It is an efficient and cost-effective analyzer designed with technicians in mind. It checks oxygen concentration, flow rate and patient outlet pressure on any oxygen concentrator. What's more, it requires no routine calibration, and there's no fuel cell to replace. The *Invacare Check O₂ Plus Oxygen Analyzer* is three analyzers rolled into one.






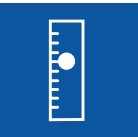

Features

- Measures three parameters, eliminating the need for multiple instruments: oxygen purity, flow rate and patient outlet pressure
- Accurately measures oxygen purity of the *Invacare XPO2* Portable Concentrator
- Features simple two-button operation
- Two AA batteries provide over 16,000 readings
- Shuts off automatically to extend battery life
- Requires no in-field calibration
- Self-diagnostics
- Includes connecting tube and protective boot



Technical data

For more comprehensive pre-sales information about this product, including the product's user manual, please see your local Invacare website.

	Weight	Quantity Per Case	Oxygen Flow	Pressure Range	O ₂ purity	Operating Temperature	Limited Warranty
Check O₂	 Weight: 185 g	 1	 0-10 l/min ± 0,2 l/min	 0,5-50 PSI ± 0,5%	 21%-95,6% ±1,5%	 15°-40°C	 2 Years

QR
Code

Tel: +41 61 487 70 70
hqeurope@invacare.com
www.invacare.eu.com

Invacare International GmbH
Benkenstrasse 260
4108 Witterswil
Switzerland



Yes, you can.®

