

**matrix**<sup>®</sup>  
SEATING SERIES

# Vi Cushion



Yes, you can.



The *Invacare* **Matrix Vi Cushion** is designed to provide skin protection and positioning for individuals at moderate to high risk of skin breakdown. The unique three part foam construction targets skin protection where it is needed most, and features excellent stability, support and comfort.



① **Outer cover**

Moisture resistant and breathable. Includes slip resistant base, hook and loop fasteners, lifting strap.

② **Inner cover**

Moisture resistant. Easily wipes clean. Zippered cover allows access for cushion modifications.

③ **Top Layer**

Soft (ILD) High Resilient (HR) foam for optimal comfort and support. Pre-ischial contour provides stability to pelvis. Abductor/adductor contour promotes midline leg positioning. Made with Ultra-Fresh™ for antibacterial and odor protection.


④ **Visco foam insert**

Molded visco-elastic foam ischial relief insert. Anatomically targeted to ischial/sacral area. Provides immersion to sensitive bony prominences.

⑤ **Bottom Layer**

Firm (ILD) High Resilient (HR) foam for support and durability. Trochanteric shelf provides stability to hips. Waterfall front edge. Rear radius for sacral support.

Features and Options

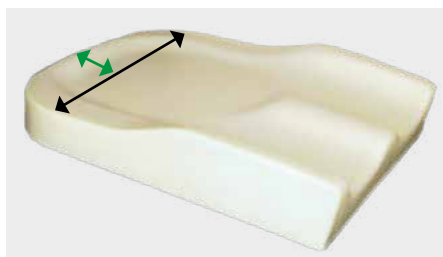

**Rear cushion radius:** +3.8 cm extra depth is designed to fill the gap between the seat and backrest/upholstery. **1.25 cm decrease in total width** allows cushion to fit both stated size and seats that are 2.5 cm narrower.



**Waterfall front end**

A diagonal cut along the front accommodates tight hamstrings, facilitates foot propulsion and standing transfers.

Technical data

For more comprehensive pre-sales information about this product, including the product's user manual, please see your local Invacare website.

|                        | Width                         | Depth      | Height at adductor/abductor | Height at leg trough | Total product weight | Max. user weight | Limited warranty |
|------------------------|-------------------------------|------------|-----------------------------|----------------------|----------------------|------------------|------------------|
| KA-Vi                  | 31 - 36 cm                    | 31 - 36 cm | 7,3 / 7 cm                  | 5 cm                 | 0,7 kg and up        | 68 kg            | 2 years          |
| MA(WWDD)-Vi (Standard) | 36, 41, 43, 46, 51, 53, 56 cm | 38 - 51 cm | 9,5 / 8,25 cm               | 6,35 cm              | 0,9 kg and up        | 136 kg           | 2 years          |

Colours and options

Black Mesthex    Kid\*ab\*ra fabric (Standard)    Pelvic Obliquity Kit

Easily installed on all Invacare **Matrix Vi Kid\*ab\*ra-Vi** cushions. Accommodates up to 25 mm obliquities

Tel: +44 61 487 70 70  
 hqeuropa@invacare.com  
 www.invacare.eu.com

Invacare International GmbH  
 Neuhofweg 51  
 4147 Aesch BL  
 Switzerland



Yes, you can.®