

*Invacare*<sup>®</sup>

## ***Propad™ Original Cushion***



Cut foam with castellations for comfort and support

Trusted Solutions,  
Passionate People



The *Invacare Propad Original Cushion* has been designed to facilitate good pressure redistribution when used with both wheelchairs and standard armchairs. This cushion is constructed to ensure improved support and prolonged longevity.

## Features and Options

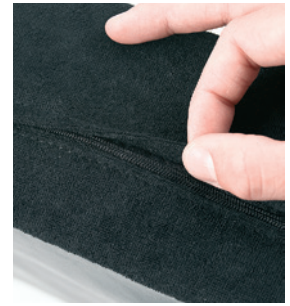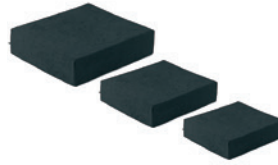

### Cover helps prevent a build up of shear and friction forces

Cover is two-way stretch, water-resistant and vapour-permeable. Cover features an anti-slip base.

### Available in a range of sizes

The *Invacare Propad Original Cushion* is available in 3 standard sizes, 410 x 410, 430 x 430, and 460 x 410 mm. (See prescription chart for details).

### Full length flap conceals zip

The full length flap concealing the zip reduces the potential for fluid ingress into the cushion and protects the foam core. The cover can be removed and laundered.



### Designed to deliver enhanced pressure care

The *Invacare Propad Original Cushion* is made from a layer of castellated foam, each segment delivers enhanced support and pressure area care to the individual. The "keyhole" cut at the base of each castellation ensures good air circulation, reducing potential for skin maceration.

### Contoured seating surface

Offers improved stability and postural support; castellations reduce shear and friction forces on tissue.

## Technical data


For more comprehensive pre-sales information about this product, including the product's user manual, please see the [Invacare website](#).

	Size range	Height	Total product weight	Max. user weight	Max. cleaning temperature	Warranty	
							Manufactured to European Standards:
<b>Propad Original Cushion</b>	406 x 406 mm 432 x 432 mm 406 x 457 mm	102 mm	1 kg	108 kg (17 stone)	80°	1 year	Interior Foam: BS EN 597-1 Cig, BS EN597-2 Match, BS 5852 (Crib 5) Cover: BS 7175 (Crib 5)

## Colours

Blue PU

Black Towelling

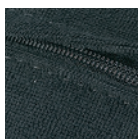

Tel: +44 1 656 776 222  
Fax: +44 1656 776 220  
E-mail: [uk@invacare.com](mailto:uk@invacare.com)  
Sales Order E-mail:  
[ordersuk@invacare.com](mailto:ordersuk@invacare.com)  
[www.invacare.co.uk](http://www.invacare.co.uk)

Invacare Limited  
Pencoed Technology Park  
Pencoed - Bridgend CF35  
5AQ - United Kingdom

Tel. +353 1 8107084  
Fax +353 1 8107085  
Email: [ireland@invacare.com](mailto:ireland@invacare.com)  
[www.invacare.ie](http://www.invacare.ie)

Invacare Ireland Ltd  
Unit 5 - Seatown  
Business Campus  
Seatown Road  
SWORDS - County  
Dublin - Ireland



**Yes, you can.®**

