

Propad™ Overlay



High specification foam
pressure redistributing overlay




Yes, you can.

The **Propad Overlay** is designed to sit on top of an existing conventional mattress, providing a high level of patient comfort with excellent pressure redistribution abilities. The overlay is lightweight, available in a range of sizes, making it ideal for the homecare environment.


Effective Pressure Ulcer Prevention

- ▶ The foam core has independent, geometrically-cut surface cells which significantly reduce shear and friction.

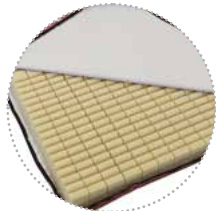

Cover helps to prevent a build up of shear and friction forces

- ▶ Cover is multi stretch, water-resistant and vapour-permeable.

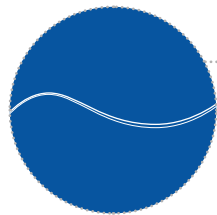

Design

- ▶ The **ProPad Overlay** has an innovative and modern cover design with the simple wave giving the mattress a fresh and current feel.



Pressure redistribution

- ▶ With larger, more supportive cells at the centre and smaller cells at each end, enabling the mattress overlay to achieve even pressure redistribution by conforming to the shape and form of the body and cradle bony prominences.

Cover has an anti-slip base

- ▶ Ensures safe use of the **Propad Overlay** and helps to prevent slippage.



Technical data ▼

For more comprehensive pre-sales information about this product, including the product's user manual, please see your local Invacare website.



Std sizes ▼



Height ▼



Total product weight ▼



Max. user weight ▼



Max. cleaning temperature ▼

	Std sizes ▼	Height ▼	Total product weight ▼	Max. user weight ▼	Max. cleaning temperature ▼
Single	1875 x 876 mm	85 mm	5.5 kg	108 kg	95°
Double	1875 x 1372 mm	85 mm	5.5 kg	108 kg	95°
King Size	1980 x 1524 mm	85 mm	5.5 kg	108 kg	95°

Additional sizes available on request

Manufactured to European Standards:
BS EN 597-1 Cig, BS EN 597-2 Match and BS 7175 (Crib 5)

** Essential nursing care is pivotal in pressure ulcer prevention. This product will positively contribute to the outcome of a pressure ulcer prevention care plan. Education, clinical judgement and action based planning based on vulnerability are fundamental factors in the prevention of pressure ulcers. A range of assessment scales can be used as a formal method of assessing risk from pressure ulcer development, and should be used in conjunction with an informal assessment (informed nursing judgement). Informal assessment is considered to be of greater importance and clinical value.*

Tel: +44 61 487 70 70
hqueurope@invacare.com
www.invacare.eu.com

Invacare International GmbH
Neuhofweg 51
4147 Aesch BL
Switzerland



Yes, you can.®

