



For safe and comfortable, single carer assisted transfers

The *Invacare Reliant 350 Stand Assist* is a compact and maneuverable mobile lifter that provides safe and comfortable assisted transfers for those with limited mobility or rehabilitation needs. The *Invacare Reliant 350* can be used to transfer a client that cannot transfer unaided, but has some level of weight bearing ability, from one seated position to another. The robust design with a weight limit of 158 kg allows for additional support whilst transferring, preventing unnecessary strain on the carer, making it a suitable option for transfers carried out by a **single carer***.

- **5 year warranty:** Frame, Actuator, Control box
- **2 year warranty:** Hand control




Adaptable and easy to use

Reliant 350 features a multi-adjustable lift boom and a selection of sling mounting points to provide an adjustable and ergonomic lift angle. Large rear and small front castors help manoeuvre the **Reliant 350** into small or confined spaces and the base has a lever adjustment mechanism for added stability. In addition, it has a generous non slip footplate providing a stable surface for the feet, as well as a height-adjustable padded knee guard.

Safe and reliable


The **Reliant 350** provides a range of safety features including an emergency stop button, emergency mechanical descent system, as well as an anti-entrapment device which stops the boom from lowering if it meets resistance. The battery recharge indicator is visual and audible. For transport and storage, the **Reliant 350** can be easily dismantled.

Features and Options

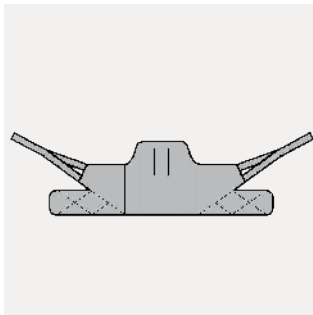

Stand Assist Sling

Designed to be used with a **Stand Assist** lifter, this sling features the latest innovative fabric for improved comfort and handling. Available in sizes; S, M, L & XL.

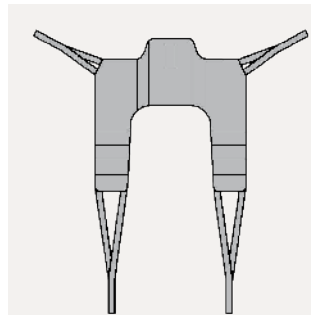

Transfer Stand Assist Sling

Designed for **Stand Assist** lifters, for those who require additional leg support, this sling provides a supportive and comfortable transfer. Available in sizes; S, M, L & XL.



Independent Charger

An off-board charger and spare battery allow for 24 hour use. Please note: The battery cannot be charged on the hoist.

Technical data

167 cm	102 cm	90 cm	94 cm	66 cm
11 cm	100-200 lifts per full charge	49 kg	Max. 158 kg	

*A risk assessment should always be carried prior to using any moving and handling product.

Frame colour



Beige

Warranty

Frame, control box & actuator
5 years
Hand control
2 years
Battery
1 year

Invacare Limited
Pencoed Technology Park -
Pencoed - Bridgend CF35 5AQ -
United Kingdom
Tel: +44 1 656 776 222
Fax: +44 1656 776 220
E-mail: uk@invacare.com
Sales Order E-mail:
ordersuk@invacare.com
www.invacare.co.uk
Invacare Ireland Ltd
Unit 5 - Seatown Business
Campus
Seatown Road - SWORDS -
County Dublin - Ireland
Tel. +353 1 8107084
Fax +353 1 8107085
Email: Ireland@invacare.com
www.invacare.ie

© 2016 Invacare International Sàrl
All rights reserved.

All information quoted is believed to be correct at time of print. Invacare reserves the right to alter product specifications without prior consultation.

Reliant 350 - UK - 01/2016
Brochure re-order code - Q99632



safepatienthandling

The *Invacare Reliant 350 Stand Assist* is part of the 'Safe Patient Handling' range.

