# **AQUATEC®**



The ultimate bathing experience



The Aquatec Orca familiy provides superior quality and durable bath lifters that are designed to offer comfort, safety and reassurance to those most in need of support as they bathe. The Orca incorporates increased leg room, ease of disassembly and an ergonomic floating hand control that contains a highly efficient, light weight battery.

- Impressive weight limit of 140 kg
- Orca XL supports users up to 170 kg
- Sits lower and further towards the back of the bath
- Increased leg room for deeper and more relaxed bathing
- Increased height of the backrest adds support whilst upright or reclined

#### Comfortable, safe and highly efficient

The Aquatec Orca range has been designed with the user in mind, giving a comfortable and straightforward bathing experience. Key features include:

- · Reclining backrest
- Integrated head support
- · Quick release hand control
- · Auto unlock for quick removal from bath
- Backrest recline up to 40°

The reinforced scissor mechanism on the Orca range combines safe lifting with a compact construction while the minimum seat height of just 60 mm allows plenty of room for a deeper and more comfortable bath.





#### Aquatec Orca XL

A powerful bath lifter with a maximum user weight of 170 kg.



### Features and options



## **Textured seating surface**

The **Orca** has a textured surface allowing use without covers if desired. The side flaps are hinged flush to the main seating surface and allow comfortable and safe side transfers, without entrapment of skin.



#### Two-piece design

When pushed forward the **Orca**'s backrest folds down and automatically unlocks. The backrest and base section can then be removed from the bath with ease. The 2 sections fit neatly into one another to give an extremely compact size for storage.



#### Easy to clean

The open, clean-cut design of the **Orca** combines with its smooth surfaces to offer effortless cleaning.

#### Accessories Aquatec bath lift



#### **Reclining special** backrest

Individually adjustable side parts for secure positioning. **1573970** White, for

**Orca** (when ordered with special backrest. All bath lift) **1471623** White, for

**Orca** (for retrofitting)



#### Height-adjustable headrest

For the safe support of the head and neck area. Recommended in conjunction with the parts for the headrest assembly included.

15461 White, for all Orca



## **Custom-built side flaps**

Facilitates transfer from the lift. **1574239** Side flaps (pair), including cover,

white, for **Orca** (when ordered with lift) **1471606** Side flap (unit), without cover, for Orca (when retrofitting)



#### Chest and pelvic belt

Soft, skin friendly material.

**1528702** for **Orca** 



#### Stopper/abduction unit

Prevents slipping off the seat surface. Available for all battery-driven bath lifters.

1525932 - Blue (7cm height)

1525933 - Blue (12cm height)




## **Height adaptors**

to suit a full range of bath heights. The 2, 4 or be quickly fitted to the bath lift base.

15170 2 cm (4 pcs.) incl. conversion kit

**14694** Height adapter 2 cm (4 pcs) Extension kit without conversion kit

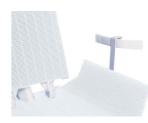

#### Side flap deflectors

The **Orca** can be adapted To be used for hand grips in the bathtub. The side flap deflector can 6cm height adapters can also be shortened when required.


**15309** for all battery Height adapter driven bath lifts (1 piece)


#### Technical data

## For more comprehensive pre-sales information about this product, including the product's user manual, please see your local Invacare website.

Dimensions	Aquatec <b>Orca</b>	Aquatec <b>Orca XL</b>
Seat width without side flaps	375 mm	375 mm
Seat width with side flaps	715 mm* (A)	715 mm (A)
Seat depth	590 mm (B)	590 mm (B)
Seat height max.	420 mm (C)	420 mm (C)
Seat height min.	60 mm	60 mm
Width of backrest	360 mm (D)	360 mm (D)
Height of backrest	665 mm (E)	665 mm (E)
Total depth (backrest upright)	650 mm (F)	650 mm (F)
Total depth (backrest reclined)	920 mm	920 mm
Total height with backrest	1075 mm (G)	1075 mm (G)
Base area	580 x 295 mm	580 x 295 mm
Cover thickness	3 mm	3 mm

More info			
Total weight without hand control	10,8 kg	11,6 kg	
Weight, seat	6,6 kg	7,4 kg	
Weight, backrest	4,2 kg	4,2 kg	
Weight, hand control	0.4 kg	0.4 kg	
Load capacity	140 kg	170 kg	
Backrest can be reclined by up to	40°	40°	
Colour	White	White	

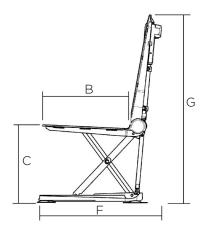

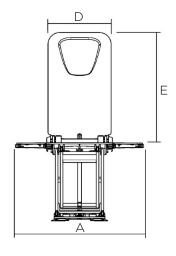

**Tel: +41 61 487 70 70** hqeurope@invacare.com www.invacare.eu.com

Invacare International GmbH Neuhofweg 51 4147 Aesch BL Switzerland

