

# Invacare® | Birdie® and Birdie® Compact



180kg weight limit  
5 year warranty

## A straightforward solution to safe and comfortable transfers

With maximum space for the client, both the *Invacare Birdie* and the *Birdie Compact* offer a truly comfortable lift and transfer to or from beds, chairs or even the floor. Ideal for use in domestic and residential care environments, the hoists are designed to fold and unfold easily without the need for tools.

Trusted Solutions, Passionate People





### Comfort and ease of use combined

The high boom maximises space for the client even when the hoist is in the highest position. This gives an ergonomic lift and transfer as the client can be easily rotated 360°, which makes positioning much easier for the carer. The optimisation of space allows increased space in front of the actuator, reducing the risk of injury to the client's knees.

### Ease of handling: simple to fold, unfold and dismantle without tools

Handling the hoist has been made as smooth and easy as possible as the hoists can easily be folded and unfolded without the need for tools. Hoists take up minimal space when stored and are easy to push and transport. If required, the mast and boom can be easily detached from the base (no tools required).

### TÜV-certification for maximum security

For total security, the folding mechanism has a snap-lock system, thus minimizing the risk of accident. The **Birdie** and **Birdie Compact** are CE-certified in line with the Medical Device Directive 93/42/ECC.

To ensure the hoists comply with EN ISO 10535 (the European standard for hoists) the hoist has been rigorously tested both by Invacare and the German TÜV institute, an accredited external body.

## Features and Options




#### Long boom - high user comfort

The client can easily be rotated and positioned. No risk of injuries to the knees.



#### Easy to fold and dismantle

Ease of transport and storage. No tools required when folding / unfolding or when dismantling the hoist.

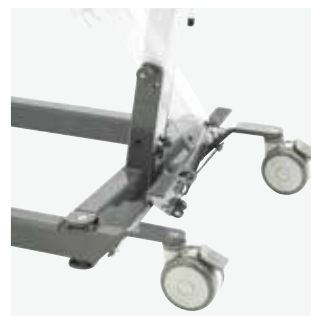

#### Snap lock

Optimal security.



#### Emergency stop, lowering and integrated charger

Manual or electric emergency lowering (compact version only electric lowering). Control box is supplied with detachable or fixed battery.

**Features and Options**



**Clip on spreader bar**

Easy to change the spreader bar - just twist the safety catch, press to open - and then detach/attach the spreader bar.



**Mounted spreader bar**

Correctly and safe positioning of the spreader bar. A locking carabiner offers exceptional security.



**Manual leg spread**

Simply step on the pedal to adjust the leg spread.



**Handle for leg spread**

Handle can easily be attached to operate the manual leg spread.



**Electrical leg spread**

Available for *Invacare Birdie* models only.



**2 point quick release spreader bars**

Easy to attach and detach. Available widths: 35 cm, 45 cm, 55 cm.



**Push bar and hand control**

Superior ergonomics for carers.




**Castors**

100 mm castors (braking castor is shown).



**Scale**

If needed, a Homecare scale is available approved, according to medical standard EN 45501.

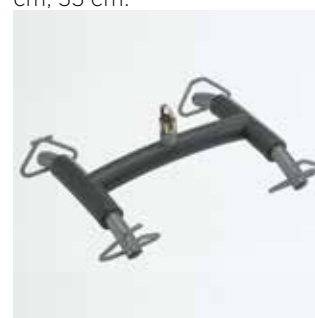

**4 point quick release spreader bars**

Easy to attach and detach. Available widths: 45 cm and 55 cm.

## Technical data

<b>Birdie</b>	124 cm	64 cm	42 kg	44,5 - 171 cm	10 cm
<b>Birdie Compact</b>	109 cm	52 cm	36 kg	52,5 - 161,5 cm	10 cm

<b>Birdie</b>	101 cm	180 kg	40 cycles	2,9 Ah	140 cm
<b>Birdie Compact</b>	84,5 cm	150 kg	40 cycles	2,9 Ah	107 cm

All measures are taken with 75 mm castors and 2 point spreaderbar.  
For 100 mm castors add 1,5 cm for height measures and 2 cm for width.

## Frame colour



Light grey /  
Dark grey

## Warranty

5 year warranty:  
Frame  
Actuator  
Control box

2 year warranty:  
Hand control

Invacare Limited  
Pencoed Technology Park -  
Pencoed - Bridgend CF35 5AQ -  
United Kingdom  
Tel: +44 1656 776 222  
Fax: +44 1656 776 220  
E-mail: [uk@invacare.com](mailto:uk@invacare.com)  
Sales Order E-mail:  
[ordersuk@invacare.com](mailto:ordersuk@invacare.com)  
[www.invacare.co.uk](http://www.invacare.co.uk)

Invacare Ireland Ltd  
Unit 5 - Seatown Business Campus  
Seatown Road - SWORDS - County  
Dublin - Ireland  
Tel. +353 1 8107084  
Fax +353 1 8107085  
Email: [ireland@invacare.com](mailto:ireland@invacare.com)  
[www.invacare.ie](http://www.invacare.ie)

© 2016 Invacare International Sàrl  
All rights reserved.

All information quoted is believed to  
be correct at time of print.  
Invacare reserves the right to alter  
product specifications without prior  
consultation.

The *Invacare Birdie* and *Birdie Compact* are part of the Safe Patient Handling range.



safepatienthandling



Yes, you can.®

Birdie - UK - 01/2016 - 1504431