

Etude[®] Plus Etude[®] Plus Low



Sleep safe and easy



Yes, you can.

The **Etude Plus** care bed, with its attractive aesthetics, flexible design and easy handling is ideal for homecare use. All the features of the popular Etude concept, such as sturdiness, ease of handling and durability, are very much present in this new version. **Etude Plus** is easy to transport and install in the home; the bed can be easily assembled and dismantled and can be used with a transport adaptor, all without compromising client comfort or safety. With the *Invacare* **Etude Plus Low**, thanks to the extra low height levels, safety is now even greater, ensuring greater protection from injury for those who are liable to fall out of bed. The low levels also make it easy for patients to get in and out of bed.



Modern and Flexible

The streamlined, discreet design perfectly complements any interior. The bed can be configured with steel or wood side rails, with or without an extra wood plate. Small details make a big difference to personal comfort: **Etude Plus** can be customised with auto-regression for increased comfort and reduction in shear and friction. Lockable hand control for increased safety. The standard, built-in double-brackets are designed to suit each client's unique transfer needs in different positions (entrance height of 400 mm, alternative of 330 mm.)

Easy handling & storage

Portability and storage are key benefits of the **Etude**. Its simplistic design allows easy assembly and dismantling without tools, and the process can be safely and comfortably managed by one person (each section weighs less than 25 kg). The bed ends can be easily attached to the transport unit making the **Etude** lightweight and manoeuvrable for handling and storage.

Highest possible safety levels

The **Etude Plus** conforms to IEC60601-2-52, the safety standard for beds. This is equal to the highest degree of safety for both client and carer. In addition to stringent internal tests, **Etude Plus** has been TÜV approved. **Etude Plus** is compatible with three different ranges of side rail. All meet the demands for increased strength required in the latest safety standard for beds.

Features and options



Built-In low version with Double Brackets

Easy and quick attachment of Mattress Support to bedends with different heights for different purposes available.



Easy Maintenance

IPX6 motors can be dismantled without tools. The cabling between the control unit and the motor is easy to access and replace.



Auto-Regression

Improved comfort when seated, reduced shear/friction.



Hand Control

Slim, ergonomic hand control. Choose with or without lockable functions.

Features and options



Etude Plus - Attractive and Homely

With its wooden plate on both sides of the bed end and full-length wooden side rail, the **Etude Plus** fits perfectly into any domestic environment, whether cosy or contemporary. The **Etude Plus** is height adjustable from 33-73 cm or 40-80 cm.



Etude Plus Low

The extra low height levels are an integral benefit of the **Etude Plus Low**, increasing safety and reducing the risk of harm for people who are prone to falling out of bed. The low levels also make it easy for patients to get in and out of bed. The **Etude Plus Low** offers a range of adjustable bed heights, including 20-60 cm and 27-67 cm.



Transport and Storage

Easy to store and to transport when placed on the transport kit.



Mattress Support Extension

For the ultimate in customisation and personal comfort, choose between an extension of 100 or 200 mm.



Cover for wooden side rail

Padded cover for wooden side rails. Cover remains in place when side rail is lowered. Available in standard, extended, short and net.



Britt V

Full-length side rails in beech. Compliant with IEC 60601-2-52.



Line/Line Extended

Full-length aluminium side rails. Compliant with IEC 60601-2-52.



Height Extender




Height of full-length side rails can be increased by 150 mm.



Transport Kit

For easy storage and enhanced portability.

Features and options

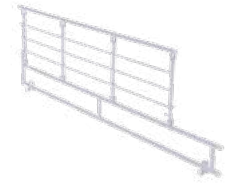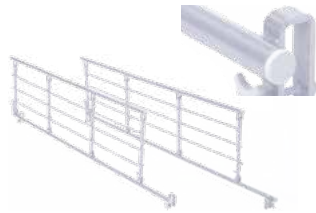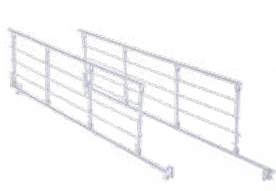

Invacare **Scala Basic 2**

Collapsible steel side rail.

Invacare **Scala Basic Plus 2**

Collapsible steel side rail with plastic inserts at the mounting brackets.

Invacare **Scala Medium 2**

Collapsible steel side rail for higher mattresses.

Invacare **Scala Decubi 2**

Collapsible steel side rail for extra high mattresses.

Compliant with the safety standard for beds IEC 60601-2-52.

Technical data

For more comprehensive pre-sales information about this product, including the product's user manual, please see your local Invacare website.

	Width outside Width inside	Length outside Length inside	Height adjustment	Mattress support dimensions	Backrest angle	Thigh angle	Legrest lift
Etude Plus	1020 mm 900 mm	2190 mm 2070 mm Mattress: 2000mm	330 - 730 mm or 400 - 800 mm	850 - 240 - 350 - 560 mm	0 - 70°	0 - 18°	0 - 25°
Etude Plus Low	1020 mm 900 mm	2190 mm 2070 mm Mattress: 2000mm	200 - 600 mm or 270 - 670 mm	850 - 240 - 350 - 560 mm	0 - 70°	0 - 18°	

	Tilt option ¹	Total weight	Max. user weight	Heaviest part of the product	
Etude Plus	11°	76,5 kg	SWL 200 kg 165 kg	24,6 kg	Complies with IEC 60601-2-52
Etude Plus Low	11°	76,5 kg	SWL 200 kg 165 kg	24,6 kg	

1. Please notice that the tilt function is only practised for users under medical supervision.

NB

The bed must not be used by patients under 12 years of age, or by patients with body size equivalent to an average 12 year old or smaller.

Frame colour and wood work

Aluminium grey

Beech



Please notice colors may vary slightly from those displayed above.

Tel: +44 61 487 70 70
hqueurope@invacare.com
www.invacare.eu.com

Invacare International GmbH
Neuhofweg 51
4147 Aesch BL
Switzerland



Yes, you can.®



© 2023 Invacare International GmbH. All rights reserved. All information quoted is believed to be correct at time of print. Invacare reserves the right to alter product specifications without prior consultation. EU - 04/23