



Kompas[™]

LIGHTWEIGHT | COMPACT | CONVENIENT



Yes, you can.™

The *Invacare Kompas* combines the freedom of a powered mobility device with the compact design of a manual wheelchair, making it ideal for individuals that want to maintain an active and independent lifestyle, without compromising on space or convenience. As one of the lightest and most compact folding power wheelchairs on the market, the *Kompas* is highly versatile and the perfect powered solution for those wanting to get the most out of life.

Easily folded and compact

Thanks to an intuitive one-step design, the *Kompas* can be quickly and easily folded to minimize time and effort. When collapsed, the power wheelchair has an impressively small footprint, comparable to the size of a manual folding wheelchair.

This compact size, coupled with the ability to stand upright when folded, makes the *Kompas* ideal for storing in tight spaces such as under the stairs or in between furniture.

Extremely lightweight and transportable

Weighing less than 21kg*, the *Kompas* is the most compact and lightweight power wheelchair within Invacare's extensive portfolio.

Optimised for ease of handling, the power wheelchair can be effortlessly lifted when necessary, for example to place inside a property or to transport in the boot of a car.

Versatile, convenient and fun!

Whether for home use, exciting daytrips or everyday outings, the highly maneuverable *Invacare Kompas* can help provide assistance for up to 11km on a variety of journeys. Simply fold and go!

*Weight without battery included.



5 KEY FEATURES

1 Fully Adjustable

- ▶ Height adjustable backrest provides additional support and comfort for taller individuals.
- ▶ Supportive flip-up anti-skid pedal has 3 height adjustments for additional convenience.

2 Removable Lithium Battery

- ▶ Airline approved removable 11.4ah battery with additional charging port and battery indicator.
- ▶ Unlike many other folding power wheelchairs, there is no need to remove the battery pack to fold.

5 Lightweight Frame

- ▶ With an impressive frame weight of less than 21kg, the compact *Kompas* can be effortlessly lifted and carried when necessary.

3 One-step Folding Mechanism

- ▶ Folding and unfolding is simple and effortless.
- ▶ Compact dimensions of 590 x 370 x 810mm when folded.

4 Anti-tippers

- ▶ Accessible when folded to allow easy maneuverability for transport and additional stability for storage.

Additional features


- ▶ **Durable UK P&G controller**
 - High quality controller indicates the remaining battery charge and allows for easy speed adjustments and nimble movement control.
 - Equipped with convenient on-board charging port.
- ▶ **Intelligent Braking System**
 - Automatically prevents unintended rolling whilst stationary, effectively locking the power wheelchair in place.
- ▶ **Rear Pneumatic Tyres**
 - Facilitates a smoother journey and are suitable for outdoor use.
- ▶ **Storage Basket**
 - Under seat, removable nylon basket enables transportation of small items.
- ▶ **Freewheel Levers**
 - Easy to operate freewheel levers to disengage the motors for manual transportation.
- ▶ **Optional 17.1Ah Lithium Battery**
 - Larger battery available as an accessory to increase travel range up to 18km.

Technical data ▼

For more comprehensive pre-sales information about this product, including the product's user manual, please see your local Invacare website.

Seat width ▼	Seat depth ▼	Seat to floor height ▼	Folding dimensions ▼	Driving unit width ▼	Total length ▼	Total height ▼
Kompas 460mm	420mm	520mm	590 x 370 x 810mm	590mm	970mm	950-1000mm
Total weight ▼	Max. user weight ▼	Battery capacity ▼	Battery weight ▼	Motor capacity ▼	Turning diameter ▼	Ground clearance ▼
Kompas 20.9kg excluding battery	115kg	24 volt/11.4ah lithium ion battery	2.1kg	DC24V,70 watt	825mm	30mm
Max. safe slope ▼	Max. driving distance ▼	Max. driving speed ▼	Tyre size ▼	Castor size ▼		
Kompas 6°	11km / 6.8miles	6.4 kph / 4 mph	8" pneumatic	7" PU		

Frame colour ▼

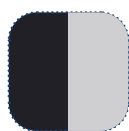

Black & Grey

Upholstery ▼



Black

Tel: +44 61 487 70 70
hqeuropa@invacare.com
www.invacare.eu.com

Invacare International GmbH
 Neuhofweg 51
 4147 Aesch BL
 Switzerland



Yes, you can.®



© 2023 Invacare International GmbH. All rights reserved. All information quoted is believed to be correct at time of print. Invacare reserves the right to alter product specifications without prior consultation. EU - 07/23