

# Medley Ergo

Providing clients comfort and support



**An advanced range of care beds**

for increased comfort and improved clinical outcomes



**Yes, you can.**

# Designed to care for the client and caregiver

Our Medley Ergo range of beds are an easy to handle solution that help deliver excellent client care and provide maximum comfort.

Designed to increase independence and improve clinical outcomes for varying needs, with a choice between our Medley Ergo or Medley Ergo Low and the availability of Medley Ergo Select for a more homely design with bed end covers.

Built-in features and functionalities ensure the client as well as the caregiver can operate the bed easily. Optimised dimensions of the mattress platform help to reduce shear and friction and in combination with one of our range of pressure redistributing mattresses, can provide optimal pressure care and micro climate management.



## 1



### Easy Handling

An optional **transportation fitting** is available so that the bed can easily be **collapsed, moved and re-installed** in a variety of locations.

## 2



### Comfort

In all resting and sitting positions with **minimised risk of slipping downwards** during profiling with optimised mattress base dimensions.

## 3



### Support

Easy repositioning of the client with hand controls that help to **increase independence** and **support the care giver**. Available with or without lockable functions.



› **Mattress support**  
Available in easy to clean steel.

› **Side rails**  
Choose from 3 options that vary in length, height and finish.  
Full length - Lisa wood finish in beech or Aria in steel  
3/4 length - Scala, collapsible in steel

› **Select design**  
A more homely design can be achieved with retro-fitted bed end covers for our Medley Ergo Select option.

› **Hand control**  
With or without lockable functions.



› **Double height positions for bed ends**  
Adjustable height ranges for  
Medley Ergo: 33 - 73 cm or 40 - 80 cm  
Medley Ergo Low: 21 - 61 cm or 28 - 68 cm

### Medley Ergo Low

- › **Positioning the bed at a lower level** increases security for clients and provides peace of mind when going to sleep. Caregivers can raise the bed to a working height when required.
- › Clients can get **in and out of bed easier** increasing independence.
- › Adjustable heights from **21 - 61 cm** and **28 - 68 cm**.

Options to choose from:



Light and easily adjustable lifting Pole



Extension kit for increasing mattress support length by 15 cm



Transport fitting for easy storage and transport



Longer side rail to fit extended bed

# Helping to deliver optimal client care

A combination of advanced functionalities, sleek features and uncompromised quality for enhanced client care.



## Shear and Friction Reduction

Optimised dimensions of the mattress base help reduce shear and friction and increase comfort, according to expert recommendations for anthropometrical measurements and pressure mapping.



## Easy to store and transport

Tool-less assembly / disassembly in combination with our optional transport fitting enables easy transportation and storage when the bed is collapsed and placed on the dedicated transporter.



## Safety and quality assured

Safety standard IEC 60601-2-52 compliant, our trusted Medley Ergo range can be relied upon for its high safety standards and optimal quality, helping to deliver reliable long-term support for clients at home or in long-term care environments.



## Technical data ▼


For more information about this product, including the user manual, please visit our website: [www.invacare.eu](http://www.invacare.eu)




Width outside  
Width inside ▼



Length outside  
Length inside ▼

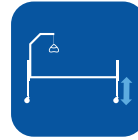


Height adjustment ▼

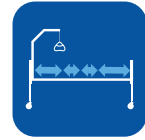


Mattress support  
dimensions ▼


|                 | Width outside<br>Width inside ▼ | Length outside<br>Length inside ▼ | Height adjustment ▼             | Mattress support<br>dimensions ▼ |
|-----------------|---------------------------------|-----------------------------------|---------------------------------|----------------------------------|
| Medley Ergo     | 1020 mm<br>900 mm               | 2240 mm<br>2090 mm                | 400 - 800 mm or<br>330 - 730 mm | 800 - 230 - 300 - 670 mm         |
| Medley Ergo Low | 1020 mm<br>900 mm               | 2240 mm<br>2090 mm                | 280 - 680 mm or<br>210 - 610 mm | 800 - 230 - 300 - 670 mm         |

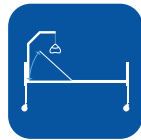


Backrest angle ▼

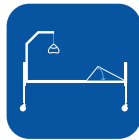

Thigh angle ▼

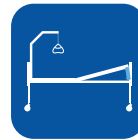

Legrest lift ▼

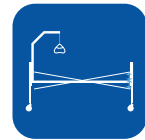

Tilt option<sup>1</sup> ▼

|                 | Backrest angle ▼ | Thigh angle ▼ | Legrest lift ▼ | Tilt option <sup>1</sup> ▼ |
|-----------------|------------------|---------------|----------------|----------------------------|
| Medley Ergo     | 0 - 70°          | 0 - 24°       | 0 - 18°        | 11°                        |
| Medley Ergo Low | 0 - 70°          | 0 - 24°       | 0 - 18°        | 11°                        |



Total weight ▼



Max. user weight ▼



Heaviest part of the  
product ▼

|                 | Total weight ▼ | Max. user weight ▼   | Heaviest part of the<br>product ▼ |
|-----------------|----------------|----------------------|-----------------------------------|
| Medley Ergo     | 64 kg          | SWL 180 kg<br>145 kg | 17 kg                             |
| Medley Ergo Low | 64 kg          | SWL 180 kg<br>145 kg | 17 kg                             |

Complies to the safety standard, IEC 60601-2-52

1. Please notice that the tilt function is only practised for users under medical supervision.

NB: The bed must not be used by clients under 12 years of age, or by clients with body size equivalent to an average 12 year old or smaller.

Silver

Beech

Please note colours may vary slightly from those displayed above.



**Tel: +44 61 487 70 70**  
[hqeuropa@invacare.com](mailto:hqeuropa@invacare.com)  
[www.invacare.eu.com](http://www.invacare.eu.com)

Invacare International GmbH  
 Neuhofweg 51  
 4147 Aesch BL  
 Switzerland



**Yes, you can.®**