

NordBed™ Ultra



ERGOMOVE™ TECHNOLOGY | FALL PREVENTION | INDEPENDENCE



Yes, you can.®

The innovative NordBed Ultra

Clients within a nursing home have a wide range of needs, in particular, those who are more dependent on others and those with dementia, therefore, the NordBed Ultra is a perfect choice as it can be **adapted easily** to suit a range of clients' requirements. This bed also has a **homely aesthetic** feel, ideal for a nursing home atmosphere. The advanced built-in auto-regression **reduces shear significantly** when bringing the client to a sitting position. **Promoting mobilisation** to ensure **independence** and meeting the needs for the future, the all-rounder NordBed Ultra is the advanced care-bed for clients and care staff.

ErgoMove technology

The Invacare NordBed Ultra includes ErgoMove technology to bring the client from a lying position to an active sitting position and back to a horizontal position **comfortably** with a push of a button. The **active sitting position** promotes client function, such as eating, drinking, reading and ensures good pressure distribution. The movement sequence ensures the pivot points at the backrest and knee break match as closely as possible the body's own natural pivot points, at the hip and knee. This guarantees that the client stays in position and does not slide down the bed. This combined with the built-in auto-regression **reduces shear significantly** when bringing the client to a sitting position and enhances pressure distribution. Supporting **single handed care**, this intuitive feature **reduces the need for care staff to reposition the client**. The optimised length of mattress support sections are set according to anthropometric measurements, which allows the profiling sections to fit to and support most adults.

PATENT PENDING



Safety, independence and mobilisation

The NordBed Ultra has a lowest height of only 25 cm :

- ▶ **Increasing safety** and reducing the risk for people who are prone to falling out of bed
- ▶ Prevents fall injuries without the use of side rails, this is a **significant benefit for those who are suffering from dementia**
- ▶ Makes it easier for clients to transfer in and out of bed
- ▶ And even allows shorter clients to place both feet safely on the floor when getting to a standing position

Up-and-Out™

- ▶ To ensure a safe and smooth transfer and **supporting independence and mobilisation**, clients are able to use the Up-and-Out feature with a simple push of a button. The profiling sections will **automatically adjust** to a suitable position for transferring in and out of bed, the movement is optimised to minimise shear

Design

- ▶ The NordBed Ultra is available in a variety of modern bed end styles, wood finishes and painted surfaces. Depending on the need and **personal taste**, NordBed Ultra can be configured with side panels, side rails and more



Perfect for clients suffering from dementia

- ▶ The low height allows a safe environment for mobile people
- ▶ The highly configurable bed can be adapted:
 - To look like the more traditional bed
 - To fit in nicely with modern homely looking beds (box mattresses) for easy recognition
 - With a range of soft covers (high/low texture) for a cosy, confident and safe environment
 - Under bed light to reduce the risk of the client being injured while trying to get in/out of bed

4

EXCLUSIVE FEATURES

1

ErgoMove technology

- ▶ Brings the client to a comfortable sitting position and back to a horizontal position without pushing the client downwards and without the need of repositioning
- ▶ Reduces shear

3

Up-and-Out

- ▶ Allows safe and smooth transfers
- ▶ Supports the ergonomic movement

2

25 cm lowest height

- ▶ With standard castors at 10 cm to allow safe transfers and increase fall prevention

4

Optimised dimensions of mattress support

- ▶ According to anthropometric measurement, for increased comfort

Additional features

- ▶ **RememberMe™:**
 - Programmable LCD handset
 - ErgoMove and Up-and-Out can be adjusted to fit individual needs
- ▶ **Integrated extendable mattress platform**
 - Additional 10 cm at the head and foot ends to accommodate clients' requirements
- ▶ **Anti-Trendelenburg and *Optional* Trendelenburg**
 - Supports care staff when repositioning the client in bed
 - Easy to redistribute pressure
- ▶ **Under bed light**
 - Easy to turn on the light
 - Fall prevention/safety



Hand control ▲

Easy recognisable colours for clients who are weak-sighted

- ▶ The key buttons, which will be used most frequently, have easily identifiable colours. These colours have been chosen in collaboration with experts
- ▶ To assist the client when selecting a feature, all buttons will light up when any button is touched for easy recognition during the day or night
- ▶ For a simple selection, the key buttons are positioned at the top of the handset, and the adjustment features are placed below the up/down buttons
- ▶ The hand control is lockable as standard

Technical data ▼


For more comprehensive pre-sales information about this product, including the product's user manual, please see your local Invacare website.




Width outside
Width inside ▼



Length outside
Length inside ▼

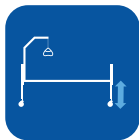

Height adjustment ▼

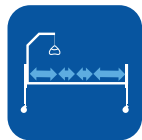

Mattress support dimensions ▼



Under bed clearance ▼




Backrest angle ▼



Thigh angle ▼

NordBed Ultra	100 / 90 cm	213 / 200 cm	25 - 85 cm	80-17-30-65 cm	15 cm	0 - 70°	0 - 49°
----------------------	-------------	--------------	------------	-------------------	-------	---------	---------

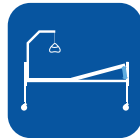

Legrest lift ▼



Anti-Trendelenburg / Trendelenburg ▼



Max. user weight ▼



Total weight (without bed ends) ▼



Heaviest part of the product ▼




Washable in washing machine ▼

NordBed Ultra	0 - 19°	+/- 15°	200 kg	105 kg	55 kg	Yes
----------------------	---------	---------	--------	--------	-------	-----

NB: The bed must not be used by patients under 12 years of age, or by patients with body size equivalent to an average 12 year old or smaller.

Frame colour ▼

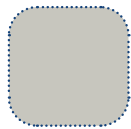

Aluminium grey

Woodwork ▼



Beech



Oak



Umber



Beech/pale blue



Pale blue



White

Tel: +41 61 487 70 70
hqeurope@invacare.com
www.invacare.eu.com

Invacare International GmbH
Neuhofweg 51
4147 Aesch BL
Switzerland



Yes, you can.®



© 2023 Invacare International GmbH. All rights reserved. All information quoted is believed to be correct at time of print. Invacare reserves the right to alter product specifications without prior consultation. EU - 04/23