

SB[®] 400



Adaptable for every purpose



The *Invacare SB 400* is an interaction between Scandinavian design and functionality. There is no better choice in modular systems, the **SB 400** has it all: A broad range of functionality and options for individual positioning, transport and storage.



Comfort, style & freedom of choice

The *Invacare SB 400* modular bed system is a byword for style and sophistication. Designed to meet every individual comfort need, it boasts a plethora of ergonomic benefits, including reclining backrest function that minimises anti-shear forces.

Central braking and an increased lifting capacity effectively optimises support for carers working with the bed, and in the transfer and transport of its users.




Handling

The design of the **SB 400** ensures trouble-free handling by the individual, without the need for tools. Snap-lock fixtures and pieces weighing less than 25 kg make the **SB 400** a light, user-friendly home care bed: the epitome of Scandinavian design-engineering excellence.

Safety

Smooth, uninterrupted use and longevity is facilitated by a battery backup which allows operation anywhere in the home. Furthermore, the **SB 400** has a proven track record in both safety and reliability. Don't just take our word for it! Listen to one of the major European test institutes (TÜV).

Features and options

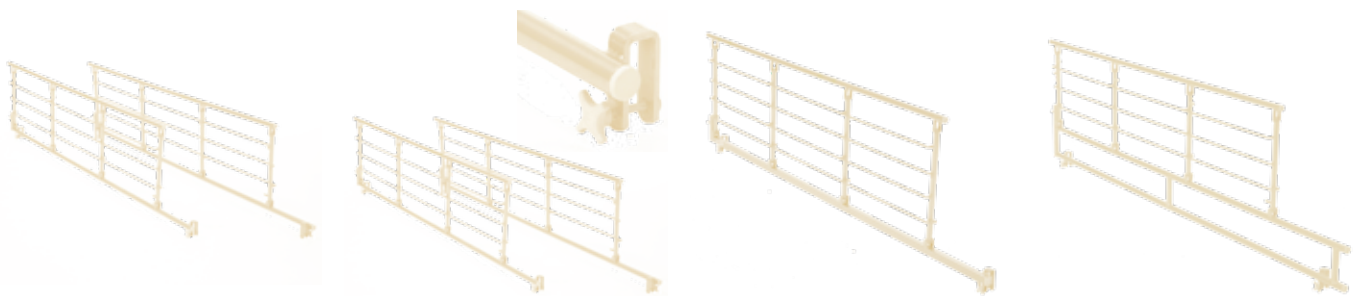

Invacare Scala Basic 2

Collapsible steel side rail. Compliant with IEC 60601-2-52.

Invacare Scala Basic Plus 2

Collapsible steel side rail with plastic inserts at the mounting brackets. Compliant with IEC 60601-2-52.

Invacare Scala Medium 2

Collapsible steel side rail for higher mattresses. Compliant with IEC 60601-2-52.

Invacare Scala Decubi 2

Collapsible steel side rail for extra high mattresses. Compliant with IEC 60601-2-52.

Features and options



Working height

Height range of 35-90 cm for optimum working position.



Hand control

The colours makes it easy to benefit from all the positioning options.




Snaplock function

Snaplock function separates mattress support from base.



Dividable base

Dividable base for easy handling of the bed.

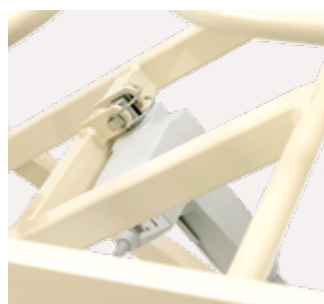

Easy maintenance Motors

Motors can be easily disengaged without the need for tools.



Bed Extension

Extension of the leg section with additional support.



Detachable bed ends

Detachable bedends in a range of designs in beech, white beech, cherry or light oak.



Bed ends

A variety of bed ends in beech, white beech, cherry and light oak wood.



Siderails

Siderails for wooden covering in standard and extended length in beech, white beech, light oak or cherry wood.



Backrest

Reclining backrest for extra comfort.




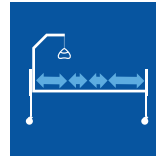
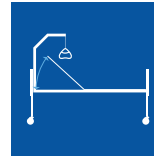

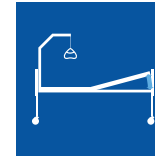




Central brake

Central braking function with steering control.

Technical data

For more comprehensive pre-sales information about this product, including the product's user manual, please see your local Invacare website.

	Width outside Width inside	Length outside Length inside	Height adjustment	Mattress support dimensions	Backrest angle	Thigh angle	Legrest lift
SB 400	 1000 mm 900 mm	 2200 mm 2070 mm Mattress: 2000 mm	 350 - 900 mm Total height: 1550 - 2100 mm	 4-sectioned: 850 - 240 - 350 - 560 mm 3-sectioned: 850 - 240 - 910 mm	 0 - 70°	 0 - 25°	 0 - 27°


	Under bed clearance	Max. user weight
SB 400	 150 mm	 SWL 180 kg 145 kg

Complies with IEC 60601-2-52

NB

The bed must not be used by patients under 12 years of age, or by patients with body size equivalent to an average 12 year old or smaller.

Frame colour and wood work

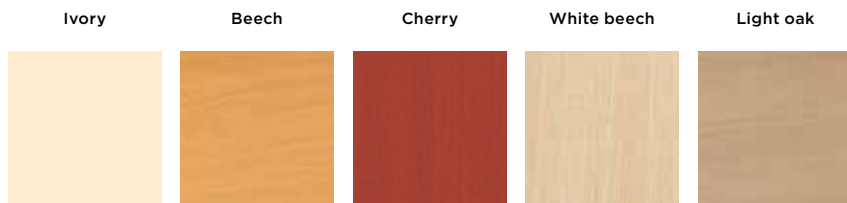

Please notice colors may vary slightly from those displayed above.

QR
Code

Tel: +41 61 487 70 70
hqeurope@invacare.com
www.invacare.eu.com

Invacare International GmbH
Benkenstrasse 260
4108 Witterswil
Switzerland

