

SB[®] 755

When comfort comes first



An adaptive support system

for varying client needs



Yes, you can.®

Comfort and support in all resting positions

Combining comfort, ease of handling and safety, our SB 755 care bed fits well into any home-care or long-term care environments.

Offering optimum comfort with an ideal seating position to minimise shear and friction. Caregivers can easily adjust the bed to meet individual client requirements.

Featuring an extendable head and foot end, the SB 755 is suitable for taller clients. Both

caregivers and clients can experience improved comfort and general well-being.

With the option to choose a bed that is specially adapted to withstand cleaning in a washing machine, the SB 755 is now even easier to maintain.

Easy to store and transport

Mounting and dismantling the SB 755 has been made simple for ease of transport and storage. Simply position the bed on its base and wheel it from the vehicle to the bedroom. Assembly can be carried out quickly and easily by one person. No part of the bed weighs more than 25 kg and dismantable scissors are available as an option.



Safety and quality assured

Safety standard IEC 60601-2-52 compliant, our trusted SB 755 range can be relied upon for its high safety standards and optimal quality, helping to deliver reliable long-term support for clients.




1 Auto-regression function

Client comfort when activating backrest, minimising shear and friction.



2 Increased comfort

Optimised dimensions of mattress platform according to anthropometrical measurement, for increased comfort.

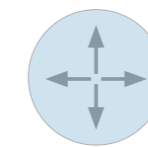

3 Extendable mattress platform

Additional 100 mm at head and foot. Bed easily extends with 50, 100, 150 or 200 mm to suit client's needs.



> **Mattress support**
Mattress support with choice of aluminium slats.

> **Side rails**
Choose from 3 options that vary in length, height and finish.
Full length - Britt V or Line
3/4 Length - Scala

> **Bed ends**
Available in a variety of designs and finishes.

> **Hand control**
Available as lockable, not lockable or as a Soft Control.

Bed ends



> **Olivia** Can be used in combination with wooden front panel. Available in beech.



> **Victoria** Detachable panel. Available in beech, light oak and white beech.



> **Sophie** Available in beech, light oak and white beech.



> **Emma** Detachable - Available in beech, light oak and white beech.



> **Susanne** Available in beech, white beech and light oak.



> **Anita** Available in beech and light oak.

Options to choose from:



Wooden side
Available in beech, white beech or light oak.



Lifting pole with adjustable handle in height and depth.



Handle to support the client reposition and get in and out of bed.



Choice of 100 or 125 mm castors.



Longer side rail to fit extended bed

SB® 755 Wide

Taller, wider, bigger...

Our popular SB 755 is available in wider variants of 1050 mm and 1200 mm and can be extended at both ends, providing even more comfort for taller or bigger clients.

These wider models provide extra space and comfort for larger clients and plenty of room to turn. The SB 755 Wide has been engineered with an improved auto-regression system to offer total comfort. The bed also offers an optimal seating position and is designed to

reduce shear and friction. The SB 755 Wide is attractive and modern; the streamlined design is easy to clean and will blend harmoniously into any homecare or long-term care environment.



Technical data ▼

For more information about this product, including the user manual, please visit our website: www.invacare.eu




Width outside
Width inside ▼

SB 755	950/850 mm 1000/900 mm
SB 755 Wide 1050 mm	1150 mm 1050 mm
SB 755 Wide 1200 mm	1300 mm 1200 mm




Length outside
Length inside ▼

2130 mm 2070 mm Mattress: 2000mm
--

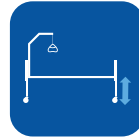

Height adjustment ▼

350 - 820 mm Total height: 1600 - 2070 mm

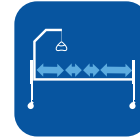


Mattress support
dimensions ▼

800 - 260 - 260 - 680 mm




Under bed clearance ▼

150 mm

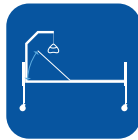


Backrest angle ▼

0 - 73°

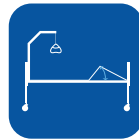

Thigh angle ▼

0 - 28°

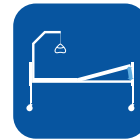

Legrest angle ▼

0 - 18°



Legrest lift ▼

0 - 18°



Angle below/above
frame ▼

below: 25° above: 14°



Max. user weight ▼

SWL 235 kg 200 kg



Heaviest part of the
product ▼

SB 755: 25 kg
SB 755 Wide 1050 mm: 40 kg
SB 755 Wide 1200 mm: 44 kg



Complies to the safety standard, IEC 60601-2-52

NB: The bed must not be used by clients under 12 years of age, or by clients with body size equivalent to an average 12 year old or smaller.

Light grey

Beech

Light oak

White
beech

Please note colours may vary slightly from those displayed above.



Tel: +41 61 487 70 70
hqeurope@invacare.com
www.invacare.eu.com

Invacare International GmbH
Neuhofweg 51
4147 Aesch BL
Switzerland



Yes, you can.®