

Action[®] Ampla[™]



Providing independence for an active lifestyle

INNOVATIVE | AGILE | AWARD WINNING



Yes, you can.™

Taking bariatric wheelchairs to the next level

The Invacare **Action® Ampla™** has been designed and developed specifically for Plus-Size users, taking into consideration body shapes, varying activity levels and weight limit requirements. But it doesn't stop there! The sleek and minimalist wheelchair, with its modern appearance, provides independence and support for a comfortable and active lifestyle.



Award winning

Already an award-winning chair, the agile and stylish **Invacare Action® Ampla™** has won a Red Dot Award for Product Design. This award is not only the largest product competition in the world, but also one of the most prestigious design awards. The **Action® Ampla™** has also won the Italian Platinum A'Design Award and the Gold German Design Award.



7 Inbuilt transportation solution

The transportation bearings on the **Action® Ampla™** have been cleverly designed to take the weight of the wheelchair, allowing it to easily go through doorways or narrow spaces:

1. Fold down the backrest
 2. Lift the legrest section of the wheelchair
 3. Tilt wheelchair onto the transportation bearings
- Whilst in this position, the wheelchair can be moved easily.



5 KEY FEATURES



2 Ergonomic backrest

- ▶ Maximum access to handrims
- ▶ No obstructions
- ▶ More freedom of movement
- ▶ Superb driving performance and manoeuvrability



3 Height adjustable transfer levers

- ▶ Support for entering and leaving the chair
- ▶ Can be turned inwards or outwards
- ▶ Can be lowered to allow table access



4 Personalised accent

- ▶ Three stylish colour choices




5 Multiple active rear wheel positions

- ▶ Very easy to propel
- ▶ Users can go further, using less energy

Innovative backrest design



The design of the wheelchair allows for a perfect fit for both Plus-Size male users, who tend to have more of an "apple" body-shape and Plus-Size female users, who tend to have a "pear" body-shape. This is thanks to the innovative backrest on the **Action® Ampla™** that has been designed to accommodate the anatomy of both body-shapes whilst offering support and comfort to the user.

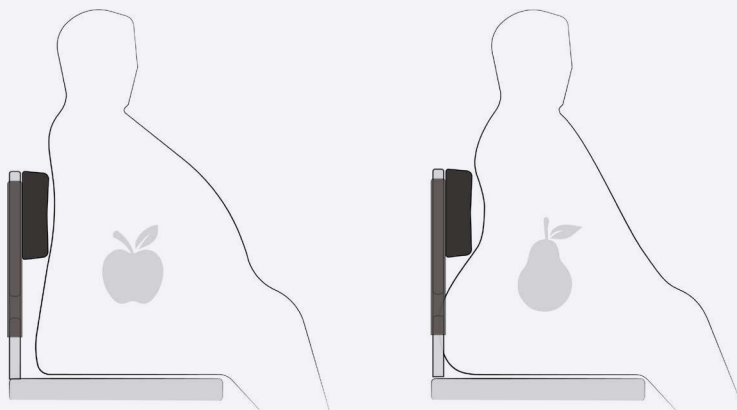

Main features

- ▶ Inbuilt transportation solution
- ▶ Ergonomic backrest
- ▶ Height adjustable transfer lever
- ▶ Personalised accent
- ▶ Multiple active rear wheel positions

Additional features


- ▶ An elastic, breathable privacy cover hides the gluteal shelf
- ▶ Push to lock brakes are nicely integrated into the side guard for accessibility

Technical data ▼

For more comprehensive pre-sales information about this product, including the product's user manual, please see your local Invacare website.

Seat width ▼	Seat depth ▼	Seat angle ▼	Seat height ▼	Backrest angle ▼	Backrest height ▼	
Action Ampla	505 / 555 / 605 / 655 / 705 / 755 mm	460 / 510 mm	1°	380 - 455 mm in steps of 25 mm	0°	450 - 525 mm
Front frame angle ▼	Length of footrests ▼	Driving unit width ▼	Total length inc. footrests ▼	Total length without footrests ▼	Total height ▼	
Action Ampla	90 - 100°	365 - 525 mm in steps of 25 mm	Seat width + 205 mm	1065 - 1145 mm	855 - 920 mm	825 - 900 mm
Max. safe slope ▼	Radius ▼	Max. user weight ▼	Total weight from ▼	Transportation weight from ▼		
Action Ampla	7°	1420 - 1660 mm	300 kg	27 kg	19.5 kg	

Frame colours ▼

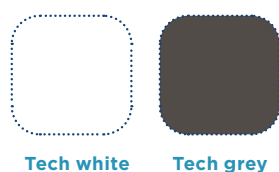

Accents ▼



Please note: colours may vary slightly from those displayed here.

Regulatory data ▼



Tel: +44 61 487 70 70
 hqeuropa@invacare.com
 www.invacare.eu.com

Invacare International GmbH
 Neuhofweg 51
 4147 Aesch BL
 Switzerland

