

# NordBed™ Optimo Family



ERGOMOVE™ TECHNOLOGY | IDEAL FOR REFURBISHMENT | INDEPENDENCE



Yes, you can.®



The NordBed Optimo family are the ideal beds for homecare environments. With the **exceptional comfort** and **ease of use** for the client, the NordBed Optimo is continuously promoting independence and ensuring mobilisation. With the ErgoMove technology, pressure care is at the heart of the NordBed Optimo. The NordBed Optimo 90's **simplicity** ensures easy handling, not only for the client but for the carer and technician. The NordBed Optimo Wide allows space for plus size clients and accommodates clients who require additional space when repositioning themselves in bed.

## ErgoMove technology

The Invacare NordBed Optimo includes ErgoMove technology to bring the client from a lying position to an **active sitting position** comfortably with a push of a button. The active position is more suitable for eating, drinking and reading and ensures good pressure distribution.

The movement sequence ensures that the client stays in position and does not slide down the bed. This combined with the built-in auto-regression **reduces shear significantly** when bringing the client to a sitting position and enhances pressure distribution.



### Up-and-Out™

To ensure a **safe and smooth transfer**, supporting independence and mobilisation, clients are able to use the Up-and-Out feature with a simple push of a button. The profiling sections will automatically adjust to a suitable position for transferring in and out of bed, the movement is **optimised to minimise shear**.



6

## EXCLUSIVE FEATURES

1

### ErgoMove technology

- ▶ Brings the client from a laying position to an active sitting position comfortably

2


### Ease of handling

- ▶ Easy dismantling for transport and storage. Ideal for refurbishment

3

### Range of sizes

- ▶ Available in 90, 105 and 120 cm wide

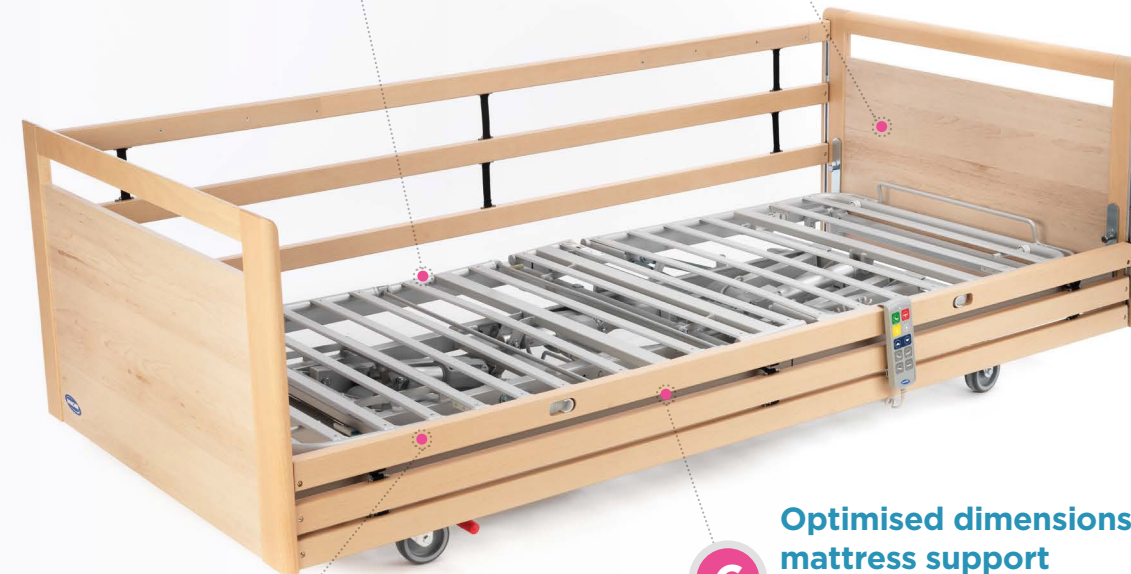

4

### Up-and-Out

- ▶ Ensures a safe and smooth transfer

5

### Design

- ▶ Available in a variety of modern bed end styles, wood finishes and painted surfaces

6

### Optimised dimensions of mattress support

- ▶ According to anthropometric measurement, for increased comfort

### Additional features

#### ▶ RememberMe™:

- Programmable LCD handset
- ErgoMove and Up-and-Out can be adjusted to fit individual needs

#### ▶ Integrated extendable mattress platform

- Additional 10 cm at the head and foot ends to accommodate clients' requirements

#### ▶ Anti-Trendelenburg and *Optional Trendelenburg*

- Supports care staff when repositioning the client in bed
- Easy to redistribute pressure

#### ▶ Under bed light

- Easy to turn on the light
- Fall prevention/safety

## Key features and benefits

### ◀ Ease of handling

- NordBed Optimo 90 cm and NordBed Optimo Wide 105 cm, dismantle into just five lightweight parts for ease of handling and easy storage and transportation. Each part weighing 24 kg or less on the 90 cm model and 26 kg or less on the 105 cm model.
- A compact design that can be easily manoeuvred through narrow spaces with the dismantled bed placed onto the base.
- The NordBed Optimo 120 cm model dismantles into three parts and can be effortlessly placed on a trolley for easy manoeuvring.



### ◀ Ideal for refurbishment

- The NordBed Optimo range has an option to be specially adapted to withstand cleaning in a washing machine.






### ◀ Hand control

- The key buttons, which will be used most frequently, have easily identifiable colours. These colours have been chosen in collaboration with experts
- To assist the client when selecting a feature, all buttons will light up when any button is touched, for easy recognition during the day or night
- For a simple selection, the key buttons are positioned at the top of the handset, and the adjustment features are placed below the up/down buttons
- The hand control is lockable as standard



# Introducing the **NordBed Optimo** family


**NordBed Optimo is available in three models with different load capacities and dimensions for all homecare environments.**

	Dismantled into:	Max user weight	Heaviest part of the product
<p><b>Nordbed Optimo 90</b></p> <p><i>Continuously promoting independence and ensuring mobilisation</i></p> 	5 parts	200 kg	24 kg
<p><b>Nordbed Optimo 105</b></p> <p><i>Wider size accommodates clients who require additional space when repositioning themselves in bed</i></p> 	5 parts	200 kg	26 kg
<p><b>Nordbed Optimo 120</b></p> <p><i>Extra comfort for the plus size</i></p> 	3 parts	250 kg	70 kg

A range of configurations is available to meet the needs of your carer and client.

## Technical data ▼

For more comprehensive pre-sales information about this product, including the product's user manual, please see your local Invacare website.

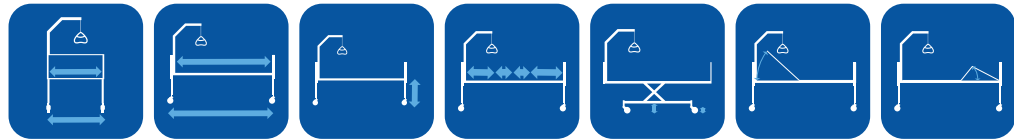


	Width outside Width inside ▼	Length outside Length inside ▼	Height adjustment ▼	Mattress support dimensions ▼	Under bed clearance ▼	Backrest angle ▼	High angle ▼
<b>NordBed Optimo 90</b>	100 / 90 cm	213 / 200 cm	32 - 81 cm	80-17-30-65 cm	15 cm	0 - 70°	0 - 49°
<b>NordBed Optimo 105</b>	115 / 105 cm	213 / 200 cm	32 - 81 cm	80-17-30-65 cm	15 cm	0 - 70°	0 - 49°
<b>NordBed Optimo 120</b>	130 / 120 cm	213 / 200 cm	32 - 81 cm	80-17-30-65 cm	15 cm	0 - 70°	0 - 49°



	Legrest lift ▼	Anti-Trendelenburg / Trendelenburg ▼	Max. user weight ▼	Total weight ▼	Heaviest part of the product ▼	Washable in washing machine ▼
<b>NordBed Optimo 90</b>	0 - 19°	+/- 15°	200 kg	100 kg	24 kg	Yes
<b>NordBed Optimo 105</b>	0 - 19°	+/- 15°	200 kg	110 kg	26 kg	Yes
<b>NordBed Optimo 120</b>	0 - 19°	+/- 15°	250 kg	130 kg	70 kg	Yes

**NB:** The bed must not be used by patients under 12 years of age, or by patients with body size equivalent to an average 12 year old or smaller.

## Frame colour ▼ Woodwork ▼

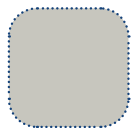

Aluminium grey



Beech



Oak



Umber



Beech/  
pale blue



Pale blue



White

Tel: +41 61 487 70 70  
hqeurope@invacare.com  
www.invacare.eu.com

Invacare International GmbH  
Neuhofweg 51  
4147 Aesch BL  
Switzerland



Yes, you can.®



© 2023 Invacare International GmbH. All rights reserved. All information quoted is believed to be correct at time of print. Invacare reserves the right to alter product specifications without prior consultation. EU - 04/23