

NEW

Invacare®

**Action<sup>®</sup>4**  
**NG**

Stylish and stable with outstanding  
functionality



The new and improved **Action<sup>4</sup><sub>NG</sub>** has been redesigned with a weight reduction of nearly 1 kg and maintains its outstanding functionality. Featuring a double cross brace system, the **Action<sup>4</sup><sub>NG</sub>** now offers a more stable, easier drive. Enhanced with lightweight options and accessories to reduce energy consumption and aids effortless propulsion, this new quality design is echoed throughout the range with its high cross-compatibility.

- **Much lower, overall product weight**
- **Improved functionality and new options**
- **Updated design with an increased colour range**



### **Significant weight reduction**

Without compromising on the functionality - or limiting the flexibility of adjusting the wheelchair - the weight of the **Action<sup>4</sup><sub>NG</sub>** has been reduced by nearly 1 kg. This is, in part, thanks to a new lightweight castor housing design and an optimised rearwheel attachment which are both made out of aluminum. This weight reduction makes the **Action<sup>4</sup><sub>NG</sub>** easier to propel and lift in and out of the car for transportation.

### **Attention to detail**

A new fork design provides a lighter look, the backrest upholstery has improved style and the new, fresh frame colours give the **Action<sup>4</sup><sub>NG</sub>** a sleek new look to complement individual personalities.

### **Improved functionality**

New, cleverly designed features make the **Action<sup>4</sup><sub>NG</sub>** even easier to use. An added lever to the swing away armrest makes it easier to release and transfer out of the wheelchair, whilst at the same time the sideguard design has been optimised with a cut-out section and foam grip to allow a smoother transfer (on the height adjustable version).

Optimum comfort can now be achieved with the new intuitive stepless angle adjustable backrest and newly designed bi-material shaped cover.

Therefore creating the perfect fit backrest angle (-15° - +15°) for the individual and ensuring all round comfort and support from the oversized backrest cover, together with a lumbar support to position where required.

### **Easy to set-up and change**

The **Action<sup>4</sup><sub>NG</sub>** remains fully cross compatible with the other models of the Action Family, meaning all new components of the wheelchair fit to previous models.



## Features and options



### Lightweight

The **Action<sup>4</sup><sub>NG</sub>** has a basic weight from 14.7 kg and a transportation weight of 8 kg (without legrests, armrests or rear wheels). With the lightweight wheels this can be reduced to 14.0 kg.



### Rearwheel support

The new aluminum rearwheel support, combines lightness with functionality.



### Castor housing and fork

Aluminum castor housing, combined new fork design, combines weight reduction with style.



### Stepless, adjustable, folding backrest

This innovative, patented, angle adjustable backrest - incorporated from the *Invacare Action<sup>5</sup>* - allows fine adjusting of the backrest angle when the individual is seated in the wheelchair. This unique folding mechanism enables the backrest to fold on top of the seat for easy storage or transportation in a car. In addition, the built-in scale allows for a more accurate backrest angle for the perfect fit.



### Flip-back Armrest No.2 and No.3

Both armrest options now feature a release lever and an open zone - with handgrip - for easy transfer out of the chair. In addition, the No.3 armrest is also fully height adjustable.



### Amputee support

Offering height, angle and depth adjustments with a durable, moisture proof design, this swing-away support is ideal for transfers.

### Front castor

Front castors available in 5", 6" and 8" pneumatic or solid tyre, with various options of tyre width also available.

### Tie-down hook

A bright, orange tie-down hook provides a secure fixation point for transportation in a motor vehicle.



### Dual handrim propelling wheels

The Dual Handrim version is ideal for hemiplegic individuals. The telescopic, spring-loaded connections can be removed so the chair remains foldable.

### Belt buckle

The Bodypoint 2-point lapbelt can be positioned where required to offer maximum postural support.

## Features and options



### **Matrix® seating**

From basic cushions, to the top of the range *Invacare Matrix Libra* cushions, Invacare offers a large variety of seat cushions to accommodate a range of individual needs. The *Invacare Matrix Elite* and *Elite Deep* backrests come complete with the EasyMount quick release system, which allows quick, simple removal when the chair is folded. This compatibility with our extensive seating and backrest range puts comfort and support at the heart of the **Action<sup>4</sup><sub>NG</sub>**.



### **Length compensating elevating legrest**

This legrest combines optimal length compensation with 0-80° elevation.



### **Reclining backrest**

The reclining backrest, both mechanical and gas spring, allow carers to make quick and easy adjustments for individual positioning.



### **Head and neck rest options**

A choice of head and neck rest options are available - including the award winning Elan headrest - to help give individuals the maximum level of comfort and support.

## Technical data

For more comprehensive information about this product, including the product's user manual, please visit [www.invacare.co.uk](http://www.invacare.co.uk)

|                        | Seat width                                   | Seat depth            | Seat angle   | Seat height                                      | Backrest angle                       | Backrest height | Front frame angle |
|------------------------|--|-----------------------|--------------|--|--------------------------------------|-----------------|-------------------|
| <b>Self Propulsion</b> | 380/405/430/<br>455/480/505/<br>555/605 mm   | 400 / 450 /<br>500 mm | 0° / 3° / 6° | 360 / 385 / 410 /<br>435 / 460 / 485<br>/ 510 mm | 0° / 7° / 0° - 30° /<br>- 15° - +15° | 325 - 510 mm    | 80° / 0° - 80°    |
| <b>Transit</b>         | 380/405/430/<br>455/480/505/<br>555 / 605 mm | 400 / 450 /<br>500 mm | 0° / 3° / 6° | 460 / 485 / 510<br>mm                            | 0° / 7° / 0° - 30° /<br>- 15° - +15° | 380 - 460 mm    | 80° / 0° - 80°    |

|                        | Length of footrests                        | Armrest height                                   | Driving unit width     | Total width, folded | Total length inc. Footrests | Total length without footrests | Total height  |
|------------------------|--|--|------------------------|---------------------|-----------------------------|--------------------------------|---------------|
| <b>Self Propulsion</b> | 200 - 400 /<br>330 - 460 /<br>440 - 510 mm | 200 - 280 /<br>240 / 170 - 300 /<br>210 - 350 mm | Seat width<br>+ 190 mm | 305 mm              | 995 - 1170 mm               | 850 - 1025 mm                  | 810 - 1040 mm |
| <b>Transit</b>         | 200 - 400 /<br>330 - 460 /<br>440 - 510 mm | 200 - 280 /<br>240 / 170 - 300 /<br>210 - 350 mm | Seat width<br>+ 160 mm | 265 mm              | 950 - 1050 mm               | 805 - 905 mm                   | 910 - 1040 mm |

|                        | Max. Safe slope | Radius | Max. User weight          | Total weight, from | Mass of the heaviest part, from |
|------------------------|-----------------|--------|---------------------------|--------------------|---------------------------------|
| <b>Self Propulsion</b> | 18°             | 850 mm | 125 / 160 kg <sup>1</sup> | 14.0 kg            | 8.0 kg                          |
| <b>Transit</b>         | 18°             | 850 mm | 125 / 160 kg <sup>1</sup> | 14.5 kg            | 11.7 kg                         |



1. Heavy duty version

## Frame colours


Hypnotic Purple

Matt Black

Blue Pop

Matt Red

Black Nylon

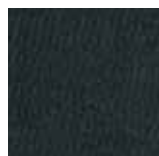


Pearl White

Water Blue

Electric Green

Anthracite Black



## Upholstery

Follow us on social media:



Please notice colours may vary slightly from those displayed above.

**Invacare Limited**  
Pencoed Technology Park - Pencoed  
Bridgend CF35 5AQ - United Kingdom  
Tel: +44 1 656 776 222  
Fax: +44 1656 776 220  
E-mail: [uk@invacare.com](mailto:uk@invacare.com)  
Sales Order E-mail: [ordersuk@invacare.com](mailto:ordersuk@invacare.com)  
[www.invacare.co.uk](http://www.invacare.co.uk)

**Invacare Ireland Ltd**  
Unit 5 - Seatown Business Campus  
Seatown Road - SWORDS  
County Dublin - Ireland  
Tel. +353 1 8107084  
Fax +353 1 8107085  
Email: [Ireland@invacare.com](mailto:Ireland@invacare.com)  
[www.invacare.ie](http://www.invacare.ie)



**Yes, you can.®**

