


# Softform<sup>®</sup> Premier

Pressure redistributing high specification mattress

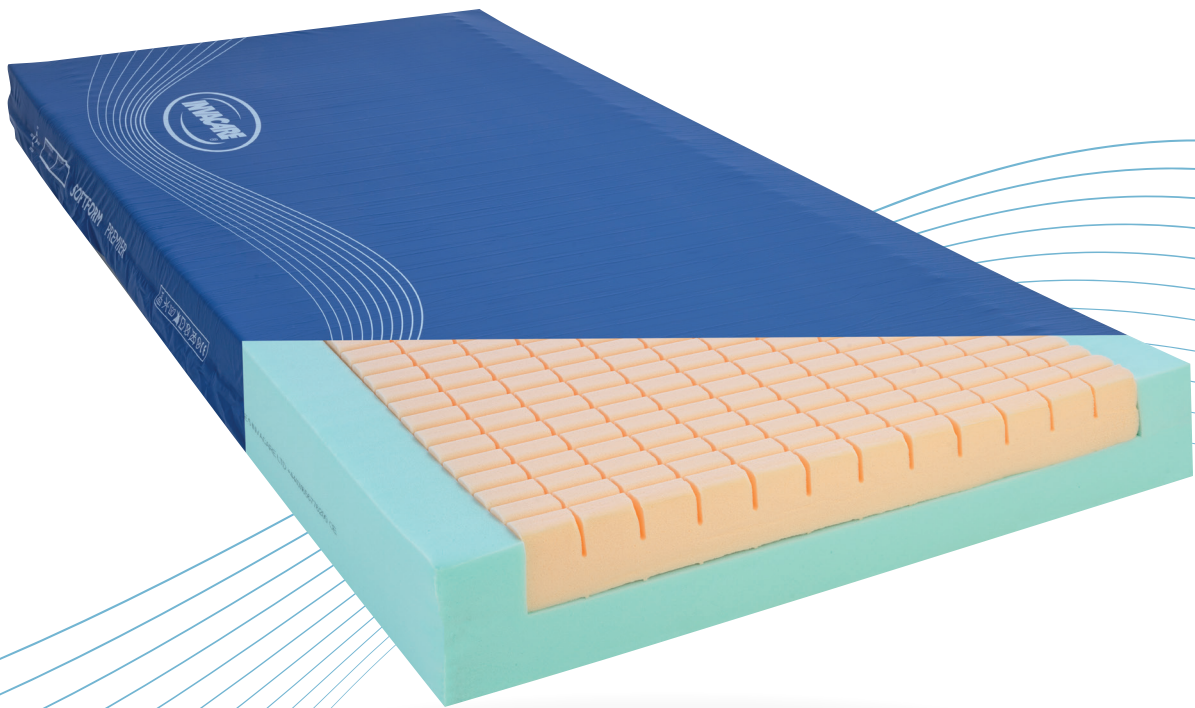

Yes, you can.<sup>®</sup>


# Effective pressure injury prevention

The Invacare Softform Premier is a high specification pressure redistributing mattress designed to meet the demands of the modern ward environment. Practical and durable, this mattress delivers exceptional levels of comfort and pressure reduction.

The Invacare Softform Premier features a unique contoured insert pad which provides excellent patient comfort and optimises pressure reduction.


The waterproof mattress has a two-way stretch that positively helps to reduce shear and friction.

The cover is a two-way stretch polyurethane material with high frequency welded seams; and a full length flap conceals the zip. These features combine to effectively reduce the potential of fluid ingress. Its toughened PU coated base also helps to prolong the longevity of the mattress.

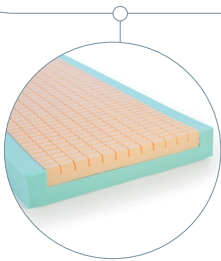

## Single piece foam "U" core

- Facilitates patient transfer by providing **effective side-wall support**.

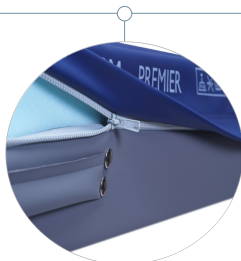


## Full length flap conceals zip

- The full length flap protects the zip. This **reduces the possible risk** of fluid ingress, preventing contamination of the foam core.



## Orientation symbol

- The orientation symbol ensures that the mattress is **positioned correctly** in the bed.

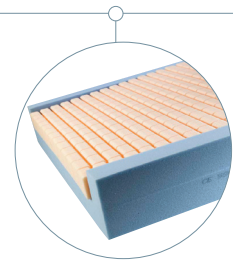

## Tapered surface in the heel zone

- The top surface of the **mattress tapers** to relieve pressure on the heels, which are particularly vulnerable to pressure injuries.

## Technical data ▼

For more information about this product, including the user manual, please visit our website: [www.invacare.co.uk](http://www.invacare.co.uk)



Width ▼



Height ▼



Length ▼



Weight ▼



User weight ▼



Clean temperature ▼



Colour

Single size:  
660 - 900 mm  
\*Double size:  
910 - 1370 mm  
\*King size:  
1380 - 1525 mm

152 mm

1410 - 2100 mm

13 kg

247.6 kg

95 °C

Manufactured to European Standards:  
BS EN 597-1 Cig,  
BS EN 597-2 Match.  
Foam Interior:  
BS 5852: 1990 Crib Source 5; Cover BS7175:  
1989 Crib Source 5.

1. Essential nursing care is pivotal in pressure injury prevention. This mattress will positively contribute to the outcome of a pressure injury prevention care plan. Education, clinical judgement and action based planning based on vulnerability are fundamental factors in the prevention of pressure injuries. A range of assessment scales can be used as a formal method of assessing risk from pressure injury development, and should be used in conjunction with an informal assessment (informed nursing judgment). Informal assessment is considered to be of greater importance and clinical value.

\*N.B. The Double and Kingsize mattresses are not available with welded seams, only sewn covers

T: +44 1656 776222  
uk@invacare.com  
www.invacare.co.uk  
www.invacare.ie

Invacare UK & Ireland  
Unit 4 Pencoed Technology Park  
Bridgend CF35 5AQ



Yes, you can.

