


Dyna-Form[®] Jupiter

Static mattress range

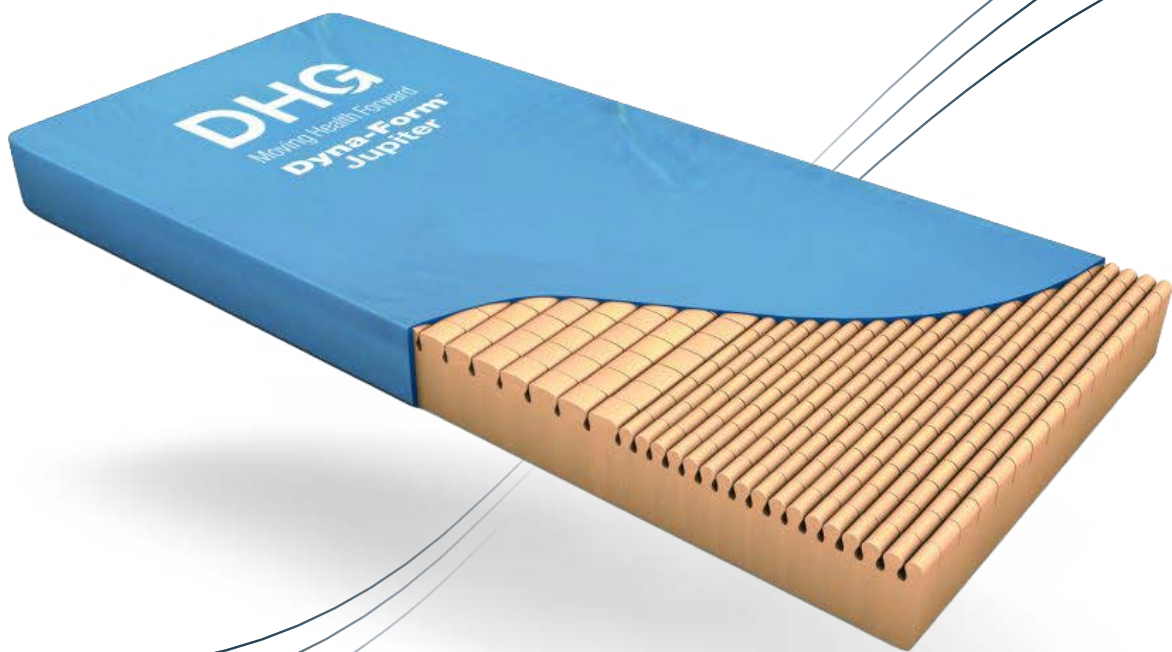


No turning required

Single-sided castellated foam



Castellated foam, no turning required

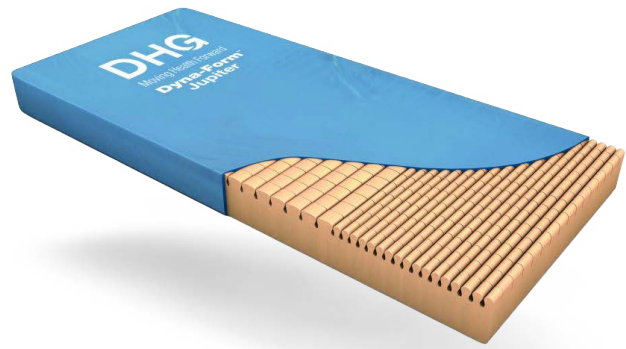

Comfort & Support

- > Castellated foam with tear duct cut: smaller cells at the heels, larger cells sacrally for targeted pressure redistribution.
- > Single-sided cut surface works without a turning programme and aids moving and handling.
- > Compatible with any electric profiling bed frame and suitable for users up to high risk of pressure ulcer development.



Multi-Directional Stretch

- > Multi-stretch, vapour permeable fabric impregnated with an anti-fungal agent reduces shear and friction forces.
- > Fully sewn or high frequency welded cover prevents fluid ingress and supports infection control and decontamination.
- > Launderable up to 95 °C.



Quality Assured

- > Finest quality CME foam and fabrics for durability and long-term value.
- > Fully replaceable parts extend product life and reduce total cost of ownership.
- > Manufactured to BS 7177:2008 (Crib 5 Composite), Very Severe fatigue class. Also available with Pro-Formance™ strikethrough resistant cover.

Technical data ▼

					
Height ▼	Width ▼	Length ▼	Max. User Weight ▼	Total weight ▼	Max. Cleaning Temperature ▼
120, 140, 150 mm	760-1200 mm	1790-2050 mm	165 kg	11 kg	95 °C

Manufactured to: BS EN 597-1 Cig, BS EN 597-2 Match and BS 7177 (Ignition Source 5)

Essential nursing care is pivotal in pressure injury management. This mattress may positively contribute to the outcome of a pressure injury prevention plan that includes education, clinical judgment and action planning based on vulnerability, which are fundamental factors in the prevention of pressure injuries. A range of assessment scales can be used as a formal method of assessing risk from pressure injury development, and should be used in conjunction with an informal assessment (informed nursing judgement) as part of a comprehensive and individualised plan.

Invacare International GmbH
Neuhofweg 51, 4147 Aesch BL, Switzerland
globalcomms@dhcare.com
www.dhcare.com

Invacare International GmbH is a member of the DHCare group of companies. Invacare International GmbH reserves the right to modify existing product models at any time. Product specifications may vary across regions. ©2026 Invacare International GmbH. All rights reserved. Dyna-Form_Jupiter_EU_20072026

