


Dyna-Form[®] Mercury[®] Advance

Powered hybrid mattress range

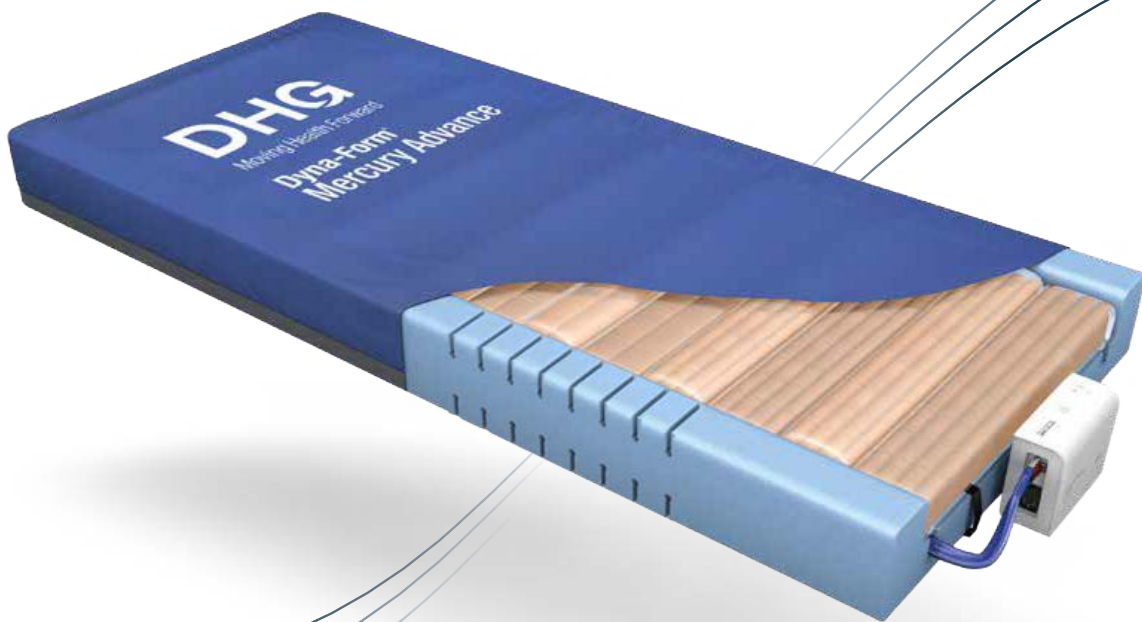

Step up / step down therapy

Static redistribution and alternating relief



Very high-risk pressure care



Hybrid pressure care

- > A 2-in-1 system combining static pressure redistribution with alternating pressure relief, able to step up to dynamic therapy when clinically required and step down as the patient improves.
- > Castellated foam construction with a CME foam base and side walls provides stability for patient transfer and prevents bottoming out.
- > Specifically designed for patients at Very High Risk of pressure ulcer development, with a weight capacity up to 254 kg.



Simple and safe operation

- > Easy-use power unit that automatically adjusts pressure to the patient's weight, with a static mode for rapid transfer and cleaning.
- > Flashing LED alert warns the carer of low cell pressure or system failure.
- > Easy-deflate CPR connector valve allows rapid deflation within 10 seconds.



Cover and sustainability

- > Multi-stretch, vapour permeable, high frequency welded cover satisfies strict infection control policies and is hypochlorite compatible.
- > Made from Orbis sustainable foam, saving around 2.6 kg of CO₂ per mattress compared with traditional polyurethane, with one post-consumer mattress recycled for every 16 produced.
- > All component parts are interchangeable and replaceable, maximising product life and reducing environmental impact.

Technical data ▼

					
Height ▼	Width ▼	Length ▼	Max. User Weight ▼	Total weight ▼	Max. Cleaning Temperature ▼
150 mm	830-1090 mm	1970, 1980, 2130 mm	254 kg	13.3 kg mattress 1.7 kg power unit	95 °C

Manufactured to: BS EN 597-1 Cig, BS EN 597-2 Match and BS 7177 (Ignition Source 5)

Essential nursing care is pivotal in pressure injury management. This mattress may positively contribute to the outcome of a pressure injury prevention plan that includes education, clinical judgment and action planning based on vulnerability, which are fundamental factors in the prevention of pressure injuries. A range of assessment scales can be used as a formal method of assessing risk from pressure injury development, and should be used in conjunction with an informal assessment (informed nursing judgement) as part of a comprehensive and individualised plan.

Invacare International GmbH
Neuhofweg 51, 4147 Aesch BL, Switzerland
globalcomms@dhcare.com
www.dhcare.com

Invacare International GmbH is a member of the DHCare group of companies. Invacare International GmbH reserves the right to modify existing product models at any time. Product specifications may vary across regions.
©2026 Invacare International GmbH. All rights reserved. Dyna-Form_Mercury Advance_EU_20072026

