


Dyna-Form[®] Mercury[®]

Static mattress

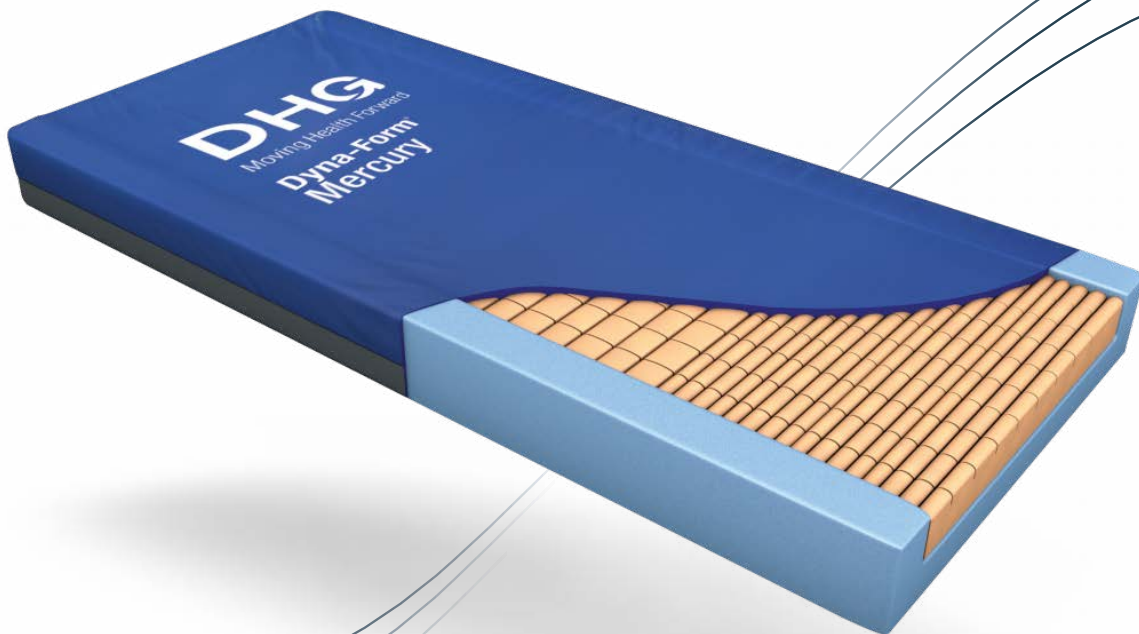


Castellated foam

Static pressure redistribution



High-risk pressure care

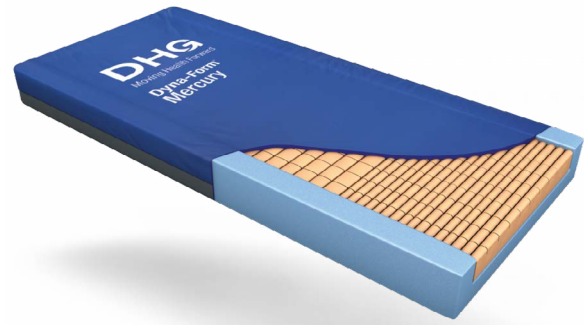

Comfort & Support

- > Castellated foam insert with tear duct cut: smaller cells at the heels, larger cells sacrally for optimal pressure redistribution.
- > CME foam base and side walls provide primary stability for patient transfer and prevent bottoming out.
- > Specifically designed for users at high risk of pressure ulcer development.



Multi-Directional Stretch

- > Multi-stretch, vapour permeable fabric impregnated with an anti-fungal agent reduces shear and friction forces.
- > Fully high frequency welded cover prevents fluid ingress and supports infection control and decontamination.
- > Launderable up to 95 °C.



Quality Assured

- > Toughened PU non-slip base provides a robust protective surface against hydrolysis and strikethrough.
- > Fully replaceable parts for long-term value and extended product life.
- > Manufactured to BS 7177:2008 (Crib 5 Composite), Very Severe fatigue class. Also available with Pro-Formance™ strikethrough resistant cover.

Technical data ▼

 Height ▼	 Width ▼	 Length ▼	 Max. User Weight ▼	 Total weight ▼	 Max. Cleaning Temperature ▼
150 mm	800-1370 mm	1980-2080 mm	254 kg	11 kg	95 °C

Manufactured to: BS EN 597-1 Cig, BS EN 597-2 Match and BS 7177 (ignition Source 5)

Essential nursing care is pivotal in pressure injury management. This mattress may positively contribute to the outcome of a pressure injury prevention plan that includes education, clinical judgment and action planning based on vulnerability, which are fundamental factors in the prevention of pressure injuries. A range of assessment scales can be used as a formal method of assessing risk from pressure injury development, and should be used in conjunction with an informal assessment (informed nursing judgement) as part of a comprehensive and individualised plan.

Invacare International GmbH
Neuhofweg 51, 4147 Aesch BL, Switzerland
globalcomms@dhcare.com
www.dhcare.com

Invacare International GmbH is a member of the DHCare group of companies. Invacare International GmbH reserves the right to modify existing product models at any time. Product specifications may vary across regions. ©2026 Invacare International GmbH. All rights reserved. Dyna-Form_Mercury_EU_20072026

