


# Dyna-Form<sup>®</sup> Mercury<sup>®</sup>

Static mattress

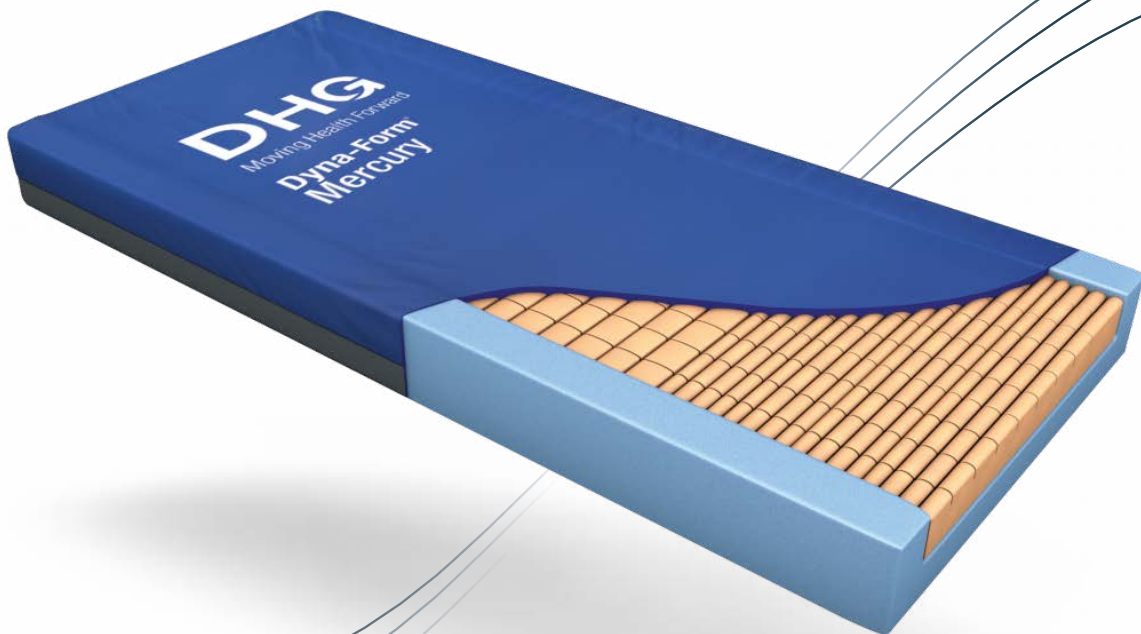


**Castellated foam**

Static pressure redistribution



# High-risk pressure care

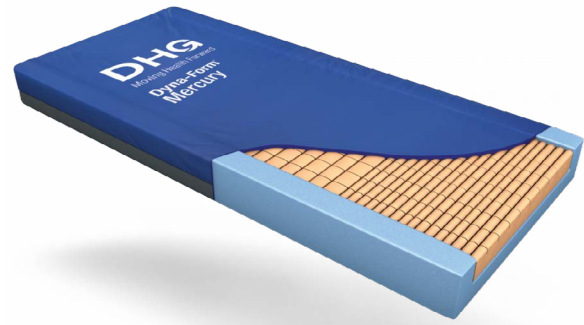

## Comfort & Support

- > Castellated foam insert with tear duct cut: smaller cells at the heels, larger cells sacrally for optimal pressure redistribution.
- > CME foam base and side walls provide primary stability for patient transfer and prevent bottoming out.
- > Specifically designed for users at high risk of pressure ulcer development.



## Multi-Directional Stretch

- > Multi-stretch, vapour permeable fabric impregnated with an anti-fungal agent reduces shear and friction forces.
- > Fully high frequency welded cover prevents fluid ingress and supports infection control and decontamination.
- > Launderable up to 80 °C.



## Quality Assured

- > Toughened PU non-slip base provides a robust protective surface against hydrolysis and strikethrough.
- > Fully replaceable parts for long-term value and extended product life.
- > Manufactured to BS 7177:2008 (Crib 5 Composite). Also available with Pro-Formance™ strikethrough resistant cover.

## Technical data ▼

| Height ▼ | Width ▼     | Length ▼     | Max. User Weight ▼ | Total weight ▼ | Max. Cleaning Temperature ▼ |
|----------|-------------|--------------|--------------------|----------------|-----------------------------|
| 150 mm   | 800-1370 mm | 1980-2080 mm | 254 kg             | 11 kg          | 80 °C                       |

Manufactured to: BS EN 597-1 Cig, BS EN 597-2 Match and BS 7177 (ignition Source 5)

Essential nursing care is pivotal in pressure injury management. This mattress may contribute to the outcome of a pressure injury prevention plan that includes education, clinical judgment and action planning based on vulnerability, which are fundamental factors in the prevention of pressure injuries.

Contact us:  
Invacare International GmbH  
Neuhofweg 51  
4147 Aesch BL, Switzerland  
globalcomms@dhcare.com  
[www.dhcare.com](http://www.dhcare.com)



Legal Manufacturer:  
DHG Ltd.  
Withey Court  
Western Industrial Estate  
Caerphilly, United Kingdom  
CF83 1BF



Invacare International GmbH is a member of the DHGcare group of companies. Invacare International GmbH reserves the right to modify existing product models at any time. Product specifications may vary across regions. ©2026 Invacare International GmbH. All rights reserved. Dyna-Form\_Mercury\_EU\_25092026