

THE OPTIMAL CUSHION FOR WHEELCHAIR USERS WITH BODY ASYMMETRY



Pre-ischial bar limits sliding forward | converts shear into normal forces

Follows the body's contours for equal pressure distribution | support IT-coccyx area | managing pelvic asymmetry



Thigh support

Mid height zip closure for easy adjustment when client is seated | cushion keeps its shape, enhancing the user's skin protection

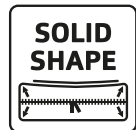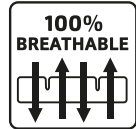

- ▲ Individual solution for asymmetry
- ▲ Minimizes sliding and shear forces
- ▲ Optimal adjustment to the body's contours for equal pressure distribution

The Vicair Adjuster 02 consists of 5 individual compartments, designed to contour under loading. Each compartment is designed to influence a specific area of the pelvis, striking a unique balance between preserving skin and underlying tissue at the ischials, coccyx and trochanters, positioning the pelvis and supporting the femurs. For specific needs the filling of each of the 5 compartments can be easily adjusted. Allowing you to change the density of each compartment benefitting users who: (1) need accommodation of a mild to severe fixed pelvic obliquity or correction of a flexible pelvic obliquity; (2) benefit from additional support of the thighs, for instance users with an above knee amputation.

AVAILABLE COVERS
Vicair 02 Top Cover, Comfair Cover, Incotec Cover and Reversible Cover.

PRODUCT PERFORMANCE

Skin protection / Anti decubitus

Pressure redistribution	▲▲▲▲▲
Shear force reduction	▲▲▲▲▲
Microclimate control	▲▲▲▲▲
Impact & vibration damping	▲▲▲▲▲

Positioning

Lateral stability	▲▲▲▲△
Frontal stability	▲▲▲▲▲
Sliding reduction	▲▲▲▲▲
Allows manual adjustments	▲▲▲▲▲
Increasing endurance	▲▲▲▲△

Maintenance & Safety

Durability	▲▲▲▲▲
Reliability	▲▲▲▲▲
Functionality when misused	▲▲▲▲△
Easy to maintain	▲▲▲▲▲
Level of hygiene	▲▲▲▲▲

Functionality & Application

Comfort	▲▲▲▲▲
Ultra light weight	▲▲▲▲▲
Allows easy transfers	▲▲▲△△
Available standard sizes	▲▲▲▲▲
Maximum user weight	▲▲▲▲△

TECHNICAL SPECIFICATIONS

Thickness	6, 10 and 12 cm
Average weight	725 grams
Max. user weight	250 kg
Machine washable	60°C max.

CUSHION SIZES

- Standard (Sizes >53 cm thickness 10 only)
- On request

