



erivo®

Sturdy & comfortable
power wheelchair



Be there:
Be free.

NEW



+ erivo. always there.

LiNX



Up and away – arrive in comfort with erivo.




Your compact companion for daily life and travel.

Being independent and self-determined, going out into the fresh air with friends, enjoying the sun in the park with the family, going on a long-awaited vacation with your partner: erivo from Alber helps keep you mobile and gives you a new quality of life. erivo folds with a simple mechanism and its compact size means it can be easily transported so you can dream of destinations again instead of worrying about the path.

Simply mobile.

erivo – powered by Alber.

As the market leader in power add-ons for wheelchairs, Alber GmbH has more than 35 years of experience in the field of electromobility. We are a recognized partner to healthcare suppliers across the world.

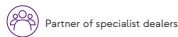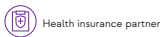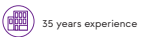

erivo.

elegant & practical

+ erivo. design



In our ever more networked world, one thing counts above all else: mobility. That's why we need means of transportation that suit our lifestyle and help us overcome obstacles. They must be reliable, flexible and user-friendly – and of course comfortable.



In top shape.

erivo offers a dynamic and contemporary drive, with its elegant shape, colour scheme and high-quality materials, it sets new standards in design gives a stylish impression.

Two-tone paint finish
in champagne gold & matte black

- ✓ **Sturdy & robust design**

+ erivo. comfort



The perfect fit.

One of the most important factors in an electric wheelchair is seat comfort. erivo boasts well-thought-out ergonomics, the best materials and top workmanship. High-quality foam ensures sitting comfort and good support at all times.

- ✓ **Ergonomic**
Comfort with firm seat & high-quality foam
- ✓ **Breathable materials**



Keep rolling.

Effective rear-wheel suspension ensures an outstanding drive comfort. Its large front wheels are puncture proof, offering good absorption and ideal support when surmounting small obstacles.

- Suspension and large pneumatic tires**
- ✓ 12 1/2" rear
- Puncture-proof PU wheels**
- ✓ 8" front



easy & flexible

+ erivo[®] function



Smart Folding Technology

With a simple mechanism, erivo folds to a compact size to save space. Hardly larger than a small suitcase, it can fit into any niche – or easily accompany you on your next trip.

- ✔ Folds quickly & easily
- ✔ Compact storage



Control unit, armrest & Footrest

erivo offers maximum flexibility with its swing-away joystick control-ler and foldable armrests. This allows you to come close to tables or change places with ease. The footrest also swivels and is height-adjustable.

- ✔ Height-adjustable footrest
- ✔ Flexible control unit & armrests



Power Unit

The powerful battery is installed under the seat and can be removed without tools. Supplied with a charger, the battery can be charged on chair. Maximum charge in only 6 hours.

- ✔ Max. charging time: 6 hours
- ✔ Range: up to (11/22 mi)*




+ erivo[®] dynamic



Intuitive Controls

The intuitive joystick offers infinitely variable and precise travel speed control. This allows erivo to navigate safely in tight spaces as well as keeping up with the grandkids in the park.

- ✔ Precise & nimble control
- ✔ Stepless speed selection (1-4 mph)

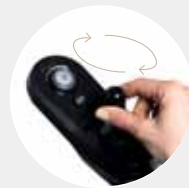

erivo - a sturdy reliable travel companion

Rediscover your favourite places.

erivo helps in everyday life. Its compact size makes it portable even in the boot of a small car. It's the ideal companion for your next trip.

Whether it's a long-planned family vacation or a spontaneous trip with friends – you and your erivo are always there for the experience. Its long-lasting battery has a range of 22 miles and is fully charged again in only 6 hours - ready for the next adventure.



© 2018 Alber Group



Be there. Be free.

Unfold experience.

erivo folds very simply and can easily be transported even in a small vehicle – and it unfolds just as easily when needed. Discover new places or enjoy your personal favourite with no worries.

During development, we placed particular emphasis on ride comfort. This is ensured by integrated suspension, and large rear wheels with pneumatic tyres. Giving a driving experience rarely seen in folding powerchairs.



alber Simply mobile.

erivo.

+ erivo® overview

Swivel arm
Allows close access to tables

Smart folding technology & compact design
Folded dimensions 25w x 34h x 21d in

Power pack
Long range up to 22 miles & rapid charging

Flip-up armrests
For easier sitting and standing up

Footrest
Swivel design & height-adjustable

Control unit
Intuitive control with LINX technology. Stepless speed control from 1 to 4 mph

Quality ride characteristics
Large, solid wheels & rear-wheel suspension

Ideal seating ergonomics
Firm, contoured & comfortable seat with breathable, high-quality materials

Innovative design
Modern design frame in champagne gold & matte black



Additional technical data

Two models available R10 & R11. The R11 is the same specification as R10 except for battery, lighting package and range*.

	R10 model	R11 model*
Shipping weight	from 60 lbs	
Battery removed:	5.3 lbs	
Standard battery pack weight:	10 lbs	
Replacement battery pack weight:	17km	35km
Range:	No	Yes
Lighting package:	11.0 Ah	30 Ah
Battery capacity		
Standard battery pack:	30.15 Ah	
Optional replacement battery pack:		
Seat width:	17.7 in	
Seat depth:	18.5 in	
Back height:	17.7 in	
Seat height:	20 in	
Total length:	38 in	
Power output:	250 watts	
Max. user weight:	23.5 st	
Protection class: (splash-protected)	IPX4	

erivo®

alber

The sturdy & comfortable
power wheelchair

NEW



10.0001.6.01.01 As at 08/2022 @Alber GmbH, Albstadt.
Subject to technical modifications.

Invacare UK & Ireland
Unit 4 Pencoed Technology Park
Bridgend CF35 5AQ

T: 01656 776222
E: uk@invacare.com
www.invacare.co.uk

Our quality standard: Your safety!

Alber is certified in accordance with the current quality management standard ISO 13485 for medical devices.

We ensure the environmentally responsible manufacture of our products with our environmental management system, certified in accordance with ISO 14001.

erivo is CE-compliant according to the European Medical Device Regulation (MDR) 2017/745 and is tested and certified on a voluntary basis by TÜV Süd, according to the current valid edition of the product standard EN 12184.

+ **erivo** always there.