

# **ESPRIT** *Action*<sup>®</sup>

DISCOVER. EXPLORE. EXPERIENCE.



**NEW** highly configurable, customisable & transportable power chair

Featuring innovative gyro technology



**Yes, you can.<sup>®</sup>**


# Advanced new design.

**Our Esprit Action® has moved forward with technical precision.**

Every part of its new design has made progressive steps.

An extensive customer-led design programme has made the difference. The new Esprit Action® is lightweight, foldable and uses integrated, high specification powered wheels on a highly customisable and configurable aluminium chassis from our popular Action® range.

The results are unprecedented. A perfect combination of a technologically advanced power chair design in a truly transportable concept, its finally time to enjoy, freedom and independence without compromise.

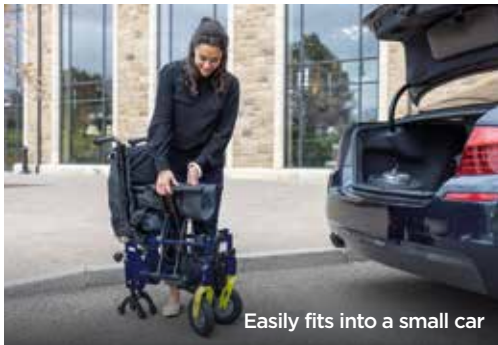

**KEY**

## FEATURES

**Negotiates narrower spaces**

Featuring new innovative gyro technology, the compact wheelbase gives enhanced manoeuvrability making the narrowest of environments easier to negotiate. The rear wheels are repositioned frontward for better traction and a smaller turning radius.

**Effortlessly handles inclines**

Inclines of up to 6° are effortless, smooth and stable. The new gyro technology reduces acceleration according to the degree of the slope and quickly detects when the wheelchair is on a tilt.



\*min. transportable weight

14" powered wheels

**Better drive performance**

14" rear wheels and pneumatic rear and front tyres give a smooth and comfortable drive.





## Easier to transport

Designed using an Action® lightweight folding frame and removable powered rear wheels and battery, the Esprit Action® can be folded for transportation without ramps or the need for adapted vehicles.



## Cross compatibility

Cross-compatibility with the Action® range allows for quick adaptations. The Esprit Action® can be customised around specific needs quickly and cost effectively.



6km/h  
speed



## Options



Lightweight Alber swivel arm



Premium anti-tipper option provides jack-up function for quick release of wheels



## Comfort & support

Optimal positioning is achieved through highly configurable seat dimensions, postural back systems as well as a vast array of options and accessories. The wheelchair can be discreetly fine-tuned to meet clinical needs.



# Discover.

A new design

# Explore.

New possibilities

# Experience.

Even better performance

**Introducing the next generation of Esprit Action®.** An exciting new design using an innovative combination of technology, performance and carefully considered refinements.





# ESPRIT Action®















## Junior

The Esprit Action® Junior enables children to enjoy an active and independent lifestyle. The new Junior encourages kids to be kids, be happy, active and ooze confidence.

Whether it's connecting with friends online indoors or meeting up and playing much more easily due to its quick and simple transportability, every child can actively participate and feel good about their Esprit Action®. Sharper visual aesthetics and a choice of three colour combinations, the Esprit Action® ensures children use their wheelchair without limitations. More possibilities, more fun.

## Technical data

For more information about this product, including the user manual, please visit our website: [www.invacare.co.uk](http://www.invacare.co.uk)

	Seat width	Seat depth	Seat height	Backrest height	Length of footrests	Seat angle	Backrest angle
<b>Esprit Action®</b>	 380 - 605 mm	 400 - 500 mm	 460 or 485 mm	 390 - 510 mm	 340 - 470 mm	 3°	 0° - 30°
<b>Esprit Action® Junior</b>	305 - 380 mm	325 - 375 mm	460 or 485 mm	390 - 510 mm	200 - 330 mm	3°	0° - 30°
	Driving unit width	Total height	Total length inc. footrests	Total length without footrests	Total weight	Max. user weight	Battery capacity
<b>Esprit Action®</b>	 Seat width + 180 mm	 945 - 1000 mm	 980 - 1115 mm	 640 - 970 mm	 40.4 - 45 kg <sup>1</sup>	 125 kg	 36 V (3 x 12V)/ 7.65 Ah
<b>Esprit Action® Junior</b>	Seat width + 180 mm	945 - 1000 mm	925 - 1075 mm	640 - 970 mm	40.4 - 45 kg <sup>1</sup>	80 kg	36 V (3 x 12V)/ 7.65 Ah
	Turning diameter	Max. climbable obstacle height	Max. safe slope	Motor capacity	Driving range	Speed	Adjustable armrests
<b>Esprit Action®</b>	 1430 - 1960 mm	 30 mm	 6°	 2 x 110 W	 15km	 6 km/h	 +50 mm
<b>Esprit Action® Junior</b>	1430 - 1960 mm	30 mm	6°	2 x 110 W	15km	6 km/h	+50 mm

<sup>1</sup>Weight of heaviest part is the frame at 15kg

### Frame colours

#### Adult



Slate Grey



Blue

#### Adult and Junior



Anthracite Black



Dark Purple

#### Junior



Water Blue

### Colour accents



Lime



Black

Tel: +44 1 656 776 222  
E-mail: [uk@invacare.com](mailto:uk@invacare.com)  
[www.invacare.co.uk](http://www.invacare.co.uk)

Invacare Limited  
Pencoed Technology Park -  
Pencoed  
Bridgend CF35 5AQ - United  
Kingdom

