

matrx[®]
SEATING SERIES

Flo-tech™ Contour Lo-Back Cushion



Allows user to sit back further into
the seat



The Invacare **Matrix Flo-tech Contour Lo-Back Cushion** has been designed with a lower rear contour allowing the user to sit further back in the seat period. The additional contouring and support under the thighs stabilizes the user further, thus reducing the damaging effects of shear.



Matrix Flo-tech Contour



Matrix Flo-tech Contour Lo-Back

Features and Options



Ingress Resistant Cover

Cover is two-way stretch, water-resistant and vapour-permeable. Cover features an anti-slip base with carry handle.



Available in 83 sizes

The *Invacare Matrix Flo-tech Contour Lo-Back Cushion* is available in a range of 83 sizes, from 360 x 330 to 600 x 550 mm. See size chart for details.



Contoured seating

Offers improved stability and postural support. Castellations reduce sheer and friction forces on tissue



Ultra-fresh™ added to all cushions for antimicrobial and odour protection

Antimicrobial treatment that has been added to all **Flo-tech** foam to help make it resistant to the growth of odour causing bacteria. Helps provide longer lasting freshness.

Contour Lo-Back Cushion optimizes support surface contact area






Comparison with standard wheelchair cushion, showing improved posture through use of the **Matrix Flo-tech Contour Lo-Back Cushion**.

Moulded Foam Cushion

Moulding allows **Matrix Flo-Tech** cushions to provide a three-dimensional form, which cannot be emulated in cut foam. Moulded foam is highly durable, increasing the product's longevity.

Technical data


For more comprehensive pre-sales information about this product, including the product's user manual, please see your local Invacare website.

Various sizes	Height	Total product weight	Max. user weight ¹	Max. cleaning temperature
				
A range of 83 sizes	103 mm	1.38 kg	152 kg	80°

Manufactured to European Standards: BS EN 1021-1 Cig, BS EN 1021-2 Match and BS 7176 (Crib 5)

1. Guide maximum user weight limits are set assuming that the cushion selected is of appropriate width and depth for the intended user/s.

Colours

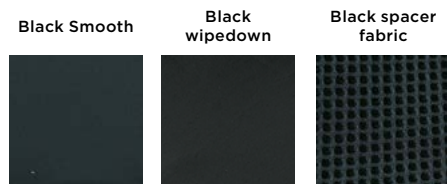

QR Code

Invacare Ltd
Pencoed Technology Park - Pencoed
Bridgend CF35 5AQ - United Kingdom
Tel +44 1 656 776 222
E-mail: uk@invacare.com
www.invacare.co.uk

Invacare Ireland Ltd
Unit 5 - Seatown Business Campus
Seatown Road - SWORDS, County Dublin - Ireland
Tel. +353 1 8107084
Email: Ireland@invacare.com
www.invacare.ie



Yes, you can.®

