


**matrx**<sup>®</sup>  
SEATING SERIES

# Flo-tech™ Lite



Slimline contoured cushion  
affording exceptional comfort and  
improved posture

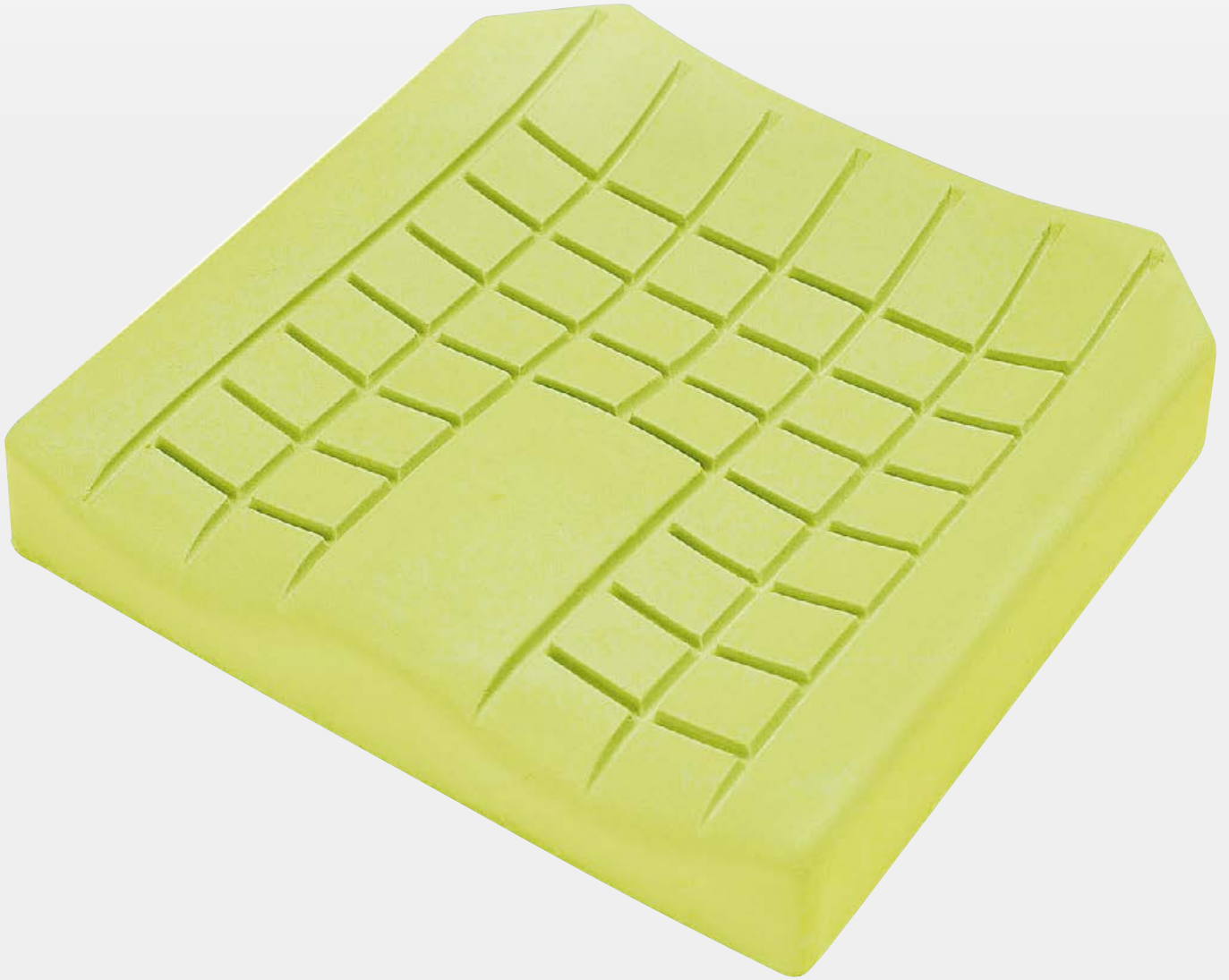

The *Invacare* **Matrix Flo-tech Lite** cushion is a slimline contoured cushion ideally suited to users at low risk of developing pressure ulcers and for comfort purposes only. The contoured seating surface offers improved stability and postural support. With sizes starting from 25 x 25 cm (250 x 250 mm) the **Matrix Flo-tech Lite** cushion is particularly suited to paediatric users.

## Features and Options



### Cover is resistant to fluid ingress

Cover is two-way stretch, water-resistant and vapour-permeable. Cover features an anti-slip base with carry handle

### Available in 61 sizes

The *Invacare* **Matrix Flo-tech Lite** is available in a range of 61 sizes, from 250 x 250 to 480 x 530 mm. (See size chart for details)

### Integral sag compensator option

Option of sag compensator to accommodate the sag in wheelchair seats, and ensure the cushion is correctly presented.

### Ultra-fresh™ added to all cushions for antimicrobial and odour protection

Antimicrobial treatment that has been added to all **Flo-tech** foam to help make it resistant to the growth of odour causing bacteria. Helps provide longer lasting freshness.

### Slimline foam base made from high quality resilient, moulded foam

The high quality resilient foam is manufactured using a unique formula which offers high levels of pressure redistribution and improves product longevity. Moulding allows the **Matrix Flo-tech Lite** cushion to provide a 3D form, which cannot be emulated in cut foam. The pommel at the front of the cushion separates the thighs, encourages good posture and allows air to circulate.

### Contoured seating surface

Offers improved stability and postural support. Also delivers a greater degree of weight distribution compared to that of a flat surface. The castellations of the cushion allow slight movement which in turn reduces shear and friction forces.

## Technical data



For more comprehensive pre-sales information about this product, including the product's user manual, please see your local *Invacare* website.

Various sizes	Height	Total product weight	Max. user weight <sup>1</sup>	Max. cleaning temperature
A range of 61 sizes	78 mm	1.5 kg	127 kg	80°

Manufactured to European Standards: BS EN 1021-1 Cig, BS EN 1021-2 Match and BS 7176 (Crib 5)

1. Guide maximum user weight limits are set assuming that the cushion selected is of appropriate width and depth for the intended user/s.

## Colours

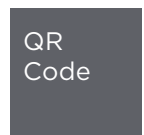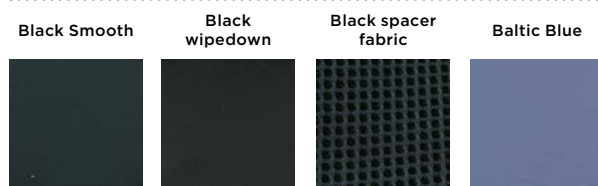

Invacare Ltd  
Pencoed Technology Park - Pencoed  
Bridgend CF35 5AQ - United Kingdom  
Tel +44 1 656 776 222  
E-mail: uk@invacare.com  
www.invacare.co.uk

Invacare Ireland Ltd  
Unit 5 - Seatown Business Campus  
Seatown Road - SWORDS, County Dublin - Ireland  
Tel. +353 1 8107084  
Email: Ireland@invacare.com  
www.invacare.ie

