


Flo-tech[®] Lite



Slimline contoured cushion
affording exceptional comfort and
improved posture

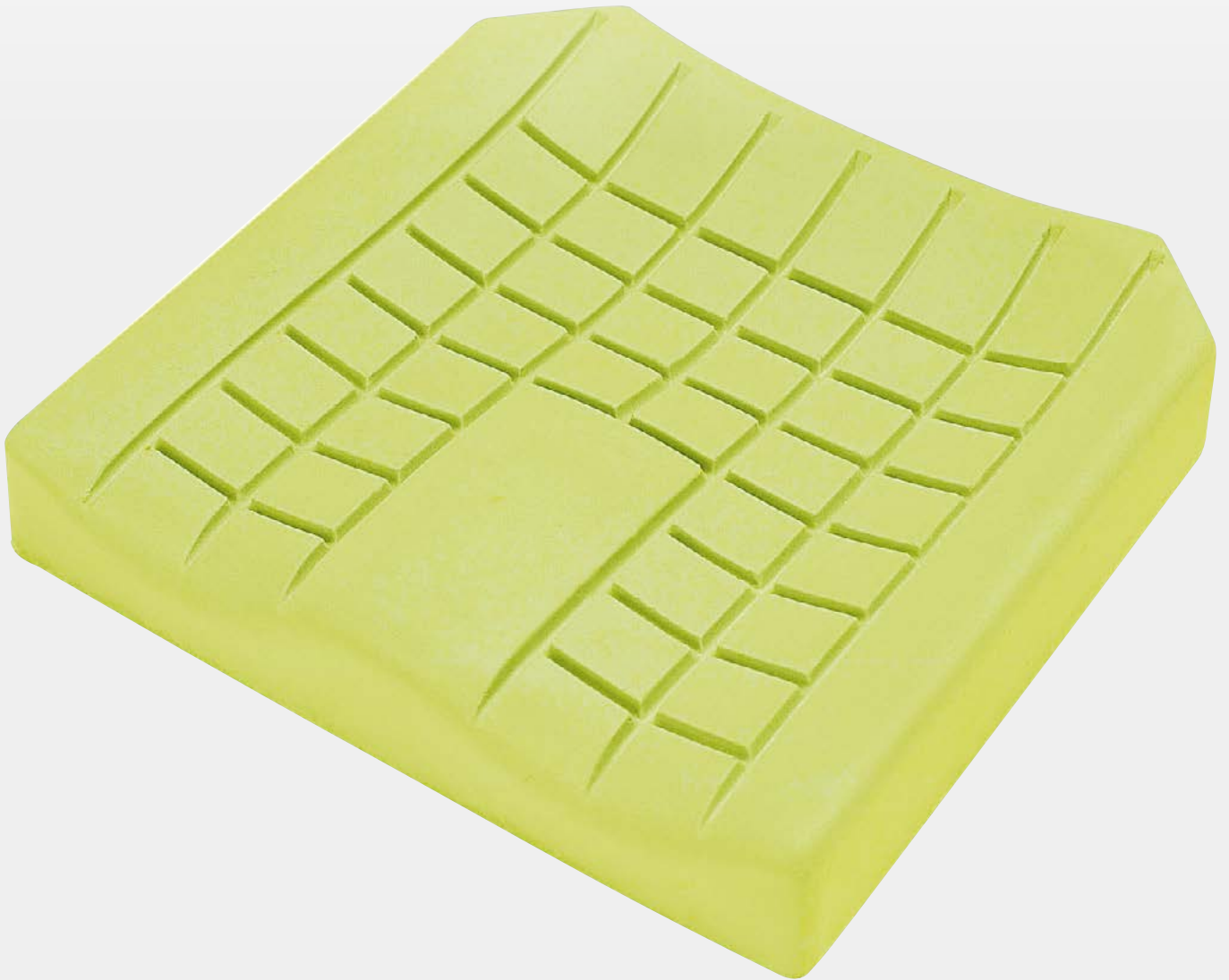

Yes, you can.

The *Invacare Flo-tech Lite* cushion is a slimline contoured cushion ideally suited to users at low risk of developing pressure ulcers and for comfort purposes only. The contoured seating surface offers improved stability and postural support. With sizes starting from 25 x 25 cm (250 x 250 mm) the **Flo-tech Lite** cushion is particularly suited to paediatric users.

Features and Options



Cover is resistant to fluid ingress

Cover is two-way stretch, water-resistant and vapour-permeable. Cover features an anti-slip base with carry handle



Available in 61 sizes

The *Invacare Flo-tech Lite* is available in a range of 61 sizes, from 250 x 250 to 480 x 530 mm. (See size chart for details)



Integral sag compensator option

Option of sag compensator to accommodate the sag in wheelchair seats, and ensure the cushion is correctly presented.



Ultra-fresh™ added to all cushions for antimicrobial and odour protection

Antimicrobial treatment that has been added to all **Flo-tech** foam to help make it resistant to the growth of odour causing bacteria. Helps provide longer lasting freshness.

Slimline foam base made from high quality resilient, moulded foam

The high quality resilient foam is manufactured using a unique formula which offers high levels of pressure redistribution and improves product longevity. Moulding allows the **Flo-tech Lite** cushion to provide a 3D form, which cannot be emulated in cut foam. Theommel at the front of the cushion separates the thighs, encourages good posture and allows air to circulate.

Contoured seating surface

Offers improved stability and postural support. Also delivers a greater degree of weight distribution compared to that of a flat surface. The castellations of the cushion allow slight movement which in turn reduces shear and friction forces.

Technical data


For more comprehensive pre-sales information about this product, including the product's user manual, please see your local Invacare website.

| Various sizes | Height | Total product weight | Max. user weight ¹ | Max. cleaning temperature |
|---|--------|----------------------|-------------------------------|---------------------------|
| | | | | |
| Flo-tech Lite A range of 61 sizes | 78 mm | 1.5 kg | 127 kg | 80° |

Manufactured to European Standards:
BS EN 1021-1 Cig, BS EN 1021-2 Match
and BS 7176 (Crib 5)

1. Guide maximum user weight limits are set assuming that the cushion selected is of appropriate width and depth for the intended user/s.

Colours

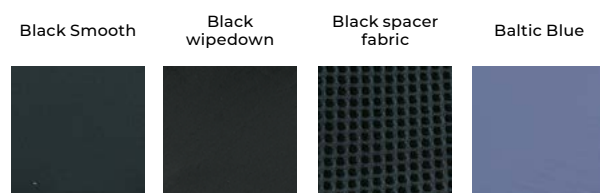

T: +44 1656 776222
uk@invacare.com
www.invacare.co.uk
www.invacare.ie

Invacare UK & Ireland
Unit 4 Pencoed
Technology Park
Bridgend CF35 5AQ



Yes, you can.®