



Küschall® **CHAMPION**

AN ORIGINAL

A BORN
CHAMPION



Küschall.com



FOLDABLE
AND SUPER
COMPACT

CHAMPION

Explore the Küschall Champion, our original and first folding, open frame designed wheelchair. Precisely engineered with a more robust hybrid folding mechanism for an advanced drive when engaged. Simple to fold and disengage in one smooth motion making for easier travel and storage.



ROBUST
DESIGN





AN ORIGINAL

Designed with passion and love,
the Kuschall Champion is like no
other. Leading in innovation and with
continuous developments to give you
a chair that has stood the test of time
for a reliable drive every time.



S U P E R **COMPACT**

Easy to handle for transporting and storage with the one step folding mechanism, folding backrest and quick release rear wheels. Upgrade to the Swiss Knife 'SK' version with a folding front frame to create a super compact Champion.



Küschall Champion SK





A B O R N ***CHAMPION***

Focused on functionality and performance,
the Kuschall champion is the perfect combination.
Advanced hydroformed technology and an
innovative folding mechanism, gives a more
rigid frame for an energy efficient drive.



K



INTRODUCING **CHAMPION CARBON**

Light when it matters. Strong when it counts. The Champion in carbon enhances portability without compromising performance. The increased resistance improves responsiveness, while the extra lightweight frame supports easier lifting, folding and transport. Look and feel confident with the beauty of carbon.



Watch
the video

A stand-out champion in
performance and aesthetics.

C



-430g*

*Lighter in total weight than the
standard aluminium Champion.



Küschall®
CHAMPION



TAKE A CLOSER LOOK

CARBOLIFE HANDRIMS

The innovative handrim profile helps improve control. The low-effort grip reduces the stresses during propelling.

FOLDING UNIT


Effortlessly folds in one movement with built in centre of gravity adjustment for tailored positioning.

CHAMPION SL

Take a look at the super light, responsive adaptation of the Champion:



**WATCH
THE VIDEO**

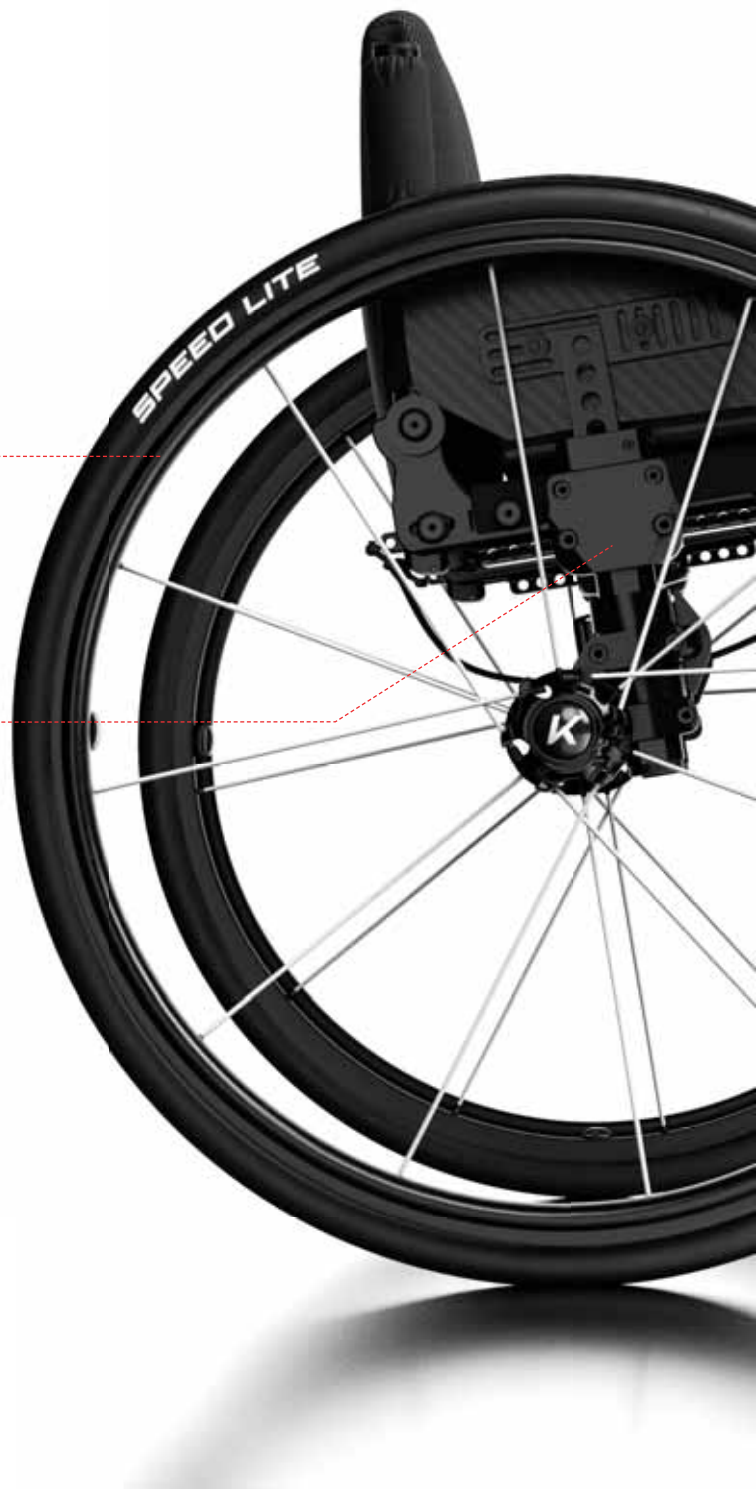



LIGHT ACTIVE BRAKE

Minimalist design built with precision even in the smallest details.

HYDROFORMED TECHNOLOGY

Strengthens the frame for increased rigidity and improved energy efficiency when propelling.

Küschall® |


LIGHTWEIGHT CARBON FRAME OPTION

Improves drive performance and gives a confident visual aesthetic.

CUSTOMISE YOUR CHAMPION

1. FRAME COLOUR

Standard frame colours

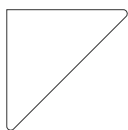

Pure White



Anthracite Grey



Space Black



Carmine Red



Azure Blue

Extra frame colours*



Flashy Yellow



Black Glossy



Silver Grey



Pearl Night Blue



Ruby Red



Dynamite Orange



Forest Green

*Extra delivery time will apply.

2. DECAL STYLE



3. DECAL COLOUR



Red



Black



Light Grey



Gold



Red



Blue



Dark Grey



Black



Pictured: Spinerger wheel with UNIQUE features. Scan the QR code on this page for more information.

4. UNIQUE SERVICE

Create your new Küschall wheelchair from the wheels up with UNIQUE. This dedicated service will create your personalized chair, designed to fit you and your lifestyle.

UNIQUE
Designed by you. Built by us.





TECHNICAL SPECIFICATION

SEAT WIDTH 360 - 480 mm In increments of 20 mm	SEAT DEPTH 340 - 480 mm In increments of 20 mm	KNEE-TO-HEEL LENGTH 320 - 500 mm Stepless adjustable	BACKREST HEIGHT 300 - 465 mm In increments of 15 mm	BACKREST ANGLE 76° / 80.5° / 85° / 89.5° / 94°
SEAT HEIGHT Front: 450 - 540 mm Rear: 390 - 490 mm	FRAME ANGLE 75° / 90° Straight or tapered	CENTRE OF GRAVITY 9 positions	CAMBER 0° / 1° / 3°	TOTAL LENGTH 75° approx. 930 mm 90° approx. 830 mm
TOTAL WIDTH SW + 160 - 240 mm	TOTAL WIDTH FOLDED Approx. 280 - 340 mm	TOTAL WEIGHT UNLOADED Approx. 8.8 kg	TRANSPORT WEIGHT Approx. 6.3 kg	MAX. USER WEIGHT 120 kg

All measurements and weights indicated are based on seat width 400 mm and seat depth 400 mm in the lightest and most compact configuration. Weight and measurements can change with the different configurations. For further information please check the user manual on our website.

/ UNLIMIT YOUR WORLD

T: +44 1656 776222
uk@invacare.com
www.invacare.co.uk
www.invacare.ie

Invacare UK & Ireland
Unit 4 Pencoed
Technology Park
Bridgend CF35 5AQ



Küschall is owned by and is a registered trademark of Invacare International GmbH. Invacare International GmbH reserves the right to modify existing models at any time. Specifications may vary across regions.
©2026 Invacare International GmbH. All rights reserved. Küschall Champion_UK_23062026