

Küschall®

K-SERIES

K

CUSTOM
POSITIONING
FOR AN
OPTIMAL
DRIVE



Küschall.com

CUSTOM
ADJUSTMENT




K-SERIES

The Küschall K-Series is an iconic modular design.

A leading trendsetter in active wheelchairs. The ingenious versatile seat module has a set-up that gives the best driving performance, wherever you go.

RIGID AND
ENERGY
EFFICIENT



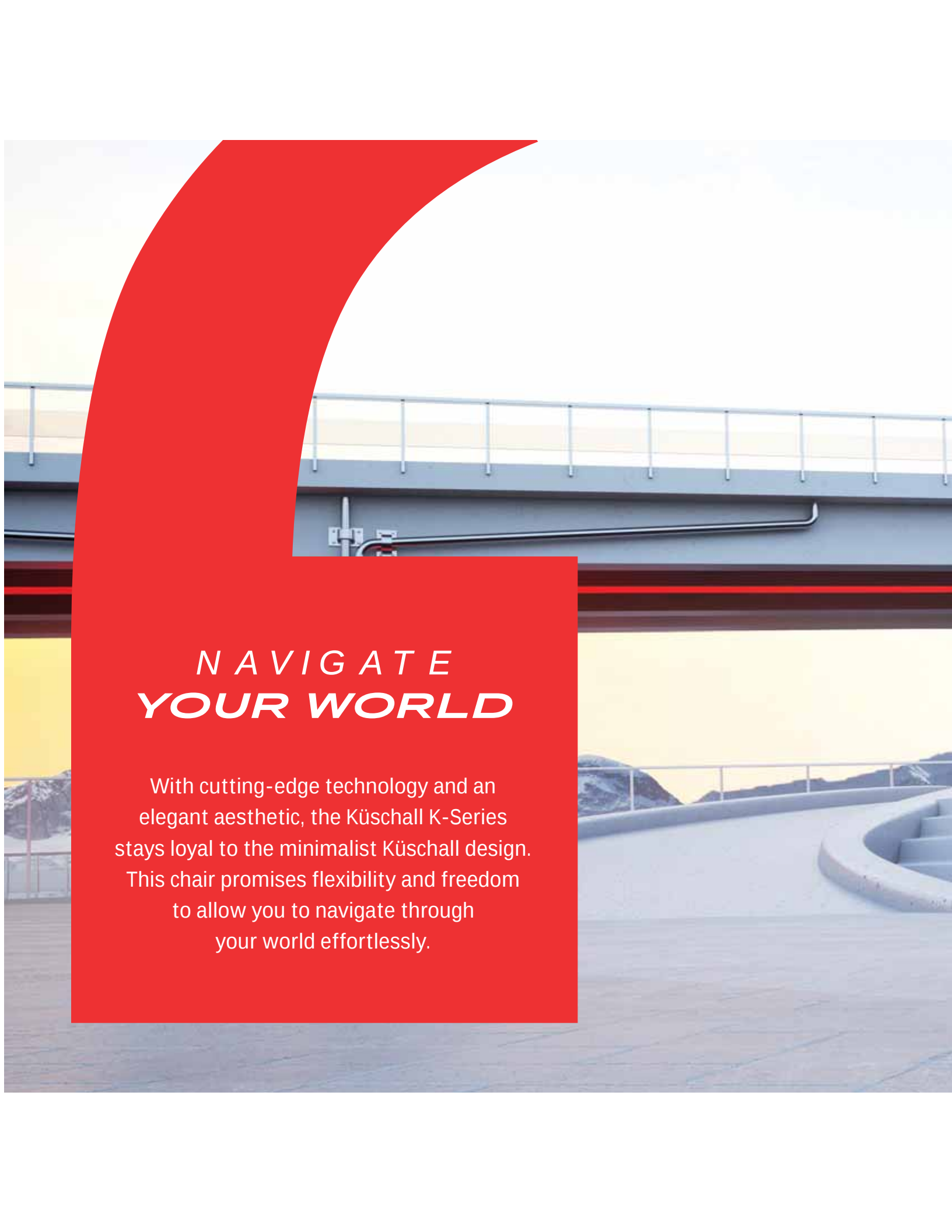

NAVIGATE YOUR WORLD

With cutting-edge technology and an elegant aesthetic, the Kuschall K-Series stays loyal to the minimalist Kuschall design. This chair promises flexibility and freedom to allow you to navigate through your world effortlessly.







MULTI-POSITIONAL S E A T

At the heart of the K-Series is a modular seat system. With five different centre of gravity positions and many seat heights, there are endless positioning possibilities – independently adjustable from the frame. This is a super ingenious design that effortlessly adapts to changing needs.



S M O O T H ***CONSISTENT*** ***DRIVE***

The K-Series has a perfect geometry to the ground. The independent seat module isolates the frame to ensure the wheels remain in the same position at all times, resulting in a smooth and consistent drive performance.



K



INTRODUCING ***K-SERIES CARBON***


Everyday performance, elevated. The K-Series in carbon brings greater efficiency to an already iconic design. With improved resistance and reduced energy loss, the chair works harder – giving a smoother, more consistent drive. The same adaptability you trust, now with a more refined and confident carbon finish.



Watch
the video

—
Enjoy the same minimalist and
modular design elevated in carbon.

C

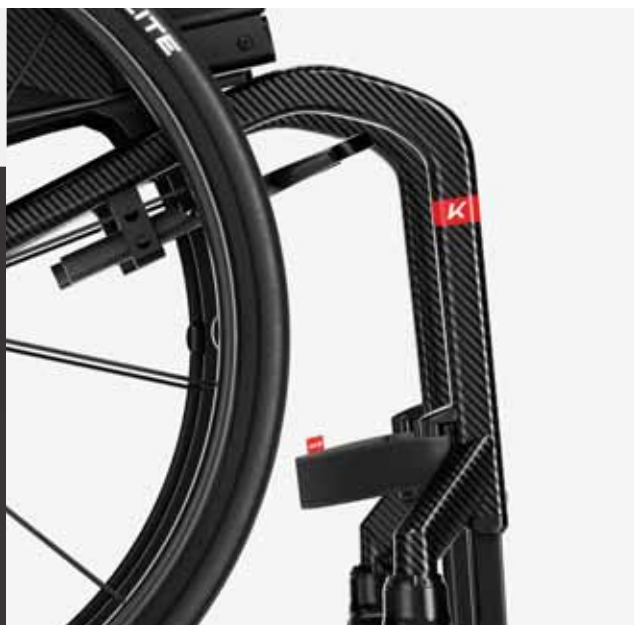

-290g*

*Lighter in total weight than the
standard aluminium K-Series.



Küschall®
K-SERIES



TAKE A CLOSER LOOK

BACKBRACE

Ergonomic handle, a crucial aspect in the chair lifting and unloading phases. Lower brace positioning, ideal for Matrx® backrests.

BACKREST

Folding backrest with six possible angle adjustments for different positioning needs.

SEAT MODULE

At the heart of the K-series. The seat module offers five different positions of the centre of gravity, making the wheelchair more active.

Küschall® | 

LIGHTWEIGHT CARBON FRAME OPTION

Improves drive performance and gives a confident visual aesthetic.





***HYDROFORMED
KNEE AREA***

Ergonomic grip, ideal for repositioning and transfers, frame angle available in both 90° and 75°.

***LIGHT
ACTIVE BRAKE***

The aluminum brake, with a minimalist design, brings lightness to unprecedented levels. Precision even in the smallest details.

CASTOR HOUSING

Aluminum profile welded to the front frame tube. Consistent design throughout the whole wheelchair.

TAPERED FRAME

Tapered frame allows the wheelchair to contour to your body for a sleek aesthetic.


FOOTBOW

Designed with a flat surface for greater comfort while driving and perfect positioning.

CUSTOMISE YOUR K-SERIES

1. FRAME COLOUR

Standard frame colours

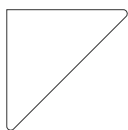

Pure White



Anthracite Grey



Space Black



Carmine Red



Azure Blue

Extra frame colours*



Flashy Yellow



Black Glossy



Silver Grey



Pearl Night Blue



Ruby Red



Dynamite Orange



Forest Green

*Extra delivery time will apply.

2. DECAL STYLE



3. DECAL COLOUR



Red



Black



Light Grey



Gold



Red



Blue



Dark Grey



Black



4. UNIQUE SERVICE

Create your new Küschall wheelchair from the wheels up with UNIQUE. This dedicated service will create your personalized chair, designed to fit you and your lifestyle.

UNIQUE
Designed by you. Built by us.





TECHNICAL SPECIFICATION

**SEAT
WIDTH**
320 – 500 mm
in increments of 20 mm



**SEAT
DEPTH**
350 – 525 mm
in increments of 25 mm



**KNEE-TO-HEEL
LENGTH**
290 – 480 mm
in increments of 10 mm



**BACKREST
HEIGHT**
270 – 480 mm
in increments of 15 mm



**BACKREST
ANGLE**
74° / 78° / 82° /
86° / 90° / 94°



**SEAT
HEIGHT**
Front: 450 – 520 mm
Rear: 390 – 490 mm
in increments of 10 mm



**FRAME
ANGLE**
75° / 90°
V-frame or tapered



**TOTAL
WIDTH**
SW + 170 mm



**TOTAL
LENGTH**
75° approx. 830 - 1055 mm
90° approx. 760 - 985 mm



**CENTER
OF GRAVITY**
80mm / 95mm / 110mm /
125mm / 140mm



CAMBER
0° / 3° / 6°



**TOTAL WEIGHT
UNLOADED**
Approx. 7.5 kg



**TRANSPORT
WEIGHT**
Approx. 4.8 kg



**MAX. USER
WEIGHT**
130 kg



CRASH TEST
Approved – Fulfills the
requirements of
ISO 7176-19



All measurements and weights indicated are based on seat width 400 mm and seat depth 400 mm in the lightest and most compact configuration. Weight and measurements can change with the different configurations. For further information please check the user manual on our website.

/ UNLIMIT YOUR WORLD

T: +44 1656 776222
uk@invacare.com
www.invacare.co.uk
www.invacare.ie

Invacare UK & Ireland
Unit 4 Pencoed
Technology Park
Bridgend CF35 5AQ

