



# Küschall® CHAMPION

AN ORIGINAL

—  
A BORN  
CHAMPION



Küschall.com



*FOLDABLE  
AND SUPER  
COMPACT*

*ROBUST  
DESIGN*

---

# *CHAMPION*

Explore the Küschall Champion, our original and first folding, open frame designed wheelchair. Precisely engineered with a more robust hybrid folding mechanism for an advanced drive when engaged. Simple to fold and disengage in one smooth motion making for easier travel and storage.





## ***AN ORIGINAL***

Designed with passion and love,  
the Küschall Champion is like no  
other. Leading in innovation and with  
continuous developments to give you  
a chair that has stood the test of time  
for a reliable drive every time.



## *SUPER COMPACT*

Easy to handle for transporting and storage with the one step folding mechanism, folding backrest and quick release rear wheels. Upgrade to the Swiss Knife 'SK' version with a folding front frame to create a super compact Champion.



Küschall Champion SK





# *A BORN CHAMPION*

Focused on functionality and performance,  
the Küschall champion is the perfect combination.  
Advanced hydroformed technology and an  
innovative folding mechanism, gives a more  
rigid frame for an energy efficient drive.



# INTRODUCING CHAMPION CARBON

**Light when it matters. Strong when it counts.** The Champion in carbon enhances portability without compromising performance. The increased resistance improves responsiveness, while the extra lightweight frame supports easier lifting, folding and transport. Look and feel confident with the beauty of carbon.



Watch  
the video

A stand-out champion in  
performance and aesthetics.

# C



-430g\*

\*Lighter in total weight than the  
standard aluminium Champion.



Küschall® | C  
CHAMPION



# TAKE A CLOSER LOOK

## **CARBOLIFE HANDRIMS**

The innovative handrim profile helps improve control. The low-effort grip reduces the stresses during propelling.

## **FOLDING UNIT**

Effortlessly folds in one movement with built in centre of gravity adjustment for tailored positioning.

## **CHAMPION SL**

Take a look at the super light, responsive adaptation of the Champion:



**WATCH  
THE VIDEO**



## **LIGHT ACTIVE BRAKE**

Minimalist design built with precision even in the smallest details.

## **HYDROFORMED TECHNOLOGY**

Strengthens the frame for increased rigidity and improved energy efficiency when propelling.

## **Küschall® |**

## **LIGHTWEIGHT CARBON FRAME OPTION**

Improves drive performance and gives a confident visual aesthetic.



# CUSTOMISE YOUR CHAMPION

## 1. FRAME COLOUR

Standard frame colours



Pure White    Anthracite Grey    Space Black



Carmine Red    Azure Blue

Extra frame colours\*



Flashy Yellow    Black Glossy    Silver Grey    Pearl Night Blue



Ruby Red    Dynamite Orange    Forest Green

\*Extra delivery time will apply.


## 2. DECAL STYLE



## 3. DECAL COLOUR



Red    Black

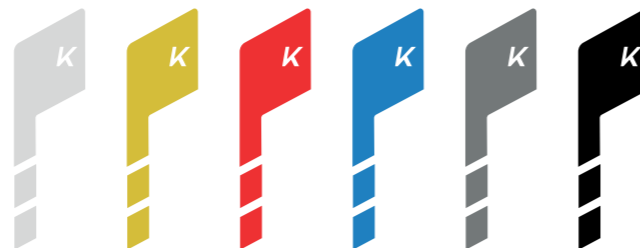

Light Grey    Gold    Red    Blue    Dark Grey    Black

## 4. UNIQUE SERVICE

Create your new Küschall wheelchair from the wheels up with UNIQUE. This dedicated service will create your personalized chair, designed to fit you and your lifestyle.

**UNIQUE**  
Designed by you. Built by us.



Pictured: Spinegy wheel with UNIQUE features. Scan the QR code on this page for more information.





# TECHNICAL SPECIFICATION

<b>SEAT WIDTH</b> 360 - 480 mm In increments of 20 mm	<b>SEAT DEPTH</b> 340 - 480 mm In increments of 20 mm	<b>KNEE-TO-HEEL LENGTH</b> 320 - 500 mm Stepless adjustable	<b>BACKREST HEIGHT</b> 300 - 465 mm In increments of 15 mm	<b>BACKREST ANGLE</b> 76° / 80.5° / 85° / 89.5° / 94°
<b>SEAT HEIGHT</b> Front: 450 - 540 mm Rear: 390 - 490 mm	<b>FRAME ANGLE</b> 75° / 90° Straight or tapered	<b>CENTRE OF GRAVITY</b> 9 positions	<b>CAMBER</b> 0° / 1° / 3°	<b>TOTAL LENGTH</b> 75° approx. 930 mm 90° approx. 830 mm
<b>TOTAL WIDTH</b> SW + 160 - 240 mm	<b>TOTAL WIDTH FOLDED</b> Approx. 280 - 340 mm	<b>TOTAL WEIGHT UNLOADED</b> Approx. 8.8 kg	<b>TRANSPORT WEIGHT</b> Approx. 6.3 kg	<b>MAX. USER WEIGHT</b> 120 kg

All measurements and weights indicated are based on seat width 400 mm and seat depth 400 mm in the lightest and most compact configuration. Weight and measurements can change with the different configurations. For further information please check the user manual on our website.

**/ UNLIMIT YOUR WORLD**

info@kuschall.com  
www.kuschall.com

Invacare International GmbH  
Neuhofweg 51  
4147 Aesch BL  
Switzerland

