



Küschall® COMPACT

VERSATILITY
& COMFORT
IN ONE CHAIR



Küschall.com



*FOLDABLE
& LIGHT*

COMPACT

Küschall Compact is the absolute standard for reliability and quality in one foldable chair.

It can easily be adapted to suit varying needs to enhance overall comfort. Available in a fixed front frame with impressive rigidity for excellent driving performance or a swing-away frame, making transfers, storage and transportation even easier.

*ENHANCED
COMFORT*





OUTSTANDING FUNCTIONALITY

A foldable backrest option makes the chair small and compact for easy transportation. Available in 70° or 80° legrest angles, in addition to the adjustable precision locking mechanism with the swing-away legrests for increased rigidity and optimal drive.



INTRODUCING COMPACT CARBON

Small detail. Noticeable difference. Optional carbon legrests bring a lighter, more refined finish to the Compact. A subtle upgrade that enhances both the feel and the overall aesthetic, without changing the reliability you depend on.



Watch
the video

Elevate your Compact and be
confident with a touch of carbon.

C

-200g*

*Lighter in total weight than the
standard aluminium Compact.



Küschall® | C
COMPACT



TAKE A CLOSER LOOK

ANGLE ADJUSTABLE BACKREST

Available as an option, the angle adjustable and foldable backrest folds neatly over the seat rails for compact storage and easy transportation.

EASY FOLDING

The unique design of the cross brace makes the Compact simple to fold. A key feature if you have reduced strength or limited hand function.

KEEP ON TRACK

Intuitive and unique castor housing adjustment, to ensure your wheelchair always performs at its most energy efficient, making propelling easier for longer.



FIXED FRAME

Increases rigidity and gives a more active aesthetic. 70° or 80° legrest angle.

SWING-AWAY LEGRESTS

Making transport, storage or transfers easier. 70° or 80° legrest angle.

Küschall® |

Optional carbon legrests for a lighter more refined finish.

FOOTPLATE

The footplate range fits all seat widths and features an anti-slip surface for increased comfort and support.



CUSTOMISE YOUR COMPACT

1. FRAME COLOUR

Standard frame colours



Pure White Anthracite Grey Space Black



Carmine Red Azure Blue

Extra frame colours*



Flashy Yellow Black Glossy Silver Grey Pearl Night Blue



Ruby Red Dynamite Orange Forest Green

*Extra delivery time will apply.


2. DECAL STYLE



3. DECAL COLOUR



Red Black

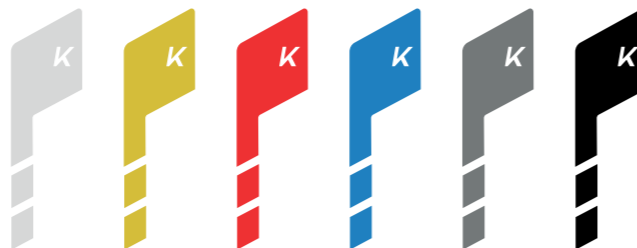

Light Grey Gold Red Blue Dark Grey Black

4. UNIQUE SERVICE

Create your new Küschall wheelchair from the wheels up with UNIQUE. This dedicated service will create your personalized chair, designed to fit you and your lifestyle.

UNIQUE
Designed by you. Built by us.





TECHNICAL SPECIFICATION

SEAT WIDTH
280 – 500 mm
in increments of 20 mm



SEAT DEPTH
280 – 500 mm
in increments of 20 mm



KNEE-TO-HEEL LENGTH
200 – 510 mm
in increments of 10 mm



BACKREST HEIGHT
300 – 510 mm
in increments of 15 mm



BACKREST ANGLE
82°/86°/90°/
94°/98°/102°



SEAT HEIGHT
Front: 380 – 530 mm
Rear: 370 – 500 mm
in increments of 10 mm



FRAME ANGLE
70° / 80°
Straight or tapered



TOTAL WIDTH
SW + 160 mm



TOTAL WIDTH FOLDED
Approx. 290 mm



TOTAL LENGTH
70° approx. 870 – 1220 mm
80° approx. 850 – 1200 mm



CAMBER
0° / 3°



TOTAL WEIGHT UNLOADED
Approx. 10.3 kg (swing away)
Approx. 9.6 kg (fixed frame)



TRANSPORT WEIGHT
Approx. 7.8 kg (swing away)
Approx. 6.7 kg (fixed frame)



MAX. USER WEIGHT
130 kg



CRASH TEST
Approved – Fulfills the requirements of ISO 7176-19



All measurements and weights indicated are based on seat width 400 mm and seat depth 400 mm in the lightest and most compact configuration. Weight and measurements can change with the different configurations. For further information please check the user manual on our website.

/ UNLIMIT YOUR WORLD

info@kuschall.com
www.kuschall.com

Invacare International GmbH
Neuhofweg 51
4147 Aesch BL
Switzerland

