

Medley[®] Ergo Low



Safe for transfer and sleep



Invacare **Medley Ergo Low**

The **Medley Ergo Low** from Invacare offers the perfect balance of economy, functionality and quality. The optimised length of the various mattress platform sections provides increased comfort for clients when sitting up in bed. It also reduces shear and friction when the bed sections are being adjusted and prevents clients sliding downwards. The **Medley Ergo Low** offers all the same benefits as the **Medley Ergo**, but boasts extra-low height levels, thus maximising safety for individuals as risk of falling out of bed. The **Medley Ergo Low** can be transformed into the *Invacare Medley Ergo Low Select*. **The Low Select** is a particularly stylish model in natural beech, creating a homely feel suited to any environment.



Optimised dimensions of the mattress platform

The mattress platform sections boast optimal dimensions designed in accordance with anthropometric measurements. The **Medley Ergo Low** offers:

- Increased comfort when sitting
- Reduced shear and friction from sliding downwards during profiling
- Suitability for shorter clients

Improved functionality and efficiency

The **Medley Ergo Low** has been specially designed for quick transportation to meet the needs of the transient working environments of Community care provision. The **Medley Ergo Low** can be supplied on a transporter kit so the bed can easily be moved and installed. In addition, any assembling can be done without the use of tools or technicians.

Features and options



Full length side rails

Available with wooden or steel full length side rails. Wooden extendable side rails (150mm) are also available.



Scala 2 range

A range of collapsible three quarter length steel side rails, available in various heights to cover standard to deep mattress heights.

Height range supports clients and carers

The flexible design of the **Medley Ergo Low** allows the bed to be mounted at a low height, or a higher position if preferred, creating a range of working heights for the carer and transfer heights for the client. This is facilitated by the unique feature of the bed end, which enables the mattress support to be mounted in both upper and lower positions.



Medley Ergo Oak



Safety in and out of bed

In the event of a client falling from the bed, a foldable mat can be placed on the floor and attached to the **Medley Ergo Low** for added safety.

Product code: 1494608

Options



Medley Ergo Low

The height of the **Medley Ergo Low** bed offers increased security for clients enabling them to sleep with real peace of mind. The **Medley Ergo Low** sits low to the floor (210mm in height) which, when combined with a floor mattress, reduces harm from falls. Carers also have the option to raise the bed to a more suitable height when necessary.

The **Medley Ergo Low** makes it easier for clients to get in and out of bed, providing more independence for those who can self transfer.

The **Medley Ergo Low** is dual height adjustable, offering a 210-610mm and 280-680mm height option to suit a wider range of needs.

Features and Options



Easy to transport

All **Medley Ergo** beds are easy to assemble and dismantle and can be transported using transporter kits (optional).



Lifting Pole

Can be adjusted in both height and depth (optional).



Select design for Medley Ergo

The homely design of **Medley Ergo Select** can also be retro-fitted (optional).



Double height positions for bed ends

Adjustable height range for **Medley Ergo Low** 210-610mm or 280-680mm.



Optimised dimensions of mattress platform

According to anthropometrical measurements, the optimised dimensions meets the body size of a wider number of clients and reduces shear and friction.

Features and options



Mattress support extension

Extension kit for increasing the length of the bed by 150mm.



Transport fitting

For easy storage and transport.



Hand controls

Available with or without lockable functions.

Technical data


For more comprehensive pre-sales information about this product, including the product's user manual, please visit www.invacare.co.uk.

	Width outside Width inside	Length outside Length inside	Height adjustment	Mattress support dimensions	Backrest angle	Thigh angle	Legrest lift
<i>Medley Ergo Low</i>	1020 mm 900 mm	2240 mm 2080 mm	280 - 680 mm or 210 - 610 mm	800 - 230 - 300 - 670 mm	0 - 70°	0 - 24°	0 - 11°
<i>Medley Ergo Low</i>	Tilt option ¹	Total weight	Max. user weight	Heaviest part of the product	Complies to the safety standard, IEC 60601-2-52.		
	11°	64 kg	SWL 180 kg 145 kg	17 kg			

1. Please note: the tilt function is only used for clients under medical supervision. If there is a risk of clients using the tilt function without supervision, please use a lockable hand control.

NB
The bed must not be used by patients under 12 years of age, or by patients with body size equivalent to an average 12 year old or smaller.

Frame colour and wood work

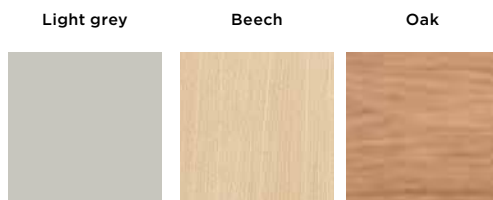

Please note colours may vary slightly from those displayed above.

Warranty

3 years frame
3 years actuator
2 years handset

The *Invacare Medley Ergo Low* is part of the Safe Patient Handling range.



safepatienthandling

Invacare Limited
Pencoed Technology Park - Pencoed
Bridgend CF35 5AQ - United Kingdom
Tel: +44 1 656 776 222
Fax: +44 1656 776 220
E-mail: uk@invacare.com
Sales Order E-mail: ordersuk@invacare.com
www.invacare.co.uk

Invacare Ireland Ltd
Unit 5 - Seatown Business Campus
Seatown Road - SWORDS - County Dublin - Ireland
Tel. +353 1 8107084
Fax +353 1 8107085
Email: Ireland@invacare.com
www.invacare.ie

