

POWER SERIES

# AVIVA RX ULM



Impressive rear wheel drive power chair  
with advanced ULM seating





| AVIVA RX ULM

# A confident and comfortable drive.

Design. Technology. Performance.

Explore the **AVIVA RX**. A high performing rear wheel drive power chair, combined with the advanced ULM seating for maximum comfort and support in any environment.





| AVIVA RX ULM

# Design. Technology. Performance.

Take a closer look at the **AVIVA RX** base. The first from the AVIVA line, the AVIVA RX sets a new standard in rear wheel drive chairs. Designed with the latest technical innovations to give a more advanced drive performance.



## CTC suspension (patented)

Designed to absorb vibrations or bumps during travel for optimum comfort and gives increased reassurance of stability with all four wheels remaining in contact with floor.



## LiNX G-Trac Gyro

Travelling long distances or at higher speeds has been made even easier with LiNX G-Trac Gyro technology. A gyro module within the electronics ensures the AVIVA RX follows control direction and stays on the right path.



## Compact design

A shorter, narrower wheelbase and low seat height makes the AVIVA RX easy to manoeuvre with precise control in all environments and even easier for use in smaller spaces with a tight turning radius.



## Sleek & stylish

Clean lines, sleek angles and a high-quality finish, sets the AVIVA RX apart when it comes to aesthetics. Modern in design to look as good as it performs and can be customised with a choice of colours.



| AVIVA RX ULM

# Designed for clinical excellence.

ULM seating has been designed by our clinical experts to offer optimal comfort and support with a personalised fit, to support greater independence and accessibility.

Built with integrated adjustment, along with the smart Ultra Rail system, ensures the ULM can be tailored to provide the correct combined seating solution.



## Extended shear reduction

Designed around the natural pivot points of the body with an optimised range of powered positioning to ensure the seating remains supportive and comfortable. Design helps minimise shear and manage pressure redistribution.



## Tailored positioning

Highly modular and adjustable to suit needs as they change over time.



## LNX legrests

Extends in both length and elevation, working with powered recline to help reduce shear.



## Award-winning Matrix

Personalised Matrix seating options to ensure an optimum fit every time.





| AVIVA RX ULM

**matrix**<sup>®</sup>  
SCIENCE OF SEATING

# Advanced seating.

Working in combination with the highly adjustable ULM seating system, the award-winning **Matrix**<sup>®</sup> seating range gives you a more supportive fit with increased protection, to ensure maximum comfort.







AVIVA RX ULM



Choose from a selection of backrests, seat cushions, headrests and a range of accessories for a personalised seating solution.



### Headrests

Our highly adjustable Matrix Elan headrest range is designed to enable effective positioning. The dual density moulded foam provides immersion, comfort and protection.



### Bodypoint positioning belts

For optimal function, mobility and stability we offer a choice of chest and positioning belts. This helps ensure both pelvis and upper body are secure for more control and improved functionality.



### Hip guides and lateral thigh/knee supports

A range of supports in various shapes and sizes. Adjustable quick release brackets for enhanced positioning, which are easy to remove for transfers.



### Backrest options

A range of backrest options including our Matrix series are available for maximum postural support, pressure relief and exceptional microclimate control. Featuring built in contoured support for increased stability.



### Seat cushions

Choose from a selection of seat cushions that support positioning and skin protection. Our advanced Matrix Libra is designed with a dual layer Flo-tech sac for enhanced shear management and can be adapted with additional positioning options.

| AVIVA RX ULM

# Featuring advanced **LiNX**<sup>®</sup> technology.

LiNX<sup>®</sup> technology helps give a personalised and confident drive experience. The innovative technology learns habits and preferences over time making it a super intuitive system.





| AVIVA RX ULM

**LINX**<sup>®</sup>

Discover more.



**New Resolve joystick mount**

- Exceptional range of achievable joystick positions for optimised comfort.
- Adjustable tension release.
- Optimised comfort and control.



**3.5" Touchscreen**

Impressive 3.5" colour touchscreen with joystick, speedometer and an odometer.



**Specialist controls**

Including compact remotes, head arrays and four-switch proximity arrays.



**Built-in Bluetooth<sup>®</sup>**

Enables connectivity to devices and environmental controls.



**Control input options**

A variety of control input options allow a personalised set up for maximum control and reassurance that the power chair will react consistently.







AVIVA RX ULM

# Technical data.

Further information about this product, including the user manual, please visit our website [www.invacare.eu.com](http://www.invacare.eu.com).

## Colour collection

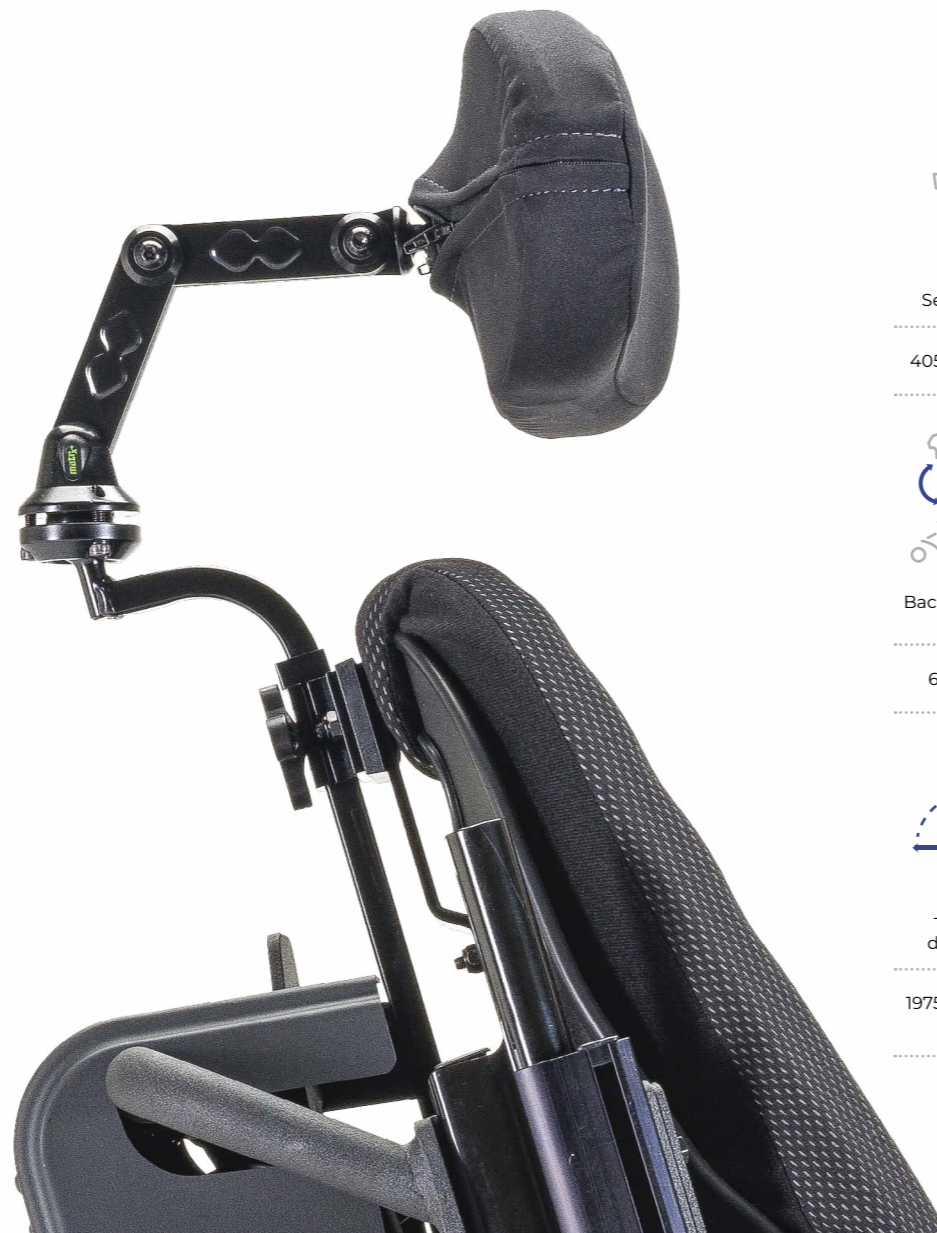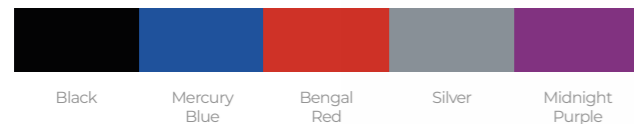

Seat width	Seat depth	Seat height	Backrest height	Armrest height	Footrests	Seat angle
405 - 555 mm	405 - 580 mm	420 - 485 mm	355 - 685 mm	230 - 405 mm	290 - 465 mm	-10° - 50°
Backrest angle	Total height	Total length inc. footrests	Total length without footrests	Driving unit width	Max. user weight	Battery capacity
60° - 170°	1101 - 1210 mm	1251 - 1300 mm	830 mm	560 / 635 mm	136 kg	60 / 73.5 Ah
Turning diameter	Max. climbable obstacle height	Max. safe slope	Speed	Driving range	Total weight	
1975 - 2105 mm	75 mm (100 mm with kerb climber)	6° / 10.5%	up to 12 km/h	up to 37 km	up to 190.1 kg	

| AVIVA RX ULM

# Share your thoughts and help shape the future of power chairs.

We invite you to **share your feedback**. Your insights will help us better serve the power chair community.

Please scan and complete the form.



Thank you for letting us join your journey.

**Invacare Power Series Team**




Share your feedback.

SUBSCRIBE



 My chair allows me to live an independent life and get outdoors everyday.


**Ian Mackay**  
Invacare Ambassador


 Joshua Evans



1 355 likes


Terry maddock Looks so peaceful!

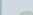 Richard Thomas

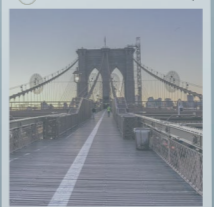

555 likes

Tyrene Stewart Wish I was there, looks nice.

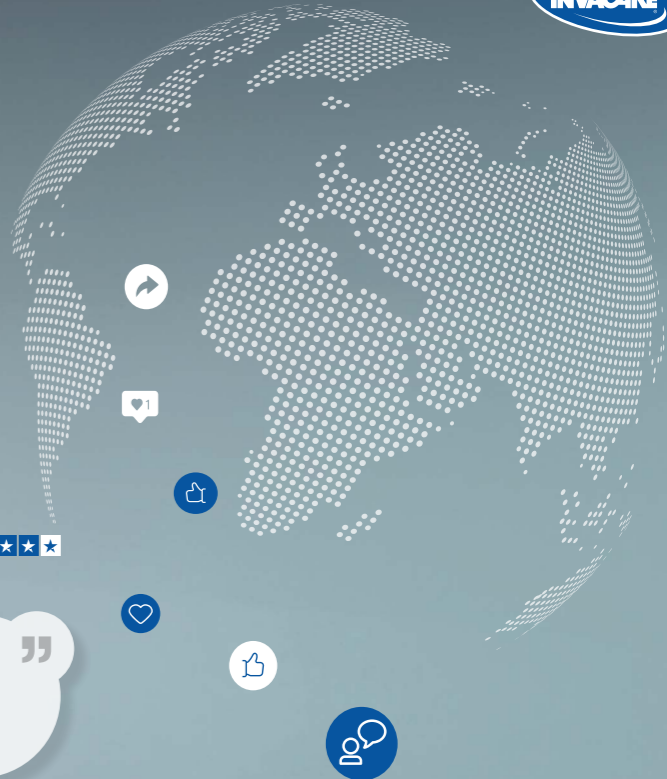



P O W E R   S E R I E S

**AVIVA RX**

AVIVA FX

TDX SP2

Esprit Action

Storm 4

**T: +44 1656 776222**

[uk@invacare.com](mailto:uk@invacare.com)

[www.invacare.co.uk](http://www.invacare.co.uk)

[www.invacare.ie](http://www.invacare.ie)

© 2024 Invacare International GmbH. All rights reserved.  
All information quoted is believed to be correct at time of print.  
Invacare reserves the right to alter product specifications without  
prior consultation. UK - AVIVA RX ULM - 11/2024

