



erivo®

The light & foldable
electric wheelchair



Be there:
Be free.

NEW



+ erivo. always there.

www.erivo.de



Simply mobile.

erivo – powered by Alber.

As the market leader in mobility aids for wheelchairs, Alber GmbH has more than 35 years of experience in the field of electromobility. We are a recognized partner to health insurance companies and work with healthcare suppliers across the world. From the initial consultation through to service: You can rely on us. Stay „simply mobile“ with Alber.

Up and away – arrive relaxed with erivo.

Your compact companion for daily life and travel.

Being independent and self-determined, going out into the fresh air with friends, enjoying the sun in the park with the family, going on a long-awaited vacation with your partner: erivo from Alber helps keep you mobile and gives you a new quality of life. erivo folds with just a flick of the wrist and thanks to its compact size can be easily carried in any form of transport. Dream of destinations again instead of worrying about the path.

erivo – introduction

Test now at your leisure at no cost!

Make your appointment for a no-obligation in-home test ride conveniently online:

www.erivo.de



35 years experience



Health insurance partner



Partner of specialist dealers



elegant & practical

+ erivo® design



In our ever more networked world, one thing counts above all else: mobility. That's why we need means of transportation that suit our lifestyle and help us overcome obstacles. They must be reliable, flexible and user-friendly – and of course comfortable.



In top shape.

erivo looks like it rides: dynamic & contemporary. With its elegant shape, color scheme and high-quality materials, it sets new standards in design and always ensures a stylish impression.

- ✓ **Two-tone paint finish**
in champagne gold & matte black
- ✓ **Sporty & rugged frame design**



+ erivo® comfort



The perfect fit.

One of the most important factors in an electric wheelchair is seat comfort. erivo boasts well-thought-out ergonomics, the best materials and top workmanship. High-quality foam ensures sitting comfort and good support at all times.

- ✓ **Ergonomic**
firm seat & high-quality foam
- ✓ **Breathable materials**



Keep rolling.

Effective rear-wheel suspension combined with pneumatic tires ensures outstanding ride comfort. The large front wheels are absolutely flat-proof, offer good absorption and ideal support when surmounting small obstacles.

- ✓ **Suspension and large pneumatic tires**
12 1/2" rear
- ✓ **Flat-proof PU wheels**
8" front



easy & flexible

+ erivo® function



Smart folding technology.

With a flick of the wrist, erivo folds to save space. Hardly larger than a small suitcase, it can fit into any niche – or easily accompany you on your next trip.

Control unit, armrests & footrest.

erivo offers maximum flexibility with its swing-away joystick controller and the foldable armrests. This allows you to come close to tables or change places with ease. The footrest also swivels and is height-adjustable.



Power pack.

The powerful battery is installed under the seat and can be removed without tools. With the supplied automatic charger, it can also be charged while installed. The maximum charge level is reached in only 6 hours.

- ✓ Folds with a flick of the wrist
- ✓ Fast stowage & space saving

- ✓ Height-adjustable footrest
- ✓ Flexible control unit & armrests

- ✓ Max. charging time: 6 hours
- ✓ Range: up to (11/22 mi)*

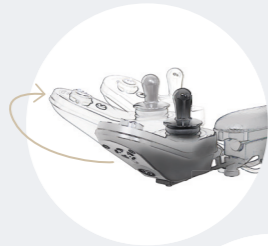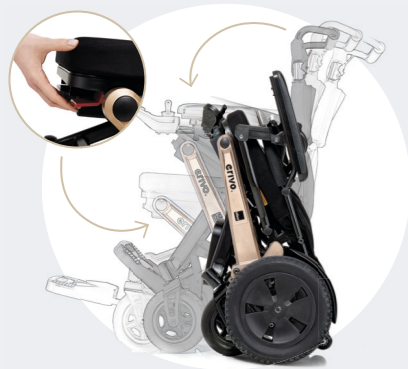

+ erivo® dynamism



Intuitive control.

The intuitive joystick offers infinitely variable and precise travel speed control. This allows erivo to navigate safely in tight spaces as well as to play tag with the grandkids in the park.

- ✓ Precise & nimble control
- ✓ Stepless speed selection (1-4 mph)



Optionally available
with attendant control
more at
www.erivo.de



* The range depends on various factors such as outside temperature, route conditions, passenger weight and tire pressure and may vary. Tested as per ISO 7176-4 TÜV Süd test procedure.

In case of wanderlust or travel fever: Take erivo daily.

Recover your favorite places.

erivo helps not only in everyday life. Its compact size makes it portable in any means of transport, and it's the ideal companion for your next trip.

Whether it's a long-planned family vacation or a spontaneous trip with friends - you and your erivo are always there for the experience. And if the long-lasting battery does run out, no need to worry: erivo is fully charged again in only 6 hours and ready for a new adventure.



car



airplane



ship



bus



caravan



train



erivo - wanderlust



Be there: Be free.

Unfold experience.

Folded with just a flick of the wrist, erivo can easily travel in any mode of transport - and unfold just as easily when needed. Discover new areas or enjoy your personal favorite place with no worries. During development, we placed particular

emphasis on your seating comfort. This is ensured by the suspension, the large rear wheels with pneumatic tires and the ergonomic seat. Take a seat and experience great moments with a little help from erivo.



+ erivo® overview

Swiveling control unit

Allows moving close to the table



Smart folding technology & compact design

Folded dimensions (w x h x d) 25 x 34 x 21 in

Power pack

Long range (11/22 mi) & rapid charging

Foldable armrests

Easy sitting down & standing up



Footrest
Swiveling & height-adjustable

Control unit

Intuitive control by joystick
Stepless speed control from 1 to 4 mph

Best ride characteristics

Large, flat-proof 8" PU front wheels
Pneumatic tires & special rear-wheel suspension

Ideal seating ergonomics

Firm seat / contoured seat
Breathable, high-quality materials

Innovative design

Elegant but sporty frame in champagne gold & matte black



Book your free test ride:



www.erivo.de

Test now for free.

Book your appointment now for a free, no-obligation erivo test ride in the familiar comfort of your home. Call us at the service number **+49 (0)7432 2006-0** or book your demonstration conveniently online at **www.erivo.de**

At erivo.de you can find even more information about erivo or about reimbursement by your health insurance company. The erivo team at Alber GmbH looks forward to meeting you and will be happy to provide you with advice and assistance.

+49 (0)7432 2006-0

Service hours
Mon–Thu: 8 a.m. to 6 p.m.
Fri: 8 a.m. to 4 p.m.

www.erivo.de

Additional technical data

Shipping weight
Battery removed: from 60 lbs

Standard battery pack weight: 5.3 lbs

Replacement battery pack weight: 10 lbs

Battery capacity
Standard battery pack: 11.0 Ah
Optional replacement battery pack: 30.15 Ah

Seat width: 17.7 in
Seat depth: 18.5 in
Back height: 17.7 in
Seat height: 20 in
Total length: 38 in
Power output: 250 watts
Max. user weight: 330 lbs

Protection class: IPX4 (splash-protected)

erivo®

alber

The light & foldable
electric wheelchair

NEW



10.0001.6.01.01 As at 08/2022 @Alber GmbH, Albstadt.
Subject to technical modifications.

Alber GmbH

Vor dem Weißen Stein 14
72461 Albstadt-Tailfingen

Phone: +49 (0)7432 2006-0
Mon-Thu: 8 a.m. to 6 p.m.
Fri: 8 a.m. to 4 p.m.
info@alber.de
www.alber.de

Our quality standard: Your safety!

Our company is certified in accordance with the current quality management standard ISO 13485 for medical devices.

We ensure the environmentally responsible manufacture of our products with our environmental management system, certified in accordance with ISO 14001.

erivo is CE-compliant according to the European Medical Device Regulation (MDR) 2017/745 and is tested and certified on a voluntary basis by TÜV Süd, according to the current valid edition of the product standard EN 12184.