

Raizer[®] II

Lifting Chair



Lifts a fallen person in 30 seconds

Safe, dignified floor-to-chair lifting – operable by one person




Feel safe and comfortable with Raizer II

Raizer II is a battery-operated mobile lifting chair that helps raise someone who has fallen from the floor to an almost standing position in just 20–30 seconds.

Raizer can easily be assembled and operated by only one person and does not require any physical effort besides a supportive hand.

The design of the Raizer makes it easy to assemble and operate, and ensures easy cleaning.

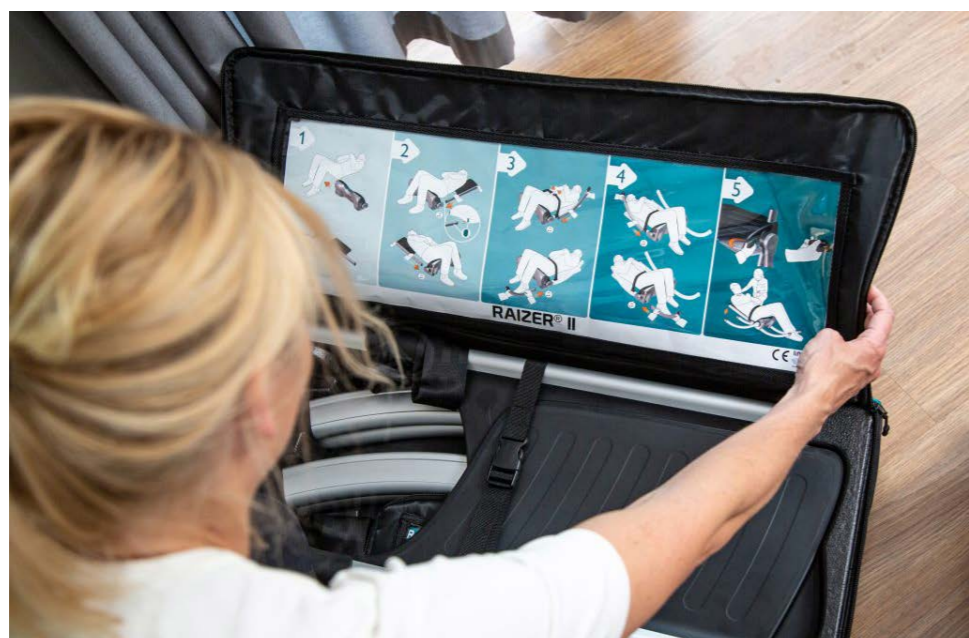

80 lifts per charge

Up to 80 lifts on a single charge – 10–15 minutes charging time for one additional lift.



4 Identical legs and 2 identical backrests

Simple assembly – all parts are interchangeable for correct setup every time.




> **All surface materials** have been specified to be **cleaned easily and disinfect**

> **Additional seat belts** can be fitted for increased stability of the fallen person if this is required

> **Intuitive sound and led indicators** for correct assembly

> **The remote control** is stored on the seat unit and assured

> **Disposable hygiene covers** are also available (optional accessory)

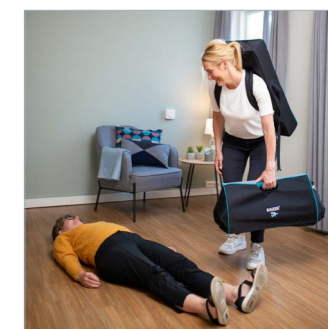

Step 1 – Place the seat unit next to the fallen person



Step 2 – Assemble legs and backrest around the person



Step 3 – Secure the person and start the lift



Step 4 – Person reaches near-standing position – transfer to chair

Technical data ▼



Lifting capacity ▼

Max. 150kg



Lifting time ▼

20-30 seconds



Weight ▼

Seat: 8.5kg
Legs and back support: 4.5kg
Total: 13kg



Number of lifts when fully charged ▼

40 at max. load.
Approx. 80 with an average load



Charging time ▼

Up to 6 hours
