



rea | ***Azalea[®] Base***

The perfect tilt in space base for custom
made seat units



Rea Azalea Base is designed for people seeking a custom-made seat combined with a first-class wheelchair base. The *Rea Azalea Base* is wider in the rear section of the seat, offering the possibility to use larger seating systems on smaller frames. As there are many balance options, the pushing is easier. Another key feature of the Base is the 'Tilt in Space' function - a weight shifting mechanism that ensures reliable stability. All the advantages of our standard *Rea Azalea* models are included. What's more, the *Azalea Base* can be used together with all the parts and accessories that come with other *Rea Azalea* models (excluding trunk supports). Available in 5 colours.

Seat adaptation

Rea Azalea Base is specifically designed to accommodate custom-made seats with extra space in the rear of the seat frame and seat adjustability. The backrest is angle adjustable including an increased recline which can be angled via a 'rigging screw' or an attendant operated gas piston.

The seat width can be adjusted from the arm and leg rests, and the depth changed via the telescopic seat frame. These features ensure that the *Rea Azalea Base* easily adapts to different-shaped seat moulds.

For other *Rea Azalea* accessories, please consult the *Rea Azalea* brochure

Features and Options



Fixing attachments

The seat offers various attachment holes, plus greater 'open space' at the back.



Harness attachment

An attachment for different kind of harnesses is also available as an option.




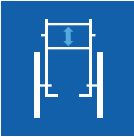














Seat plate with turnable wedge

An additional +5° seat angle can be obtained, using the optional seat plate with a wedge.

Technical data


For more comprehensive pre-sales information about this product, including the product's user manual, please visit www.invacare.co.uk

	Seat width	Seat depth	Seat height	Backrest height	Armrest height	Length of footrests	Seat angle	Backrest angle
<i>Rea Azalea Base</i>	 340 / 390 / 440 / 490 mm (+ 50 / 100 mm)	 380 - 450 mm ¹ 430 - 500 mm	 400 / 450 mm	 600 / 700 mm	 240 - 360 mm	 330 - 500 mm	 -1° to 23° (22 / 24°) (+1° to 25°) (16°)	 0° - +30°
	Driving unit width	Total height	Total length inc. Footrests	Total weight	Max. User weight	Transport weight	Crash tested	Maximum slope for use of parking brake
<i>Rea Azalea Base</i>	 SW + 220 mm (SW + 230 mm)	 1300 / 1400 mm	 930 - 1040 / 950 - 1070 mm	 20 kg	 Max 135 kg	 15 kg	 ISO 7176-19	 7°

1. Crash tested in specific configurations

Frame colour

Pearl Grey Azurite Blue Sand Happy Red Electric Green

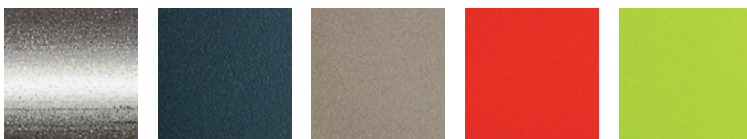

Invacare Limited
Pencoed Technology Park - Pencoed
Bridgend CF35 5AQ - United Kingdom
Tel: +44 1 656 776 222
Fax: +44 1656 776 220
E-mail: uk@invacare.com
Sales Order E-mail: ordersuk@invacare.com
www.invacare.co.uk

Invacare Ireland Ltd
Unit 5 - Seatown Business Campus
Seatown Road - SWORDS - County Dublin - Ireland
Tel. +353 1 8107084
Fax +353 1 8107085
Email: Ireland@invacare.com
www.invacare.ie

