

Spare Parts Catalogue

Version: 2017/03/02a NO



**Invacare®
TDX-SP2 / TDX-SP2 NB**



Yes, you can.®

General Information

For any order of spare parts the following information is required:

- Product type*
- Product name & model*
- Serial Number*
- Part-Number*
- Part-Description*

Exploded view drawings and code numbers are not contractual and subject to modification, in relation to Invacare® continuous product improvement policy to meet end users needs.

The reproduction or use of this work, or of any substantial part of this work, in any form or by any electronic, mechanical or other means, including photocopying, scanning and recording, and in any information storage and retrieval system, including electronic transmission, is forbidden without the prior written permission of Invacare® International Sarl.



Yes, you can.®

TABLE OF CONTENTS

Informasjon & brukermanual

(User Manuals TDX-SP2)	15
(User Manuals Ultra Low Maxx Seat)	16
(User Manuals Shark2)	17
(User Manuals REM A & B)	18
(User Manuals REM 24 SD)	19
(User Manuals REM 550)	20
(User Manuals G90A)	21
(User Manuals G91S)	22
(User Manuals Chin Control (2015-10-01 =>))	23
(User Manuals REM 34 & REM 41 Tray/Centre Control)	24
(User Manuals IDC (Intuitive Dual Control))	25
(User Manuals Compact Dual Control)	26
(User Manuals ACS Handprogrammer)	27
(User Manual Proximity Head Remote DSL-105)	28
(User Manuals Headcontrol)	29
(User Manuals iPortal)	30

Chassis

Chassis	31
(Chassis Standard Base)	32
(Battery Box Hardware)	33
(Walking Beam)	34
(Rear Pivot Standard Base)	35
(Shrouds Standard Base)	36
(Rear Shrouds - Standard Base)	37
(Bezel Shroud)	38
(Front Shroud - Standard Base)	39
(Fender)	40
(Rear Caps)	41
(Low Rider Option)	42
(Chassis Narrow Base)	43
(Battery Box Hardware)	44
(Walking Beam (Narrow Base))	45
(Rear Pivot Narrow Base)	46
(Shrouds Narrow Base)	47
(Rear Shrouds Narrow Base)	48
(Bezel Shroud)	38
(Front Shroud Narrow Base)	49
(Fender)	40
(Rear Caps)	41

Sete

Sete	50
(Modulite Seat)	51
(man. settable Tilt)	52
(powered simple Tilt 20°)	53
(powered Tilt 0°-30°)	54
(Powered seat lifter with tilt 300mm/30°)	55
(Anterior Tilt)	56



Yes, you can.®

TABLE OF CONTENTS

(One Piece Seat Plate)	57
(Telescopic Seat Frame and Plate)	58
(Telescopic Seat Frame for man. /pow. 20° Tilt)	59
(Telescopic Seat Frame for Tilt 30° / Lifter)	60
(Armrest Bracket)	61
(Legrest Support)	62
(Seat Plate)	63
(Telescopic Seat Frame and Sling Seat)	64
(Telescopic Seat Frame for man. / pow. 20° Tilt)	65
(Telescopic Seat Frame for Tilt 30° / Lifter)	66
(Armrest Bracket)	61
(Legrest Support)	67
(Sling Seat)	68
(Telescopic Seat Frame and Universal Adapter)	69
(Telescopic Seat Frame for man. /pow. 20° Tilt)	59
(Telescopic Seat Frame for Tilt 30° / Lifter)	60
(Armrest Bracket)	61
(Legrest Support)	62
(Adapter Kit for Universal Seat)	70
(Seat Cushions "Standard" & "Matrx Flo-Tech")	71
(Seat Cushions "Matrx Libra")	74
(Seat Cushion "Matrx Vi")	75
(Ultra Low Maxx Seat)	76
(Seat Tilt)	77
(Seat Lifter with Tilt)	78
(Anterior- / Posterior-Tilt)	79
(Seat Plates & Side Rails)	80
(Seat Cushions "Matrx Vi")	81
(Seat Cushions "Matrx Flo-Tech")	82
(Seat Cushions "Matrx Libra")	83
(Max Seat)	84
(man. settable Tilt)	85
(powered simple Tilt)	86
(Seat Frame and Plate)	87
(Seat Cushion)	88

Rygg

Rygg	89
(Backrest for Modulite Seat)	90
(Support for Backrest Frame)	91
(Backrest Frame)	92
(Sling Back, tension adjustable)	93
(Backrest support)	94
(Backposts)	95
(Connection tube)	96
(Back Recliner / Stabilization bracket)	97
(Backrest Cushions)	98
(Flex 3 / Laguna)	99
(Matrx Posture Back)	100
(Matrx Elite T/R Back)	101
(Matrx Elite Deep Back)	102
(Matrx Elite (Full) Back)	103



Yes, you can.®

TABLE OF CONTENTS

(Mounting Hardware for Matrix backrests)	104
(Flex3 back, tension adjustable)	105
(Backrest for Ultra Low Maxx Seat)	106
(Backrest Recliner)	107
(High Back)	108
(Backrest Hardware)	109
(Backrest Cushions)	110
(Matrix Back)	112
(Backrest Hardware)	113
(Posture Back)	114
(Elite T/R Back)	115
(Elite Deep Back)	116
(Privacy Flaps)	117
(Backrest for Max Seat)	118
(Manual backrest adjustment)	119
(Powered backrest adjustment)	120
(Seat depth adapter)	121
(Backrest)	122

Benstøtte

Benstøtte	123
(Legrests - Standard Base)	124
(Legrests for Modulite Seat)	125
(Legrests 80° & man. adjustable)	126
(Legrest Support, 80° & man. elevating - Standard & Modulite Seat)	127
(Hanger, Std. 80° (01.07.15 =>))	128
(Hanger, Std. 80° (=> 30.06.15))	129
(Hanger, manually height adjustable)	130
(Footplate, fixed 80°)	131
(Footplate, depth-, angleadjustable, 80°)	132
(Amputee Support)	133
(Legrest Support, Vari series - Modulite Seat)	134
(Hanger Vari A (04/2009 =>))	135
(Hanger Vari F)	136
(Hangers ADM LH/RH manual)	137
(Hanger Vari ADE/ADP electric (07/2005 =>))	138
(Footrests Vari, Flip-Up (01.08.10 =>))	139
(Footrests Vari, Flip-Up & angle adjustable (01.08.10 =>))	140
(Footrests Vari, angle & depth adjustable)	141
(Footboard Vari, pre-assembled)	142
(Central Legrests - Modulite Seat)	143
(Support Centr. Legrest man.)	144
(Support Centr. Legrest pow. LNX (=> 2016-06-30))	145
(Support Centr. Legrest pow. LNX (2016-07-01 =>))	146
(Central Legrest man.)	147
(Central Legrest pow. LNX (=> 2016-06-30))	148
(Central Legrest pow. LNX (2016-07-01 =>))	149
(Legrests for Ultra Low Maxx Seat)	150
(Central Legrests)	151
(Supports, Central Legrests)	152
(Central Legrest pow. LNX)	153
(Central Legrest man.)	147



Yes, you can.®

TABLE OF CONTENTS

(Vari Legrests)	154
(Supports, Vari)	155
(Hanger Vari A (04/2009 =>))	135
(Hanger Vari F)	136
(Hangers ADM LH/RH manual)	137
(Hanger Vari ADE/ADP electric (07/2005 =>))	138
(Footrests Vari, Flip-Up (01.08.10 =>))	139
(Footrests Vari, Flip-Up & angle adjustable (01.08.10 =>))	140
(Footrests Vari, angle & depth adjustable)	141
(Pivot Plus Legrests)	156
(Cross Brace, Pivot Plus)	157
(Legrests)	158
(Legrests for Max Seat)	159
(Legrest Support, Vari series - Max Seat)	160
(Hanger Vari A (04/2009 =>))	135
(Hanger Vari F)	136
(Hangers ADM LH/RH manual)	137
(Hanger Vari ADE/ADP electric (07/2005 =>))	138
(Footrests Vari, Flip-Up (01.08.10 =>))	139
(Footrests Vari, Flip-Up & angle adjustable (01.08.10 =>))	140
(Footrests Vari, angle & depth adjustable)	141
(Footboard Vari, pre-assembled)	142
(Support Centr. Legrest man. LNX - Max Seat)	161
(Support Centr. Legrest LNX el. - Max Seat)	162
(Central Legrest man.)	147
(Central Legrest pow. LNX (=> 2016-06-30))	148
(Fixed footboard)	163
(Calf plates)	164
(Legrests - Narrow Base)	165
(Legrests for Modulite Seat)	166
(Legrests 80° & man. adjustable)	126
(Legrest Support, 80° & man. elevating - Standard & Modulite Seat)	127
(Hanger, Std. 80° (01.07.15 =>))	128
(Hanger, Std. 80° (=> 30.06.15))	129
(Hanger, manually height adjustable)	130
(Footplate, fixed 80°)	131
(Footplate, depth-, angleadjustable, 80°)	132
(Amputee Support)	133
(Legrest Support, Vari series - Modulite Seat)	134
(Hanger Vari A (04/2009 =>))	135
(Hanger Vari F)	136
(Hangers ADM LH/RH manual)	137
(Hanger Vari ADE/ADP electric (07/2005 =>))	138
(Footrests Vari, Flip-Up (01.08.10 =>))	139
(Footrests Vari, Flip-Up & angle adjustable (01.08.10 =>))	140
(Footrests Vari, angle & depth adjustable)	141
(Central Legrests - Modulite Seat)	143
(Support Centr. Legrest man.)	144
(Support Centr. Legrest pow. LNX (=> 2016-06-30))	145
(Support Centr. Legrest pow. LNX (2016-07-01 =>))	146
(Central Legrest man.)	147
(Central Legrest pow. LNX (=> 2016-06-30))	148



Yes, you can.®

TABLE OF CONTENTS

(Central Legrest pow. LNX (2016-07-01 =>))	149
(Legrests for Ultra Low Maxx Seat)	150
(Central Legrests)	151
(Supports, Central Legrests)	152
(Central Legrest pow. LNX)	153
(Central Legrest man.)	147
(Vari Legrests)	154
(Supports, Vari)	155
(Hanger Vari A (04/2009 =>))	135
(Hanger Vari F)	136
(Hangers ADM LH/RH manual)	137
(Hanger Vari ADE/ADP electric (07/2005 =>))	138
(Footrests Vari, Flip-Up (01.08.10 =>))	139
(Footrests Vari, Flip-Up & angle adjustable (01.08.10 =>))	140
(Footrests Vari, angle & depth adjustable)	141
(Pivot Plus Legrests)	156
(Cross Brace, Pivot Plus)	157
(Legrests)	158

Armlen

Armlen	167
(Armrests for Modulite Seat)	168
(T-Armrest)	169
(Padded Armrest)	170
(Inrigged Armrest)	171
(Flip-Up Armrest)	172
(Following Armrest)	173
(Hip Support)	174
(Armrests for Ultra Low Maxx Seat)	175
(Recline Flip-Up Armrest)	176
(Cantilever Flip-Up Armrest)	177
(Armpads, Desk & Full)	178
(Armpad, Ergonomic)	179
(Elbow Block)	180
(MACES Multi Axis Armpad Holder)	181
(Hip Support)	174
(Armrests for Max Seat)	182

Svinghjul

(Castors)	183
-----------	-----

Drivhjul

(Drive Wheels 14")	184
(Wheel Hub 5-hole)	185

Fremdrift

(Overview Propulsion SSD)	186
(Propulsion SSD - Standard Base)	189
(Propulsion SSD - Narrow Base)	190
(Wheel Hub 5-hole)	185



Yes, you can.®

TABLE OF CONTENTS

Elektronikk

Elektronikk	191
(Electronic - Standard Base)	192
(ACS 2 Electronic)	193
(ACS Electronic Modules)	194
(Power Modul - ACS)	195
(Bracket for Power Modul - Standard Base)	196
(ACTs)	197
(Bracket ACT - Modulite Seat)	198
(Bracket ACT - Ultra Low Maxx Seat)	199
(Bracket ACT - Max Seat)	200
(Light Board)	201
(G-Trac Sensor/Gyro)	202
(ACS Remotes & Brackets)	203
(REM 24 SD with brackets)	204
(Remote REM 24SD)	205
(Remote Bracket Standard)	206
(Remote Bracket Standard - ULM and Soft Armpad on Modulite)	207
(Remote Bracket Swing away & height adjust.)	208
(Remote Bracket Swing away & height adjust. - ULM and Soft Armpad on Modulite)	209
(Remote Bracket Attendant)	210
(Remote Bracket Attendant Ultra Low Maxx Seat)	211
(REM A & B with brackets)	212
(Remote Bracket Standard)	206
(Remote Bracket Standard - ULM and Soft Armpad on Modulite)	207
(Remote Bracket Swing away & height adjust.)	213
(Remote Bracket Swing away & height adjust. - ULM and Soft Armpad on Modulite)	209
(Remote Bracket Attendant)	210
(Remote Bracket Attendant Ultra Low Maxx Seat)	211
(REM 550 with brackets)	214
(Remote Bracket Standard)	206
(Remote Bracket Standard - ULM and Soft Armpad on Modulite)	207
(Remote Bracket Attendant)	210
(Remote Bracket Attendant Ultra Low Maxx Seat)	211
(Remote Bracket Swing away & height adjust.)	208
(Remote Bracket Swing away & height adjust. - ULM and Soft Armpad on Modulite)	209
(Remote Bracket Swing away Quad Link Ultra Low Maxx Seat)	215
(REM 34 with brackets)	216
(Remote Bracket Standard)	217
(Remote Bracket Swing away & height adjust.)	218
(G90A with brackets)	219
(Remote Bracket Standard)	220
(Remote Bracket Standard - ULM and Soft Armpad on Modulite)	221
(Remote Bracket Swing away & height adjust.)	222
(Remote Bracket Swing away & height adjust. - ULM and soft Armpad on Modulite)	223
(Remote Bracket Attendant)	224
(Remote Bracket Attendant Ultra Low Maxx Seat)	225
(G91S with brackets)	226



Yes, you can.®

TABLE OF CONTENTS

(Remote Bracket Standard)	227
(Centre control REM 41)	228
(Chin Control, small (=> 16/12/31))	229
(Small Chincontrol)	230
(Framework for chincontrol)	231
(Chin Control - Modulite Seat (15/10/01 =>))	232
(User Manuals Chin Control (2015-10-01 =>))	23
(Chin Control DX-RJM)	233
(Chin Control powered - Details)	234
(Chin Control manual - Details)	235
(Mounting Hardware - Chin Control)	236
(Dual Control with brackets)	237
(Remote Bracket Dual Control for Modulite Seat)	238
(Remote Bracket Dual Control for Max Seat)	239
(Compact Dual Control with brackets)	240
(Compact Dual Control)	241
(Bracket, Compact Dual Control - Ultra Low Maxx Seat)	242
(Bracket, Compact Dual Control - short)	243
(Bracket, Compact Dual Control - long)	244
(Intuitive Dual Control (IDC) with brackets)	245
(IDC with bracket)	246
(User Manuals IDC (Intuitive Dual Control))	25
(IDC with bracket (Ultra Low Maxx Seat))	247
(User Manuals IDC (Intuitive Dual Control))	25
(Height adjustable bracket IDC)	248
(ACS Cables & Accessories)	249
(Bus-Cables ACS)	250
(Light Cable ACS)	251
(T-Adapter (Modulite & Max Seat))	252
(T-Adapter (Ultra Low Maxx Seat))	253
(T-Adapter (Low Rider))	254
(Tilt Sensors - Ultra Low Maxx Seat)	255
(Programming Kits)	256
(Magnetkey for Remotes)	257
(ACS Components)	258
(iPortal Modules and Accessories)	259
(Remote, Joystick, Control)	260
(Joystick Knob)	261
(Buddy Button)	262
(Piko Button)	263
(Interface, Adapter, Other)	264
(LiNX Electronic)	265
(Modules)	266
(Power Modules)	267
(Mounting Hardware, Power Module & Light Board)	268
(Actuator Modules (ACTs))	269
(Bracket ACT - Modulite Seat)	198
(Bracket ACT - Max Seat)	200
(G-Trac Sensor (Gyro))	270
(Remotes & Brackets)	271
(User Manuals LiNX)	272
(Remotes & Brackets)	273



Yes, you can.®

TABLE OF CONTENTS

(Remote REM110)	274
(Remote REM211)	275
(Remote REM216)	276
(Remote Bracket, Standard)	277
(Remote Bracket, Swing away)	278
(Remote Bracket, Attendant)	279
(Dual Control & Bracket)	280
(Dual Control)	281
(Remote Bracket, Attendant)	279
(Cables & Accessories)	282
(Bus Cables)	283
(Light Cable)	284
(Light Board)	201
(Lifter/Tilt Module Cables)	285
(Actuator Cable)	286
(4-way Expansion Block)	287
(4-way Expansion Block (with Tilt module))	288
(4-way Expansion Block (with Lifter/Tilt module))	289
(4-way Expansion Block (Max Seat))	290
(4-way Expansion Block (Low Rider))	291
(4-way Expansion Block (without any module))	292
(Components)	293
(Accessory cables)	294

(10-way switch & brackets)

(10-way switch & brackets)	295
(10-way switch)	296
(Bracket, front)	297
(Bracket, opposite side)	298
(Bracket, rear)	299
(Buddy Button)	262
(Piko Button)	263
(Electronic - Narrow Base)	300
(ACS 2 Electronic)	301
(ACS Electronic Modules)	302
(Power Modul - ACS)	195
(Bracket for Power Modul - Narrow Base)	303
(ACTs)	197
(Bracket ACT - Modulite Seat)	198
(Bracket ACT - Ultra Low Maxx Seat)	199
(Bracket ACT at Chassis - Narrow Base)	304
(Light Board)	305
(G-Trac Sensor/Gyro)	202
(ACS Remotes & Brackets)	306
(REM 24 SD with brackets)	204
(Remote REM 24SD)	205
(Remote Bracket Standard)	206
(Remote Bracket Standard - ULM and Soft Armpad on Modulite)	207
(Remote Bracket Swing away & height adjust.)	208
(Remote Bracket Swing away & height adjust. - ULM and Soft Armpad on Modulite)	209
(Remote Bracket Attendant)	210



Yes, you can.®

TABLE OF CONTENTS

(Remote Bracket Attendant Ultra Low Maxx Seat)	211
(REM A & B with brackets)	212
(Remote Bracket Standard)	206
(Remote Bracket Standard - ULM and Soft Armpad on Modulite)	207
(Remote Bracket Swing away & height adjust.)	213
(Remote Bracket Swing away & height adjust. - ULM and Soft Armpad on Modulite)	209
(Remote Bracket Attendant)	210
(Remote Bracket Attendant Ultra Low Maxx Seat)	211
(REM 550 with brackets)	214
(Remote Bracket Standard)	206
(Remote Bracket Standard - ULM and Soft Armpad on Modulite)	207
(Remote Bracket Attendant)	210
(Remote Bracket Attendant Ultra Low Maxx Seat)	211
(Remote Bracket Swing away & height adjust.)	208
(Remote Bracket Swing away & height adjust. - ULM and Soft Armpad on Modulite)	209
(Remote Bracket Swing away Quad Link Ultra Low Maxx Seat)	215
(REM 34 with brackets)	216
(Remote Bracket Standard)	217
(Remote Bracket Swing away & height adjust.)	218
(G90A with brackets)	219
(Remote Bracket Standard)	220
(Remote Bracket Standard - ULM and Soft Armpad on Modulite)	221
(Remote Bracket Swing away & height adjust.)	222
(Remote Bracket Swing away & height adjust. - ULM and soft Armpad on Modulite)	223
(Remote Bracket Attendant)	224
(Remote Bracket Attendant Ultra Low Maxx Seat)	225
(G91S with brackets)	226
(Remote Bracket Standard)	227
(Centre control REM 41)	228
(Chin Control, small (=> 16/12/31))	229
(Small Chincontrol)	230
(Framework for chincontrol)	231
(Chin Control - Modulite Seat (15/10/01 =>))	232
(User Manuals Chin Control (2015-10-01 =>))	23
(Chin Control DX-RJM)	233
(Chin Control powered - Details)	234
(Chin Control manual - Details)	235
(Mounting Hardware - Chin Control)	236
(Dual Control with brackets - Modulite Seat)	307
(Remote Bracket Dual Control for Modulite Seat)	238
(Compact Dual Control with brackets)	240
(Compact Dual Control)	241
(Bracket, Compact Dual Control - Ultra Low Maxx Seat)	242
(Bracket, Compact Dual Control - short)	243
(Bracket, Compact Dual Control - long)	244
(Intuitive Dual Control (IDC) with brackets)	245
(IDC with bracket)	246
(User Manuals IDC (Intuitive Dual Control))	25
(IDC with bracket (Ultra Low Maxx Seat))	247
(User Manuals IDC (Intuitive Dual Control))	25



Yes, you can.®

TABLE OF CONTENTS

(Height adjustable bracket IDC)	248
(ACS Cables & Accessories)	308
(Bus-Cables ACS)	250
(Light Cable ACS)	251
(T-Adapter - Modulite Seat)	309
(T-Adapter (Ultra Low Maxx Seat))	253
(Tilt Sensors - Ultra Low Maxx Seat)	255
(Programming Kits)	256
(Magnetkey for Remotes)	257
(ACS Components)	258
(iPortal Modules and Accessories)	259
(Remote, Joystick, Control)	260
(Joystick Knob)	261
(Buddy Button)	262
(Piko Button)	263
(Interface, Adapter, Other)	264
(Shark Electronic)	310
(Shark Electronic Modules)	311
(Power Modul)	312
(Bracket for Power Modul)	313
(Light Board)	201
(Shark Remote & Brackets)	314
(Remote Shark with Brackets)	315
(Remote Bracket Standard)	206
(Remote Bracket Swing away & height adjust.)	213
(Remote Dual Control with Brackets)	316
(Remote Bracket Dual Control)	317
(Shark Cables & Accessories)	318
(Bus-Cables Shark)	319
(Cables Shark)	320
(Drive Lockout)	321
(Drive Lockout with bracket)	322
(Drive Lockout Cables - Shark)	323
(Programming Kits)	256
(LiNX Electronic)	324
(Modules)	325
(Power Modules)	267
(Mounting Hardware, Power Module & Light Board)	326
(Actuator Modules (ACTs))	269
(Bracket ACT - Modulite Seat)	198
(Bracket ACT at Chassis - Narrow Base)	304
(G-Trac Sensor (Gyro))	270
(Remotes & Brackets)	271
(User Manuals LiNX)	272
(Remotes & Brackets)	273
(Remote REM110)	274
(Remote REM211)	275
(Remote REM216)	276
(Remote Bracket, Standard)	277
(Remote Bracket, Swing away)	278
(Remote Bracket, Attendant)	279
(Dual Control & Bracket)	280



Yes, you can.®

TABLE OF CONTENTS

(Dual Control)	281
(Remote Bracket, Attendant)	279
(Cables & Accessories)	327
(Bus Cables)	283
(Light Cable)	284
(Light Board)	201
(Lifter/Tilt Module Cables)	285
(Actuator Cable)	286
(4-way Expansion Block)	328
(4-way Expansion Block (with Tilt module))	288
(4-way Expansion Block (with Lifter/Tilt module))	289
(4-way Expansion Block (without any module))	292
(Components)	293
(Accessory cables)	294

(10-way switch & brackets)

(10-way switch & brackets)	295
(10-way switch)	296
(Bracket, front)	297
(Bracket, opposite side)	298
(Bracket, rear)	299
(Buddy Button)	262
(Piko Button)	263

Batterier, kabler og ladere

Batterier, kabler og ladere	329
(Battery, Tray & Wiring - Standard Base)	330
(Battery, Tray & Wiring - Standard Base)	331
(Battery, Tray & Wire Harness - Narrow Base)	332
(Chargers)	333

Sikkerhedsudstyr og tilbehør

(Lightsystem)	334
(Lightsystem LED)	335
(Lightsystem german Homologation)	336
(Rear Light Protector Kit)	337
(Reflective Sign)	338
(Reflective Sign)	339
(Reflective Sign (Ultra Low Maxx Seat))	340
(Adapter-Kit DAHL Docking Station - Standard Base)	341
(Adapter-Kit DAHL Docking Station - Narrow Base)	342
(Belts)	343
(Lap Belts (Modulite Seat))	344
(Lap Belt (Ultra Low Maxx))	345
(Lap Belt (Max Seat))	346
(H-Belts (Modulite Seat))	347
(Rearview Mirror)	348

Setetilbehør

(Headrests)	349
-------------------	-----



Yes, you can.®

TABLE OF CONTENTS

(Head- / Neckrest)	350
(Headrest sidewise adjustable)	351
(Headrest with cheek support)	352
(Headrests - Ultra Low Maxx Seat)	353
(Lateral Supports)	354
(Lateral Supports - Modulite Seat)	355
(Lateral supports, Flex3 backrest)	356
(Lateral supports, Flex3 backrest - multi adjustable)	357
(Lateral supports, Flex3 backrest (tension adjustable) - multi adjustable)	358
(Lateral Support - Ultra Low Maxx Seat)	359

Tilbehør

(Cane holder)	360
(Cane holder (Modulite: backrest frame))	361
(Cane holder (Modulite: back posts man recliner))	362
(Cane holder (Modulite: back posts pow recliner))	363
(Cane holder (MAX Seat))	364
(Trays)	365
(Trays)	366
(Trays (Ultra Low Maxx Seat))	367
(Hour Counter)	368
(T-Nuts for Accessories Modulite & ULM Seat)	369

(User Manuals TDX-SP2)

languages.iso

BG	Български	FR	Français	NL	Nederlands	SL	Slovenščina
CS	Čeština	HE	עברית	NO	Norsk	SV	Svenska
DA	Dansk	HR	Hrvatski	PL	Polski	UK	Українська
DE	Deutsch	HU	Magyar	PT	Português		
EN	English	IT	Italiano	RU	Русский		
ES	Español	JA	日本語	RO	Română		
ET	Eesti	LV	Latviešu	SA	العربية		
FI	Suomi	LT	Lietuvių	SK	Slovenčina		

CS	SP1580120	(User Manual TDX-SP2 - Czech)	1
DA	SP1580070	(User Manual TDX-SP2 - Danish)	1
DE	SP1578996	(User Manual TDX-SP2 - German)	1
EN	SP1578997	(User Manual TDX-SP2 - English)	1
ES	SP1580066	(User Manual TDX-SP2 - Spanish)	1
FI	SP1580069	(User Manual TDX-SP2 - Finnish)	1
FR	SP1580064	(User Manual TDX-SP2 - French)	1
IT	SP1580065	(User Manual TDX-SP2 - Italian)	1
LT	SP1580121	(User Manual TDX-SP2 - Lithuanian)	1
NL	SP1580067	(User Manual TDX-SP2 - Dutch)	1
NO	SP1580072	(User Manual TDX-SP2 - Norwegian)	1
PL	SP1580122	(User Manual TDX-SP2 - Polish)	1
PT	SP1580068	(User Manual TDX-SP2 - Portuguese)	1
RU	SP1580123	(User Manual TDX-SP2 - Russian)	1
SL	SP1580124	(User Manual TDX-SP2 - Slovenian)	1
SV	SP1580071	(User Manual TDX-SP2 - Swedish)	1

(User Manuals Ultra Low Maxx Seat)

languages.iso

BG	Български	FR	Français	NL	Nederlands	SL	Slovenščina
CS	Čeština	HE	עברית	NO	Norsk	SV	Svenska
DA	Dansk	HR	Hrvatski	PL	Polski	UK	Українська
DE	Deutsch	HU	Magyar	PT	Português		
EN	English	IT	Italiano	RU	Русский		
ES	Español	JA	日本語	RO	Română		
ET	Eesti	LV	Latviešu	SA	العربية		
FI	Suomi	LT	Lietuvių	SK	Slovenčina		

DA	SP1585726	(User Manual Ultra Low Maxx Seat - Danish)	1
DE	SP1585727	(User Manual Ultra Low Maxx Seat - German)	1
EN	SP1585725	(User Manual Ultra Low Maxx Seat - English)	1
ES	SP1601887	(User Manual Ultra Low Maxx Seat - Spanish)	1
FI	SP1585728	(User Manual Ultra Low Maxx Seat - Finnish)	1
FR	SP1585729	(User Manual Ultra Low Maxx Seat - French)	1
IT	SP1588266	(User Manual Ultra Low Maxx Seat - Italian)	1
NL	SP1585730	(User Manual Ultra Low Maxx Seat - Dutch)	1
NO	SP1585731	(User Manual Ultra Low Maxx Seat - Norwegian)	1
PT	SP1585732	(User Manual Ultra Low Maxx Seat - Portuguese)	1
SV	SP1585733	(User Manual Ultra Low Maxx Seat - Swedish)	1

(User Manuals Shark2)

languages.iso

BG	Български	FR	Français	NL	Nederlands	SL	Slovenščina
CS	Čeština	HE	עברית	NO	Norsk	SV	Svenska
DA	Dansk	HR	Hrvatski	PL	Polski	UK	Українська
DE	Deutsch	HU	Magyar	PT	Português		
EN	English	IT	Italiano	RU	Русский		
ES	Español	JA	日本語	RO	Română		
ET	Eesti	LV	Latviešu	SA	العربية		
FI	Suomi	LT	Lietuvių	SK	Slovenčina		

DA	1451349	(User Manual Shark II Danish)	1
PT	1451353	(User Manual Shark II Portuguese)	1
DE	1451346	(User Manual Shark II German)	1
NO	1451350	(User Manual Shark II Norwegian)	1
EN	1451356	(User Manual Shark II English)	1
FR	1451347	(User Manual Shark II French)	1
IT	1451354	(User Manual Shark II Italian)	1
SV	1451348	(User Manual Shark II Swedish)	1
FI	1451355	(User Manual Shark II Finnish)	1
NL	1451351	(User Manual Shark II Dutch)	1
SL	1451584	(User Manual Shark II Slovenian)	1
ES	1451352	(User Manual Shark II Spanish)	1
BG	SP1546061	(User Manual Shark II Bulgarian)	1
CS	SP1533297	(User Manual Shark2 CS)	1
ET	SP1533298	(User Manual Shark2 EE)	1
HR	SP1574028	(User Manual Shark2 Croatian)	1
HU	SP1533299	(User Manual Shark2 HU)	1
LV	SP1576075	(User Manual Shark2 Latvian)	1
LT	SP1533300	(User Manual Shark2 LT)	1
PL	SP1533301	(User Manual Shark2 PL)	1
RU	SP1570186	(User Manual Shark2 Russian)	1

(User Manuals REM A & B)

languages.iso

BG	Български	FR	Français	NL	Nederlands	SL	Slovenščina
CS	Čeština	HE	עברית	NO	Norsk	SV	Svenska
DA	Dansk	HR	Hrvatski	PL	Polski	UK	Українська
DE	Deutsch	HU	Magyar	PT	Português		
EN	English	IT	Italiano	RU	Русский		
ES	Español	JA	日本語	RO	Română		
ET	Eesti	LV	Latviešu	SA	العربية		
FI	Suomi	LT	Lietuvių	SK	Slovenčina		

IT	1529230	(User Manual REM A & B - Italian)	1
NO	1529232	(User Manual REM A & B - Norwegian)	1
ES	1529174	(User Manual REM A & B - Spanish)	1
FI	1529228	(User Manual REM A & B - Finnish)	1
EN	1529173	(User Manual REM A & B - English)	1
DA	1529171	(User Manual REM A & B - Danish)	1
NL	1529231	(User Manual REM A & B - Dutch)	1
DE	1529172	(User Manual REM A & B - German)	1
PT	1529233	(User Manual REM A & B - Portuguese)	1
SV	1529234	(User Manual REM A & B Swedish)	1
FR	1529229	(User Manual REM A & B - French)	1
CS	SP1533259	(User Manual REM A & B - Czech)	1
ET	SP1533265	(User Manual REM A & B - Estonian)	1
HR	SP1533258	(User Manual REM A & B - Croatian)	1
HU	SP1533266	(User Manual REM A & B - Hungarian)	1
LT	SP1533267	(User Manual REM A & B - Lithuanian)	1
PL	SP1533268	(User Manual REM A & B - Polish)	1
RO	SP1533270	(User Manual REM A & B - Romanian)	1
SL	SP1533269	(User Manual REM A & B - Slovenian)	1

(User Manuals REM 24 SD)

languages.iso

BG	Български	FR	Français	NL	Nederlands	SL	Slovenščina
CS	Čeština	HE	עברית	NO	Norsk	SV	Svenska
DA	Dansk	HR	Hrvatski	PL	Polski	UK	Українська
DE	Deutsch	HU	Magyar	PT	Português		
EN	English	IT	Italiano	RU	Русский		
ES	Español	JA	日本語	RO	Română		
ET	Eesti	LV	Latviešu	SA	العربية		
FI	Suomi	LT	Lietuvių	SK	Slovenčina		

CS	SP1533252	(User Manual REM24SD Czech)	1
DA	1436146	(User Manual REM24SD Danish)	1
DE	1436142	(User Manual REM24SD German)	1
EN	1436144	(User Manual REM24SD English)	1
ES	1436149	(User Manual REM24SD Spanish)	1
ET	SP1533253	(User Manual REM24SD Estonian)	1
FI	1436386	(User Manual REM24SD Finnish)	1
FR	1436143	(User Manual REM24SD French)	1
HR	SP1533251	(User Manual REM24SD Croatian)	1
HU	SP1533254	(User Manual REM24SD Hungarian)	1
IT	1436385	(User Manual REM24SD Italian)	1
LV	SP1543152	(User Manual REM24SD Latvian)	1
LT	SP1533255	(User Manual REM24SD Lithuanian)	1
NL	1436148	(User Manual REM24SD Dutch)	1
NO	1436147	(User Manual REM24SD Norwegian)	1
PL	SP1533256	(User Manual REM24SD Polish)	1
PT	1436150	(User Manual REM24SD Portuguese)	1
SL	SP1533257	(User Manual REM24SD Slovenian)	1
SV	1436145	(User Manual REM24SD Swedish)	1

(User Manuals REM 550)

languages.iso

BG	Български	FR	Français	NL	Nederlands	SL	Slovenščina
CS	Čeština	HE	עברית	NO	Norsk	SV	Svenska
DA	Dansk	HR	Hrvatski	PL	Polski	UK	Українська
DE	Deutsch	HU	Magyar	PT	Português		
EN	English	IT	Italiano	RU	Русский		
ES	Español	JA	日本語	RO	Română		
ET	Eesti	LV	Latviešu	SA	العربية		
FI	Suomi	LT	Lietuvių	SK	Slovenčina		

IT	1528094	(User Manual REM550 Italiano)	1
NO	1528096	BRUKSANVISNING DX REM550 MANUAL NORSK	1
ES	1528091	(User Manual REM550 Spanish)	1
FI	1528092	(User Manual REM550 Finsk)	1
EN	1528090	BRUKSANVISNING DX-REM550 MANUAL	1
DA	1528088	(User Manual REM550 Danish)	1
NL	1528095	(User Manual REM550 Dutch)	1
DE	1528089	(User Manual REM550 German)	1
PT	1528097	(User Manual REM550 Portuguese)	1
SV	1528098	(User Manual REM550 Swedish)	1
FR	1528093	(User Manual REM550 French)	1
CS	SP1533280	(User Manual REM550 Czech)	1
ET	SP1533281	(User Manual REM550 Estonian)	1
HR	SP1533279	(User Manual REM550 Croatian)	1
HU	SP1533282	(User Manual REM550 Hungarian)	1
LT	SP1533283	(User Manual REM550 Lithuanian)	1
PL	SP1533284	(User Manual REM550 Polish)	1
SL	SP1533285	(User Manual REM550 Slovenian)	1

(User Manuals G90A)

languages.iso

BG	Български	FR	Français	NL	Nederlands	SL	Slovenščina
CS	Čeština	HE	עברית	NO	Norsk	SV	Svenska
DA	Dansk	HR	Hrvatski	PL	Polski	UK	Українська
DE	Deutsch	HU	Magyar	PT	Português		
EN	English	IT	Italiano	RU	Русский		
ES	Español	JA	日本語	RO	Română		
ET	Eesti	LV	Latviešu	SA	العربية		
FI	Suomi	LT	Lietuvių	SK	Slovenčina		

EN	SP1524441	(User Manual G90A - English)	1
DE	SP1536742	(User Manual G90A - German)	1
ES	SP1638870	(User Manual G90A - Spanish)	1
FR	SP1638871	(User Manual G90A - French)	1
NO	SP1638872	(User Manual G90A - Norwegian)	1
PT	SP1638873	(User Manual G90A - Portuguese)	1
SV	SP1638874	(User Manual G90A - Swedish)	1
DA	SP1638875	(User Manual G90A - Danish)	1
NL	SP1638876	(User Manual G90A - Dutch)	1

(User Manuals G91S)

languages.iso

BG	Български	FR	Français	NL	Nederlands	SL	Slovenščina
CS	Čeština	HE	עברית	NO	Norsk	SV	Svenska
DA	Dansk	HR	Hrvatski	PL	Polski	UK	Українська
DE	Deutsch	HU	Magyar	PT	Português		
EN	English	IT	Italiano	RU	Русский		
ES	Español	JA	日本語	RO	Română		
ET	Eesti	LV	Latviešu	SA	العربية		
FI	Suomi	LT	Lietuvių	SK	Slovenčina		

(User Manuals Chin Control (2015-10-01 =>))

languages.iso

BG	Български	FR	Français	NL	Nederlands	SL	Slovenščina
CS	Čeština	HE	עברית	NO	Norsk	SV	Svenska
DA	Dansk	HR	Hrvatski	PL	Polski	UK	Українська
DE	Deutsch	HU	Magyar	PT	Português		
EN	English	IT	Italiano	RU	Русский		
ES	Español	JA	日本語	RO	Română		
ET	Eesti	LV	Latviešu	SA	العربية		
FI	Suomi	LT	Lietuvių	SK	Slovenčina		

CS	SP1589326	(User Manual Chin Control - Czech)	1
DA	SP1589394	(User Manual Chin Control - Danish)	1
DE	SP1589395	(User Manual Chin Control - German)	1
EN	SP1589325	(User Manual Chin Control - English)	1
ES	SP1589535	(User Manual Chin Control - Spanish)	1
FI	SP1589436	(User Manual Chin Control - Finnish)	1
FR	SP1589437	(User Manual Chin Control - French)	1
HR	SP1589438	(User Manual Chin Control - Croatian)	1
IT	SP1589439	(User Manual Chin Control - Italian)	1
LT	SP1589440	(User Manual Chin Control - Lithuanian)	1
NL	SP1589441	(User Manual Chin Control - Dutch)	1
NO	SP1589442	(User Manual Chin Control - Norwegian)	1
PL	SP1589443	(User Manual Chin Control - Polish)	1
PT	SP1589444	(User Manual Chin Control - Portuguese)	1
SL	SP1589445	(User Manual Chin Control - Slovenian)	1
SV	SP1589465	(User Manual Chin Control - Swedish)	1

(User Manuals REM 34 & REM 41 Tray/Centre Control)

languages.iso

BG	Български	FR	Français	NL	Nederlands	SL	Slovenščina
CS	Čeština	HE	עברית	NO	Norsk	SV	Svenska
DA	Dansk	HR	Hrvatski	PL	Polski	UK	Українська
DE	Deutsch	HU	Magyar	PT	Português		
EN	English	IT	Italiano	RU	Русский		
ES	Español	JA	日本語	RO	Română		
ET	Eesti	LV	Latviešu	SA	العربية		
FI	Suomi	LT	Lietuvių	SK	Slovenčina		

CS	SP1538031	(User Manual REM41 Tray Control Czech)	1
DA	SP1538022	(User Manual REM41 Tray Control Danish)	1
DE	SP1541402	(User Manual REM41 Tray Control German)	1
EN	SP1538023	(User Manual REM41 Tray Control English)	1
ES	SP1538024	(User Manual REM41 Tray Control Spanish)	1
FR	SP1538025	(User Manual REM41 Tray Control French)	1
IT	SP1538026	(User Manual REM41 Tray Control Italian)	1
NL	SP1538027	(User Manual REM41 Tray Control Dutch)	1
NO	SP1538028	(User Manual REM41 Tray Control Norwegian)	1
PT	SP1538029	(User Manual REM41 Tray Control Portuguese)	1
SL	SP1538032	(User Manual REM41 Tray Control Slovenian)	1
SV	SP1538030	(User Manual REM41 Tray Control Swedish)	1

(User Manuals IDC (Intuitive Dual Control))

languages.iso

BG	Български	FR	Français	NL	Nederlands	SL	Slovenščina
CS	Čeština	HE	עברית	NO	Norsk	SV	Svenska
DA	Dansk	HR	Hrvatski	PL	Polski	UK	Українська
DE	Deutsch	HU	Magyar	PT	Português		
EN	English	IT	Italiano	RU	Русский		
ES	Español	JA	日本語	RO	Română		
ET	Eesti	LV	Latviešu	SA	العربية		
FI	Suomi	LT	Lietuvių	SK	Slovenčina		

DE	1521480	(User Manual IDC - German)	1
EN	1521481	(User Manual IDC - English)	1
FR	1521482	(User Manual IDC - French)	1
NL	1521483	(User Manual IDC - Dutch)	1
NO	1521484	BRUKSANVISNING INTUITIV LEDSAGERSTYRING	1
SV	1521485	(User Manual IDC - Swedish)	1
DA	1521486	(User Manual IDC - Danish)	1
IT	1522252	(User Manual IDC - Italian)	1
PT	1522254	(User Manual IDC - Portuguese)	1
ES	1522255	(User Manual IDC - Spanish)	1

(User Manuals Compact Dual Control)

languages.iso

BG	Български	FR	Français	NL	Nederlands	SL	Slovenščina
CS	Čeština	HE	עברית	NO	Norsk	SV	Svenska
DA	Dansk	HR	Hrvatski	PL	Polski	UK	Українська
DE	Deutsch	HU	Magyar	PT	Português		
EN	English	IT	Italiano	RU	Русский		
ES	Español	JA	日本語	RO	Română		
ET	Eesti	LV	Latviešu	SA	العربية		
FI	Suomi	LT	Lietuvių	SK	Slovenčina		

BG	SP1589186	(User Manual Compact Dual Control - Bulgarian)	1
CS	SP1589313	(User Manual Compact Dual Control - Czech)	1
DA	SP1588488	(User Manual Compact Dual Control - Danish)	1
DE	SP1588489	(User Manual Compact Dual Control - German)	1
EN	SP1588462	(User Manual Compact Dual Control - English)	1
ES	SP1588496	(User Manual Compact Dual Control - Spanish)	1
ET	SP1589311	(User Manual Compact Dual Control - Estonian)	1
FI	SP1588490	(User Manual Compact Dual Control - Finnish)	1
FR	SP1588491	(User Manual Compact Dual Control - French)	1
HR	SP1589314	(User Manual Compact Dual Control - Croatian)	1
HU	SP1589315	(User Manual Compact Dual Control - Hungarian)	1
IT	SP1588492	(User Manual Compact Dual Control - Italian)	1
LV	SP1589317	(User Manual Compact Dual Control - Latvian)	1
LT	SP1589316	(User Manual Compact Dual Control - Lithuanian)	1
NL	SP1588493	(User Manual Compact Dual Control - Dutch)	1
NO	SP1588494	(User Manual Compact Dual Control - Norwegian)	1
PL	SP1589318	(User Manual Compact Dual Control - Polish)	1
PT	SP1588495	(User Manual Compact Dual Control - Portuguese)	1
RU	SP1589320	(User Manual Compact Dual Control - Russian)	1
RO	SP1589319	(User Manual Compact Dual Control - Romanian)	1
SK	SP1589321	(User Manual Compact Dual Control - Slovak)	1
SL	SP1589322	(User Manual Compact Dual Control - Slovenian)	1
SV	SP1588497	(User Manual Compact Dual Control - Swedish)	1
UK	SP1589323	(User Manual Compact Dual Control - Ukrainian)	1

(User Manuals ACS Handprogrammer)

languages.iso

BG	Български	FR	Français	NL	Nederlands	SL	Slovenščina
CS	Čeština	HE	עברית	NO	Norsk	SV	Svenska
DA	Dansk	HR	Hrvatski	PL	Polski	UK	Українська
DE	Deutsch	HU	Magyar	PT	Português		
EN	English	IT	Italiano	RU	Русский		
ES	Español	JA	日本語	RO	Română		
ET	Eesti	LV	Latviešu	SA	العربية		
FI	Suomi	LT	Lietuvių	SK	Slovenčina		

DE	1419316	(User Manual Handheld Programmer German)	1
EN	1419317	(User Manual Handheld Programmer English)	1

(User Manual Proximity Head Remote DSL-105)

languages.iso

BG	Български	FR	Français	NL	Nederlands	SL	Slovenščina
CS	Čeština	HE	עברית	NO	Norsk	SV	Svenska
DA	Dansk	HR	Hrvatski	PL	Polski	UK	Українська
DE	Deutsch	HU	Magyar	PT	Português		
EN	English	IT	Italiano	RU	Русский		
ES	Español	JA	日本語	RO	Română		
ET	Eesti	LV	Latviešu	SA	العربية		
FI	Suomi	LT	Lietuvių	SK	Slovenčina		

EN	SP1555731	(User Manual Proximity Head Remote DSL-105 English)	1
DE	SP1555732	(User Manual Proximity Head Remote DSL-105 German)	1
FR	SP1555733	(User Manual Proximity Head Remote DSL-105 French)	1
DA	SP1555839	(User Manual Proximity Head Remote DSL-105 Danish)	1
ES	SP1555840	(User Manual Proximity Head Remote DSL-105 Spanish)	1
FI	SP1555841	(User Manual Proximity Head Remote DSL-105 Finnish)	1
IT	SP1555842	(User Manual Proximity Head Remote DSL-105 Italian)	1
NL	SP1555843	(User Manual Proximity Head Remote DSL-105 Dutch)	1
NO	SP1555844	(User Manual Proximity Head Remote DSL-105 Norwegian)	1
PT	SP1555845	(User Manual Proximity Head Remote DSL-105 Portuguese)	1
SV	SP1555846	(User Manual Proximity Head Remote DSL-105 Swedish)	1

(User Manuals Headcontrol)

languages.iso

BG	Български	FR	Français	NL	Nederlands	SL	Slovenščina
CS	Čeština	HE	עברית	NO	Norsk	SV	Svenska
DA	Dansk	HR	Hrvatski	PL	Polski	UK	Українська
DE	Deutsch	HU	Magyar	PT	Português		
EN	English	IT	Italiano	RU	Русский		
ES	Español	JA	日本語	RO	Română		
ET	Eesti	LV	Latviešu	SA	العربية		
FI	Suomi	LT	Lietuvių	SK	Slovenčina		

DE	SP1493028	(User Manual Headrest Headset German)	1
EN	SP1493029	(User Manual Headrest Headset English)	1
NL	SP1493030	(User Manual Headrest Headset Dutch)	1

(User Manuals iPortal)

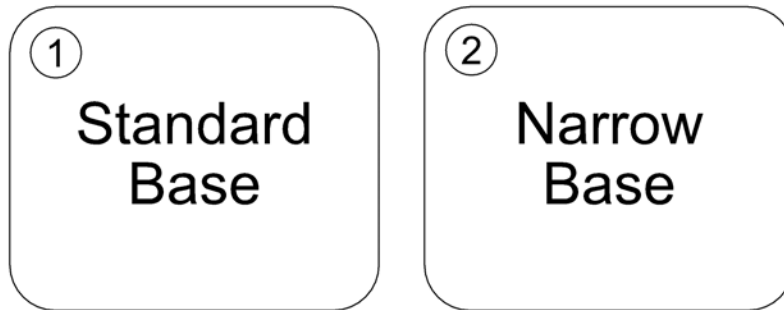
languages.iso

BG	Български	FR	Français	NL	Nederlands	SL	Slovenščina
CS	Čeština	HE	עברית	NO	Norsk	SV	Svenska
DA	Dansk	HR	Hrvatski	PL	Polski	UK	Українська
DE	Deutsch	HU	Magyar	PT	Português		
EN	English	IT	Italiano	RU	Русский		
ES	Español	JA	日本語	RO	Română		
ET	Eesti	LV	Latviešu	SA	العربية		
FI	Suomi	LT	Lietuvių	SK	Slovenčina		

DA	SP1552493	(User Manual iPortal danisch)	1
DE	SP1552535	(User Manual iPortal german)	1
EN	SP1548231	(User Manual iPortal english)	1
ES	SP1568026	(User Manual iPortal Spanish)	1
FI	SP1552512	(User Manual iPortal finnish)	1
FR	SP1552494	(User Manual iPortal french)	1
IT	SP1552534	(User Manual iPortal Italian)	1
NL	SP1559807	(User Manual iPortal Dutch)	1
NO	SP1552513	(User Manual iPortal norwegian)	1
PT	SP1567967	(User Manual iPortal Portuguese)	1
SL	SP1552703	(User Manual iPortal slovenian)	1
SV	SP1552533	(User Manual iPortal swedish)	1

Chassis

SPCTDX_SP210.iso




1	-	(Chassis Standard Base)	1
2	-	(Chassis Narrow Base)	1

(Chassis Standard Base)

-	(Battery Box Hardware)	1
-	(Walking Beam)	1
-	(Rear Pivot Standard Base)	1
-	(Shrouds Standard Base)	1
-	(Low Rider Option)	1

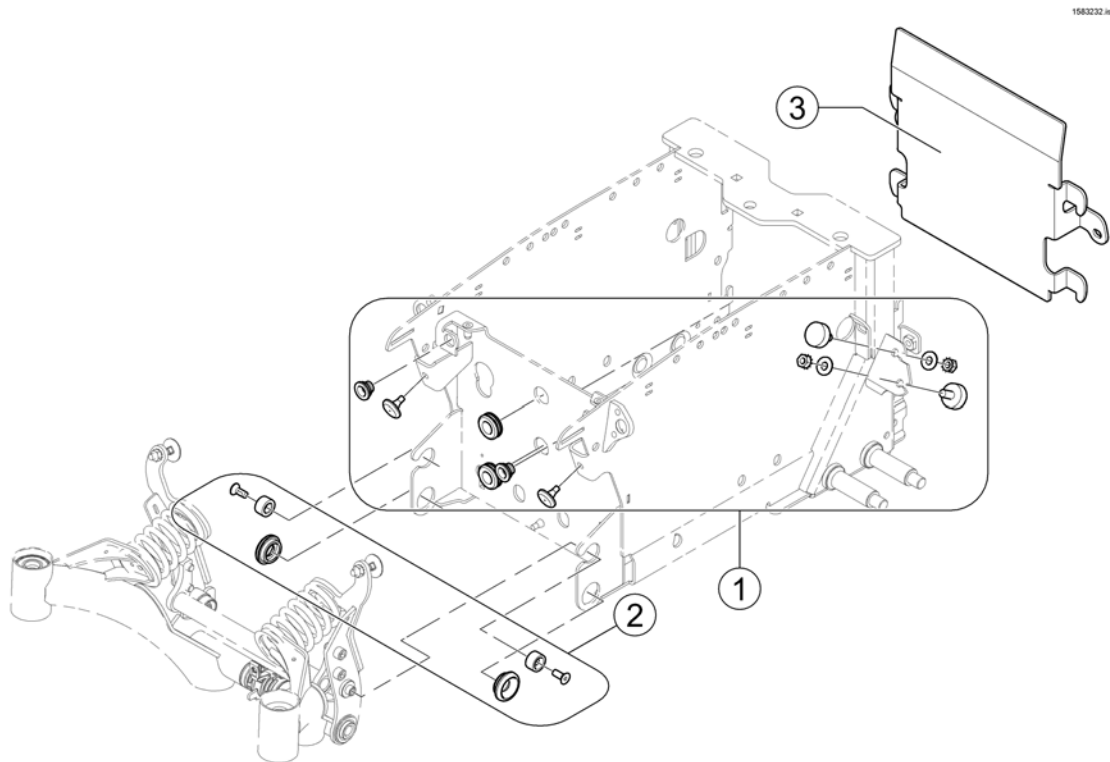
(Battery Box Hardware)



1	SP1523179	BATTERIBOKS TDX EXTRA GUMMIPAD OG FORINGER	1
2	SP1523193	FESTE HODESTØTTE BAKRE TDX	1
3	SPTRX8891	(Recessed Battery door for Anterior Tilt mounted Ultra Low Maxx Seat)	1
3	SP1589636	(Recessed Battery Door for Anterior Tilt mounted Ultra Low Maxx Seat & Low Rider Option)	1

(Walking Beam)

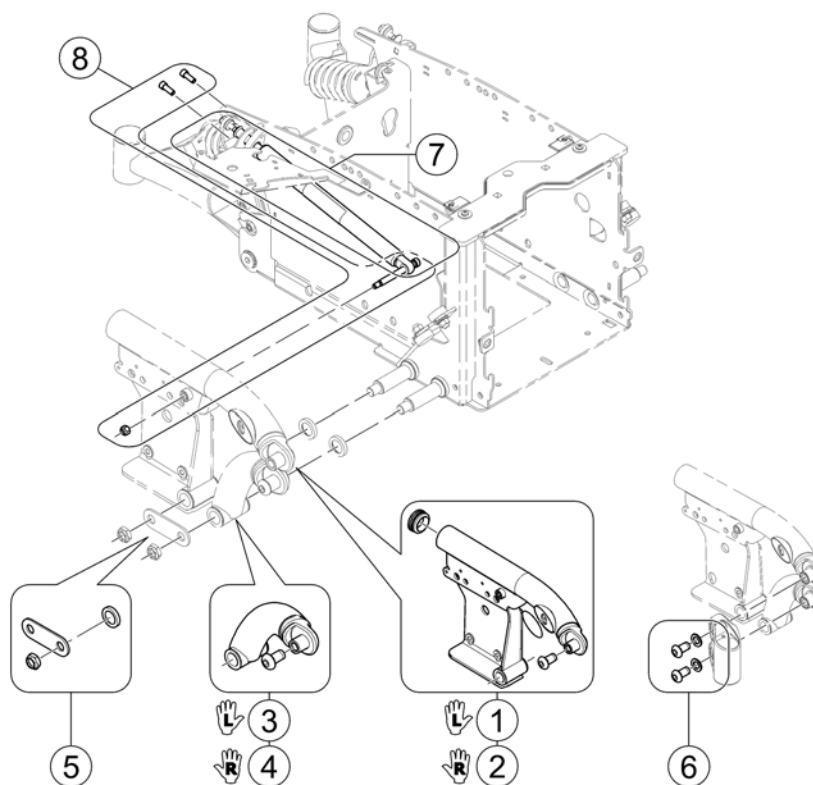
1583299.no



-	SP1142270	(Slide Bushing Kit (12 pcs.)) <i>It's recommended to replace all bearings at once</i>	1
1	SP1179092	(Kit, Walking Beam Assy Black, SSD LH)	1
2	SP1179091	(Kit, Walking Beam Assy Black, SSD RH)	1
3	SP1178814	(Link Assy, Lower Pivot Black LH)	1
4	SP1178813	()	1
5	SP1523168	(Kit, Lower Pivot Link / 4-Pole Motor Mounting Hardware)	1
6	SP1523176	SKRUER GAFFELHUS TDX	1
7	SP1523158	GASSFJÆR 180 KG TDX	1
8	SP1179253	MONTERINGSKIT GASSFJÆR TDX-SP2	1

(Rear Pivot Standard Base)

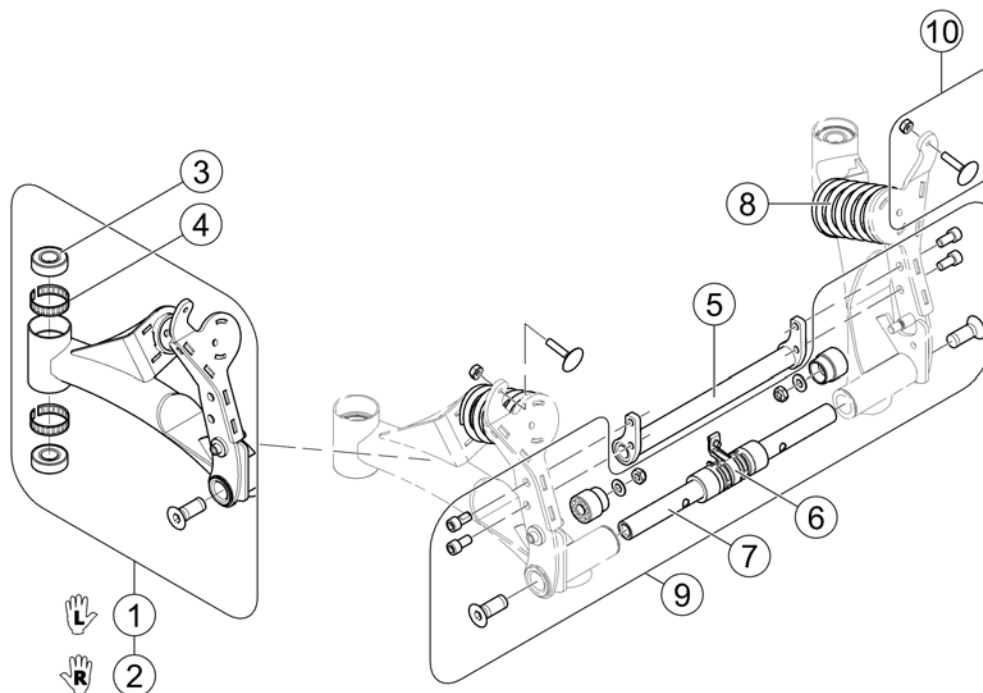
1583286.iss




1	SP1179094	(Kit, Rear Swing Arm Assy w/ Hardware, black, LH)	1
2	SP1179093	(Kit, Rear Swing Arm Assy w/ Hardware, black, RH)	1
3	SP40117X006	KULELAGER TDX SP2 BAK	1
4	SP1110717	(Ring, Tolerance)	1
5	SP1176365	(Cross strut, black - Standard Base)	1
6	SP1523166	FESTE DEKSEL TDX TDX	1
7	SP1155020	(Bar, Pivot - Standard Base)	1
8	SP1579643	(Spring - silver)	1
9	SP1523177	FESTESETT ROTASJONSAKSEL TDX	1
10	SP1151949	(Kit, Elevator Bolt w/ Keps Nut (1/4-20))	1

(Shrouds Standard Base)

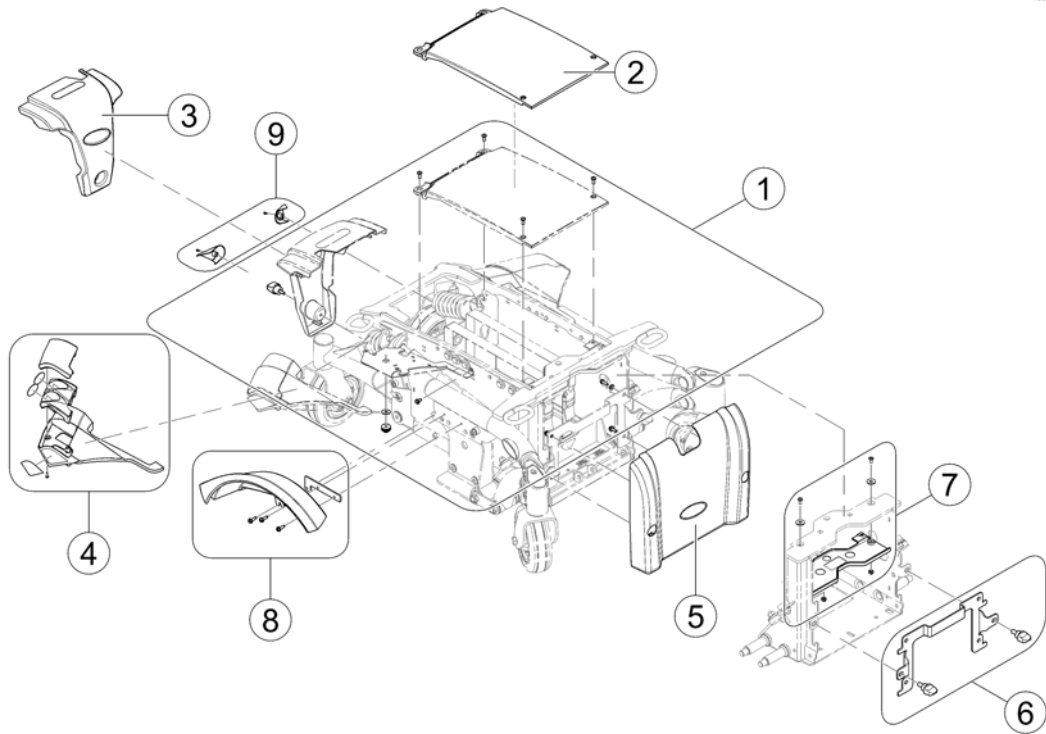1583238.00




1	SP1179326	(Kit, Shroud Hardware)	1
2	SP1174636	(Shroud, Center, black - Standard Base)	1
3	-	(Rear Shrouds - Standard Base)	1
4	-	(Bezel Shroud)	1
5	-	(Front Shroud - Standard Base)	1
6	SP1179898	(Front Shroud, Bracket & Hardware - Standard Base)	1
7	SP1179739	(Isolation Shroud Assy - Standard Base)	1
8	-	(Fender)	1
9	-	(Rear Caps)	1

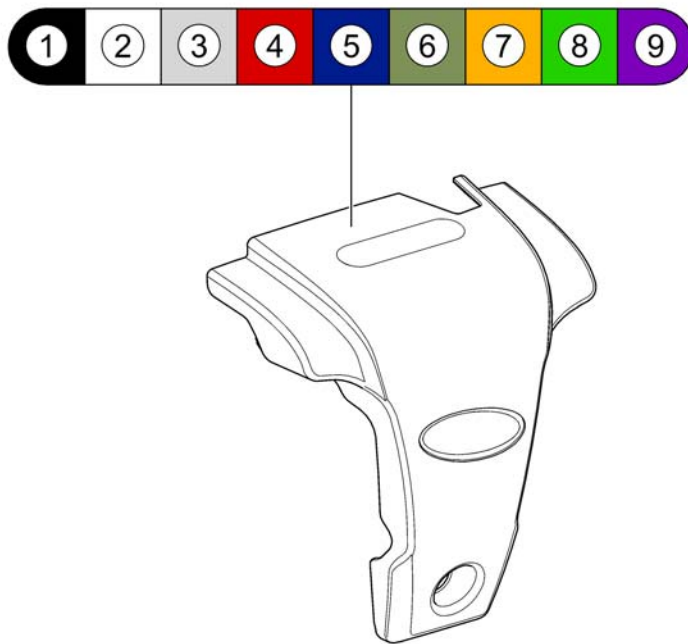


1	SP1176568	(Rear Shroud, Wet Black)		1
2	SP1176569	(Rear Shroud, Arctic White)		1
3	SP1176570	(Rear Shroud, Electric Silver)		1
4	SP1176571	(Rear Shroud, Electric Red)		1
5	SP1176572	(Rear Shroud, Electric Blue)		1
6	SP1176573	(Rear Shroud, Aviator Green)	→ 30.06.2015	1
7	SP1176574	(Rear Shroud, Desert Orange)	→ 30.06.2015	1
8	SP1588787	(Rear Shroud, Jade Green)	01.07.2015 →	1
9	SP1588850	(Rear Shroud, Transparent Purple)	01.07.2015 →	1

(Bezel Shroud)

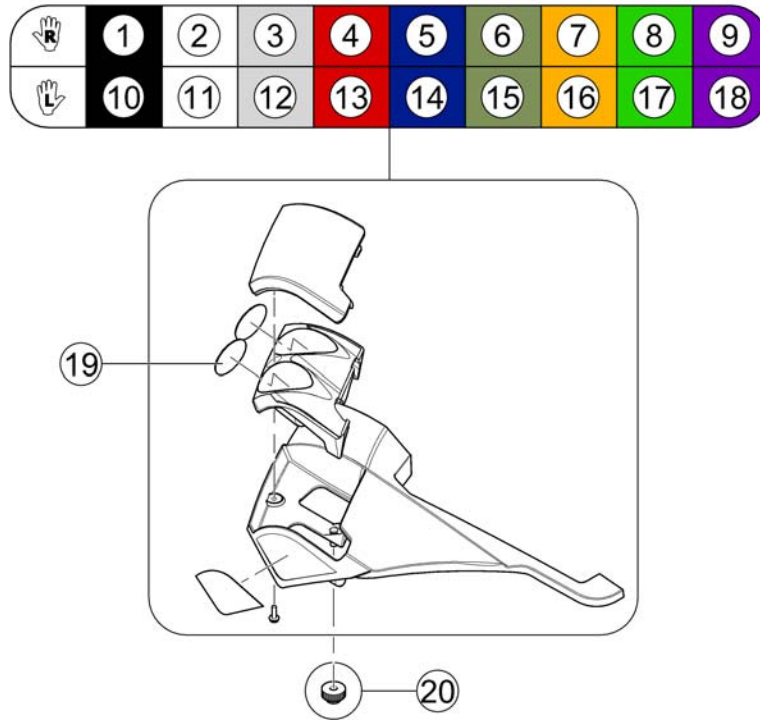1583256 Rev




1	SP1179237	(Coloured Bezel Shroud RH w/ Mount Wet Black)		1
2	SP1179238	(Coloured Bezel Shroud RH w/ Mount Artic White)		1
3	SP1179239	(Coloured Bezel Shroud RH w/ Mount Electric Silver)		1
4	SP1179240	(Coloured Bezel Shroud RH w/ Mount Electric Red)		1
5	SP1179241	(Coloured Bezel Shroud RH w/ Mount Electric Blue)		1
6	SP1179242	(Coloured Bezel Shroud RH w/ Mount Aviator Green)	→ 30.06.2015	1
7	SP1179243	(Coloured Bezel Shroud RH w/ Mount Desert Orange)	→ 30.06.2015	1
8	SP1588779	(Coloured Bezel Shroud RH w/ Mount Jade Green)	01.07.2015 →	1
9	SP1588780	(Coloured Bezel Shroud RH w/ Mount Transparent Purple)	01.07.2015 →	1
10	SP1179244	(Coloured Bezel Shroud LH w/ Mount Wet Black)		1
11	SP1179245	(Coloured Bezel Shroud LH w/ Mount Artic White)		1
12	SP1179246	(Coloured Bezel Shroud LH w/ Mount Electric Silver)		1
13	SP1179247	(Coloured Bezel Shroud LH w/ Mount Electric Red)		1
14	SP1179248	(Coloured Bezel Shroud LH w/ Mount Electric Blue)		1
15	SP1179249	(Coloured Bezel Shroud LH w/ Mount Aviator Green)	→ 30.06.2015	1
16	SP1179250	(Coloured Bezel Shroud LH w/ Mount Desert Orange)	→ 30.06.2015	1
17	SP1588781	(Coloured Bezel Shroud LH w/ Mount Jade Green)	01.07.2015 →	1
18	SP1588782	(Coloured Bezel Shroud LH w/ Mount Transparent Purple)	01.07.2015 →	1
19	SP1540286	(Reflector-Kit (4pcs), red, Dia. 28mm)		1
20	SP1607356	(Knurled Nut, plastic M5)		1

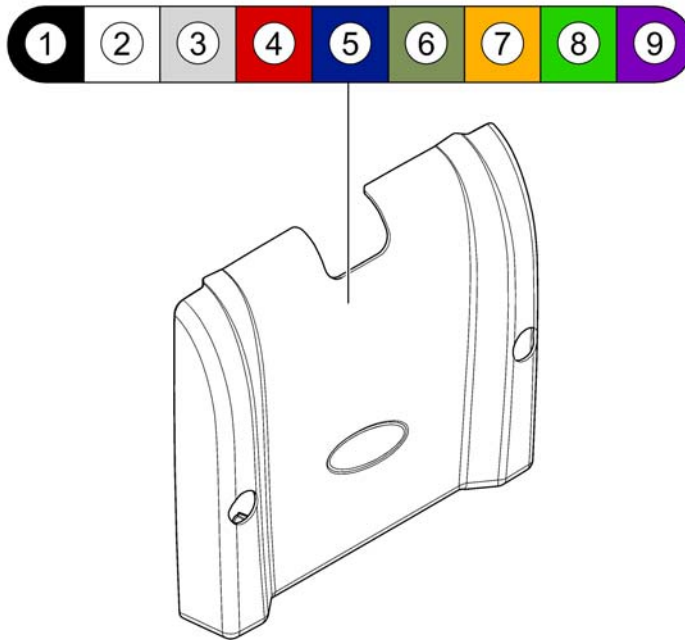

1	SP1179717	(Front Shroud Assy, Wet Black)		1
2	SP1189505	(Front Shroud Assy, Artic White)		1
3	SP1189506	(Front Shroud Assy, Electric Silver)		1
4	SP1179718	(Front Shroud Assy, Electric Red)		1
5	SP1189507	(Front Shroud Assy, Electric Blue)		1
6	SP1179719	(Front Shroud Assy, Aviator Green)	→ 30.06.2015	1
7	SP1189508	(Front Shroud Assy, Desert Orange)	→ 30.06.2015	1
8	SP1588784	(Front Shroud Assy, Jade Green)	01.07.2015 →	1
9	SP1588786	(Front Shroud Assy, Transparent Purple)	01.07.2015 →	1

(Fender)

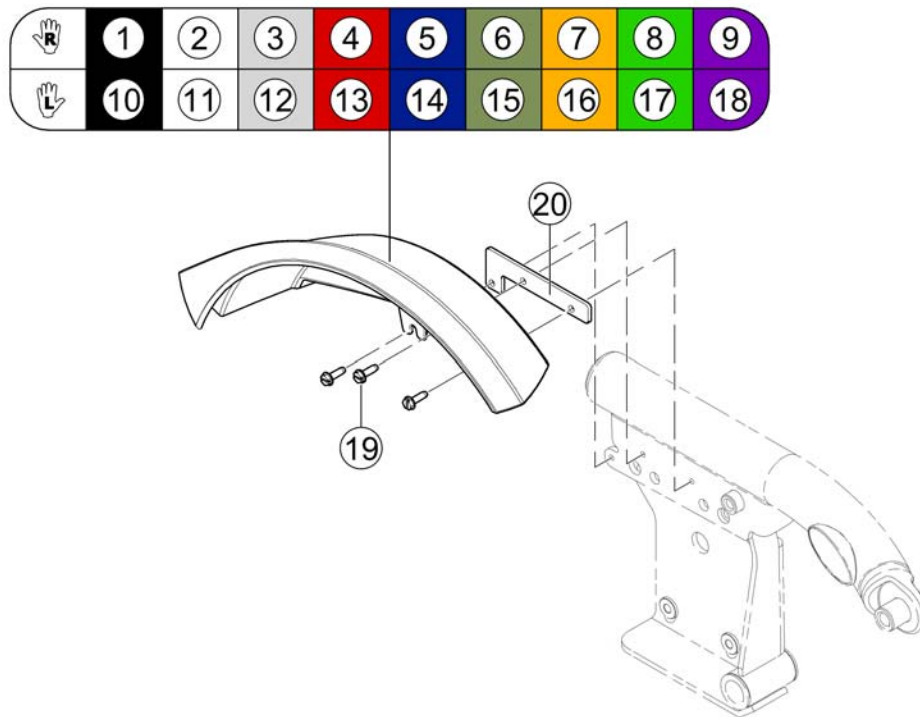
1583236 Rev



1	SP1176696	(Fender Shroud Assy, Wet Black, RH)		1
2	SP1176697	(Fender Shroud Assy, Artic White, RH)		1
3	SP1176698	SKJERM HØ TDX KPL SØLV TDX SP2		1
4	SP1176699	(Fender Shroud Assy, Electric Red, RH)		1
5	SP1176700	(Fender Shroud Assy, Electric Blue, RH)		1
6	SP1176701	(Fender Shroud Assy, Aviator Green, RH)	→ 30.06.2015	1
7	SP1176702	(Fender Shroud Assy, Desert Orange, RH)	→ 30.06.2015	1
8	SP1588737	(Fender Shroud Assy, Jade Green, RH)	01.07.2015 →	1
9	SP1588738	(Fender Shroud Assy, Transparent Purple, RH)	01.07.2015 →	1
10	SP1176703	(Fender Shroud Assy, Wet Black, LH)		1
11	SP1176704	(Fender Shroud Assy, Artic White, LH)		1
12	SP1176705	SKJERM VE TDX KPL SØLV TDX SP2		1
13	SP1176706	(Fender Shroud Assy, Electric Red, LH)		1
14	SP1176707	(Fender Shroud Assy, Electric Blue, LH)		1
15	SP1176708	(Fender Shroud Assy, Aviator Green, LH)	→ 30.06.2015	1
16	SP1176709	(Fender Shroud Assy, Desert Orange, LH)	→ 30.06.2015	1
17	SP1588739	(Fender Shroud Assy, Jade Green, LH)	01.07.2015 →	1
18	SP1588740	(Fender Shroud Assy, Transparent Purple, LH)	01.07.2015 →	1
19	SP1025451	(Act.Tech Screw HexHead10-32)		3
20	SP1179913	(Spacer, Shroud, Black)		1

(Rear Caps)

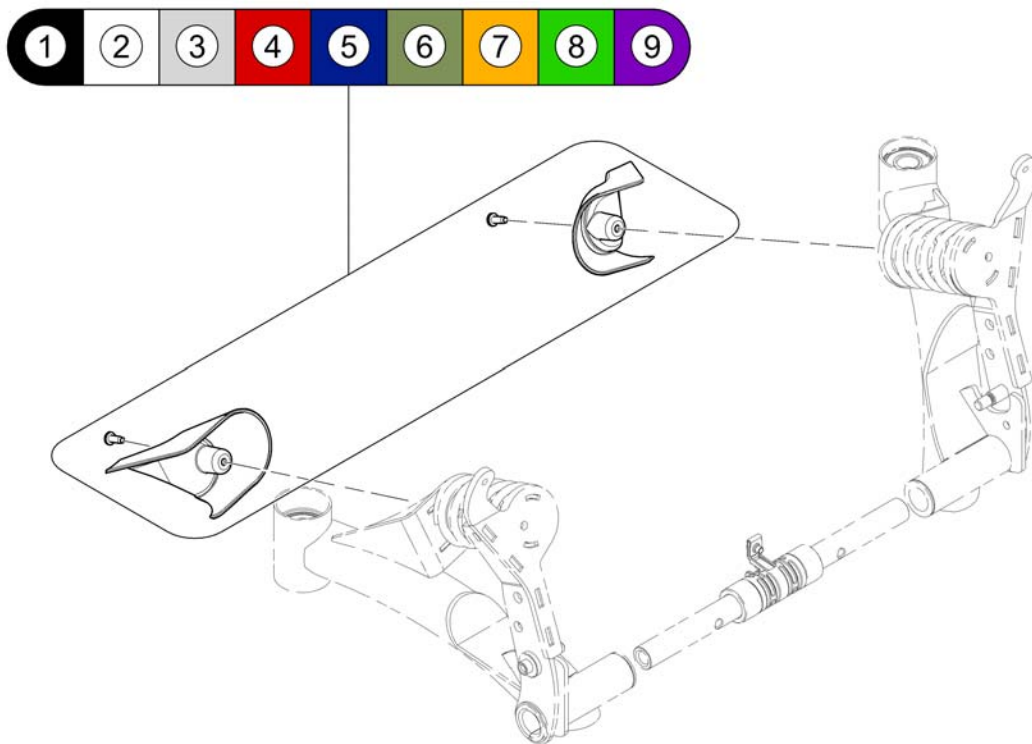
1583287 Rev



1	SP1179084	(Kit, Coloured Rear Caps w/ Mount Wet Black)		1
2	SP1179085	(Kit, Coloured Rear Caps w/ Mount Artic White)		1
3	SP1179086	(Kit, Coloured Rear Caps w/ Mount Electric Silver)		1
4	SP1179087	(Kit, Coloured Rear Caps w/ Mount Electric Red)		1
5	SP1179088	(Kit, Coloured Rear Caps w/ Mount Electric Blue)		1
6	SP1179089	(Kit, Coloured Rear Caps w/ Mount Aviator Green)	→ 30.06.2015	1
7	SP1179090	(Kit, Coloured Rear Caps w/ Mount Desert Orange)	→ 30.06.2015	1
8	SP1588741	(Kit, Coloured Rear Caps w/ Mount Jade Green)	01.07.2015 →	1
9	SP1588778	(Kit, Coloured Rear Caps w/ Mount Transparent Purple)	01.07.2015 →	1

(Low Rider Option)

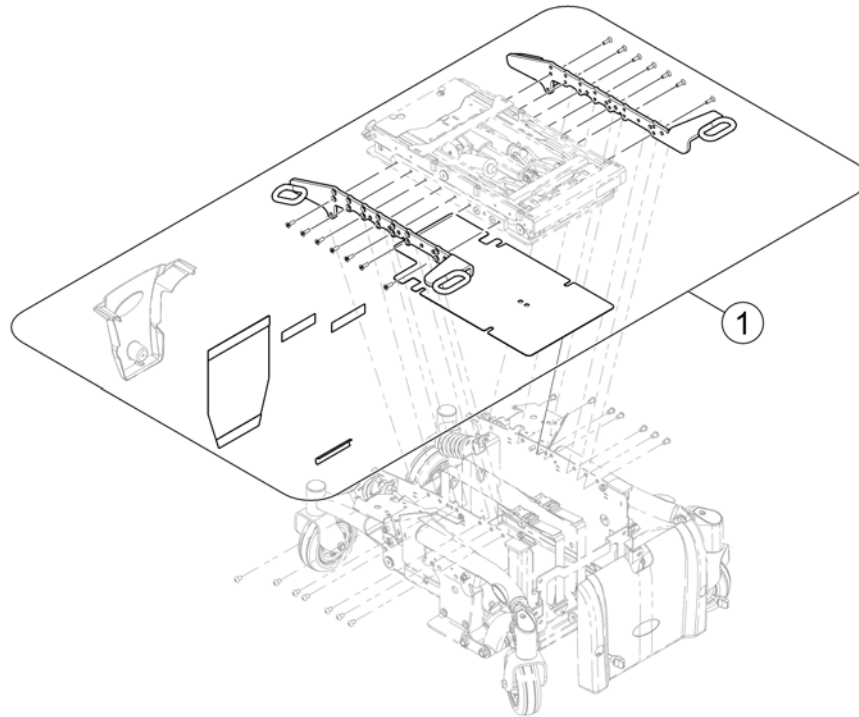1587261 Rev




-	(The Low Rider Kit may not be used for retrofitting!)	0
1	SP1587260 LOW-RIDER KIT KPL TDX-SP2	1
1	SP1590597 (Low Rider Kit - for Ultra Low Maxx Seat)	1

(Chassis Narrow Base)

-	(Battery Box Hardware)	0
-	(Walking Beam (Narrow Base))	1
-	(Rear Pivot Narrow Base)	1
-	(Shrouds Narrow Base)	1

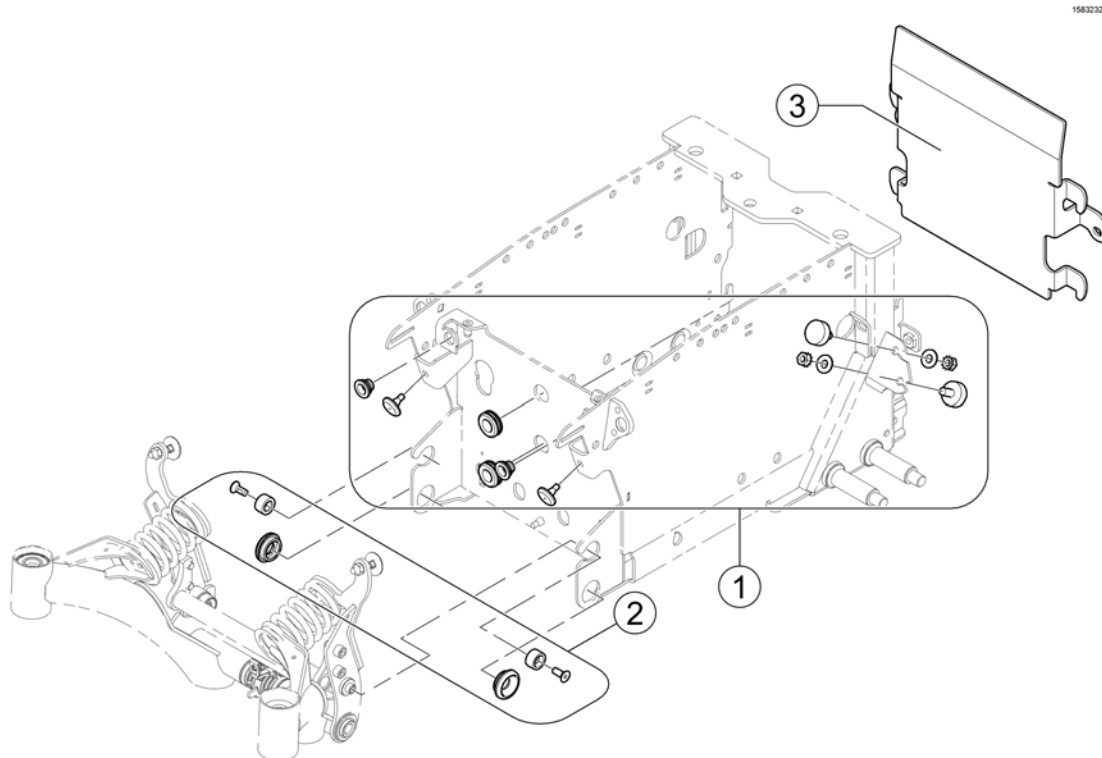(Battery Box Hardware)




1583232.000

1	SP1523179	BATTERIBOKS TDX EXTRA GUMMIPAD OG FORINGER	1
2	SP1523193	FESTE HODESTØTTE BAKRE TDX	1
3	SPTRX8890	(Recessed Battery door for Anterior Tilt mounted Ultra Low Maxx Seat)	1

(Walking Beam (Narrow Base))

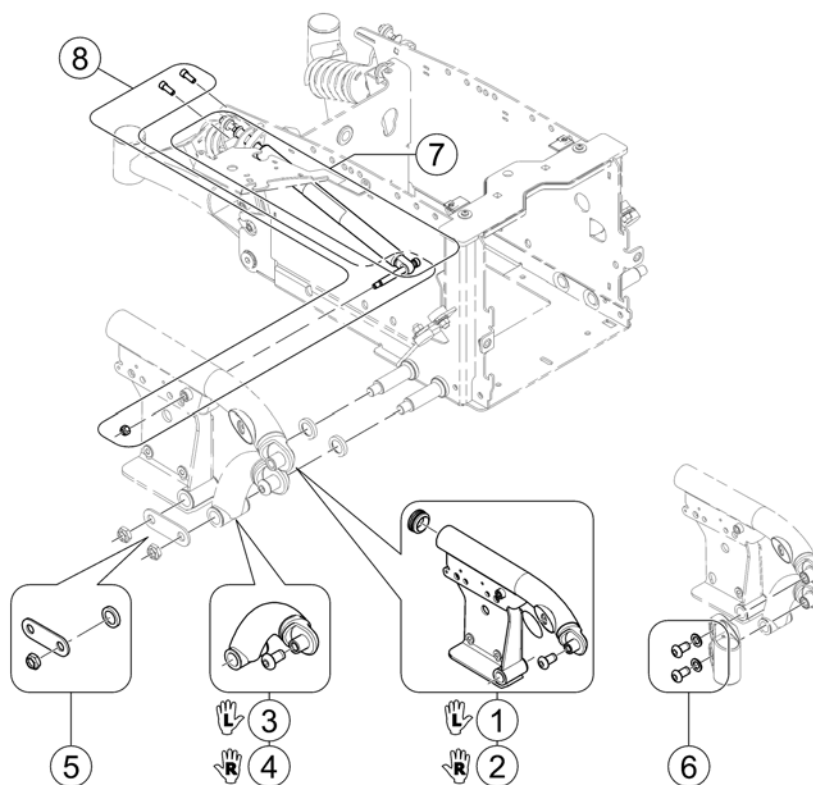
1583259.000



-	SP1142270	(Slide Bushing Kit (12 pcs.)) <i>It's recommended to replace all bearings at once</i>	1
1	SP1179092	(Kit, Walking Beam Assy Black, SSD LH)	1
2	SP1179091	(Kit, Walking Beam Assy Black, SSD RH)	1
3	SP1178814	(Link Assy, Lower Pivot Black LH)	1
4	SP1178813	()	1
5	SP1523168	(Kit, Lower Pivot Link / 4-Pole Motor Mounting Hardware)	1
6	SP1523176	SKRUER GAFFELHUS TDX	1
7	SP1523162	(Cylinder, Locking Gas (300lb / 135kg user weight limit))	1
8	SP1179253	MONTERINGSKIT GASSFJÆR TDX-SP2	1

(Rear Pivot Narrow Base)

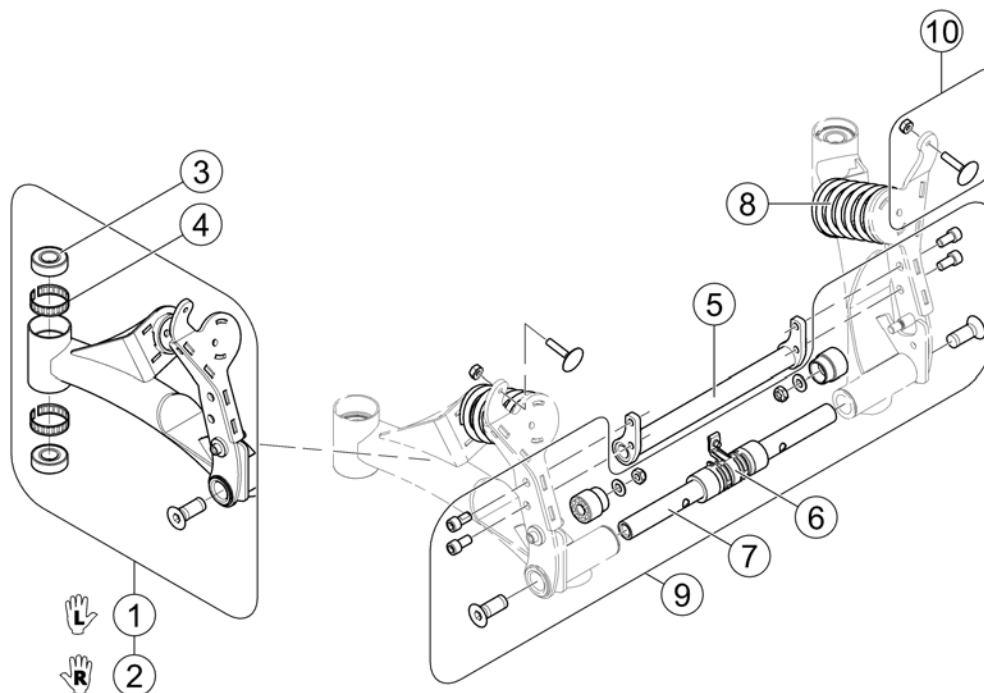
1563286.iss



1	SP1179094	(Kit, Rear Swing Arm Assy w/ Hardware, black, LH)	1
2	SP1179093	(Kit, Rear Swing Arm Assy w/ Hardware, black, RH)	1
3	SP40117X006	KULELAGER TDX SP2 BAK	1
4	SP1110717	(Ring, Tolerance)	1
5	SP1176604	(Crossmember, Pivot, black - Narrow Base)	1
6	SP1523166	FESTE DEKSEL TDX TDX	1
7	SP1155018	()	1
8	SP1579643	(Spring - silver)	1
9	SP1523177	FESTESETT ROTASJONSAKSEL TDX	1
10	SP1151949	(Kit, Elevator Bolt w/ Keps Nut (1/4-20))	1

(Shrouds Narrow Base)

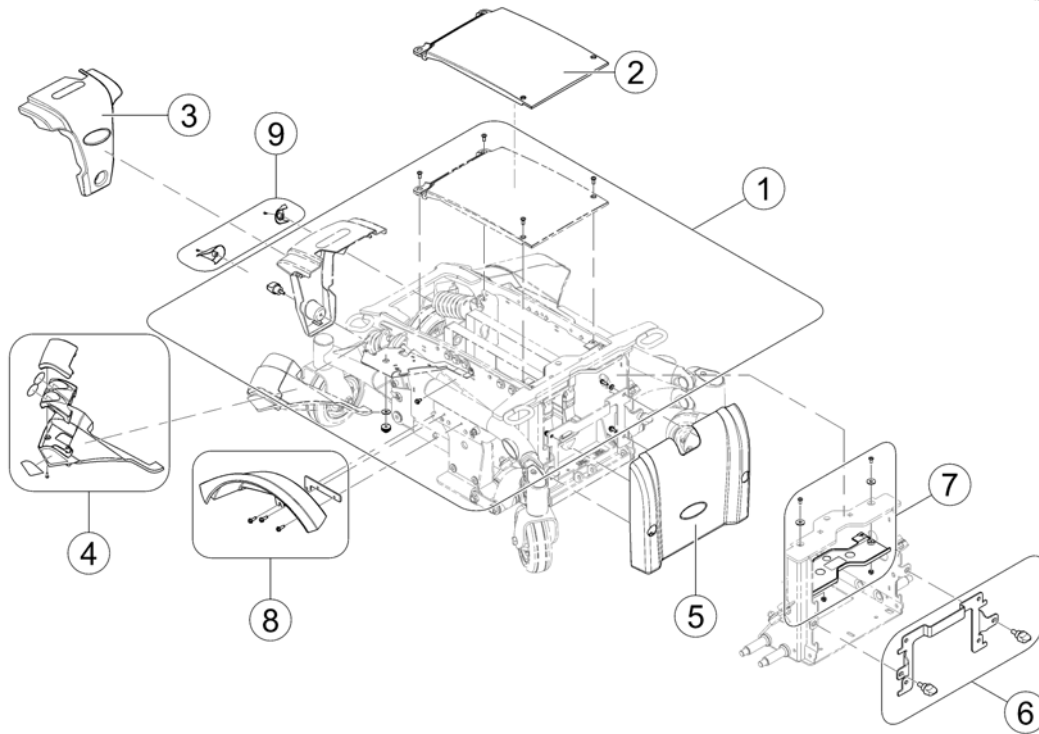
1583238.00



1	SP1179326	(Kit, Shroud Hardware)	1
2	SP1174634	(Shroud, Center, black - Narrow Base)	1
3	-	(Rear Shrouds Narrow Base)	1
4	-	(Bezel Shroud)	1
5	-	(Front Shroud Narrow Base)	1
6	SP1179897	(Front Shroud, Bracket & Hardware - Narrow Base)	1
7	SP1179738	(Isolation Shroud Assy - Narrow Base)	1
8	-	(Fender)	1
9	-	(Rear Caps)	1

(Rear Shrouds Narrow Base)

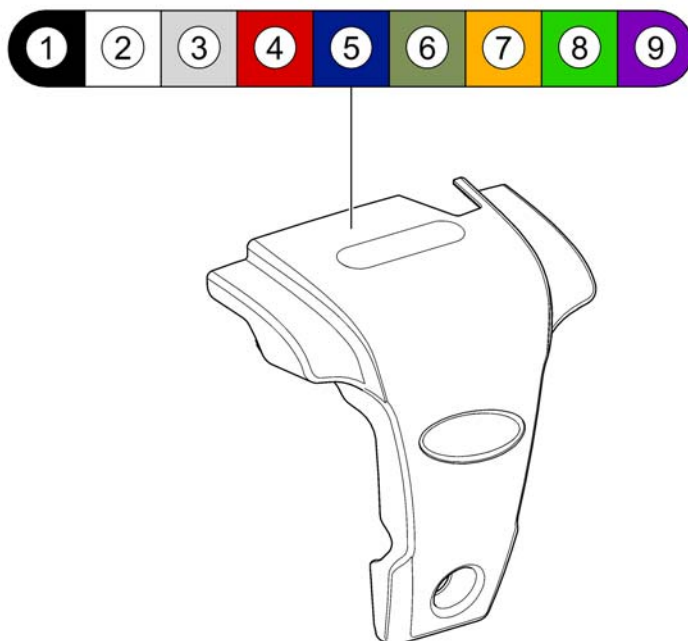
1583254.no



1	SP1176547	(Rear Shroud w/ Decals, Wet Black)		1
2	SP1176548	(Rear Shroud w/ Decals, Artic White)		1
3	SP1176549	(Rear Shroud w/ Decals, Electric Silver)		1
4	SP1176550	(Rear Shroud w/ Decals, Electric Red)		1
5	SP1176551	(Rear Shroud w/ Decals, Electric Blue)		1
6	SP1176552	(Rear Shroud w/ Decals, Aviator Green)	→ 30.06.2015	1
7	SP1176553	(Rear Shroud w/ Decals, Desert Orange)	→ 30.06.2015	1
8	SP1588735	(Rear Shroud w/ Decals, Jade Green)	01.07.2015 →	1
9	SP1588736	(Rear Shroud w/ Decals, Transparent Purple)	01.07.2015 →	1

(Front Shroud Narrow Base)

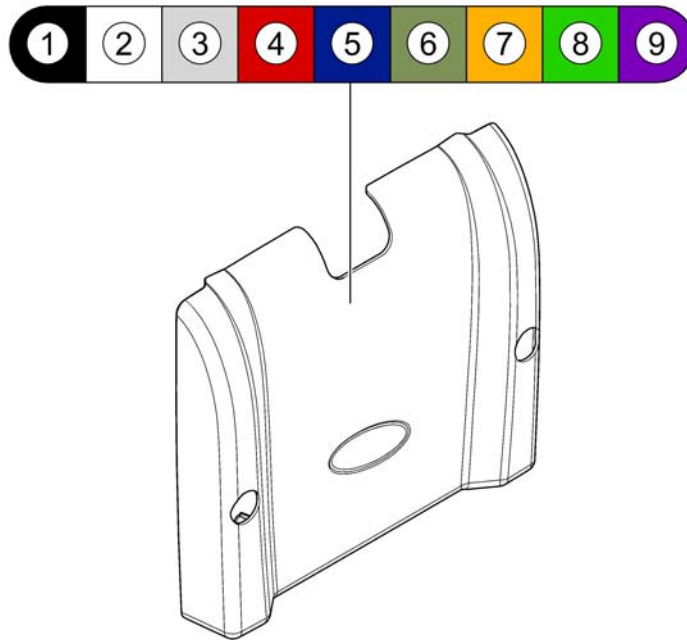1583258.iss






1	SP1179714	(Front Shroud Assy, Wet Black)		1
2	SP1189214	(Front Shroud Assy, Artic White)		1
3	SP1189215	(Front Shroud Assy, Electric Silver)		1
4	SP1179715	(Front Shroud Assy, Electric Red)		1
5	SP1189216	(Front Shroud Assy, Electric Blue)		1
6	SP1179716	(Front Shroud Assy, Aviator Green)	→ 30.06.2015	1
7	SP1189217	(Front Shroud Assy, Desert Orange)	→ 30.06.2015	1
8	SP1588783	(Front Shroud Assy, Jade Green)	01.07.2015 →	1
9	SP1588785	(Front Shroud Assy, Transparent Purple)	01.07.2015 →	1

Sete

SPCTDX_SP220.kso




1	-	(Modulite Seat)	1
2	-	(Ultra Low Maxx Seat)	1
3	-	(Max Seat)	1

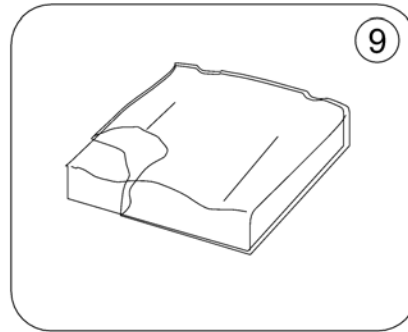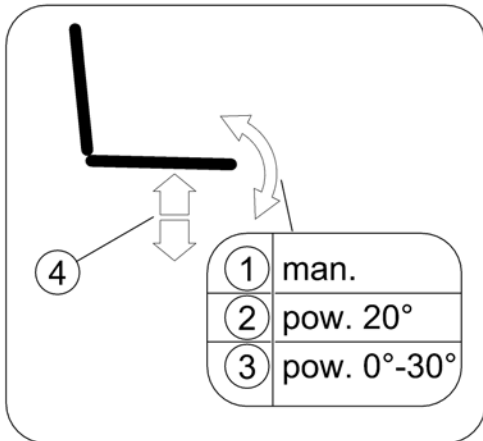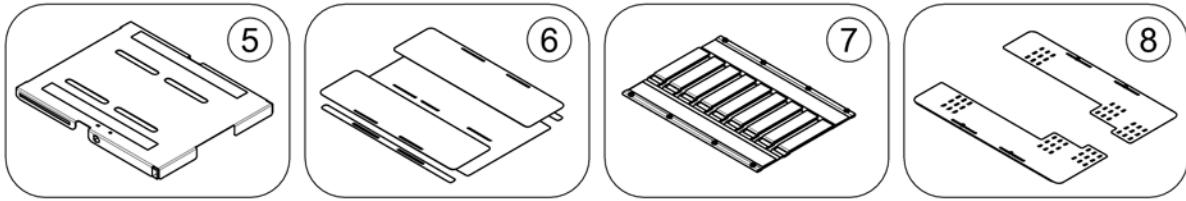

1	-	(man. settable Tilt) [1091]	1
2	-	(powered simple Tilt 20°) [1092]	1
3	-	(powered Tilt 0°-30°) [1093]	1
4	-	(Powered seat lifter with tilt 300mm/30°) [1094]	1
-	-	(Anterior Tilt) [2054]	1
5	-	(One Piece Seat Plate) [2110]	1
6	-	(Telescopic Seat Frame and Plate) [2111]	1
7	-	(Telescopic Seat Frame and Sling Seat) [2112]	1
8	-	(Telescopic Seat Frame and Universal Adapter) [2109]	1
9	-	(Seat Cushions "Standard" & "Matrx Flo-Tech")	1
9	-	(Seat Cushions "Matrx Libra") [2147]	1
9	-	(Seat Cushion "Matrx Vi") [2144]	1

(man. settable Tilt)

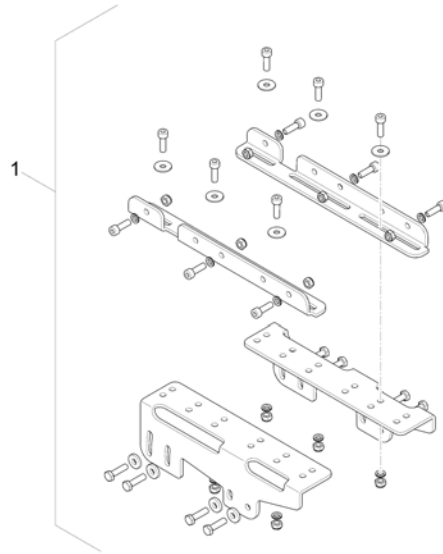
1577301 iso



1 SP1563229

(Mounting Rails)

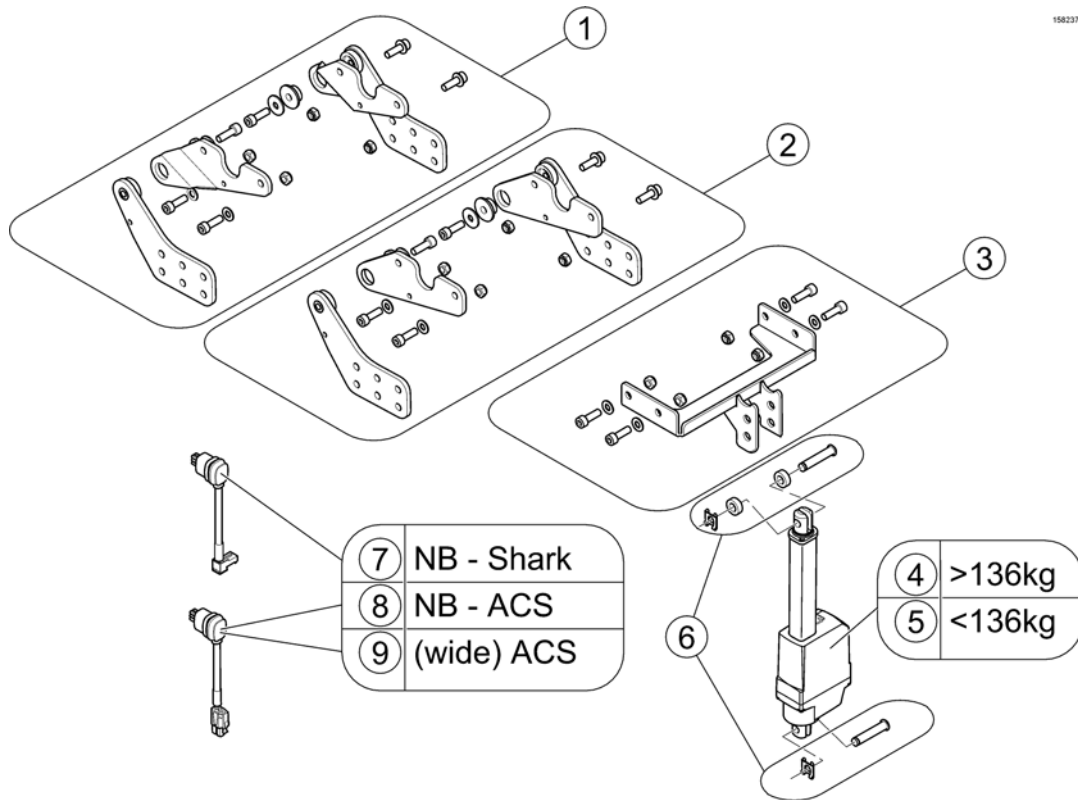
1



1	SP1582738	(Pivot Bearing Back - Wide Base)	1
2	SP1582736	(Pivot Bearing Back - Narrow Base)	1
3	SP1582731	(Support Plate Front for Tilt)	1
4	SP1556487	(Actuator LA23 - stroke 135mm / 2500N)	1
5	SP1581020	(Actuator LA23 - stroke 135mm / 1800N)	1
6	SP1582737	(Mounting Hardware for Tilt)	1
7	SP1552297	(Actuator Cable LA23 Shark (~500mm))	1
8	SP1556449	(Actuator Cable LA23 ACS2 (~400mm))	1
9	SP1572166	KABEL ACTUATOR MODULITE	1

(powered Tilt 0°-30°)

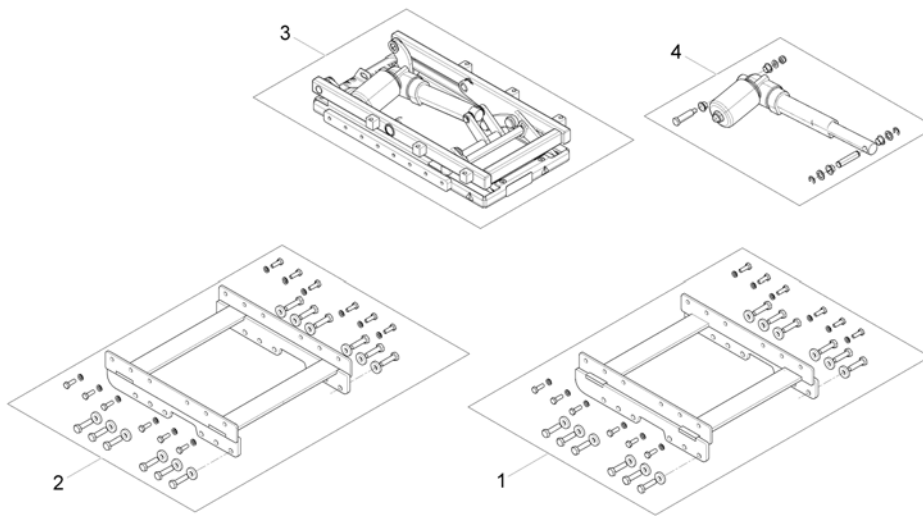
1577302.no



1	SP1568426	(Lifter- /Tilt Bracket for wide Base)	1
2	SP1568424	(Lifter- /Tilt Bracket for narrow Base)	1
3	SP1558714	(Tilt unit)	1
4	SP1571809	(Actuator)	1
-	SP1570291	(Adapter Cable)	1

(Powered seat lifter with tilt 300mm/30°)

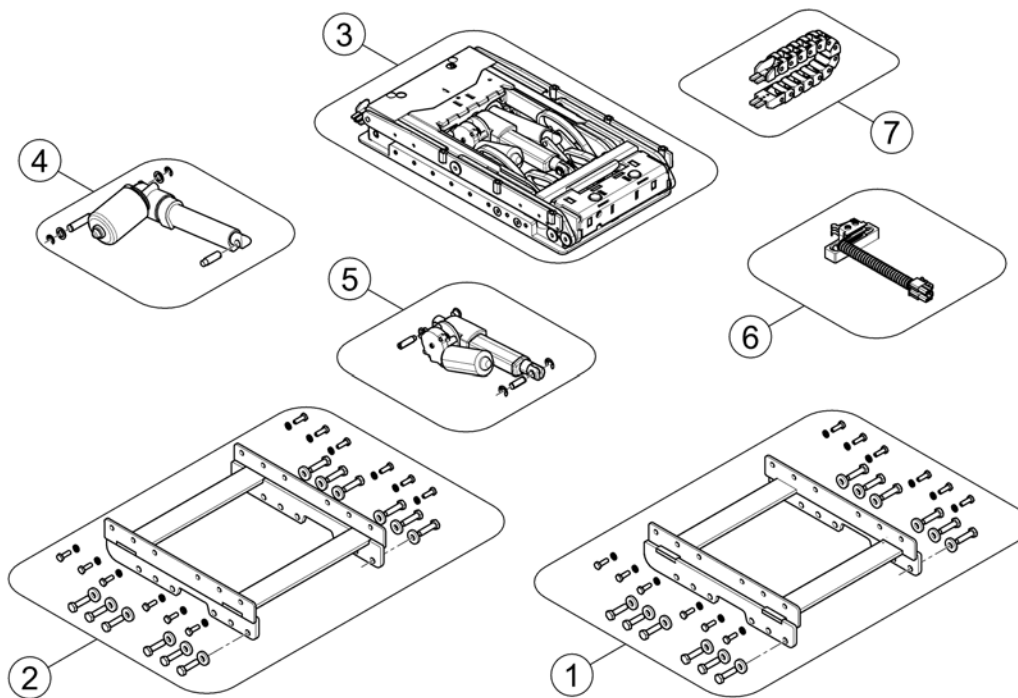
1577303 iso



-	SP1570291	(Adapter Cable) <i>(2x for wide Base / 1x for Narrow Base)</i>		1
-	SP1576062	(Extension Cable - 1100mm) <i>(only for Narrow Base)</i>		1
-	SP1570294	(Extension cable A (700mm))		1
-	SP1533561	(Light Cable w. DCI)		1
1	SP1568426	(Lifter- /Tilt Bracket for wide Base)		1
2	SP1568424	(Lifter- /Tilt Bracket for narrow Base)		1
3	SP1558715	SETELØFT MED TILT KPL TDX SP2 LOW-RIDER		1
4	SP1571823	(Lifter actuator, cpl. with mounting hardware)	→ 30.06.2016	1
4	SP1588442	(Lifter actuator, cpl. with hardware)	01.07.2016 →	1
5	SP1571822	AKTUATOR FOR TILT M. LØFT TDX SP2		1
6	SP1577762	(Slow-Stop Switchcable with Mounting (LS18KIT))		1
7	SP1604199	(Energie chain)		1

(Anterior Tilt)

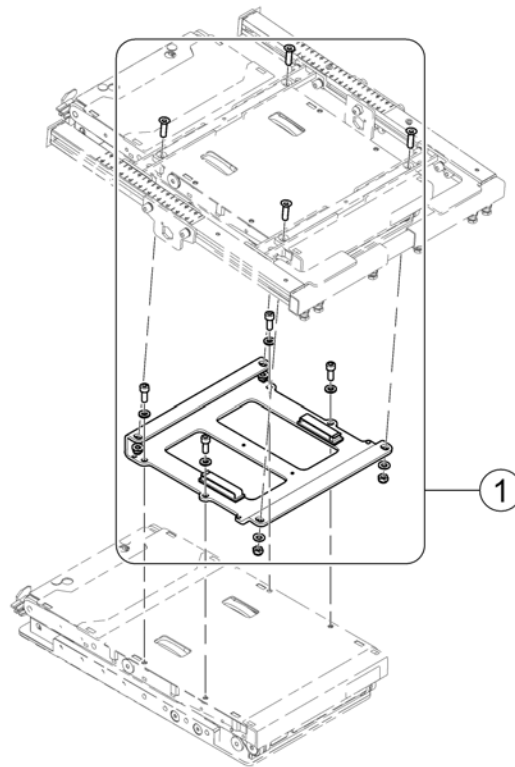1602720.kim




1 SP1590707

(Anterior Tilt, cpl.)

1

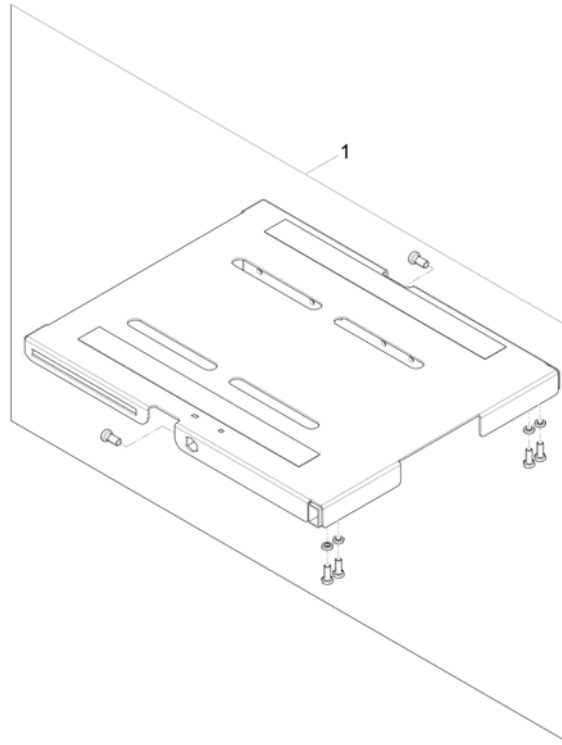1571064.000




1 SP1571297

(One Piece Seat Plate)

1


(Telescopic Seat Frame and Plate)

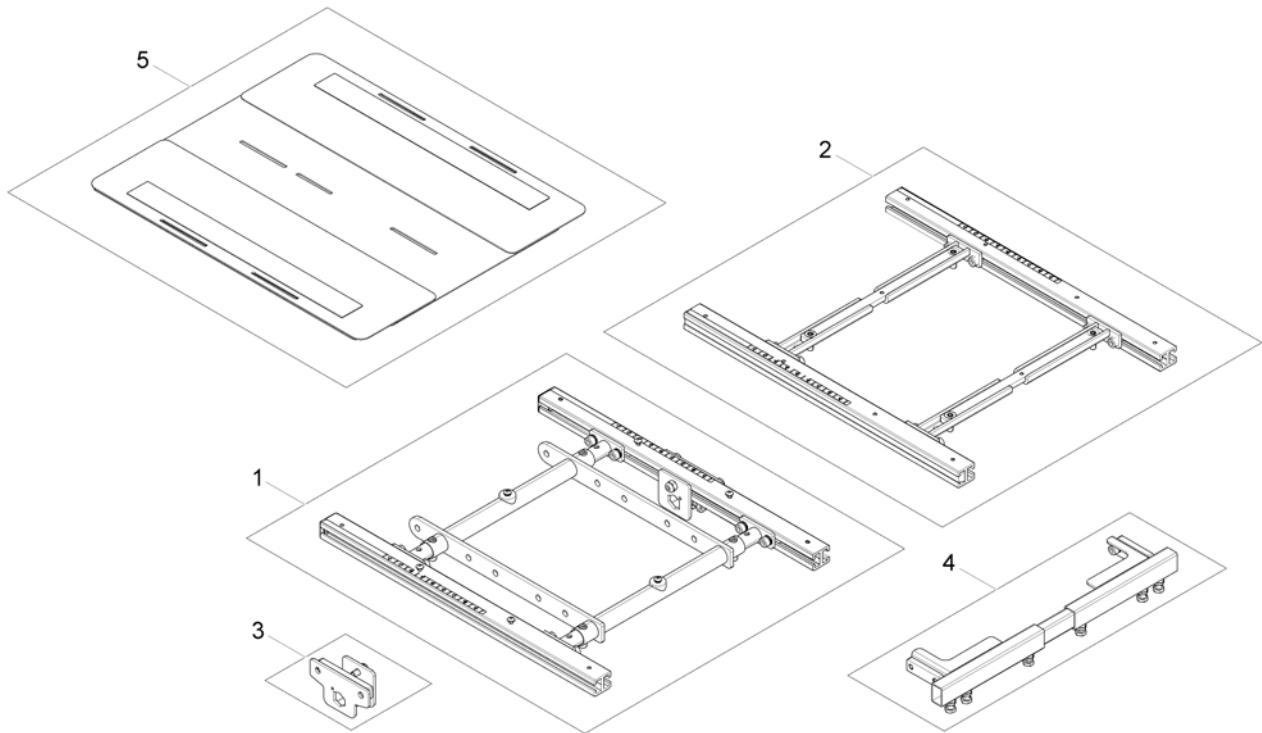1572312 iso



1	-	(Telescopic Seat Frame for man. /pow. 20° Tilt)	1
2	-	(Telescopic Seat Frame for Tilt 30° / Lifter)	1
3	-	(Armrest Bracket)	1
4	-	(Legrest Support)	1
5	-	(Seat Plate)	1

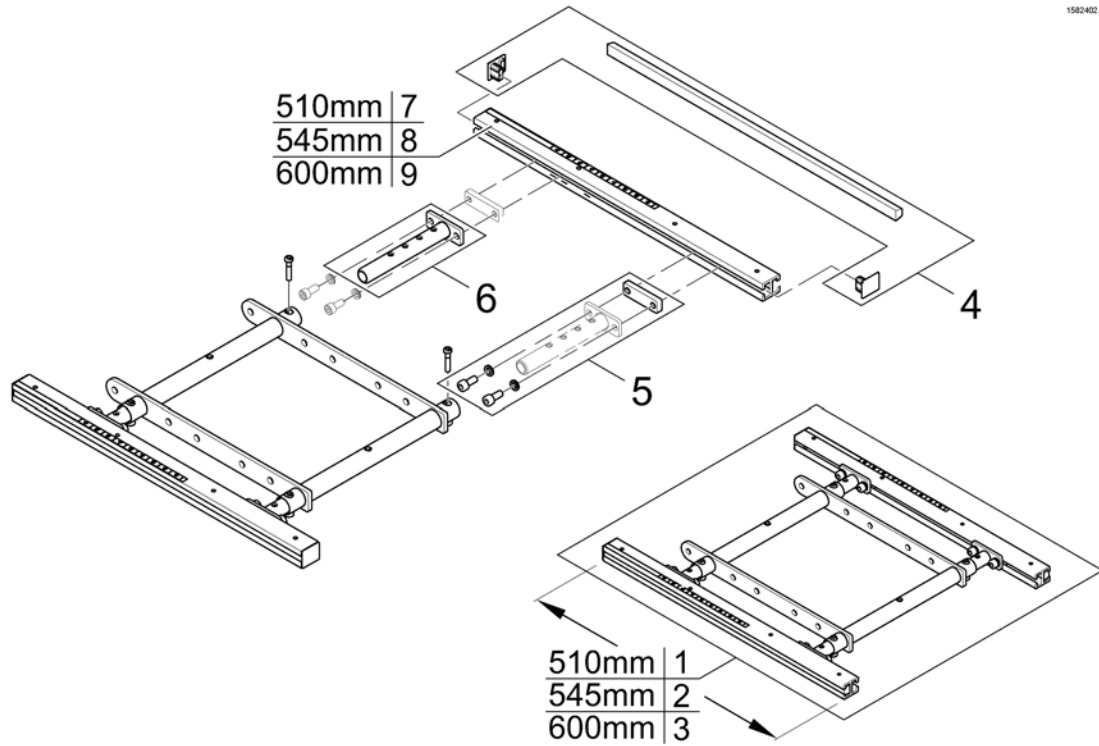
(Telescopic Seat Frame for man. /pow. 20° Tilt)



1	SP1572977	(Seatframe 510mm, cpl., pre set up to SW 43cm)	1
2	SP1582216	SETERAMME KPL MODULITE MODULITE	1
3	SP1582739	(Seatframe 600mm, cpl., pre set up to SW 43cm)	1
4	SP1571810	(Sealing strap with plugs)	1
5	SP1571811	(Groove stone with mounting parts)	1
6	SP1554646	RØRBRAKETT SETEBUNN MODULITE	1
7	SP1588443	C-PROFIL SETE 51 CM LANG MODULITE	1
8	SP1588444	C-PROFIL SETE 44 CM LANG MODULITE	1
9	SP1588445	C-PROFIL SETE 60CM LANG MODULITE	1


(Telescopic Seat Frame for Tilt 30° / Lifter)

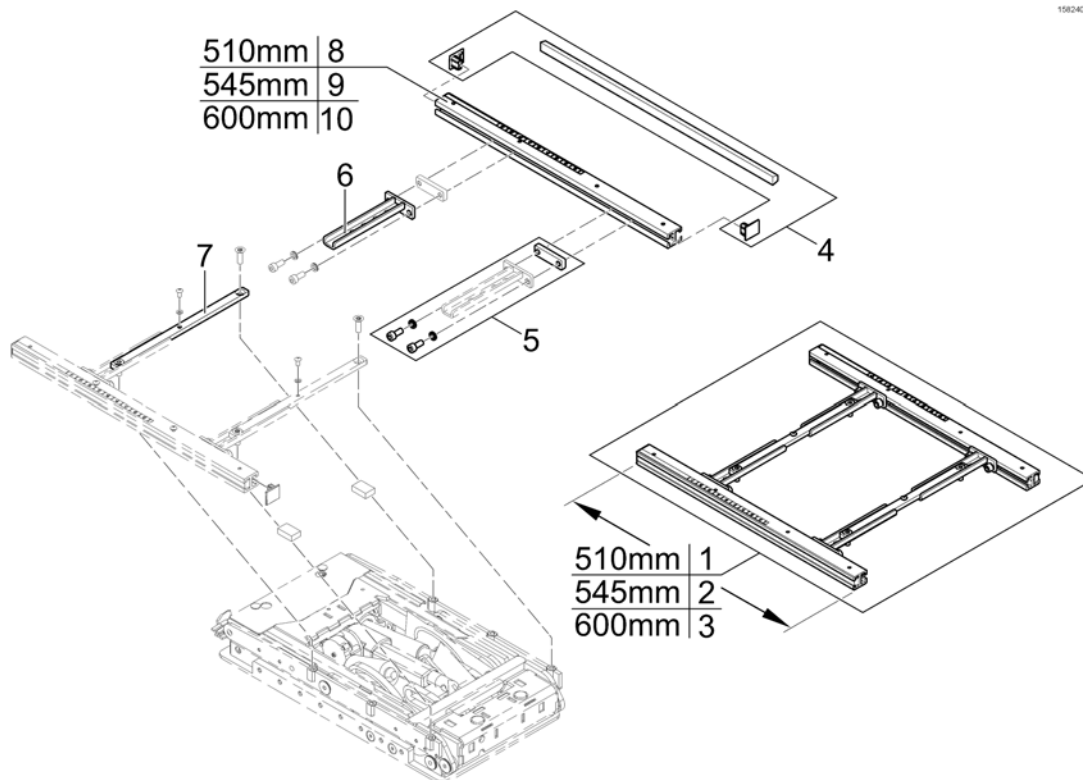1582403 Rev



1	SP1572992	(Seatframe 510mm, pre set up to SW 430mm, cpl.)	1
2	SP1582219	(Seatframe 545mm, pre set up to SW 430mm, cpl.)	1
3	SP1582740	(Seatframe 600mm, pre set up to SW 430mm, cpl.)	1
4	SP1571810	(Sealing strap with plugs)	1
5	SP1571811	(Groove stone with mounting parts)	1
6	SP1567007	(Seat Frame Adapter)	1
7	SP1567008	(Middle Rail)	2
8	SP1588443	C-PROFIL SETE 51 CM LANG MODULITE	1
9	SP1588444	C-PROFIL SETE 44 CM LANG MODULITE	1
10	SP1588445	C-PROFIL SETE 60CM LANG MODULITE	1

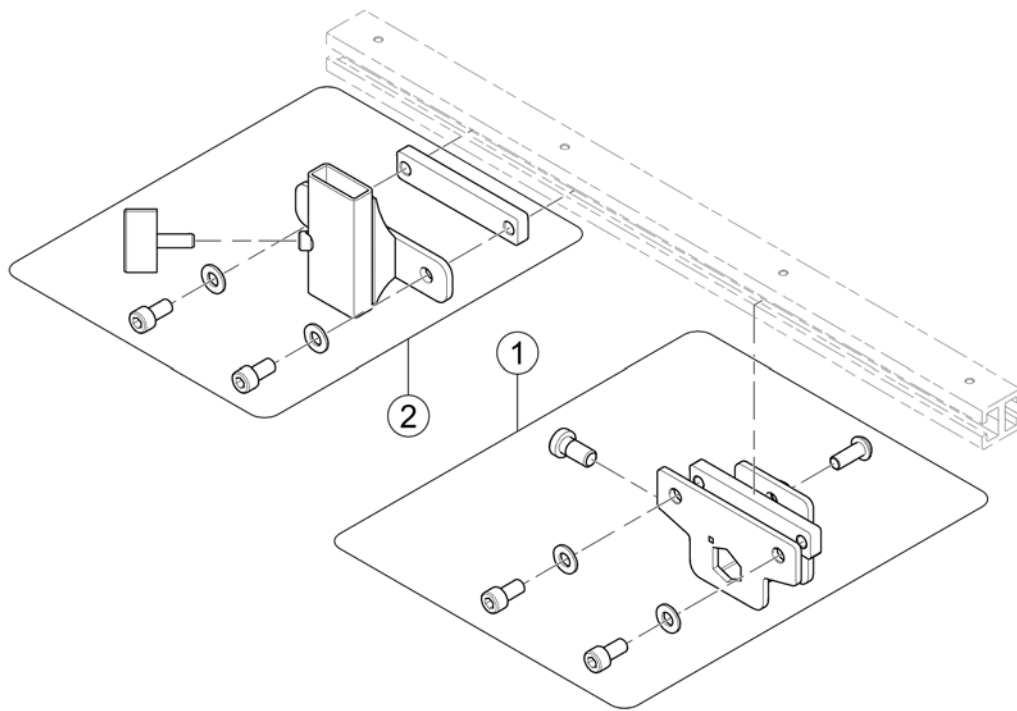
1587804 Rev




1	SP1570189	ARMLENSBRAKETT MODULITE	1
2	SP1585469	HOLDER ARMLEN LOW-RIDER TDX SP2	1

(Legrest Support)

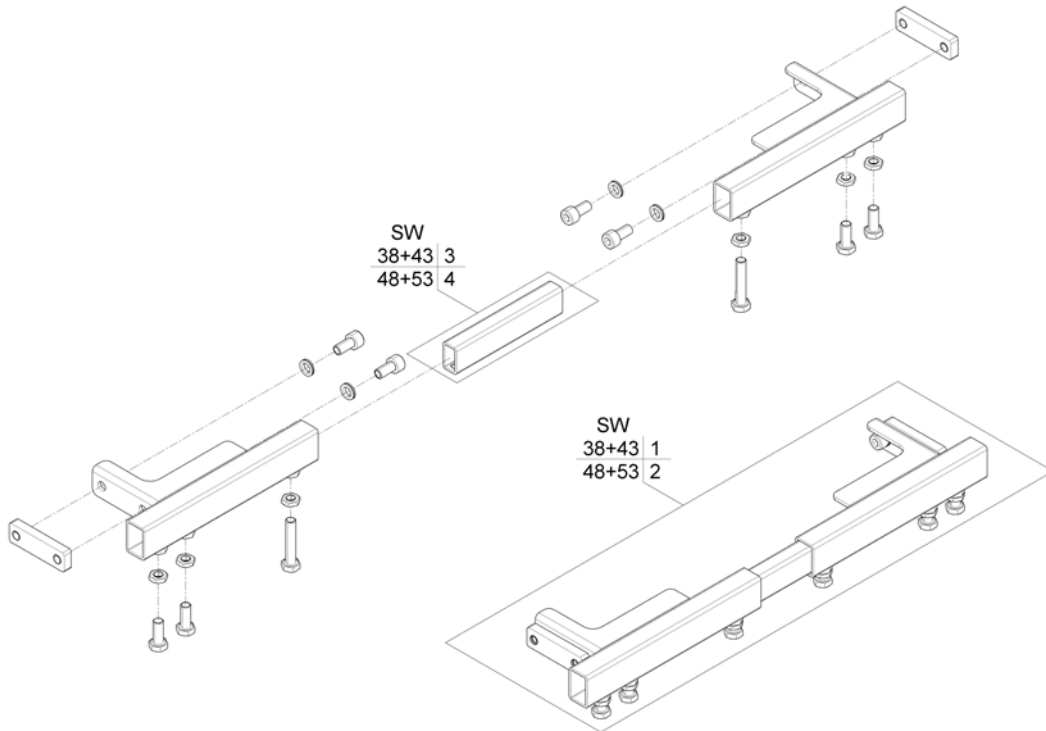1572021.000




1	SP1558935	BENSTØTTE FESTE SB38+43 MODULITE	1
2	SP1559294	BENSTØTTE FESTE SB48+53 MODULITE	1
3	SP1558134	MELLOMSTYKKE FESTE BENST MODULITE	1
4	SP1558135	BRAKETT SETEPLATE BENST MODULITE SB48-53	1

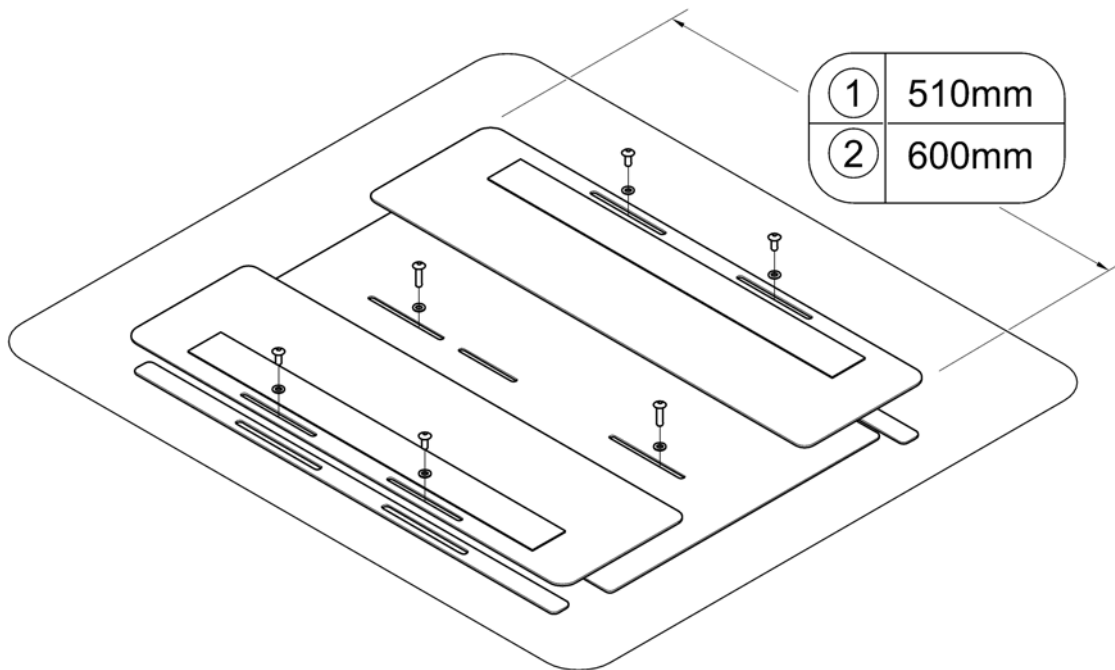

1	SP1572199	SETEPLATE MODULITE	1
2	SP1580756	SETEPLATE 60CM KPL LANG MODULITE	1

(Telescopic Seat Frame and Sling Seat)

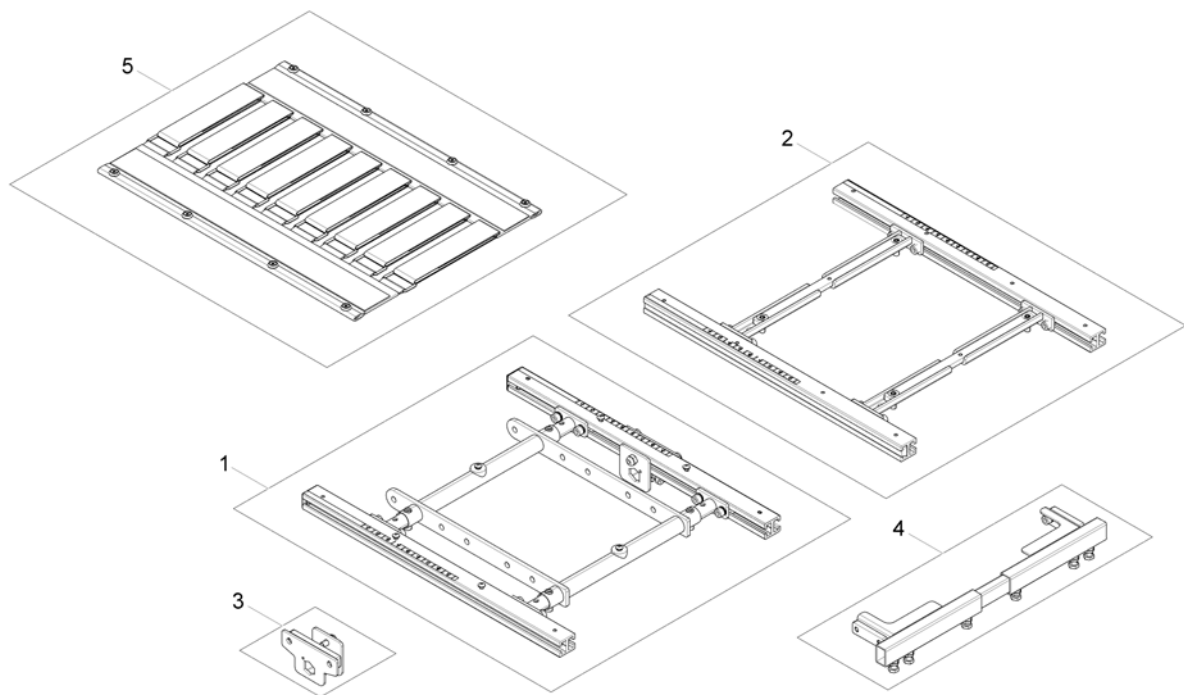
1572311.00



1	-	(Telescopic Seat Frame for man. / pow. 20° Tilt)	1
2	-	(Telescopic Seat Frame for Tilt 30° / Lifter)	1
3	-	(Armrest Bracket)	1
4	-	(Legrest Support)	1
5	-	(Sling Seat)	1


(Telescopic Seat Frame for man. / pow. 20° Tilt)

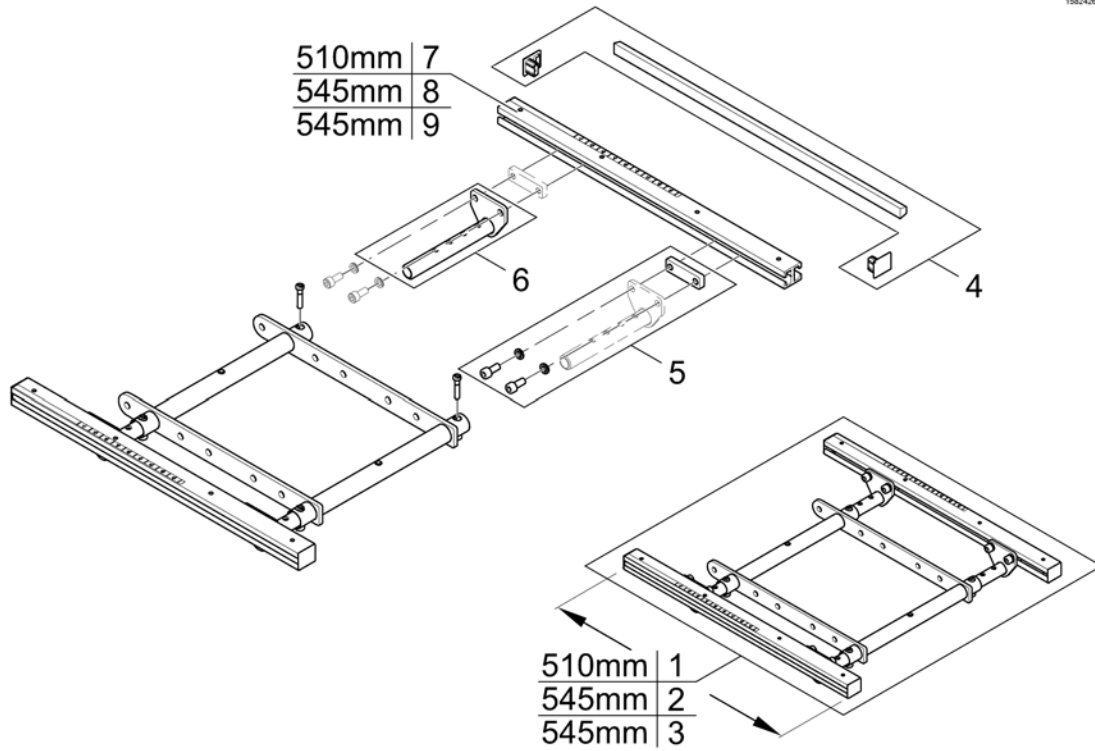1582426.iss



1	SP1572978	(Seatframe 510mm, cpl., pre set up to SW 43cm)	1
2	SP1582217	(Seatframe 545mm, cpl., pre set up to SW 43cm)	1
3	SP1582741	(Seatframe 600mm, cpl., pre set up to SW 43cm)	1
4	SP1571810	(Sealing strap with plugs)	1
5	SP1571811	(Groove stone with mounting parts)	1
6	SP1555400	(Adjustment Tube)	1
7	SP1588443	C-PROFIL SETE 51 CM LANG MODULITE	1
8	SP1588444	C-PROFIL SETE 44 CM LANG MODULITE	1
9	SP1588445	C-PROFIL SETE 60CM LANG MODULITE	1

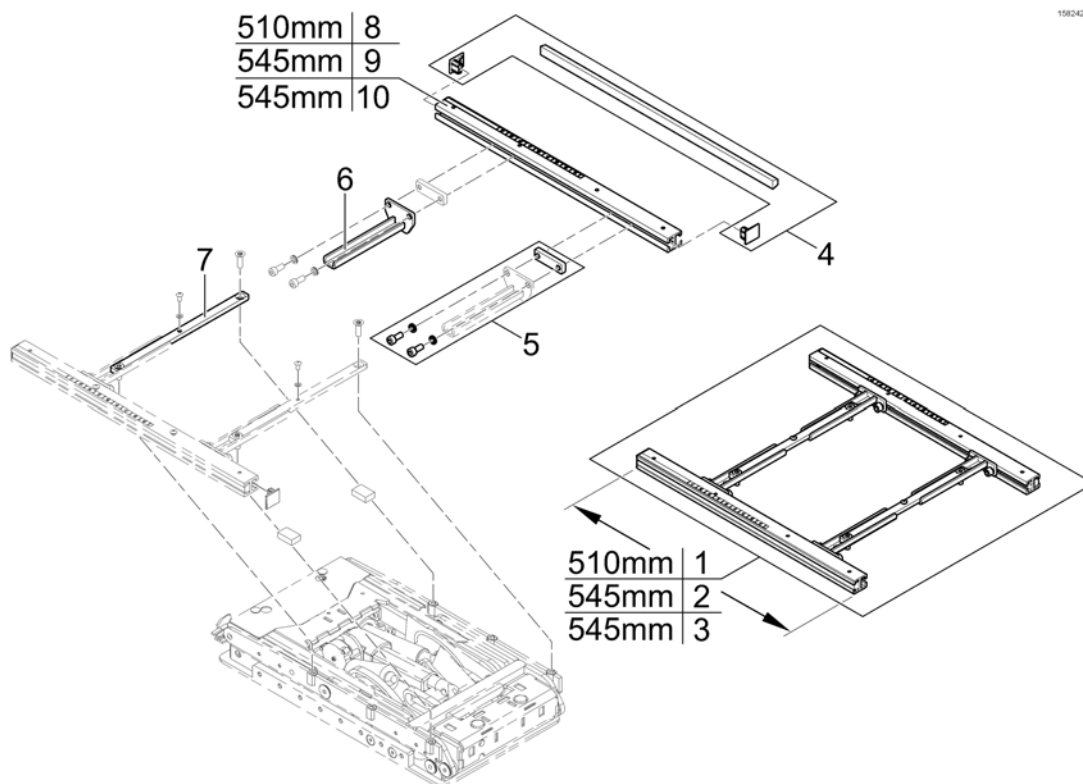
(Telescopic Seat Frame for Tilt 30° / Lifter)



1	SP1572993	(Seatframe 510mm, pre set up to SW 430mm, cpl.)	1
2	SP1582230	(Seatframe 545mm, pre set up to SW 430mm, cpl.)	1
3	SP1582742	(Seatframe 600mm, pre set up to SW 430mm, cpl.)	1
4	SP1571810	(Sealing strap with plugs)	1
5	SP1571811	(Groove stone with mounting parts)	1
6	SP1567720	VELCRO SETEADAPTER MODULITE	1
7	SP1567008	(Middle Rail)	1
8	SP1588443	C-PROFIL SETE 51 CM LANG MODULITE	1
9	SP1588444	C-PROFIL SETE 44 CM LANG MODULITE	1
10	SP1588445	C-PROFIL SETE 60CM LANG MODULITE	1

(Legrest Support)

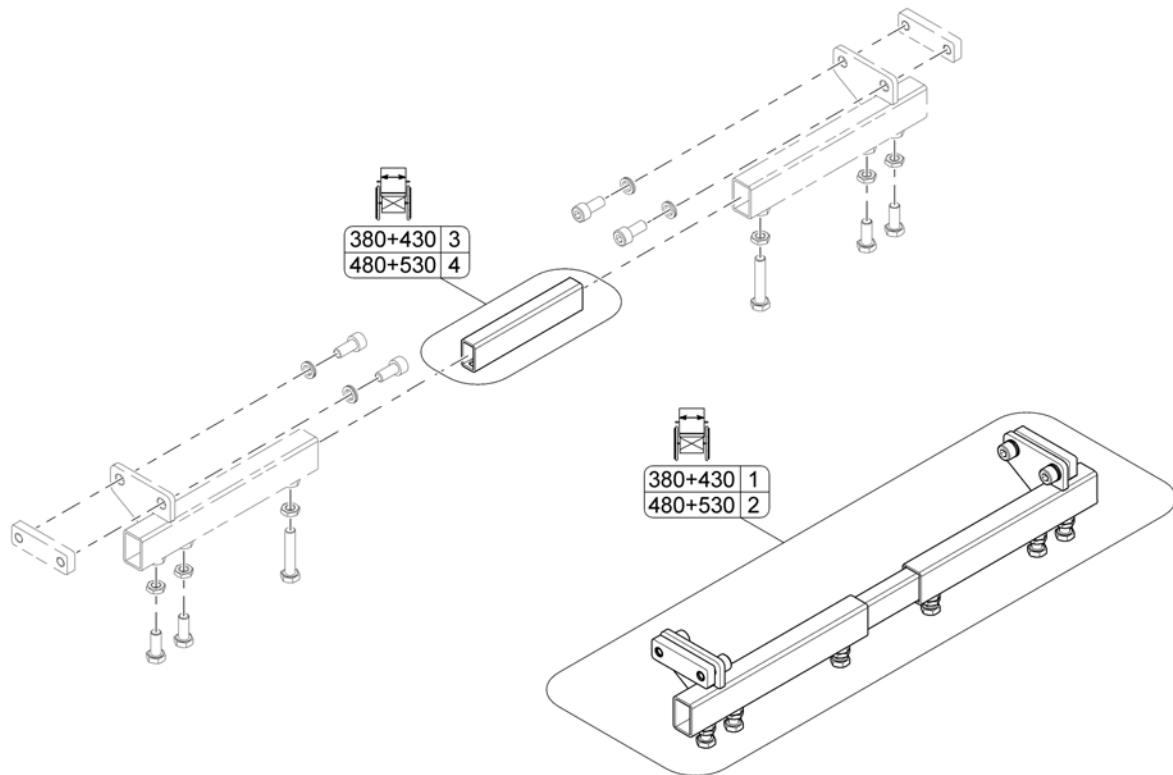1572030.060




1	SP1559305	FIRKANTRØR FESTE BENST MODULITE	1
2	SP1559307	(Legrest Support SW480+530mm, cpl.)	1
3	SP1558134	MELLOMSTYKKE FESTE BENST MODULITE	1
4	SP1558135	BRAKETT SETEPLATE BENST MODULITE SB48-53	1

(Sling Seat)

1582429 iso


	SW	X
①	380/430	510mm
②	480/530	510mm
③	380/430	600mm
④	480/530	600mm

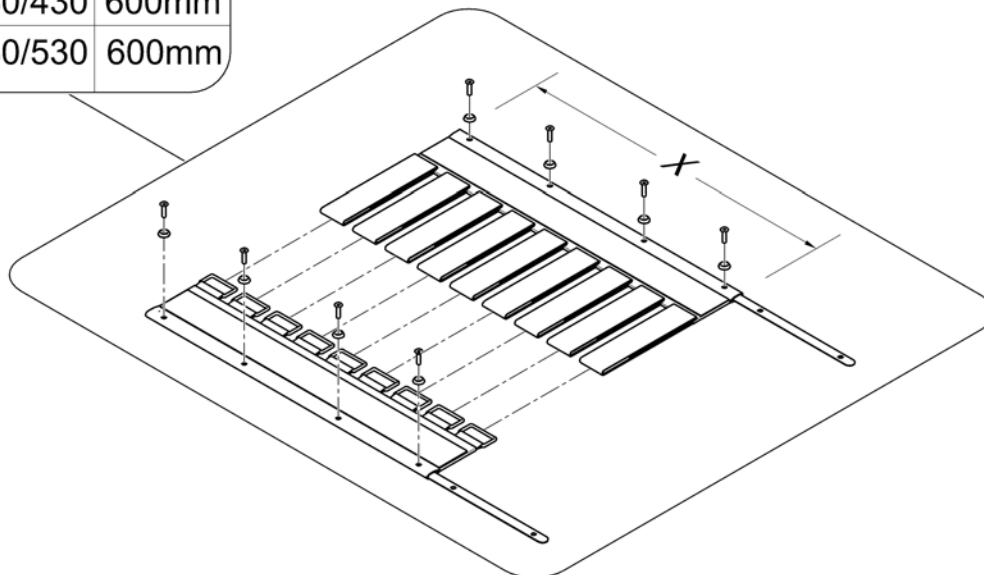


1	SP1557877	VELCRO SETE SB38-44 STORM4 TDX SP2	1
2	SP1555911	VELCRO SETEBUNN SB46 53 STORM4 TDX SP2	1
3	SP1580757	(Sling Seat SW380/430, 600mm)	1
4	SP1580758	(Sling Seat SW480/530, 600mm)	1

(Telescopic Seat Frame and Universal Adapter)

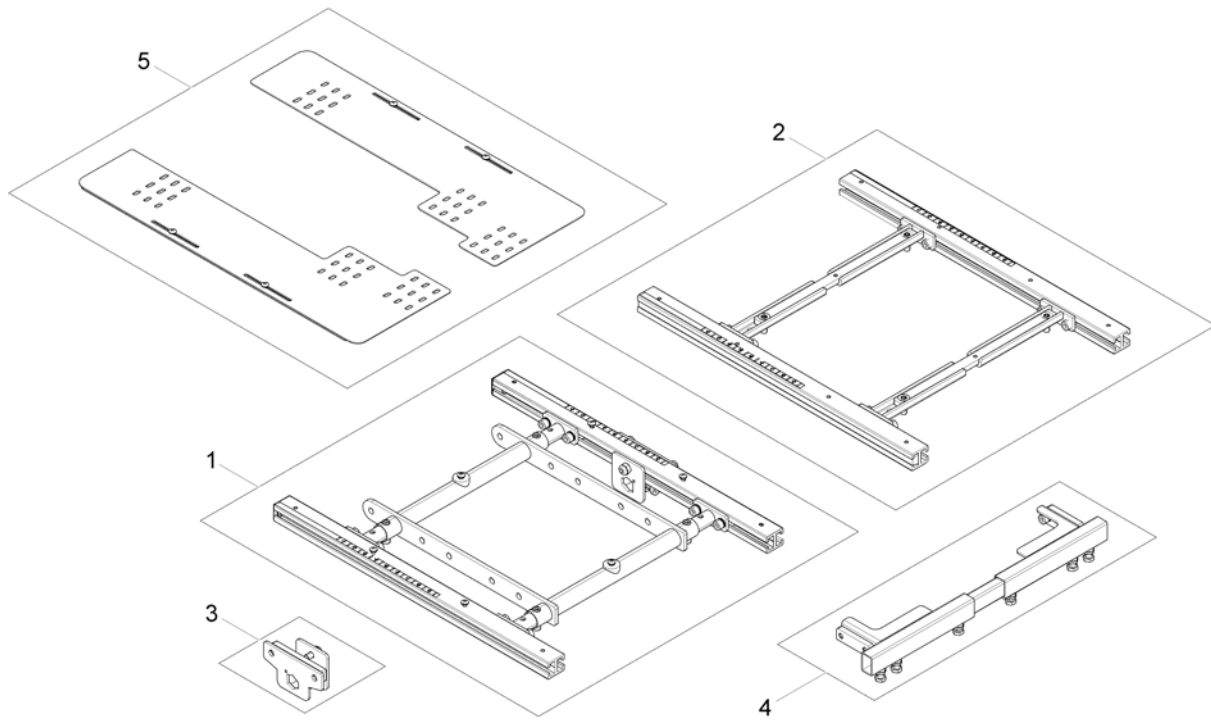1581864 iso



1	-	(Telescopic Seat Frame for man. /pow. 20° Tilt)	1
2	-	(Telescopic Seat Frame for Tilt 30° / Lifter)	1
3	-	(Armrest Bracket)	1
4	-	(Legrest Support)	1
5	-	(Adapter Kit for Universal Seat)	1

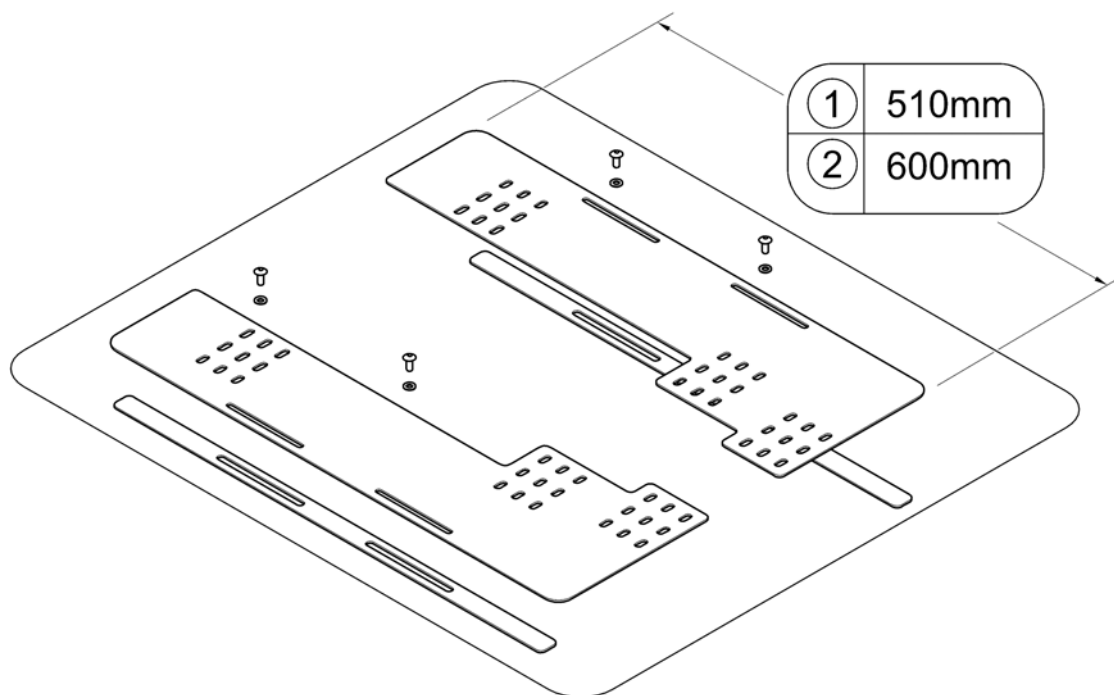



1	SP1580750	(Adapter Kit for Universal Seat 510mm)	1
2	SP1580759	(Adapter Kit for Universal Seat 600mm)	1

(Seat Cushions "Standard" & "Matrix Flo-Tech")



1582643 iso

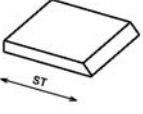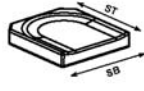Standard


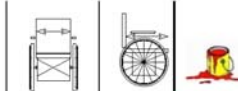
1	41cm	black
2	46cm	
3	51cm	

**Matrix Flo-tech
Lite / Contour**

4	38cm	41cm	black smooth	
5	38cm	46cm		
6	38cm	51cm		
7	43cm	41cm		
8	43cm	46cm		
9	43cm	51cm		
10	48cm	41cm		
11	48cm	46cm		
12	48cm	51cm		
13	53cm	43cm		
14	53cm	46cm		
15	53cm	51cm		
16	53cm	56cm		
17	38cm	41cm		Dartex black
18	38cm	46cm		
19	38cm	51cm		
20	43cm	41cm		
21	43cm	46cm		
22	43cm	51cm		
23	48cm	41cm		
24	48cm	46cm		
25	48cm	51cm		
26	53cm	43cm		
27	53cm	46cm		
28	53cm	51cm		

**Matrix Flo-tech
Visco**

55	43cm	43cm	black
----	------	------	-------

1	1420460	SETEPUTE SD41 STD DRAGON VERTIC <i>[2141]</i>	1
2	1420461	SETEPUTE SD46 STD DRAGON VERTIC <i>[2141]</i>	1
3	1420462	SETEPUTE SD51 STD DRAGON VERTIC <i>[2141]</i>	1
4	SPFLB15/16	SETEPUTE SB38 SD41 SORT STORM4 <i>[2142+2151]</i>	1
5	SPFLB15/18	SETEPUTE SB38 SD46 SORT STORM4 <i>[2142+2151]</i>	1
6	SPFLB15/20	SETEPUTE SB38 SD51 SORT STORM4 <i>[2142+2151]</i>	1
7	SPFLB17/16	SETEPUTE SB43 SD41 SORT STORM4 <i>[2142+2151]</i>	1
8	SPFLB17/18	SETEPUTE SB43 SD46 STD STORM4 <i>[2142+2151]</i>	1
9	SPFLB17/20	SETEPUTE SB43 SD51 SORT STORM4 <i>[2142+2151]</i>	1
10	SPFLB19/16	SETEPUTE SB48 SD41 SORT STORM4 <i>[2142+2151]</i>	1
11	SPFLB19/18	SETEPUTE SB48 SD46 SORT STORM4 <i>[2142+2151]</i>	1
12	SPFLB19/20	SETEPUTE SB48 SD51 SORT STORM4 <i>[2142+2151]</i>	1
13	SPFCB21/17	SETEPUTE SB53 SD43 SORT STORM4 <i>[2142+2151]</i>	1
14	SPFCB21/18	SETEPUTE SB53 SD46 SORT STORM4 <i>[2142+2151]</i>	1
15	SPFCB21/20	SETEPUTE SB53 SD51 SORT STORM4 <i>[2142+2151]</i>	1
16	SP1578042	(Seatcushion Flo-tech Contour SW53 x SD56cm black) <i>[2142+2151]</i>	1
17	SP1529546-TR26	(seat cushion Flo-tech Lite SW 38 x SD 41 Dartex black) <i>[2142+2152]</i>	1
18	SP1529547-TR26	(seat cushion Flo-tech Lite SW 38 x SD 46 Dartex black) <i>[2142+2152]</i>	1

(Seat Cushions "Standard" & "Matrx Flo-Tech")


19	SP1529548-TR26	(seat cushion Flo-tech Lite SW 38 x SD 51 Dartex black) [2142+2152]	1
20	SP1529549-TR26	(seat cushion Flo-tech Lite SW 43 x SD 41 Dartex black) [2142+2152]	1
21	SP1529550-TR26	SETEPUTE SB43SD46 DARTEX STORM4 [2142+2152]	1
22	SP1529551-TR26	SETEPUTE SB43SD51 DARTEX STORM4 [2142+2152]	1
23	SP1529552-TR26	(seat cushion Flo-tech Lite SW 48 x SD 41 Dartex black) [2142+2152]	1
24	SP1529553-TR26	(seat cushion Flo-tech Lite SW 48 x SD 46 Dartex black) [2142+2152]	1
25	SP1529554-TR26	SETEPUTE SB48SD51 DARTEX STORM4 [2142+2152]	1
26	SP1529555-TR26	(seat cushion Flo-tech Contour SW 53 x SD 43 Dartex black) [2142+2152]	1
27	SP1529556-TR26	(seat cushion Flo-tech Contour SW 53 x SD 46 Dartex black) [2142+2152]	1
28	SP1529570-TR26	(seat cushion Flo-tech Contour SW 53 x SD 51 Dartex black) [2142+2152]	1
29	SP1529546-TR35	(seat cushion Flo-tech Lite SW 38 x SD 41 black chenille fabric) [2142+2154]	1
30	SP1529547-TR35	(seat cushion Flo-tech Lite SW 38 x SD 46 black chenille fabric) [2142+2154]	1
31	SP1529548-TR35	(seat cushion Flo-tech Lite SW 38 x SD 51 black chenille fabric) [2142+2154]	1
32	SP1529549-TR35	(seat cushion Flo-tech Lite SW 43 x SD 41 black chenille fabric) [2142+2154]	1
33	SP1529550-TR35	(seat cushion Flo-tech Lite SW 43 x SD 46 black chenille fabric) [2142+2154]	1
34	SP1529551-TR35	(seat cushion Flo-tech Lite SW 43 x SD 51 black chenille fabric) [2142+2154]	1
35	SP1529552-TR35	(Seat cushion Flo-tech Lite SW 48 x SD 41 black chenille fabric) [2142+2154]	1
36	SP1529553-TR35	SETEPUTE SB48 SD46 SORT STORM4 [2142+2154]	1
37	SP1529554-TR35	SETEPUTE SB48SD51 PLYSJ STORM4 [2142+2154]	1
38	SP1529555-TR35	(seat cushion Flo-tech Contour SW 53 X SD 43 black chenille fabric) [2142+2154]	1
39	SP1529556-TR35	(seat cushion Flo-tech Contour SW 53 x SD 46 black chenille fabric) [2142+2154]	1
40	SP1529570-TR35	(seat cushion Flo-tech Contour SW 53 x SD 51 black chenille fabric) [2142+2154]	1
41	SP1529546	SETEPUTE EASY-A SB38SD41 SORT M ORANG SØM [2142+2153]	1
42	SP1529547	SETEPUTE EASY-A SB38SD46 SORT M ORANG SØM [2142+2153]	1
43	SP1529548	SETEPUTE EASY-A SB38SD51 SORT M ORANG SØM [2142+2153]	1
44	SP1529549	SETEPUTE EASY-A SB43SD41 SORT M ORANG SØM [2142+2153]	1
45	SP1529550	SETEPUTE EASY-A SB43SD46 SORT M ORANG SØM [2142+2153]	1
46	SP1529551	SETEPUTE EASY-A SB43SD51 SORT M ORANG SØM [2142+2153]	1
47	SP1529552	SETEPUTE EASY-A SB48SD41 SORT M ORANG SØM [2142+2153]	1
48	SP1529553	SETEPUTE EASY-A SB48SD46 SORT M ORANG SØM [2142+2153]	1
49	SP1529554	SETEPUTE EASY-A SB48SD51 SORT M ORANG SØM [2142+2153]	1
50	SP1529555	SETEPUTE EASY-A SB53SD43 SORT M ORANG SØM [2142+2153]	1

(Seat Cushions "Standard" & "Matrx Flo-Tech")

51	SP1529556	SETEPUTE EASY-A SB53SD46 SORT M ORANG SØM <i>[2142+2153]</i>	1
52	SP1529570	SETEPUTE EASY-A SB53 SD5 SORT M ORANG SØM <i>[2142+2153]</i>	1
55	SP1501113	SETEPUTE 43*43 TDX FLOTEC CONTUR VISCO <i>[2143]</i>	1

(Seat Cushions "Matrx Libra")

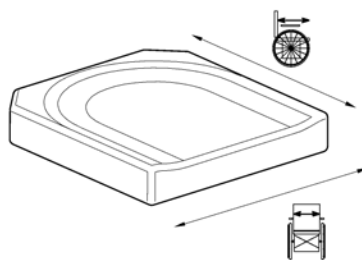1603356.iss




Matrx Libra

		mm			
		380	430	480	530
mm	410	1	4	7	-
	450	2	5	8	10
	510	3	6	9	11
	560	-	-	-	12



1	LC1516	(Matrx Libra cushions 380x410 mm)	1
2	LC1518	(Matrx Libra cushions 380x460 mm)	1
3	LC1520	(Matrx Libra cushions 380x510 mm)	1
4	LC1716	(Matrx Libra cushions 430x410 mm)	1
5	LC1718	(Matrx Libra cushions 430x460 mm)	1
6	LC1720	(Matrx Libra cushions 430x510 mm)	1
7	LC1916	(Matrx Libra cushions 480x410 mm)	1
8	LC1918	(Matrx Libra cushions 480x460 mm)	1
9	LC1920	(Matrx Libra cushions 480x510 mm)	1
10	LC2118	(Seat Cushion Matrx Libra 21x18")	1
11	LC2120	(Seat Cushion Matrx Libra 21x20")	1
12	LC2122	(Seat Cushion Matrx Libra 21x22")	1

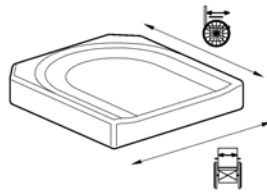(Seat Cushion "Matrx Vi")

1606756.kim



Matrx Vi

		 mm			
		380	430	480	530
 mm	410	1	4	8	-
	460	2	5	9	12
	510	3	6	10	13
	560	-	7	11	14


1	SP1579786	(Matrx Vi 16x16")	1
2	SP1579787	(Matrx Vi 16x18")	1
3	SP1579788	(Matrx Vi 16x20")	1
4	SP1579920	(Matrx Vi 18x16")	1
5	SP1579921	(Matrx Vi 18x18")	1
6	SP1579922	(Matrx Vi 18x20")	1
7	SP1584183	(Matrx Vi 18x22")	1
8	SP1579923	(Matrx Vi 20x16")	1
9	SP1579924	(Matrx Vi 20x18")	1
10	SP1579925	(Matrx Vi 20x20")	1
11	SP1584184	(Matrx Vi 20x22")	1
12	SP1579926	(Matrx Vi 22x18")	1
13	SP1579927	(Matrx Vi 22x20")	1
14	SP1584247	(Matrx Vi 22x22")	1

(Ultra Low Maxx Seat)

1	-	(Seat Tilt)	1
2	-	(Seat Lifter with Tilt)	1
3	-	(Anterior- / Posterior-Tilt)	1
4	-	(Seat Plates & Side Rails)	1
5	-	(Seat Cushions "Matrx Vi")	1
6	-	(Seat Cushions "Matrx Flo-Tech")	1
7	-	(Seat Cushions "Matrx Libra")	1

(Seat Tilt)

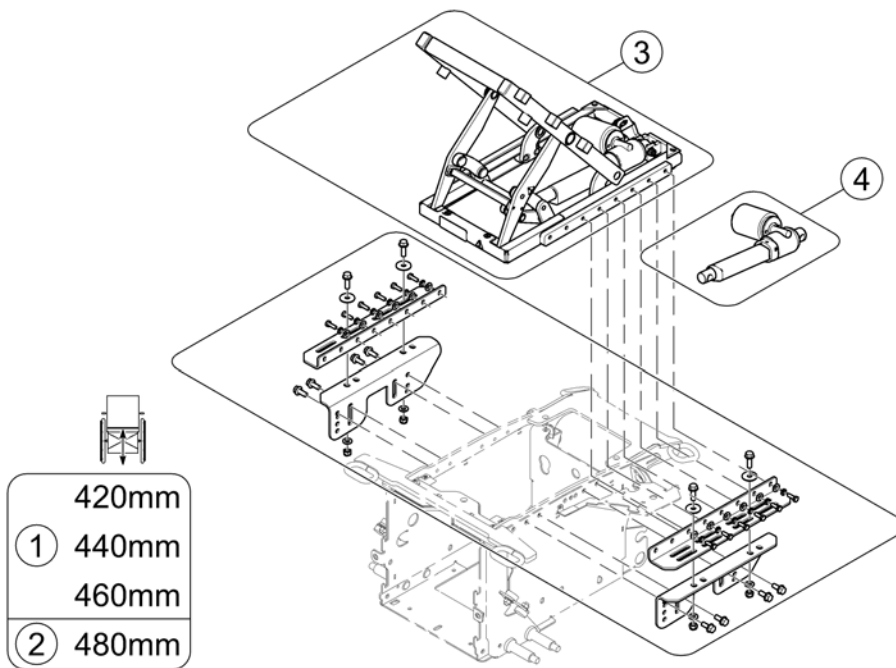
1588100 Rev




1	SPTRC2038	(Base Interface Kit STF 420/440/460mm)	1
2	SPTRC2103	(Base Interface Kit STF 480mm)	1
3	SP1585523	(Tilt Module, cpl.)	1
4	SP1588440	(Tilt Actuator, cpl. with hardware)	1
-	SPTRE0755	(Cable, Tilt Actuator)	1

(Seat Lifter with Tilt)

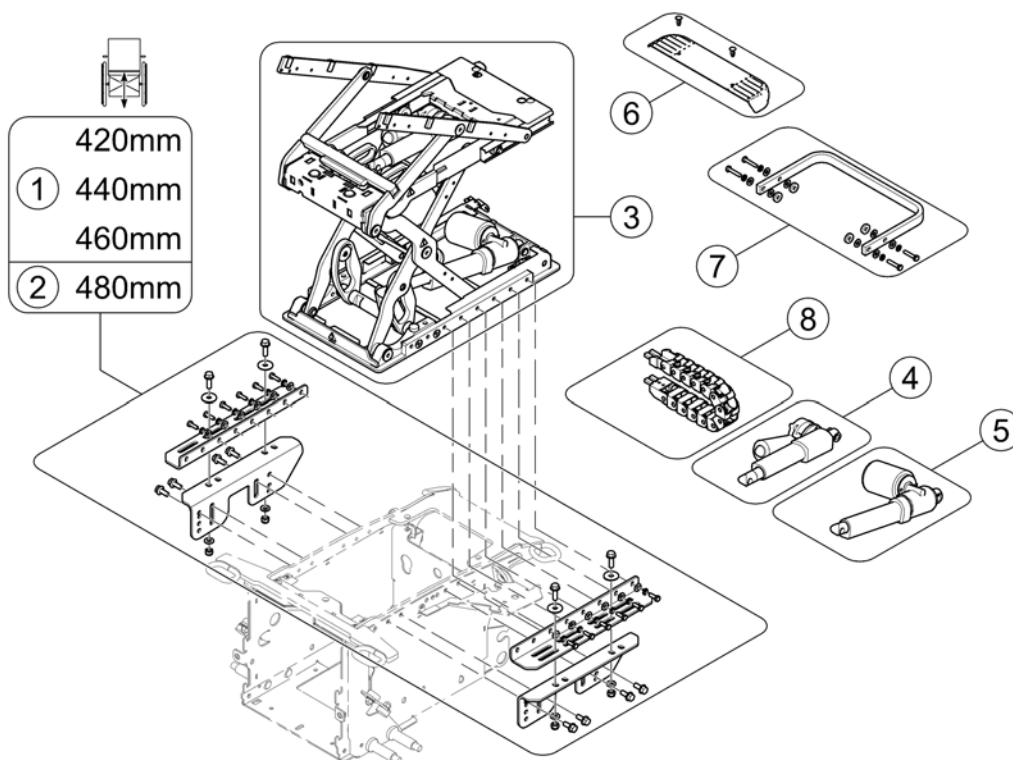1588151.00




-	SP1602307	(Wiring Kit Lifter) <i>replaces item SPTRW0289</i>	1
1	SPTRC2038	(Base Interface Kit STF 420/440/460mm)	1
2	SPTRC2103	(Base Interface Kit STF 480mm)	1
3	SP1585524	(Lifter/Tilt Module, cpl.)	1
4	SP1588442	(Lifter actuator, cpl. with hardware)	1
5	SP1588441	(Tilt actuator, cpl. with hardware)	1
6	SPTRC2071	(Cover, Rear Wiring Lifter, cpl.)	1
7	SPTRC2111	(Bumper with mounting hardware)	1
8	SP1604199	(Energie chain)	1

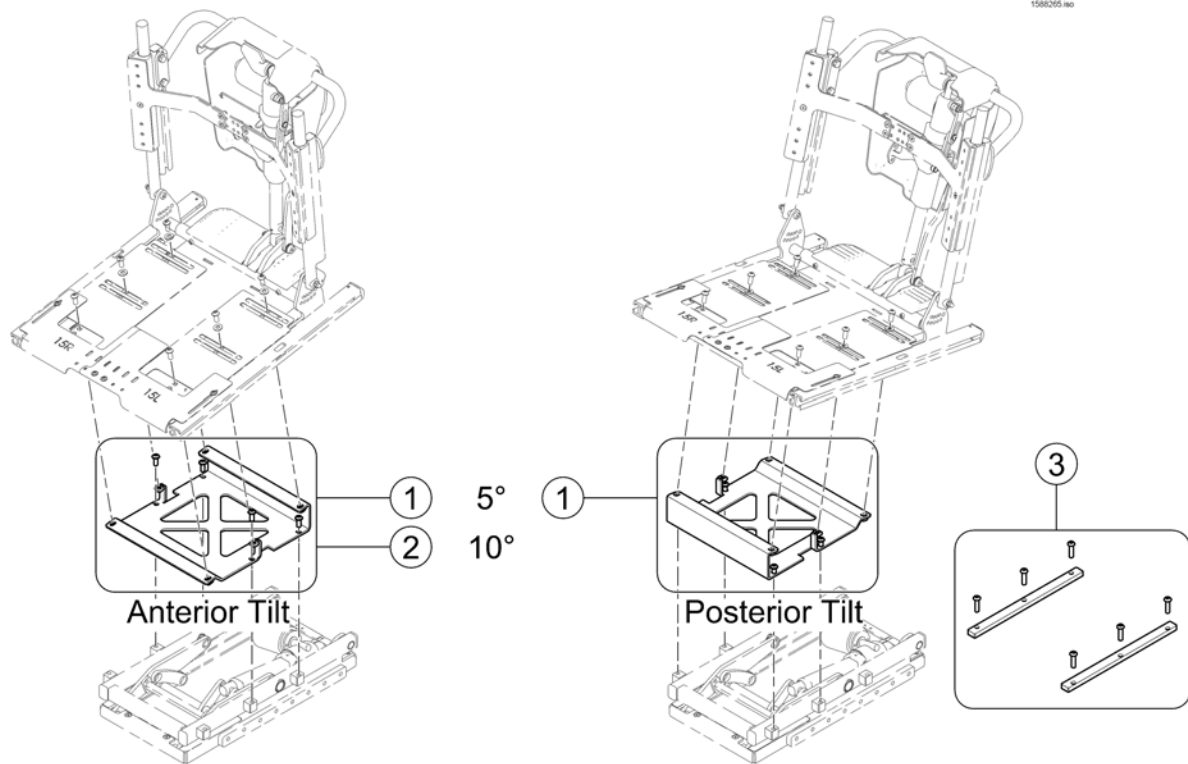

1	SPTRC1923	(Pre-Tilt Mounting Kit 5°) <i>Anterior-/Posterior-Tilt</i>	1
2	SPTRC1924	(Pre-Tilt Mounting Kit 10°) <i>Anterior-Tilt</i>	1
3	SP1589797	(Seat Adapter, cpl. (only for Low Rider option))	1

(Seat Plates & Side Rails)

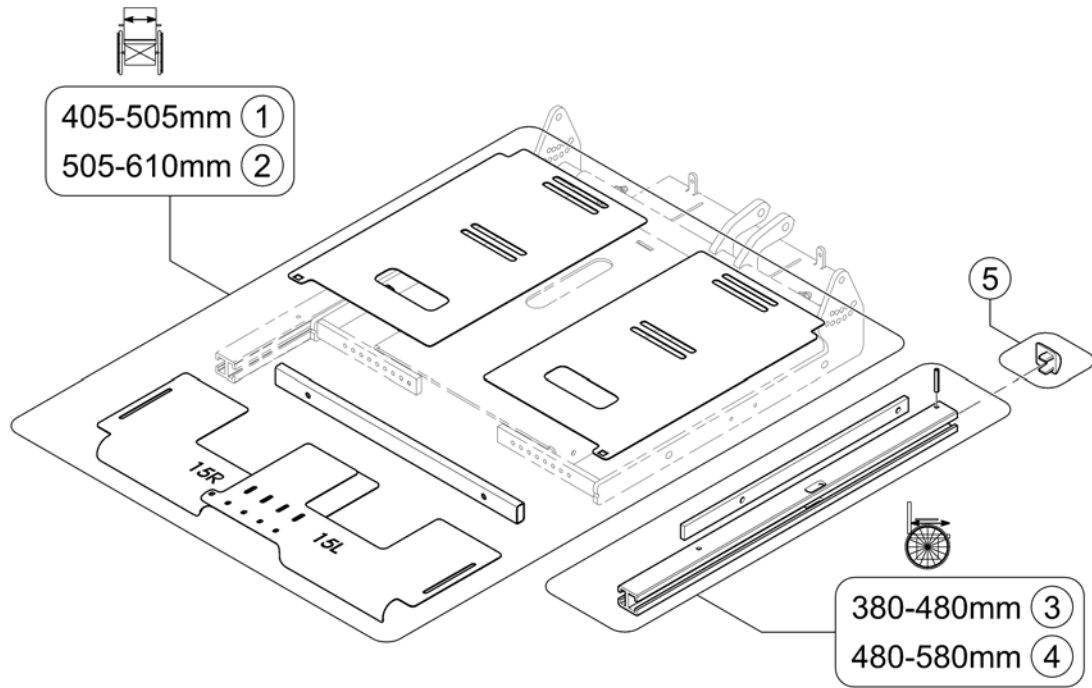
1588152 Rev



1	SPTRC2187	(Seat Pans - SW 405-505mm (16-20"))	1
2	SPTRC2189	(Seat Pans - SW 505-610mm (19-24"))	1
3	SPTRC2141	(Seat Rail Kit SD 380-480mm, cpl.)	1
4	SPTRC2142	(Seat Rail Kit SD 480-580mm, cpl.)	1
5	SPTRC2090	(End cap kit (2 pcs))	1

(Seat Cushions "Matrx Vi")

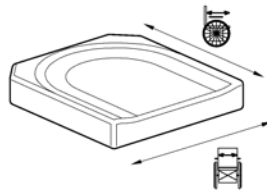1588153.iss




Matrx Vi


		mm									
		405	430	455	480	505	530	555	580	610	
mm	380	1	4	4	8	8	-	-	-	-	
	405	1	4	4	8	8	-	-	-	-	
	430	1	4	4	8	8	-	-	-	-	
	455	2	5	5	9	9	12	12	15	15	
	480	2	5	5	9	9	12	12	15	15	
	505	3	6	6	10	10	13	13	16	16	
	530	3	6	6	10	10	13	13	16	16	
	555	-	7	7	11	11	14	14	17	17	
	580	-	7	7	11	11	14	14	17	17	


1	SP1579786	(Matrx Vi 16x16")	1
2	SP1579787	(Matrx Vi 16x18")	1
3	SP1579788	(Matrx Vi 16x20")	1
4	SP1579920	(Matrx Vi 18x16")	1
5	SP1579921	(Matrx Vi 18x18")	1
6	SP1579922	(Matrx Vi 18x20")	1
7	SP1584183	(Matrx Vi 18x22")	1
8	SP1579923	(Matrx Vi 20x16")	1
9	SP1579924	(Matrx Vi 20x18")	1
10	SP1579925	(Matrx Vi 20x20")	1
11	SP1584184	(Matrx Vi 20x22")	1
12	SP1579926	(Matrx Vi 22x18")	1
13	SP1579927	(Matrx Vi 22x20")	1
14	SP1584247	(Matrx Vi 22x22")	1
15	SP1584384	(Matrx Vi 24x18")	1
16	SP1584417	(Matrx Vi 24x20")	1
17	SP1584418	(Matrx Vi 24x22")	1

(Seat Cushions "Matrix Flo-Tech")


1588508.iss

						
		405	430	455	480	505
Matrix Flo-Tech Lite	380	1	1	5	5	9
	405	2	2	6	6	9
	430	2	2	6	6	10
	455	3	3	7	7	10
	480	3	3	7	7	-
	505	4	4	8	8	-
	530	4	4	8	8	-

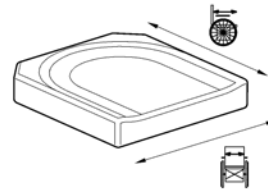

											
		405	430	455	480	505	530	555	580	610	
Matrix Flo-Tech Contour	380	11	11	15	15	19	-	-	-	-	
	405	12	12	16	16	19	-	-	-	-	
	430	12	12	16	16	20	20	24	24	27	
	455	13	13	17	17	21	21	24	24	27	
	480	13	13	17	17	21	21	25	25	28	
	505	14	14	18	18	22	22	25	25	28	
	530	14	14	18	18	22	22	26	26	29	
	555	-	-	-	-	23	23	26	26	29	

1	FLB17/15	(Flo-Tech Lite 17X15")	1
2	FLB17/16	(Flo-Tech Lite 17X16")	1
3	FLB17/18	(Flo-Tech Lite 17X18")	1
4	FLB17/20	(Flo-Tech Lite 17X20")	1
5	FLB19/15	(Flo-Tech Lite 19x15")	1
6	FLB19/16	(Flo-Tech Lite 19x16")	1
7	FLB19/18	(Flo-Tech Lite 19x18")	1
8	SPFLB19/20	SETEPUTE SB48 SD51 SORT STORM4	1
9	SP1584714	(Flo-Tech Lite 20x15")	1
10	SP1584713	(Flo-Tech Lite 20x17")	1
11	FCB17/15	(Flo-Tech Contour 17x15")	1
12	SPFCB17_16	(Flo-Tech Contour 17x16")	1
13	SPFCB17_18	(Flo-Tech Contour 17x18")	1
14	SPFCB17_20	(Flo-Tech Contour 17x20")	1
15	SPFCB19_15	(Flo-Tech Contour 19x15")	1
16	SPFCB19_16	(Flo-Tech Contour 19x16")	1
17	SPFCB19_18	(Flo-Tech Contour 19x18")	1
18	SPFCB19_20	(Flo-Tech Contour 19x20")	1
19	SPFCB20_15	(Flo-Tech Contour 20x15")	1
20	SPFCB21/17	SETEPUTE SB53 SD43 SORT STORM4	1
21	SPFCB21/18	SETEPUTE SB53 SD46 SORT STORM4	1
22	SPFCB21/20	SETEPUTE SB53 SD51 SORT STORM4	1
23	SP1584739	(Flo-Tech Contour 21x22")	1
24	SP1584740	(Flo-Tech Contour 23x17")	1
25	SP1584741	(Flo-Tech Contour 23x19")	1
26	SP1584742	(Flo-Tech Contour 23x21")	1
27	SP1584743	(Flo-Tech Contour 24x17")	1
28	SP1584770	(Flo-Tech Contour 24x19")	1
29	SP1584771	(Flo-Tech Contour 24x21")	1

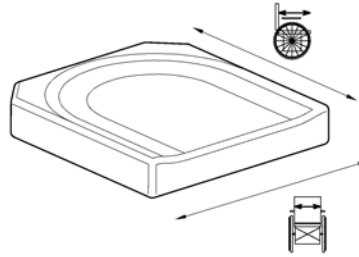


Matrx Libra

		mm						
		405	430	455	480	505	530	555
mm	405	1	4	7	10	13	-	-
	430	1	4	7	10	13	-	-
	455	2	5	8	11	14	16	19
	480	2	5	8	11	14	16	19
	505	3	6	9	12	15	17	20
	530	3	6	9	12	15	17	20
	555	-	-	-	-	-	18	21
	580	-	-	-	-	-	18	21


1	LC1616	(Matrx Libra cushions 410x410 mm)	1
2	LC1618	(Matrx Libra cushions 410x460 mm)	1
3	LC1620	(Matrx Libra cushions 410x510 mm)	1
4	LC1716	(Matrx Libra cushions 430x410 mm)	1
5	LC1718	(Matrx Libra cushions 430x460 mm)	1
6	LC1720	(Matrx Libra cushions 430x510 mm)	1
7	LC1816	(Matrx Libra cushions 460x410 mm)	1
8	LC1818	(Matrx Libra cushions 460x460 mm)	1
9	LC1820	(Matrx Libra cushions 380x510 mm)	1
10	LC1916	(Matrx Libra cushions 480x410 mm)	1
11	LC1918	(Matrx Libra cushions 480x460 mm)	1
12	LC1920	(Matrx Libra cushions 480x510 mm)	1
13	LC2016	(Matrx Libra cushions 510x410 mm)	1
14	LC2018	(Matrx Libra cushions 510x460 mm)	1
15	LC2020	(Matrx Libra cushions 510x510 mm)	1
16	LC2118	(Seat Cushion Matrx Libra 21x18")	1
17	LC2120	(Seat Cushion Matrx Libra 21x20")	1
18	LC2122	(Seat Cushion Matrx Libra 21x22")	1
19	LC2218	(Matrx Libra cushions 560x460 mm)	1
20	LC2220	(Matrx Libra cushions 560x510 mm)	1
21	LC2222	(Seat Cushion Matrx Libra 22x22")	1

(Max Seat)

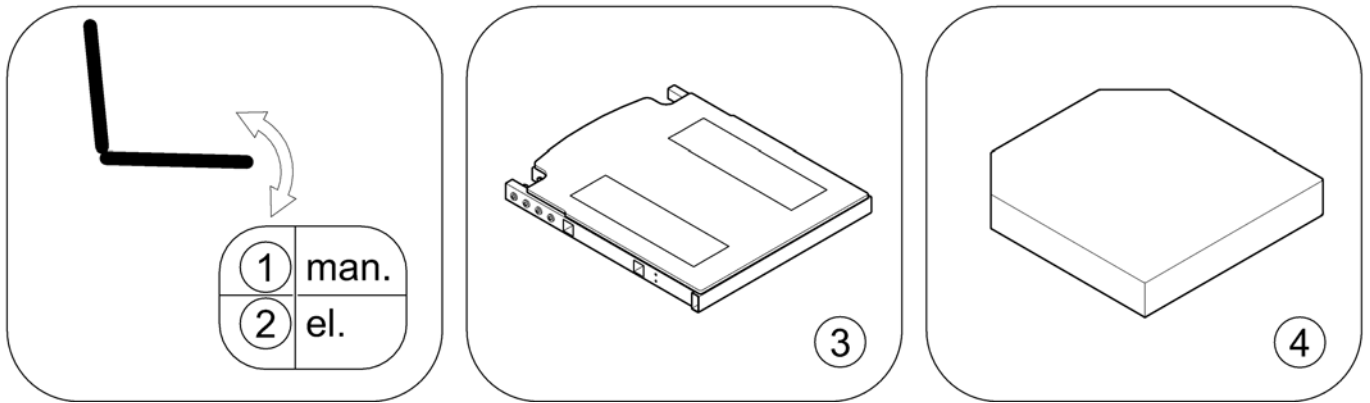1582428.no




1	-	(man. settable Tilt)	1
2	-	(powered simple Tilt)	1
3	-	(Seat Frame and Plate)	1
4	-	(Seat Cushion)	1

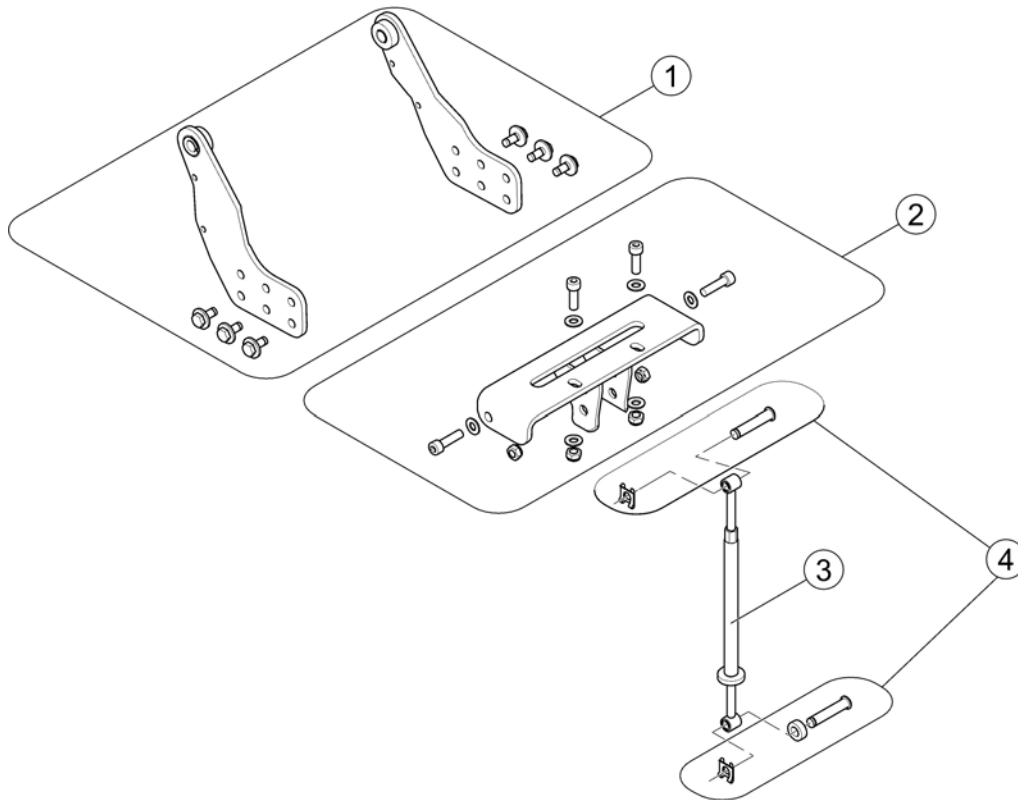

1	SP1582746	(Pivot Bearing Back - Max Seat)	1
2	SP1582747	(Support Plate Tilt - Max Seat)	1
3	1420003	SETELØFT MAN SPINDEL STORM3	1
4	SP1582737	(Mounting Hardware for Tilt)	1

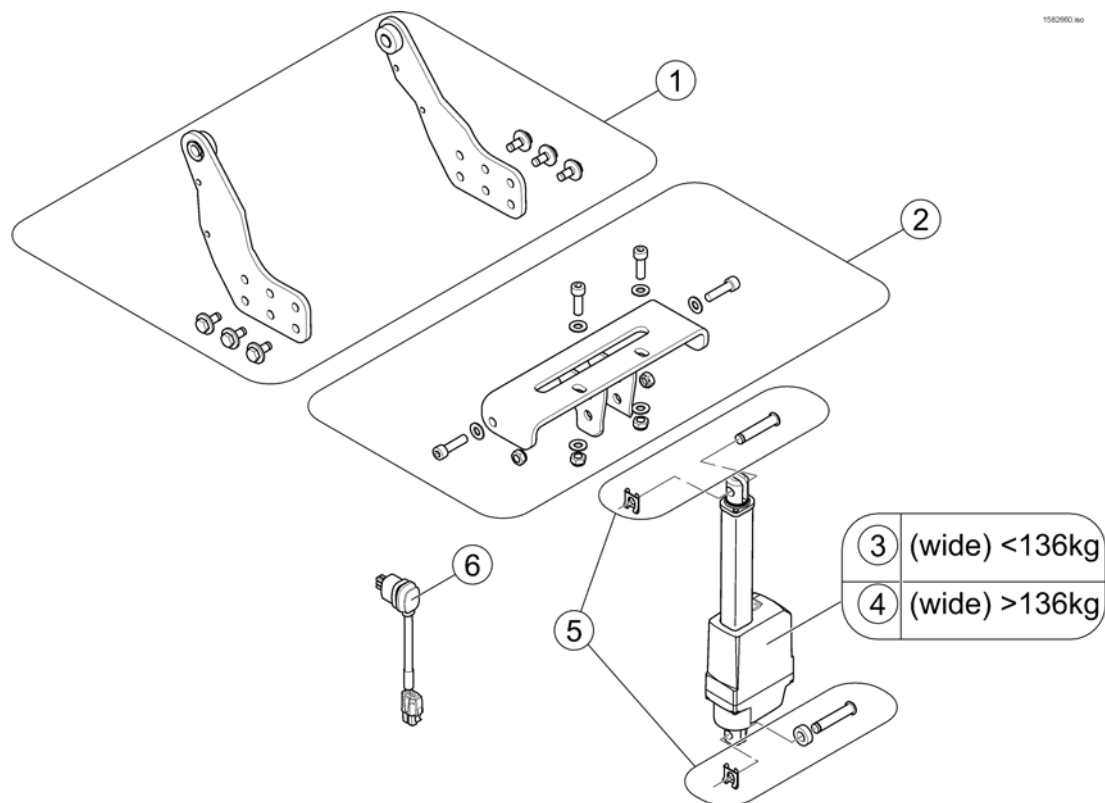
(powered simple Tilt)



1	SP1582746	(Pivot Bearing Back - Max Seat)	1
2	SP1582747	(Support Plate Tilt - Max Seat)	1
3	SP1581020	(Actuator LA23 - stroke 135mm / 1800N)	1
4	SP1556487	(Actuator LA23 - stroke 135mm / 2500N)	1
5	SP1582737	(Mounting Hardware for Tilt)	1
6	SP1572166	KABEL ACTUATOR MODULITE	1

(Seat Frame and Plate)

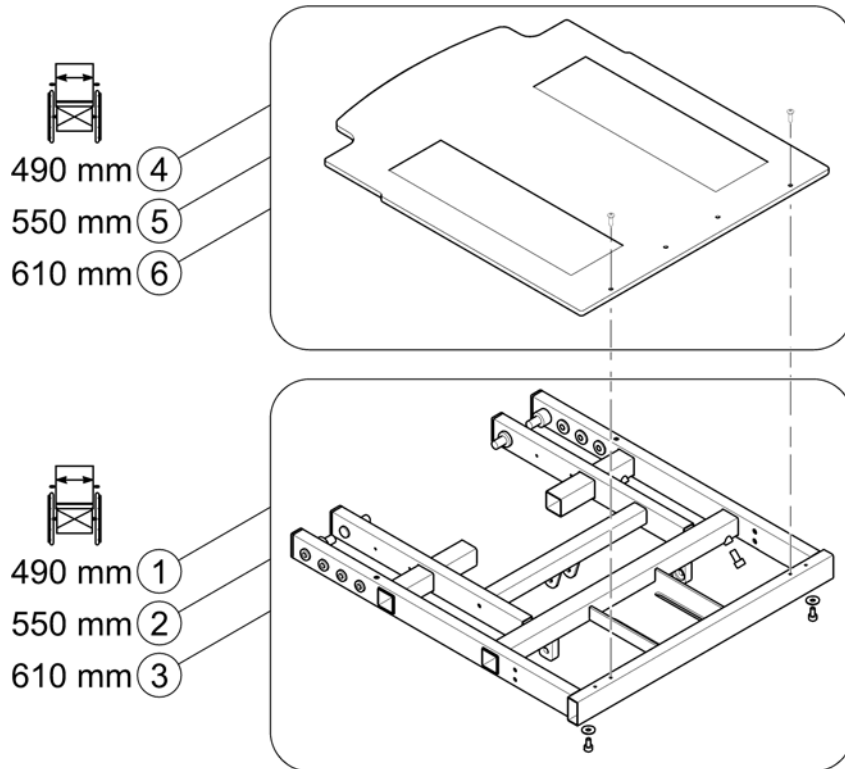
1582061.000



1	SP1582748	(Seat Frame SW 490mm)	1
2	SP1582749	(Seat Frame SW 550mm)	1
3	SP1582772	(Seat Frame SW 610mm)	1
4	SP1508137	(Seat plate SW 490mm)	1
5	SP1525083	(Seat plate SW 550mm)	1
6	SP1525082	(Seat plate SW 610mm)	1

(Seat Cushion)

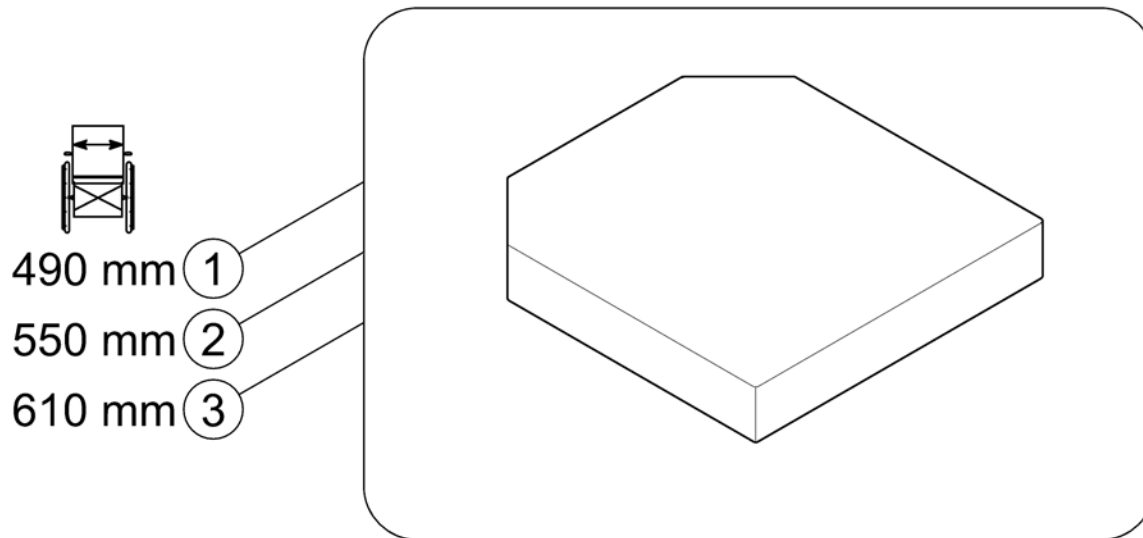1582666.iso




1	SP1508059	SETEPUTE FLOSHAPE SB50 STORM4 MAX DARTEX SD53	1
2	SP1536595	SETEPUTE SB55 SD65 MAX STORM4 MAX	1
3	SP1536596	SETEPUTE SB61 SD65 MAX STORM4 MAX	1



1	-	(Backrest for Modulite Seat)	1
2	-	(Backrest for Ultra Low Maxx Seat)	1
3	-	(Backrest for Max Seat)	1


(Backrest for Modulite Seat)

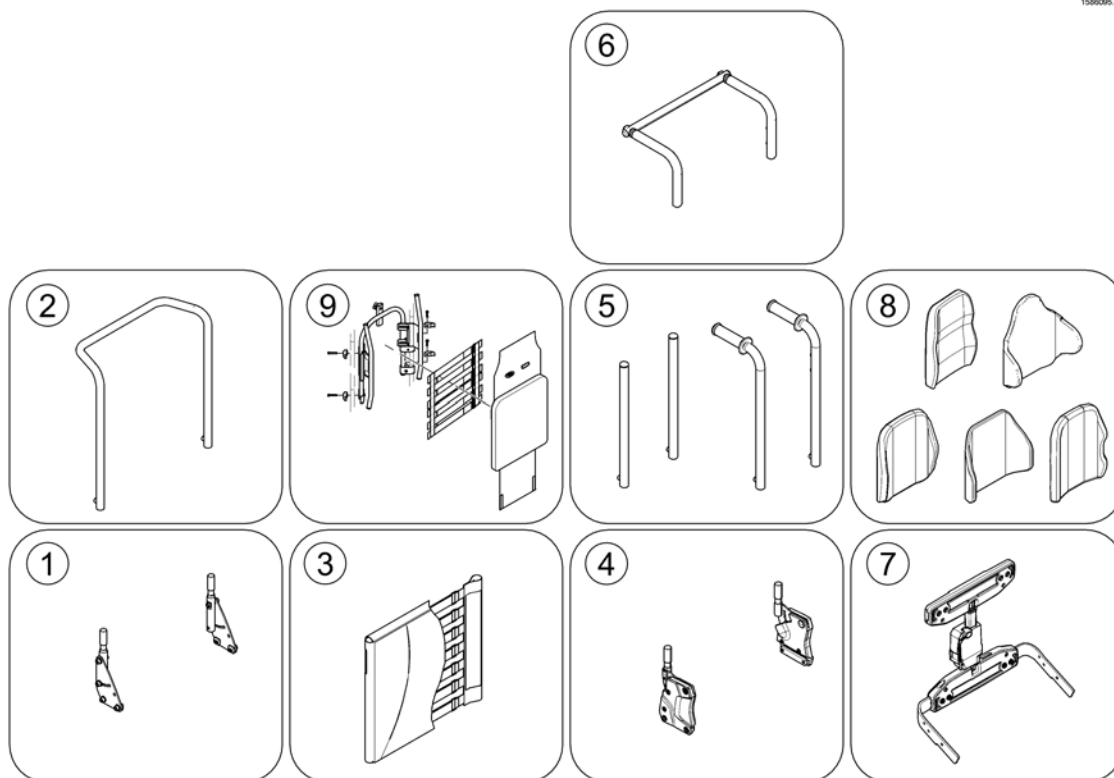1566095.no



1	-	(Support for Backrest Frame) [0975+2161]	1
2	-	(Backrest Frame) [2161]	1
3	-	(Sling Back, tension adjustable) [2171]	1
4	-	(Backrest support)	1
5	-	(Backposts) [2162/2163]	1
6	-	(Connection tube) [1307]	1
7	-	(Back Recliner / Stabilization bracket) [0977]	1
8	-	(Backrest Cushions)	1
9	-	(Flex3 back, tension adjustable) [2179]	1

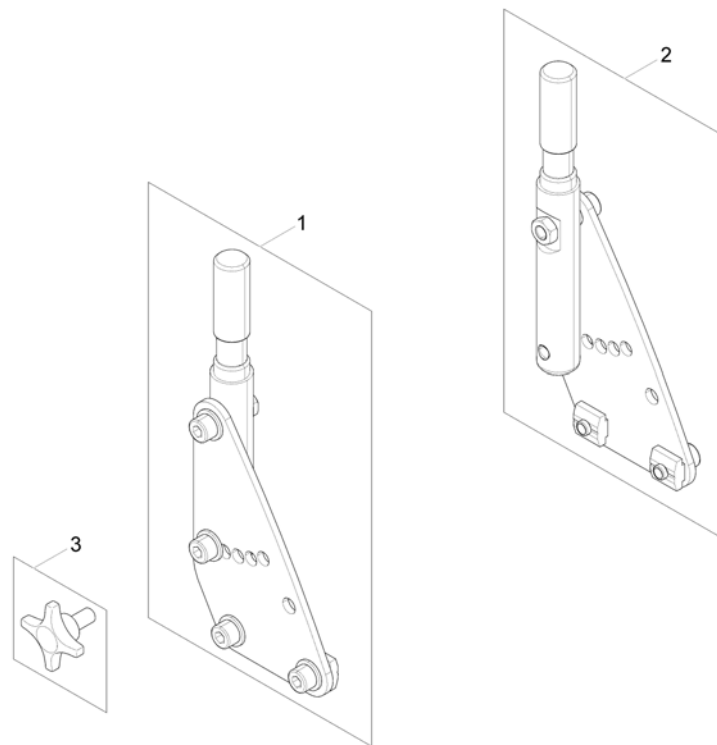
1571553.no



1	SP1571542	(Backrest Support RH)	1
2	SP1571553	(Backrest Support LH)	1
3	5824117	STJERNEGREP FESTE NY MANØVERBOKS MOOVER <i>(only swedish homologation)</i>	2


(Backrest Frame)

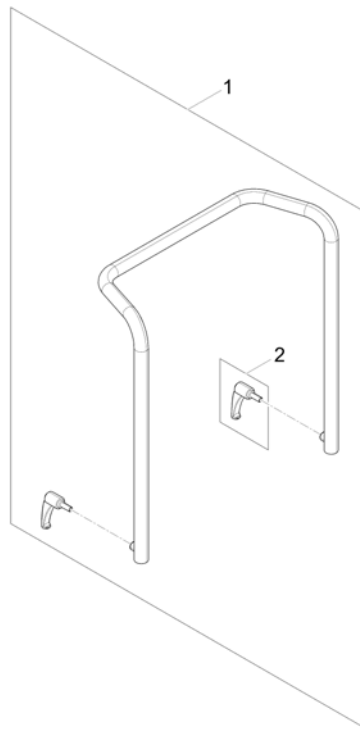1571831.no



1	SP1567810	RYGGRØR TDX MODULITE	1
2	0021-70002	FINGERSKRUE M8 ITALIA RYGGRØR	1

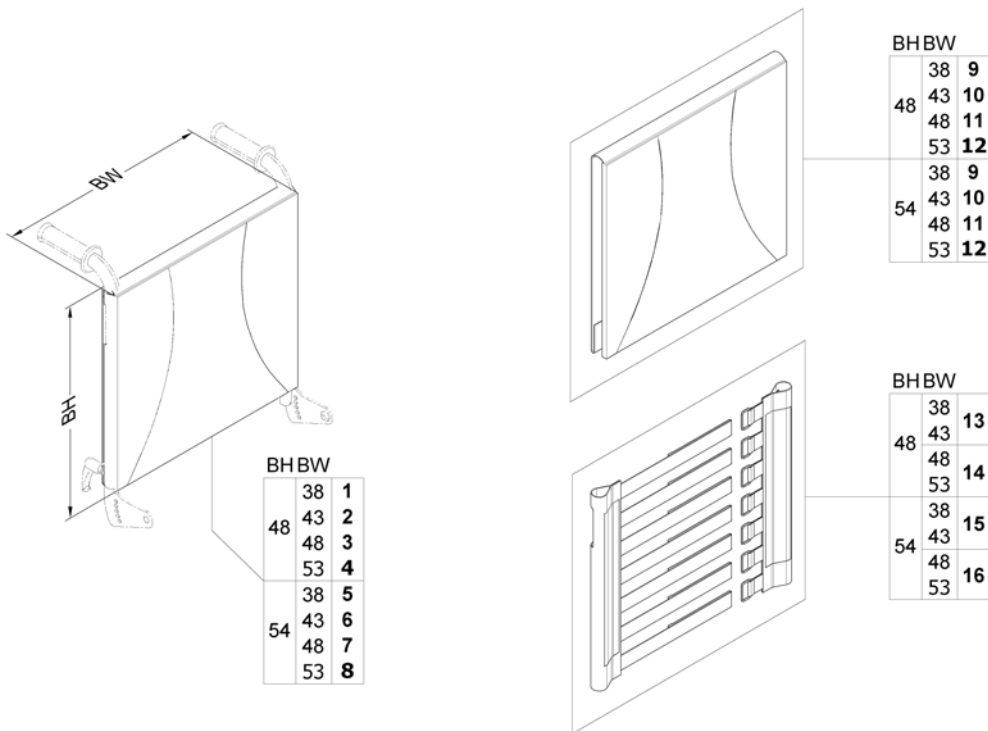
1571829.no



1	SP1567529	(Sling Back, tension adjustable & Upholstery, BW 380mm / BH 480mm, cpl.)	1
2	SP1567465	VELCRO POLSTER HØ48 SB43 MODULITE	1
3	SP1567530	RYGG VELCRO SB48 RH48 MODULITE	1
4	SP1567531	(Sling Back, tension adjustable & upholstery, BW 530mm / BH 480mm, cpl.)	1
5	SP1567627	(Sling Back, tension adjustable & upholstery, BW 380mm / BH 540mm. cpl.)	1
6	SP1567746	RYGG VELCRO SB43 RH43 MODULITE	1
7	SP1567747	RYGG VELCRO SB48 RH54 MODULITE	1
8	SP1567748	VELCRORYGG KPL RH54 RB53 MODULITE	1
9	SP1563563	(Backrest upholstery)	1
10	SP1563653	(Backrest upholstery)	1
11	SP1563654	(Backrest upholstery)	1
12	SP1563655	(Backrest upholstery)	1
13	SP1571492	(Sling back)	1
14	SP1571493	(Sling back)	1
15	SP1571494	(Sling back)	1
16	SP1571510	(Sling back)	1


(Backrest support)

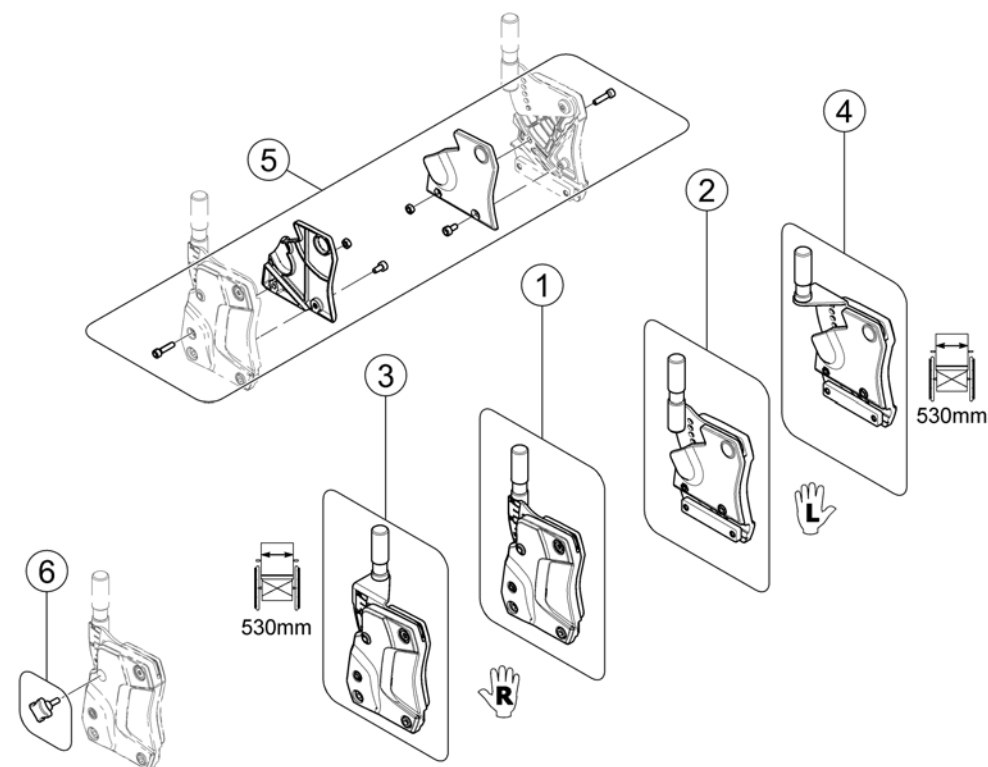1572928 Rev



-	-	(Design change on 02. May 2016 - we suggest to change older supports on both sides.)		1
1	SP1582709	FESTE RYGGRØR HØ MODULITE <i>Compatible with Flip-Up & Following Armrest [2193/2192/2194] - replaces SP1571560</i>		1
2	SP1582708	FESTE RYGGRØR VE MODULITE <i>Compatible with Flip-Up & Following Armrest [2193/2192/2194] - replaces SP1571561</i>		1
3	SP1582952	RYGGFESTE 53 HØ MODULITE <i>Compatible with Flip-Up & Following Armrest [2193/2192/2194] - replaces SP1571568</i>		1
4	SP1582951	RYGGFESTE 53 VE MODULITE <i>Compatible with Flip-Up & Following Armrest [2193/2192/2194] - replaces SP1571569</i>		1
5	SP1571676	(Covers)		1
6	SP1579877	(Cross Handle M6) <i>(Swedish Homologation only!)</i>	→ 31.07.2015	2
6	SP1589324	(Cross handle M8) <i>(Swedish Homologation only!)</i>	01.08.2015 →	2
7	SP1590656	(Hole plug (for powered backrest only))	01.07.2016 →	1

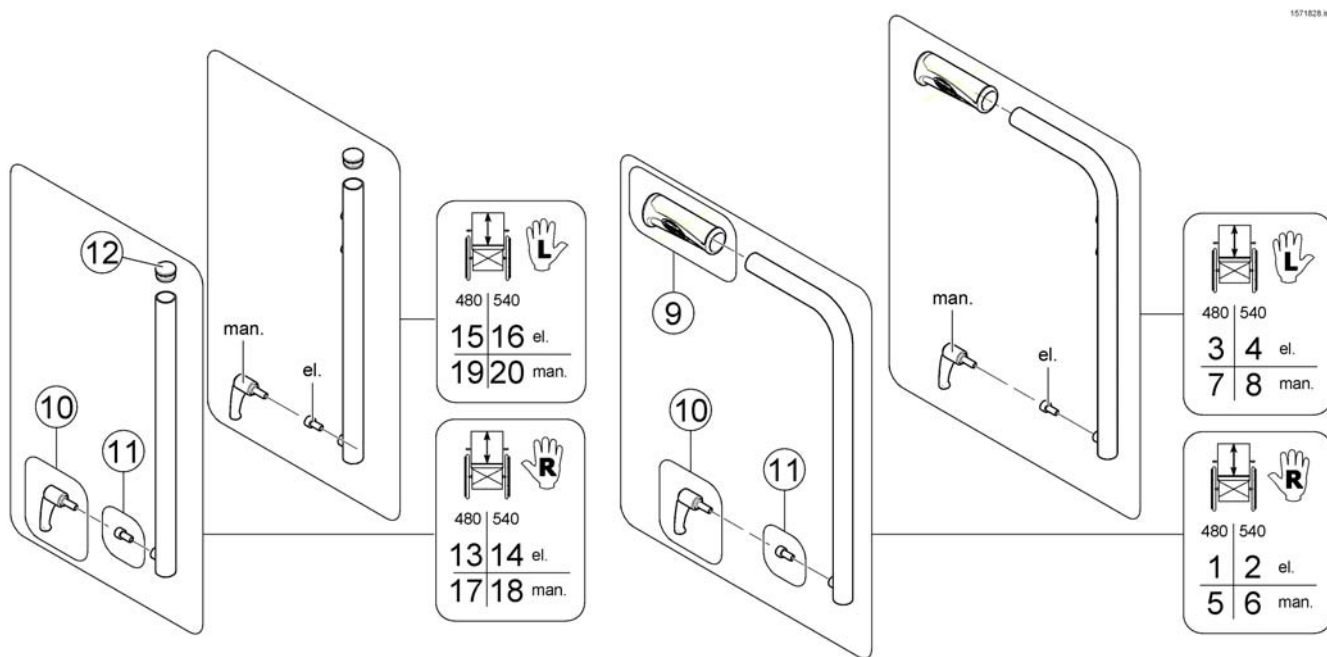
(Backposts)



1	SP1587177	RYGGRØR HØYDE 48 HØ MODULITE <i>Compatible with Flip-Up & Following Armrest [2193/2192/2194] - replaces SP1563102</i>	1
2	SP1587184	RYGGRØR HØ 54 MODULITE <i>Compatible with Flip-Up & Following Armrest [2193/2192/2194] - replaces SP1563103</i>	1
3	SP1587176	RYGGRØR HØYDE 48 VE MODULITE <i>Compatible with Flip-Up & Following Armrest [2193/2192/2194] - replaces SP1560576</i>	1
4	SP1587183	RYGGRØR VE 54 MODULITE <i>Compatible with Flip-Up & Following Armrest [2193/2192/2194] - replaces SP1559578</i>	1
5	SP1587163	(Backpost man. BH 480mm - RH) <i>Compatible with Flip-Up & Following Armrest [2193/2192/2194] - replaces SP1572489</i>	1
6	SP1587182	(Backpost man. BH 540mm - RH) <i>Compatible with Flip-Up & Following Armrest [2193/2192/2194] - replaces SP1572490</i>	1
7	SP1587162	(Backpost man. BH 480mm - LH) <i>Compatible with Flip-Up & Following Armrest [2193/2192/2194] - replaces SP1572491</i>	1
8	SP1587178	(Backpost man. BH 540mm - LH) <i>Compatible with Flip-Up & Following Armrest [2193/2192/2194] - replaces SP1572492</i>	1
9	SP1586122	(Handles (pair)) <i>replaces item SP1563651</i>	1
10	0021-70002	FINGERSKRUE M8 ITALIA RYGGRØR	1
11	1567999	(Bolt M8x16 8.8 black DIN912)	1
13	SP1606879	(Backpost el. BH 480mm - RH w/o Pushhandle)	1
14	SP1606893	(Backpost el. BH 540mm - RH w/o Pushhandle)	1
15	SP1606880	(Backpost el. BH 480mm - LH w/o Pushhandle)	1
16	SP1606894	(Backpost el. BH 540mm - LH w/o Pushhandle)	1
17	SP1606877	(Backpost man. BH 480mm - RH w/o Pushhandle)	1
18	SP1606881	(Backpost man. BH 540mm - RH w/o Pushhandle)	1
19	SP1606878	(Backpost man. BH 480mm - LH w/o Pushhandle)	1
20	SP1606882	(Backpost man. BH 540mm - LH w/o Pushhandle)	1

(Connection tube)

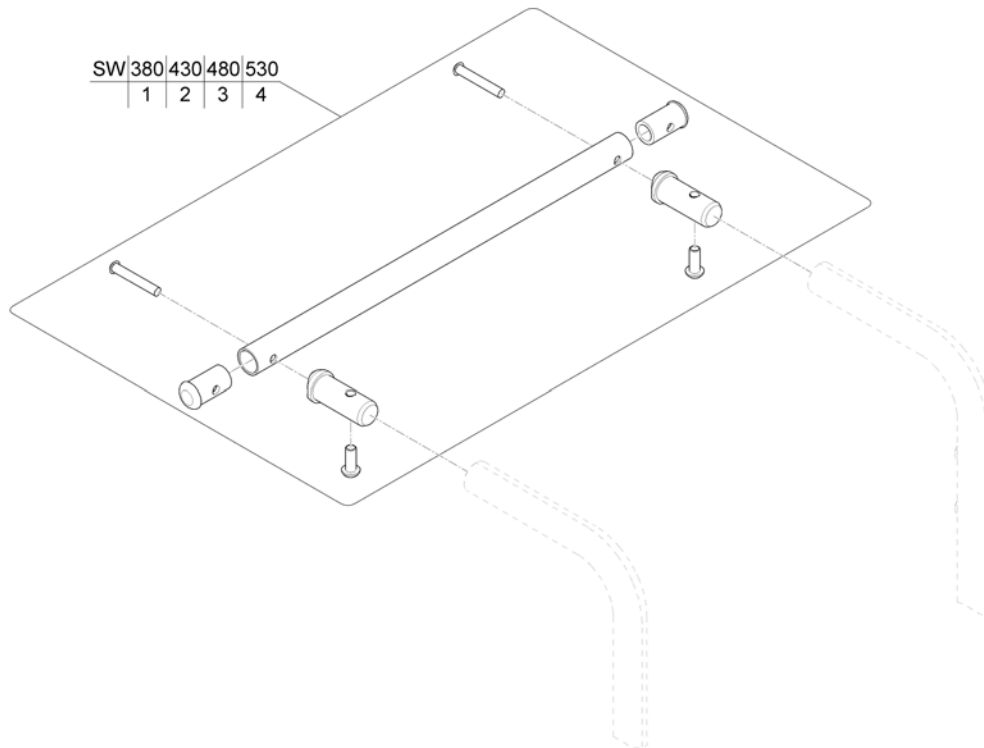
1571835 no



1	SP1570370	KJØREBØYLE MODULITE SB38 MODULITE	1
2	SP1570371	KJØREBØYLE SB43 MODULITE STORM4	1
3	SP1570372	BRAKETT RYGG SB48 LS1 MODULITE	1
4	SP1570373	KJØREBØYLE SB53 MODULITE STORM4 TDXSP2	1

Invacare® TDX-SP2 Lowrider, TDX-SP2 NB | Rygg
(Back Recliner / Stabilization bracket)

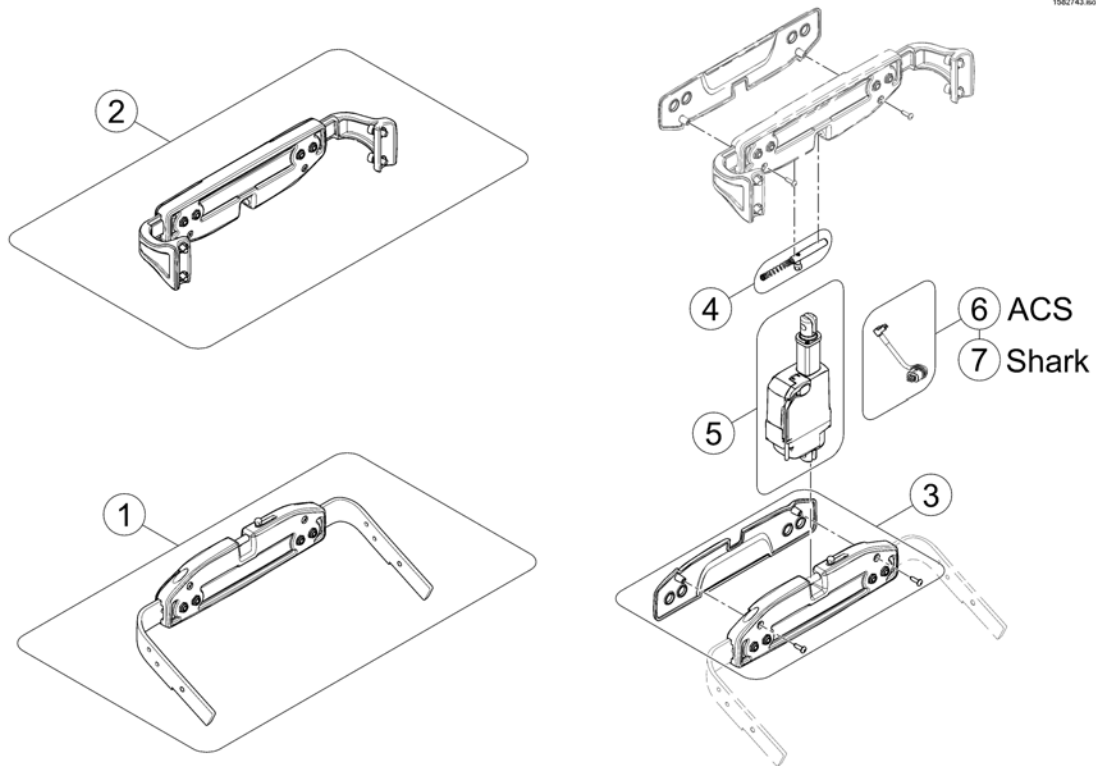
1582743.kno




1	SP1567169	ACTUATOR BRAKETT MODULITE	1
2	SP1567168	BRAKETT STABILISATOR RYG STORM4 MODULITE	1
3	SP1571702	(Cover)	1
4	SP1571704	(Locking bolt)	1
5	SP1557820	ACTUATOR RYGG MODULITE	1
6	SP1572166	KABEL ACTUATOR MODULITE	1
7	SP1552297	(Actuator Cable LA23 Shark (~500mm))	1

(Backrest Cushions)

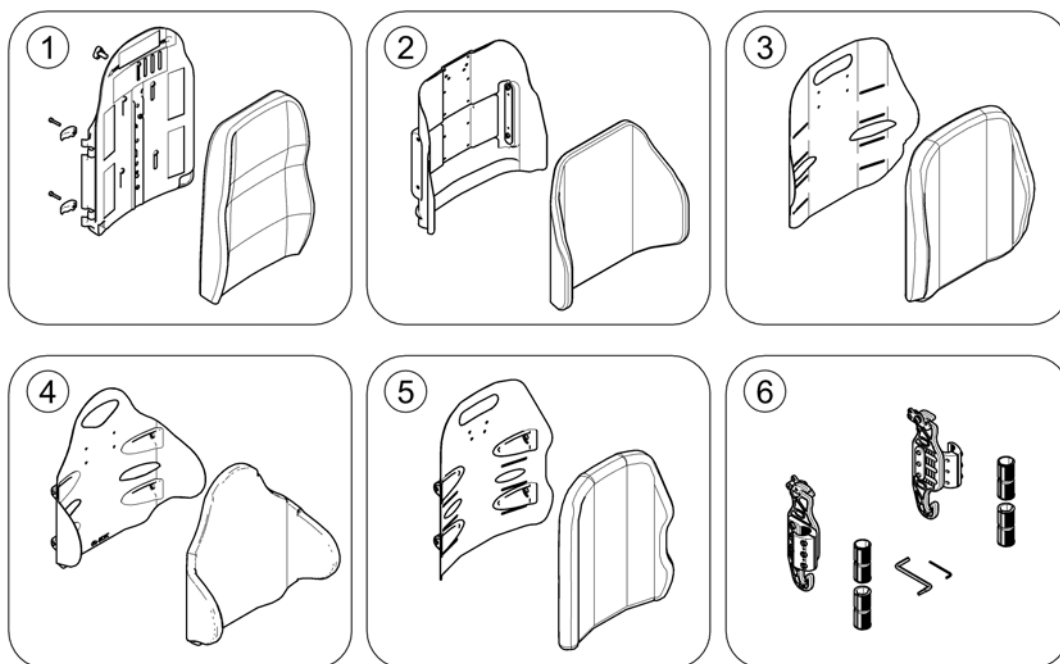1588915.no




1	-	(Flex 3 / Laguna)		1
2	-	(Matrx Posture Back)		1
3	-	(Matrx Elite T/R Back)	01.07.2015 →	1
4	-	(Matrx Elite Deep Back)	01.07.2015 →	1
5	-	(Matrx Elite (Full) Back)		1
6	-	(Mounting Hardware for Matrx backrests)		1

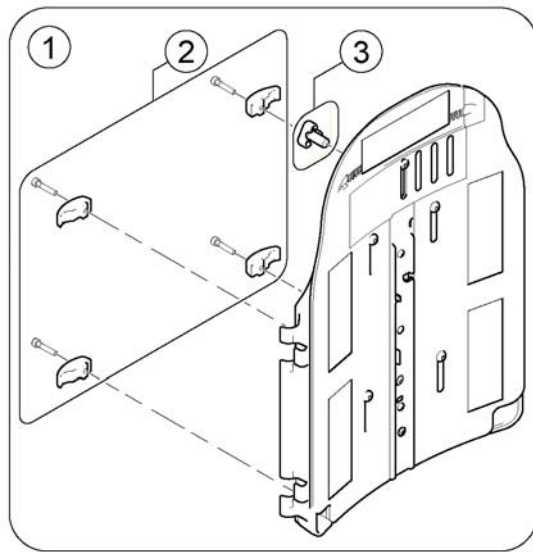

Flex 3

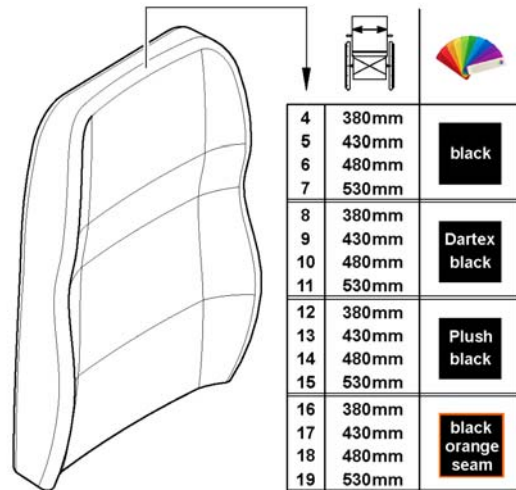
Laguna



1	SP1568489	FLEX3 RYGG TIL MODULITE STORM4 TDX [2172]	1
2	SP1572321	(Clamps)	1
3	1497411	RATT M12 LANG KPL NAKKE/HODESTØTTE	1
4	SP1529585	RYGGPUTE EASYADAPT SB38 STORM4 M SORT [2173+2151]	1
5	SP1529586	RYGGPUTE EASYADAPT SB43 STORM4 STD SORT [2173+2151]	1
6	SP1529587	RYGGPUTE EASYADAPT SB48 STORM4 SORT [2173+2151]	1
7	SP1529596	RYGGPUTE EASYADAPT SB53 STORM4 M SORT [2173+2151]	1
8	SP1529538-TR26	(Backrest Cushion Laguna / EasyAdapt, Dartex black - SW 380mm) [2173+2152]	1
9	SP1529539-TR26	(Backrest Cushion Laguna / EasyAdapt, Dartex black - SW 430mm) [2173+2152]	1
10	SP1529540-TR26	(Backrest Cushion Laguna / EasyAdapt, Dartex black - SW 480mm) [2173+2152]	1
11	SP1529545-TR26	(Backrest Cushion Laguna / EasyAdapt, Dartex black - SW 530mm) [2173+2152]	1
12	SP1529538-TR35	(Backrest Cushion Laguna, Plush black - SW 380mm) [2173+2154]	1
13	SP1529539-TR35	(Backrest Cushion Laguna, Plush black - SW 430mm) [2173+2154]	1
14	SP1529540-TR35	(Backrest Cushion Laguna, Plush black - SW 480mm) [2173+2154]	1
15	SP1529545-TR35	(Backrest Cushion Laguna, Plush black - SW 530mm) [2173+2154]	1
16	SP1529538	RYGGPUTE LAGUNA 37CM SORT M ORANG SØM [2173+2153]	1
17	SP1529539	RYGGPUTE LAGUNA 43CM SORT M ORANG SØM [2173+2153]	1
18	SP1529540	RYGGPUTE LAGUNA 48CM SORT M ORANG SØM [2173+2153]	1
19	SP1529545	RYGGPUTE LAGUNA 50CM SORT M ORANG SØM [2173+2153]	1

(Matrx Posture Back)

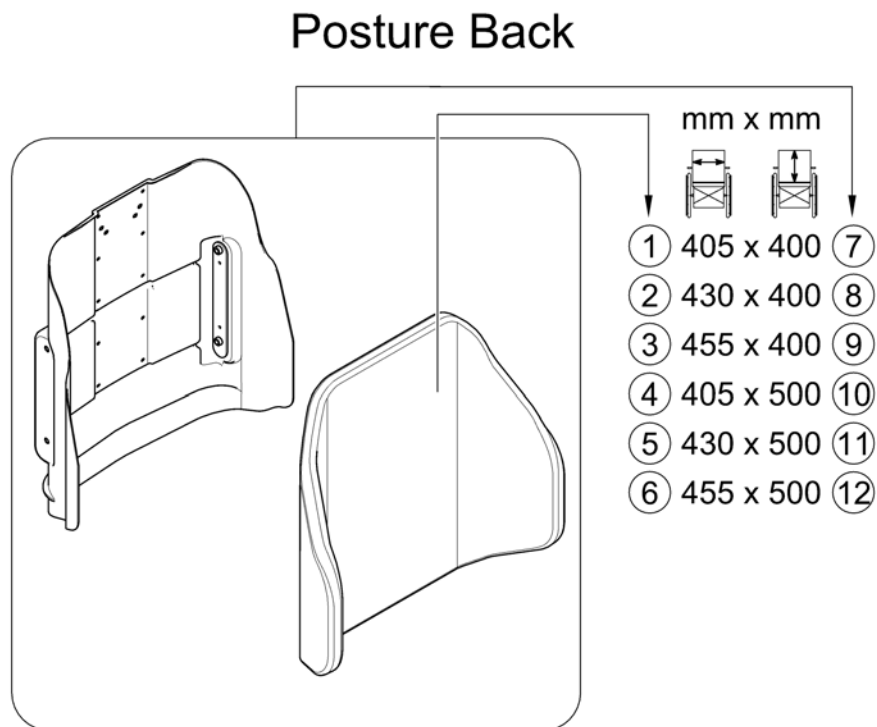158885 no




1	1544896	Ryggtrekk Posture Back s b40,5 rh40,5	1
2	1544899	Ryggtrekk Posture Back t rekk sb43 rh40,5	1
3	1544995	Ryggtrekk Posture Back s b45,5 rh40,5	1
4	1544897	Ryggtrekk Posture Back s b40,5 rh50,5	1
5	1544982	Ryggtrekk Posture Back t rekk sb43 rh50,5	1
6	1544996	Ryggtrekk Posture Back s b45,5 rh50,5	1
7	SP1515727	(Posture Back 405x400mm (16x16"), cpl.) [2174]	1
8	SP1515729	(Posture Back 430x400mm (17x16"), cpl.) [2174]	1
9	SP1515731	(Posture Back 455x400mm (18x16"), cpl.) [2174]	1
10	SP1515728	Elite PB Back 405x500mm MATRIX [2175]	1
11	SP1515730	(Posture Back 430x500mm (17x20"), cpl.) [2175]	1
12	SP1515732	(Posture Back 455x500mm (18x20"), cpl.) [2175]	1

(Matrx Elite T/R Back)

1580006 no

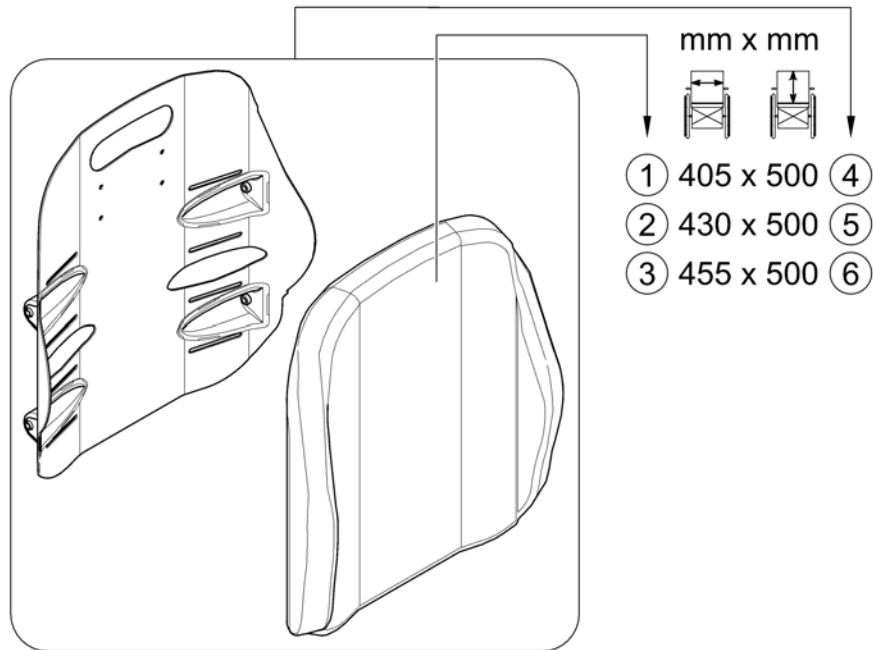Elite T/R Back




1	1580110	Ryggtrekk PB Elite TiltR ecline sb48 rh46	1
2	1580297	Ryggtrekk PB Elite TiltR ecline sb48 rh51	1
3	1580298	Ryggtrekk PB Elite TiltR ecline sb51 rh41	1
4	SP1580029	(Elite T/R Back 405x500mm (16x20"), cpl.) [2177+2121]	1
5	SP1580030	(Elite T/R Back 430x500mm (17x20"), cpl.) [2177+2122]	1
6	SP1580031	(Elite T/R Back 455x500mm (18x20"), cpl.) [2177+2123]	1

(Matrx Elite Deep Back)

158913.no

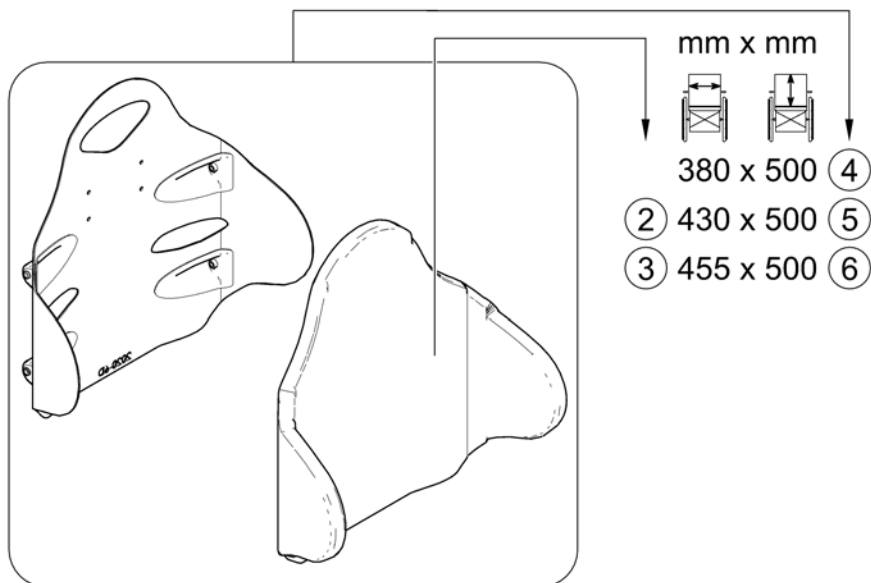Elite Deep Back




2	1544077	Ryggtrekk Elite dyp Back sb43 rh50,5	1
3	1544081	Ryggtrekk Elite dyp Back sb45,5 rh50,5	1
4	SP1543511	(Elite Deep Back 380x500mm (15x20"), cpl.) [2178+2121]	1
5	SP1543519	(Elite Deep Back 430x500mm (17x20"), cpl.) [2178+2122]	1
6	SP1543523	(Elite Deep Back 455x500mm (18x20"), cpl.) [2178+2123]	1

(Matrx Elite (Full) Back)

158914.no

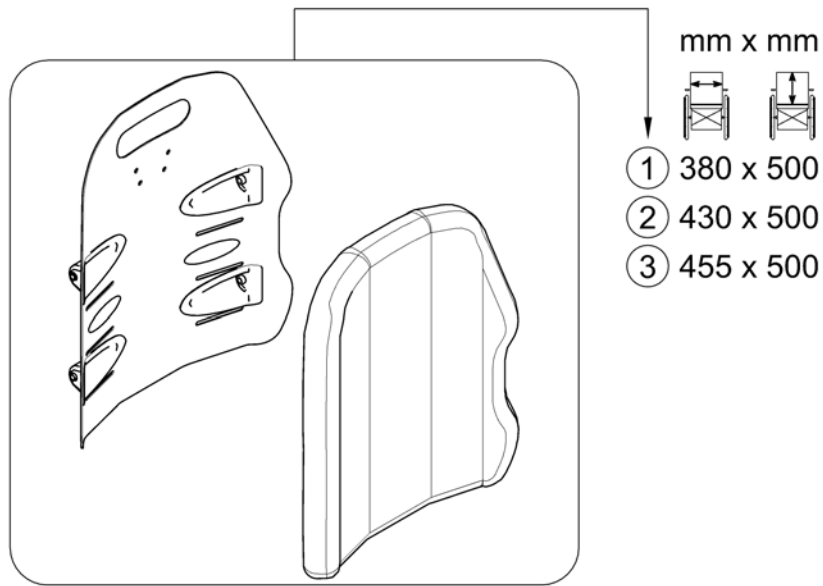
Elite Full Back



1	SP1532094	(Matrx Elite Back 15x20") [2176]	1
2	SP1517755	(Matrx Elite Back 17x20") [2176]	1
3	SP1517756	(Matrx Elite Back 18x20") [2176]	1

(Mounting Hardware for Matrix backrests)

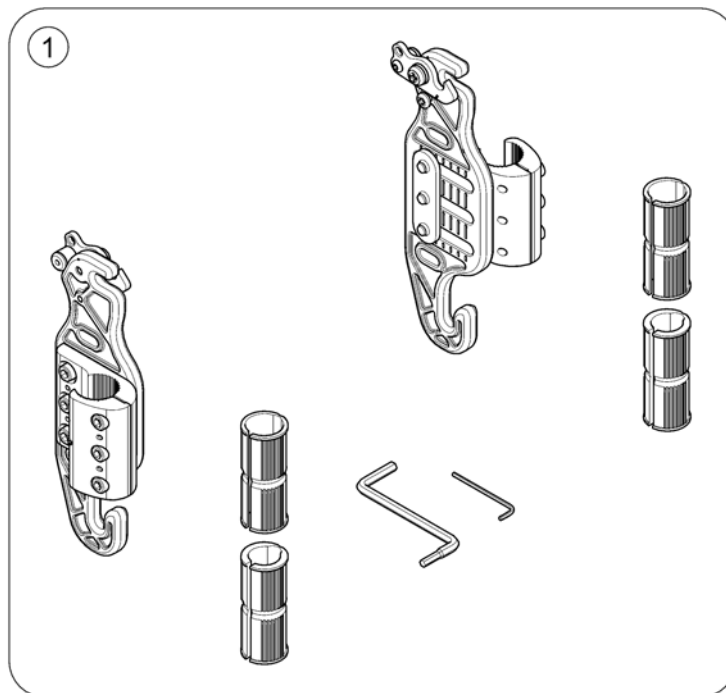
1637607/00



1 SP1545275

(Mounting hardware, Matrix Backrest)

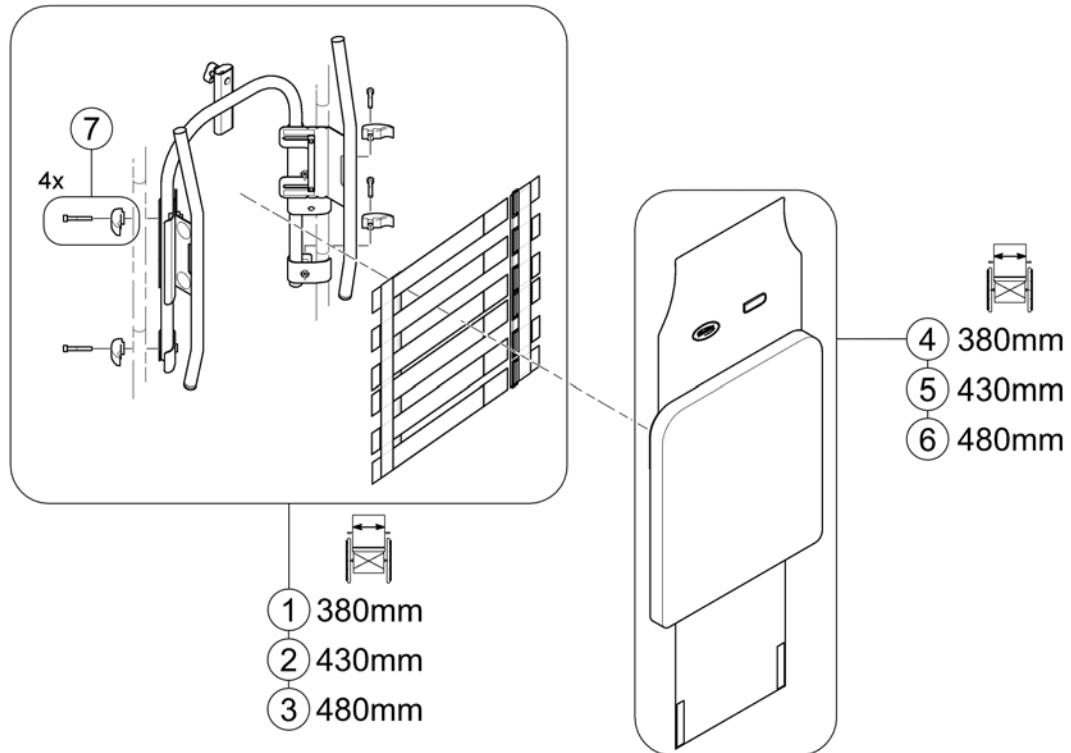
1




1	SP1579675	RAMME STILLBAR RYGG SB39 AZALEA	1
2	SP1579692	RAMME STILLBAR RYGG SB44 AZALEA	1
3	SP1579693	RAMME STILLBAR RYGG SB49 AZALEA	1
4	SP1543667	Ryggtrekk Lateral sb44 D artex sort	1
5	SP1543669	Ryggtrekk Lateral sb49 D artex sort	1
6	SP1543671	Ryggtrekk LATERAL sb55 D artex sort	1
7	SP1580557	Festesett feste av rygg mot ryggamme (sett à 4)	1

Invacare® TDX-SP2 Lowrider, TDX-SP2 NB | Rygg
(Backrest for Ultra Low Maxx Seat)

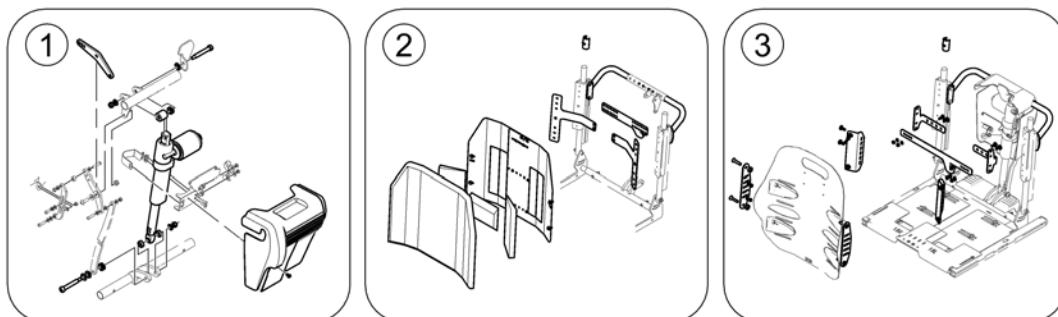156815.iss



-	(Backrest Recliner)	1
-	(High Back)	1
-	(Matrix Back)	1


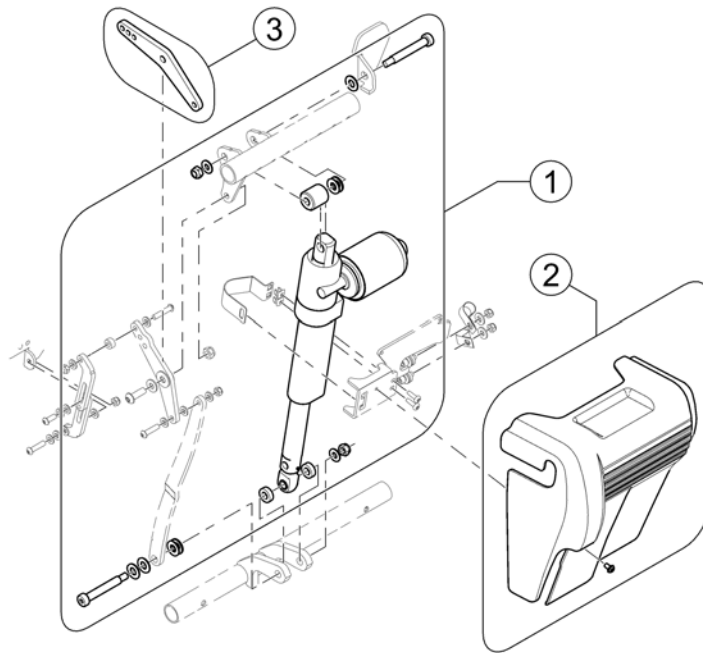


1	SPTRC2191	(Backrest Recliner with Hardware)	1
2	SPTRC2138	(Backrest shroud with mounting hardware)	1
3	SPTRC1926	(Backrest Precline Crankarm 30°)	1
-	SPTRE0322	(Actuator Cable, Backrest Recliner)	1


(High Back)

1588155.lin

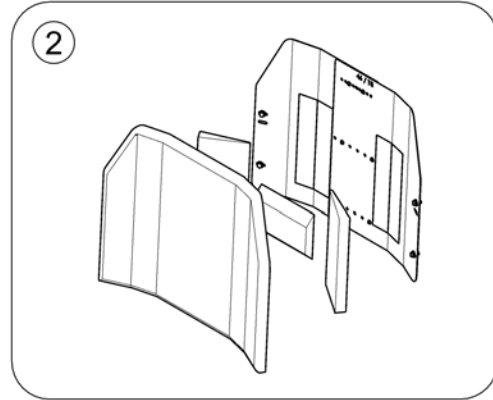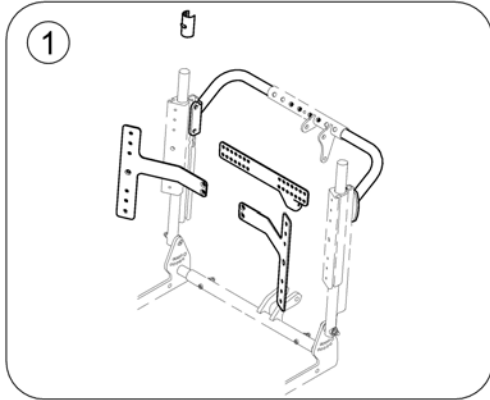High Back




1	-	(Backrest Hardware)	1
2	-	(Backrest Cushions)	1

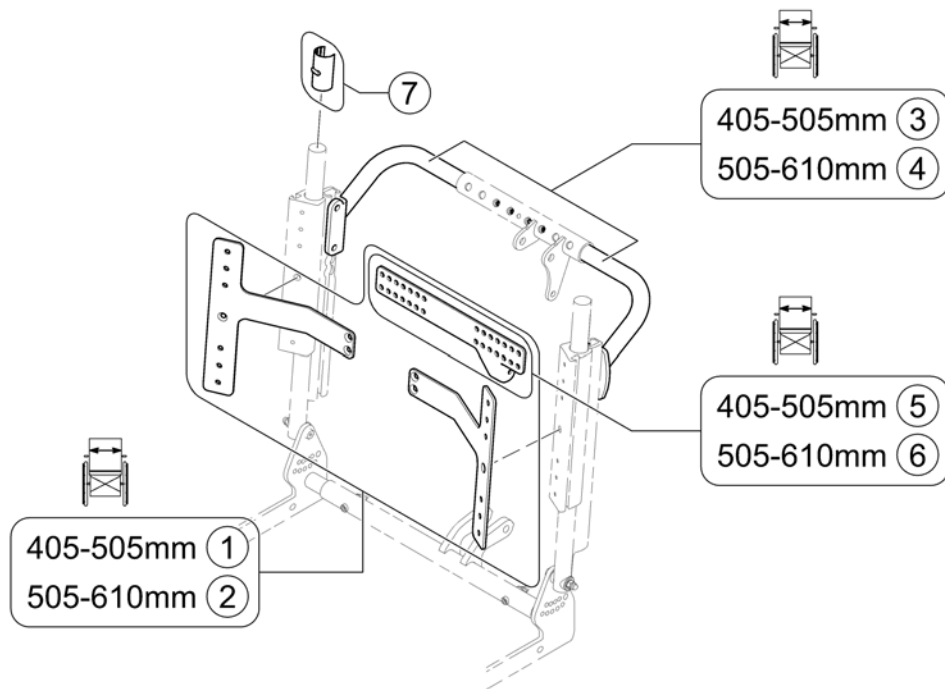

1	SPTRC2136	(High Back Interface Kit SW 405-505mm (16-20"))	1
2	SPTRC2137	(High Back Interface Kit SW 505-610mm (19-24"))	1
3	SPTRC2188	(Push Bars SW 405-505mm (16-20"))	1
4	SPTRC2190	(Push Bars SW 505-610mm (19-24"))	1
5	SPTRX8497	(High Back Center Plate SW 405-505mm (16-20"))	1
6	SPTRX8655	(High Back Center Plate SW 505-610mm (19-24"))	1
7	SPTRX7883	(Linear Bearing)	4

(Backrest Cushions)

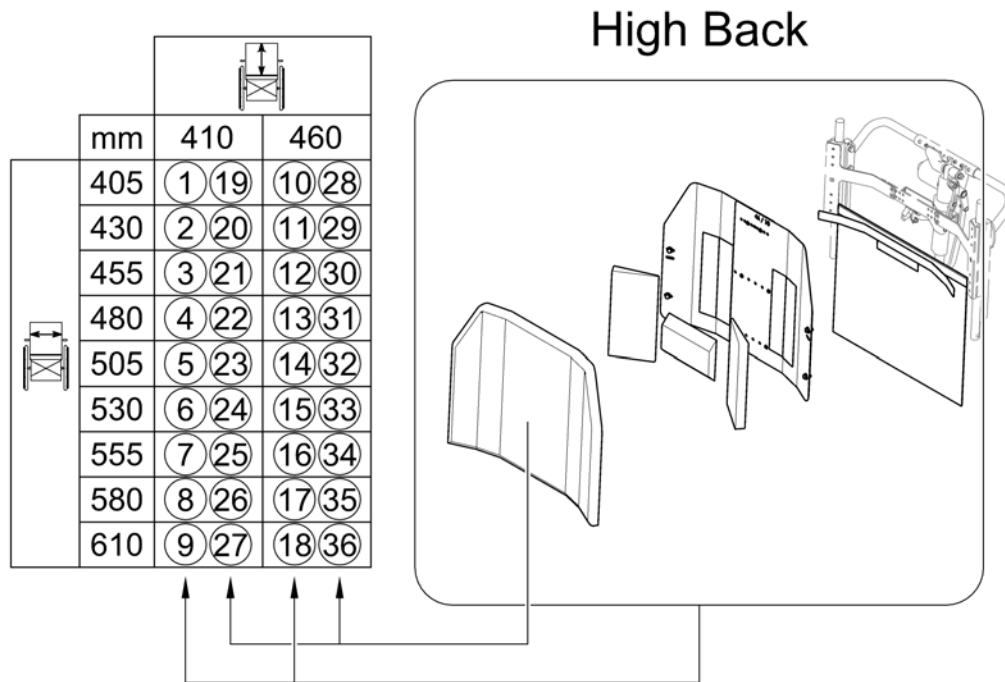158201 no



1	SP1583769	(High Back 405x410mm (16x16"), cpl.)	1
2	SP1583772	(High Back 430x410mm (17x16"), cpl.)	1
3	SP1583774	(High Back 405x410mm (18x16"), cpl.)	1
4	SP1583857	(High Back 480x410mm (19x16"), cpl.)	1
5	SP1583930	(High Back 505x410mm (20x16"), cpl.)	1
6	SP1583932	(High Back 530x410mm (21x16"), cpl.)	1
7	SP1584053	(High Back 555x410mm (22x16"), cpl.)	1
8	SP1584055	(High Back 580x410mm (23x16"), cpl.)	1
9	SP1584081	(High Back 610x410mm (24x16"), cpl.)	1
10	SP1583771	(High Back 405x460mm (16x18"), cpl.)	1
11	SP1583773	(High Back 430x460mm (17x18"), cpl.)	1
12	SP1583775	(High Back 455x460mm (18x18"), cpl.)	1
13	SP1583858	(High Back 480x460mm (19x18"), cpl.)	1
14	SP1583931	(High Back 505x460mm (20x18"), cpl.)	1
15	SP1583934	(High Back 530x460mm (21x18"), cpl.)	1
16	SP1584054	(High Back 555x460mm (22x18"), cpl.)	1
17	SP1584056	(High Back 580x460mm (23x18"), cpl.)	1
18	SP1584118	(High Back 610x460mm (24x18"), cpl.)	1
19	SPTRS0425-1616	(High Back Cover and foam 405x410mm (16x16"))	1
20	SPTRS0425-1716	(High Back Cover and foam 430x410mm (17x16"))	1
21	SPTRS0425-1816	(High Back Cover and foam 455x410mm (18x16"))	1
22	SPTRS0425-1916	(High Back Cover and foam 480x410mm (19x16"))	1
23	SPTRS0425-2016	(High Back Cover and foam 505x410mm (20x16"))	1
24	SPTRS0425-2116	(High Back Cover and foam 530x410mm (21x16"))	1
25	SPTRS0425-2216	(High Back Cover and foam 555x410mm (22x16"))	1
26	SPTRS0425-2316	(High Back Cover and foam 580x410mm (23x16"))	1
27	SPTRS0425-2416	(High Back Cover and foam 610x410mm (24x16"))	1
28	SPTRS0425-1618	(High Back Cover and foam 405x460mm (16x18"))	1
29	SPTRS0425-1718	(High Back Cover and foam 430x460mm (17x18"))	1
30	SPTRS0425-1818	(High Back Cover and foam 455x460mm (18x18"))	1
31	SPTRS0425-1918	(High Back Cover and foam 480x460mm (19x18"))	1


(Backrest Cushions)

32	SPTRS0425-2018	(High Back Cover and foam 505x460mm (20x18"))	1
33	SPTRS0425-2118	(High Back Cover and foam 430x460mm (21x18"))	1
34	SPTRS0425-2218	(High Back Cover and foam 455x460mm (22x18"))	1
35	SPTRS0425-2318	(High Back Cover and foam 480x460mm (23x18"))	1
36	SPTRS0425-2418	(High Back Cover and foam 610x460mm (24x18"))	1


(Matrx Back)

1588242 Rev

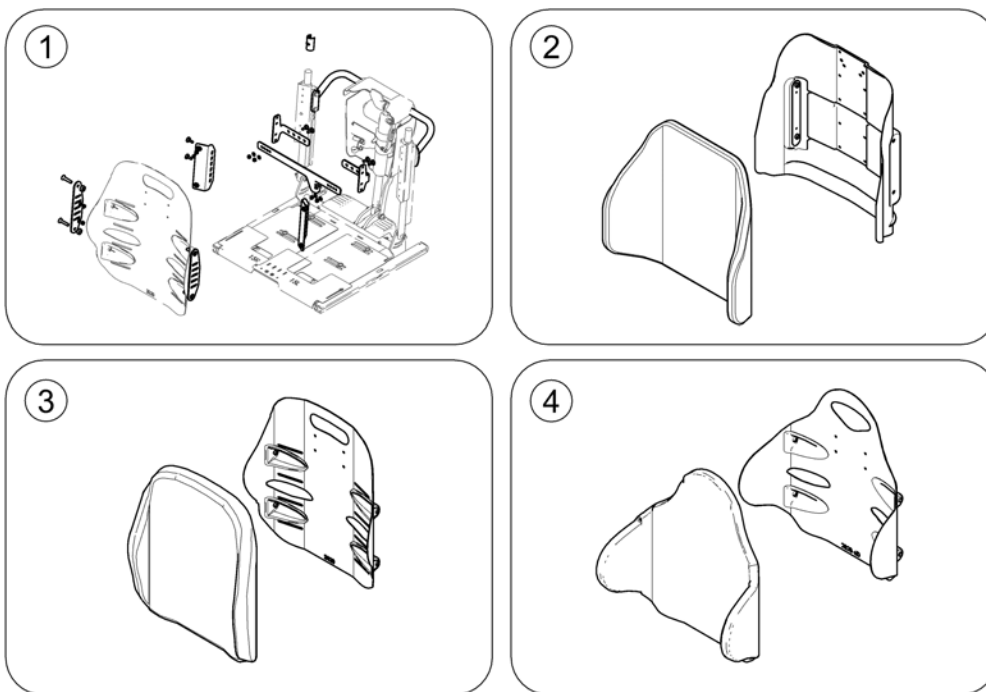Matrx Back



1	-	(Backrest Hardware)	1
2	-	(Posture Back)	1
3	-	(Elite T/R Back)	1
4	-	(Elite Deep Back)	1
-	-	(Privacy Flaps)	1

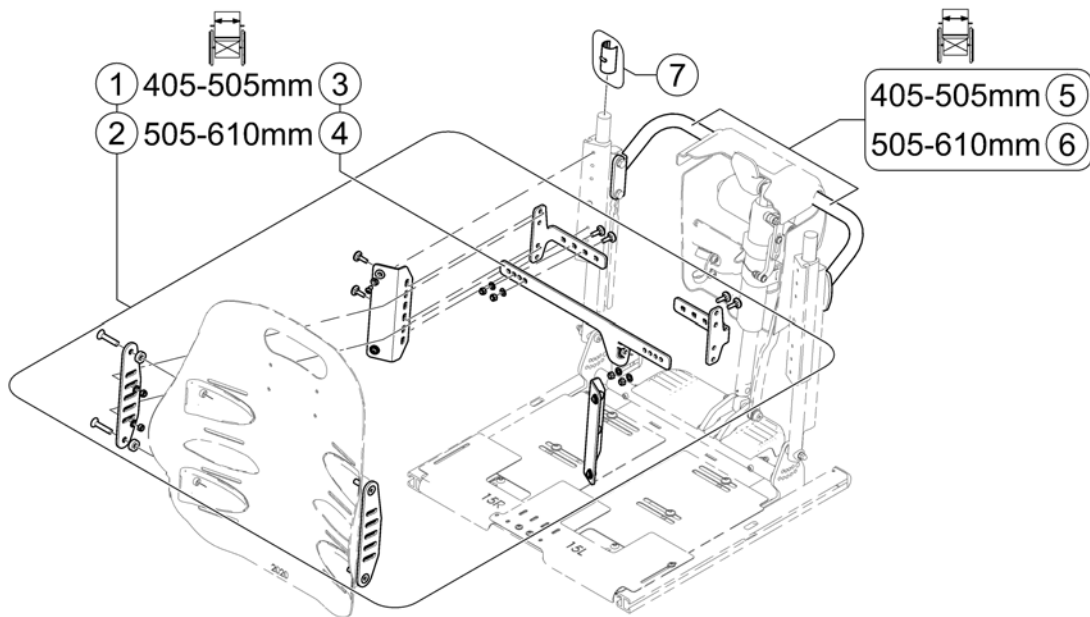


1	SPTRC1951	(Matrx Back Interface Kit SW 405-505mm (16-20"))	1
2	SPTRC2095	(Matrx Back Interface Kit SW 505-610mm (19-24"))	1
3	SPTRX8299	(Matrx Center Plate SW 405-505mm (16-20"))	1
4	SPTRX8779	(Matrx Center Plate SW 505-610mm (19-24"))	1
5	SPTRC2188	(Push Bars SW 405-505mm (16-20"))	1
6	SPTRC2190	(Push Bars SW 505-610mm (19-24"))	1
7	SPTRX7883	(Linear Bearing)	1


(Posture Back)

158204 no

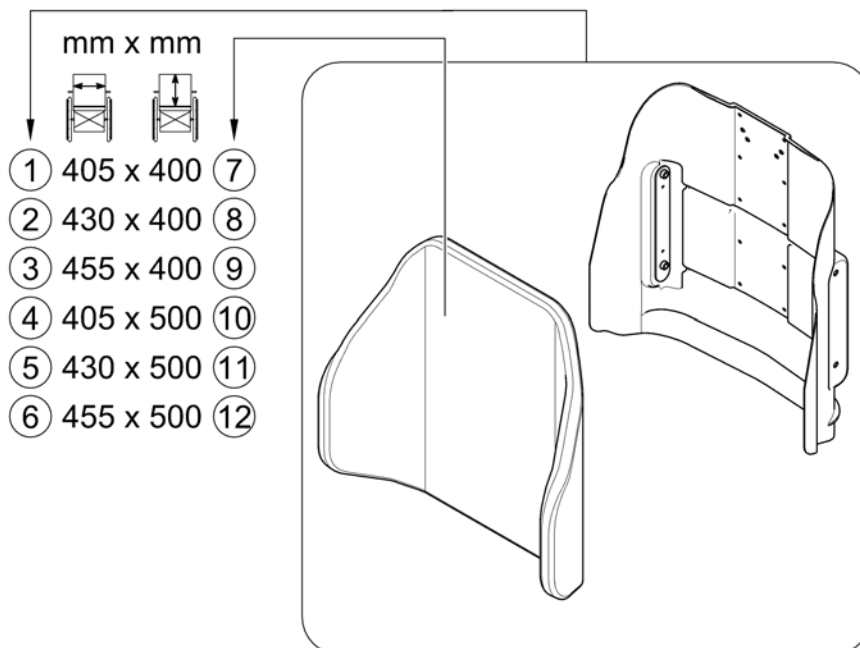Posture Back



1	SP1515727	(Posture Back 405x400mm (16x16"), cpl.)	1
2	SP1515729	(Posture Back 430x400mm (17x16"), cpl.)	1
3	SP1515731	(Posture Back 455x400mm (18x16"), cpl.)	1
4	SP1515728	Elite PB Back 405x500mm MATRX	1
5	SP1515730	(Posture Back 430x500mm (17x20"), cpl.)	1
6	SP1515732	(Posture Back 455x500mm (18x20"), cpl.)	1
7	1544896	Ryggtrekk Posture Back s b40,5 rh40,5	1
8	1544899	Ryggtrekk Posture Back t rekk sb43 rh40,5	1
9	1544995	Ryggtrekk Posture Back s b45,5 rh40,5	1
10	1544897	Ryggtrekk Posture Back s b40,5 rh50,5	1
11	1544982	Ryggtrekk Posture Back t rekk sb43 rh50,5	1
12	1544996	Ryggtrekk Posture Back s b45,5 rh50,5	1

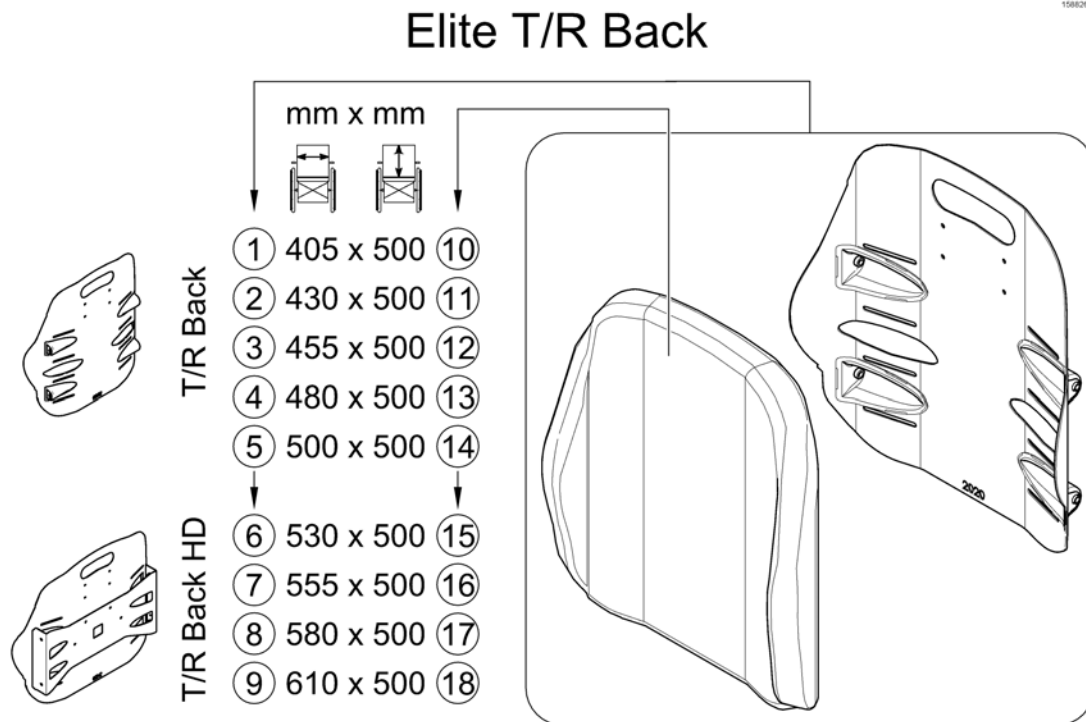


1	SP1580029	(Elite T/R Back 405x500mm (16x20"), cpl.)	1
2	SP1580030	(Elite T/R Back 430x500mm (17x20"), cpl.)	1
3	SP1580031	(Elite T/R Back 455x500mm (18x20"), cpl.)	1
4	SP1580032	(Elite T/R Back 480x500mm (19x20"), cpl.)	1
5	SP1580033	ELITE TR RYGG 505x500mm 20x20"	1
6	SP1584160	(Elite T/R Back HD 530x500mm (21x20"), cpl.)	1
7	SP1584178	(Elite T/R Back HD 555x500mm (22x20"), cpl.)	1
8	SP1584180	(Elite T/R Back HD 580x500mm (23x20"), cpl.)	1
9	SP1584182	(Elite T/R Back HD 610x500mm (24x20"), cpl.)	1
10	1580110	Ryggtrekk PB Elite TiltR ecline sb48 rh46	1
11	1580297	Ryggtrekk PB Elite TiltR ecline sb48 rh51	1
12	1580298	Ryggtrekk PB Elite TiltR ecline sb51 rh41	1
13	1580299	Ryggtrekk PB Elite TiltR ecline sb51 rh46	1
14	1580300	Ryggtrekk PB Elite TiltR ecline sb51 rh51	1
15	SPTRS0375-2120	(Backrest Cushion - Elite T/R Back HD 530x500mm (21x20"))	1
16	SPTRS0375-2220	(Backrest Cushion - Elite T/R Back HD 555x500mm (22x20"))	1
17	SPTRS0375-2320	(Backrest Cushion - Elite T/R Back HD 580x500mm (23x20"))	1
18	SPTRS0375-2420	(Backrest Cushion - Elite T/R Back HD 610x500mm (24x20"))	1


(Elite Deep Back)

158208 no

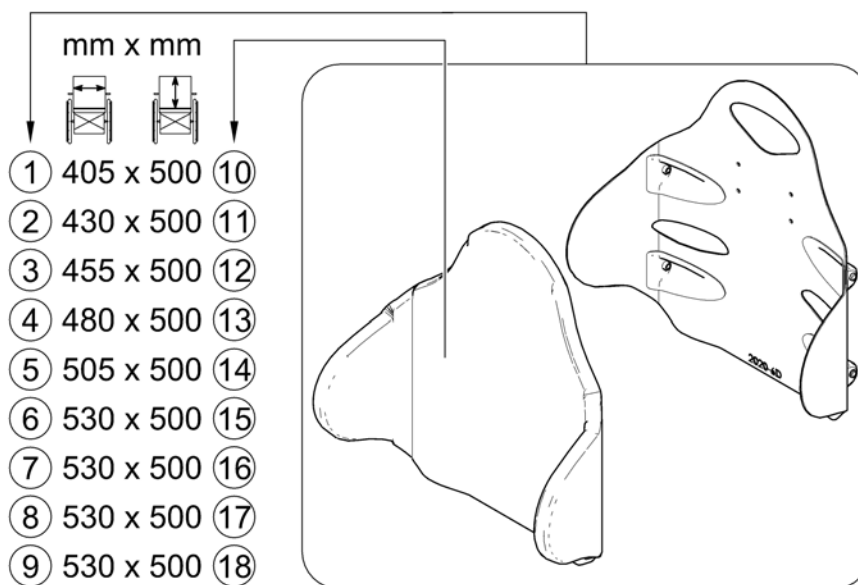Elite Deep Back



1	SP1543515	(Elite Deep Back 405x500mm (16x20"), cpl.)	1
2	SP1543519	(Elite Deep Back 430x500mm (17x20"), cpl.)	1
3	SP1543523	(Elite Deep Back 455x500mm (18x20"), cpl.)	1
4	SP1543987	(Elite Deep Back 480x500mm (19x20"), cpl.)	1
5	SP1543991	(Elite Deep Back 505x500mm (20x20"), cpl.)	1
6	SP1584177	(Elite Deep Back 530x500mm (21x20"), cpl.)	1
7	SP1584179	(Elite Deep Back 555x500mm (22x20"), cpl.)	1
8	SP1584181	(Elite Deep Back 580x500mm (23x20"), cpl.)	1
9	SP1583933	(Elite Deep Back 610x500mm (24x20"), cpl.)	1
10	1544073	Ryggtrekk Elite dyp Back sb40,5 rh50,5	1
11	1544077	Ryggtrekk Elite dyp Back sb43 rh50,5	1
12	1544081	Ryggtrekk Elite dyp Back sb45,5 rh50,5	1
13	1544085	Ryggtrekk Elite dyp Back sb48 rh50,5	1
14	1544124	Ryggtrekk Elite dyp Back sb50,5 rh50,5	1
15	SPTRS0344-2120	(Backrest Cushion - Elite Deep Back 530x500mm (21x20"))	1
16	SPTRS0344-2220	(Backrest Cushion - Elite Deep Back 555x500mm (22x20"))	1
17	SPTRS0344-2320	(Backrest Cushion - Elite Deep Back 580x500mm (23x20"))	1
18	SPTRS0344-2420	(Backrest Cushion - Elite Deep Back 610x500mm (24x20"))	1

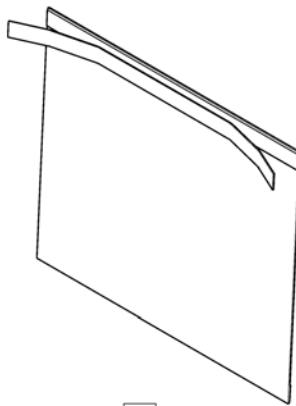


405/430mm	①
455/480mm	②
505/530mm	③

1	SP1551800	(Privacy Flap SW 405/430mm)	1
2	SP1551903	(Privacy Flap SW 455/480mm)	1
3	SP1551904	(Privacy Flap SW 505/530mm)	1



(Backrest for Max Seat)

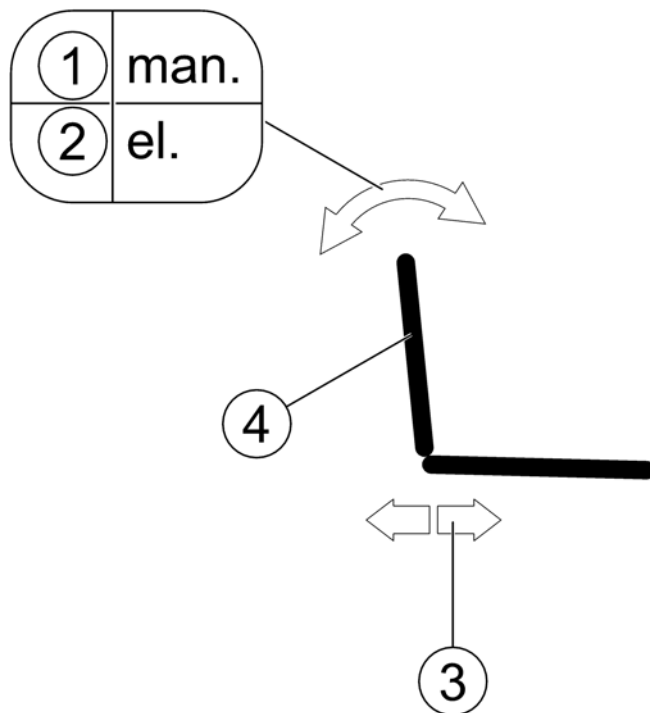1582669.iso




1	-	(Manual backrest adjustment)	1
2	-	(Powered backrest adjustment)	1
3	-	(Seat depth adapter)	1
4	-	(Backrest)	1

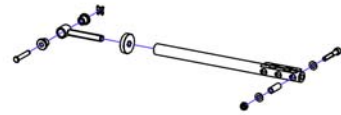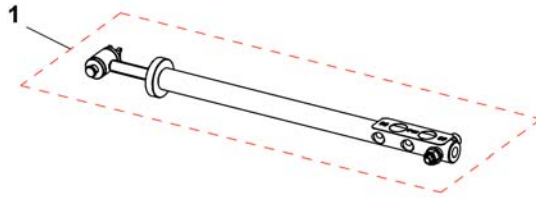1545185.jpg



1 SP1545185 (Spindle, backrest adjustment, cpl.)

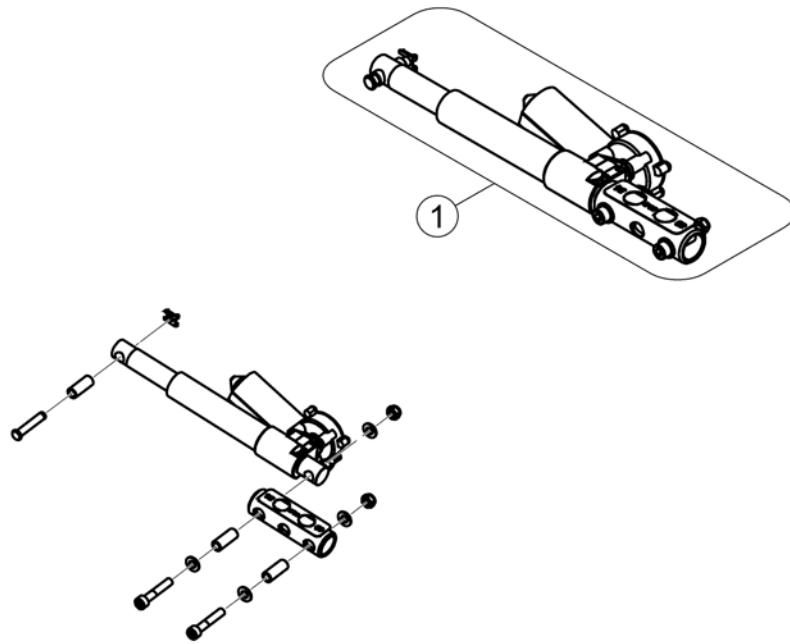
1



1 SP1544669 (Actuator, cpl.)

1

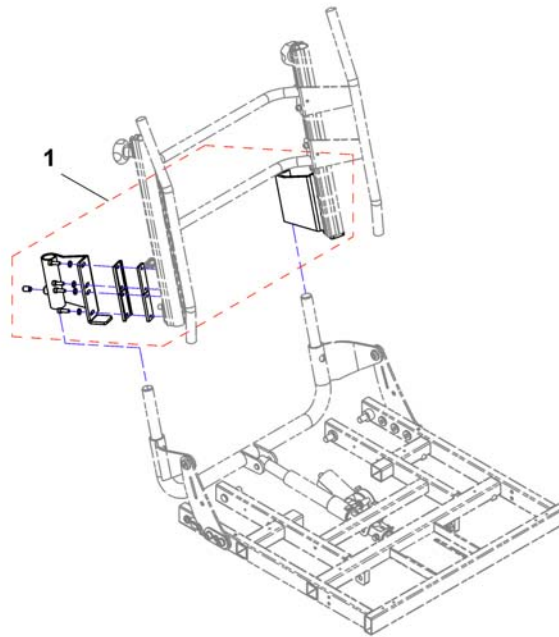1551724K.dwg




1 SP1551724


(Seat depth adapter, cpl.)

1

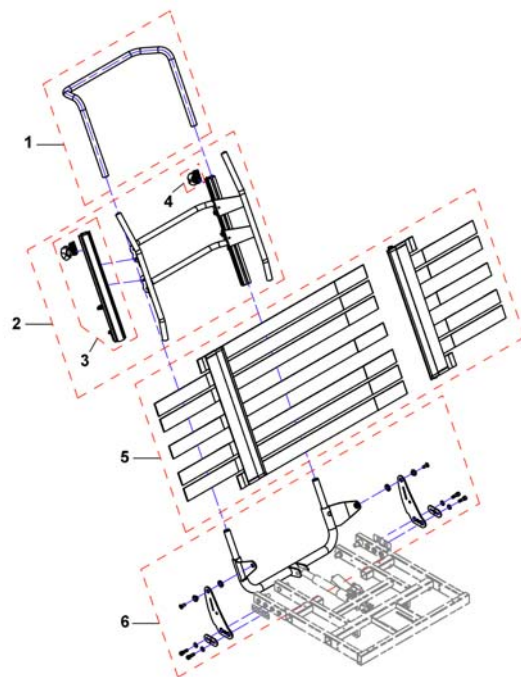(Backrest)



1	SP1543082	(Support tube, SW49)	1
1	SP1543081	(Support tube, SW55)	1
1	SP1543080	(Support tube, SW61)	1
2	SP1551199	(Backrest frame, cpl., SW49)	1
2	SP1526824	(Backrest frame, cpl., SW55)	1
2	SP1526823	(Backrest frame, cpl., SW61)	1
3	1445848	RYGGRØR TIL VELCRO KPL FLEX 3	1
4	5253088	LÅSKLOSS KPL F RYGGRØR REA COMFORT	1
-	1508867	(Cover SW49, dartex grey) <i>not shown</i>	1
-	SP1526459	(Backrest upholstery SW 55 cm) <i>not shown</i>	1
-	1526424	(Cover, Lateral 700X500, Dartex Grey) <i>not shown</i>	1
5	1445694	RYGGTREKK BORRELÅS STILL BAR RYGG BELLIS AZALEA	1
5	SP1526308	(Backrest, tension adjustable, SW 55 - 61cm)	1
6	SP1544310	(Backrest bracket, cpl., SW49)	1
6	SP1544311	(Backrest bracket, cpl., SW55)	1
6	SP1544312	(Backrest bracket, cpl., SW61)	1

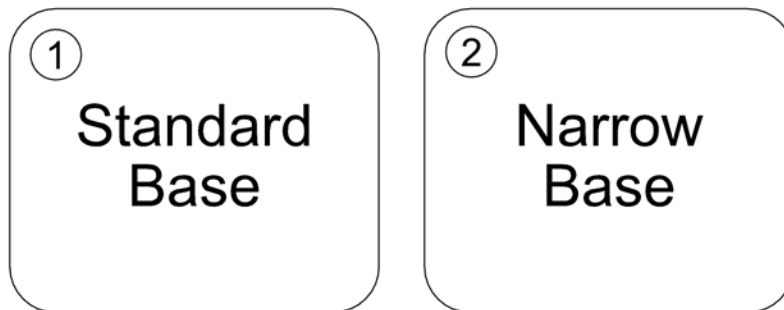

1	-	(Legrests - Standard Base)	0
2	-	(Legrests - Narrow Base)	0

(Legrests - Standard Base)

SPCTDX_SP220.kso



1	-	(Legrests for Modulite Seat)	0
2	-	(Legrests for Ultra Low Maxx Seat)	0
3	-	(Legrests for Max Seat)	0

(Legrests for Modulite Seat)


-	(Legrests 80° & man. adjustable)	0
-	(Legrest Support, Vari series - Modulite Seat)	0
-	(Hanger Vari A (04/2009 =>))	0
-	(Hanger Vari F)	0
-	(Hangers ADM LH/RH manual)	0
-	(Hanger Vari ADE/ADP electric (07/2005 =>))	0
-	(Footrests Vari, Flip-Up (01.08.10 =>))	0
-	(Footrests Vari, Flip-Up & angle adjustable (01.08.10 =>))	0
-	(Footrests Vari, angle & depth adjustable)	0
-	(Footboard Vari, pre-assembled)	0
-	<i>Football Chair only!</i>	0
-	(Central Legrests - Modulite Seat)	0

(Legrests 80° & man. adjustable)

-	(Legrest Support, 80° & man. elevating - Standard & Modulite Seat)	1
-	(Hanger, Std. 80° (01.07.15 =>))	1
-	(Hanger, Std. 80° (=> 30.06.15))	1
-	(Hanger, manually height adjustable)	1
-	(Footplate, fixed 80°)	1
-	(Footplate, depth-, angleadjustable, 80°)	1
-	(Amputee Support)	1

(Legrest Support, 80° & man. elevating - Standard & Modulite Seat)

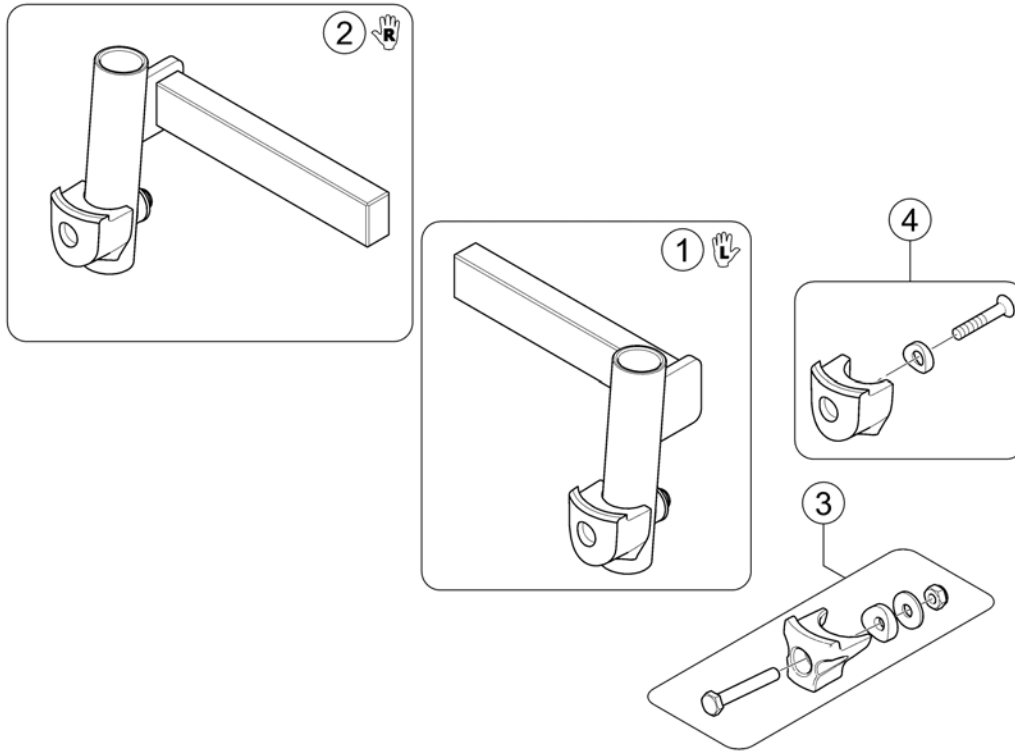
1441526PAC.00




1	SP1531195	ADAPTER BENST KPL VE BORA <i>We recommend to exchange the hangers by pairs for products built before 01.01.12, because of new design.</i>	01.07.2015 →	1
2	SP1531196	ADAPTER BENST KPL HØ BORA <i>We recommend to exchange the hangers by pairs for products built before 01.01.12, because of new design.</i>	01.07.2015 →	1
3	SP1589136	(Kit, Support Nose)	→ 30.06.2015	1
4	SP1588726	FESTE BENSTØTTE MANUELL BORA	01.07.2015 →	1

(Hanger, Std. 80° (01.07.15 =>))

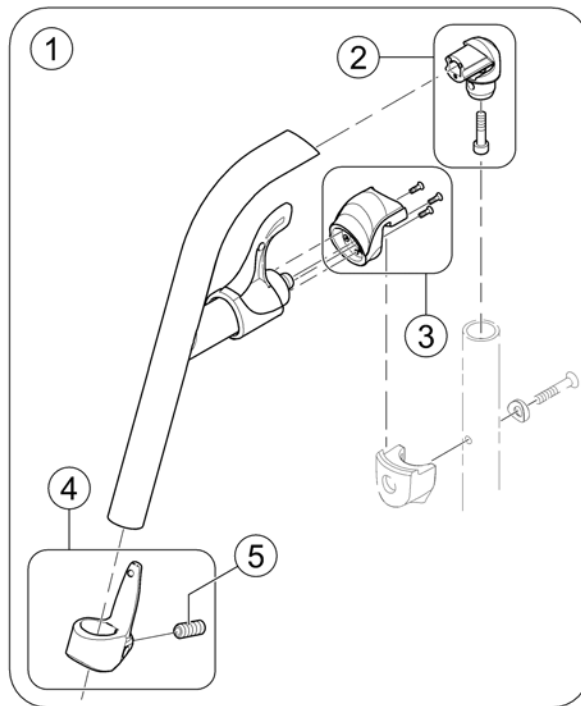1588852.000



1	SP1574349-9005	(Hanger Standard 80° cpl. - black)	1
2	SP1588727	(Kit, Pivot Std. 80° Hanger)	1
3	SP1588731	(Nose locking Kit)	1
4	SP1588730	(Calf Strap Holder, cpl.)	1
5	1424571	FJÆRENDE TRYKSKRUE M10x23	1

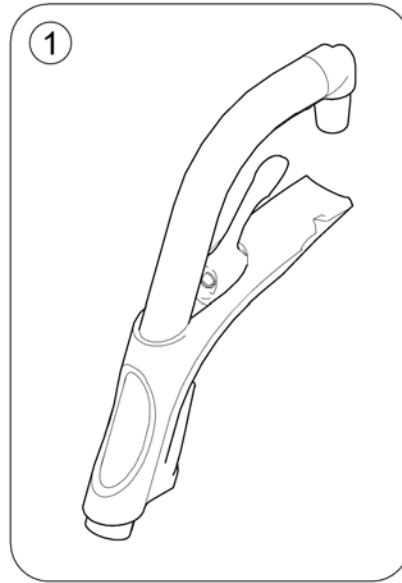


1 SP1493055

BENSTØTTEOVERDEL FAST BORA

1


(Hanger, manually height adjustable)

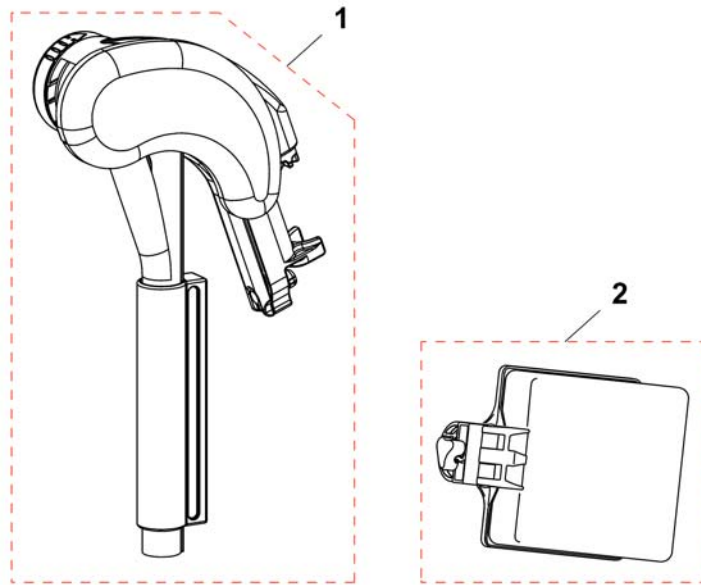1513006.dwg



1	1513006	BENST OVERDEL VREG HØ M/ KNEPUTE ALU 10 2008	1
1	1513007	BENST OVERDEL VREG VE M/ KNEPUTE ALU 10 2008	1
2	1516314	LEGGSTØTTE sb44 ,5X17X8 TR26 KPL	1
2	1518532	LEGGSTØTTE KPL 19" FOCUS,SPIREA	1

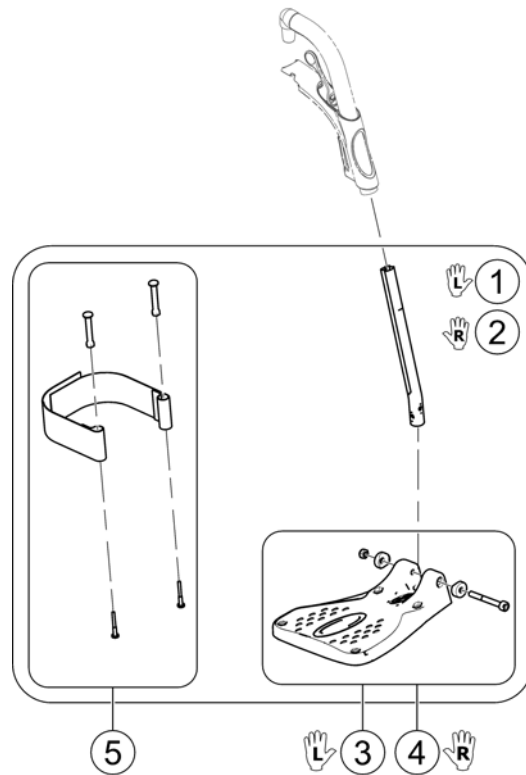
1588087.joo




1	SP1537943	(Footplate LH fixed angle with Support Tube, cpl.38-40,5cm)	1
2	SP1537941	(Footplate RH fixed angle with Support Tube, cpl.38-40,5cm)	1
3	SP1529643	(Footplate LH)	1
4	SP1529642	(Footplate RH)	1
5	SP1442467	LEGGBÅND HELKAPPE PAR MIRAGE, STORM4, TDX	1

(Footplate, depth-, angleadjustable, 80°)

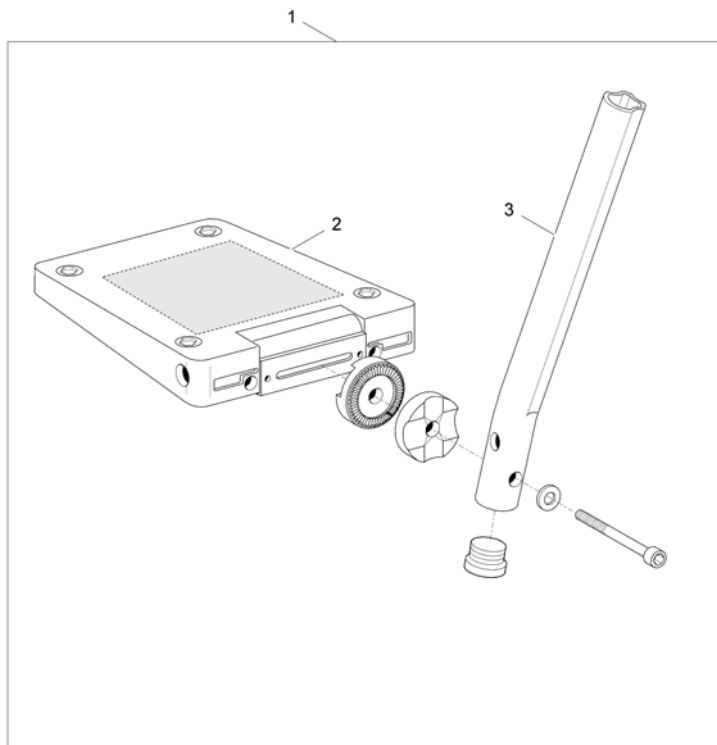1428654X.iso




1	1428660	FOTPLATE M RØR NEDRE VE V.REG ACTION 3 SB45,5	1
1	1428654	FOTPLATE M RØR NEDRE HØ V.REG ACTION 3 SB45,5	1
2	5106594	FOTPLATE 48 BELLIS AZAL/FOC 166	1
3	1434561-0010	INNERRØR BENSTØTTE KORTE AZALEA/BELLIS	1

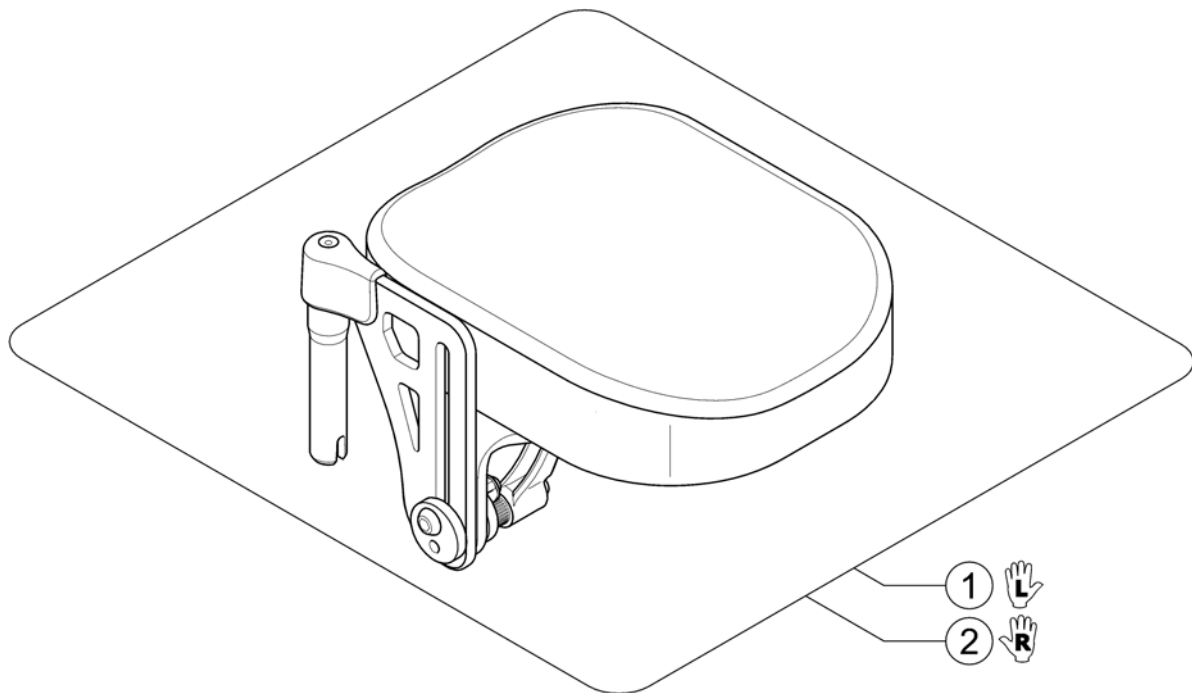


1	1442991	AMP STØTTE VE KPL 11" [0555]	1
2	1442990	AMP STØTTE HØ KPL 11" [0554]	1

(Legrest Support, Vari series - Modulite Seat)

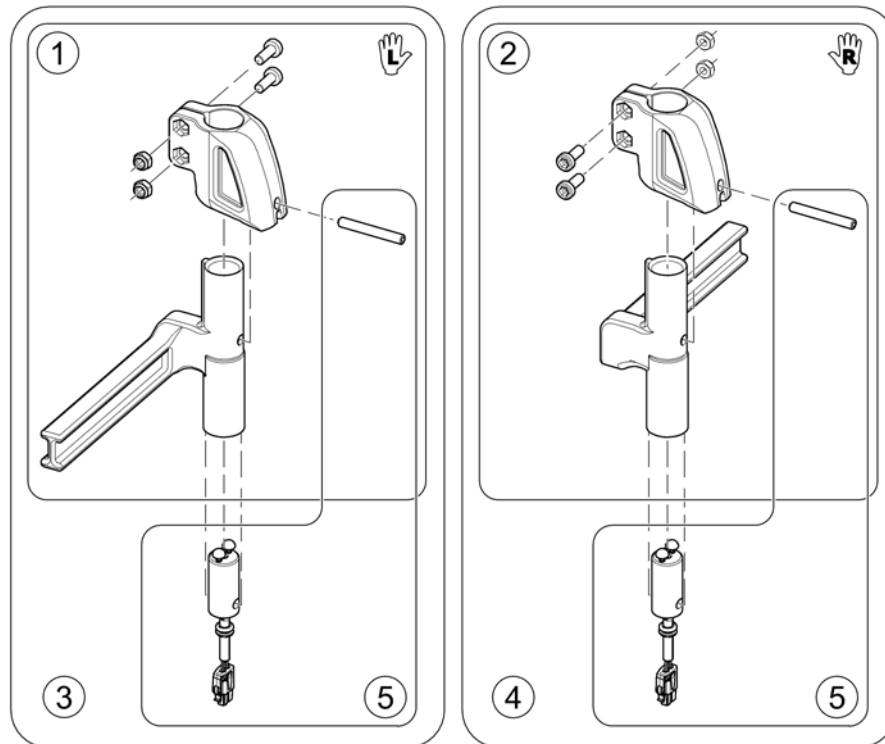1573651.000




1	SP1449733	BENSTØTTEHOLDER MAN VE 5CM STORM	1
2	SP1449732	BENSTØTTEHOLDER MAN HØ 5CM STORM	1
3	SP1529951	BENSTHOLD EL 5CM VE TDX	1
4	SP1529949	BENSTHOLDER 5CM HØ TDX	1
5	SP1527775	KABEL HOLDER BENST 65MM TDX	1

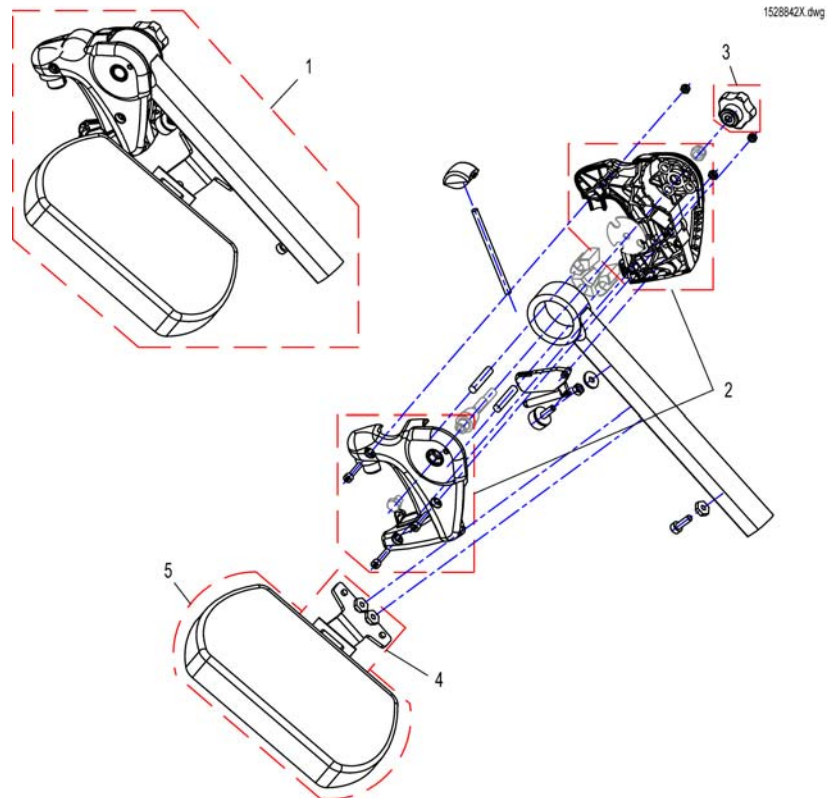


-	SP1527820	(Rework-Kit Vari A Legrest Version 1 => Version 2) <i>Vari A Version 1 => Version 2 Upgrade Kit</i>	1
1	SP1527498	BENSTØTTE OVERDEL VE MAN VREG LKOMP FRA 07/2005	1
1	SP1527492	BENST OVERDEL MAN HØ STORM4	1
2	SP1434955	(Cover Legrest cpl.)	1
3	L34011	(Handknob M8 Tristarfemale)	1
4	SP1450447	LEGGPLATEHOLDER KPL VE NY MAI 05	1
4	SP1450446	LEGGPLATEHOLDER KPL HØ NY MAI 05	1
5	6301333	LEGGPLATE EL OG MAN VREG BENST STORM/TWISTER	1

(Hanger Vari F)

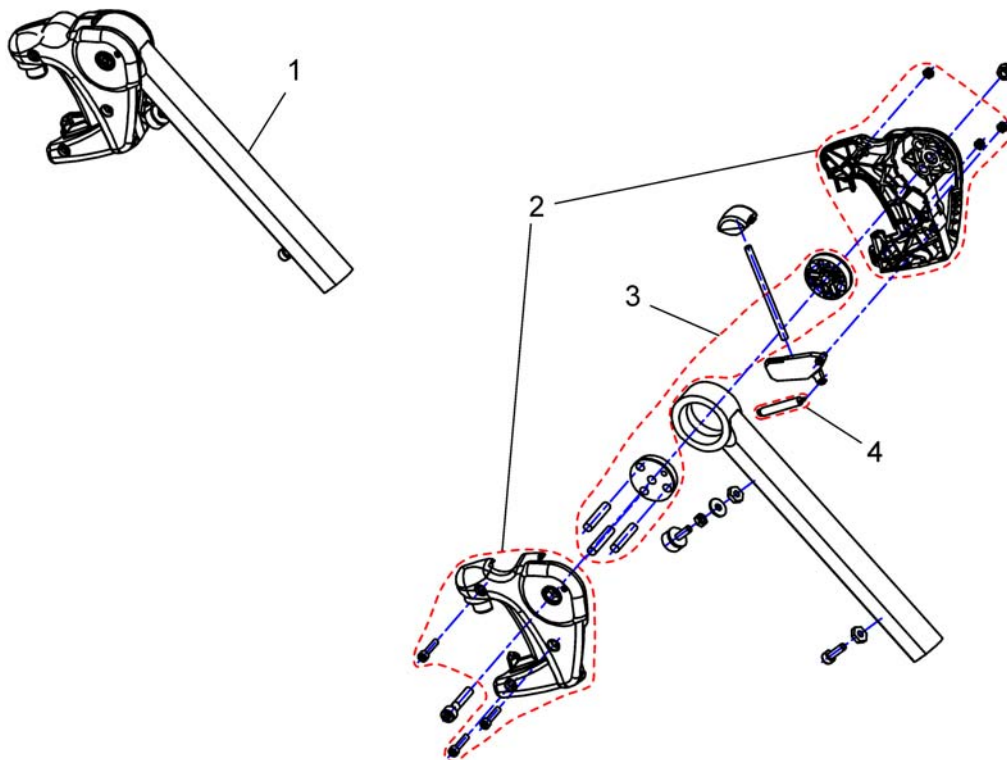144934 / X.dwg




1	SP1449347	BENSTØTTE OVERDEL VARI F STORM TDX BORA VE	1
1	SP1449346	BENSTØTTE OVERDEL VARI F STORM TDX BORA HØ	1
2	SP1434955	(Cover Legrest cpl.)	1
3	1434956	(Mechanic)	1
4	SPZ173178	FJÆR STORM4 BORA	1

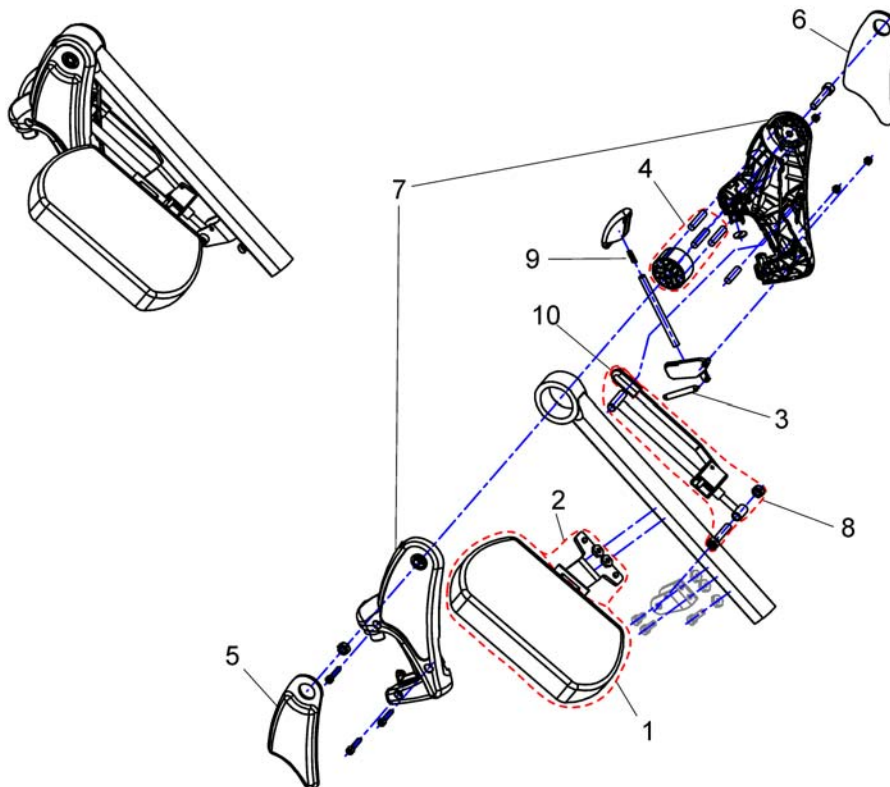

1	6301333	LEGGPLATE EL OG MAN VREG BENST STORM/TWISTER	1
2	SP1450447	LEGGPLATEHOLDER KPL VE NY MAI 05	1
2	SP1450446	LEGGPLATEHOLDER KPL HØ NY MAI 05	1
3	SPZ173178	FJÆR STORM4 BORA	1
4	SPZ174134	(SLIDING BUSHING)	1
5	SP1450410	KNEPOLSTER BENST EL VE 07.2005->	1
6	SP1450448	KNEPOLSTER BENSTG EL HØ 07.2005->	1
7	SP1449293	BENST HOLDER OVERDEL PAR EL VREG 07.2005->	1
8	SPZ174130	BENST ERGO REGULLERINGSMEK KPL	1
9	SPN12196	FJÆR STORM4 BORA	1
10	SP1006490	(Brake Lever Rubber -LT2000)	1
11	SP1449410	BENSTØTTE OVERDEL ADM HØ STORM4 MAN VREG	1
11	SP1449409	BENSTØTTE OVERDEL ADM VE STORM4 MAN VREG	1

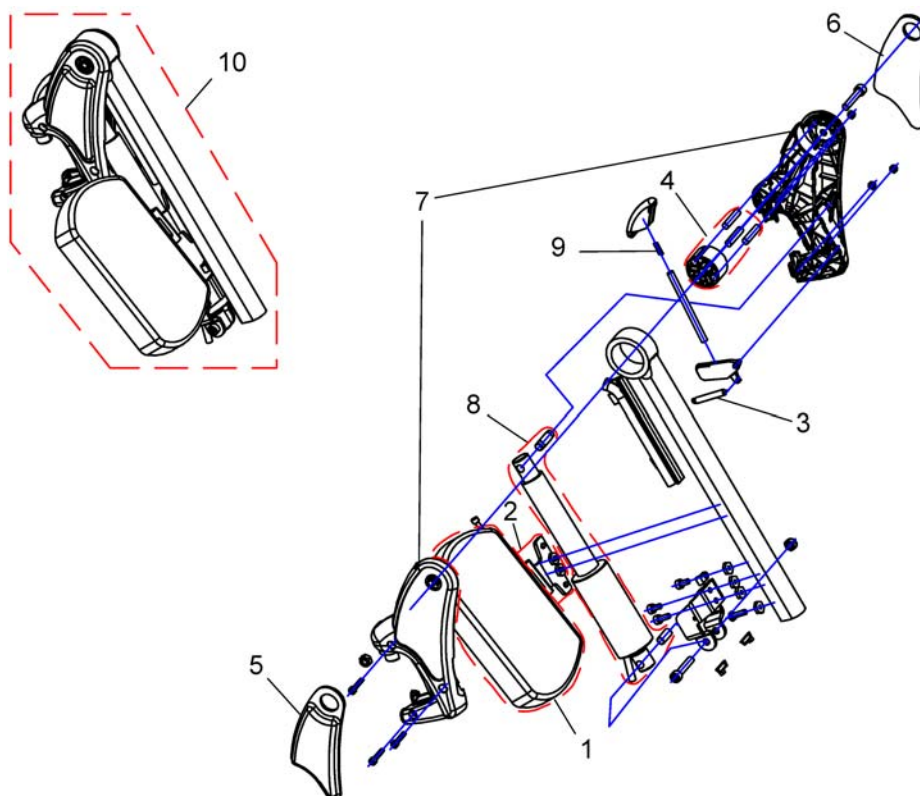(Hanger Vari ADE/ADP electric (07/2005 =>))




1449407X.dwg

1	6301333	LEGGPLATE EL OG MAN VREG BENST STORM/TWISTER	1
2	SP1450447	LEGGPLATEHOLDER KPL VE NY MAI 05	1
2	SP1450446	LEGGPLATEHOLDER KPL HØ NY MAI 05	1
3	SPZ173178	FJÆR STORM4 BORA	1
4	SPZ174134	(SLIDING BUSHING)	1
5	SP1450448	KNEPOLSTER BENSTG EL HØ 07.2005->	1
6	SP1450410	KNEPOLSTER BENST EL VE 07.2005->	1
7	SP1449293	BENST HOLDER OVERDEL PAR EL VREG 07.2005->	1
8	SPZ151006	MOTOR EL BENST STORM4	1
9	SPN12196	FJÆR STORM4 BORA	1
10	SP1449408	BENST OVERDEL ADP HØ STORM4	1
10	SP1449407	BENST OVERDEL ADP VE STORM4	1

Invacare® TDX-SP2 Lowrider, TDX-SP2 NB | Benstøtte
(Footrests Vari, Flip-Up (01.08.10 =>))

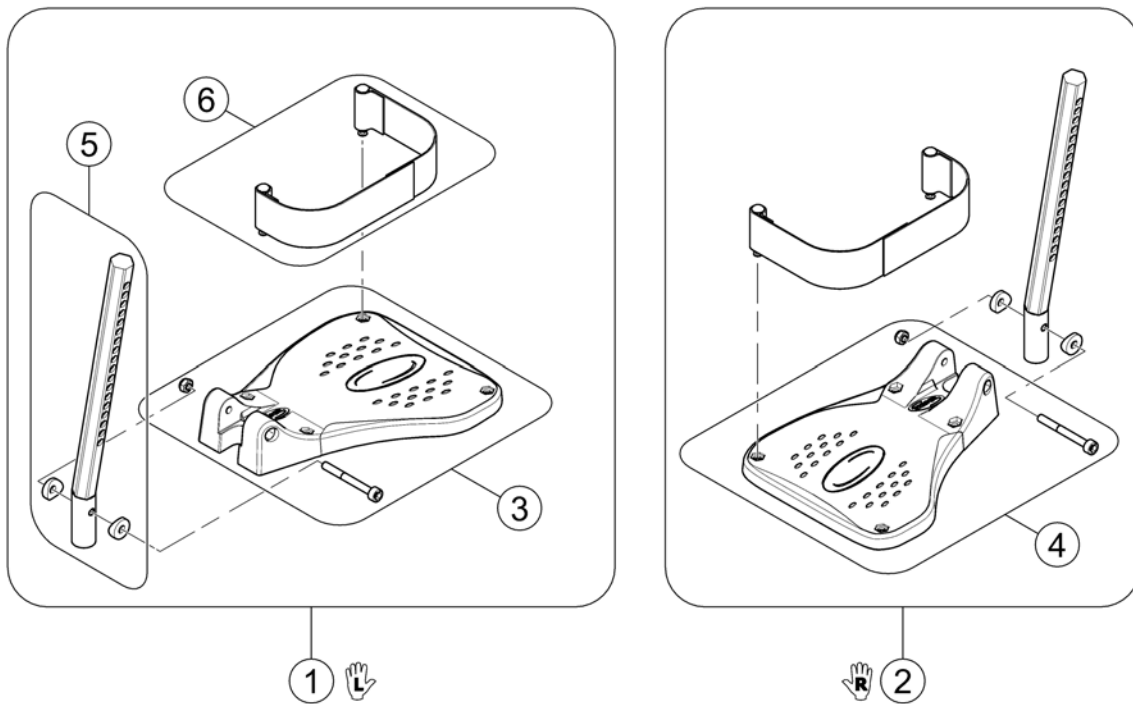
152964X.lin



1	SP1541575	FORPLATE OPPFELLBAR KPL VENSTRE	1
2	SP1541574	FOTPLATE OPPFELLBAR KPL HØYRE	1
3	SP1529645	FOTPLATE KLP SB43-45 VE STORM4	1
4	SP1529644	FOTPLATE KLP SB43-45 HØ STORM4	1
5	SP1543045	(Plug-in tube LH/RH)	1
6	SP1442467	LEGGBÅND HELKAPPE PAR MIRAGE, STORM4, TDX	1

(Footrests Vari, Flip-Up & angle adjustable (01.08.10 =>))

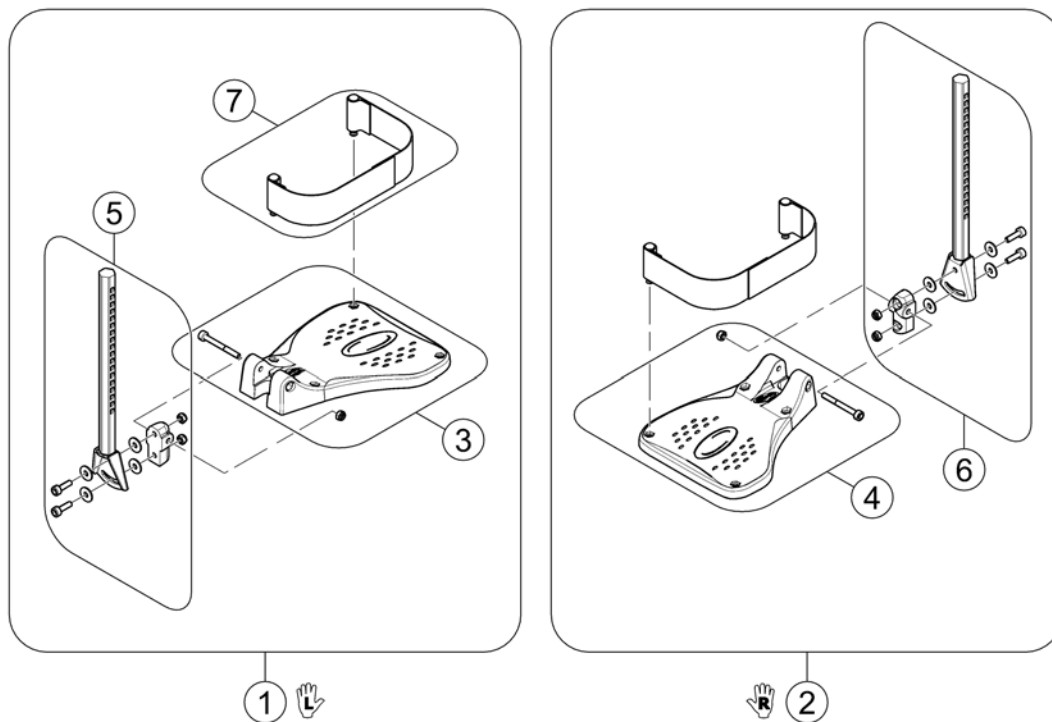
152964EX.no



1	SP1541580	FOTPLATE M/INNERRØR KPL VE VREG	1
2	SP1541579	FOTPLATE M/INNERRØR KPL	1
3	SP1529645	FOTPLATE KLP SB43-45 VE STORM4	1
4	SP1529644	FOTPLATE KLP SB43-45 HØ STORM4	1
5	SP1423204	(FOOT PLATE TUBE, ANGLE ADJUSTMENT LH CPL.)	1
6	SP1423203	(FOOT PLATE TUBE, ANGLE ADJUSTMENT RH CPL.)	1
7	SP1442467	LEGGBÅND HELKAPPE PAR MIRAGE, STORM4, TDX	1

(Footrests Vari, angle & depth adjustable)

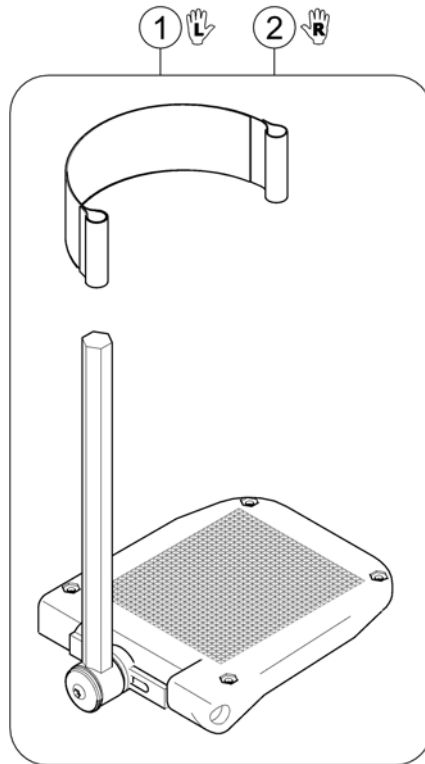
1583898.no



1	SP1419344	FOTPLATE VE VREG OG DREG KPL MED INNERRØR	1
2	SP1419345	FOTPLATE HØ VREG OG DREG KPL MED INNERRØR	1

(Footboard Vari, pre-assembled)

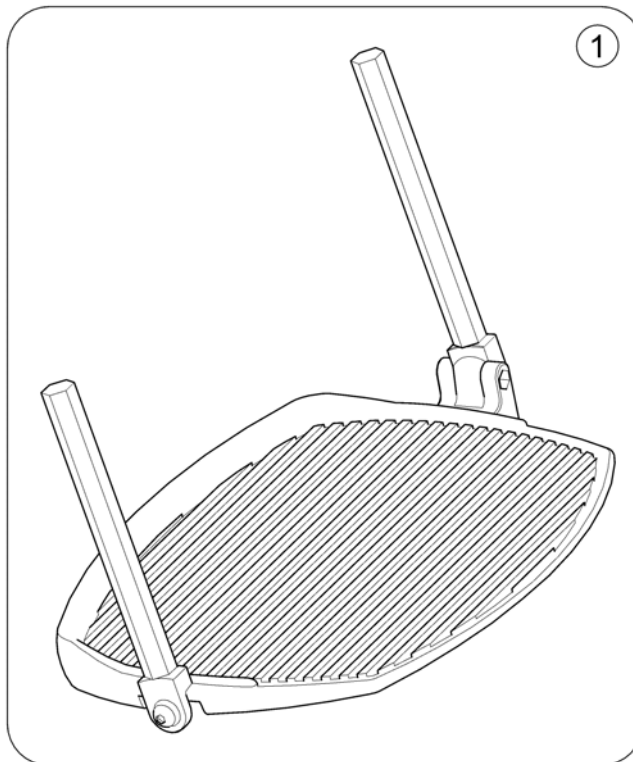1583899-iso



1 SP1423156


FOTPLATE HEL KPL STORM3

1

-	(Support Centr. Legrest man.)	1
-	(Support Centr. Legrest pow. LNX (=> 2016-06-30))	1
-	(Support Centr. Legrest pow. LNX (2016-07-01 =>))	1
-	(Central Legrest man.)	1
-	(Central Legrest pow. LNX (=> 2016-06-30))	1
-	(Central Legrest pow. LNX (2016-07-01 =>))	1

(Support Centr. Legrest man.)

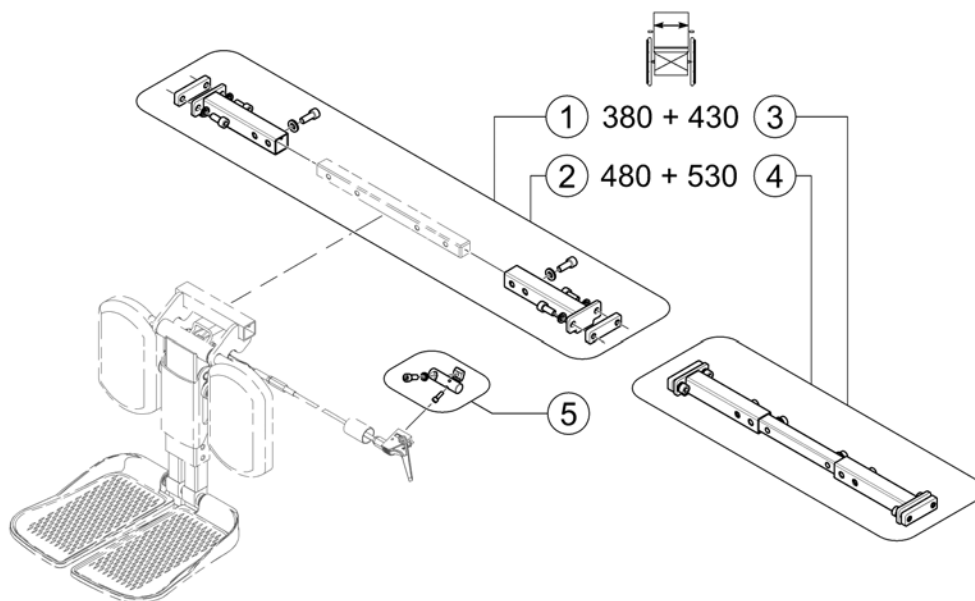
1563891.jpg



1	SP1579989	(Support, Centr. Legrest man. - SW 380+430mm)	1
2	SP1579990	(Support, Centr. Legrest man. LNX - SW 480+530mm)	1
3	SP1579991	(Support, Centr. Legrest man. LNX - SW 380+430mm, cpl.)	1
4	SP1579992	(Support, Centr. Legrest man. LNX - SW 480+530mm, cpl.)	1
5	SP1579988	(Support, Lever, cpl.)	1

Invacare® TDX-SP2 Lowrider, TDX-SP2 NB | Benstøtte
(Support Centr. Legrest pow. LNX (=> 2016-06-30))

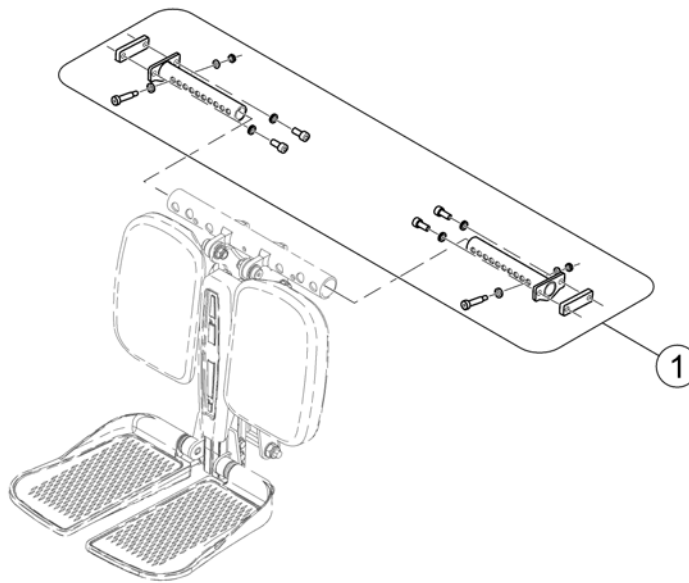1583893.kim



1 SP1575144

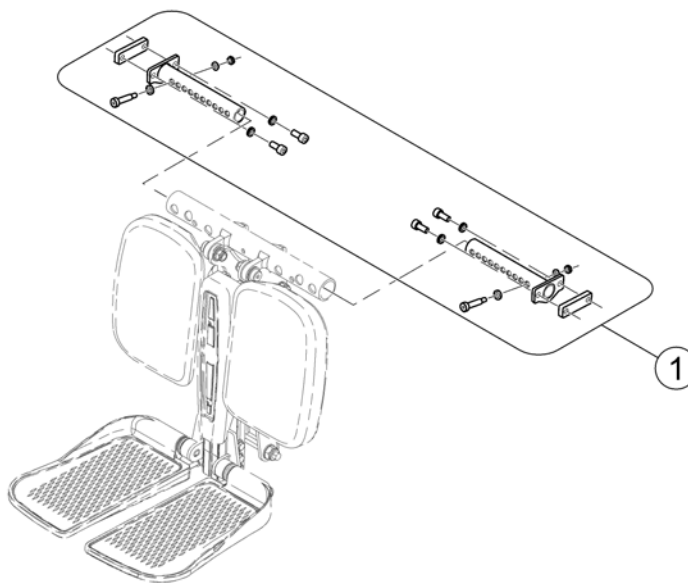
(Support, Centr. Legrest el. LNX, cpl.)

→ 30.06.2016

1

(Support Centr. Legrest pow. LNX (2016-07-01 =>))

1583893.kim




1 SP1602090

(Support, Central Legrest pow. LNX)

01.07.2016 →

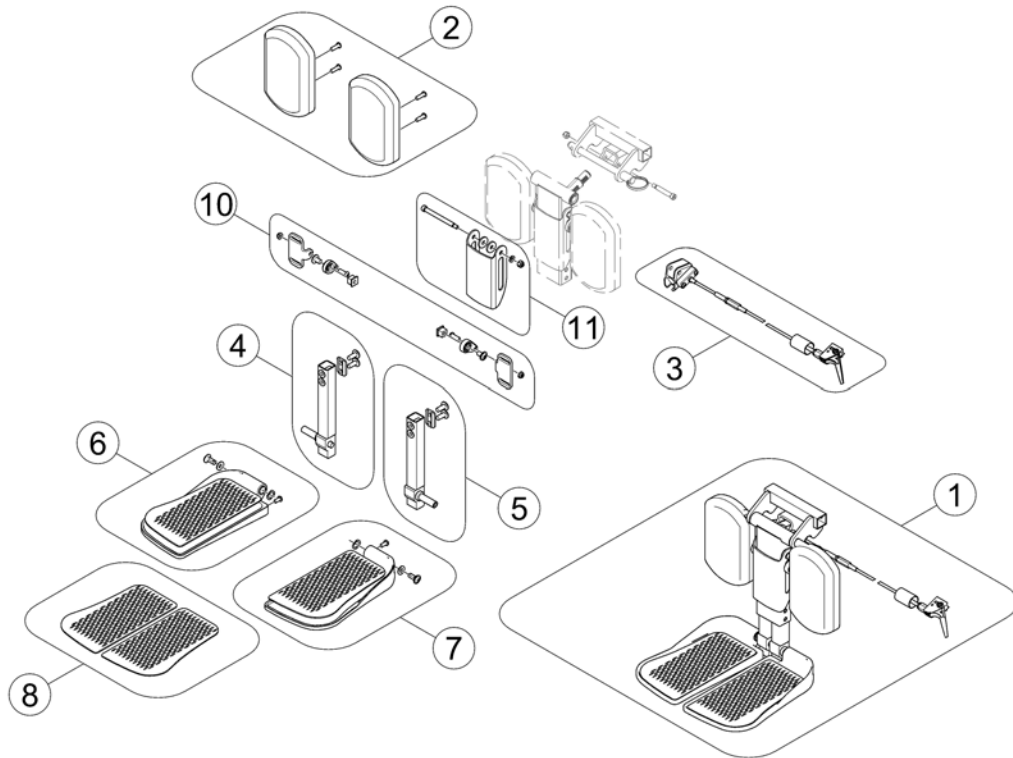
1



1	SP1537926	SENER BENSTØTTE LNX MAN TDX SP2 MODULITE	1
2	SP1163373	LEGGPUTER TDX (PAR)	1
3	SP1121079	(Bowden cable central legrest, manual)	1
4	SP1161986	(Kit, extension tube, cpl. RH)	1
5	SP1161987	(Kit, extension tube, cpl. LH)	1
6	SP1161990	(Footplate, cpl. RH)	1
7	SP1161991	(Footplate, cpl. LH)	1
8	SP1163375	MATTE FOTPLATE LNX BENSTØTTE LNX	1
10	SP1165505	(Kit, Adjustment Hardware Calf Pads)	1
11	SP1165506	(Kit, Calf Pad mounting bracket w/ hardware)	1


(Central Legrest pow. LNX (=> 2016-06-30))

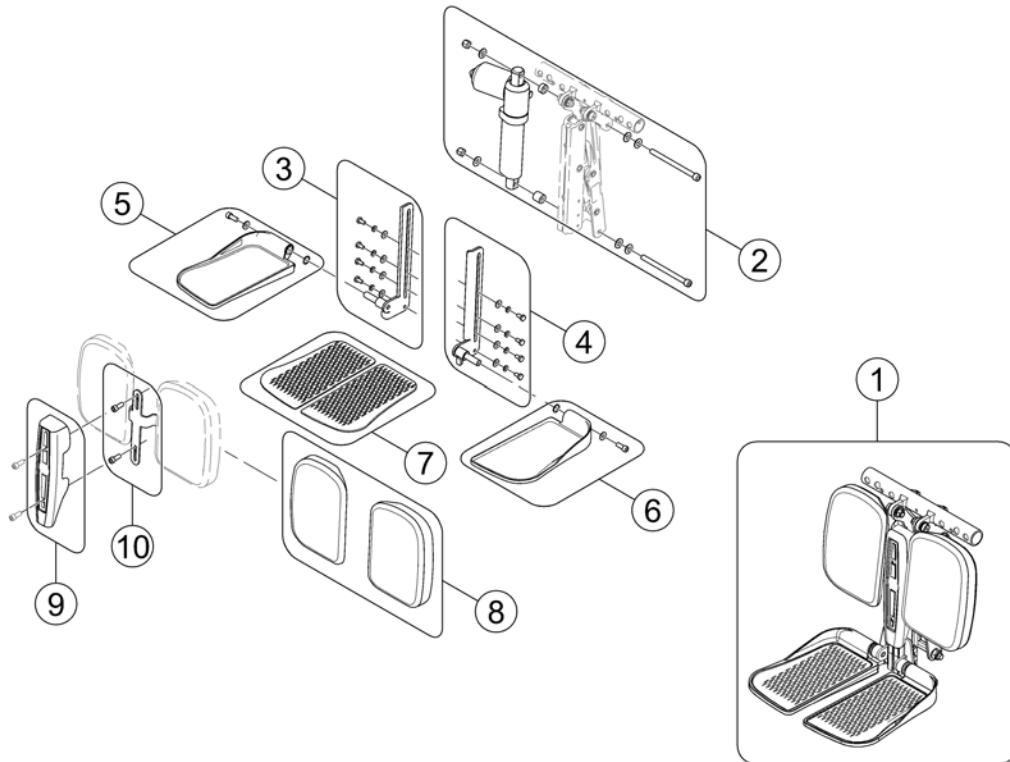1583894.kso



-	SP1541677	(Adaptercable 300mm) <i>For chairs with lifter</i>	1
-	SP1543012	(Adaptercable 1200mm) <i>For chairs without lifter</i>	1
1	-	(No longer available - please refer to Replacement Info for which part is to be used instead.) <i>SP1571959 & SP1571882 obsolete - Replacement: SPTRC1946</i>	1
2	SP1543031	(Legrest actuator, cpl. with mounting hardware)	1
3	SP1543035	(Footplate mount, RH)	1
4	SP1543034	(Footplate mount, LH)	1
5	SP1543037	FOTPLATE MEDIUM LNX VE FOTSTØTTE LNX	1
6	SP1543036	FOTPLATE MEDIUM LNX HØ BENSTØTTE LNX	1
7	SP1163375	MATTE FOTPLATE LNX BENSTØTTE LNX	1
8	SP1543033	(Calfpads (pair))	1
9	SP1543032	DEKSEL LNX BENSTØTTE LNX BENSTØTTE	1
10	SP1579985	(Calf Plate Bracket)	1

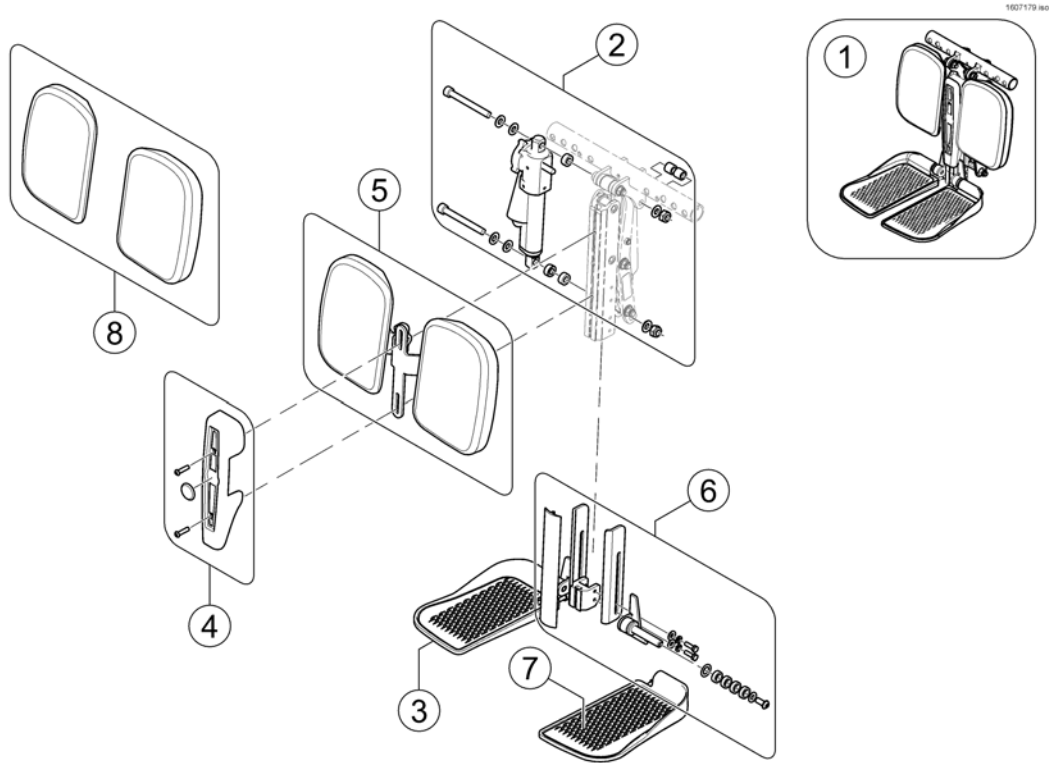(Central Legrest pow. LNX (2016-07-01 =>))



-	SP1573204	(Adater cable - 500mm)	1
1	SPTRC1946	(Central Legrest powered LNX, cpl. (with 240mm footplates))	1
1	SP1602089	(Central Legrest powered LNX, cpl. (with 240mm footplates) - SW 380mm only!)	1
2	SPTRC2124	(Legrest actuator, cpl. with mounting hardware)	1
3	SPTRC2125	(Foot plates 240mm depth (pair), cpl.)	1
4	SPTRC2126	(Cover, cpl. with mounting hardware)	1
5	SPTRC2035	(Calf Pads (pair), cpl. with mounting hardware)	1
6	SPTRC1664-C01	(Mounting hardware, foot plates)	1
7	SP1163375	MATTE FOTPLATE LNX BENSTØTTE LNX	1
8	SP1543033	(Calfpads (pair))	1

(Legrests for Ultra Low Maxx Seat)


-	(Central Legrests)	1
-	(Vari Legrests)	1
-	(Pivot Plus Legrests)	1

(Central Legrests)

-	(Supports, Central Legrests)	1
-	(Central Legrest pow. LNX)	1
-	(Central Legrest man.)	1

(Supports, Central Legrests)

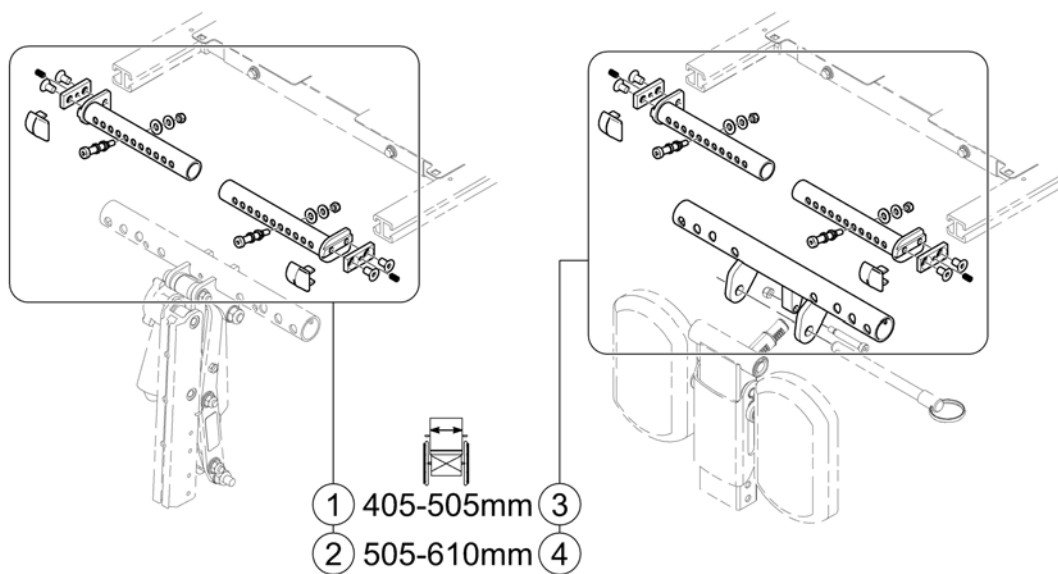
158402 Rev



1	SPTRC1945	(Legrest Support, el. - SW 405-505mm (16-20"))	1
2	SPTRC2096	(Legrest Support, el. - SW 505-610mm (19-24"))	1
3	SPTRC1947	(Legrest Support, man. - SW 405-505mm (16-20"))	1
4	SPTRC2097	(Legrest Support, man. - SW 505-610mm (19-24"))	1

(Central Legrest pow. LNX)

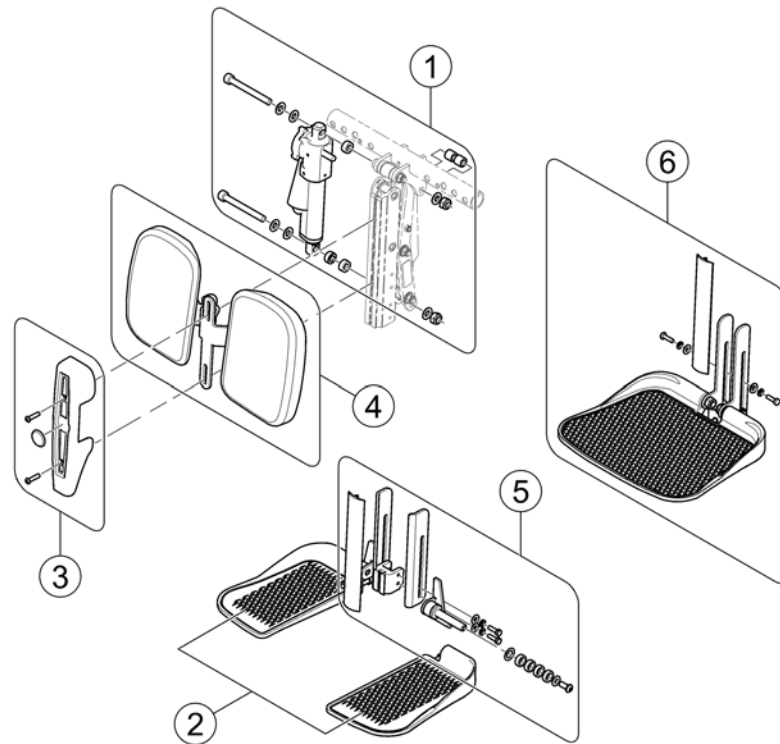
158403 no



-	SPTRC1946	(Central Legrest powered LNX, cpl. (with 240mm footplates))		1
-	SPTRE0756	(Actuator cable)		1
1	SPTRC2124	(Legrest actuator, cpl. with mounting hardware)		1
2	SPTRC2125	(Foot plates 240mm depth (pair), cpl.)		1
2	SPTRC2125_C03	(Foot plates 290mm depth (pair), cpl.)	01.01.2016 →	1
3	SPTRC2126	(Cover, cpl. with mounting hardware)		1
4	SPTRC2035	(Calf Pads (pair), cpl. with mounting hardware)		1
5	SPTRC1664-C01	(Mounting hardware, foot plates)		1
6	SPTRC2222	(Footboard 305x240mm, cpl. with mounting hardware)	01.01.2016 →	1

(Vari Legrests)


-	(Supports, Vari)	1
-	(Hanger Vari A (04/2009 =>))	1
-	(Hanger Vari F)	1
-	(Hangers ADM LH/RH manual)	1
-	(Hanger Vari ADE/ADP electric (07/2005 =>))	1
-	(Footrests Vari, Flip-Up (01.08.10 =>))	1
-	(Footrests Vari, Flip-Up & angle adjustable (01.08.10 =>))	1
-	(Footrests Vari, angle & depth adjustable)	1

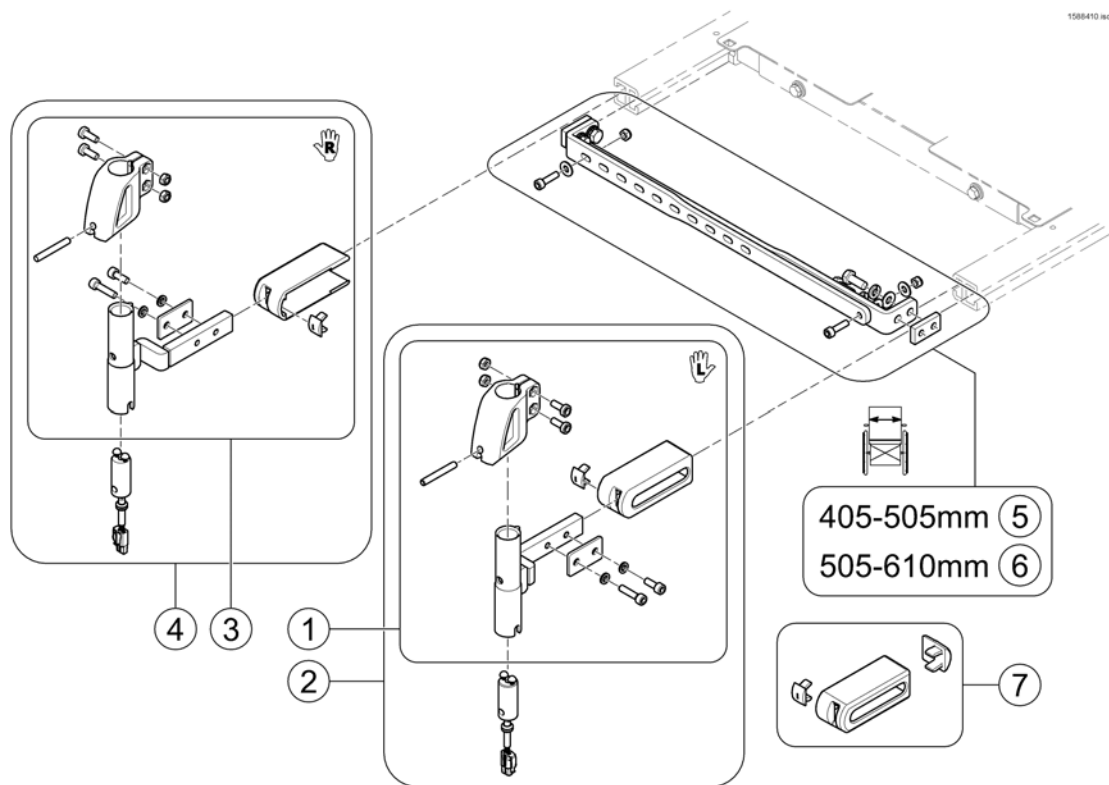

1	SP1586651	BENSTØTTE OVERDEL MAN KPL VE	1
2	SP1587187	BENSTØTTE OVERDEL EL KPL VE	1
3	SP1586706	BENSTØTTE OVERDEL MAN KPL HØ	1
4	SP1587186	(Legrest Support Vari RH el. cpl.)	1
5	SPTRC2098	(Cross brace SW 405-505mm (16-20"))	1
6	SPTRC2099	(Cross brace SW 505-610mm (19-24"))	1
7	SPTRC2061	(End cap kit)	1

(Pivot Plus Legrests)


-	(Cross Brace, Pivot Plus)	1
-	(Legrests)	1

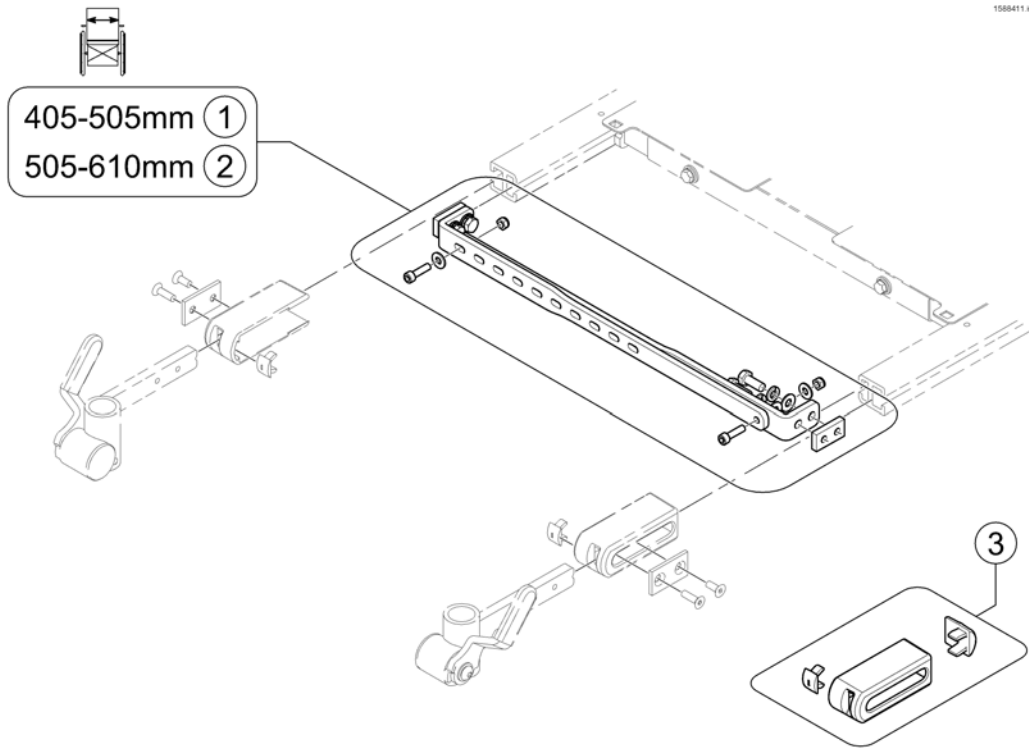

1	SPTRC2098	(Cross brace SW 405-505mm (16-20"))	1
2	SPTRC2099	(Cross brace SW 505-610mm (19-24"))	1
3	SPTRC2061	(End cap kit)	1

(Legrests)

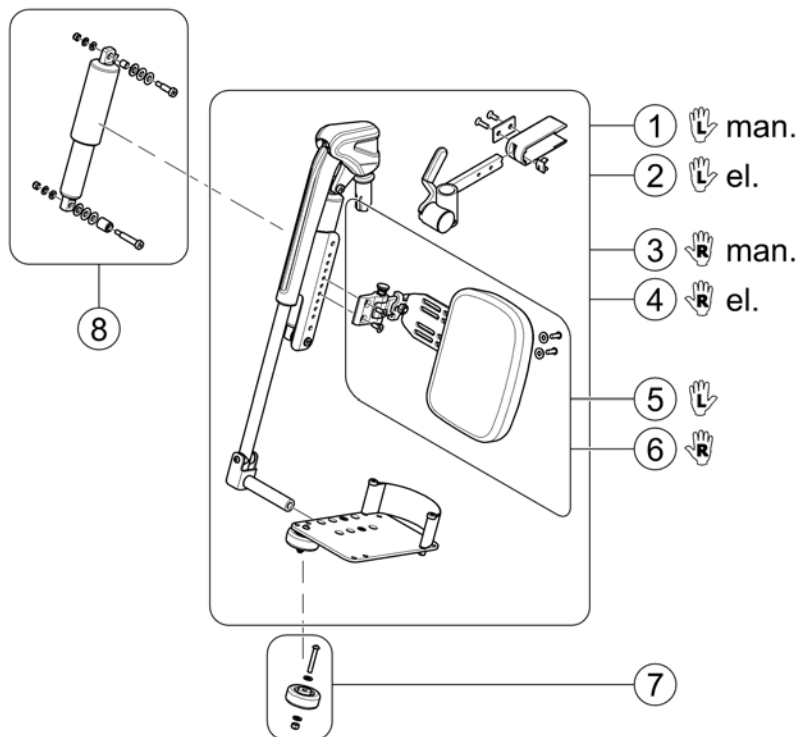1588425.iss




1	SPTRC2070	(Pivot Plus Legrest man. LH, cpl.)	1
2	SPTRC2068	(Pivot Plus Legrest el. LH, cpl.)	1
3	SPTRC2069	(Pivot Plus Legrest man. RH, cpl.)	1
4	SPTRC2067	(Pivot Plus Legrest el. RH, cpl.)	1
5	SPTRC2127-C02	(Calf Pad LH with mounting hardware)	1
6	SPTRC2127-C01	(Calf Pad RH with mounting hardware)	1
7	SPTRC2129	(Wheel, Footplate - cpl.)	1
8	SPTRC2128	(Legrest Actuator, cpl.)	1
-	SPTRE0756	(Actuator cable) <i>LH / links</i>	1
-	SPTRE0754	(Actuator Cable) <i>RH / rechts</i>	1

(Legrests for Max Seat)

-	(Legrest Support, Vari series - Max Seat)	0
-	(Hanger Vari A (04/2009 =>))	0
-	(Hanger Vari F)	0
-	(Hangers ADM LH/RH manual)	0
-	(Hanger Vari ADE/ADP electric (07/2005 =>))	0
-	(Footrests Vari, Flip-Up (01.08.10 =>))	0
-	(Footrests Vari, Flip-Up & angle adjustable (01.08.10 =>))	0
-	(Footrests Vari, angle & depth adjustable)	0
-	(Footboard Vari, pre-assembled)	0
	<i>[0522]</i>	
-	(Support Centr. Legrest man. LNX - Max Seat)	0
-	(Support Centr. Legrest LNX el. - Max Seat)	0
-	(Central Legrest man.)	0
-	(Central Legrest pow. LNX (=> 2016-06-30))	0
-	(Fixed footboard)	0
	<i>[0599]</i>	
-	(Calf plates)	0

(Legrest Support, Vari series - Max Seat)

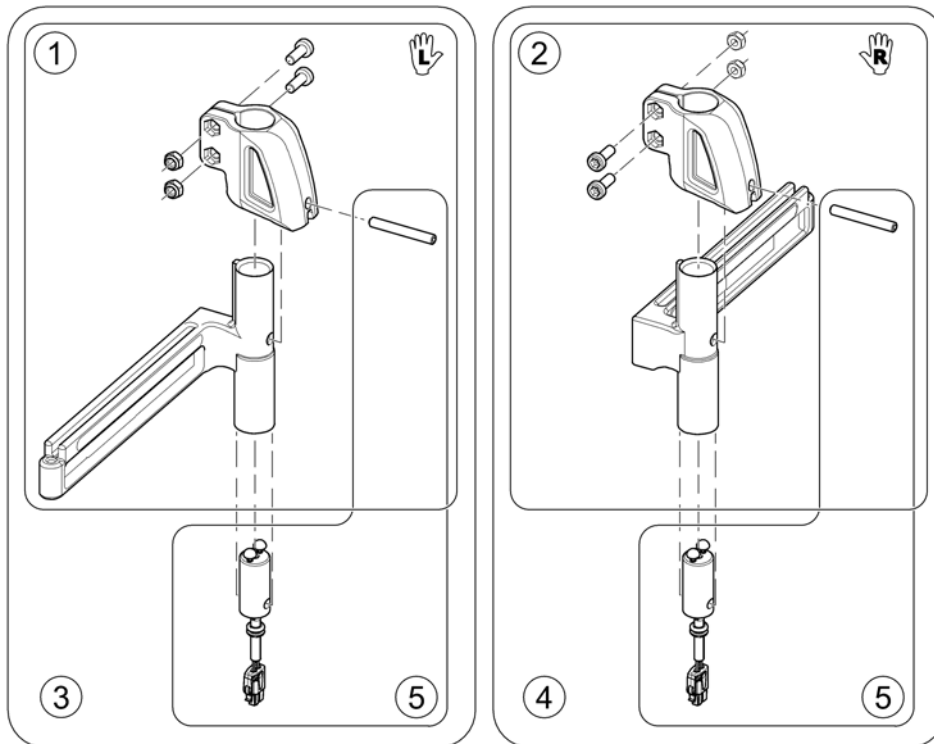
1506819_000



1	SP1526909	FESTE BENST MAN BORA	1
2	SP1526910	FESTE BENST MAN BORA	1
3	SP1555996	(Footrest/Legrest receptacle, powered, LH, 5cm)	1
4	SP1555995	(Footrest/Legrest receptacle, powered, RH, 5cm)	1
5	SP1530795	KABEL I HOLDER F BENST STORM4 BORA TDX	1

(Support Centr. Legrest man. LNX - Max Seat)

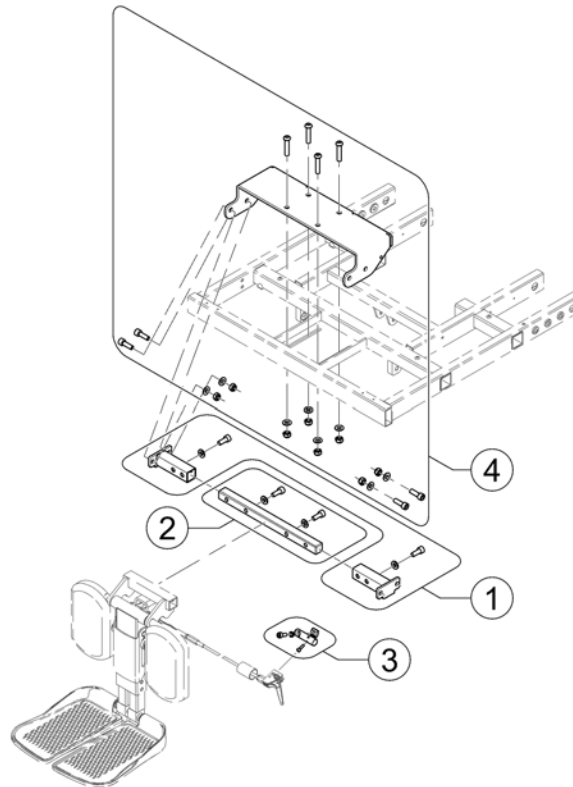
1583889.jpg



1	SP1583887	(Mounting Kit, Central Legrest LNX man. - Max Seat)	1
2	SP1576712	(Insert, Centr. Legrest LNX)	1
3	SP1579988	(Support, Lever, cpl.)	1
4	SP1580689	(Upper Support, Centr. Legrest LNX - Max Seat)	1


(Support Centr. Legrest LNX el. - Max Seat)

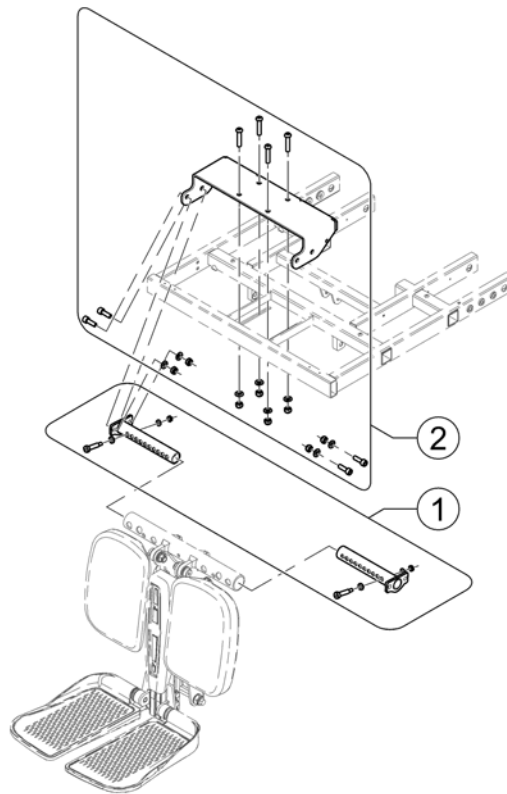1583890.kso



1	SP1583888	(Mounting Kit, Central Legrest LNX el. - Max Seat)	1
2	SP1580689	(Upper Support, Centr. Legrest LNX - Max Seat)	1

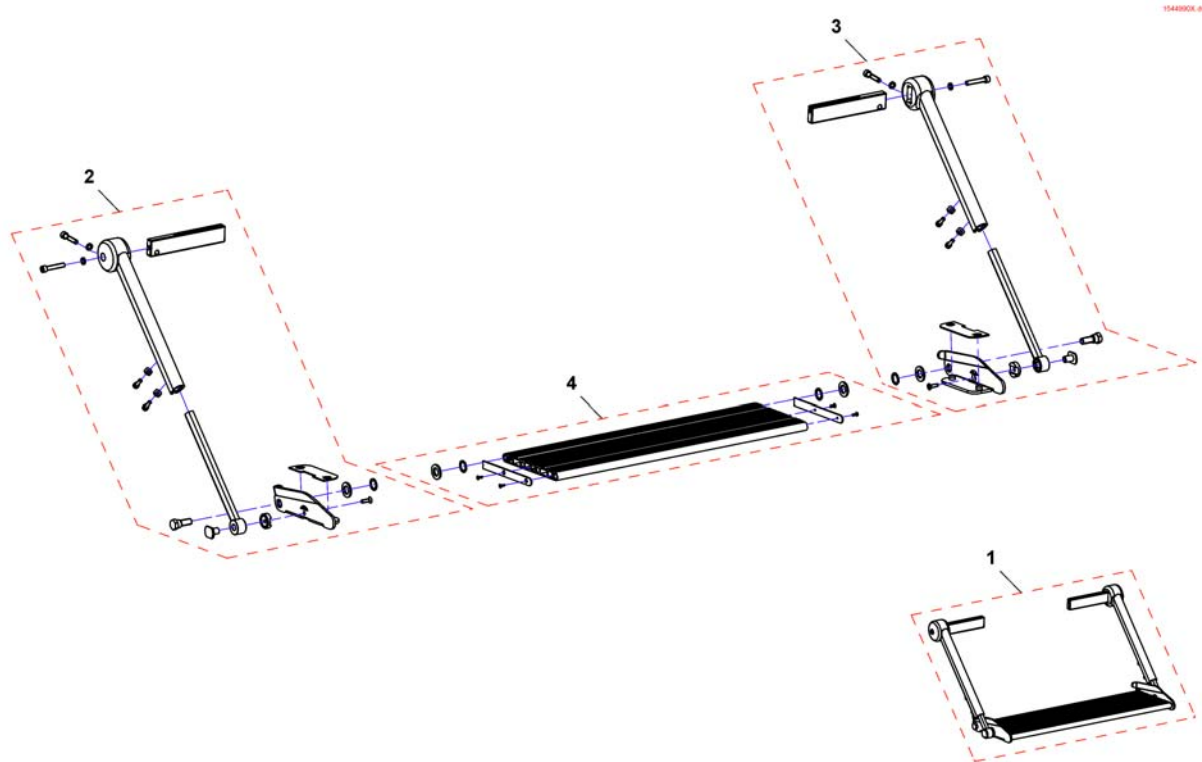
(Fixed footboard)




1	SP1544990	FOTBRETT 55 FAST STORM4 MAX	1
1	SP1545454	FOTBRETT 55 FAST STORM4 MAX	1
1	SP1545455	FOTBRETT FAST MAX SB49 STORM4 MAX	1
2	SP1552445	(Footboard bracket RH)	1
3	SP1552446	FESTE FOTBRETT KPL VE STORM4 MAX	1
4	SP1552466	(Footboard SW49)	1
4	SP1552465	(Footboard SW55)	1
4	SP1552447	(Footboard SW61)	1

(Calf plates)

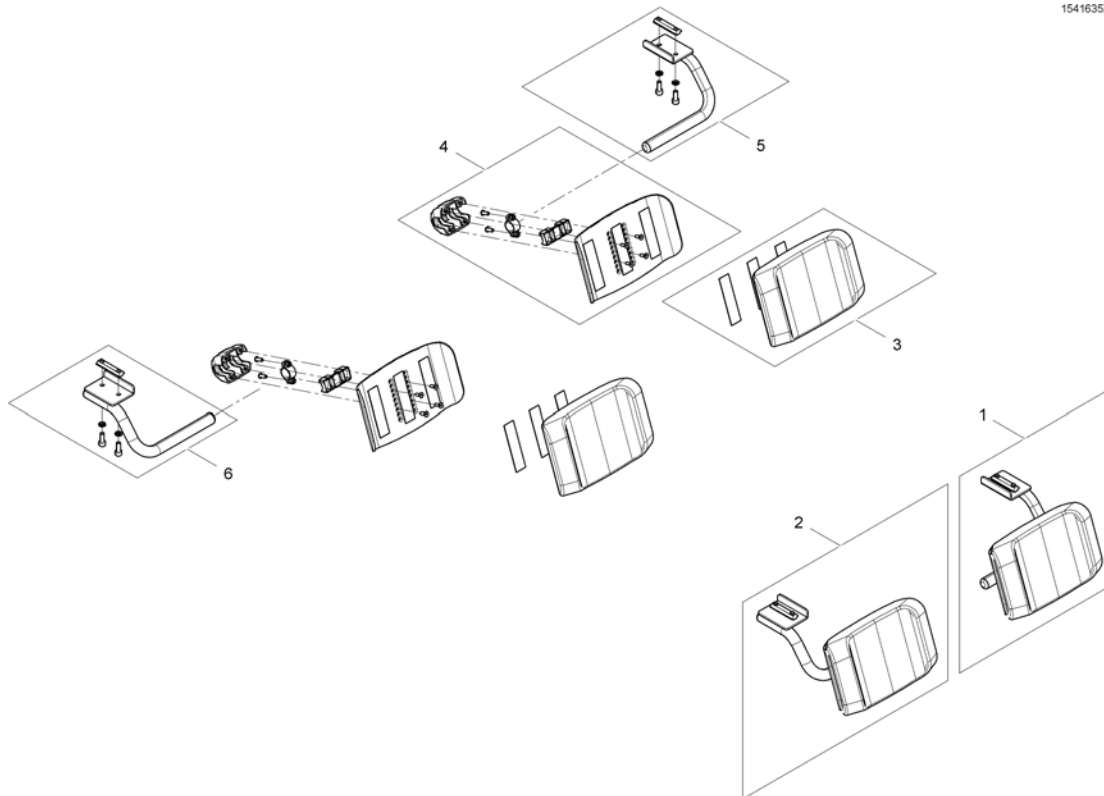1541635X.iso



1	SP1544875	(Calf pad, SW 490mm, LH, cpl.)	1
1	SP1541635	(Calf pad, SW 550mm & 610mm, LH, cpl.)	1
2	SP1544876	(Calf pad, SW 490mm, LH, cpl.)	1
2	SP1541636	(Calf pad, SW 550mm & 610mm, RH, cpl.)	1
3	SP1541672	(Calf Pad SW 490mm)	1
3	SP1541491	(Calf Pad SW 550mm & 610mm)	1
4	SP1541669	(Calf pad bracket SW 490mm, LH/RH, cpl.)	1
4	SP1541489	(Calf pad bracket SW 550mm & 610mm, LH/RH, cpl.)	1
5	SP1541633	(Calf pad bracket, LH, cpl.)	1
6	SP1541634	(Calf pad bracket, RH, cpl.)	1

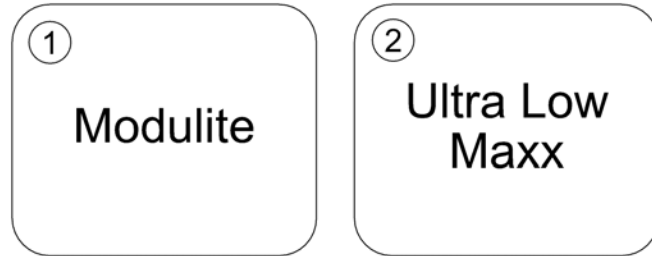

1	-	(Legrests for Modulite Seat)	0
2	-	(Legrests for Ultra Low Maxx Seat)	0

(Legrests for Modulite Seat)


-	(Legrests 80° & man. adjustable)	0
-	(Legrest Support, Vari series - Modulite Seat)	0
-	(Hanger Vari A (04/2009 =>))	0
-	(Hanger Vari F)	0
-	(Hangers ADM LH/RH manual)	0
-	(Hanger Vari ADE/ADP electric (07/2005 =>))	0
-	(Footrests Vari, Flip-Up (01.08.10 =>))	0
-	(Footrests Vari, Flip-Up & angle adjustable (01.08.10 =>))	0
-	(Footrests Vari, angle & depth adjustable)	0
-	(Central Legrests - Modulite Seat)	0



1	-	(Armrests for Modulite Seat)	0
2	-	(Armrests for Ultra Low Maxx Seat)	0
3	-	(Armrests for Max Seat)	0

(Armrests for Modulite Seat)

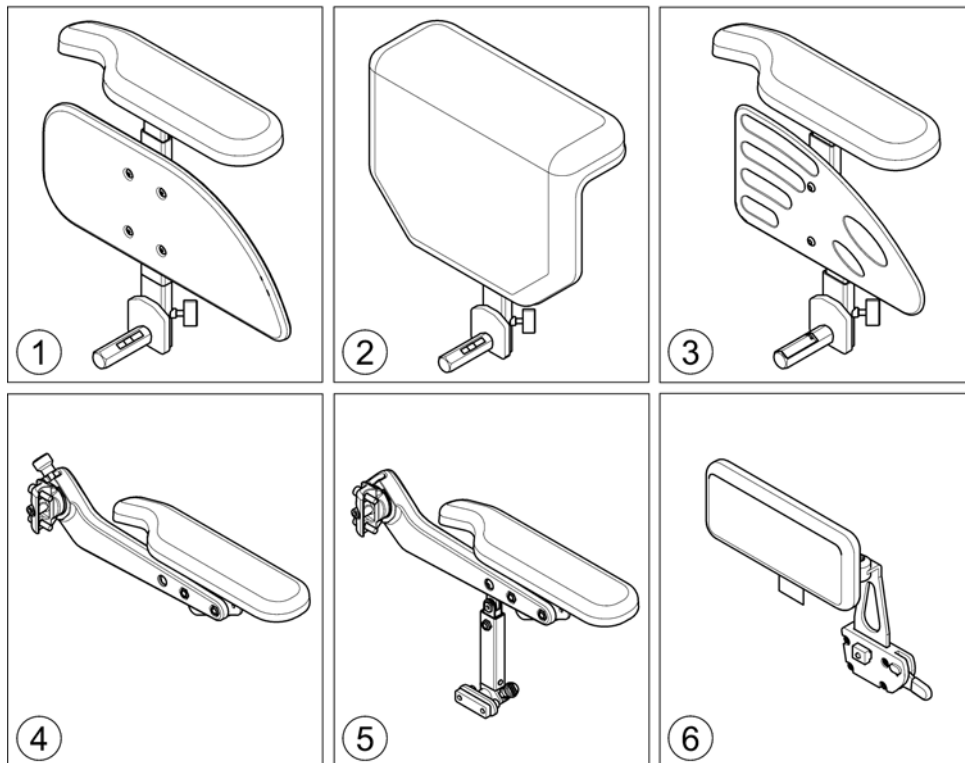
1571738.ilo



1	-	(T-Armrest) [2191+2195/2196/2197]	1
2	-	(Padded Armrest) [2198]	1
3	-	(Inrigged Armrest) [2183+2195/2196/2197/2199]	0
4	-	(Flip-Up Armrest) [2193+2195/2196/2197/2199]	0
5	-	(Following Armrest) [2192/2194+2195/2196/2197/2199]	0
6	-	(Hip Support) [2184/2185/2186/2187/2188/2189]	0

(T-Armrest)

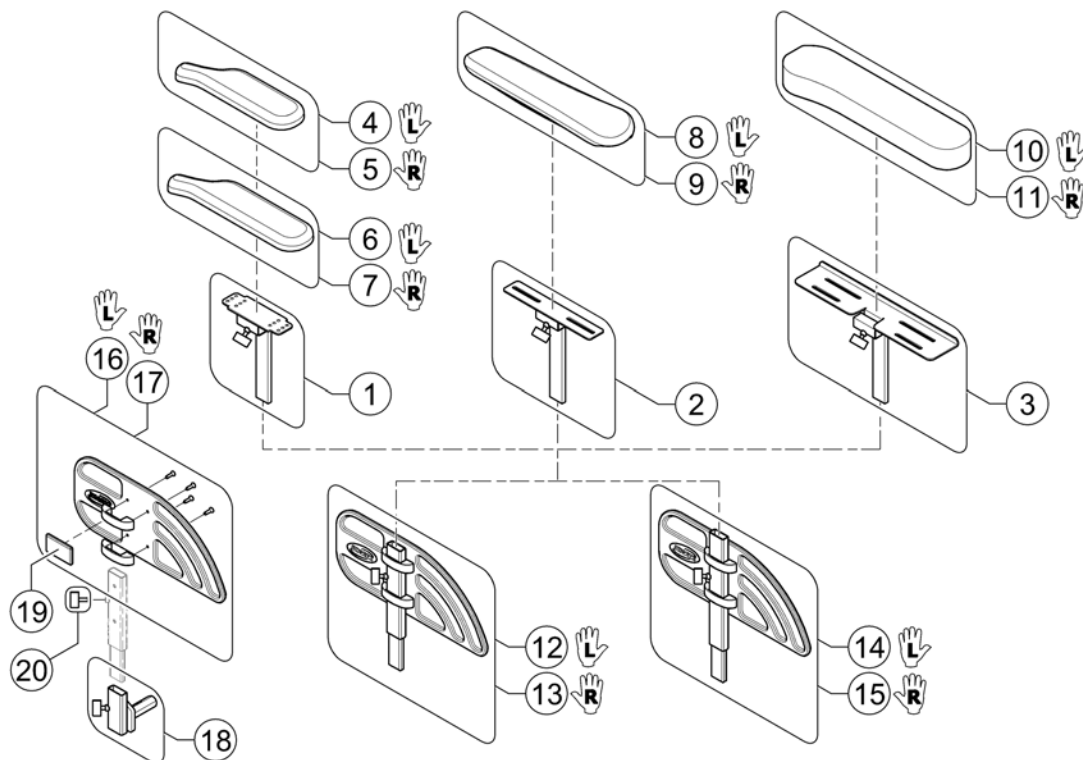
1572570 iso




1	SP1567750	FESTE POSTER ARMLEN TDX	1
2	SP1570319	FESTE ARMLEN ERGONOMISK MODULITE	1
3	SP1585088	(Armrest Bracket for Soft Armpad)	1
4	SP1443004	ARMLENSPOLSTER VE STORM STORM TDX BORA	1
5	SP1449043	ARMLENSPOLSTER HØ STORM STORM TDX BORA	1
6	SP1528846	ARMLENSPOLSTER VE LANG STORM TDX BORA	1
7	SP1528845	ARMLENSPOLSTER HØ LANG STORM4 TDX BORA	1
8	SP1571589	ARMLENSPOLSTER ERGONO VE STORM4 TDX	1
9	SP1571588	ARMLENSPOLSTER ERGONO HØ STORM4 TDX	1
10	SP1585068	ARMLENSPOLSTER VE MYKT TDX SP2	1
11	SP1585069	ARMLENSPOLSTER HØ MYKT TDX SP2	1
12	SP1588374	ARMLEN U POLSTER VE TDX SP2 MODULITE	1
13	SP1588375	ARMLEN U POLSTER HØ TDX SP2 MODULITE	1
14	SP1601882	SIDESTØTTE HØY VE UTEN FESTE	1
15	SP1601883	SIDESTØTTE HØY HØ UTEN FESTE	1
16	SP1572940	KLESBESKYTTER VE MODULITE	1
17	SP1572948	KLESBESKYTTER HØ MODULITE	1
18	SP1567749	BRAKETT KLESBESKYTTER MODULITE	1
19	SPP31298	(Reflector, yellow)	1
20	SPL10024	(Wing knobs M6x16)	1

(Padded Armrest)

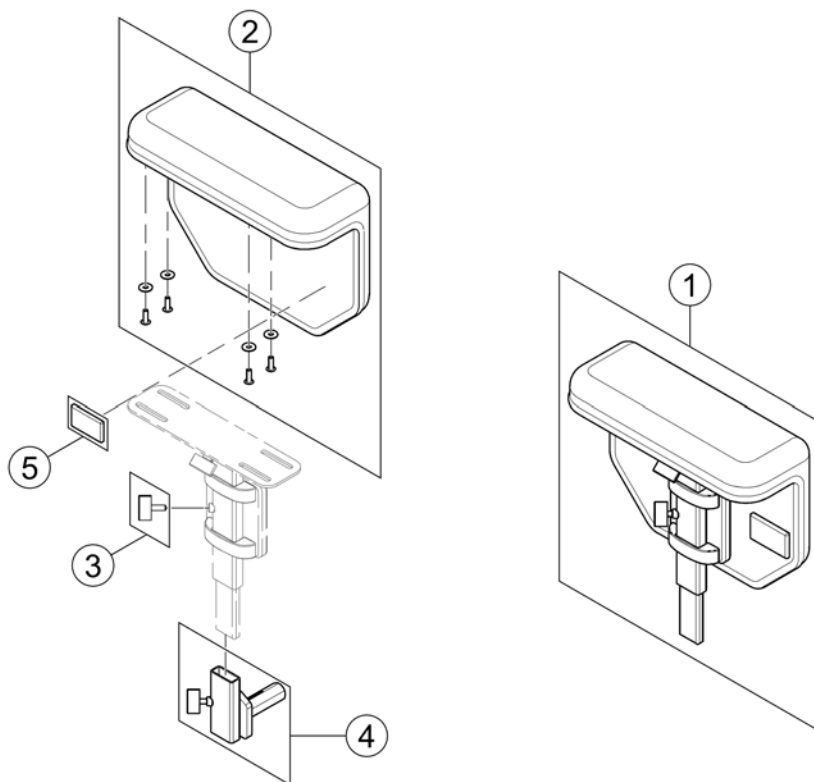1573214.no




1	SP1571987	ARMLEN PADDED KPL U/FEST MODULITE HØ/VE	1
2	SP1571962	ARMLENSPOLSTER PADDED MODULITE	1
3	SPL10024	(Wing knobs M6x16)	1
4	SP1567749	BRAKETT KLESBESKYTTER MODULITE	1
5	SPP31298	(Reflector, yellow)	1

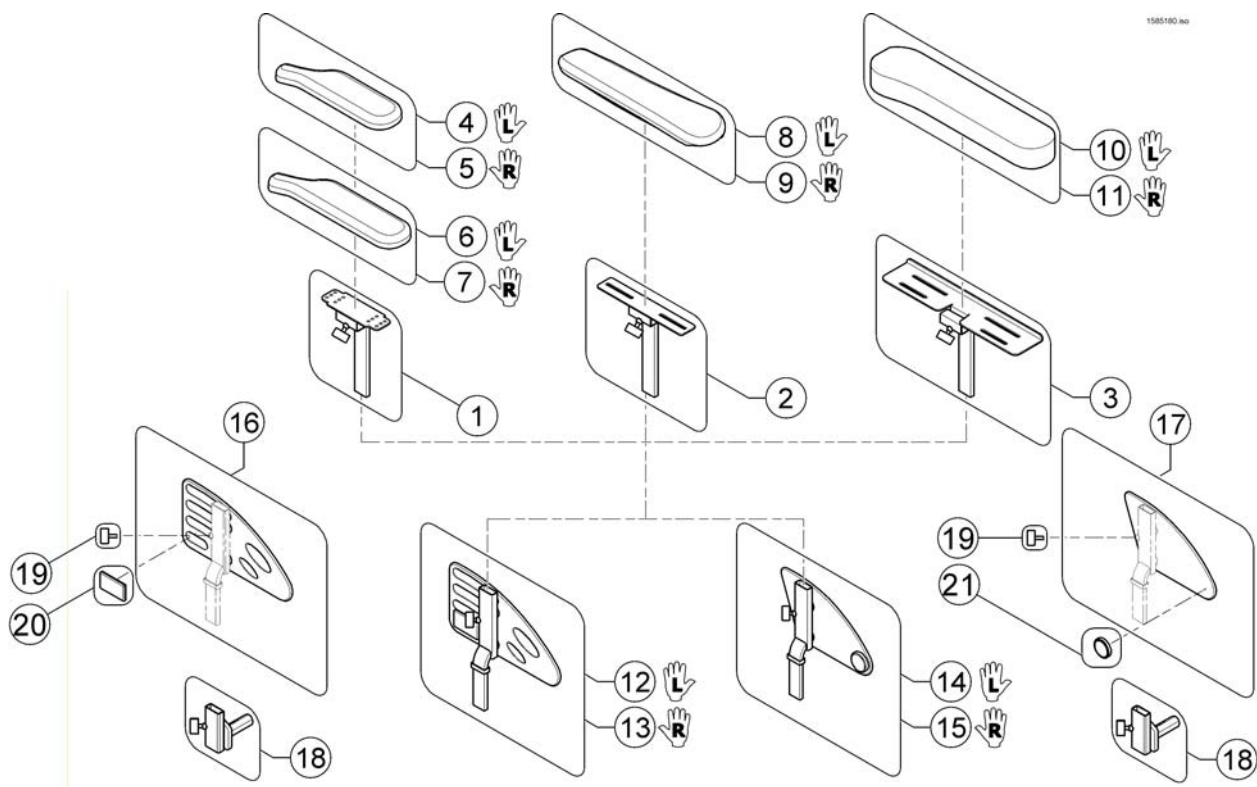(Inrigged Armrest)

1	SP1567750	FESTE POSTER ARMLEN TDX	1
2	SP1570319	FESTE ARMLEN ERGONOMISK MODULITE	1
3	SP1585088	(Armrest Bracket for Soft Ampad)	1
4	SP1443004	ARMLENSPOLSTER VE STORM STORM TDX BORA	1
5	SP1449043	ARMLENSPOLSTER HØ STORM STORM TDX BORA	1
6	SP1528846	ARMLENSPOLSTER VE LANG STORM TDX BORA	1
7	SP1528845	ARMLENSPOLSTER HØ LANG STORM4 TDX BORA	1
8	SP1571589	ARMLENSPOLSTER ERGONO VE STORM4 TDX	1
9	SP1571588	ARMLENSPOLSTER ERGONO HØ STORM4 TDX	1
10	SP1585068	ARMLENSPOLSTER VE MYKT TDX SP2	1
11	SP1585069	ARMLENSPOLSTER HØ MYKT TDX SP2	1
12	SP1585037	(Siderest inrigged LH)	1
13	SP1585038	(Siderest inrigged RH)	1
14	SP1605654	(Siderest inrigged LH (SD410mm))	1
15	SP1605653	(Siderest inrigged RH (SD410mm))	1
16	SP1578119	KLESBESKYTTER HØ/VE TDX	1
17	SP1603118	(Skirt Guard LH/RH (SD410mm))	1
18	SP1567749	BRAKETT KLESBESKYTTER MODULITE	1
19	SPL10024	(Wing knobs M6x16)	1
20	SPP31298	(Reflector, yellow)	1
21	SP1582989	(Reflector yellow D=43/38mm)	1

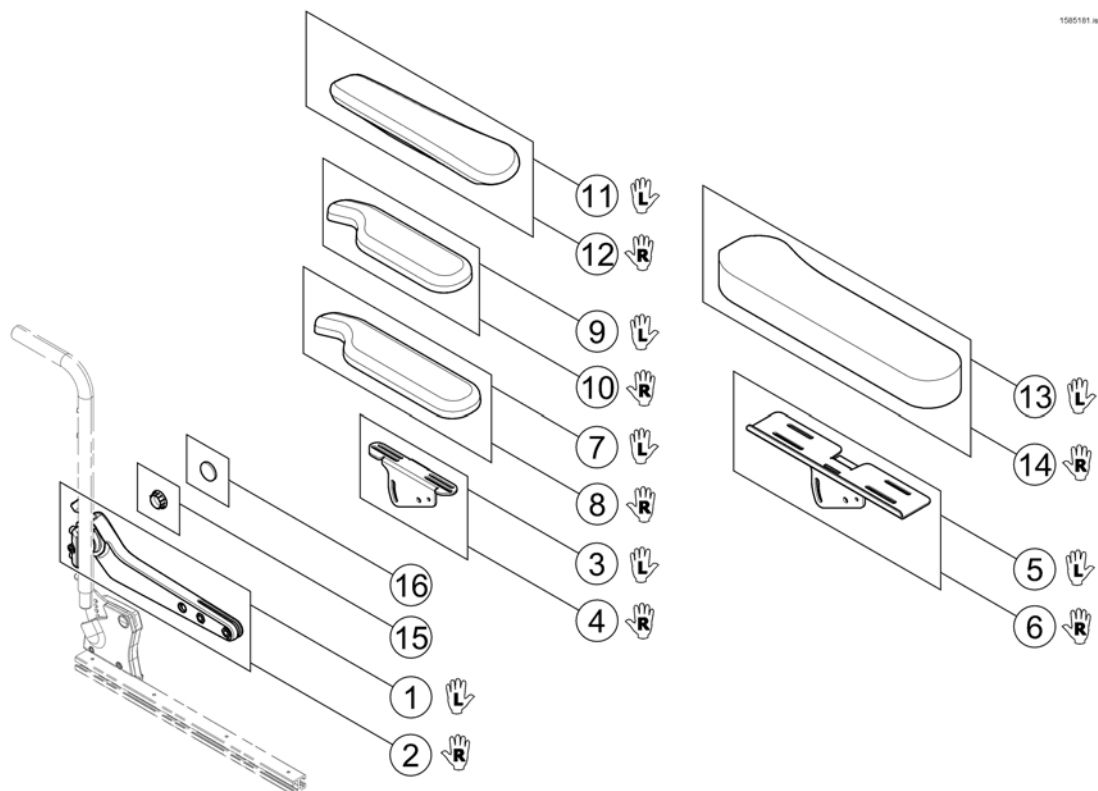
(Flip-Up Armrest)



1580181 no

		(Retrofitting the armrests on chairs built before 15-03-27 also requires new Pos. A & B and Sling Back with respective cut-outs!) !!!		0
1	SP1582668	ARMLLEN FLIP UP VE KPL MUDULITE	19.01.2015 →	1
2	SP1582644	ARMLLEN FLIP UP HØ KPL MODULITE	19.01.2015 →	1
3	SP1585412	HOLDER VE POLSTER OPPFEL MODULITE		1
4	SP1585413	HOLDR HØ POLSTER OPPFELL MODULITE		1
5	SP1585414	FESTE MYKT POLSTER HØ MODULITE		1
6	SP1585415	FESTE MYKT POLSTER VE MODULITE		1
7	SP1528846	ARMLLENSPOLSTER VE LANG STORM TDX BORA		1
8	SP1528845	ARMLLENSPOLSTER HØ LANG STORM4 TDX BORA		1
9	SP1443004	ARMLLENSPOLSTER VE STORM STORM TDX BORA		1
10	SP1449043	ARMLLENSPOLSTER HØ STORM STORM TDX BORA		1
11	SP1571589	ARMLLENSPOLSTER ERGONO VE STORM4 TDX		1
12	SP1571588	ARMLLENSPOLSTER ERGONO HØ STORM4 TDX		1
13	SP1585068	ARMLLENSPOLSTER VE MYKT TDX SP2		1
14	SP1585069	ARMLLENSPOLSTER HØ MYKT TDX SP2		1
15	SP1590599	(Cap)		1
16	SP1582989	(Reflector yellow D=43/38mm)		1

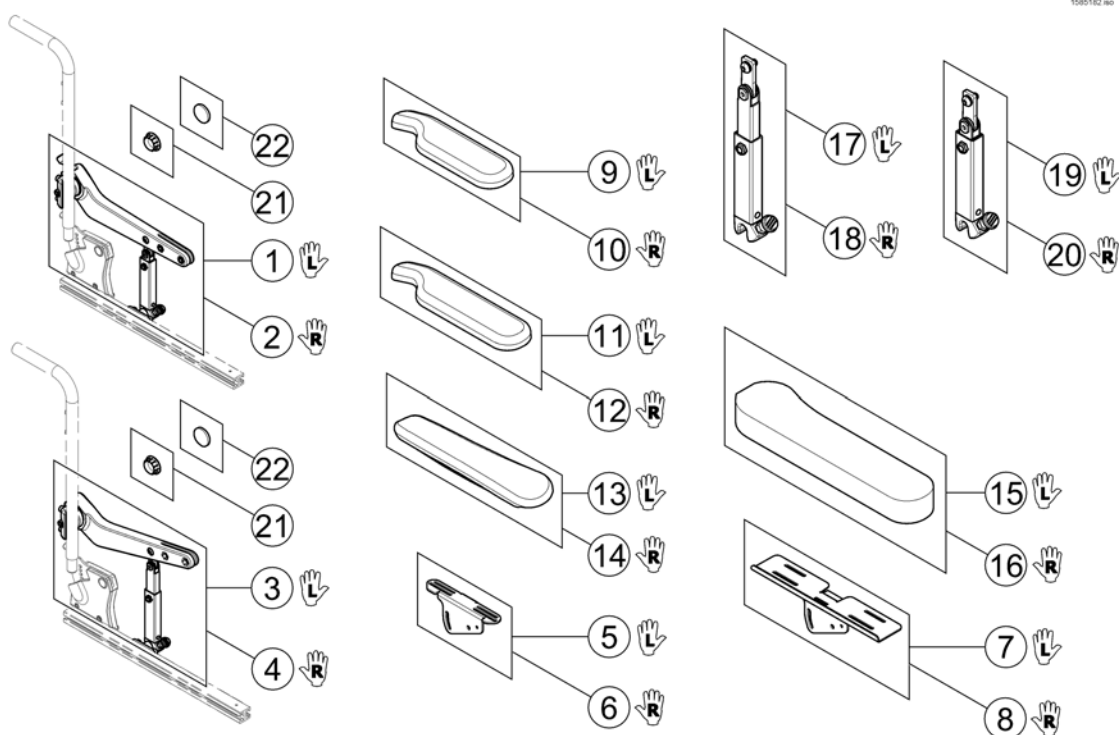
(Following Armrest)



-	-	(Retrofitting the armrests on chairs built before 15-03-27 also requires new Pos. A & B and Sling Back with respective cut-outs!) !!!		0
1	SP1585263	ARMLÉN OPPFELLBAR VE MODULITE KORT	19.01.2015 →	1
2	SP1585264	ARMLÉN OPPFELLBAR HØ MODULITE KORT	19.01.2015 →	1
3	SP1585267	ARMLÉN OPPFELL VE LANG MODULITE FOLLOWING	19.01.2015 →	1
4	SP1585268	ARMLÉN OPPFELLB HØ LANG MODULITE FOLLOWING	19.01.2015 →	1
5	SP1585412	HOLDER VE POLSTER OPPFELL MODULITE		1
6	SP1585413	HOLDR HØ POLSTER OPPFELL MODULITE		1
7	SP1585414	FESTE MYKT POLSTER HØ MODULITE		1
8	SP1585415	FESTE MYKT POLSTER VE MODULITE		1
9	SP1443004	ARMLÉNSPOLSTER VE STORM STORM TDX BORA		1
10	SP1449043	ARMLÉNSPOLSTER HØ STORM STORM TDX BORA		1
11	SP1528846	ARMLÉNSPOLSTER VE LANG STORM TDX BORA		1
12	SP1528845	ARMLÉNSPOLSTER HØ LANG STORM4 TDX BORA		1
13	SP1571589	ARMLÉNSPOLSTER ERGONO VE STORM4 TDX		1
14	SP1571588	ARMLÉNSPOLSTER ERGONO HØ STORM4 TDX		1
15	SP1585068	ARMLÉNSPOLSTER VE MYKT TDX SP2		1
16	SP1585069	ARMLÉNSPOLSTER HØ MYKT TDX SP2		1
17	SP1583156	(Amrest Support, long, LH)		1
18	SP1583155	(Amrest Support, long, RH)		1
19	SP1583154	STØTTE ARMLÉN VE KORT		1
20	SP1582639	(Amrest Support, short, RH)		1
21	SP1590599	(Cap)		1
22	SP1582989	(Reflector yellow D=43/38mm)		1

(Hip Support)

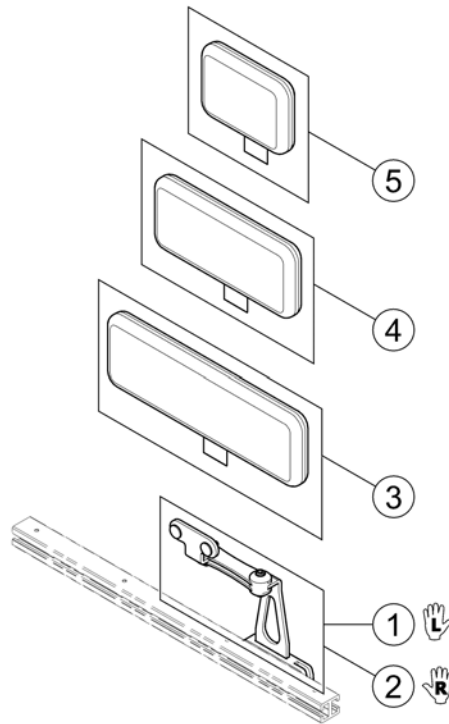1580183 no




1	SP1584411	FESTE BRACKET LÅRSTØTTE TDX SP2 MODULITE VE	1
2	SP1584410	FESTE SIDESTØTTE MODULIT STORM4 TDX LOW-RIDER HØ	1
3	SP1584414	LÅRSTØTTE PUTE 33X10 MODULITE	1
4	SP1584413	LÅRSTØTTE PUTE 23*10 CM MODULITE	1
5	SP1584412	LÅRSTØTTE PUTE 14*10 CM MODULITE	1


1	-	(Recline Flip-Up Armrest)	0
2	-	(Cantilever Flip-Up Armrest)	0
3	-	(Armpads, Desk & Full)	0
4	-	(Armpad, Ergonomic)	0
5	-	(Elbow Block)	0
6	-	(MACES Multi Axis Armpad Holder)	0
7	-	(Hip Support)	0

(Recline Flip-Up Armrest)

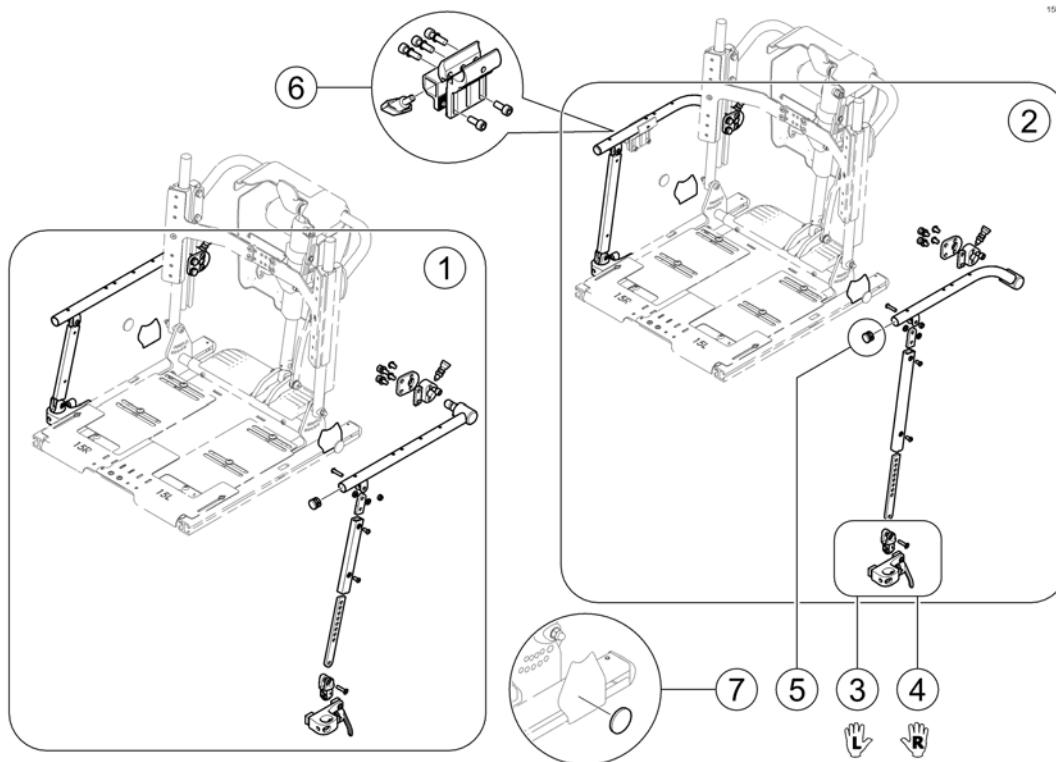1588503 Rev




1	SPTRC1933	(Reclining Arm Flip-Back 240-330mm)	1
2	SPTRC1934	(Reclining Arm Flip-Back 330-405mm)	1
3	SPTRC2135	(Flip Back Lever with connector LH)	1
4	SPTRC2134	(Flip Back Lever with connector RH)	1
5	SPTRM0130	(Armrest-pipe cap)	1
6	SPTRC2060	(Clamp with square tube adapter)	1
7	SP1582989	(Reflector yellow D=43/38mm)	1

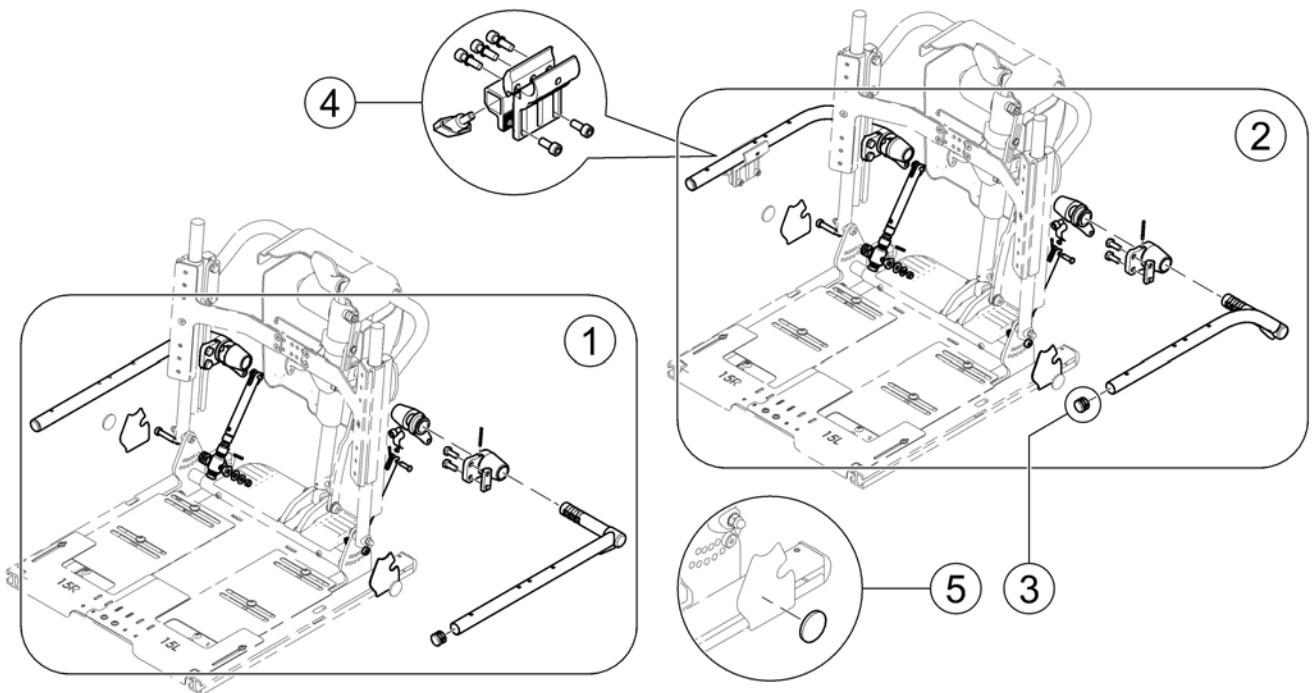


1	SPTRC1936	(Cantilever Flip-Back Desk 240-330mm)	1
2	SPTRC1937	(Cantilever Flip-Back Desk 320-405mm)	1
3	SPTRM0130	(Armrest-pipe cap)	1
4	SPTRC2060	(Clamp with square tube adapter)	1
5	SP1582989	(Reflector yellow D=43/38mm)	1

(Armpads, Desk & Full)

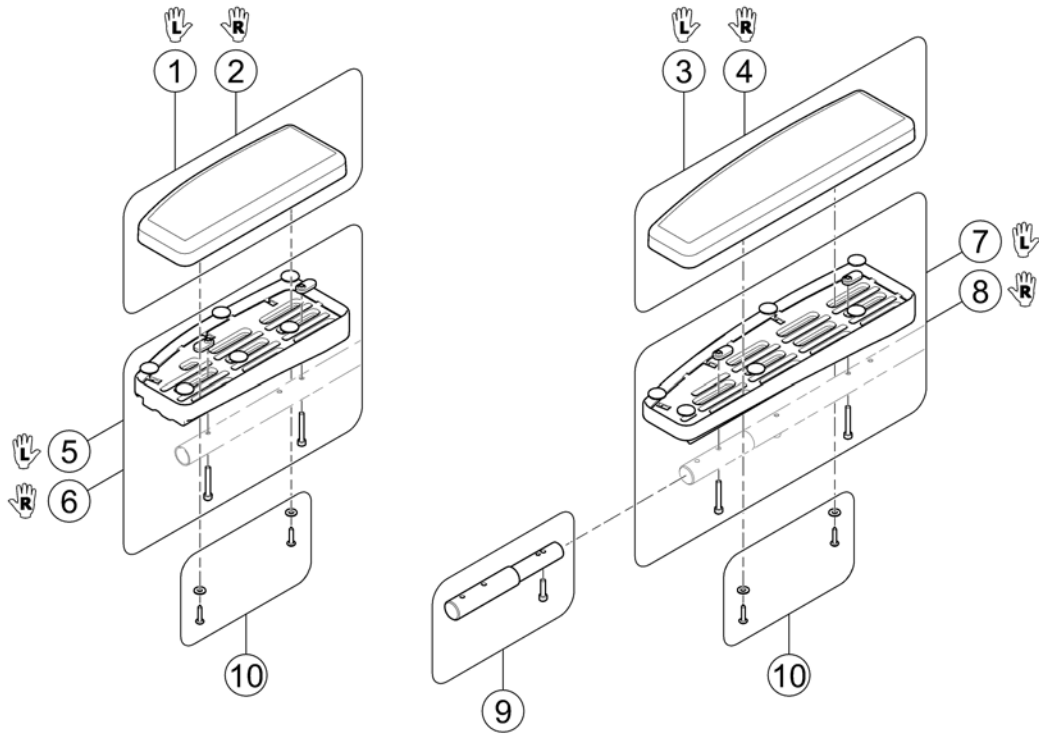158805 iso




1	SPTRX8038	(Armpad, Desk Amrest, LH)	1
2	SPTRX8037	(Armpad, Desk Amrest, RH)	1
3	SPTRX8036	(Armpad, Full Amrest, LH)	1
4	SPTRX8035	(Armpad, Full Amrest, RH)	1
5	SPTRC2133	(Desk Armrest Body LH with mounting hardware)	1
6	SPTRC2132	(Desk Armrest Body RH with mounting hardware)	1
7	SPTRC2131	(Full Armrest Body LH with mounting hardware)	1
8	SPTRC2130	(Full Armrest Body RH with mounting hardware)	1
9	SPTRC1935	(Desk to Full Armrest Extension Kit) <i>Needed for Pos. 7 & 8!</i>	1
10	SP1602292	(Screw set for armpad fixation) <i>[1494]</i>	1

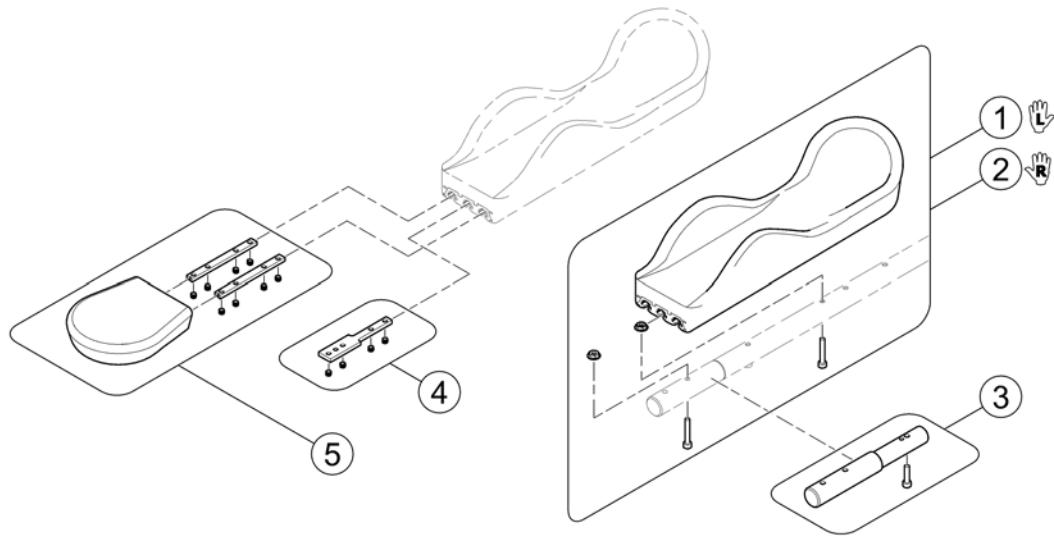


1	SPTRC2100-C02	(Ergonomic Arm Trough LH)	1
2	SPTRC2100-C01	(Ergonomic Arm Trough RH)	1
3	SPTRC1935	(Desk to Full Armrest Extension Kit) <i>Needed for Pos. 1 & 2!</i>	1
4	SPTRC2102	(Mounting Hardware for Otto Bock)	1
5	SPTRC2101	(Flat Hand Pad)	1

(Elbow Block)

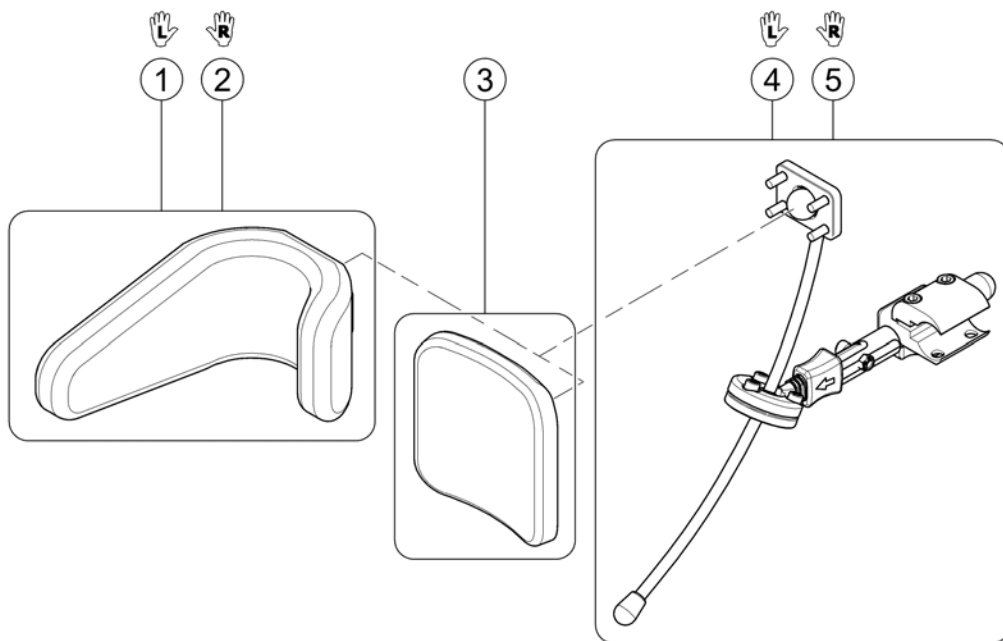1603428.iss



1	SPTRM1059	(Elbow Block pad 3x10 LH)	1
2	SPTRM1060	(Elbow Block pad 3x10 RH)	1
3	SPTRM1058	(Elbow Block pad 4x4)	1
4	SPTRM1055	(Elbow Block Hardware, cpl. LH)	1
5	SPTRM1056	(Elbow Block Hardware, cpl. RH)	1

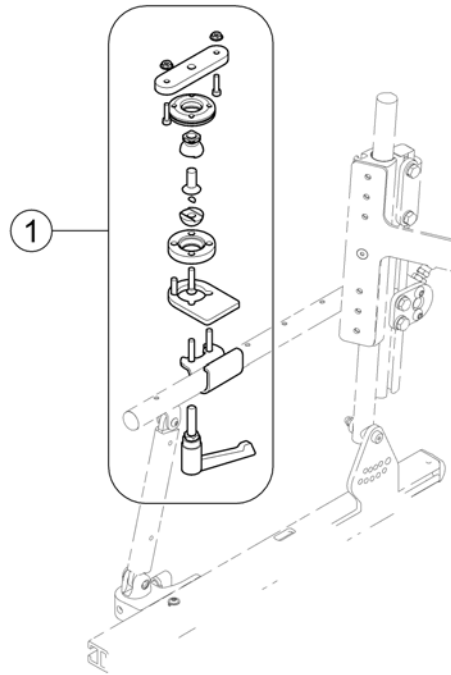


1 SPTRC1938

(MACES Armpad Holder Kit, Multi Axis)

1

(Armrests for Max Seat)

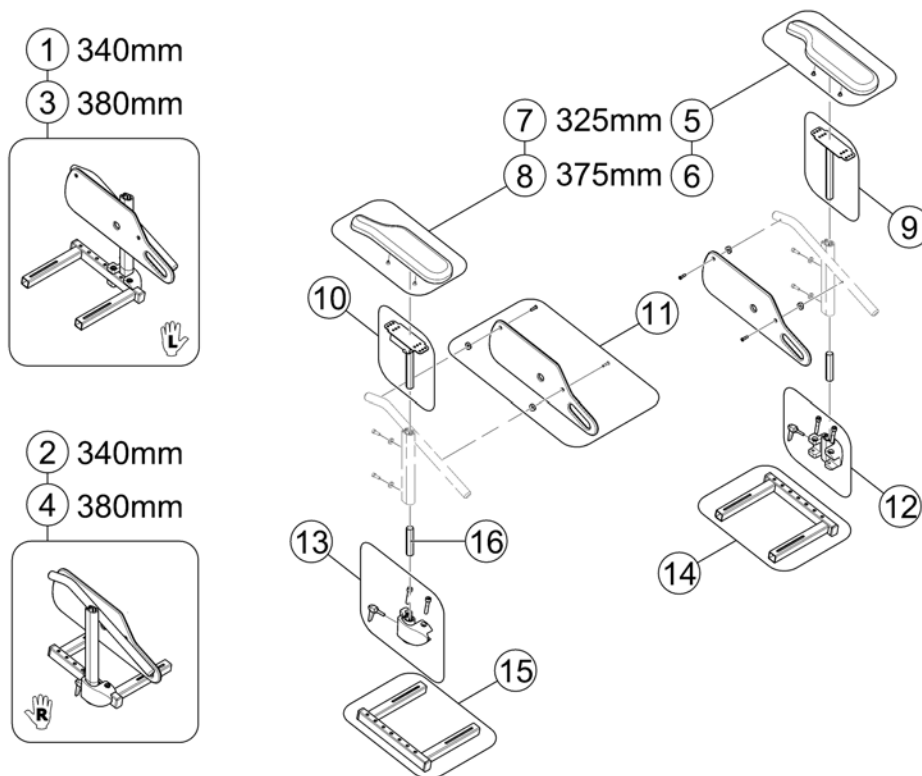
1583900.kim



1	SP1583906	(Siderest 340mm LH, without Armrest - Max Seat)	1
2	SP1583907	(Siderest 340mm RH, without Armrest - Max Seat)	1
3	SP1583908	ARMLEN VE H 38 KPL MAX STORM4 MAX	1
4	SP1583909	ARMLEN HØ H38 KPL MAX STORM4 MAX	1
5	SP1443004	ARMLENSPOLSTER VE STORM STORM TDX BORA	1
6	SP1528846	ARMLENSPOLSTER VE LANG STORM TDX BORA	1
7	SP1449043	ARMLENSPOLSTER HØ STORM STORM TDX BORA	1
8	SP1528845	ARMLENSPOLSTER HØ LANG STORM4 TDX BORA	1
9	SP1510857	ARMLENSPOLSTER HOLDER VE STORM3/4 2009->	1
10	SP1510856	ARMLENSPOLSTER HOLDER HØ STORM3/4 2009->	1
11	SP1428078	KLESBESKYTTER STORM4 TDX	1
12	SP1541647	(Clamp, LH, cpl.)	1
13	SP1541646	(Clamp, RH, cpl.)	1
14	SP1537209	(Bracket, cpl., LH)	1
15	SP1537208	(Bracket, cpl., RH)	1
16	1422092-0002	(CONNECTION PIECE)	1


(Castors)

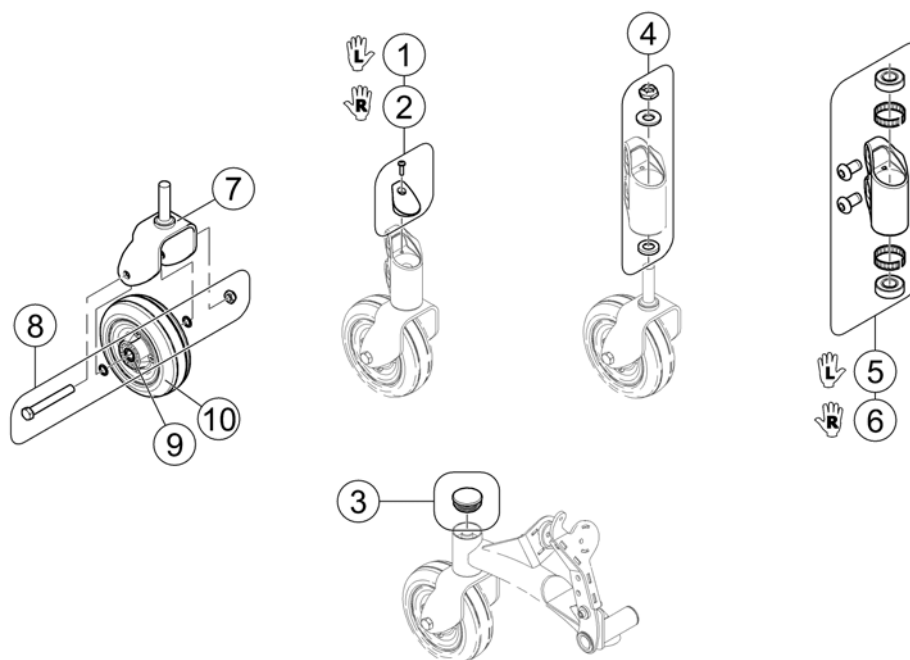1583268.no



1	SP1523141	TOPPHETTE FORAN VE TDX	1
2	SP1523140	TOPPHETTE FORAN HØ TDX	1
3	SP1110709	(Plug Button, Head Tube Black)	1
4	SP1522999	GAFFELHUS TDX	1
5	SP1523139	GAFFELHUS M FESTE VE TDX	1
6	SP1523094	GAFFELHUS M FESTE HØ TDX	1
7	SP1523016	GAFFEL TDX	1
8	SP1537431	AKSLINGSBOLT SVINGHJUL TDX	1
9	SP1522935	DISTANSELAGER TDX	1
10	SP1175881	SVINGHJUL TDX SP2	1

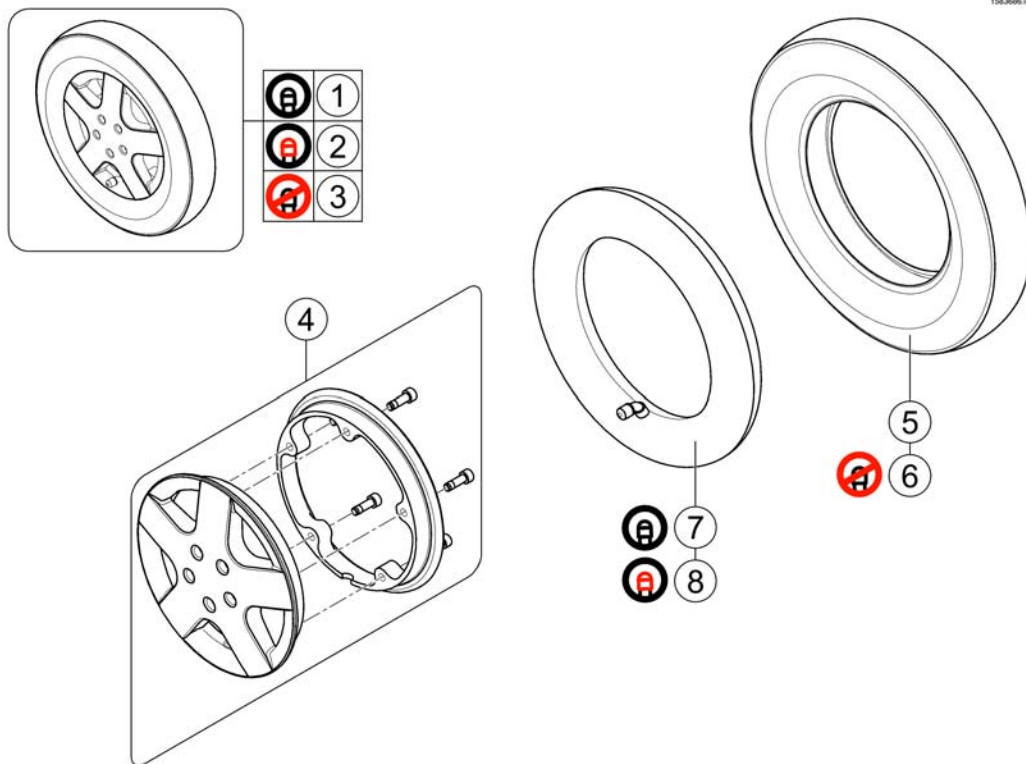
(Drive Wheels 14")



1	SP1532414	DRIVHJUL KPL M LUFT STORM4	1
2	SP1532415	DRIV HJUL PUNKT SIKKERT STORM4 KPL	1
3	SP1532416	HJUL KOMPAKT 3.00-8 GRÅ STORM4	1
4	SP1541373	FELG KPL DRIVHJUL STORM4	1
5	SP1522940	DEKK LUFT 14X 3 (3,00-8) STORM4	1
6	SP1141342	DEKK KOMPACT GB2 14X3 STORM4	1
7	SP1039651	(Tube 14")	1
8	SP1541374	(Tube 14" puncture protected)	1

(Wheel Hub 5-hole)

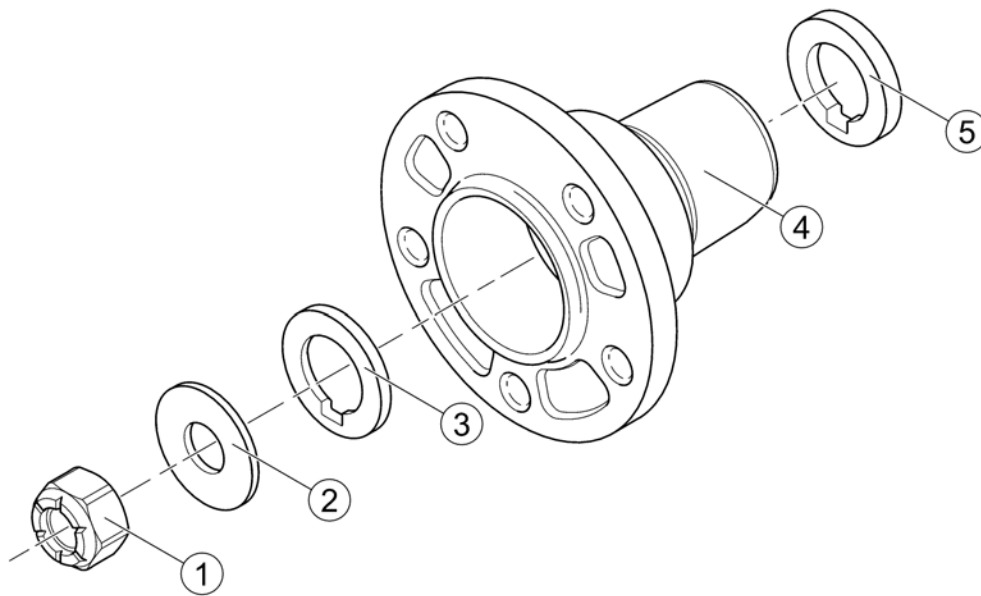
1583687 Rev



1	F00108	(Hexagon nut DIN985-M12x1,5 zc)	1
2	G00104	SKIVE 13X24X2,5 G40	1
3	SP1544342	(Washer with feather nut, 5mm thick, d=19mm)	1
4	SP1163728	(Wheel hub 5-hole)	1
5	SP1163729	(Spacer, Hub, 5mm thick)	1

(Overview Propulsion SSD)

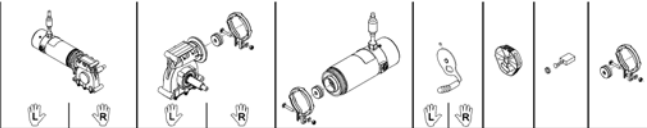1163090X.iso



Model	Speed	Gearbox	L		R		Cable		Other		
			1	2	33	34	37	---	---	51	45
Bora/ SpectraXTR	6km/h	ACS/LINX /Shark	1	2	33	34	37	---	---	51	45
		VR2	3	4	33	34	43	---	---	51	45
		R-Net	5	6	33	34	41	---	---	51	45
	10km/h	ACS/LINX /Shark	7	8	31	32	36	---	---	51	45
		VR2	9	10	31	32	40	---	---	51	45
		R-Net	11	12	31	32	42	---	---	51	45
BoraHD/ Spectra XTR2HD	6km/h	ACS/LINX /Shark	54	55	33	34	39	---	---	51	46
	10km/h	ACS/LINX /Shark	56	57	31	32	38	---	---	51	46
	6km/h	R-Net	58	59	33	34	62	---	---	51	46
	10km/h	R-Net	60	61	31	32	63	---	---	51	46
FDX	6km/h ^c	ACS/ Shark	13	14	33	34	37	52	53	---	45
	6km/h ^d		54	55	33	34	39	52	53	---	46
	10km/h		15	16	31	32	38	52	53	---	46
Kite	6km/h	ACS/LINX /Shark	17	18	33	34	37	---	---	51	45
	10km/h	ACS/LINX /Shark	19	20	31	32	36	---	---	51	45
Kite HD	6km/h	ACS/LINX /Shark	21	22	33	34	39	---	---	51	46
	10km/h	ACS/LINX /Shark	23	24	31	32	38	---	---	51	46
Storm4/ Storm4 X-plore	6km/h	ACS/ Shark	47	48	34	33	39	---	---	51	46
	10km/h		27	28	32	31	38	---	---	51	46
TDX-SP TDX-SP2	6km/h ^a	ACS/LINX /Shark	25	26	34	33	39	52	53	---	46
	8km/h ^a		29	30	34	33	38	52	53	---	46
	10km/h ^b		49	50	32	31	38	52	53	---	46
TDX-SP TDX-SP2 Narrow	6km/h	ACS/LINX /Shark	25	26	34	33	39	52	53	---	46
	10km/h		49	50	32	31	38	52	53	---	46

a) used from 28.09.11 / b) used from 07.09.11 / c) used till 31.05.13 / d) used from 01.06.13 (S/N 13F06000078)


-	N23004	KIL KRYSS DRIVAXEL (MO	1
1	SP1524195	DRIVENHET 6KMT VE BORA	1
2	SP1524194	DRIVENHET 6KMT HØ BORA	1
3	SP1541408	(Propulsion cpl. LH, Gearbox i = 23:1, Motor, 2-pole, 250W, slow, 600mm cable, VR2)	1
4	SP1541410	(Propulsion cpl. RH, Gearbox i = 23:1, Motor, 2-pole, 250W, slow, 600mm cable, VR2)	1
5	SP1541409	(Propulsion cpl. LH, Gearbox i = 23:1, Motor, 2-pole, 250W, slow, 600mm cable, R-Net)	1
6	SP1541411	(Propulsion cpl. RH, Gearbox i = 23:1, Motor, 2-pole, 250W, slow, 600mm cable, R-Net)	1
7	SP1524197	(Propulsion cpl. LH, Gearbox i = 18:1, Motor, 2-pole, 250W, fast, 600mm cable)	1
8	SP1524196	(Propulsion cpl. RH, Gearbox i = 18:1, Motor, 2-pole, 250W, fast, 600mm cable)	1
9	SP1541404	(Propulsion cpl. LH, Gearbox i = 18:1, Motor, 2-pole, 250W, fast, 600mm cable, VR2)	1
10	SP1541406	(Propulsion cpl. RH, Gearbox i = 18:1, Motor, 2-pole, 250W, fast, 600mm cable, VR2)	1
11	SP1541405	(Propulsion cpl. LH, Gearbox i = 18:1, Motor, 2-pole, 250W, fast, 600mm cable, R-Net)	1
12	SP1541407	(Propulsion cpl. RH, Gearbox i = 18:1, Motor, 2-pole, 250W, fast, 600mm cable, R-Net)	1
13	SP1168983	(Propulsion cpl. LH, Gearbox i = 23:1, Motor, 2-pole, 250W, slow, 600mm cable)	1
14	SP1168984	(Propulsion cpl. RH, Gearbox i = 23:1, Motor, 2-pole, 250W, slow, 600mm cable)	1
15	SP1168989	(Propulsion cpl. LH, Gearbox i = 18:1, Motor, 4-pole, 340W, fast, 760mm cable)	1
16	SP1168988	(Propulsion cpl. RH, Gearbox i = 18:1, Motor, 4-pole, 340W, fast, 760mm cable)	1
17	SP1568147	(Propulsion cpl. LH, Gearbox i = 23:1, Motor, 2-pole, 250W, slow, 600mm cable)	1
18	SP1568148	(Propulsion cpl. RH, Gearbox i = 23:1, Motor, 2-pole, 250W, slow, 600mm cable)	1

(Overview Propulsion SSD)

19	SP1568149	(Propulsion cpl. LH, Gearbox i = 18:1, Motor, 2-pole, 250W, fast, 600mm cable)	1
20	SP1568150	(Propulsion cpl. RH, Gearbox i = 18:1, Motor, 2-pole, 250W, fast, 600mm cable)	1
21	SP1568151	(Propulsion cpl. LH, Gearbox i = 23:1, Motor, 4-pole, 340W, slow, 760mm cable)	1
22	SP1568152	(Propulsion cpl. RH, Gearbox i = 23:1, Motor, 4-pole, 340W, slow, 760mm cable)	1
23	SP1568153	(Propulsion cpl. LH, Gearbox i = 18:1, Motor, 4-pole, 340W, fast, 760mm cable)	1
24	SP1568208	(Propulsion cpl. RH, Gearbox i = 18:1, Motor, 4-pole, 340W, fast, 760mm cable)	1
25	SP1602974	(Propulsion cpl. LH, Gearbox i = 23:1, Motor 4-pole 340W, slow, 760mm cable (incl. declutching lever))	1
26	SP1602975	(Propulsion cpl. RH, Gearbox i = 23:1, Motor 4-pole 340W, slow, 760mm cable (incl. declutching lever))	1
27	SP1541366	DRIVMOTOR KPL VE STORM4 SSD <i>S/N 061122610151 =></i>	1
28	SP1541365	DRIVMOTOR KPL HØ STORM4 SSD <i>S/N 061122610151 =></i>	1
29	SP1170307	DRIVMOTOR KPL VE TDX TDX SSD TDX SP2	1
30	SP1170306	DRIVMOTOR KPL HØ TDX TDX SSD TDX SP2	1
31	SP1163413	(Gearbox, LH, i = 18:1)	1
32	SP1163412	(Gearbox, RH, i = 18:1)	1
33	SP1163415	(Gearbox, LH, i = 23:1)	1
34	SP1163414	(Gearbox, RH, i = 23:1)	1
35	SP1542953	(Clamp band with clutch)	1
36	SP1163099	(Motor, 2-pole, 250W, fast, 600mm cable, ACS)	1
37	SP1164046	(Motor, 2-pole, 250W, slow, 600mm cable, ACS)	1
38	SP1165863	(Motor, 4-pole, 340W, fast, 760mm cable, ACS)	1
39	SP1165865	(Motor, 4-pole, 340W, slow, 760mm cable, ACS)	1
40	SP1544840	(Motor, 2-pole, 250W, fast, 600mm cable, VR2)	1
41	SP1544796	(Motor, 2-pole, 250W, slow, 600mm cable, R-Net)	1
42	SP1544841	(Motor, 2-pole, 250W, fast, 600mm cable, R-Net)	1
43	SP1544795	(Motor, 2-pole, 250W, slow, 600mm cable, VR2)	1
45	SP1540358	(Kit, brushes with caps, 2-pole (Set of 2))	1
46	SP1522967	BØRSTER 4 STK KIT TDX STORM 4	1
47	SP1536022	(Propulsion cpl. LH, Gearbox i = 23:1, Motor, 4-pole, 340W, slow, 760mm cable (incl. declutching knob)) <i>S/N 061122812868 =></i>	1
48	SP1536021	(Propulsion cpl. RH, Gearbox i = 23:1, Motor, 4-pole, 340W, slow, 760mm cable (incl. declutching knob)) <i>S/N 061122812868 =></i>	1
49	SP1602976	(Propulsion cpl. LH, Gearbox i = 18:1, Motor 4-pole 340W, fast, 760mm cable (incl. declutching lever))	1
50	SP1602998	(Propulsion cpl. RH, Gearbox i = 18:1, Motor 4-pole 340W, fast, 760mm cable (incl. declutching lever))	1
51	SP1579595	FRIKOPLINGS HÅNDTAK STORM4 OG ALLE SSD MOTOR <i>Single part. Consider to buy another one for the other side if required. - Replaces SP1168668, SP1168669 & SP1169005</i>	1
52	SP1168667	(Declutching lever, cpl. - LH)	1
53	SP1168666	(Declutching lever, cpl. - RH)	1
54	SP1552782	(Propulsion cpl. LH, Gearbox i = 23:1, Motor 4-pole, 340W, slow, 760mm cable, ACS)	1
55	SP1552783	(Propulsion cpl. RH, Gearbox i = 23:1, Motor 4-pole, 340W, slow, 760mm cable, ACS)	1
56	SP1552872	(Propulsion cpl. LH, Gearbox i = 18:1, Motor 4-pole, 340W, fast, 760mm cable, ACS)	1
57	SP1552873	(Propulsion cpl. RH, Gearbox i = 18:1, Motor 4-pole, 340W, fast, 760mm cable, ACS)	1

(Overview Propulsion SSD)


58	SP1552832	(Propulsion cpl. LH, Gearbox i = 23:1, Motor 4-pole, 340W, slow, 760mm cable, R-Net)	1
59	SP1552833	(Propulsion cpl. RH, Gearbox i = 23:1, Motor 4-pole, 340W, slow, 760mm cable, R-Net)	1
60	SP1552836	(Propulsion cpl. LH, Gearbox i = 18:1, Motor 4-pole, 340W, fast, 760mm cable, R-Net)	1
61	SP1552837	(Propulsion cpl. RH, Gearbox i = 18:1, Motor 4-pole, 340W, fast, 760mm cable, R-Net)	1
62	SP1562498	(Motor, 4-Pol, 340W, slow, 760mm cable, R-Net)	1
63	SP1562499	(Motor, 4-Pol, 340W, fast, 760mm cabel, R-Net)	1

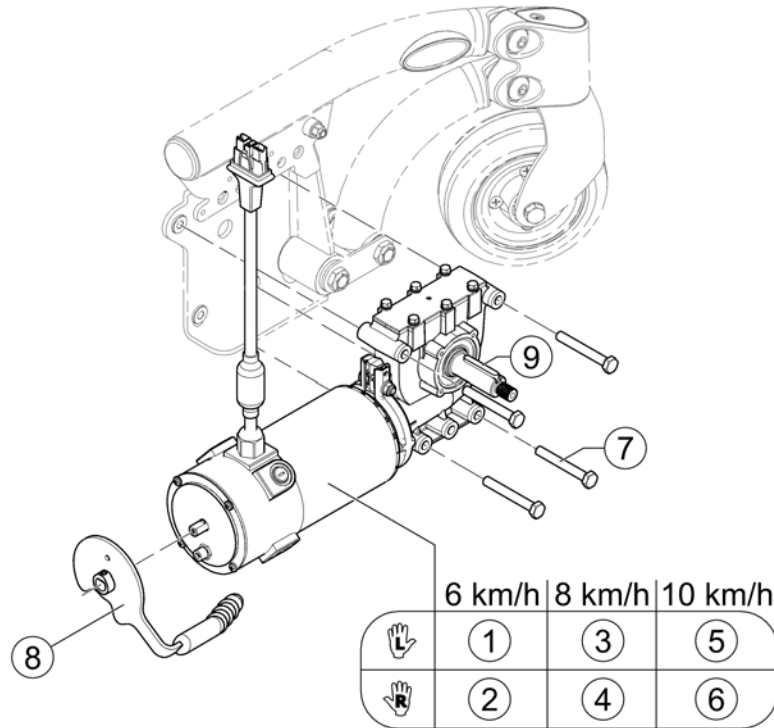


1	SP1602974	(Propulsion cpl. LH, Gearbox i = 23:1, Motor 4-pole 340W, slow, 760mm cable (incl. declutching lever))	1
2	SP1602975	(Propulsion cpl. RH, Gearbox i = 23:1, Motor 4-pole 340W, slow, 760mm cable (incl. declutching lever))	1
3	SP1170307	DRIVMOTOR KPL VE TDX TDX SSD TDX SP2	1
4	SP1170306	DRIVMOTOR KPL HØ TDX TDX SSD TDX SP2	1
5	SP1602976	(Propulsion cpl. LH, Gearbox i = 18:1, Motor 4-pole 340W, fast, 760mm cable (incl. declutching lever))	1
6	SP1602998	(Propulsion cpl. RH, Gearbox i = 18:1, Motor 4-pole 340W, fast, 760mm cable (incl. declutching lever))	1
7	D00100	(Hexagon bolt DIN931-M08x60-8.8 bk)	4
8	SP1168667	(Declutching lever, cpl. - LH)	1
8	SP1168666	(Declutching lever, cpl. - RH)	1
9	N23004	KIL KRYSS DRIVAXEL (MO)	1

(Propulsion SSD - Narrow Base)

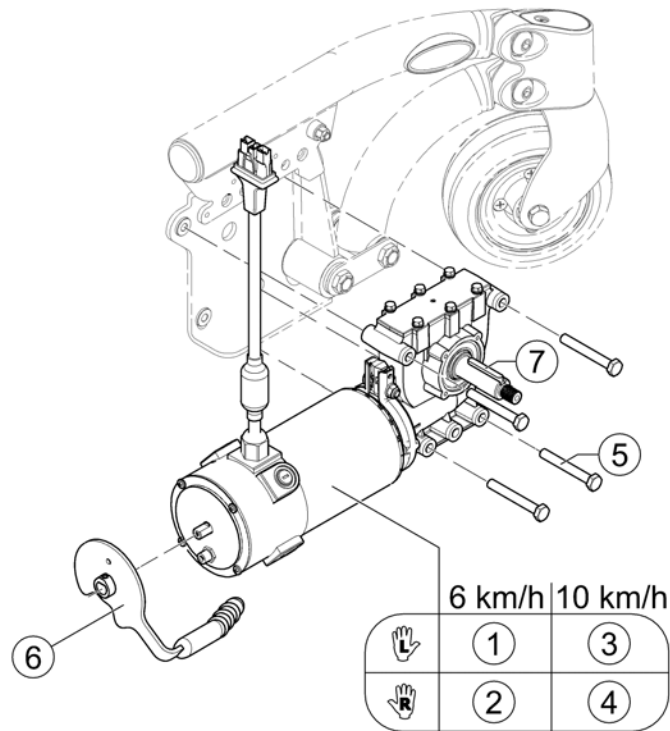
1563700.no







1	SP1602974	(Propulsion cpl. LH, Gearbox i = 23:1, Motor 4-pole 340W, slow, 760mm cable (incl. declutching lever))	1
2	SP1602975	(Propulsion cpl. RH, Gearbox i = 23:1, Motor 4-pole 340W, slow, 760mm cable (incl. declutching lever))	1
3	SP1602976	(Propulsion cpl. LH, Gearbox i = 18:1, Motor 4-pole 340W, fast, 760mm cable (incl. declutching lever))	1
4	SP1602998	(Propulsion cpl. RH, Gearbox i = 18:1, Motor 4-pole 340W, fast, 760mm cable (incl. declutching lever))	1
5	D00100	(Hexagon bolt DIN931-M08x60-8.8 bk)	4
6	SP1168667	(Declutching lever, cpl. - LH)	1
6	SP1168666	(Declutching lever, cpl. - RH)	1
7	N23004	KIL KRYSS DRIVAXEL (MO)	1

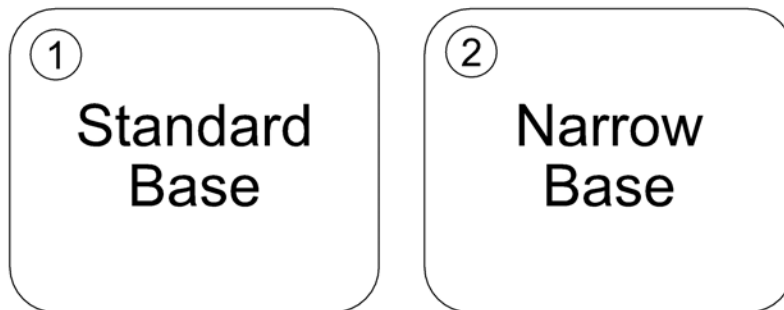

1	-	(Electronic - Standard Base)	1
2	-	(Electronic - Narrow Base)	1

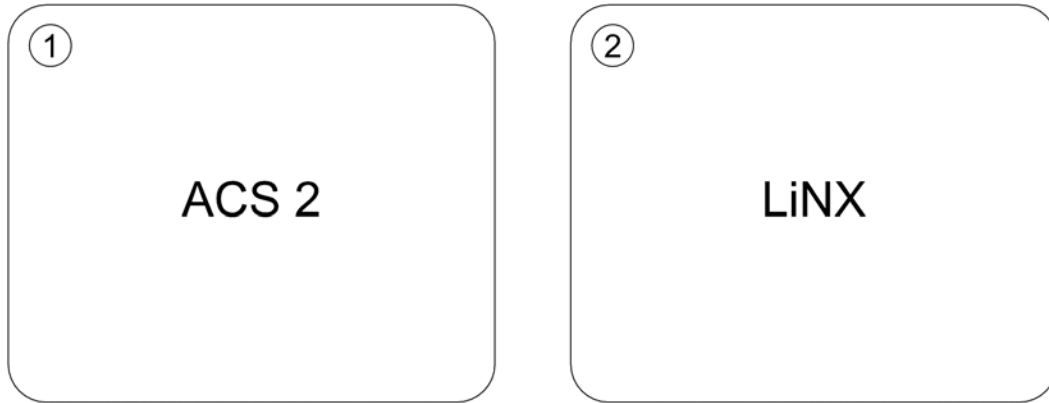

1	-	(ACS 2 Electronic)	1
2	-	(LiNX Electronic)	1

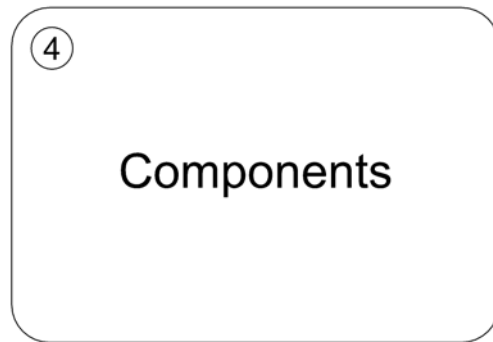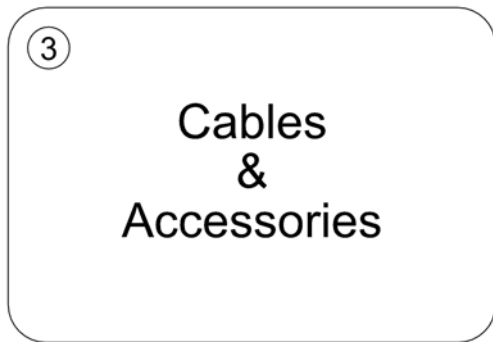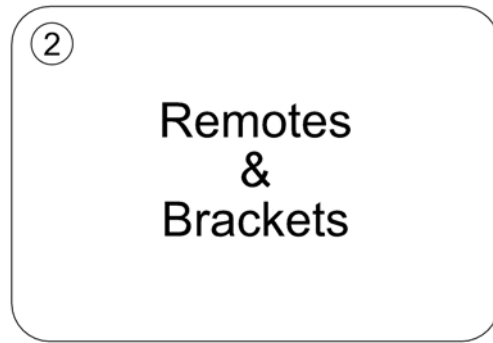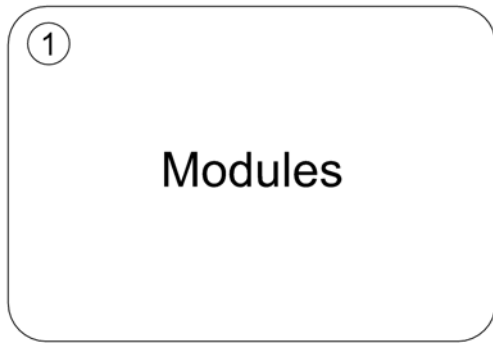


1	-	(ACS Electronic Modules)	1
2	-	(ACS Remotes & Brackets)	1
3	-	(ACS Cables & Accessories)	1
4	-	(ACS Components)	1

(ACS Electronic Modules)

-	(Power Modul - ACS)	1
-	(Bracket for Power Modul - Standard Base)	1
-	(ACTs)	1
-	(Bracket ACT - Modulite Seat)	1
-	(Bracket ACT - Ultra Low Maxx Seat)	1
-	(Bracket ACT - Max Seat)	1
-	(Light Board)	1
-	(G-Trac Sensor/Gyro)	1

(Power Modul - ACS)

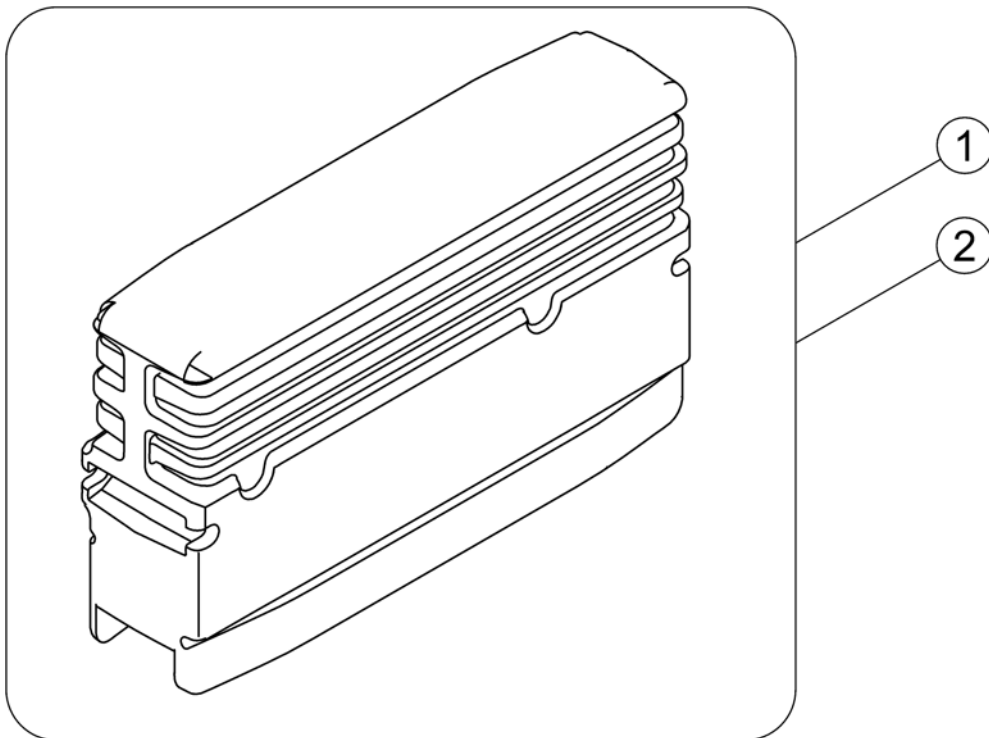
1583335.iso



1	SP1524244	POWERMODUL DX2 PMA90L STORM4/CHASSWHEEL <i>Lock washer between module and bracket!</i>	1
2	SP1537975	POWERMODUL ACS2 FOR GYRO ACS2 <i>Lock washer between module and bracket!</i>	1

(Bracket for Power Modul - Standard Base)

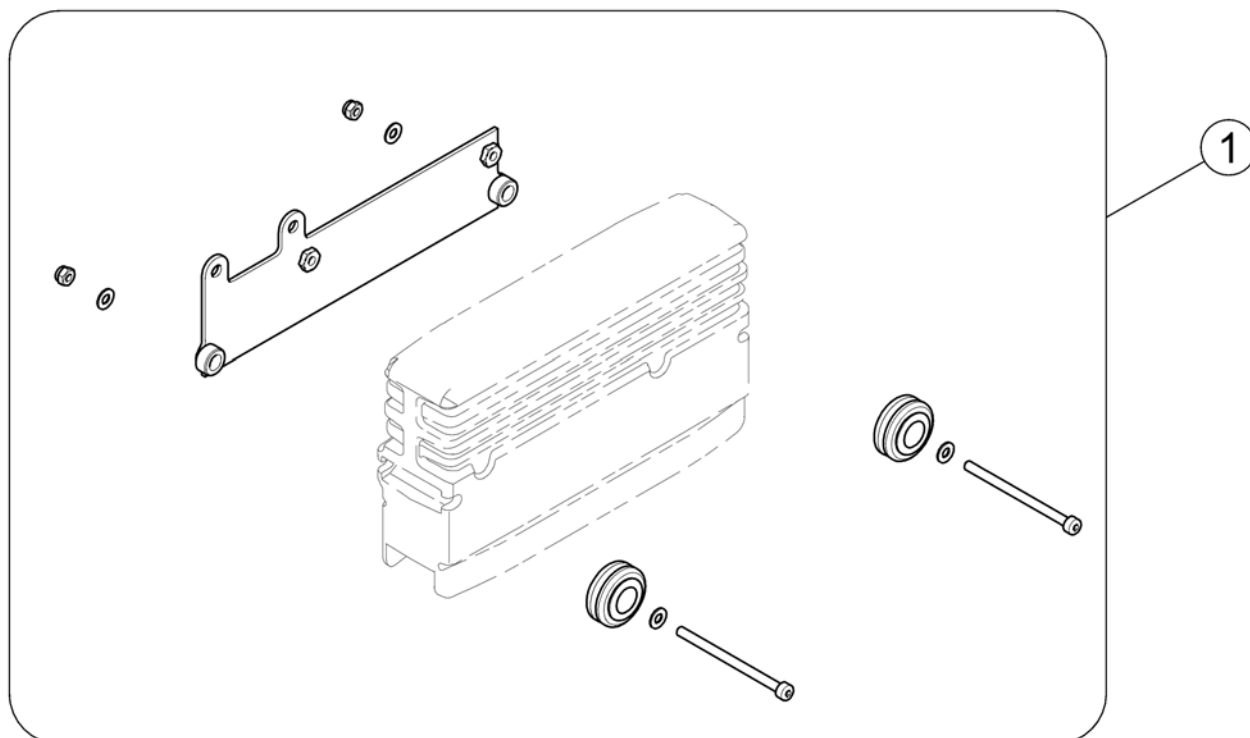
1583336.iso



1 SP1583340

(Bracket Power Module)

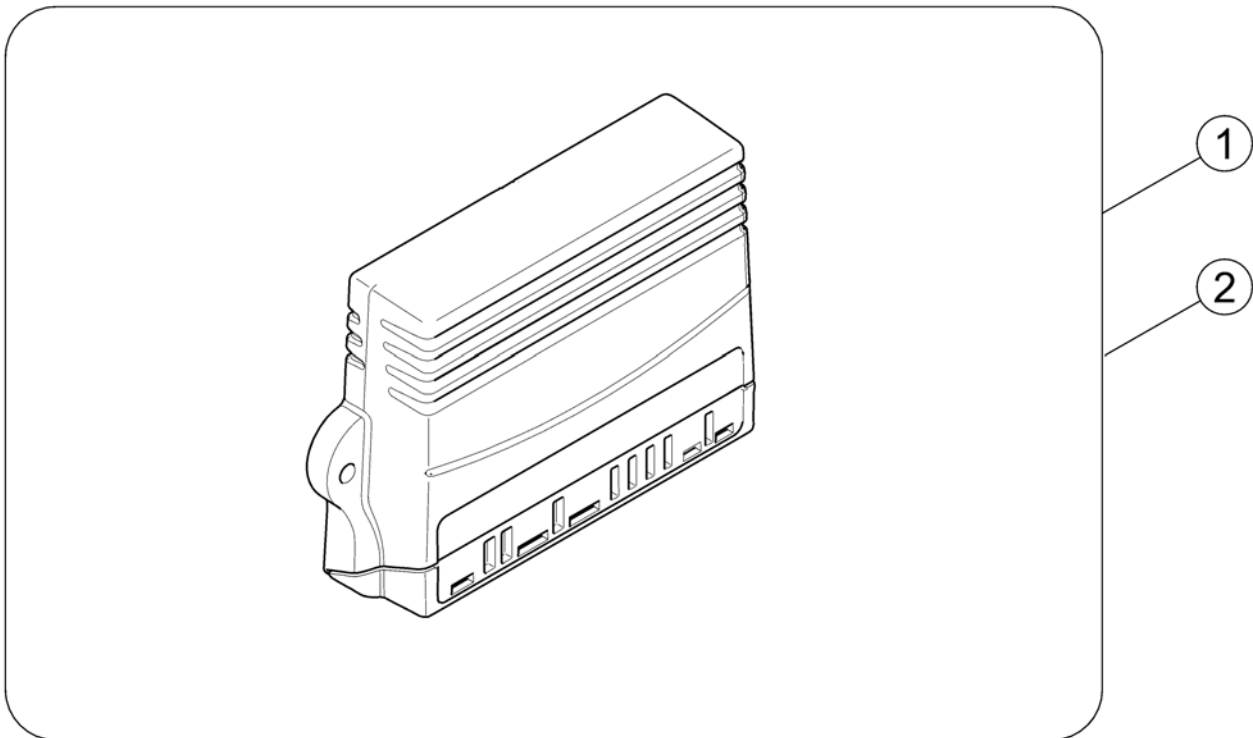
1




1	SP1528542	ACTUATORMODUL ACT2 DX2 STORM4 TDX BORA	1
2	SP1529382	ACTUATORMODUL ACT4 DX2 STORM4 TDX BORA	1

(Bracket ACT - Modulite Seat)

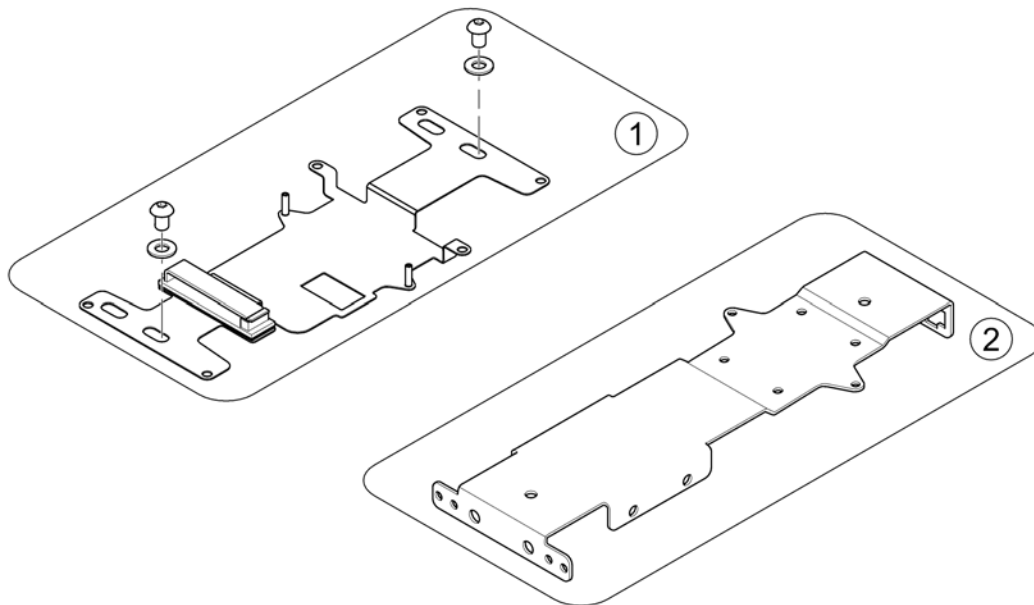1583343.no




1	SP1569960	(ACT bracket)	1
2	SP1564650	(Multi support bracket)	1

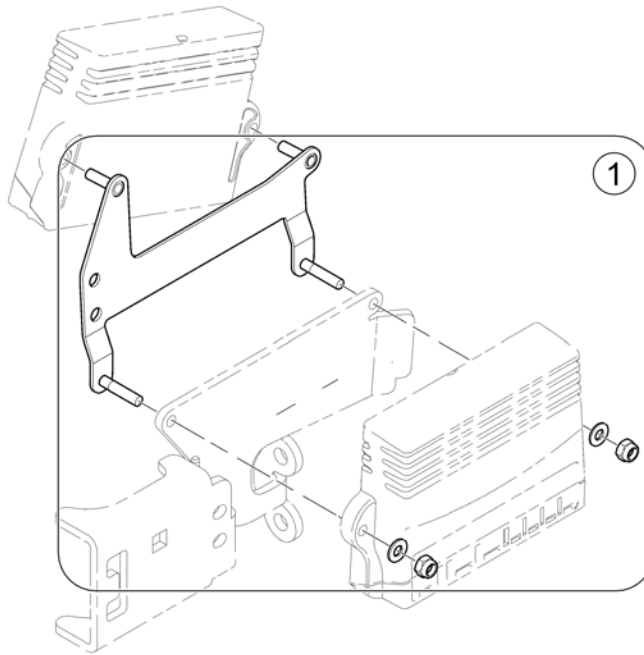

1 SP1584343

(Bracket ACTs (Ultra Low Maxx Seat))

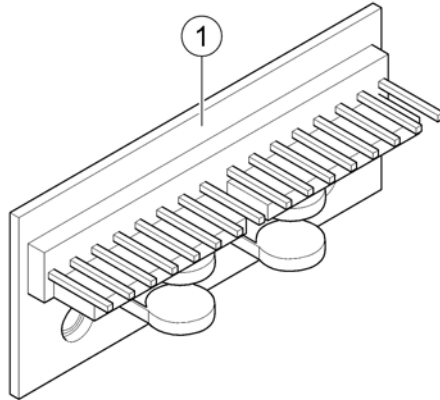1




1 1415004

LYSFORDELER PRINT G40

1

(G-Trac Sensor/Gyro)

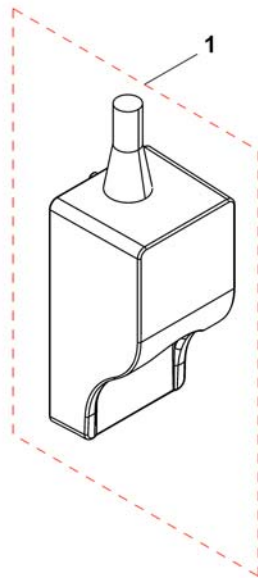1537848X.dwg




1 SP1541535

(G-Trac (Gyro) Sensor, cpl.)

1

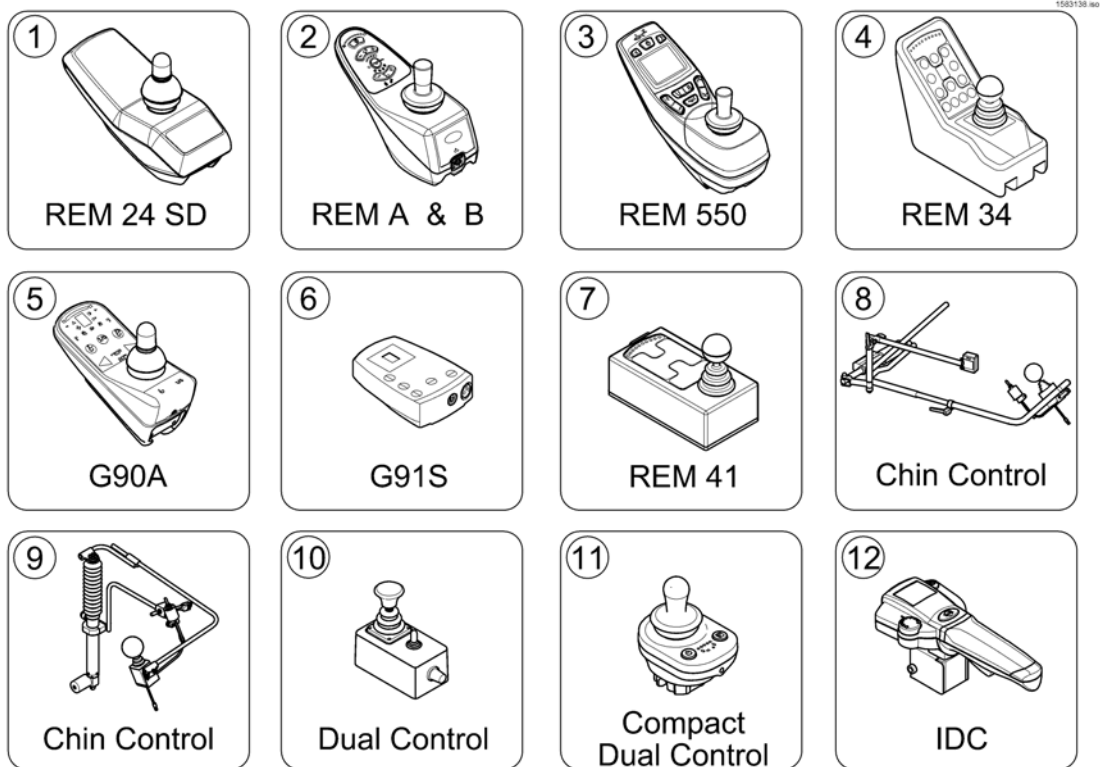
(ACS Remotes & Brackets)



1	-	(REM 24 SD with brackets) [1601/1602]		1
2	-	(REM A & B with brackets) [1638/1639/1643/1644]		1
3	-	(REM 550 with brackets) [1617/1618]		1
4	-	(REM 34 with brackets) [1611]		1
5	-	(G90A with brackets) [1614/1616]		1
6	-	(G91S with brackets) [1615]		1
7	-	(Centre control REM 41) [1607/1608]		1
8	-	(Chin Control, small (=> 16/12/31)) [1609]		1
9	-	(Chin Control - Modulite Seat (15/10/01 =>)) [1719]		1
10	-	(Dual Control with brackets) [1291/1292]	→ 30.06.2015	1
11	-	(Compact Dual Control with brackets) [1720]		1
12	-	(Intuitive Dual Control (IDC) with brackets) [1645]		1

(REM 24 SD with brackets)

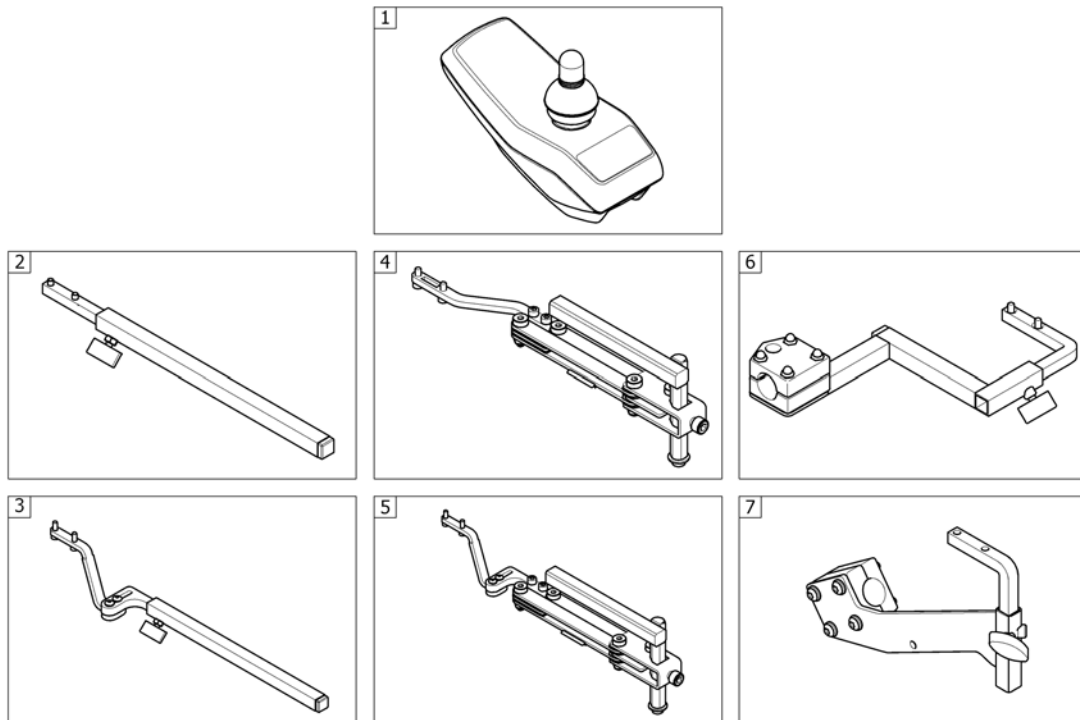
158433 no



1	-	(Remote REM 24SD) [1601/1602]	0
2	-	(Remote Bracket Standard) [1601+1207/1208]	1
3	-	(Remote Bracket Standard - ULM and Soft Armpad on Modulite) [1601+1207/1208 +2020 or +2100+2199]	1
4	-	(Remote Bracket Swing away & height adjust.) [1601+1707/1708]	1
5	-	(Remote Bracket Swing away & height adjust. - ULM and Soft Armpad on Modulite) [1601+1707/1708 +2020 or +2100+2199]	1
6	-	(Remote Bracket Attendant) [1602+1207/1208]	1
7	-	(Remote Bracket Attendant Ultra Low Maxx Seat) [1602+1207/1208+2020]	1

(Remote REM 24SD)

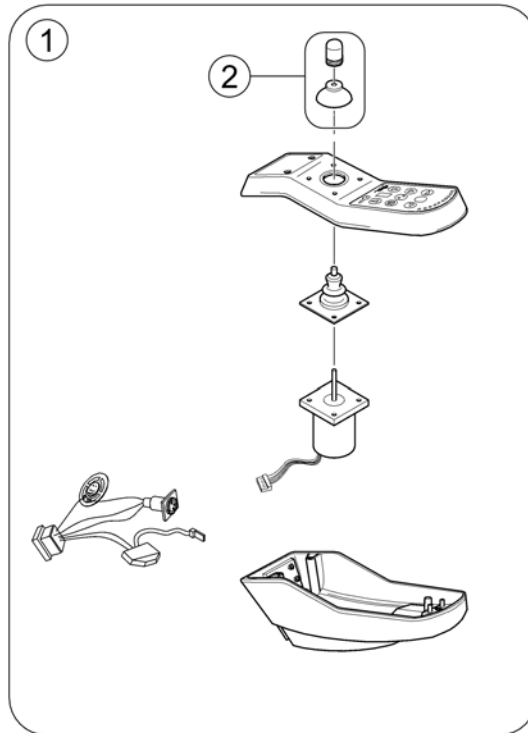
142802 860



1	1436541	(Remote REM24SD)	1
2	SPZ183229	JOYSTIKK KOMPLETT	1

(Remote Bracket Standard)

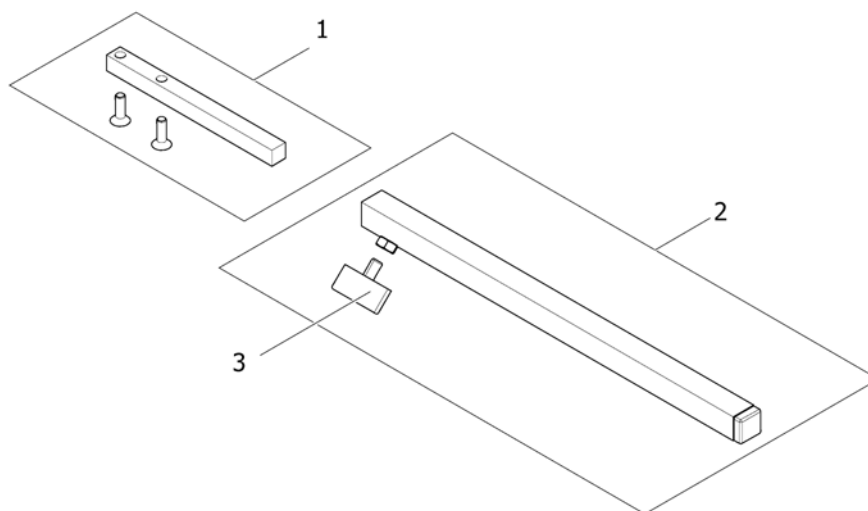
1559851.iso



1	SP1436542	HOLDER DX-REM 24 SD RETT STORM SPORT/3 OPPFL.ARML	1
2	SPS183040	STYREBOKSHOLDER STD KPL STORM3 G40	1
3	L31001	VINGSKRUV M 6X10 GARANT G24 S	1

(Remote Bracket Standard - ULM and Soft Armpad on Modulite)

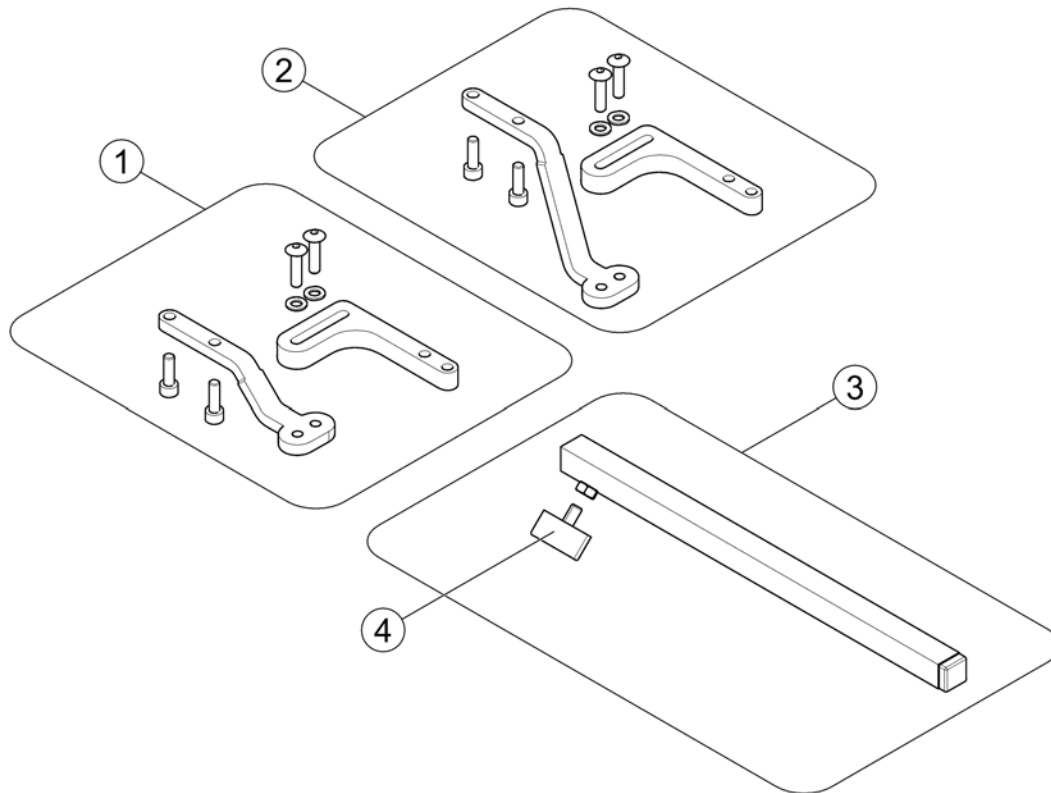
1588513.no



1	SP1587309	(Remote bracket Std. REM 24SD/A/B550)	1
2	SP1587693	(Remote bracket MACES Std - REM 24SD/A/B/550)	1
3	SPS183040	STYREBOKSHOLDER STD KPL STORM3 G40	1
4	L31001	VINGSKRUV M 6X10 GARANT G24 S	1

(Remote Bracket Swing away & height adjust.)

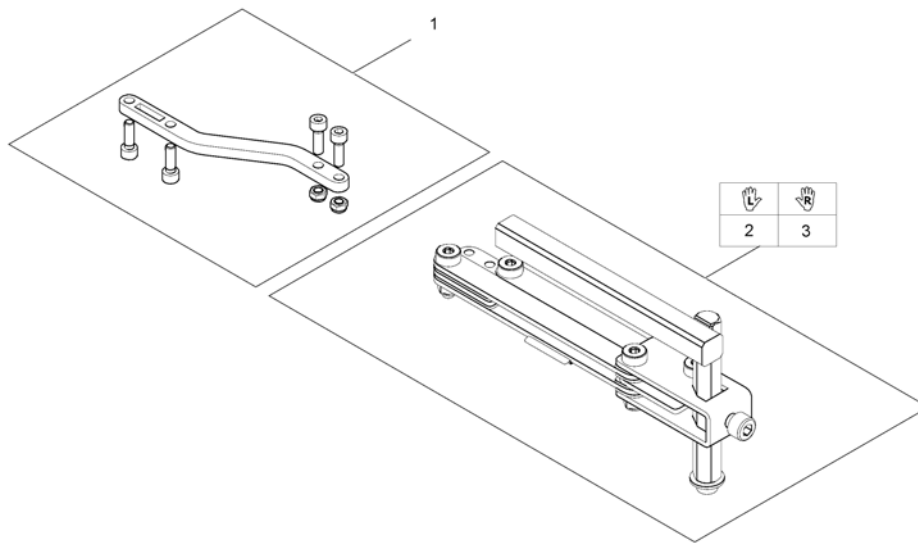
1559852.iso



1	SP1498222	FESTE FOR STYREBOKST TIL SWING AWAY	1
2	SP1496249	SWING AWAY VE NY MOD2008 KPL	1
3	SP1496248	SWING AWAY HØ NY MOD2008 KPL	1

(Remote Bracket Swing away & height adjust. - ULM and Soft Armpad on Modulite)

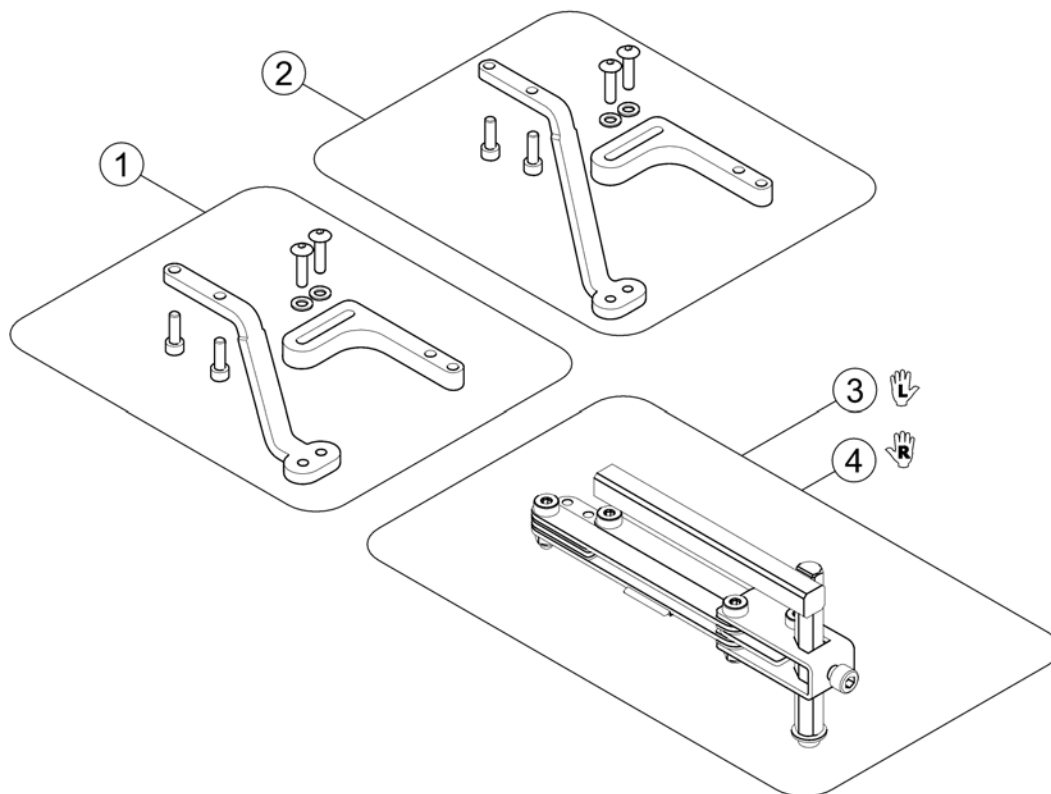
1588514.sop



1	SP1587329	(Remote bracket swingaway REM 24SD/A/B/550)	1
2	SP1587695	(Remote bracket MACES swingaway - REM 24SD/A/B/550)	1
3	SP1496249	SWING AWAY VE NY MOD2008 KPL	1
4	SP1496248	SWING AWAY HØ NY MOD2008 KPL	1

(Remote Bracket Attendant)

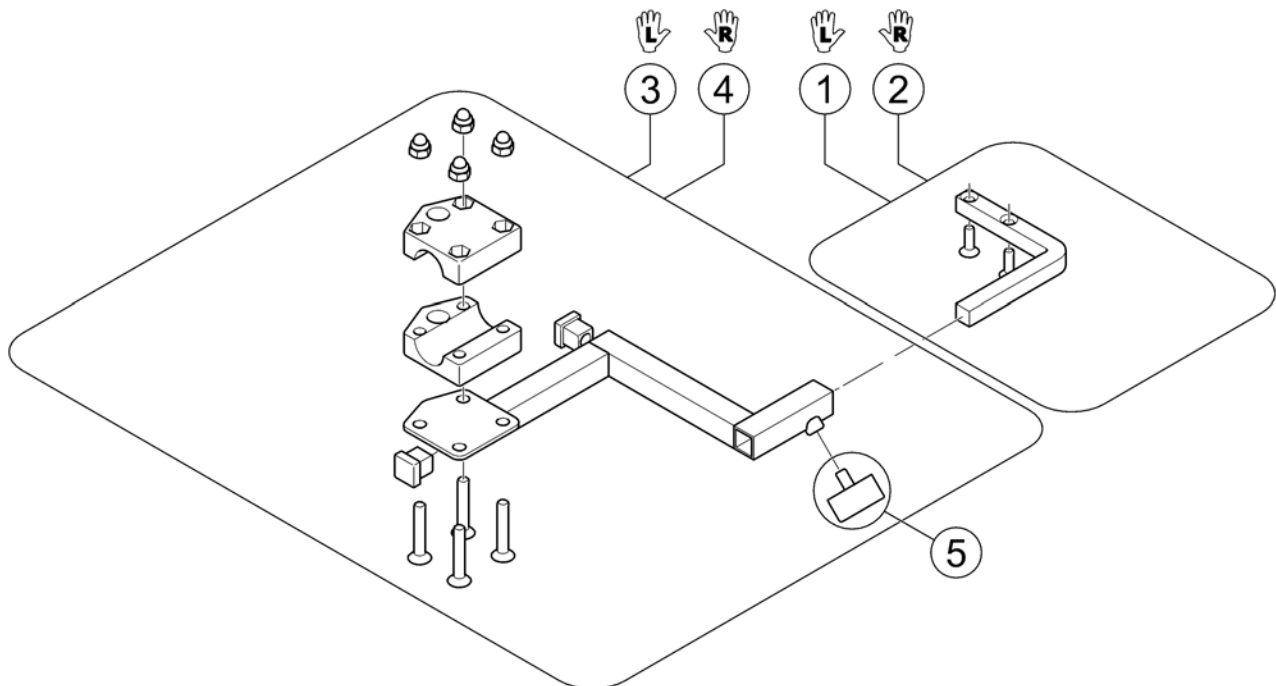
1560837.no



1	1422806	(Remote Support LH)	1
2	1422810	(Remote Support RH)	1
3	SP1422937	(Bracket for attendant control LH)	1
4	SP1422936	(Bracket for attendant control RH)	1
5	L31001	VINGSKRUV M 6X10 GARANT G24 S	1

(Remote Bracket Attendant Ultra Low Maxx Seat)

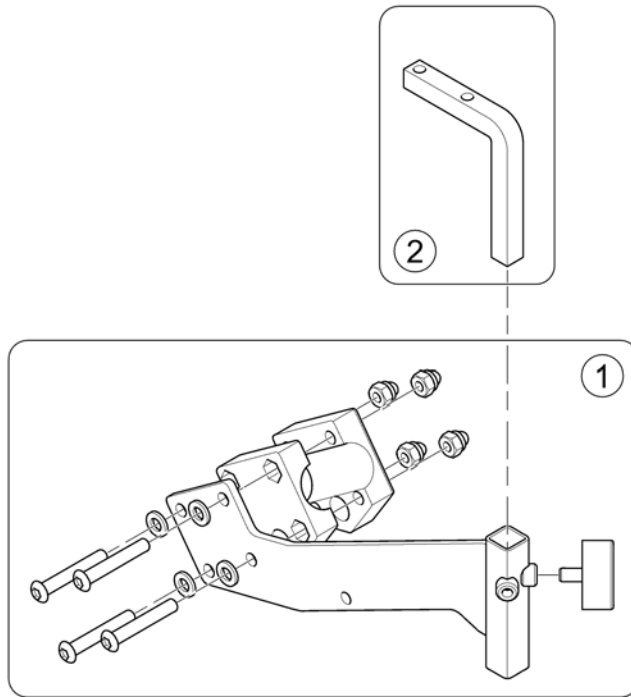
1588437 no



1	SP1588427	BRAKETT LEDSAGERSTYRING STYREBOKS BAK	1
2	SP1587228	HOLDER STYREBOKS LEDS REM 24/A/B/550 ULM	1

(REM A & B with brackets)

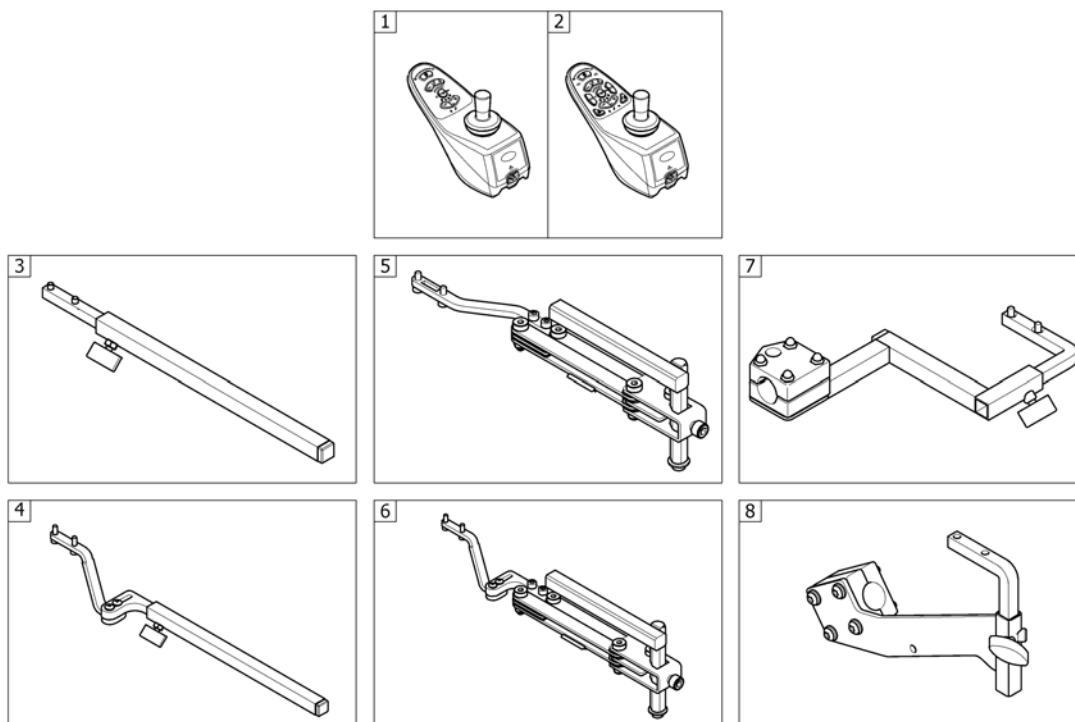
158434 No



1	SP1523440	(Remote REM A) <i>[1638/1643]</i>	1
2	SP1523445	(Remote REM B) <i>[1639/1644]</i>	1
3	-	(Remote Bracket Standard) <i>[1638/1639+1207/1208]</i>	1
4	-	(Remote Bracket Standard - ULM and Soft Armpad on Modulite) <i>[1638/1639+1207/1208 +2020 or +2100+2199]</i>	1
5	-	(Remote Bracket Swing away & height adjust.) <i>[1638/1639+1707/1708]</i>	1
6	-	(Remote Bracket Swing away & height adjust. - ULM and Soft Armpad on Modulite) <i>[1638/1639+1707/1708 +2020 or +2100+2199]</i>	1
7	-	(Remote Bracket Attendant) <i>[1643/1644+1207/1208]</i>	1
8	-	(Remote Bracket Attendant Ultra Low Maxx Seat) <i>[1643/1644+1207/1208+2020]</i>	1

(Remote Bracket Swing away & height adjust.)

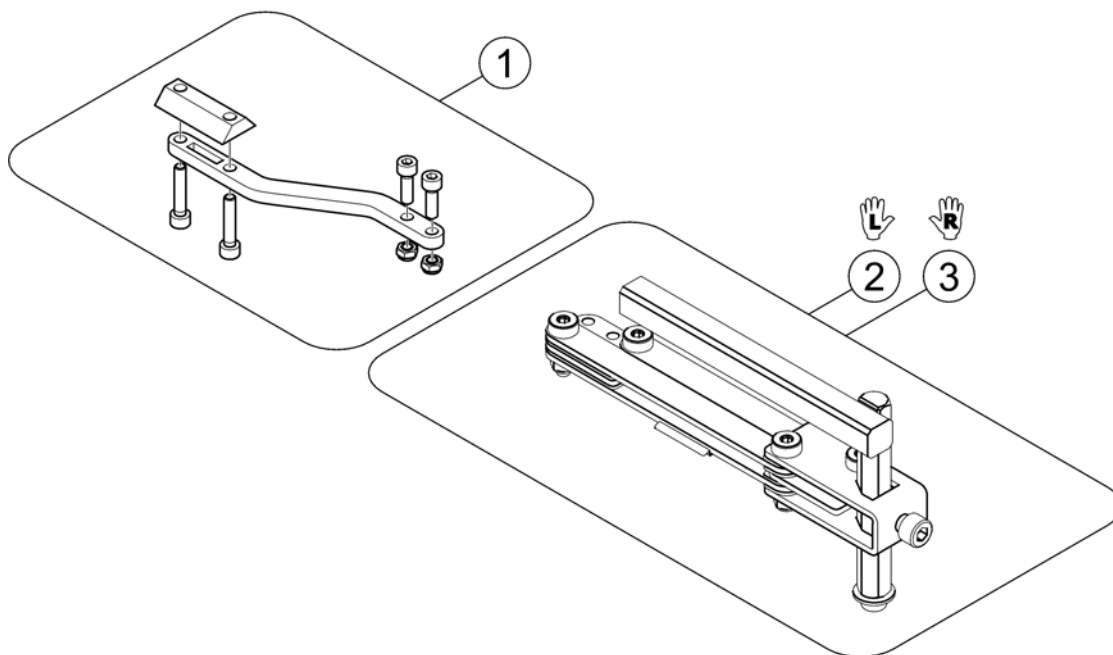
1559857.iso



1	SP1536180	HOLDER STYREBOKS SHARK/REM A/REM B	1
2	SP1496249	SWING AWAY VE NY MOD2008 KPL	1
3	SP1496248	SWING AWAY HØ NY MOD2008 KPL	1


(REM 550 with brackets)

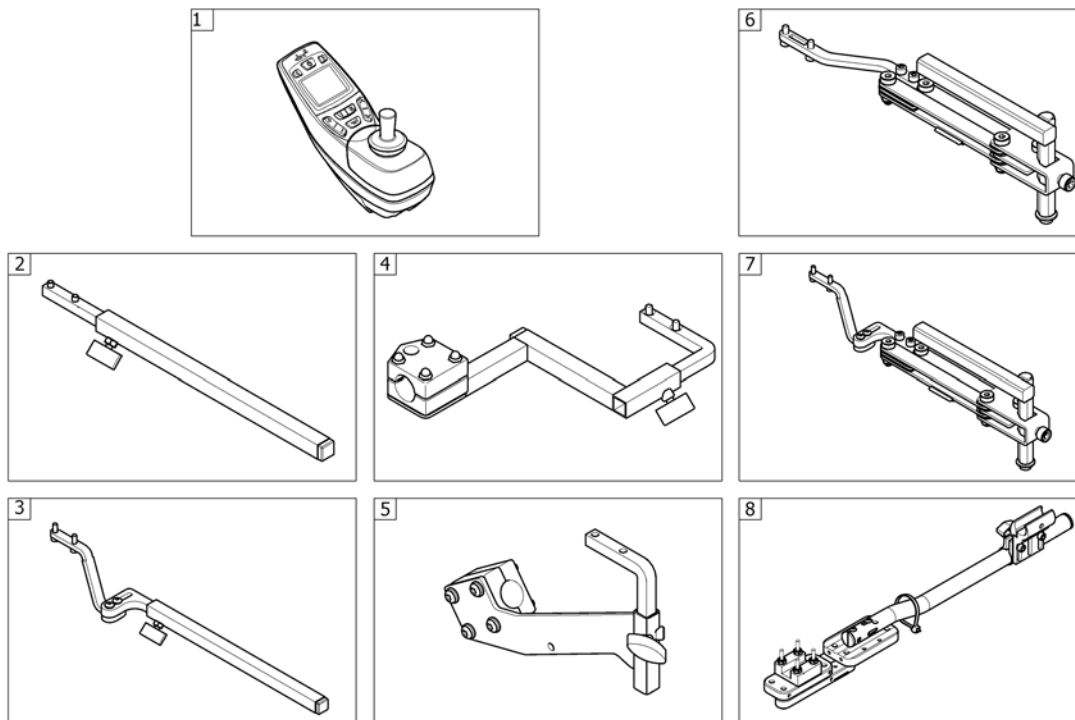1588435 Rev



1	SP1523444	STYREBOKS REM550 STORM4 TDX BORA <i>[1617/1618]</i>	1
2	-	(Remote Bracket Standard) <i>[1617+1207/1208]</i>	1
3	-	(Remote Bracket Standard - ULM and Soft Armpad on Modulite) <i>[1617+1207/1208 +2020 or +2100+2199]</i>	1
4	-	(Remote Bracket Attendant) <i>[1618+1207/1208]</i>	1
5	-	(Remote Bracket Attendant Ultra Low Maxx Seat) <i>[1618+1207/1208+2020]</i>	1
6	-	(Remote Bracket Swing away & height adjust.) <i>[1617+1707/1708]</i>	1
7	-	(Remote Bracket Swing away & height adjust. - ULM and Soft Armpad on Modulite) <i>[1617+1707/1708 +2020 or +2100+2199]</i>	1
8	-	(Remote Bracket Swing away Quad Link Ultra Low Maxx Seat) <i>[1617+2075/2076+2020]</i>	1


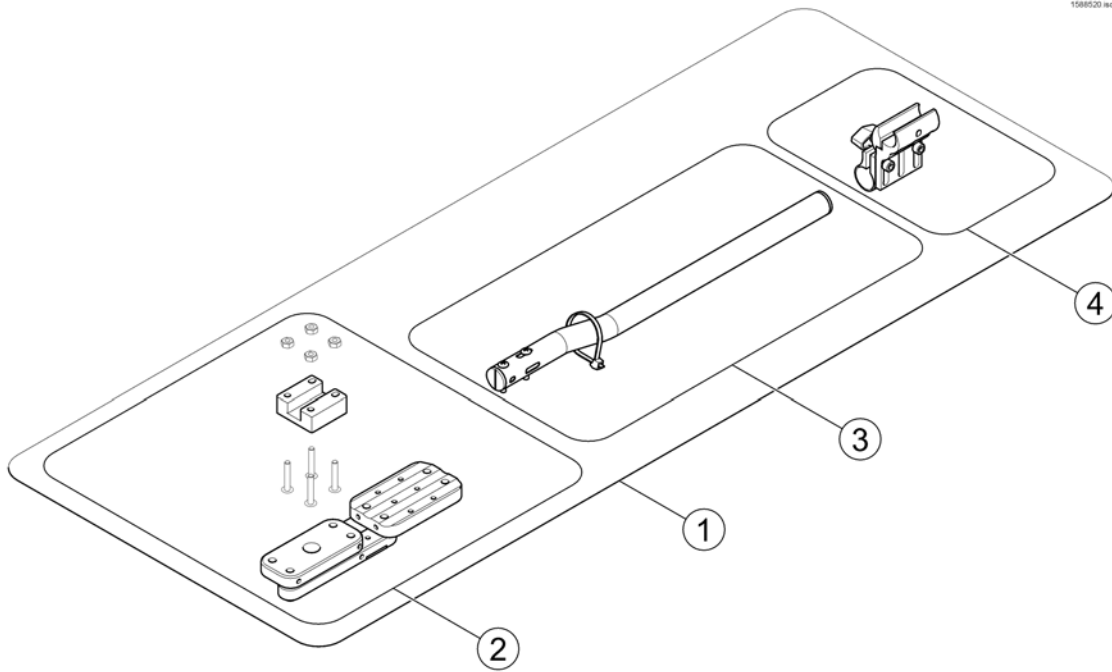
(Remote Bracket Swing away Quad Link Ultra Low Maxx Seat)




1	SP1584346	(Quad Link, cpl.)	1
2	SP1140896	(Quad Link)	1
3	SP1139416	(Quad Link bracket)	1
4	SPTRC2089	(Joystick Clamp Assy Quad Link)	1

(REM 34 with brackets)

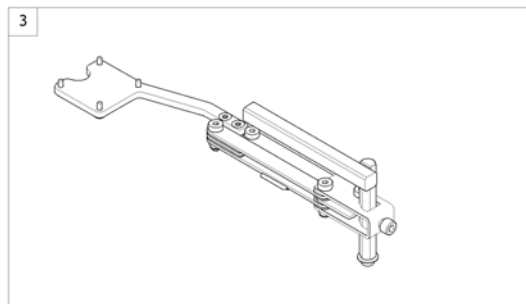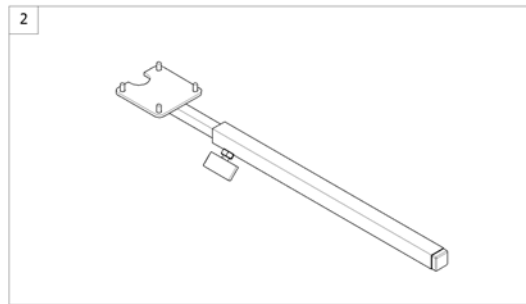1560846 iso




1	SP1440677	(Remote REM 34) [1611]	1
2	-	(Remote Bracket Standard) [1611+1207/1208]	1
3	-	(Remote Bracket Swing away & height adjust.) [1611+1707/1708]	1

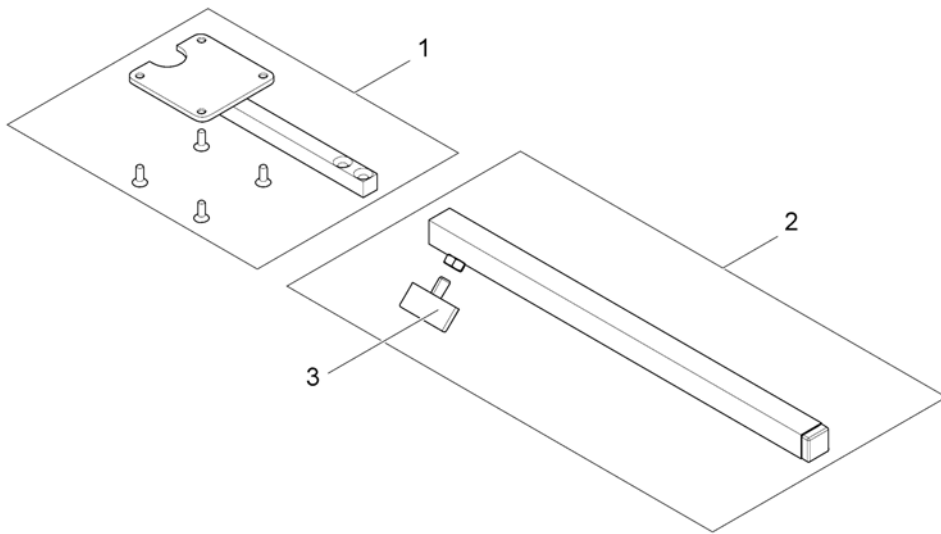


1	SP1434921	HOLDER FOR REM 34 STORM4	1
2	SPS183040	STYREBOKSHOLDER STD KPL STORM3 G40	1
3	L31001	VINGSKRUV M 6X10 GARANT G24 S	1

(Remote Bracket Swing away & height adjust.)

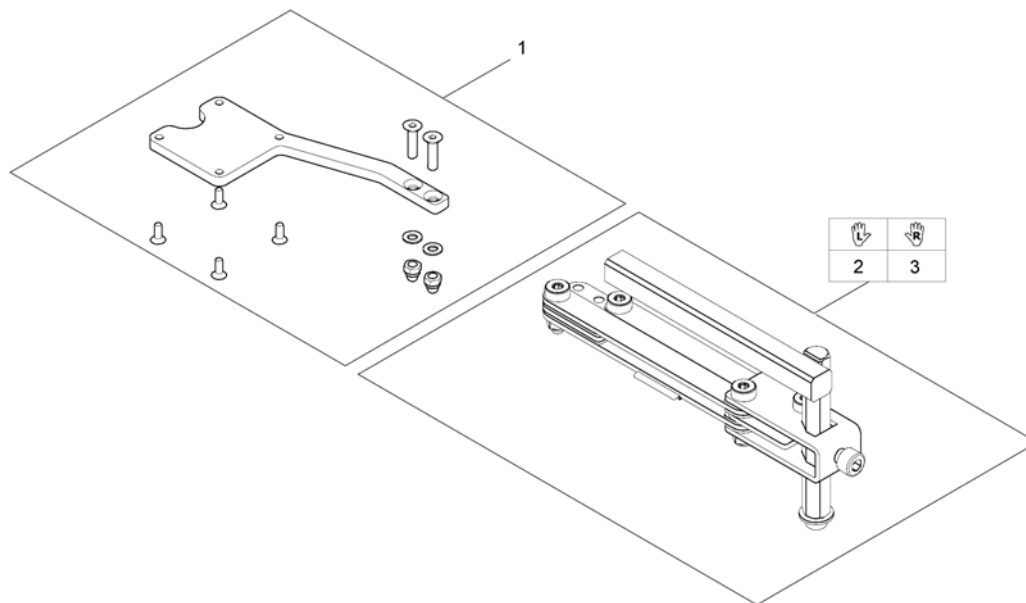1560845.iso




1	SP1527125	FESTE FOR REM34 STORM4	1
2	SP1496249	SWING AWAY VE NY MOD2008 KPL	1
3	SP1496248	SWING AWAY HØ NY MOD2008 KPL	1

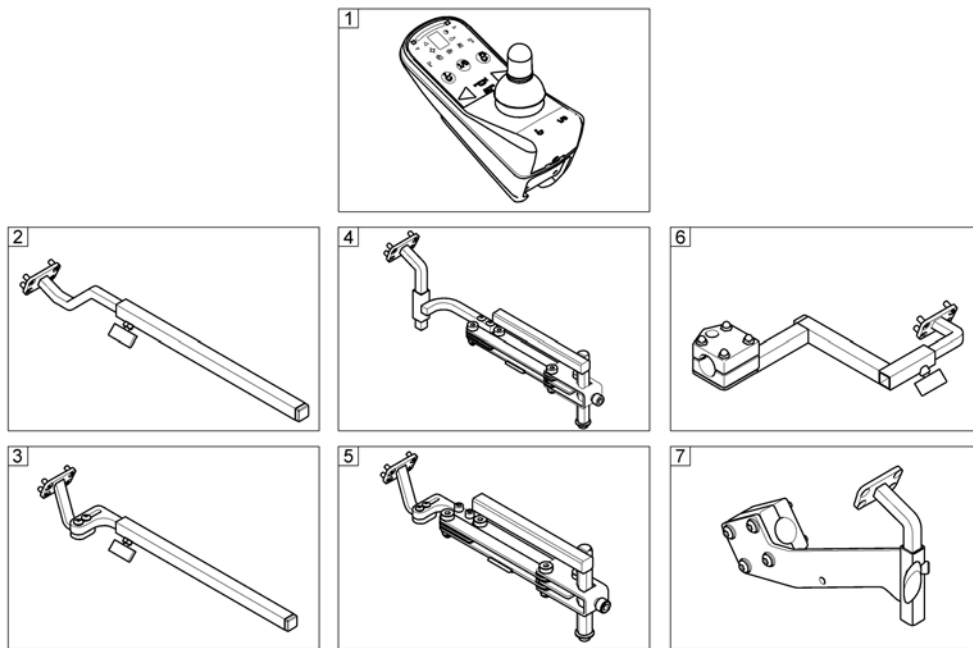

1	SP1504720	STYREBOKS DX-REM G90A STORM4 [1614/1616]	1
2	-	(Remote Bracket Standard) [1614+1207/1208]	1
3	-	(Remote Bracket Standard - ULM and Soft Armpad on Modulite) [1614+1207/1208 +2020 or +2100+2199]	1
4	-	(Remote Bracket Swing away & height adjust.) [1614+1707/1708]	1
5	-	(Remote Bracket Swing away & height adjust. - ULM and soft Armpad on Modulite) [1614+1707/1708 +2020 or +2100+2199]	1
6	-	(Remote Bracket Attendant) [1616+1207/1208]	1
7	-	(Remote Bracket Attendant Ultra Low Maxx Seat) [1616+1207/1208+2020]	1

(Remote Bracket Standard)

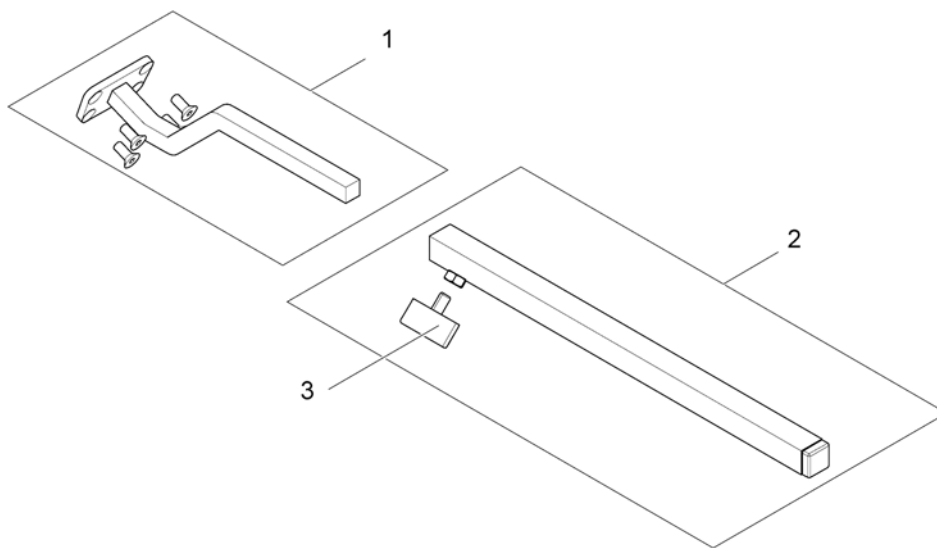
1562939.iso



1	SP1504822	(Remote Support Standard)	1
2	SPS183040	STYREBOKSHOLDER STD KPL STORM3 G40	1
3	L31001	VINGSKRUV M 6X10 GARANT G24 S	1

(Remote Bracket Standard - ULM and Soft Armpad on Modulite)

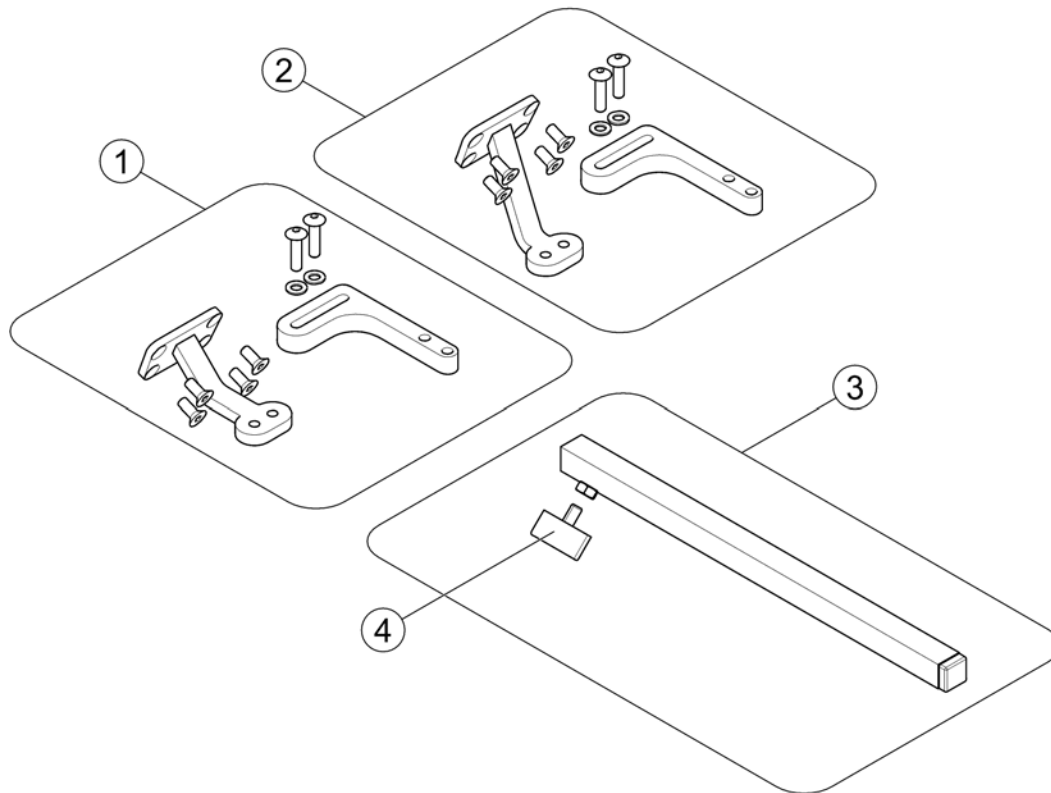
1588518.no



1	SP1587310	(Remote bracket Std G90)	1
2	SP1587694	(Remote bracket MACES Std - G90)	1
3	SPS183040	STYREBOKSHOLDER STD KPL STORM3 G40	1
4	L31001	VINGSKRUV M 6X10 GARANT G24 S	1

(Remote Bracket Swing away & height adjust.)

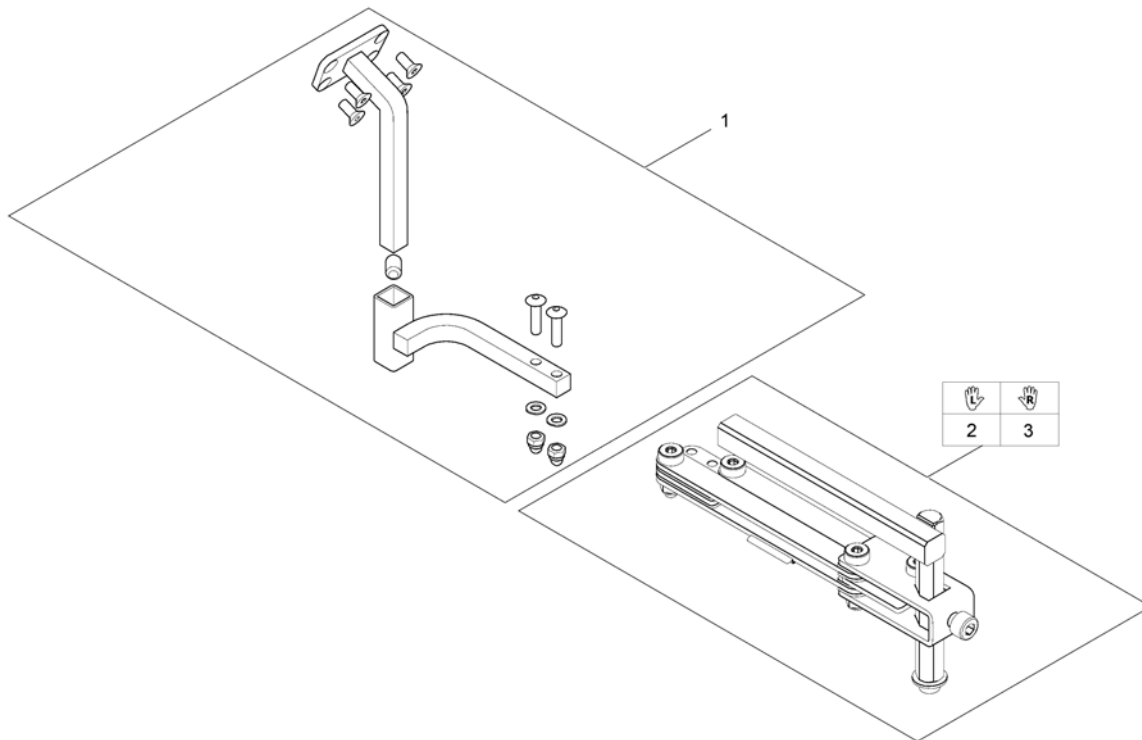
1560310.iso



1	SP1517352	FESTE FOR G90 SWINGAWAY STORM4 TDX BORA	1
2	SP1496249	SWING AWAY VE NY MOD2008 KPL	1
3	SP1496248	SWING AWAY HØ NY MOD2008 KPL	1

(Remote Bracket Swing away & height adjust. - ULM and soft Armpad on Modulite)

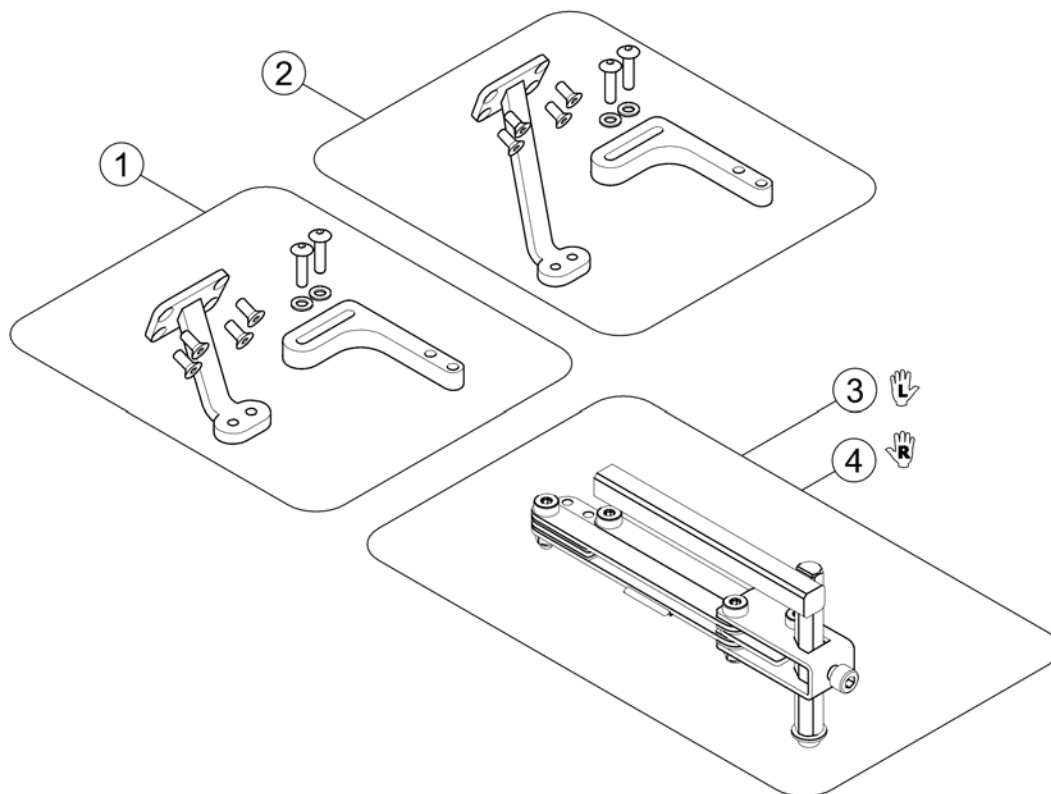
1508519.no



1	SP1587331	(Remote bracket swingaway G90)	1
2	SP1587696	(Remote bracket MACES, swingaway - G90)	1
3	SP1496249	SWING AWAY VE NY MOD2008 KPL	1
4	SP1496248	SWING AWAY HØ NY MOD2008 KPL	1

(Remote Bracket Attendant)

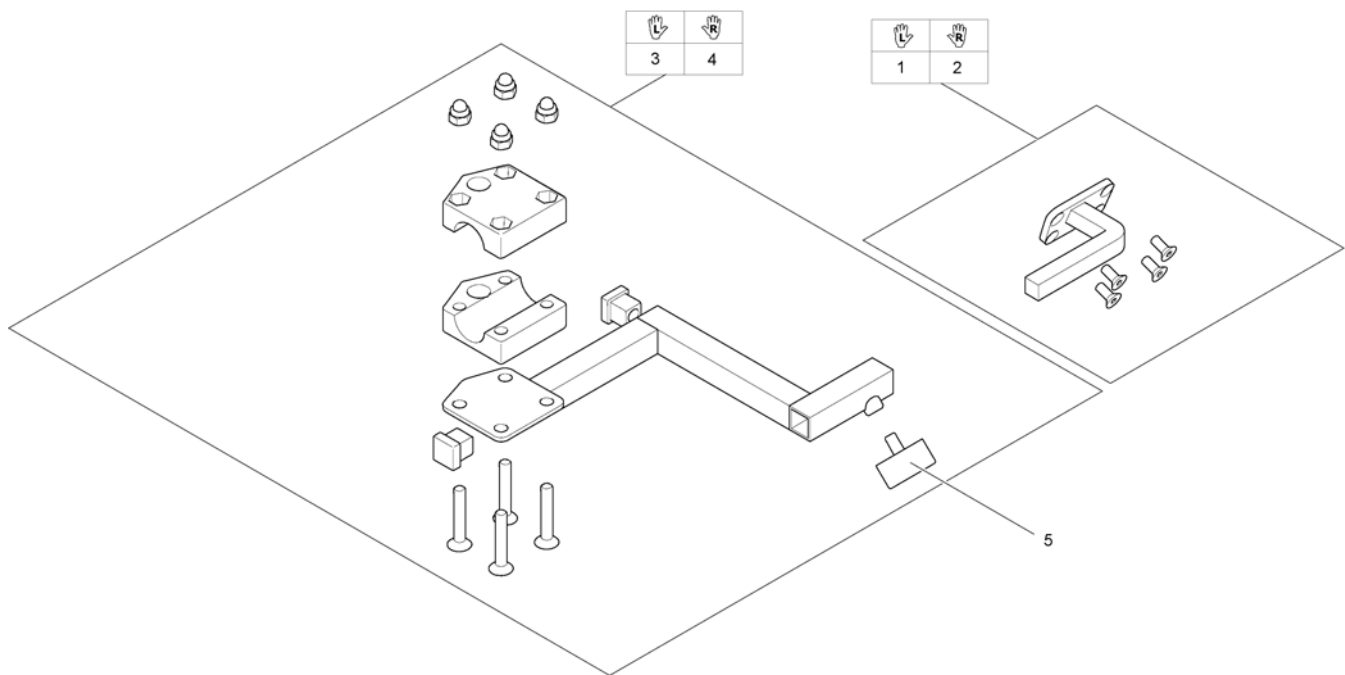
1560853.iso





1	SP1509435	(Remote Support attendant control G90A, LH)	1
2	SP1509439	(Remote Support attendant control G90A, RH)	1
3	SP1422937	(Bracket for attendant control LH)	1
4	SP1422936	(Bracket for attendant control RH)	1
5	L31001	VINGSKRUV M 6X10 GARANT G24 S	1

(Remote Bracket Attendant Ultra Low Maxx Seat)

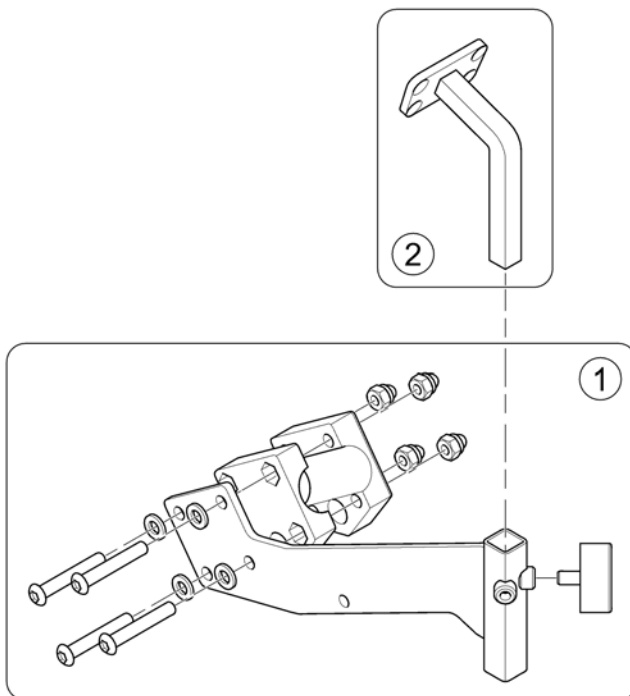
1588438 no



1	SP1588427	BRAKETT LEDSAGERSTYRING STYREBOKS BAK	1
2	SP1587229	(Adapter G90 (with screws))	1

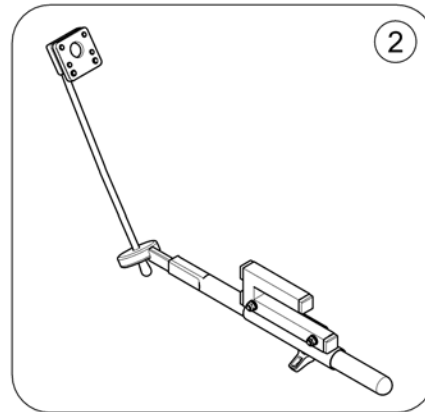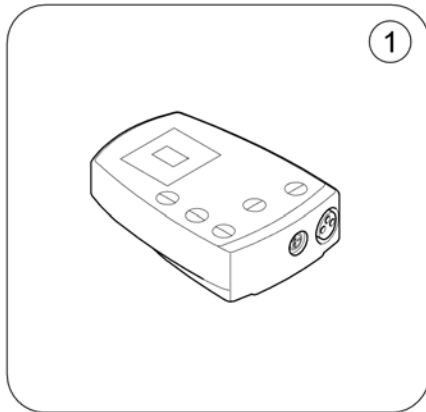

1	SP1505187	(Remote Control G91S) [1615]	1
2	-	(Remote Bracket Standard) [1615+1207/1208]	1

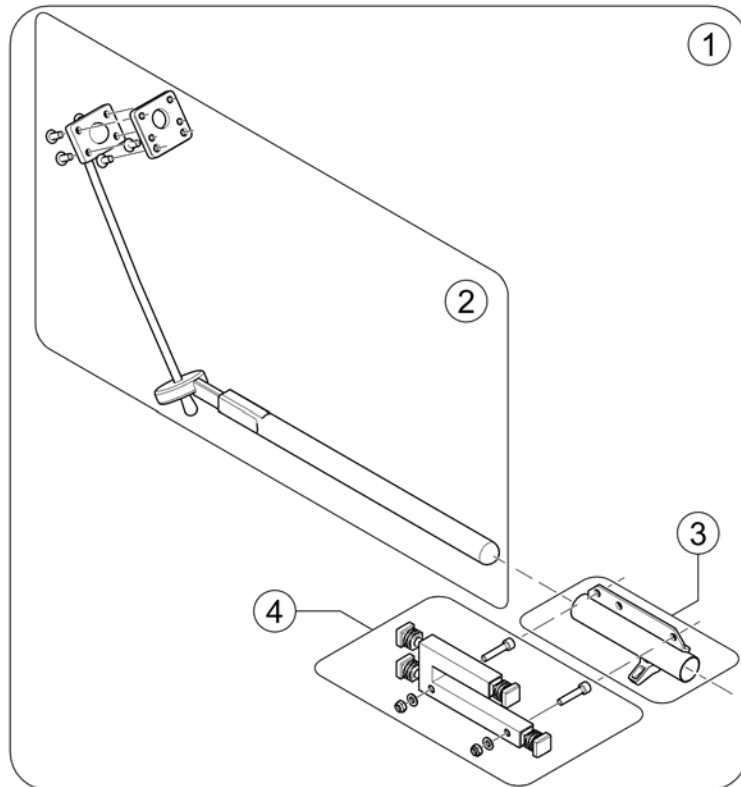


1560654.lin

1	SP1509690	HOLDER FOR DX-REMG90 HØ JUSTERBAR HØYDE/DYBDE		1
2	SP1527504	(Swan neck remote bracket with adapter plate)		1
3	SP1603801	(Mounting bracket, cpl.) <i>Old bracket no longer available!</i>	05.08.2013 →	1
4	SP1603847	(Mounting Adapter G91) <i>Old bracket no longer available!</i>	05.08.2013 →	1

(Centre control REM 41)


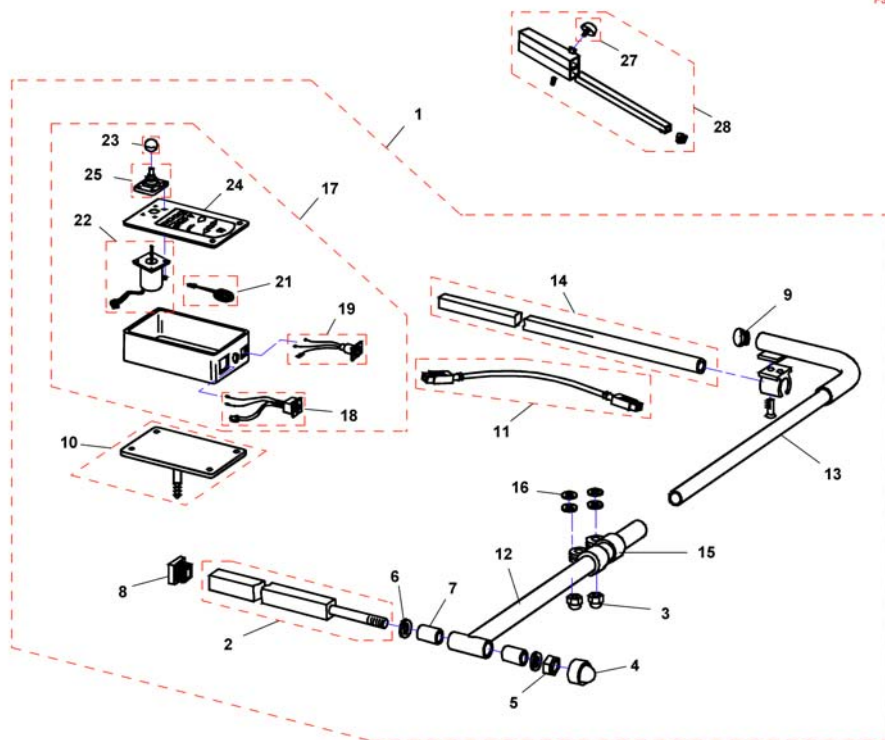
F383057.dwg



1	SP1530588	DX-REM41 MED SENTERFESTE TDX	1
1	SP1530589	(CENTRAL STEERING CPL., SWING-AWAY LH, REM41-D)	1
2	Z177024.00	FESTE F MITTSTYRNING,	1
3	F60208	(Cap nut DIN1587-M08 zc)	1
4	L01037	(Sks Nut Cap M12)	1
5	F20106	MUTTER SEKSKANT SELVLÅSENDE G40	1
6	G00104	SKIVE 13X24X2,5 G40	1
7	Z17032A	(Bush Rd. 19 X 3.5 X 25)	1
8	L01055	ENDESTYKKE FESTEBR.	1
9	L01008	(Tube End Plug SI 22)	1
10	Z183216.P40	FESTE F MANBOX MITTST,	1
11	1559862	KABEL BUS 1,7 M	1
12	Z183089.00	BREDD JUST ARM MITTST,	1
13	Z183086.00	ÄNDRÖR HÖ, STORM	1
13	Z183092.00	ÄNDRÖR VÅ, STORM	1
14	H281005A.00	(Support rail, chrome plated)	1
15	5006720	KLEMME 22MM TDX	1
16	G00106	SKIVE M8 DRIVHJUL MIRAGE	1
17	SP1527930	BORDSTYRING DX-REM41-D DOLPHIN MASTER	1
18	P41108	KONTAKT BUS KABEL DX-REM41	1
19	P41106	LADEKONTAKT REM41	1
21	P30493	(HORN)	1
22	GSM55120	(ACS-remote WITH 180MM CABLE)	1
23	L36999	(Joystick Knopf Standard Ball)	1
24	-	(No longer available.)	0
25	Z183001A	JOYSTICK GUMMI ACS	1
26	Z183289	SENTERSTYRT STØTTE RØR KPL	1
27	L31001	VINGSKRUV M 6X10 GARANT G24 S	1
28	Z183289	SENTERSTYRT STØTTE RØR KPL	1


-	(Small Chincontrol)	1
-	(Framework for chincontrol)	1

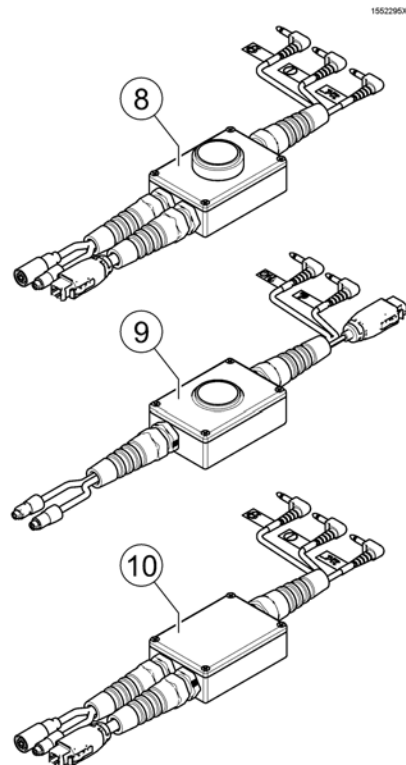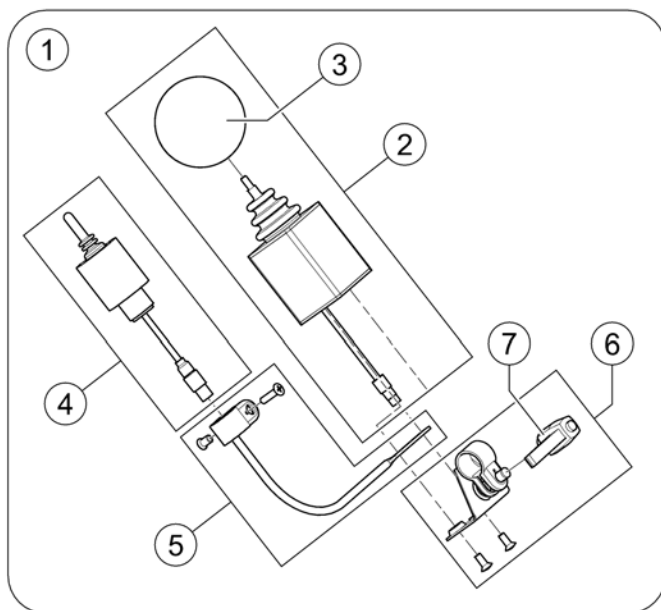(Small Chincontrol)



1	1448705	\$HAKESTYRING KPL MONT EKSTRA JOYSTICK	1
1	1421770	(Chin control G80/G90, cpl.)	1
1	SP1544255	(Chin control G91, cpl.)	1
2	Z183282	(Joystick, DX-RJM (with Binder plug, 4-pole))	1
3	Z183265	BRUK GJS26826 DYNAMICS	1
4	1442943	(Switchjoystick for REM 24 SD, cpl.)	1
4	Z183283	(Switchjoystick G80/G90, cpl.)	1
4	SP1544250	(Switchjoystick for G91S, cpl.)	1
5	SP1552294	(Switchjoystick bracket)	1
6	SP1552295	(Joystick clamp)	1
7	L32901	SPAK M8 STORM4	1
8	SP1442944	(Harness Chin Control for REM 24 SD)	1
8	1442963	\$KABELSTAMME VERTIKAL SWINGAWAY HAKESTYRING	1
9	Z183275	(Harness Chin Control for G80/G90)	1
10	SP1541715	(Harness Chin Control for G91S)	1


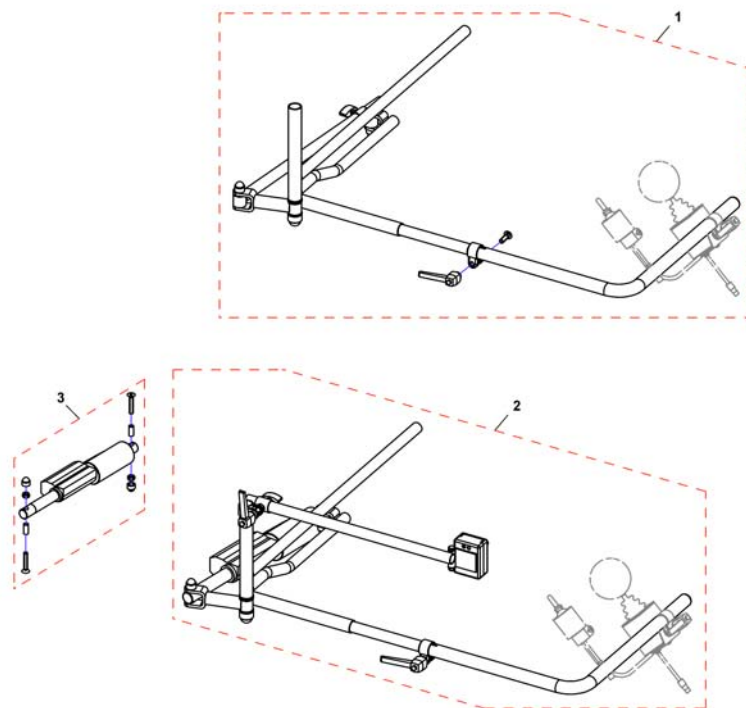
(Framework for chincontrol)





-	Z167994	(CHIN CONTROL SUPPORT . CPL.) <i>Not shown / Nicht abgebildet / Non illustré</i>	1
1	SP1422943	HAKESTYRING RAMMEVERK HORIZONTAL SWING AWAY VE	1
1	SP1422942	(Framework for chin control, manual, horizontal swing away, RH, cpl.)	1
2	SP1422945	HOLDER HAKEBETJ EL HORIS SVINGAWAY VE KPL <i>incl. Actuator</i>	1
2	SP1422944	(Framework for chin control, powered, horizontal swing away, RH, cpl.) <i>incl. Actuator</i>	1
-	SP1552335	(Hardware, chin control, powered, horizontal swing away, LH) <i>Not shown / Nicht abgebildet / Non illustré</i>	1
-	SP1552336	(Hardware, chin control, powered, horizontal swing away, RH) <i>Not shown / Nicht abgebildet / Non illustré</i>	1
-	SP1422946	(Framework for chincontrol, powered, vertikal swing away) <i>Not shown / Nicht abgebildet / Non illustré</i>	1
3	SP1552337	(Actuator for horizontal swing away chin control, cpl.)	1
3	SP1552418	(Gestell für Kinnbedienung, elektrisch, vertikal abschenkbar) <i>Not shown / Nicht abgebildet / Non illustré</i>	1

(Chin Control - Modulite Seat (15/10/01 =>))

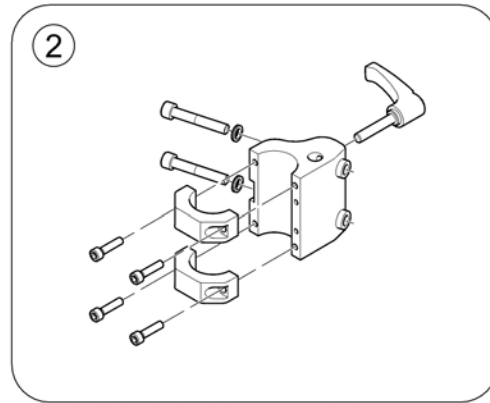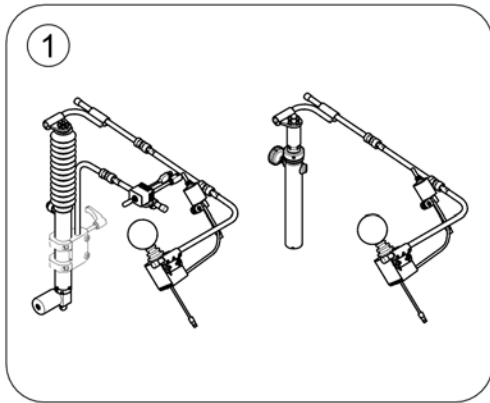1591326.iss




-	-	(If retrofitting, the user manual must be ordered separately!)	1
-	-	(User Manuals Chin Control (2015-10-01 =>))	1
1	-	(Chin Control DX-RJM)	1
2	-	(Mounting Hardware - Chin Control)	1

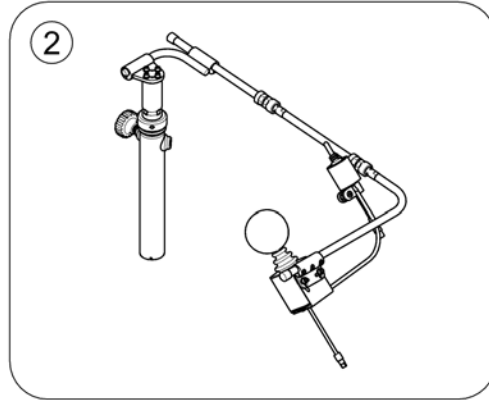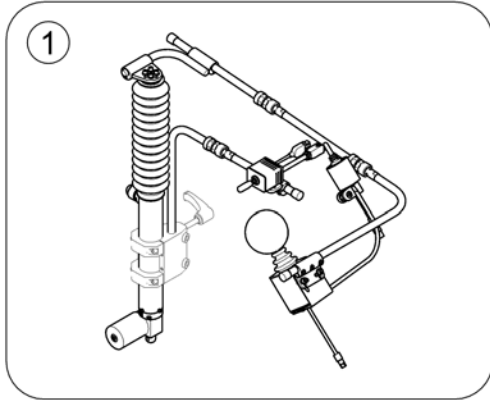

1	-	(Chin Control powered - Details)	1
2	-	(Chin Control manual - Details)	1
-	SP1601389	(Kit, Swing-away mounting arm, powered, LH, cpl. with stereo-switch (no joystick nor harness included)) <i>Please order respective harness & joystick separately!</i>	1
-	SP1601390	(Kit, Swing-away mounting arm, powered, RH, cpl. with stereo-switch (no joystick nor harness included)) <i>Please order respective harness & joystick separately!</i>	1
-	SP1601391	(Kit, Swing-away mounting arm, manual, LH, cpl. (no joystick nor harness included)) <i>Please order respective harness & joystick separately!</i>	1
-	SP1601408	(Kit, Swing-away mounting arm, manual, RH, cpl. (no joystick nor harness included)) <i>Please order respective harness & joystick separately!</i>	1
-	SP1601528	(Kit, Joysticks & Harness for REM 24 SD)	1
-	SP1601529	(Kit, Joysticks & Harness for G91S)	1


(Chin Control powered - Details)

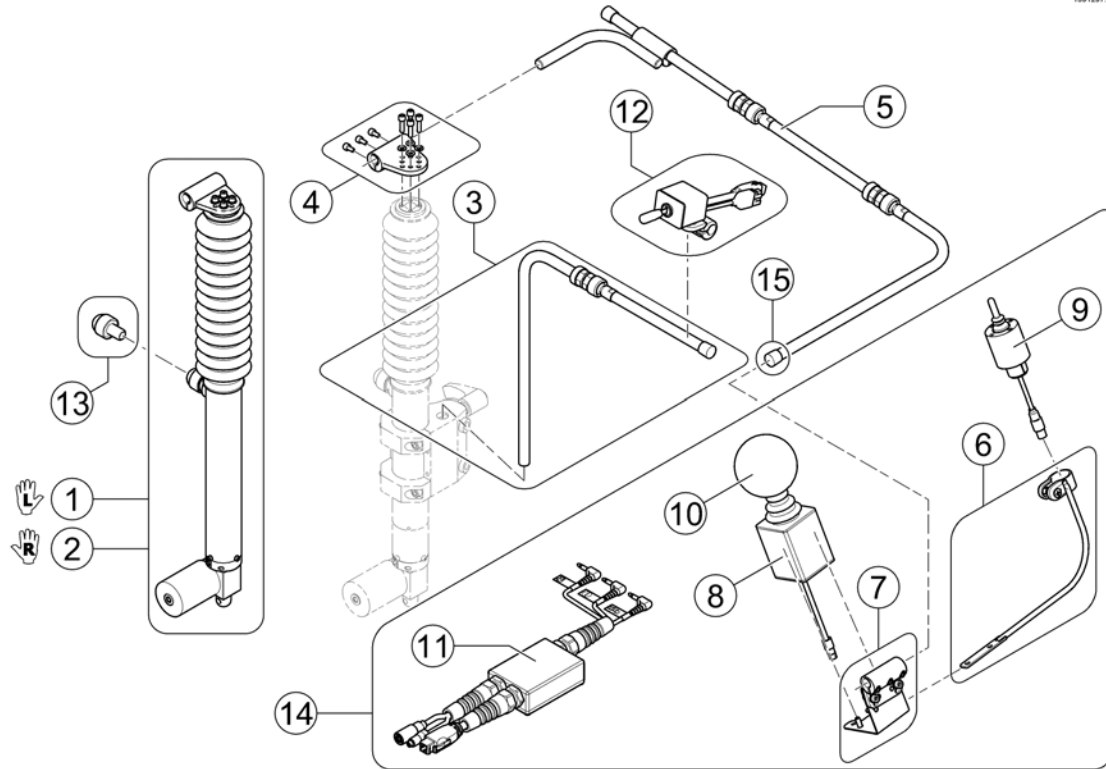1591297.jpg




-	SP1605661	(Fixation Kit (cable ties, spiral coiled tube & end caps))	01.04.2016 →	1
1	SP1580094	(Actuator cpl. LH)		1
2	SP1591292	(Actuator cpl. RH)		1
3	SP1605659	(Remote Rod Assy) <i>New designed bracket (2016-04-01 =>) for Modulite Seat needed!</i>	01.04.2016 →	1
4	SP1589144	(Rod cylinder connection chin, cpl.)		1
5	SP1605660	(Joystick Rod Assy)	01.04.2016 →	1
6	SP1589146	(Mounting arm for Switch Joystick)		1
7	SP1589142	(Mounting hardware for DX-RJM, cpl.)		1
8	Z183282	(Joystick, DX-RJM (with Binder plug, 4-pole))		1
9	1442943	(Switchjoystick for REM 24 SD, cpl.)		1
9	SP1544250	(Switchjoystick for G91S, cpl.)		1
10	Z183265	BRUK GJS26826 DYNAMICS		1
11	SP1442944	(Harness Chin Control for REM 24 SD)		1
11	SP1541715	(Harness Chin Control for G91S)		1
12	SP1580701	(Stereo-switch box, cpl. with mounting hardware)		1
13	SP1589145	(Adjustment Screw)		1
14	SP1601528	(Kit, Joysticks & Harness for REM 24 SD)		1
14	SP1601529	(Kit, Joysticks & Harness for G91S)		1
15	SP1602721	(End cap A12, black)	01.04.2016 →	3



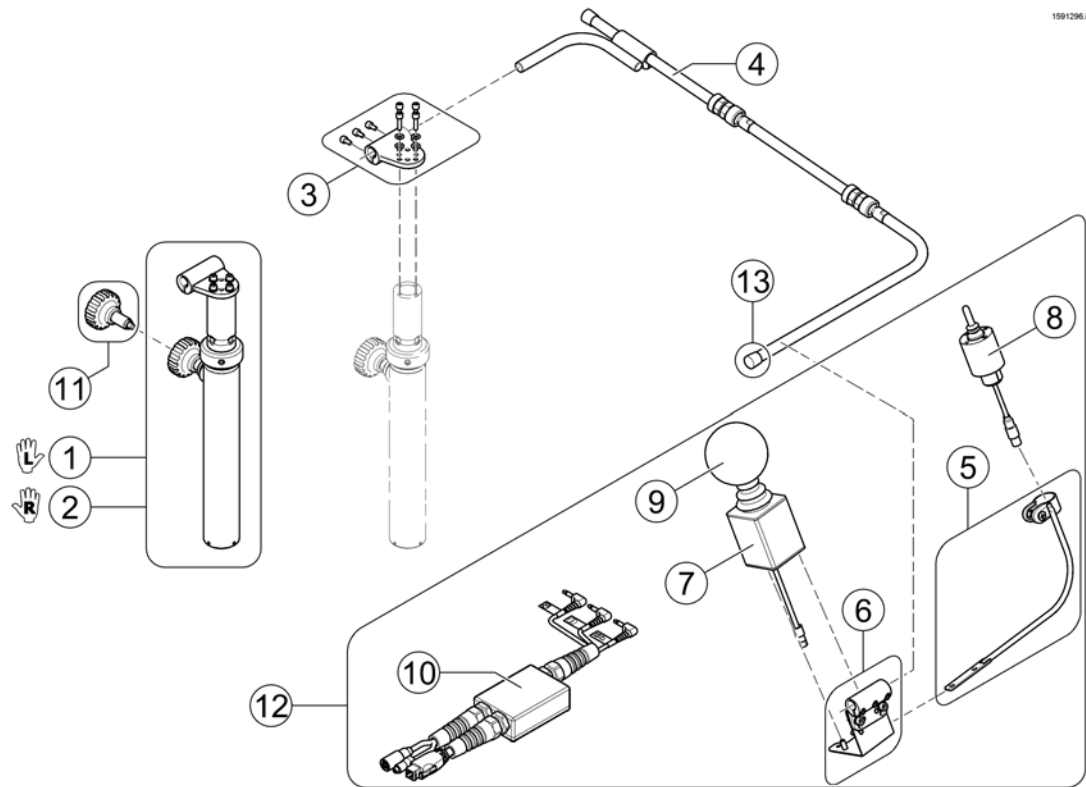
(Chin Control manual - Details)



-	SP1605661	(Fixation Kit (cable ties, spiral coiled tube & end caps))	01.04.2016 →	1
1	SP1577456	(ASM Chin Control man. LH, cpl.)		1
2	SP1591293	(ASM Chin Control man. RH, cpl.)		1
3	SP1589144	(Rod cylinder connection chin, cpl.)		1
4	SP1605660	(Joystick Rod Assy)	01.04.2016 →	1
5	SP1589146	(Mounting arm for Switch Joystick)		1
6	SP1589142	(Mounting hardware for DX-RJM, cpl.)		1
7	Z183282	(Joystick, DX-RJM (with Binder plug, 4-pole))		1
8	1442943	(Switchjoystick for REM 24 SD, cpl.)		1
8	SP1544250	(Switchjoystick for G91S, cpl.)		1
9	Z183265	BRUK GJS26826 DYNAMICS		1
10	SP1442944	(Harness Chin Control for REM 24 SD)		1
10	SP1541715	(Harness Chin Control for G91S)		1
11	SP1580666	(Adjustment Handle)		1
12	SP1601528	(Kit, Joysticks & Harness for REM 24 SD)		1
12	SP1601529	(Kit, Joysticks & Harness for G91S)		1
13	SP1602721	(End cap A12, black)	01.04.2016 →	2

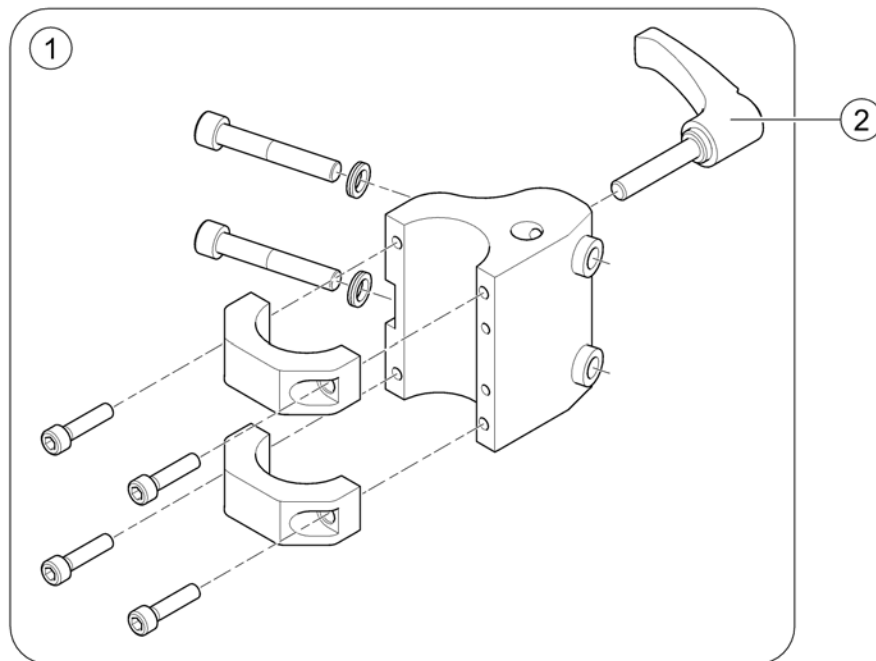


1	SP1577173	(Mounting Hardware, Chin Control - Modulite Seat) <i>Design change 2016-04-01 =></i>	1
2	SP1604355	(Clamp Lever M8x40)	1

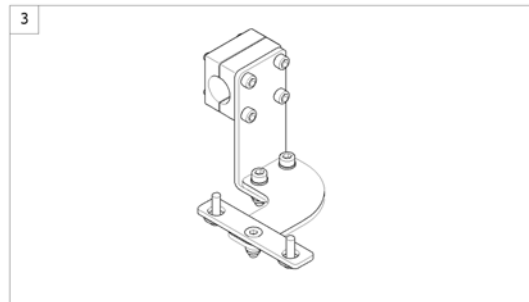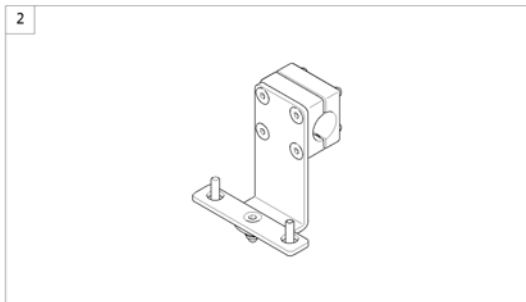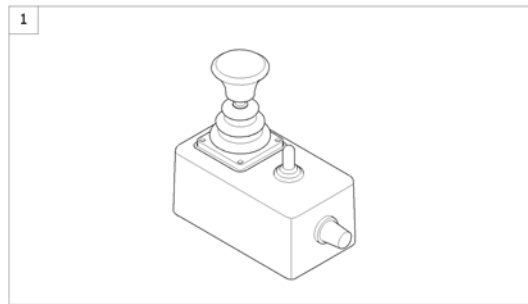

1	6001425	(ACS dual control) <i>[1291/1292]</i>	1
2	-	(Remote Bracket Dual Control for Modulite Seat)	1
3	-	(Remote Bracket Dual Control for Max Seat)	1

(Remote Bracket Dual Control for Modulite Seat)

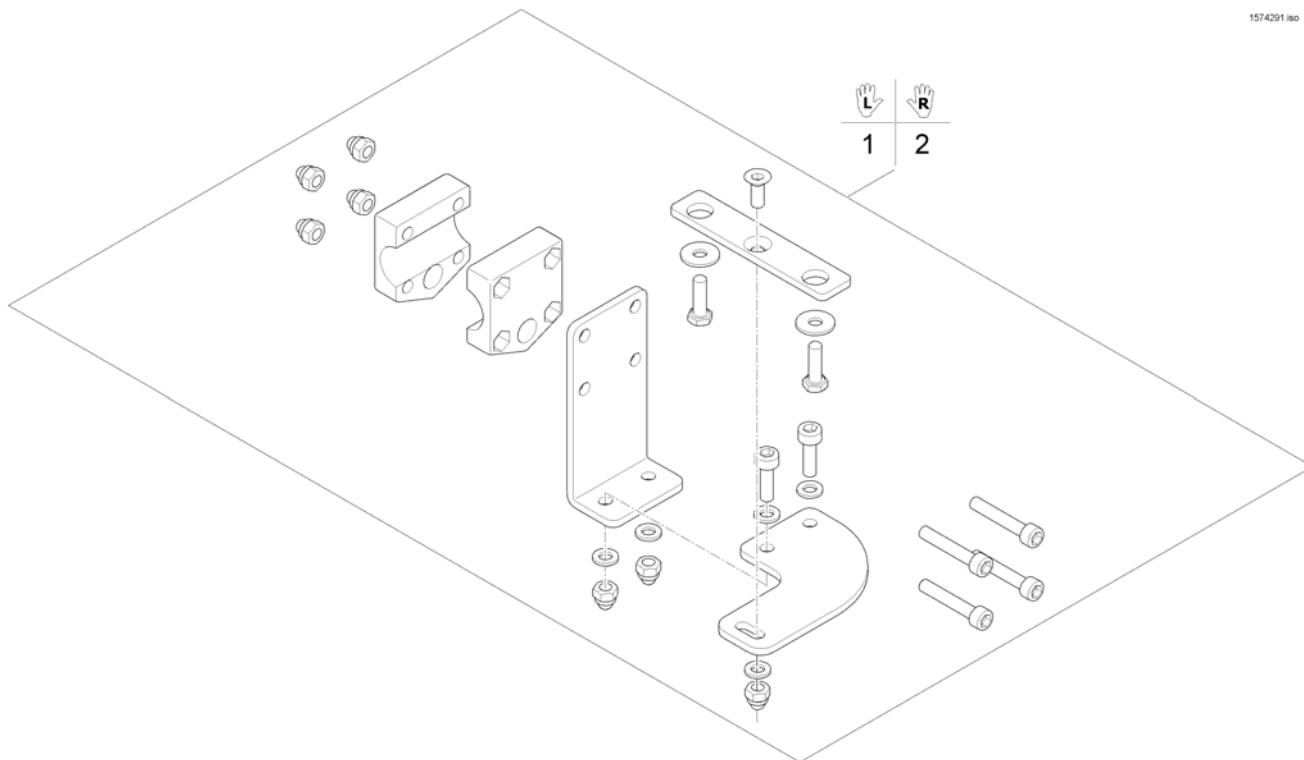
1574291 iso



1	SP1571799	(Remote Bracket Dual Control LH)	1
2	SP1571800	(Remote Bracket Dual Control RH)	1

Invacare® TDX-SP2 Lowrider, TDX-SP2 NB | Elektronikk
(Remote Bracket Dual Control for Max Seat)

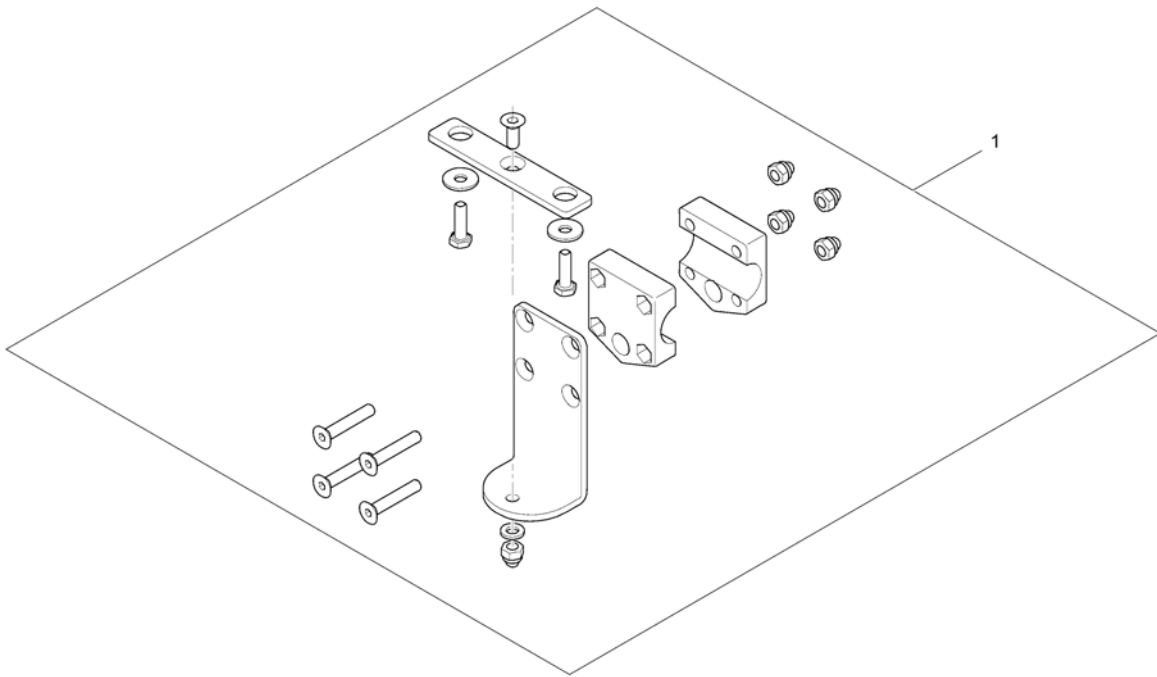1560885.iso




1 SP1422818

FESTEBRAKETT LEDSAGERSTY RYGGMONT ENEKL


1

(Compact Dual Control with brackets)

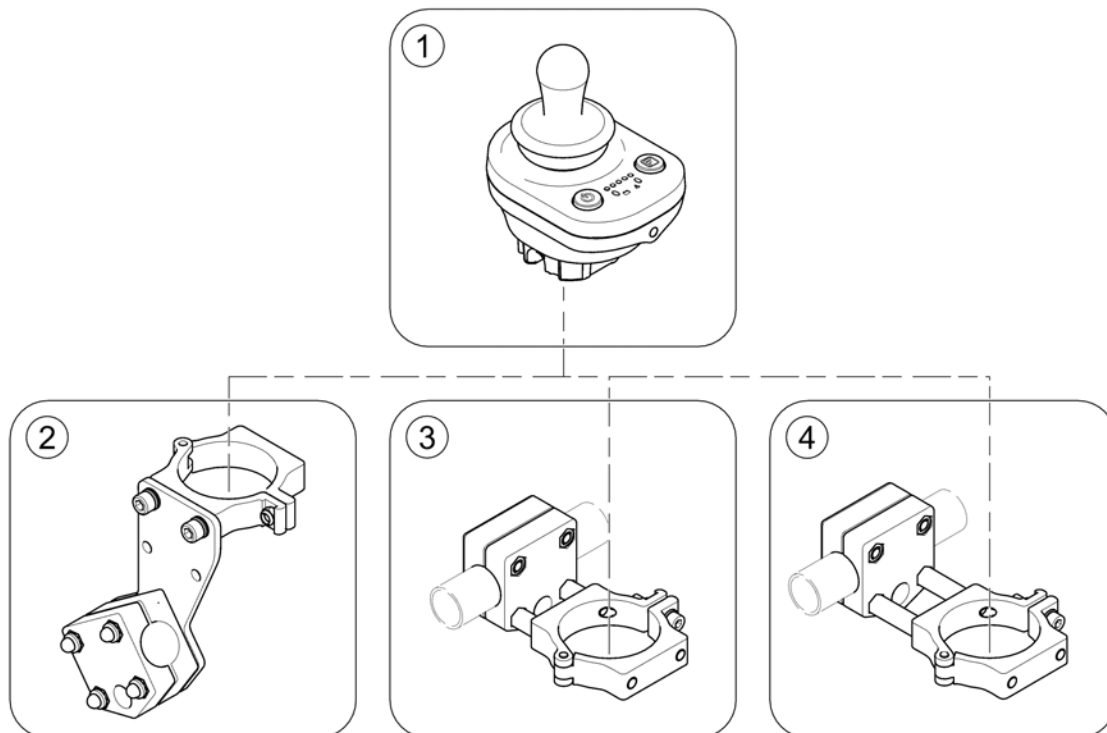1588511.000



1	-	(Compact Dual Control) [1720]		1
2	-	(Bracket, Compact Dual Control - Ultra Low Maxx Seat) [1207/1208 + 1720 for 2020]	03.06.2015 →	1
3	-	(Bracket, Compact Dual Control - short)	01.07.2015 →	1
4	-	(Bracket, Compact Dual Control - long)	01.07.2015 →	1

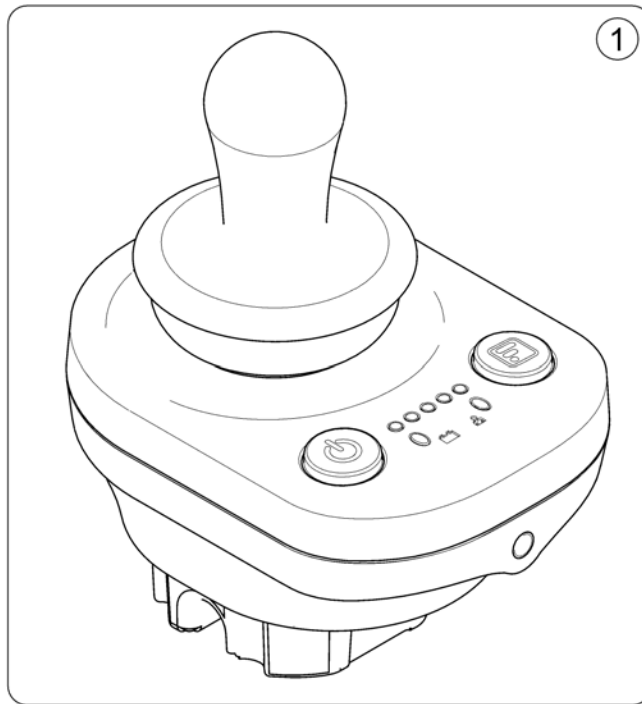


1 SP1586988

(Remote DX2-ACU)

1

(Bracket, Compact Dual Control - Ultra Low Maxx Seat)

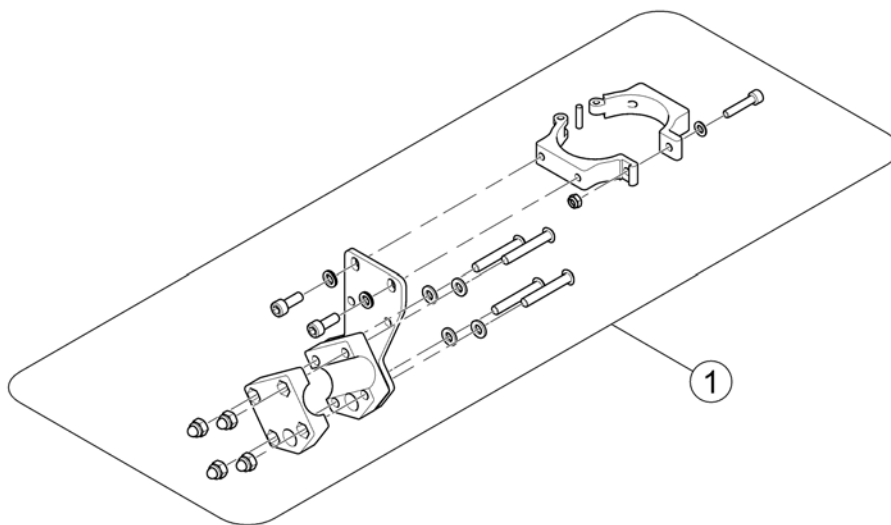
1588512.no



1 SP1586795

HOLDER STYREBOKS DUAL CONTROL

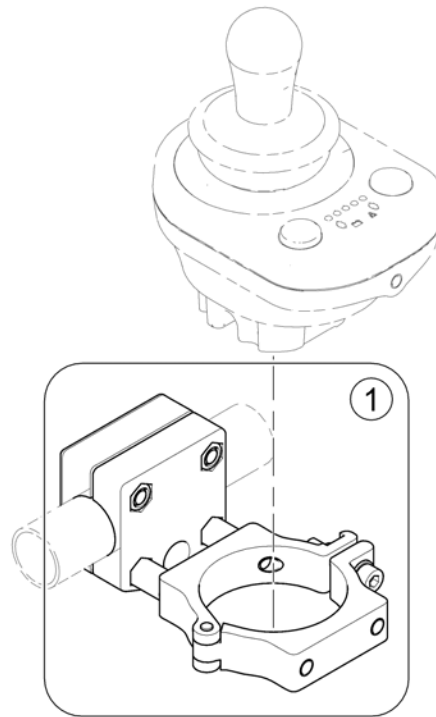1




1 SP1587904

(Remote Attachment - short)

1

(Bracket, Compact Dual Control - long)

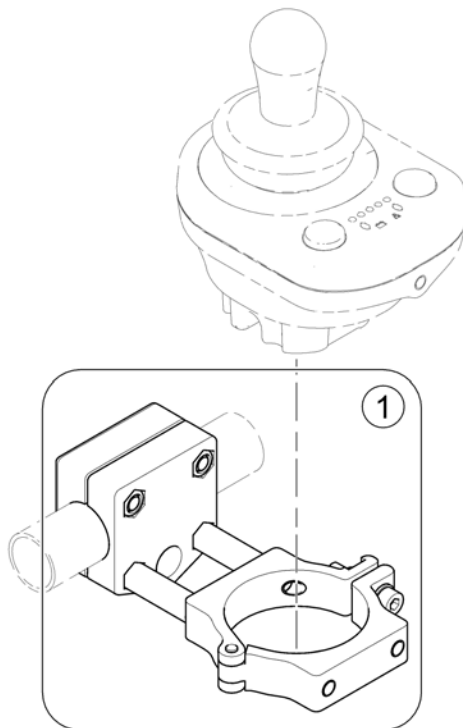158856 iso



1 SP1587119

(Remote Attachment - long)


1

(Intuitive Dual Control (IDC) with brackets)

-	(IDC with bracket)	1
-	(IDC with bracket (Ultra Low Maxx Seat))	1
-	(Height adjustable bracket IDC)	1

(IDC with bracket)

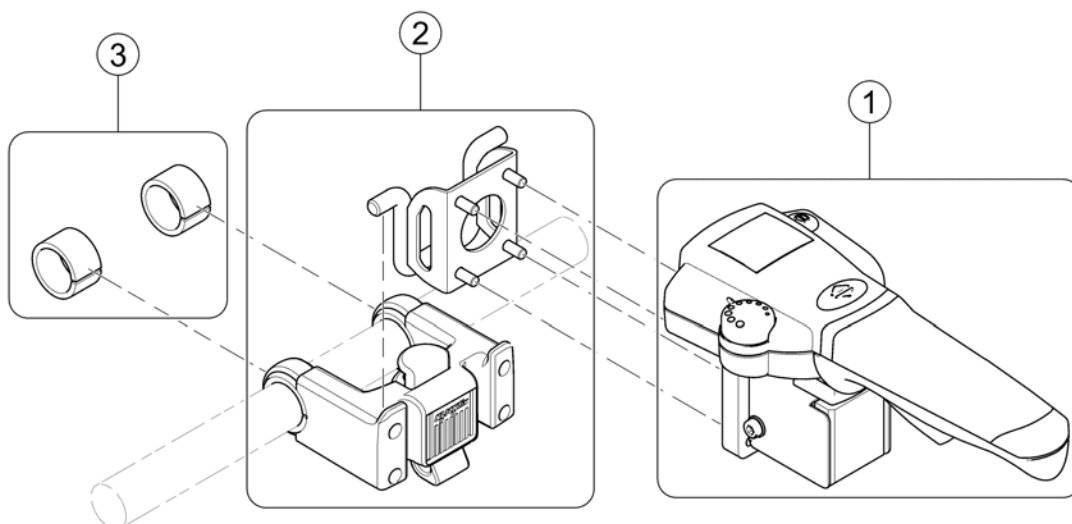
1602413.iss



-	(The IDC incl. mounting hardware has been revised, the old hardware is no longer available.)	01.01.2016 →	0
-	(If retrofitting, the user manual must be ordered separately!)		1
-	(User Manuals IDC (Intuitive Dual Control))		1
1	SP1603425 INTUITIV LEDSAGERSTYRING ALLE STOLER MED ACS		1
2	SP1603426 HOLDER INTUITIV LEDSAGERSTYRING NY TYPE		1
3	SP1601919 (Rubber cuts)		1

Invacare® TDX-SP2 Lowrider, TDX-SP2 NB | Elektronik
(IDC with bracket (Ultra Low Maxx Seat))

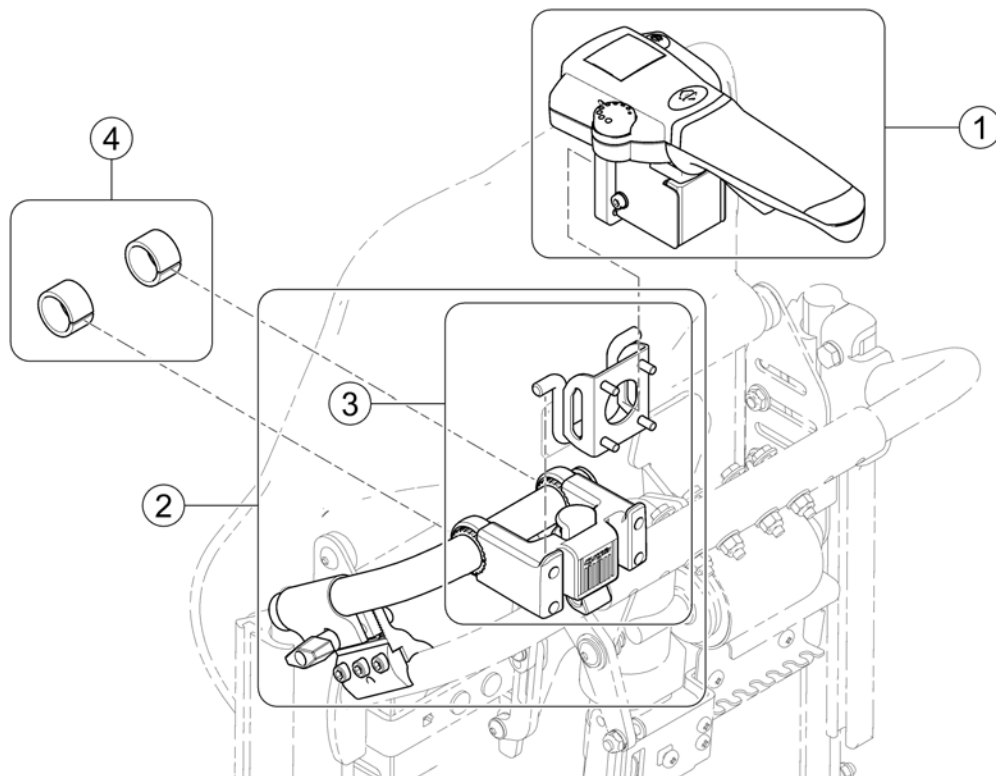
1602412.kim



-	(The IDC incl. mounting hardware has been revised, the old hardware is no longer available.)	01.01.2016 →	0
- -	(If retrofitting, the user manual must be ordered separately!)		1
- -	(User Manuals IDC (Intuitive Dual Control))		1
1	SP1603425 INTUITIV LEDSAGERSTYRING ALLE STOLER MED ACS		1
2	SP1603427 (IDC Bracket, cpl. with mounting hardware for Ultra Low Maxx Seat)		1
3	SP1603426 HOLDER INTUITIV LEDSAGERSTYRING NY TYPE		1
4	SP1601919 (Rubber cuts)		1

(Height adjustable bracket IDC)

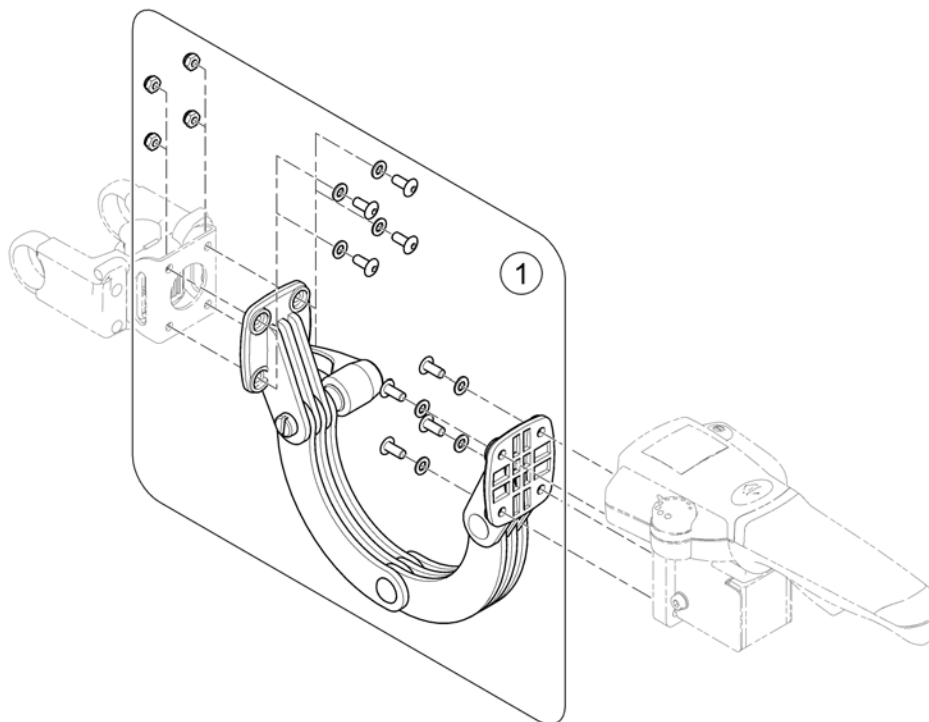1605484 iso



1 SP1604195

(Height adjustable bracket for IDC, cpl.)


1

(ACS Cables & Accessories)

-	(Bus-Cables ACS)	0
-	(Light Cable ACS)	0
-	(T-Adapter (Modulite & Max Seat))	0
-	(T-Adapter (Ultra Low Maxx Seat))	0
-	(T-Adapter (Low Rider))	0
-	(Tilt Sensors - Ultra Low Maxx Seat)	0
-	(Programming Kits)	0
-	(Magnetkey for Remotes)	0


(Bus-Cables ACS)

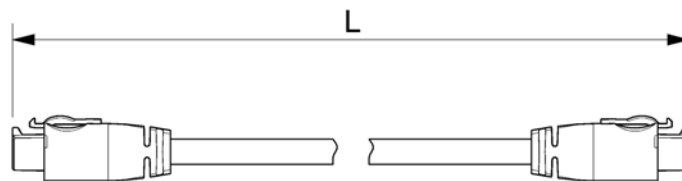ACSBUSKABEL.no



Position	L [mm]
1	300
2	500
3	700
4	1100
5	1700
6	2000
7	2500
8	2700
9	1060

1	GSM63003	KABEL BUS 0,3M ACS,	1
2	SPGSM63005	(Bus cable ACS - 500mm) <i>Replace GSM63004</i>	1
3	GSM63007	ACS-BUSKABEL L=700 MM	1
4	GSM63011	KABEL-BUS 1,1M ACS, ST utgått bruk 63010	1
5	1559862	KABEL BUS 1,7 M <i>Replace GSM63016</i>	1
6	GSM63020	DX BUS CABLE 2.0M STORM <i>Replace 1419386</i>	1
7	1440424	(Bus cable ACS, 2500mm)	1
8	SP1559873	(Bus cable ACS - 2700mm) <i>Replace GSM63026</i>	1
9	SP1550402	(Extension cable ACS - 1060mm)	1
-	GSM9887	DX-BUS GURLY CABLE	1

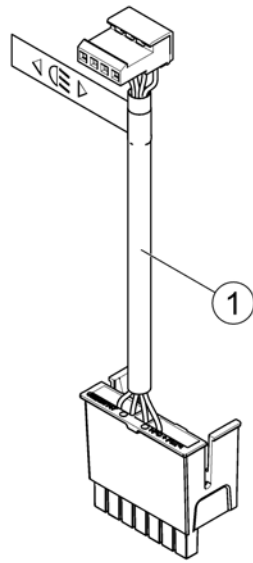


1 SP1524237

KABEL LYS TDX


1

(T-Adapter (Modulite & Max Seat))

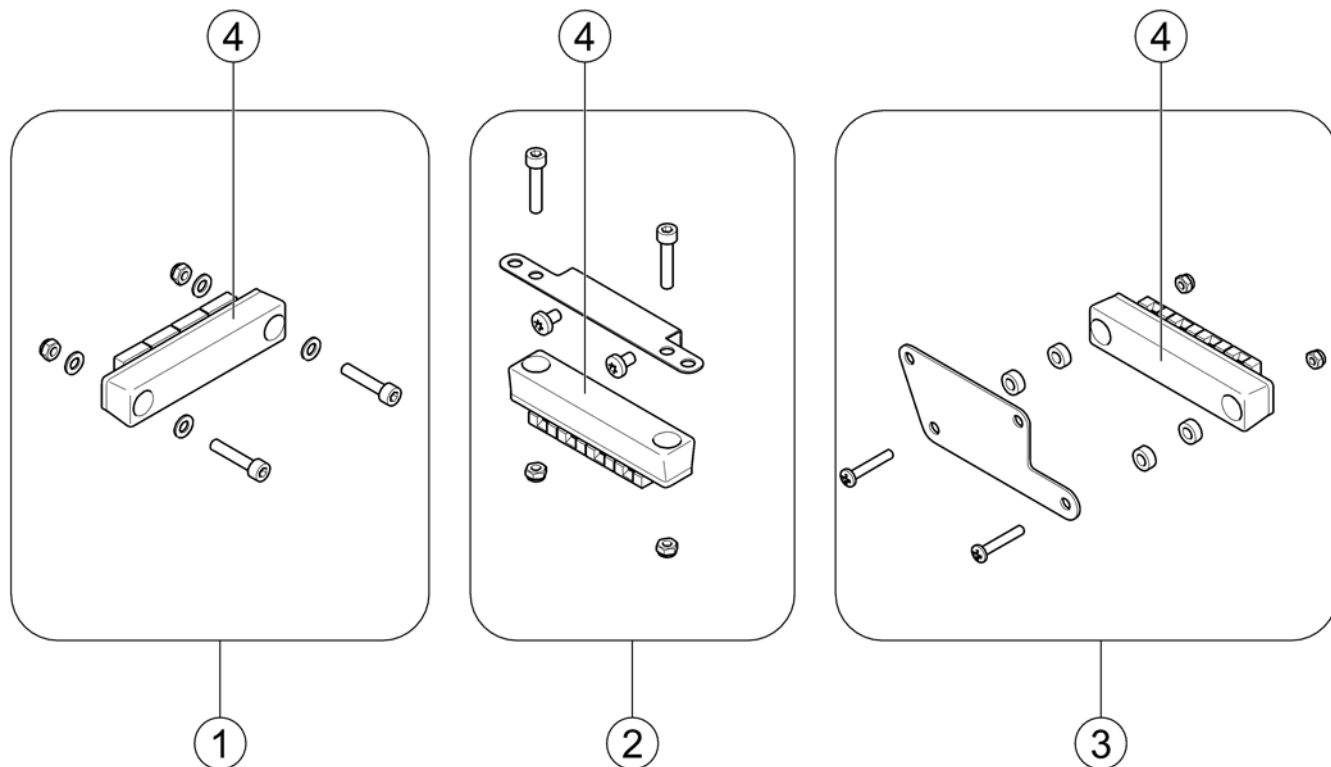1583359.iso




1	SP1571826	(T-Adapter - Modulite Seat)	1
2	SP1569958	(Bracket T-Adapter - Modulite Seat with MBL-Modul)	1
3	SP1580691	(T-Adapter - Max Seat)	1
4	P10072	(Bus 4-way socket (T-Adapter))	1

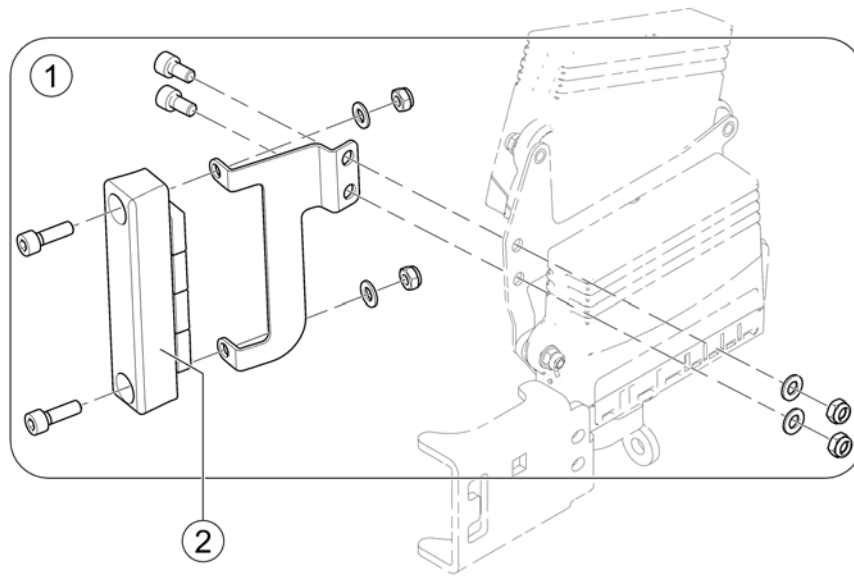

1	SP1586798	(T-Adapter with mounting hardware)	1
2	P10072	(Bus 4-way socket (T-Adapter))	1

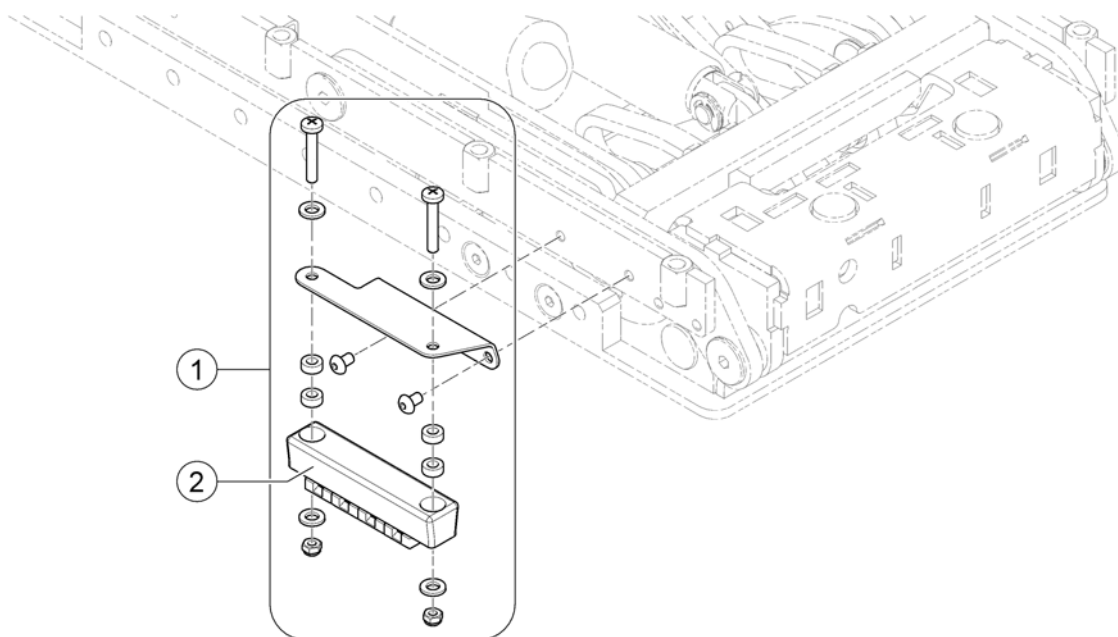(T-Adapter (Low Rider))




1039745.no

1	SP1608410	(T-Adapter, cpl. (Low Rider only))	1
2	P10072	(Bus 4-way socket (T-Adapter))	1

Invacare® TDX-SP2 Lowrider, TDX-SP2 NB | Elektronik
(Tilt Sensors - Ultra Low Maxx Seat)

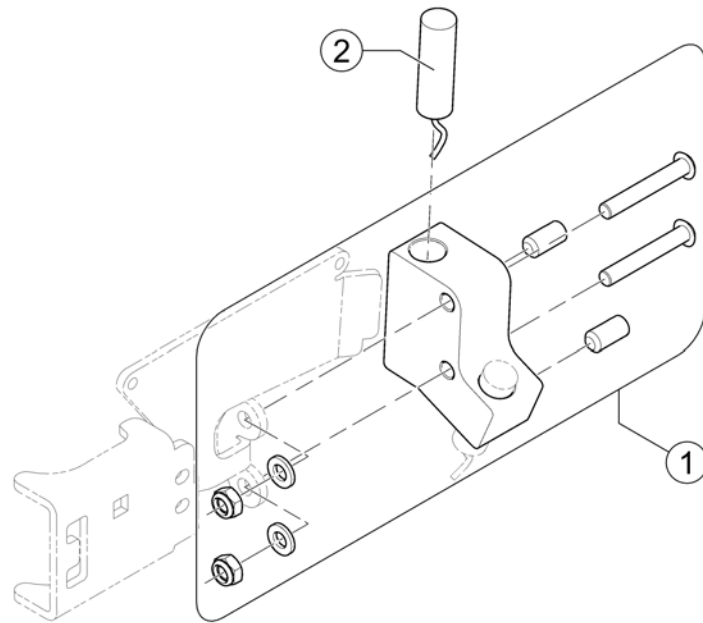1586648 Rev




1	SP1586648	(Switch Bracket)	1
2	SPTRE0753	(Tilt Switch)	2
-	SPTRE0784	(Tilt Switch Box)	1


Not shown / Non illustré / Nicht abgebildet

(Programming Kits)

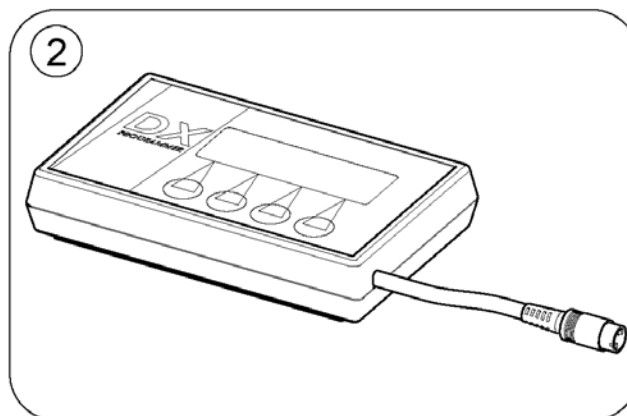1584805 iso



1	1449280	(USW wizard dongle kit, EDL-USW, cpl.)	1
2	6001450	(Hand Programmer, Instrucions In French, German, Spanish)	1
	SP1441551	(Programming Adapter Shark DK-ADAPT (necessary to program Shark))	1

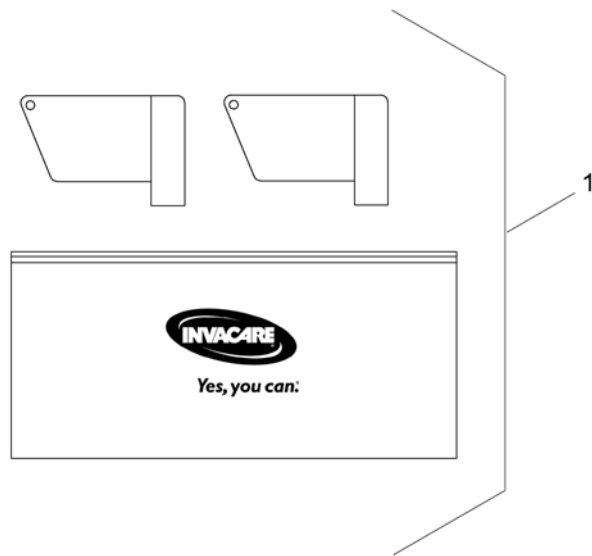

1 SP1419695


(Magnet key with bag)

1

(ACS Components)

-	(iPortal Modules and Accessories)	1
-	(Remote, Joystick, Control)	1
-	(Joystick Knob)	1
-	(Buddy Button)	1
-	(Piko Button)	1
-	(Interface, Adapter, Other)	1

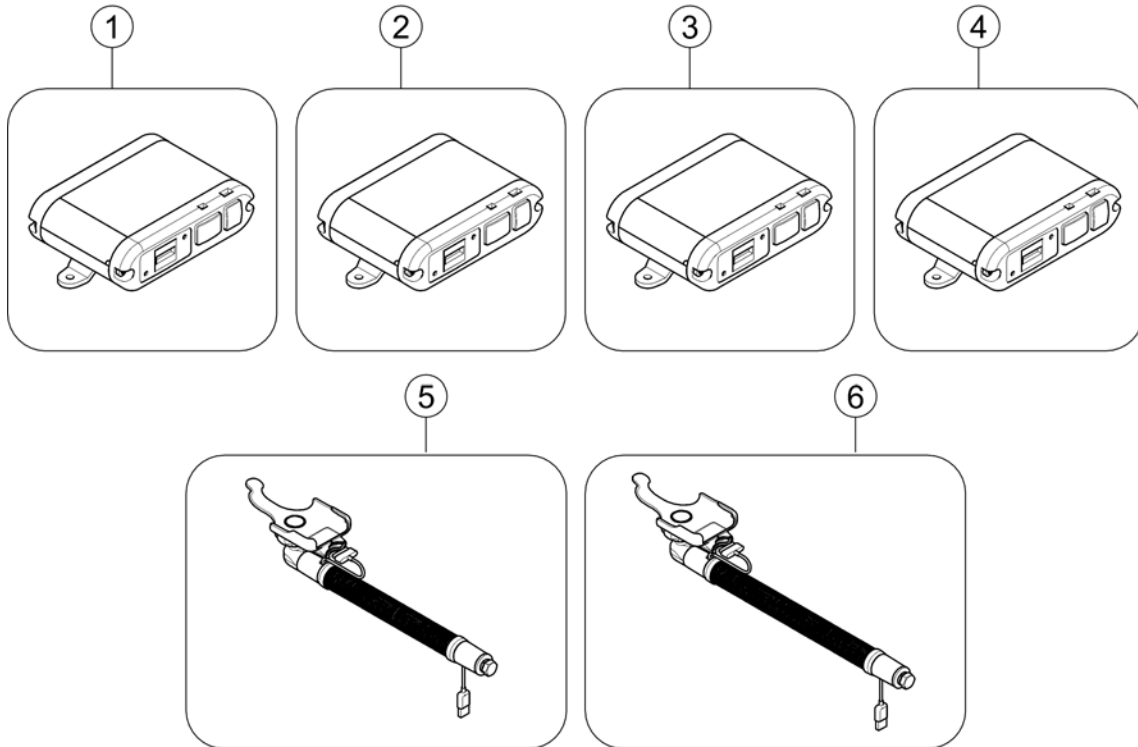
1603-426 iso



1	SP1553373	(iPortal Bluetooth Module Accessibility) [1629]	1
2	SP1570392	IPORTAL 2 MUS CTR IPORTAL 2 [1624]	1
3	SP1570393	(iPortal Bluetooth Module Mouse Mover & Accessibility) [1662]	1
4	SP1545530	(iPortal Module without Accessibility) [1634]	1
5	SP1545531	(Gooseneck short for iPod touch & iPhone) [1635]	1
6	SP1545532	(Gooseneck long for iPod touch & iPhone) [1636]	1

(Remote, Joystick, Control)

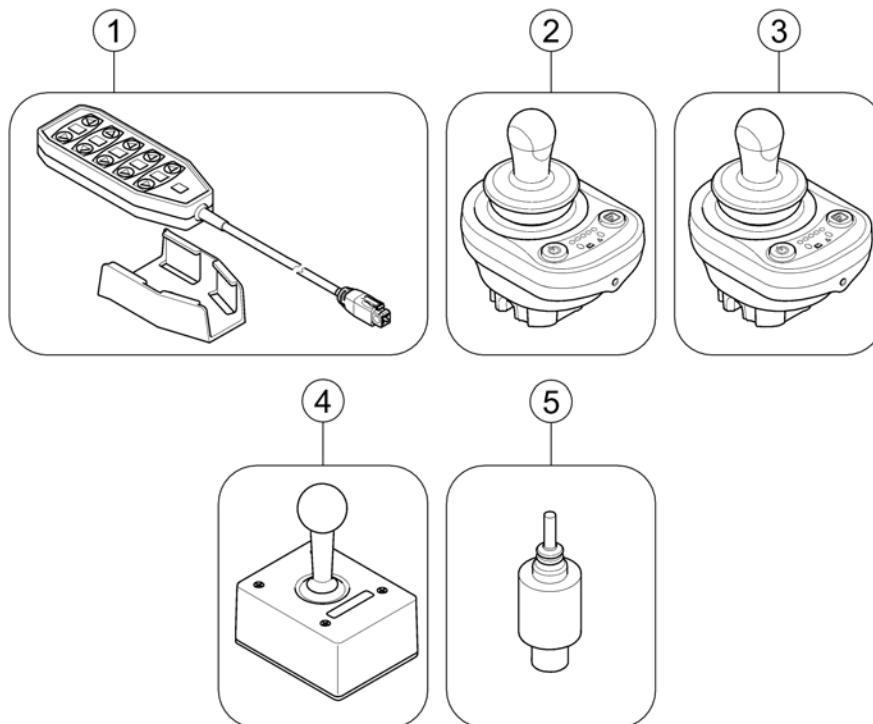1603436.lin



1	P10079	(Pendant control ARC 5) [1641/1642]	1
2	SP1586989	(Remote Joystick) [1621]	1
3	SP1586990	STYREBOKS LEDSAGER DX2 LOW-FORCE [1722]	1
4	P10171	(Heavy Duty Switch Joystick) [1620]	1
5	P10175	(Mini-Joystick with Push function) [1622]	1
-	SP1442679	(Head-Remote) [1613]	1
-	SP1603442	(Cable for Head-Remote ultrasonic porximity sensor)	1

(Joystick Knob)


1603437.000



1	L36999	(Joystick Knopf Standard Ball) [1650]	1
2	P10397	PROPP HAKESTYRING [1651]	1
3	P39302	JOYSTICK KNOTT CARROT ACS [1652]	1
4	P39303	(Joystick Knopf Flat T-Plate) [1653]	1
5	P39305	(Joystick Knob Flat T-Plate) [1654]	1
6	P10306	\$PROGRAMMERINGSSENHET DM4 GARANT 63E PRO [1655]	1
7	Z183265	BRUK GJS26826 DYNAMICS [1656]	1
8	P39306	(Joystick Knob Straight Bar) [1657]	1
9	P39307	(Joystick Knob Finger Tube) [1658]	1
10	P39304	(Joystick Knob Rount T-Bar) [1659]	1

(Buddy Button)

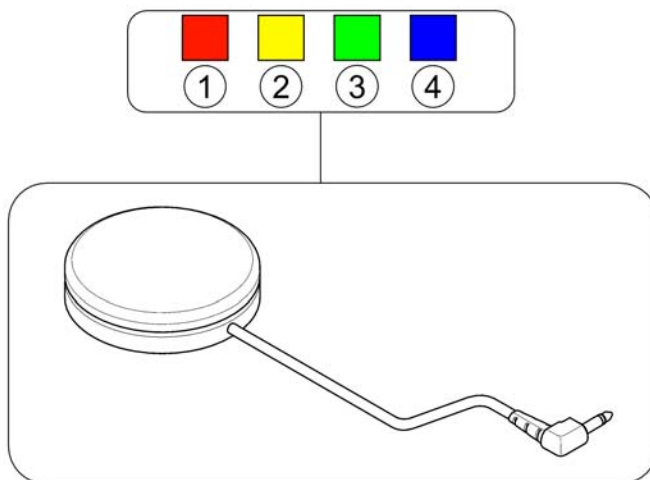
1603438.sas



1	5813964	BUDDY BUTTONS RØD MOOVER MED PILOT+ <i>[1685]</i>	1
2	P30093	(Buddy Button yellow) <i>[1686]</i>	1
3	P30092	(Buddy Button green) <i>[1687]</i>	1
4	P30089	(Buddy Button blue) <i>[1690]</i>	1

(Piko Button)

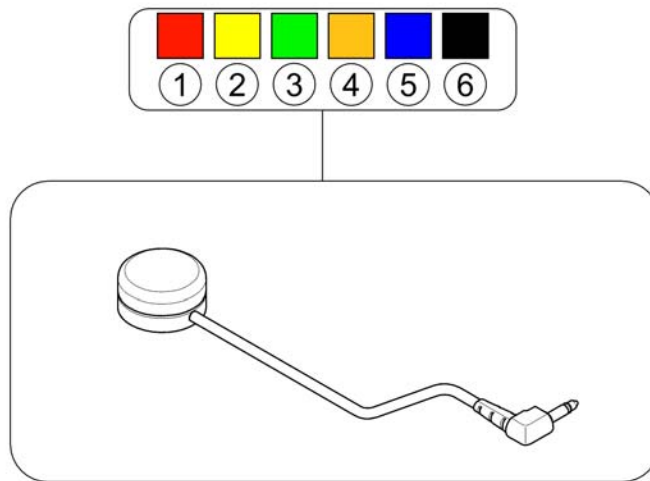
1603436.iss




1	SP1535717	PIKO BUTTON RØD DYNAMIC <i>[1711]</i>	1
2	SP1535718	(Piko Button yellow) <i>[1712]</i>	1
3	SP1535723	PIKO BUTTON GRØNN DYNAMIC <i>[1713]</i>	1
4	SP1535724	(Piko Button orange) <i>[1714]</i>	1
5	SP1535725	(Piko Button blue) <i>[1715]</i>	1
6	SP1535726	PIKO BUTTON SORT DYNAMIC <i>[1716]</i>	1

(Interface, Adapter, Other)

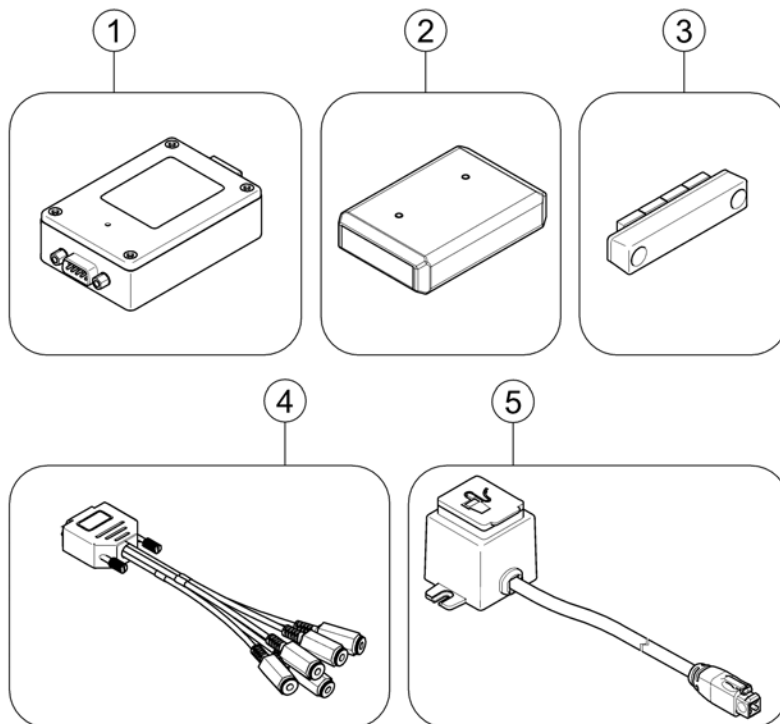1603403.iss



1	SP1544036	(4 Switch Interface) [1630]	1
2	P10163	(Sip & Puff Module(DX-SNP)) [1633]	1
-	Z183267	(Bracket for Sipp and Puff) [1668]	1
3	P10072	(Bus 4-way socket (T-Adapter)) [1680]	1
4	P30080	(5-1 Switch Adaptor) [1681]	1
5	P40107	(Charger Socket Modul) [1661]	1

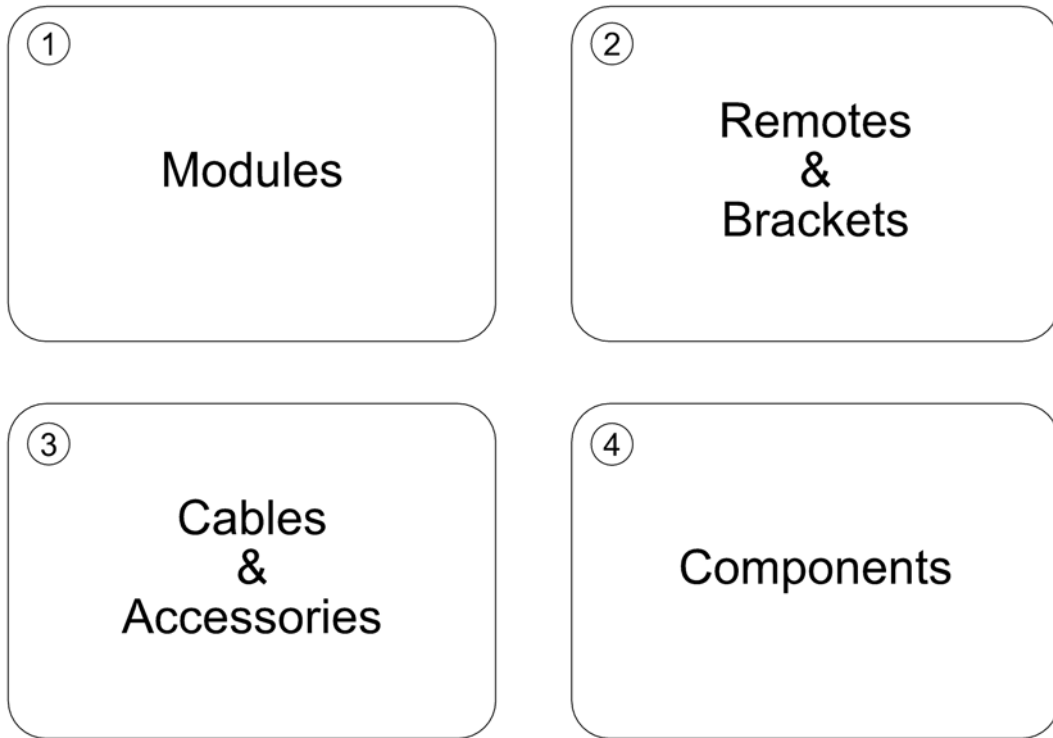


1	-	(Modules)	1
2	-	(Remotes & Brackets)	1
3	-	(Cables & Accessories)	1
4	-	(Components)	1

(Modules)

-	(Power Modules)	1
-	(Mounting Hardware, Power Module & Light Board)	1
-	(Actuator Modules (ACTs))	1
-	(Bracket ACT - Modulite Seat)	1
-	(Bracket ACT - Max Seat)	1
-	(G-Trac Sensor (Gyro))	1

(Power Modules)

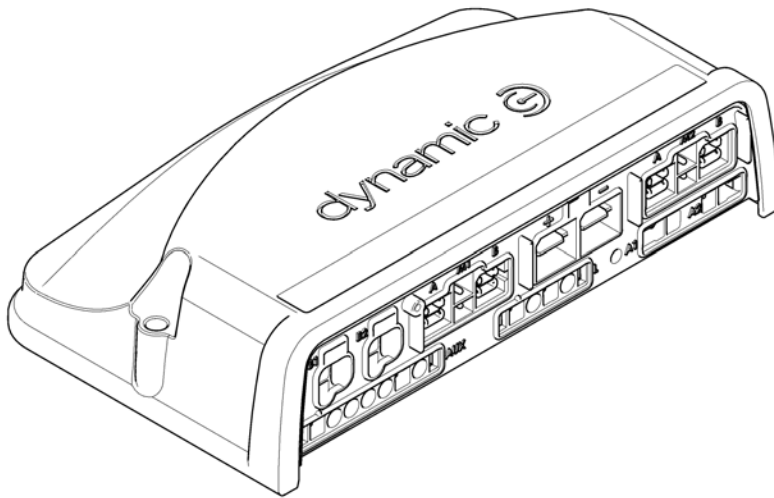160624 iso




TDX-SP2 Profiles

1	C02G06A0A00X-00 C02G08A0A00X-00 C02G10A0A00X-00
2	C02H06A1A00L-00 C02H06A2A00L-00 C02H06A3A00L-00 C02H06A1A00T-00 C02H06A2A00T-00 C02H06A3A00T-00
3	C02H08A1A00L-00 C02H08A2A00L-00 C02H08A3A00L-00 C02H08A1A00T-00 C02H08A2A00T-00 C02H08A3A00T-00
4	C02H10A1A00L-00 C02H10A2A00L-00 C02H10A3A00L-00 C02H10A1A00T-00 C02H10A2A00T-00 C02H10A3A00T-00

1	SP1605799	(Power Module LiNX 120A Spectra XTR2 HD, Kite HD, TDX-SP2 - 6 & 8km/h, Drive Only)	1
2	SP1605802	(Power Module LiNX 120A TDX-SP2 - 6km/h, Actuators & Light)	1
3	SP1605803	(Power Module LiNX 120A TDX-SP2 - 8km/h, Actuators & Light)	1
4	SP1605806	(Power Module LiNX 120A TDX-SP2 - 10km/h, Actuators & Light)	1

(Mounting Hardware, Power Module & Light Board)


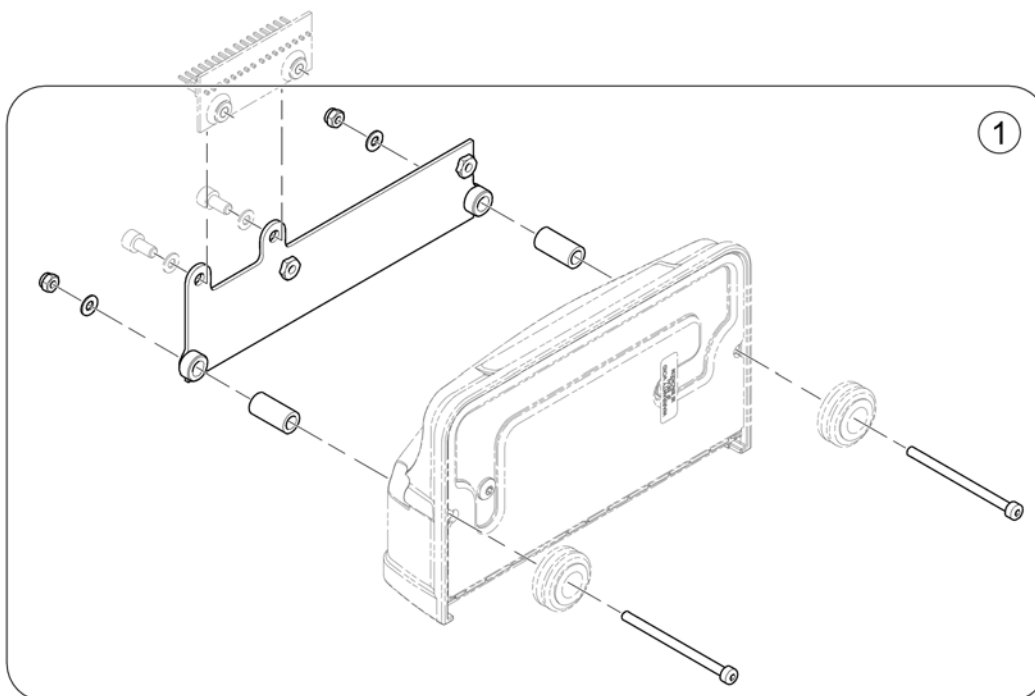
160664 iso



1 SP160665

(Mounting Hardware, Powermodule & Light Board)

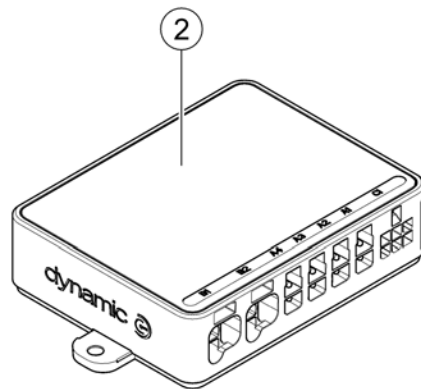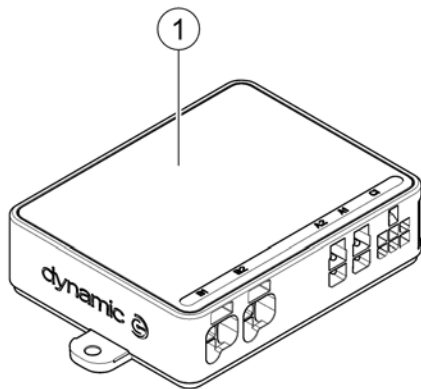
1



1	SP1602883	(LiNX 2-Channel Actuator Module (DLX-ACT200))	1
2	SP1602884	(LiNX 4-Channel Actuator Module (DLX-ACT400))	1

(G-Trac Sensor (Gyro))

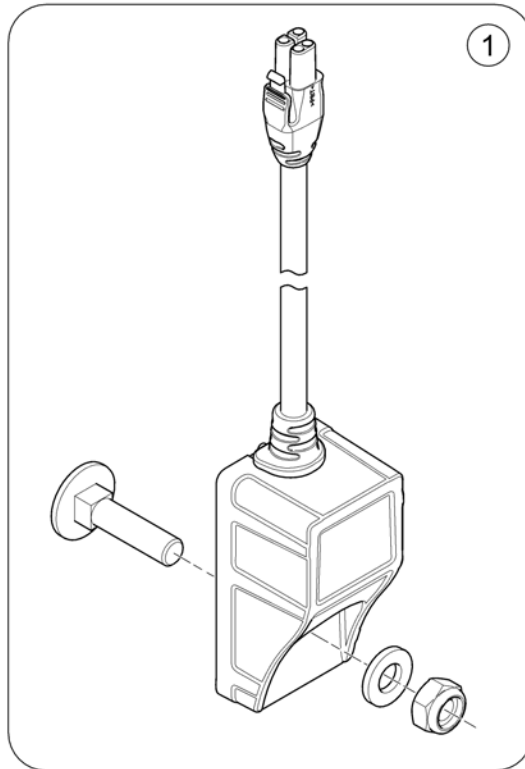1605414 iso



1 SP1603486

(LiNX G-Trac Sensor (Gyro) DLX-GYR100)

1

-	(User Manuals LiNX)	1
-	(Remotes & Brackets)	1
-	(Dual Control & Bracket)	1

(User Manuals LiNX)


languages.iso

BG	Български	FR	Français	NL	Nederlands	SL	Slovenščina
CS	Čeština	HE	עברית	NO	Norsk	SV	Svenska
DA	Dansk	HR	Hrvatski	PL	Polski	UK	Українська
DE	Deutsch	HU	Magyar	PT	Português		
EN	English	IT	Italiano	RU	Русский		
ES	Español	JA	日本語	RO	Română		
ET	Eesti	LV	Latviešu	SA	العربية		
FI	Suomi	LT	Lietuvių	SK	Slovenčina		

BG	SP1603635	(User Manual LiNX - Bulgarian)	1
CS	SP1603637	(User Manual LiNX - Czech)	1
DA	SP1603638	(User Manual LiNX - Danish)	1
DE	SP1603644	(User Manual LiNX - German)	1
EN	SP1602969	(User Manual LiNX - English)	1
ES	SP1603711	(User Manual LiNX - Spanish)	1
ET	SP1603641	(User Manual LiNX - Estonian)	1
FI	SP1603642	(User Manual LiNX - Finnish)	1
FR	SP1603643	(User Manual LiNX - French)	1
HR	SP1603636	(User Manual LiNX - Croatian)	1
HU	SP1603700	(User Manual LiNX - Hungarian)	1
IT	SP1603701	(User Manual LiNX - Italian)	1
LV	SP1603702	(User Manual LiNX - Latvian)	1
LT	SP1603703	(User Manual LiNX - Lithuanian)	1
NL	SP1603639	(User Manual LiNX - Dutch)	1
NO	SP1603704	(User Manual LiNX - Norwegian)	1
PL	SP1603705	(User Manual LiNX - Polish)	1
PT	SP1603706	(User Manual LiNX - Portuguese)	1
RU	SP1603707	(User Manual LiNX - Russian)	1
SK	SP1603799	(User Manual LiNX - Slovak)	1
SL	SP1603708	(User Manual LiNX - Slovenian)	1
SV	SP1603712	(User Manual LiNX - Swedisch)	1

(Remotes & Brackets)

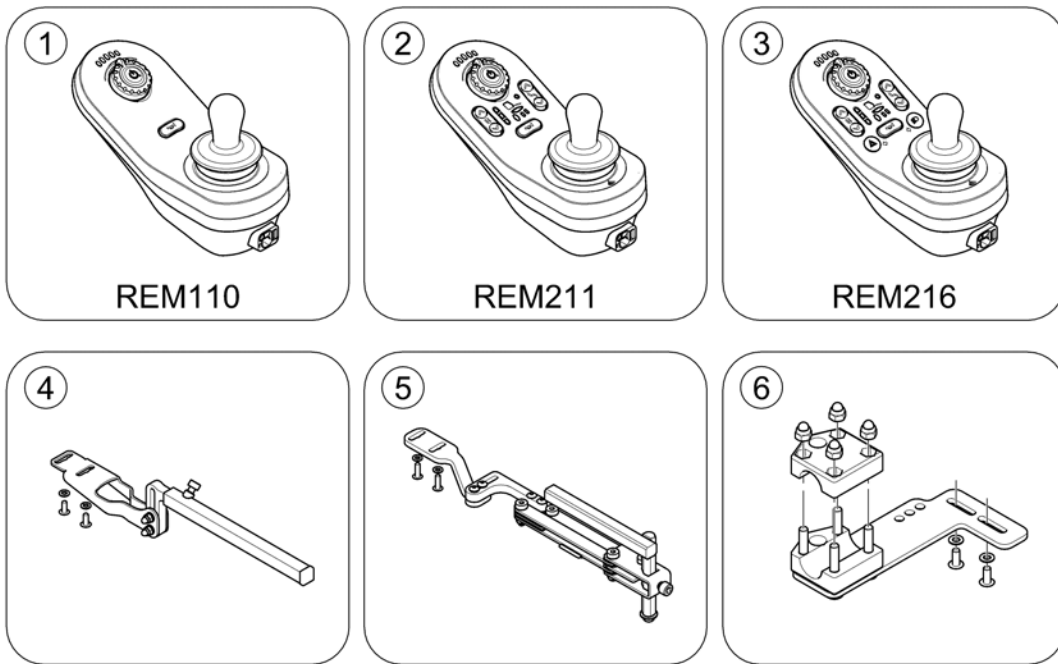
1605187.jpg




1	-	(Remote REM110)	1
2	-	(Remote REM211)	1
3	-	(Remote REM216)	1
4	-	(Remote Bracket, Standard)	1
5	-	(Remote Bracket, Swing away)	1
6	-	(Remote Bracket, Attendant)	1

(Remote REM110)

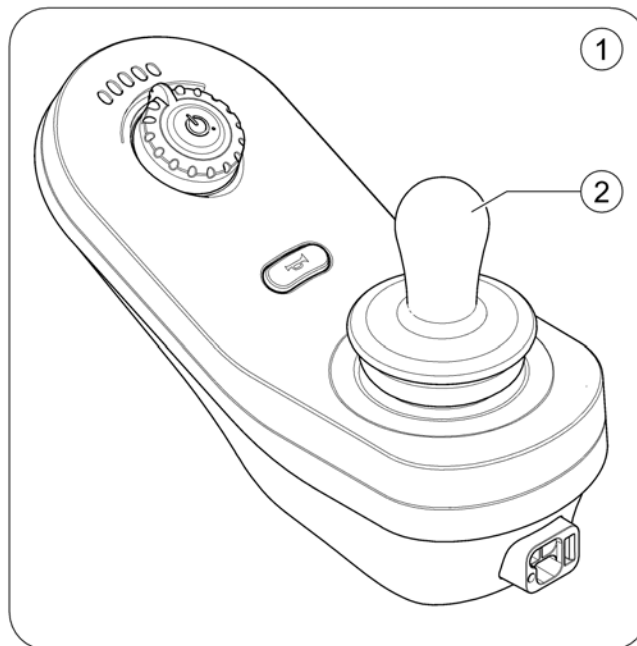1605135 iso




1	SP1602806	(LiNX Remote DLX-REM110 - Drive Only)	1
2	SP1607776	(Joystick Knob)	1

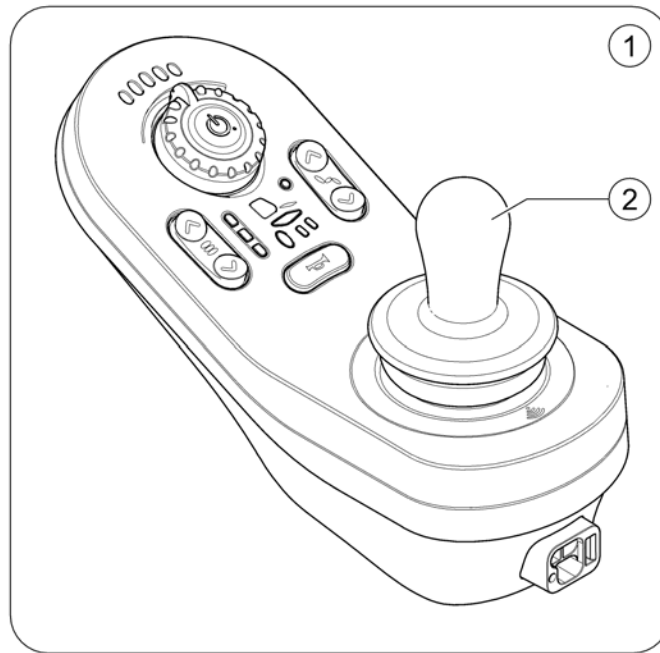


1	SP1602825	(LiNX Remote DLX-REM211 - Drive and Actuators)	1
2	SP1607776	(Joystick Knob)	1

(Remote REM216)

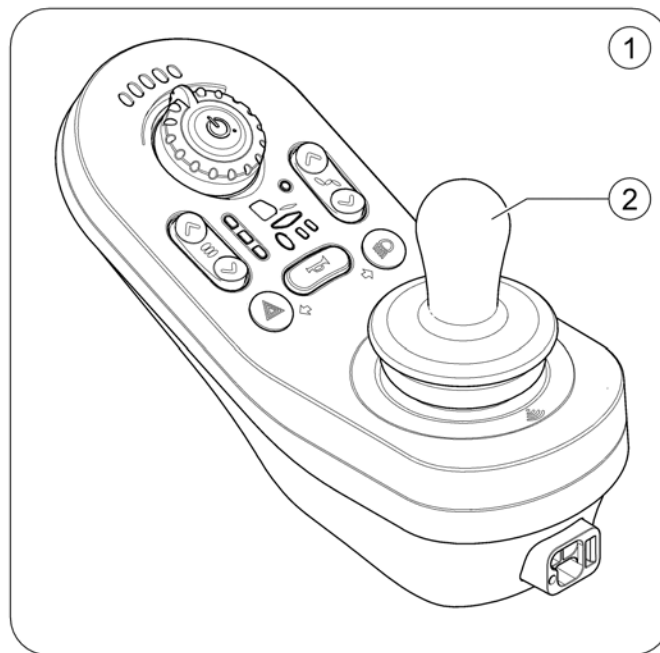1605137.indd




1	SP1602826	(LiNX Remote DLX-REM216 - Drive, Actuators and Lights)	1
2	SP1607776	(Joystick Knob)	1

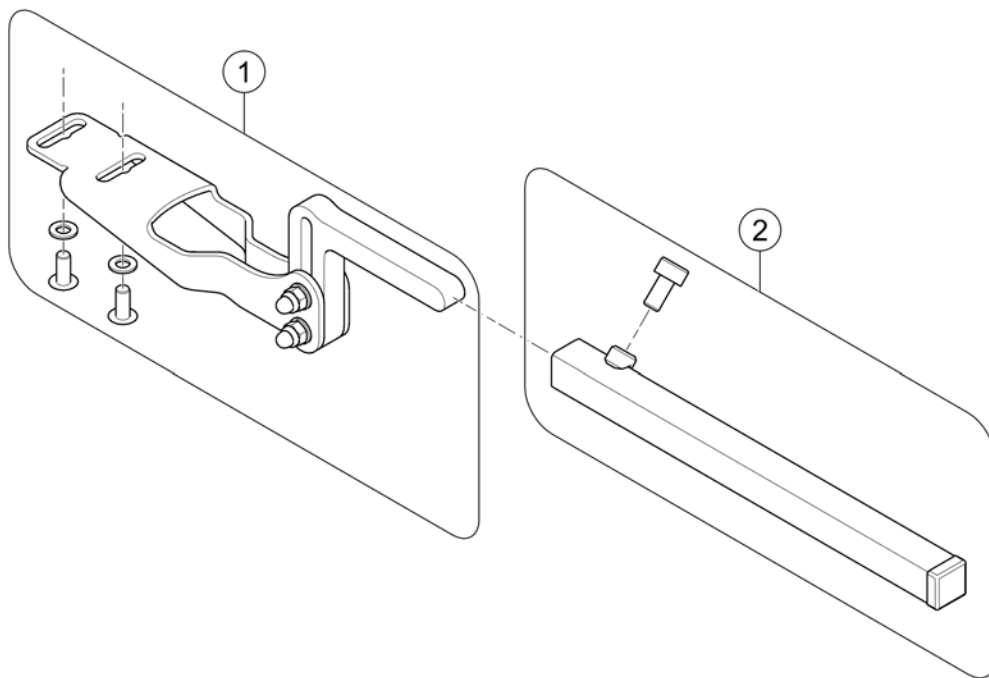


1	SP1603417	(LiNX Remote bracket, Standard)	1
2	SP1603419	(Remote support tube)	1

(Remote Bracket, Swing away)

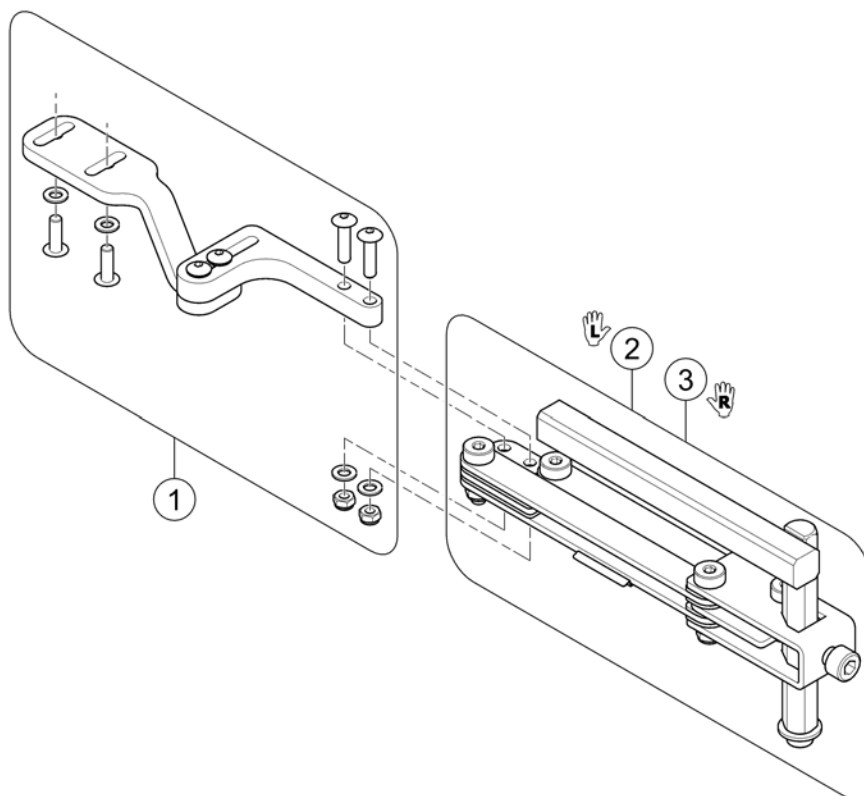1605139 iso





1	SP1603418	(LiNX Remote support for height adjustable & swing-away bracket)	1
2	SP1496249	SWING AWAY VE NY MOD2008 KPL	1
3	SP1496248	SWING AWAY HØ NY MOD2008 KPL	1

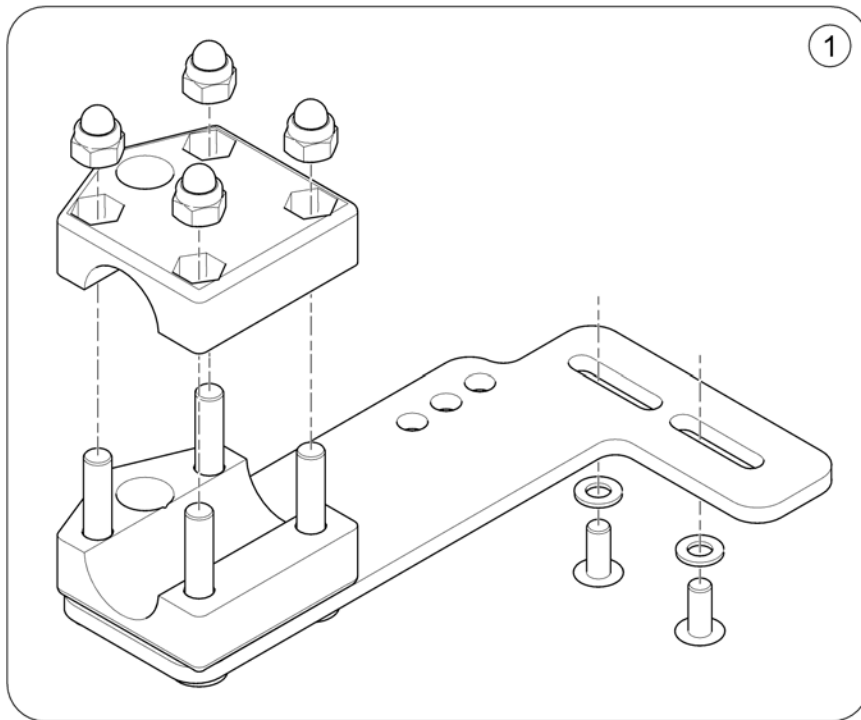


1 SP1603363

(LiNX Remote bracket, Attendant)

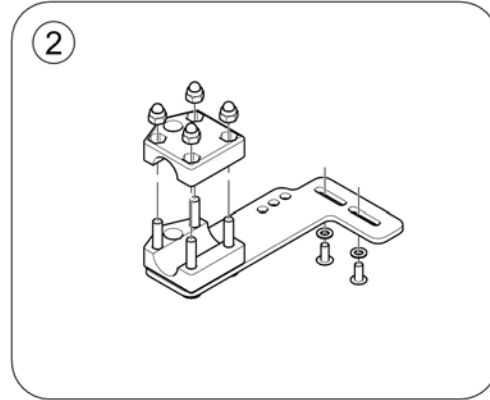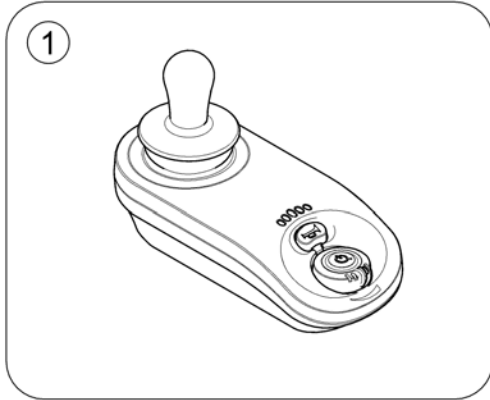1



1	-	(Dual Control)	1
2	-	(Remote Bracket, Attendant)	1

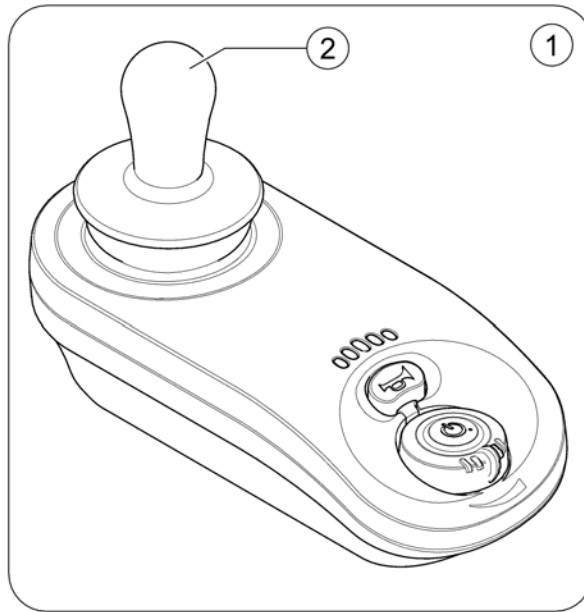

1 SP1602792

(LiNX Dual Control Remote DLX-REM050)

1



(Cables & Accessories)

-	(Bus Cables)	1
-	(Light Cable)	1
-	(Light Board)	1
-	(Lifter/Tilt Module Cables)	1
-	(Actuator Cable)	1
-	(4-way Expansion Block)	1

(Bus Cables)

LiNX

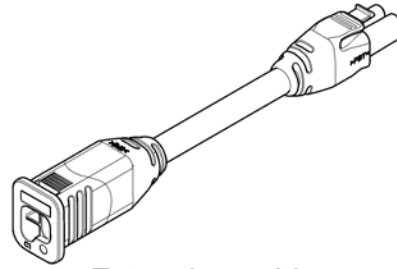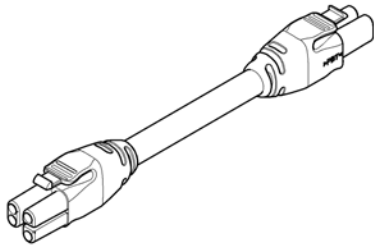1605413.jpg



Bus cable

Lenght [mm]

- ① 300
- ② 500
- ③ 700
- ④ 1000
- ⑤ 1200
- ⑥ 1500
- ⑦ 1700
- ⑧ 2000
- ⑨ 2500

Extension cable


Lenght [mm]

- ⑩ 640
- ⑪ 900

1	SP1602888	(LiNX Bus cable 300mm (GLM-BUS030))	1
2	SP1602889	(LiNX Bus cable 500mm (GLM-BUS050))	1
3	SP1602890	(LiNX Bus cable 700mm (GLM-BUS070))	1
4	SP1602891	(LiNX Bus cable 1000mm (GLM-BUS100))	1
5	SP1603489	(LiNX Bus cable 1200mm (GLM-BUS120))	1
6	SP1602892	(LiNX Bus cable 1500mm (GLM-BUS150))	1
7	SP1602893	(LiNX Bus cable 1700mm (GLM-BUS170))	1
8	SP1602894	(LiNX Bus cable 2000mm (GLM-BUS200))	1
9	SP1602904	(LiNX Bus cable 2500mm (GLM-BUS250))	1
10	SP1602905	(LiNX Extension cable 640mm (GLM-EXT064))	1
11	SP1602906	(LiNX Extension cable 900mm (GLM-EXT090))	1

(Light Cable)

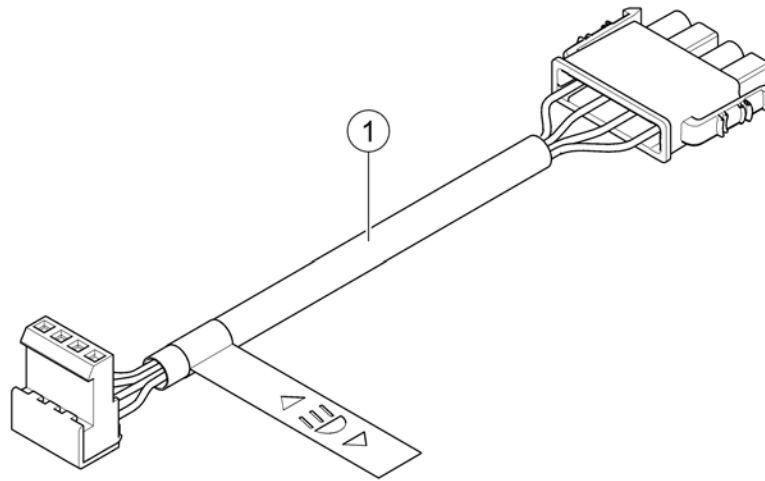
160623 iso



1 SP1603336

(LiNX Light cable)

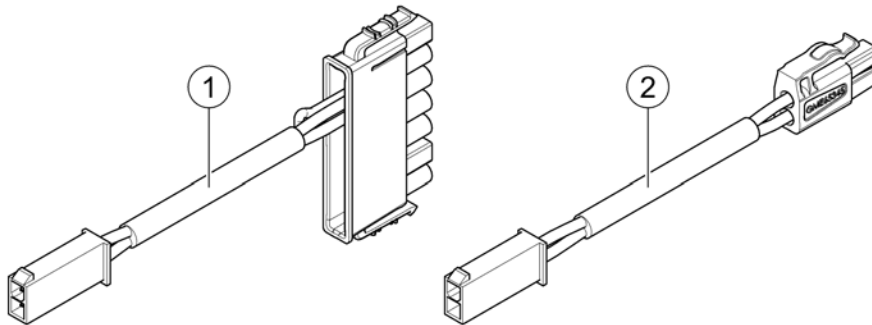
1



1	SP1604507	(LiNX Sensor Cable Lifter/Tilt Module - 300mm)	1
2	SP1576062	(Extension Cable - 1100mm)	1

(Actuator Cable)

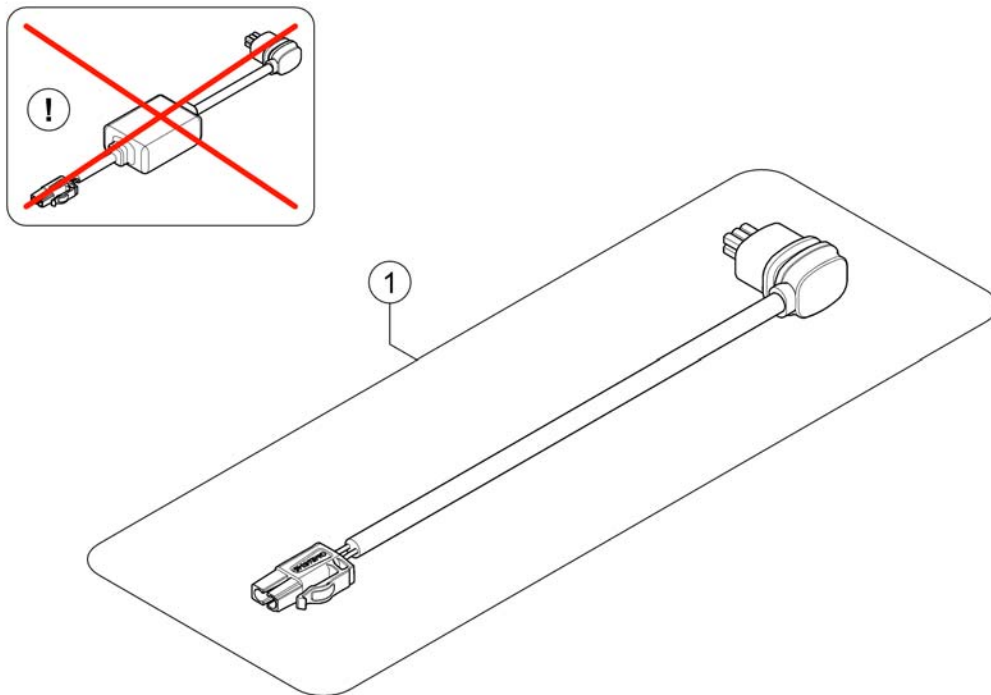
1505463 iso



!	-	(No longer available - please refer to Replacement Info for which part is to be used instead.)	0
1	SP1607105	(LiNX Actuator cable LA23 - 950mm) <i>Replaces SP1605680</i>	1

(4-way Expansion Block)

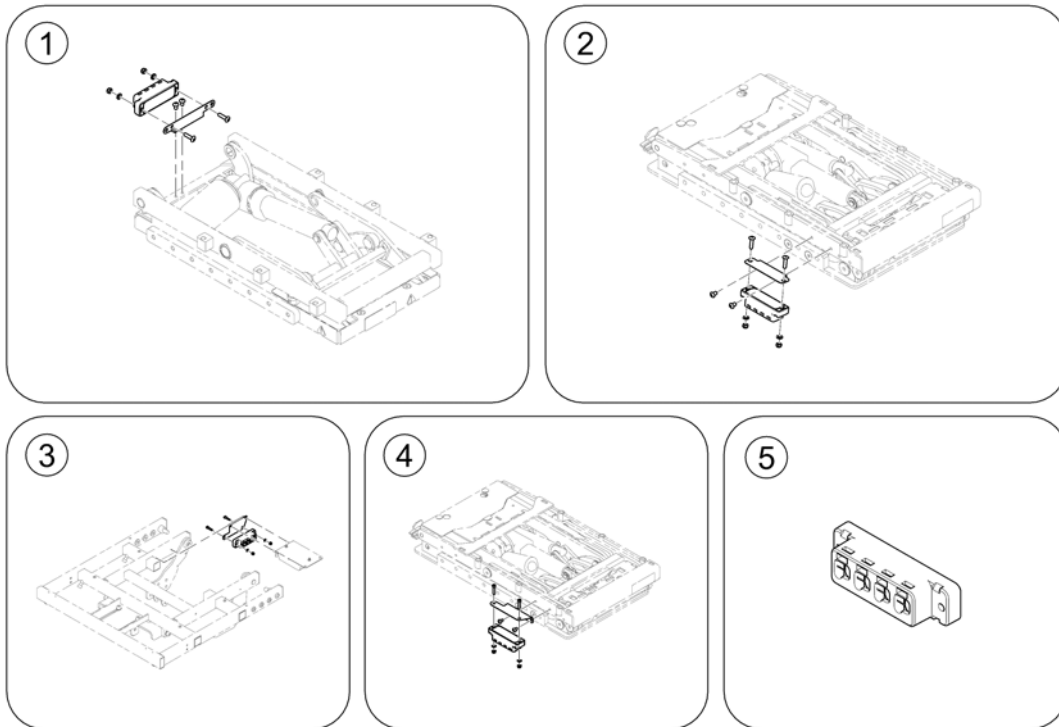
1506721 no



1	-	(4-way Expansion Block (with Tilt module))	1
2	-	(4-way Expansion Block (with Lifter/Tilt module))	1
3	-	(4-way Expansion Block (Max Seat))	1
4	-	(4-way Expansion Block (Low Rider))	1
5	-	(4-way Expansion Block (without any module))	1

(4-way Expansion Block (with Tilt module))

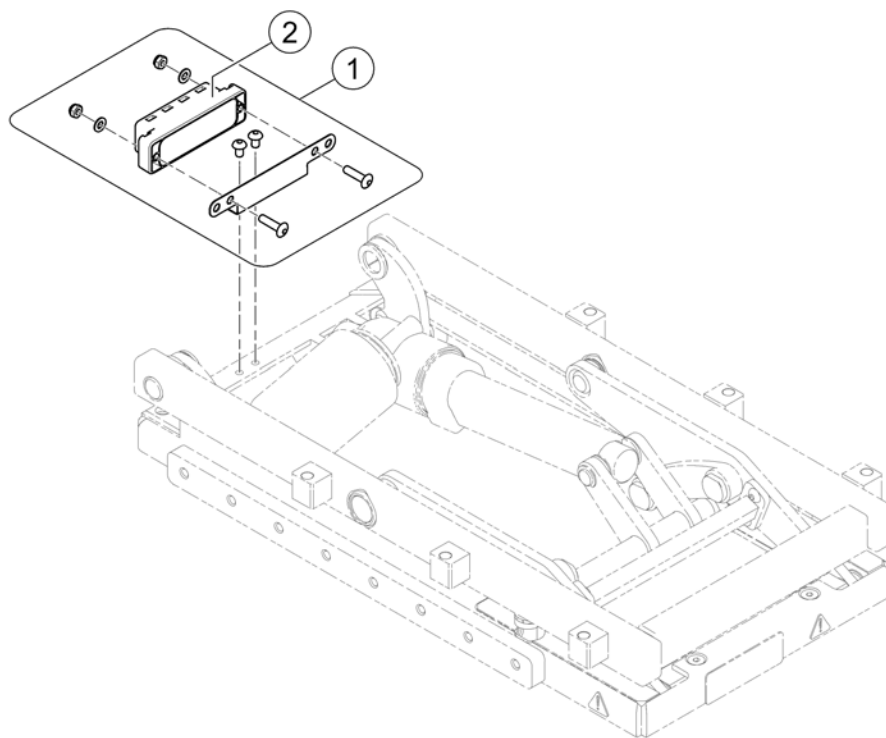
160673 iso



1	SP1604608	(LiNX 4-way Expansion Block with mounting hardware for Tilt module)	1
2	SP1602907	(4-way Expansion Block, LiNX, GLM-CONX4 (single part))	1

(4-way Expansion Block (with Lifter/Tilt module))

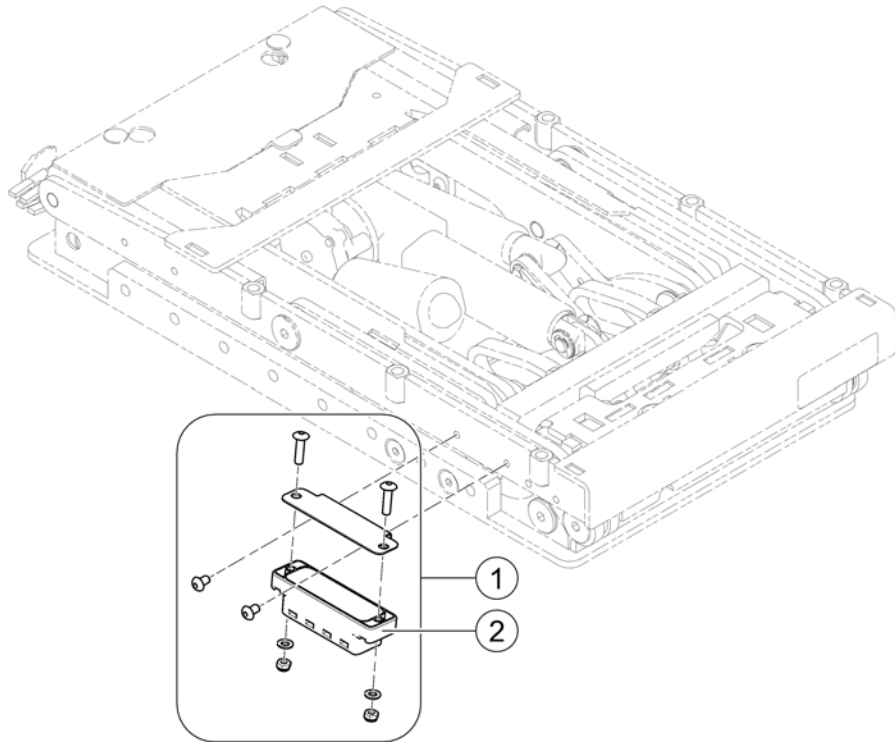
1604612.jpg



1	SP1603350	(LiNX 4-way Expansion Block with mounting hardware for Lifter/Tilt module)	1
2	SP1602907	(4-way Expansion Block, LiNX, GLM-CONX4 (single part))	1


(4-way Expansion Block (Max Seat))

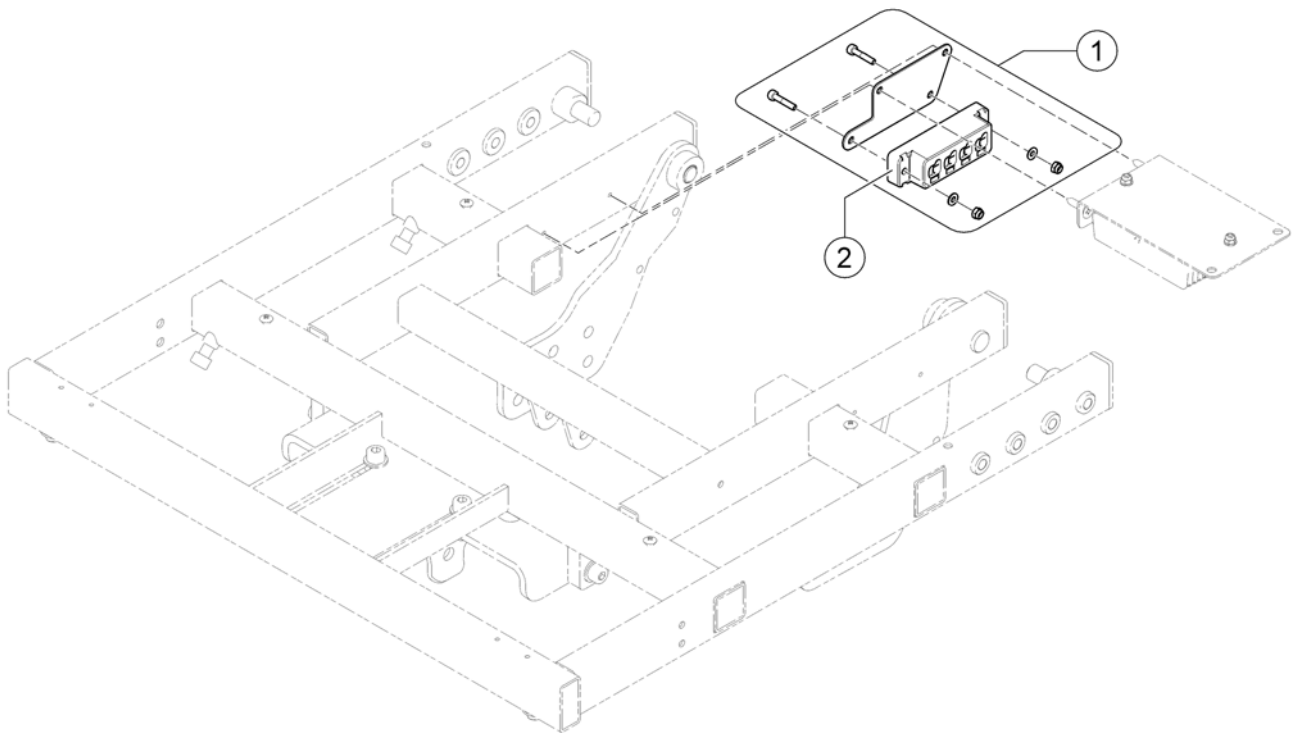1606/20 no



1	SP1605658	(4-way Expansion Block, Max Seat, cpl.)	1
2	SP1602907	(4-way Expansion Block, LiNX, GLM-CONX4 (single part))	1


Invacare® TDX-SP2 Lowrider, TDX-SP2 NB | Elektronik
(4-way Expansion Block (Low Rider))

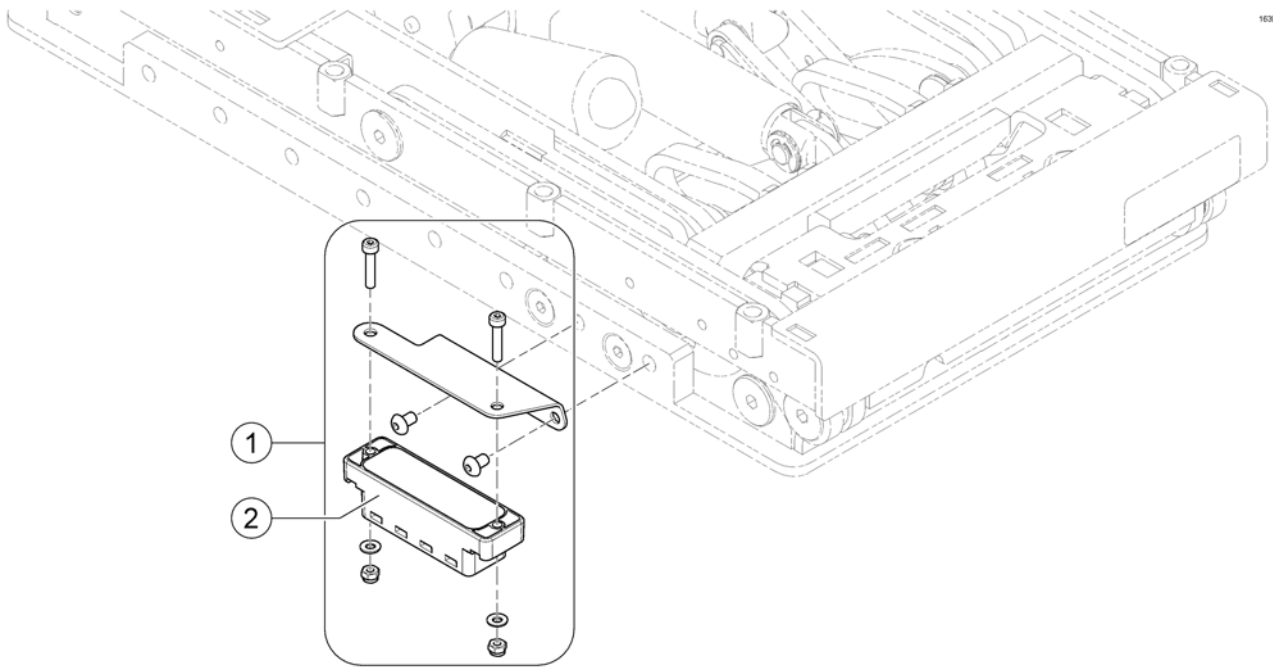

1	SP1605134	(4-way Expansion Block, LiNX, cpl. (Low Rider only!))	1
2	SP1602907	(4-way Expansion Block, LiNX, GLM-CONX4 (single part))	1

(4-way Expansion Block (without any module))

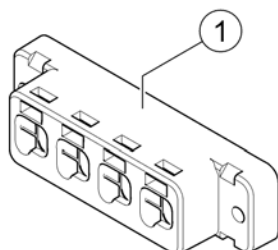160655 iso







1 SP1602907

(4-way Expansion Block, LiNX, GLM-CONX4 (single part))

1

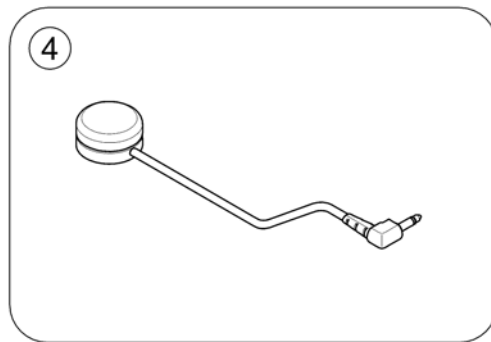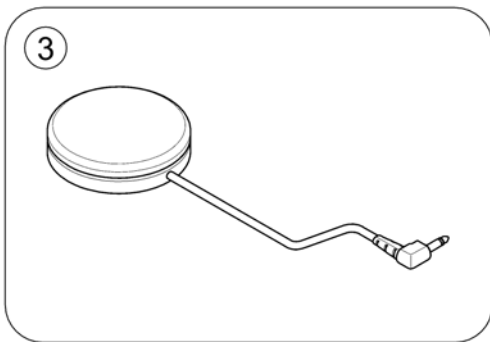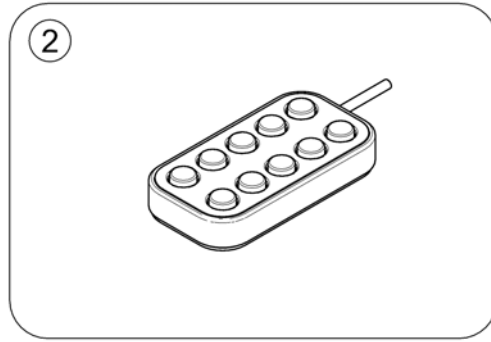
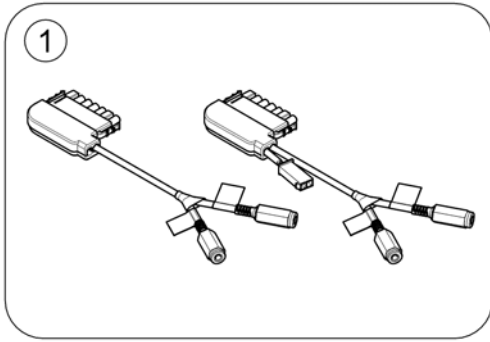
1028154.no



1	-	(Accessory cables)	1
2	-	(10-way switch & brackets)	1
3	-	(Buddy Button)	1
4	-	(Piko Button)	1





(Accessory cables)

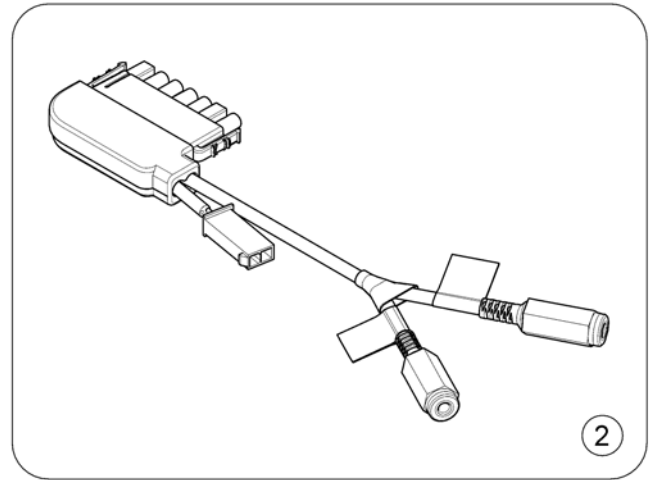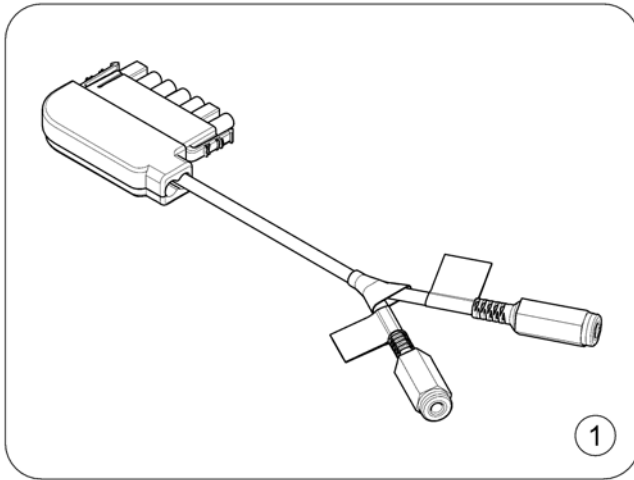1608109-00




1	SP1637339	(Accessory cable LiNX (w/o lifter))	1
2	SP1608269	(Accessory cable LiNX (with lifter))	1

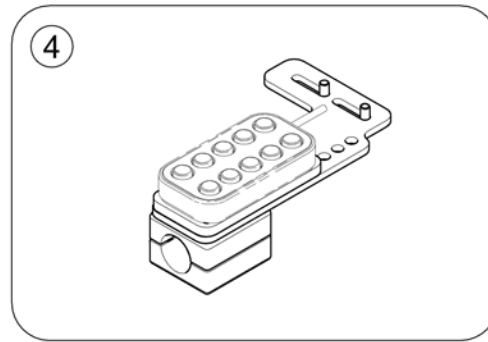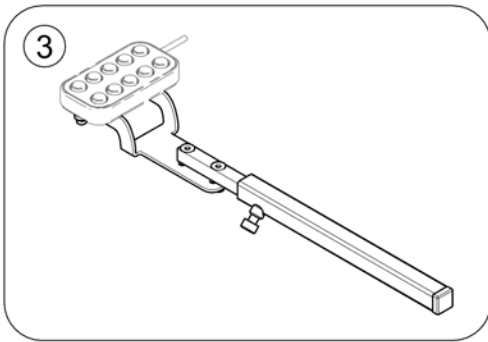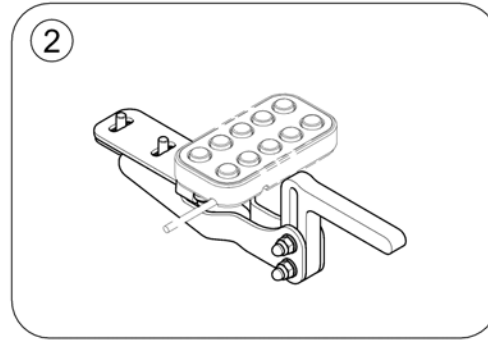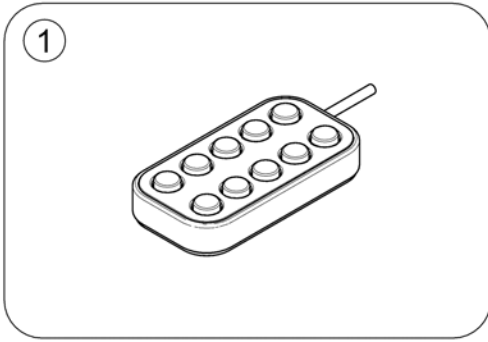1608155.kim



1	-	(10-way switch)	1
2	-	(Bracket, front)	1
3	-	(Bracket, opposite side)	1
4	-	(Bracket, rear)	1

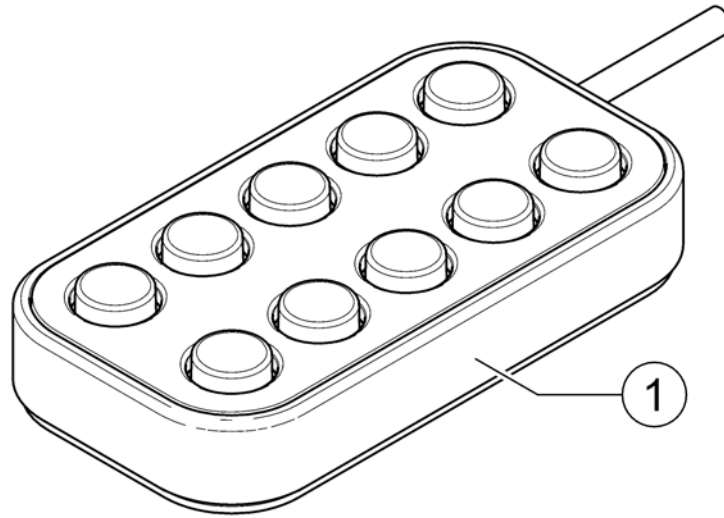(10-way switch)




1 SP1602838

(LiNX 10-way switch (DLX-SW100))

1

(Bracket, front)

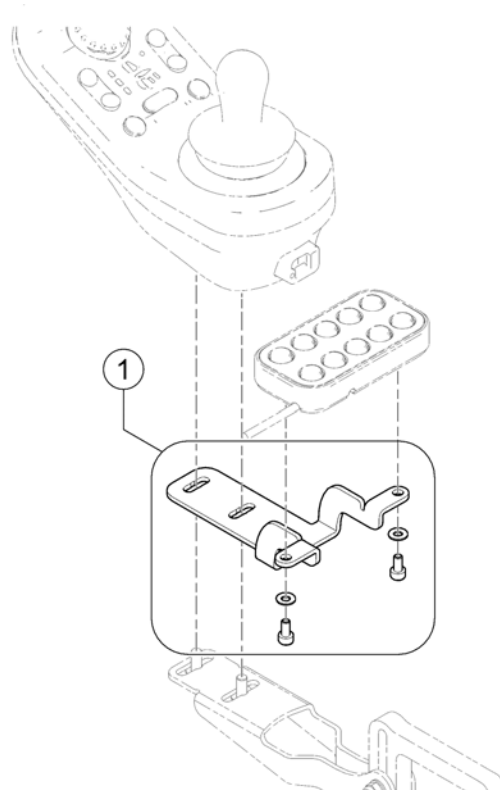1008157.no




1 SP1606388

(Bracket 10-way switch - front)

1

(Bracket, opposite side)

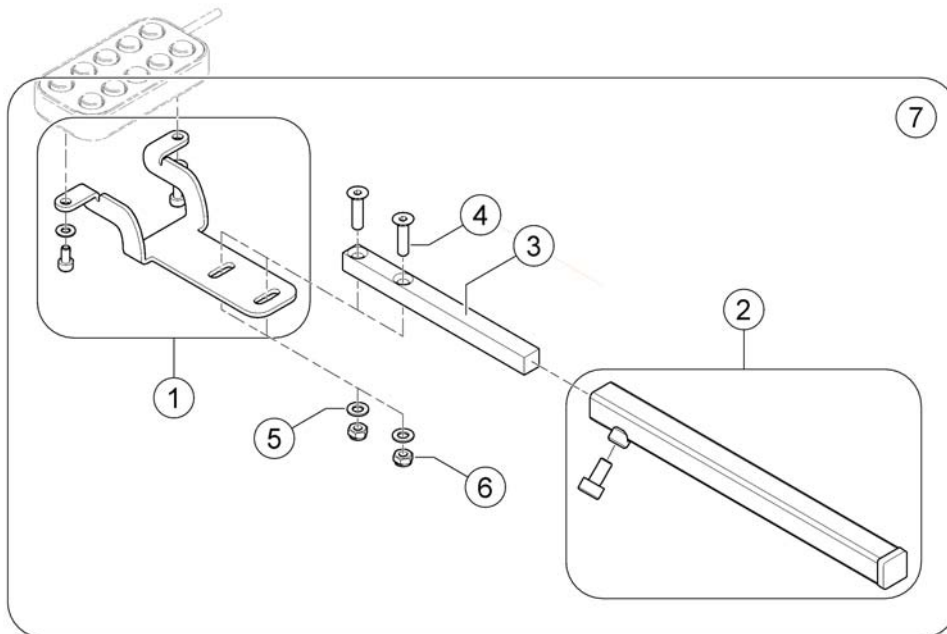
1639748-00



1	SP1606388	(Bracket 10-way switch - front)	1
2	SP1603419	(Remote support tube)	1
3	1440437-9005	(Joystick support Standard)	1
4	1417938	(Countersunk bolt DIN7991-M05x20-8.8 zc)	2
5	1422890	(Washer DIN125-Ø5,3 bk)	2
6	F20111	(Hexagon nut DIN985-M05 bk)	2
7	SP1606371	(Bracket, opposite side, cpl.)	1

(Bracket, rear)

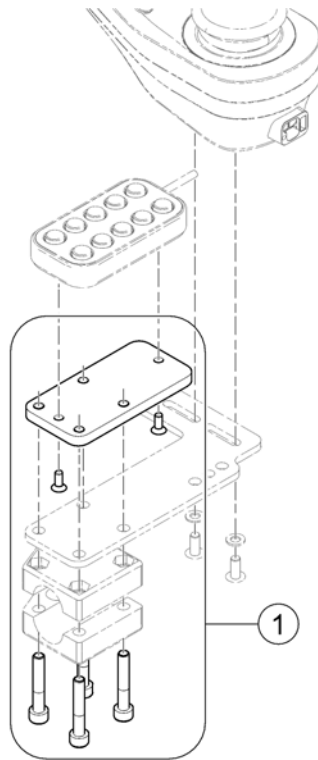1608165.00




1 SP1606389

(Bracket 10-way switch - rear)





1

(Electronic - Narrow Base)

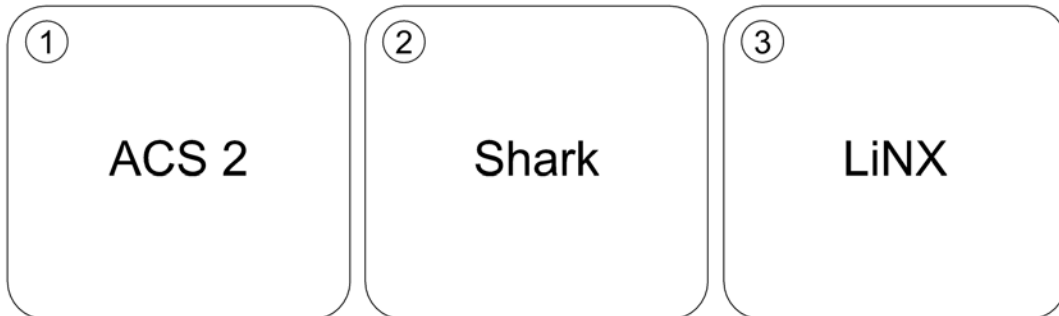1583351.no



1	-	(ACS 2 Electronic)	1
2	-	(Shark Electronic)	1
3	-	(LiNX Electronic)	1

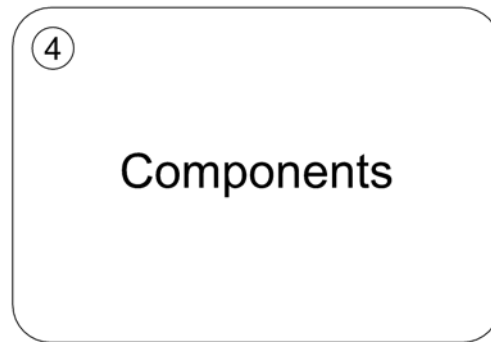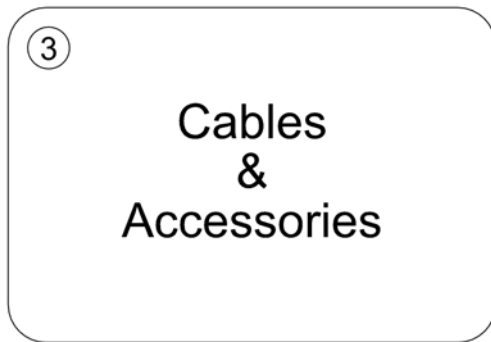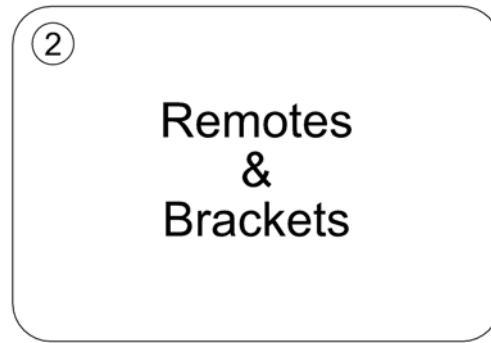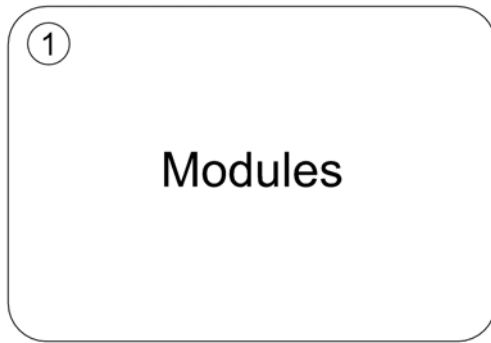


1	-	(ACS Electronic Modules)	1
2	-	(ACS Remotes & Brackets)	1
3	-	(ACS Cables & Accessories)	1
4	-	(ACS Components)	1

(ACS Electronic Modules)

-	(Power Modul - ACS)	1
-	(Bracket for Power Modul - Narrow Base)	1
-	(ACTs)	1
-	(Bracket ACT - Modulite Seat)	1
-	(Bracket ACT - Ultra Low Maxx Seat)	1
-	(Bracket ACT at Chassis - Narrow Base)	1
-	(Light Board)	1
-	(G-Trac Sensor/Gyro)	1

Invacare® TDX-SP2 Lowrider, TDX-SP2 NB | Elektronik
(Bracket for Power Modul - Narrow Base)

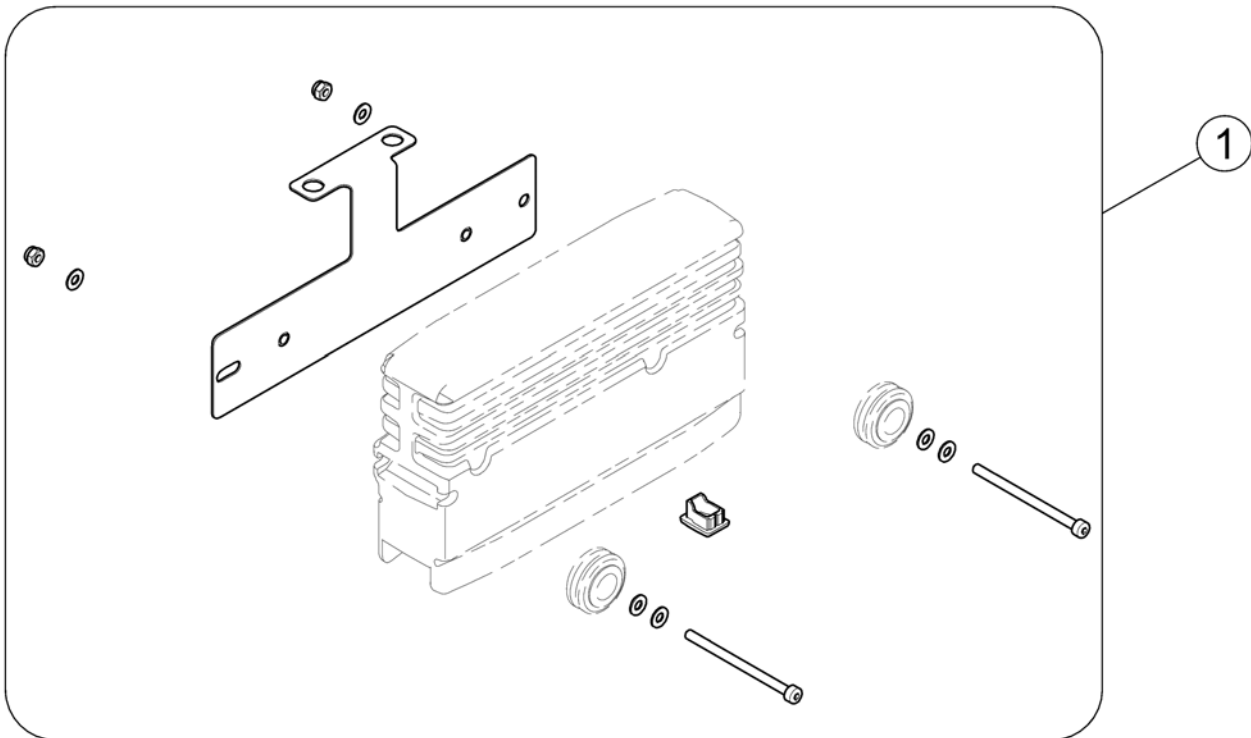1583338.iso




1 SP1583339

(Bracket Power Modul)

1

(Bracket ACT at Chassis - Narrow Base)

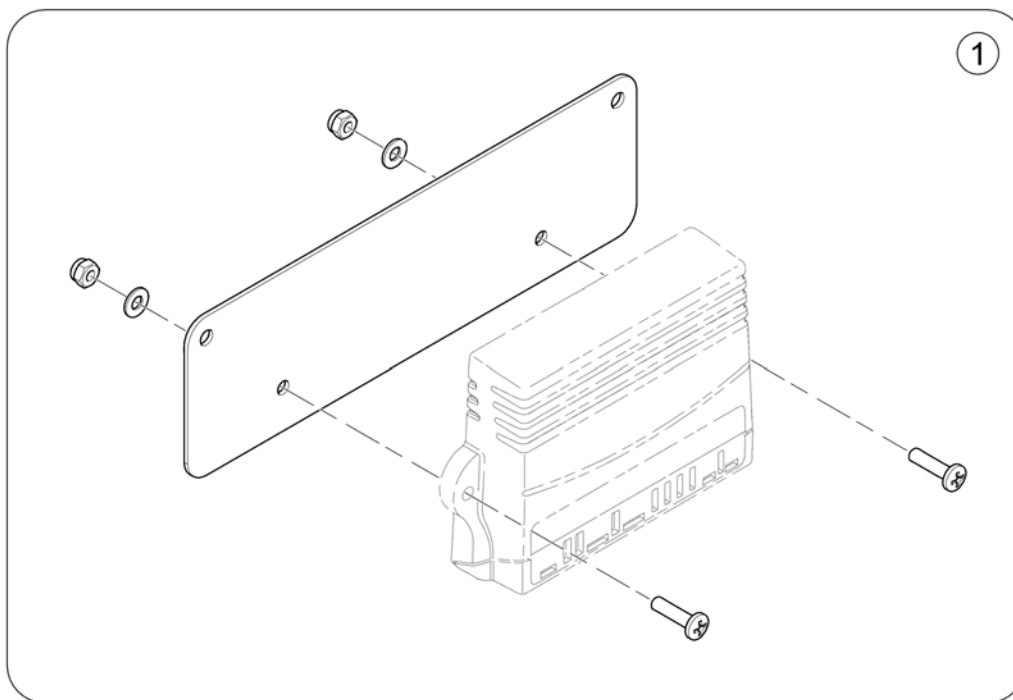
1583341.lin



1 SP1583342


(Bracket for ACT on Chassis (Narrow Base))

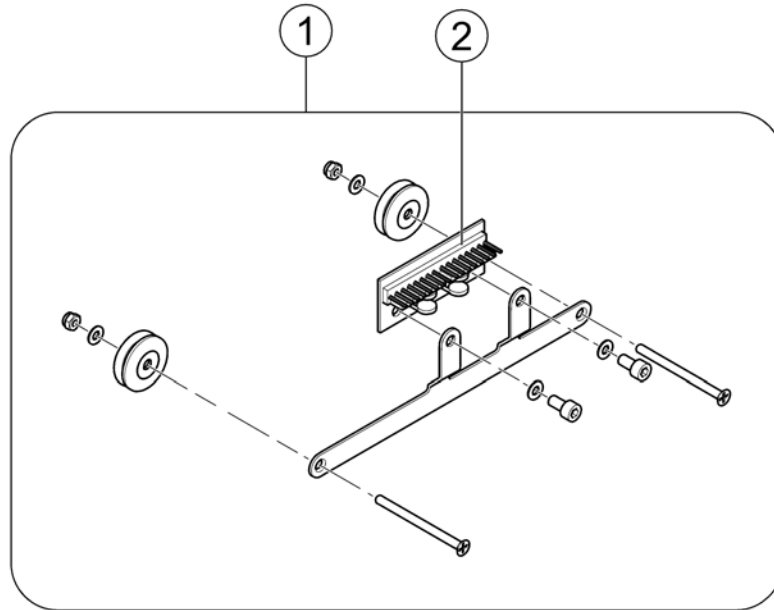1



1	SP1574115	(Light board with bracket)	1
2	1415004	LYSFORDELER PRINT G40	1


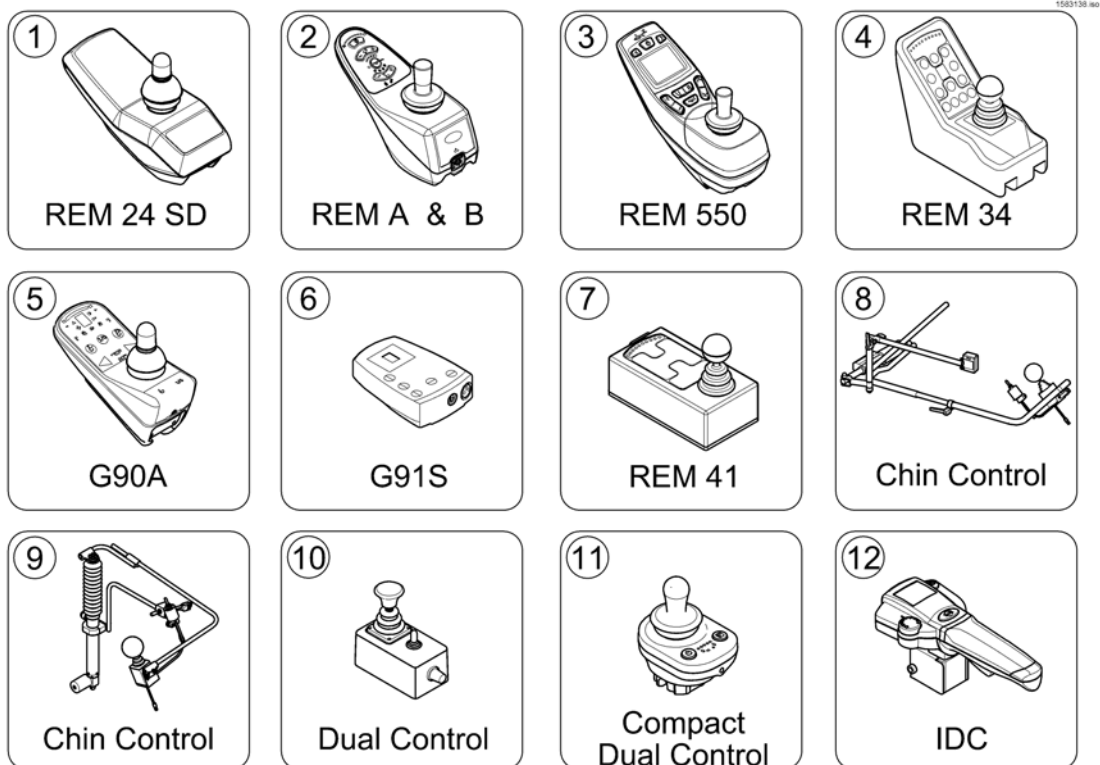
(ACS Remotes & Brackets)



1	-	(REM 24 SD with brackets) [1601/1602]	1
2	-	(REM A & B with brackets) [1638/1639/1643/1644]	1
3	-	(REM 550 with brackets) [1617/1618]	1
4	-	(REM 34 with brackets) [1611]	1
5	-	(G90A with brackets) [1614/1616]	1
6	-	(G91S with brackets) [1615]	1
7	-	(Centre control REM 41) [1607/1608]	1
8	-	(Chin Control, small (=> 16/12/31)) [1609]	1
9	-	(Chin Control - Modulite Seat (15/10/01 =>)) [1719]	1
10	-	(Dual Control with brackets - Modulite Seat) [1291/1292]	1
11	-	(Compact Dual Control with brackets) [1720]	1
12	-	(Intuitive Dual Control (IDC) with brackets) [1645]	1

(Dual Control with brackets - Modulite Seat)

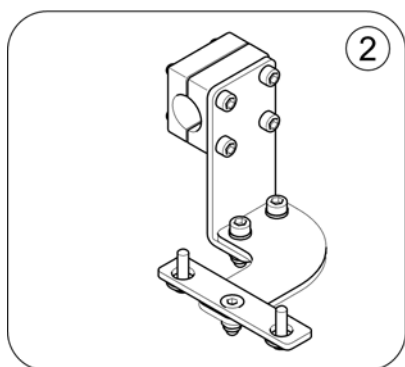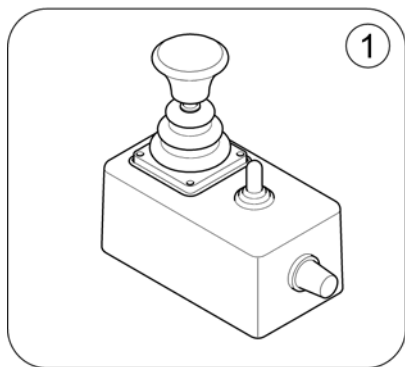
1583353.no



1	6001425	(ACS dual control) [1291/1292]	1
2	-	(Remote Bracket Dual Control for Modulite Seat)	1

(ACS Cables & Accessories)




-	(Bus-Cables ACS)	0
-	(Light Cable ACS)	0
-	(T-Adapter - Modulite Seat)	0
-	(T-Adapter (Ultra Low Maxx Seat))	0
-	(Tilt Sensors - Ultra Low Maxx Seat)	0
-	(Programming Kits)	0
-	(Magnetkey for Remotes)	0

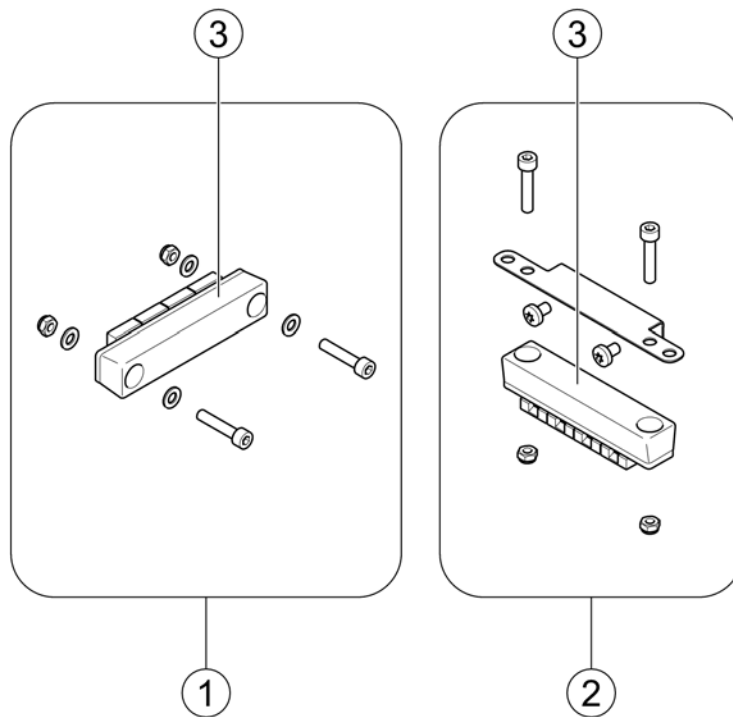

1	SP1571826	(T-Adapter - Modulite Seat)	1
2	SP1569958	(Bracket T-Adapter - Modulite Seat with MBL-Modul)	1
3	P10072	(Bus 4-way socket (T-Adapter))	1

(Shark Electronic)

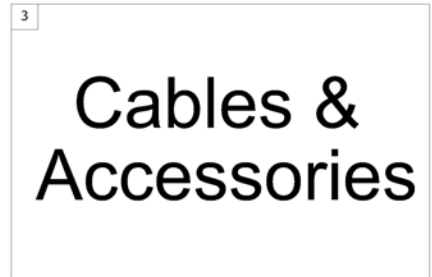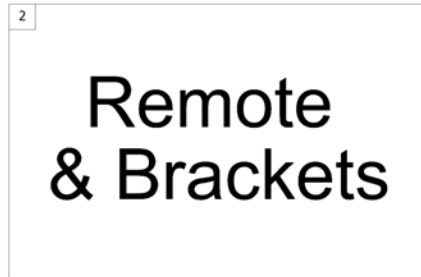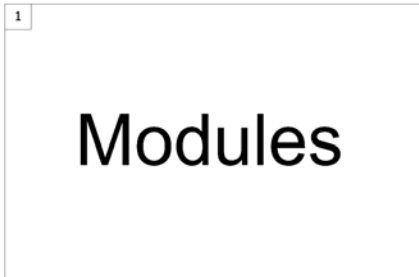1500784.no




1	-	(Shark Electronic Modules)	1
2	-	(Shark Remote & Brackets)	1
3	-	(Shark Cables & Accessories)	1

-	(Power Modul)	1
-	(Bracket for Power Modul)	1
-	(Light Board)	1

(Power Modul)

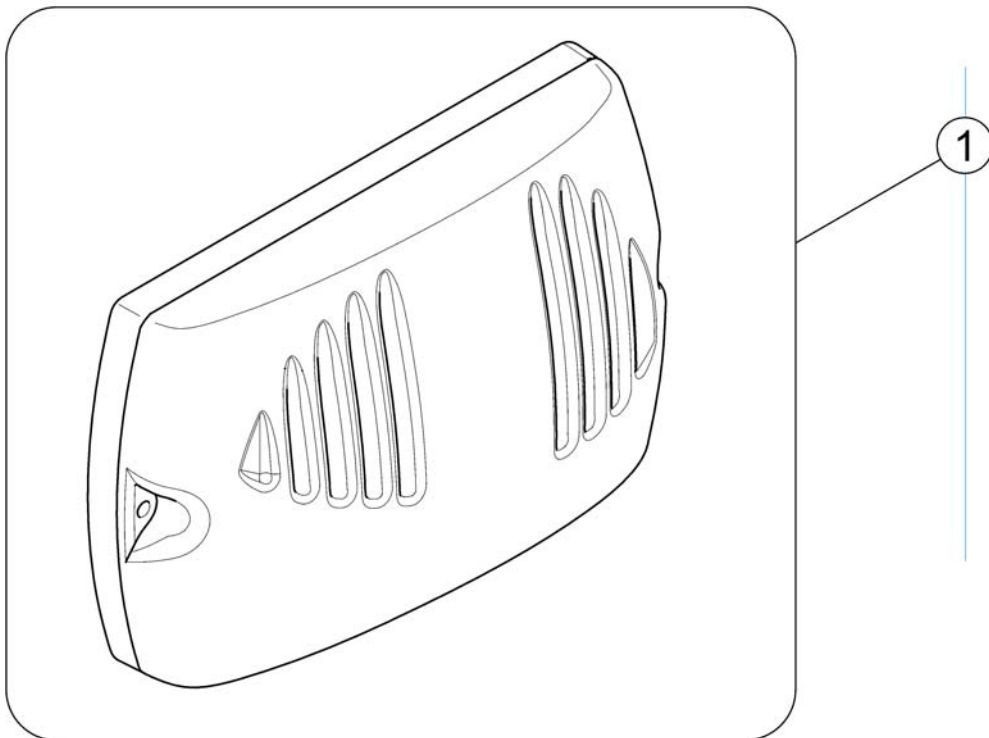
1583356.iso



1 1136898

(Powermodule, Shark, 90A, 2 actuators, DK-PMC28)

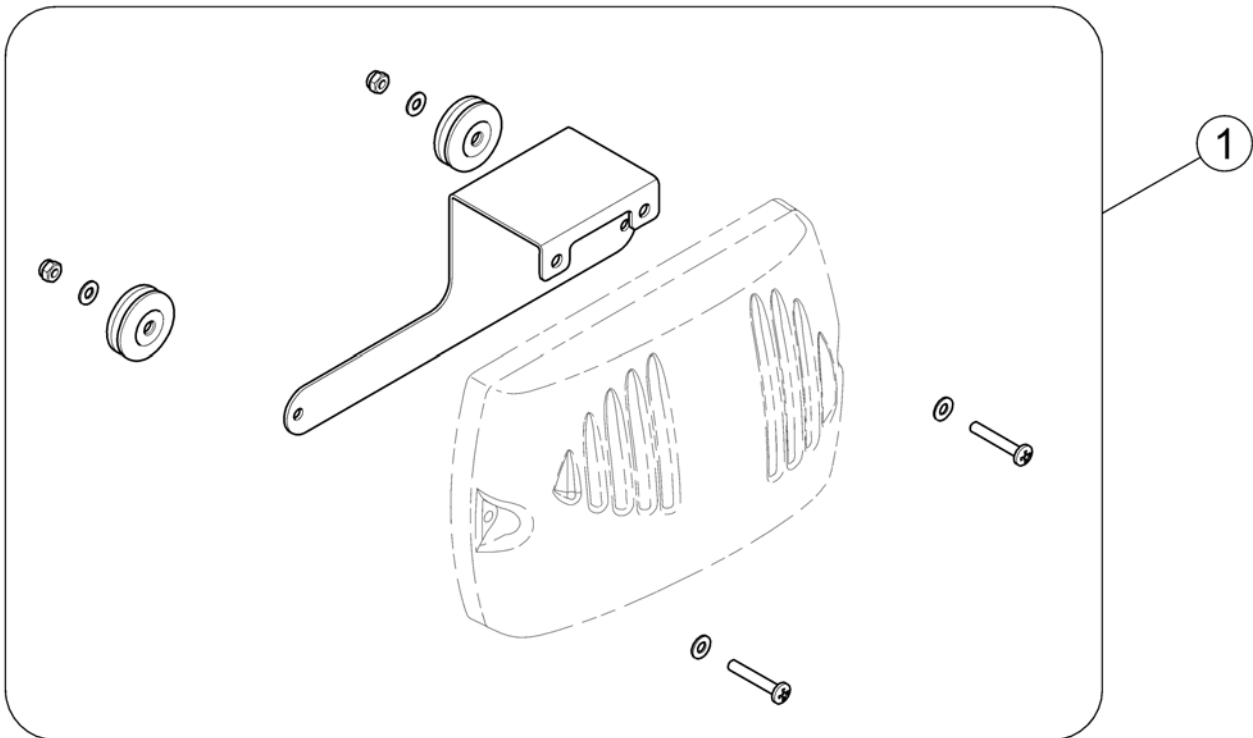1




1 SP1583358

(Bracket for Shark Power Modul)


1

(Shark Remote & Brackets)

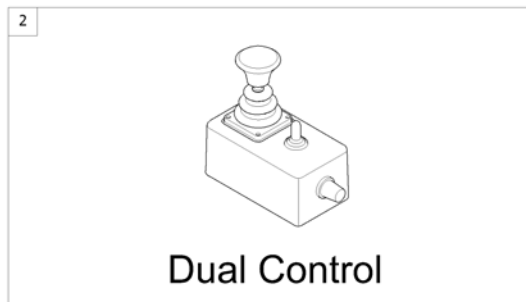156025 no





1	-	(Remote Shark with Brackets) <i>[1801]</i>	1
2	-	(Remote Dual Control with Brackets) <i>[1891/1892]</i>	1

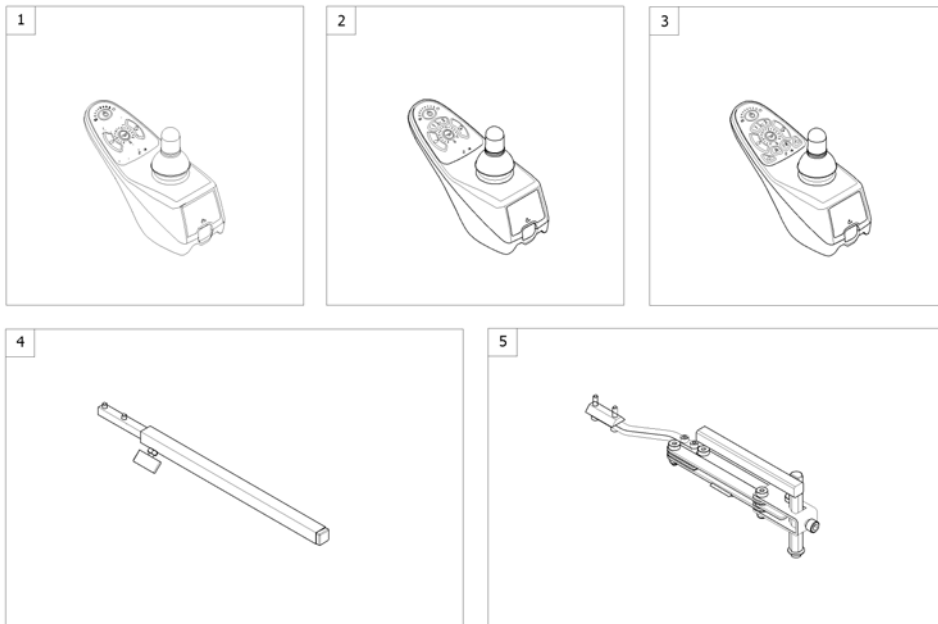


1	1492567	STYREBOKS SHARK DKREMD01 DRAGON [1801]	1
2	1449177	STYREBOKS SHARK II DX-REMD11 DRAGON [1801]	1
3	1484904	STYREBOKS SHARK II DX-REMD11 LYS OG EL DRAG [1801]	1
4	-	(Remote Bracket Standard) [1801+1207/1208]	1
5	-	(Remote Bracket Swing away & height adjust.) [1801+1707/1708]	1

(Remote Dual Control with Brackets)

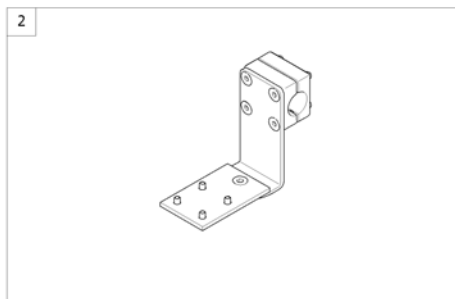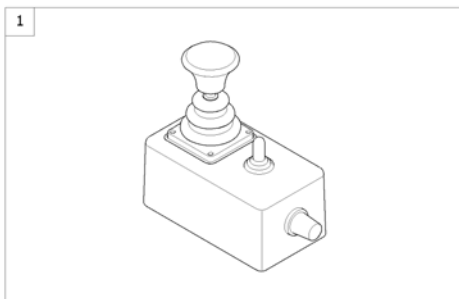1563329.iss



1	1444717	(Shark Dual Control, switchable) <i>[1891/1892]</i>	1
2	-	(Remote Bracket Dual Control) <i>[1891/1892]</i>	1

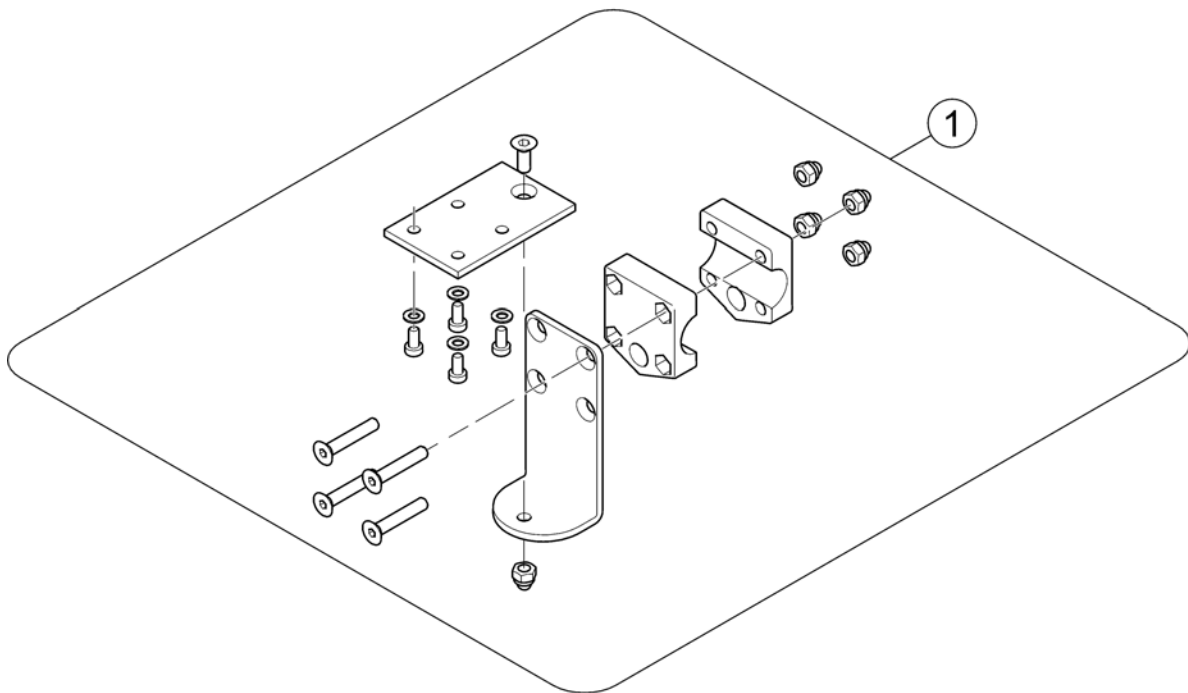

1 1450307

(Remote Bracket Attendent Control Shark)


1

(Shark Cables & Accessories)

-	(Bus-Cables Shark)	0
-	(Cables Shark)	0
-	(Drive Lockout)	0
-	(Programming Kits)	0

(Bus-Cables Shark)

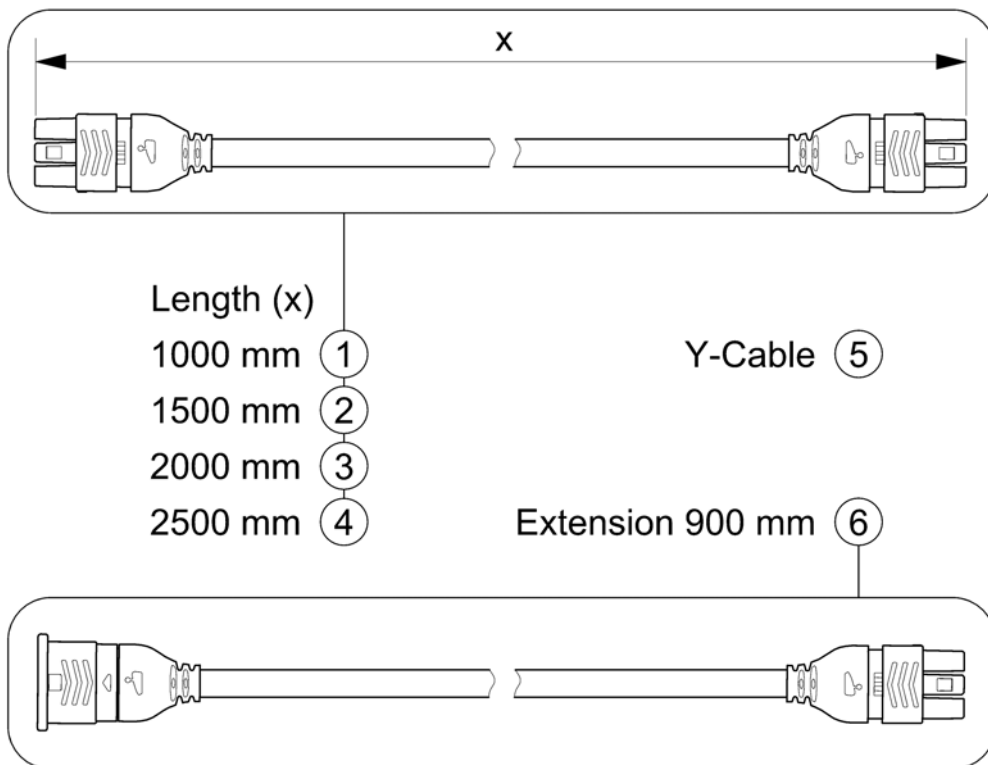
SHARKBUSKABEL.00



1	1448922	(Bus Cable Shark 1000mm)	1
2	1450313	BUSKABEL 1500 MM SHARK	1
3	1484911	BUSKABEL SHARK 2M SHARK	1
4	SP1510073	(Buscable Shark 2500mm)	1
5	1450314	BUSKABEL Y SHARK	1
6	1116404	(Extension Cable Shark 900mm)	1

(Cables Shark)


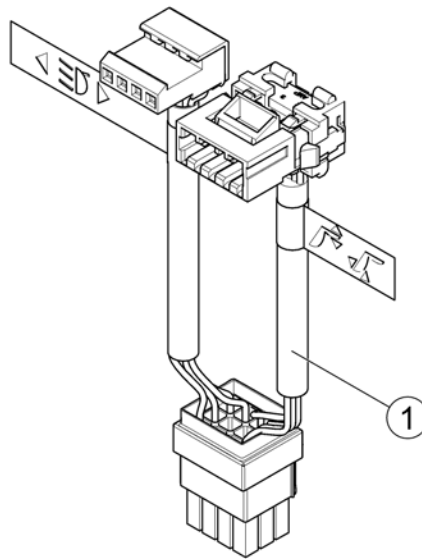
1583705.lin



1 1488951

KABELSETT KPL FOR LYS OG AKTUATOR DRAGON

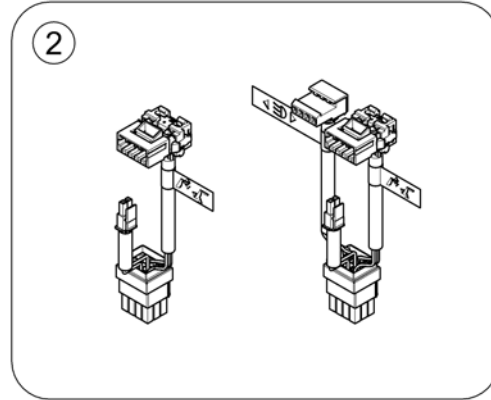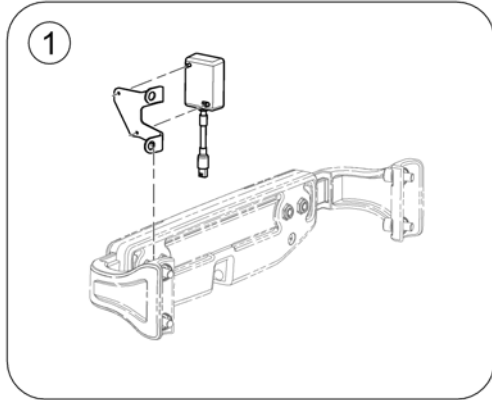
1



1	-	(Drive Lockout with bracket)	1
2	-	(Drive Lockout Cables - Shark)	1

(Drive Lockout with bracket)

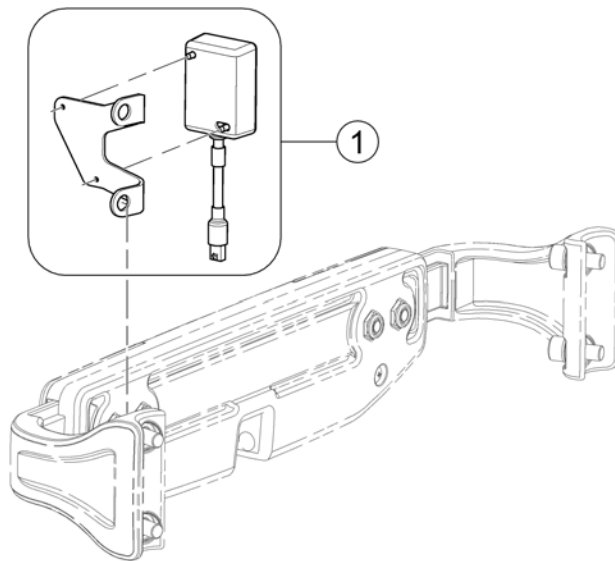
1583877.00







1 SP157707

(Drive Lockout with bracket)
Shark+[0977] / ACS+[2100]+[1092/1093/1094]

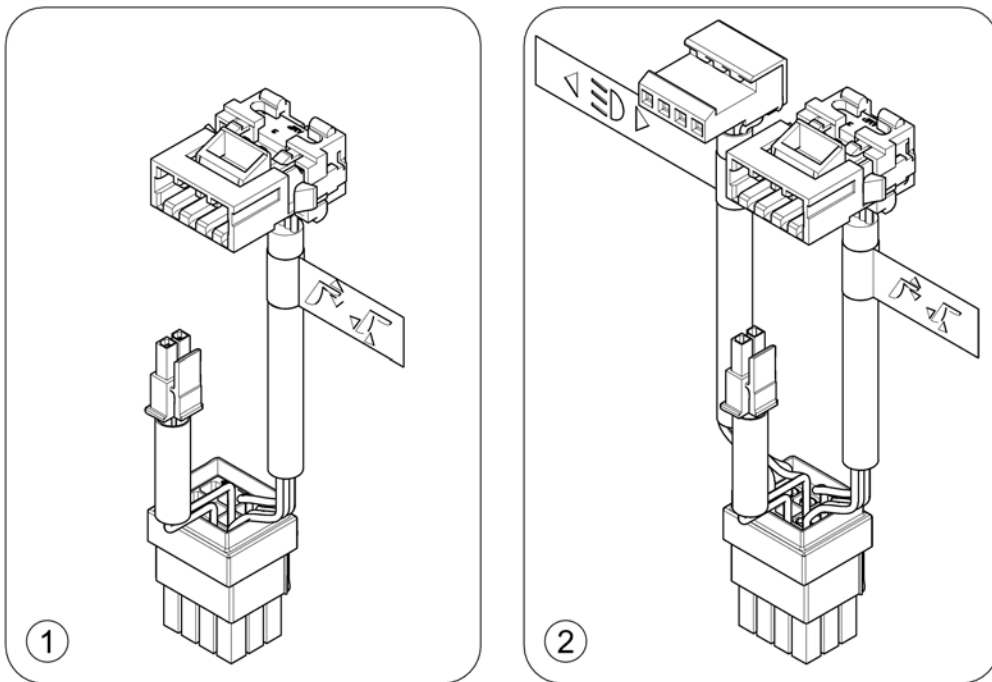1



1	SP1577326	(Wire Harness - Drive Lockout & el. Adjustments) [0977]+[1254]	1
2	SP1577325	(Wire Harness - Drive Lockout, el. Adjustments & Light) [0977]+[1367/1437]	1

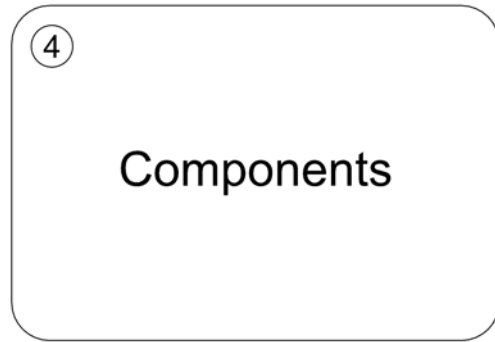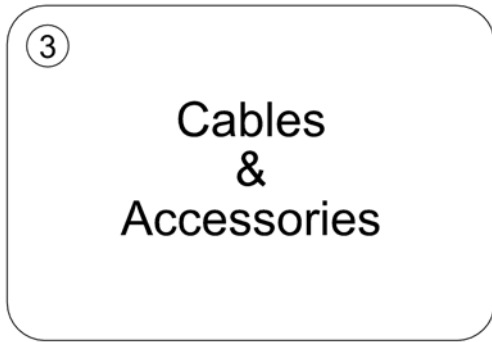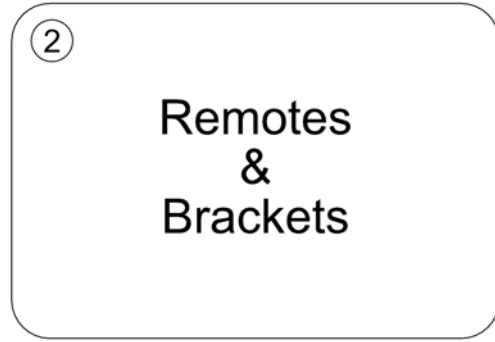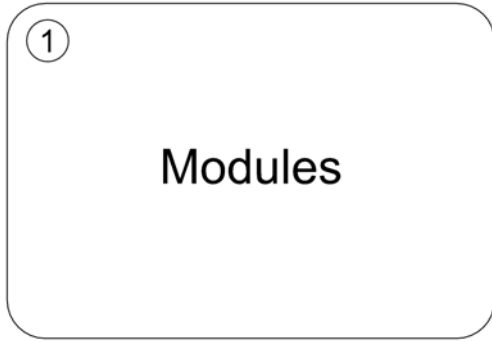


1	-	(Modules)	1
2	-	(Remotes & Brackets)	1
3	-	(Cables & Accessories)	1
4	-	(Components)	1

(Modules)

-	(Power Modules)	1
-	(Mounting Hardware, Power Module & Light Board)	1
-	(Actuator Modules (ACTs))	1
-	(Bracket ACT - Modulite Seat)	1
-	(Bracket ACT at Chassis - Narrow Base)	1
-	(G-Trac Sensor (Gyro))	1

(Mounting Hardware, Power Module & Light Board)

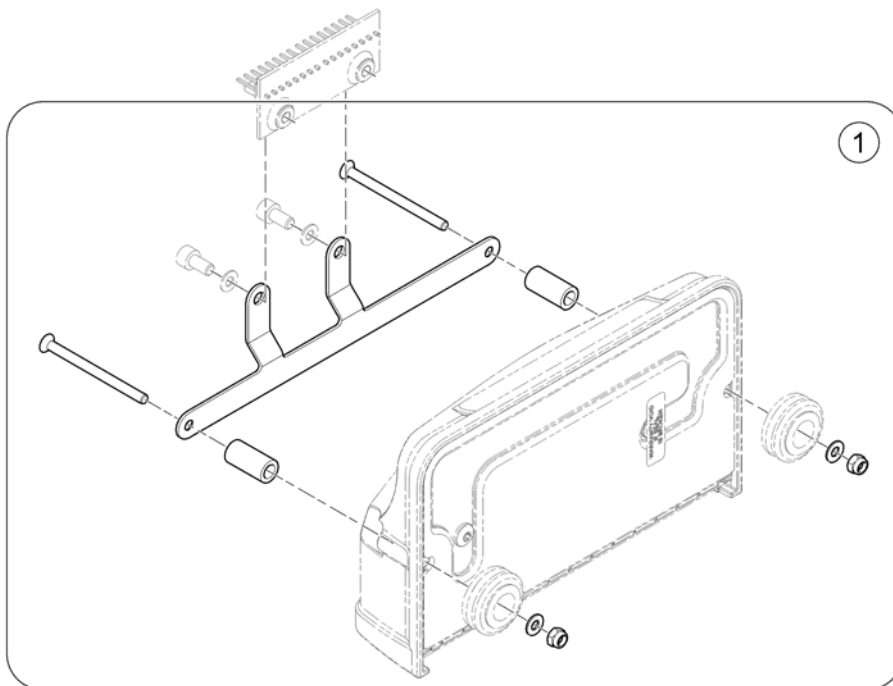160612.no



1 SP1606710

(Hardware, Power Module & Light Board)


1

(Cables & Accessories)

-	(Bus Cables)	1
-	(Light Cable)	1
-	(Light Board)	1
-	(Lifter/Tilt Module Cables)	1
-	(Actuator Cable)	1
-	(4-way Expansion Block)	1

(4-way Expansion Block)


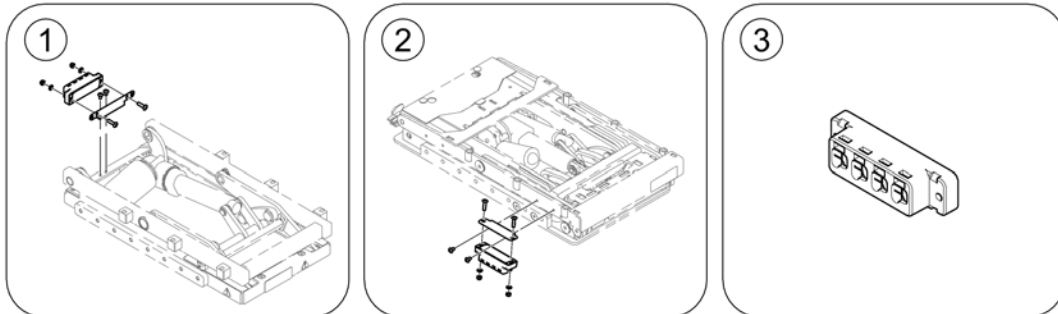
160054.lin



1	-	(4-way Expansion Block (with Tilt module))	1
2	-	(4-way Expansion Block (with Lifter/Tilt module))	1
3	-	(4-way Expansion Block (without any module))	1


Batterier, kabler og ladere

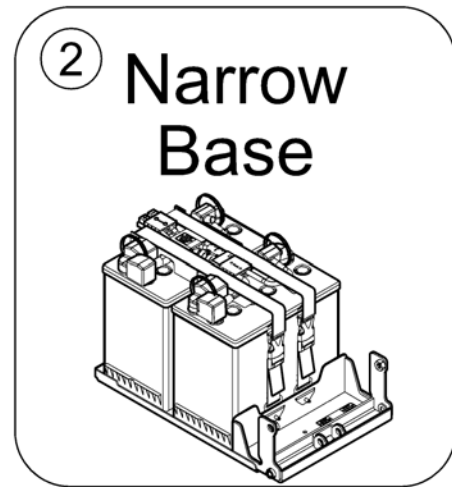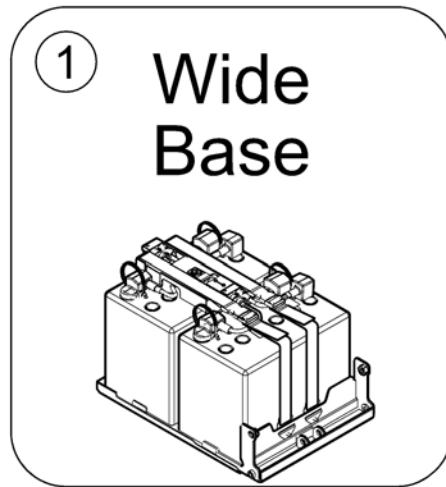SPCTDX_SP280.iso



1	-	(Battery, Tray & Wiring - Standard Base)	0
1	-	(Battery, Tray & Wiring - Standard Base)	1
2	-	(Battery, Tray & Wire Harness - Narrow Base)	0
3	-	(Chargers)	1

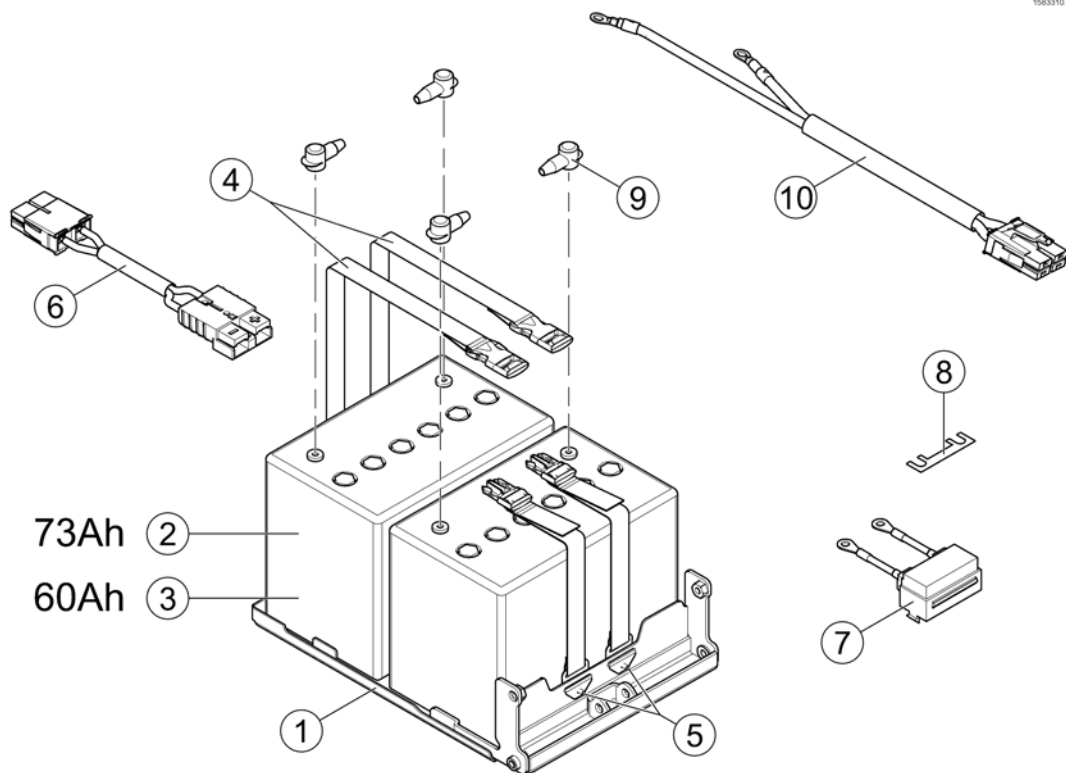
(Battery, Tray & Wiring - Standard Base)



-	(Important advice! Replace batteries in pairs. Never mix: manufactures, date codes or technologies (Gel with AGM batteries))	1
-	SP1523150 (Strip, Self Adhesive Wear (2-1/2'Long x1/2"Wide))	1
1	SP1179723 (Kit, Battery Tray Assy - Standard Base)	1
2	P50016 BATTERI 70AH MK12-73 TROJAN	2
3	1432944 BATTERI 60AH MK12-60 TROJAN <i>Lowrider Option only!</i>	2
4	SP1179296 (Battery Strap Assy EU, GP24)	2
5	SP1523266 (Rivet (3/16 x7/16"))	4
6	SP1527780 BATTERIKABEL TDX	1
7	1432561 SIKRINGSHOLDER STORM3 KOMP.M 80A SIKRING	1
8	1433734 SIKRING 80A TIL BATTERI	1
9	P52996 POLSKO DEKSEL	4
10	- (Not available so far.)	0

Invacare® TDX-SP2 Lowrider, TDX-SP2 NB | Batterier, kabler og ladere
(Battery, Tray & Wiring - Standard Base)

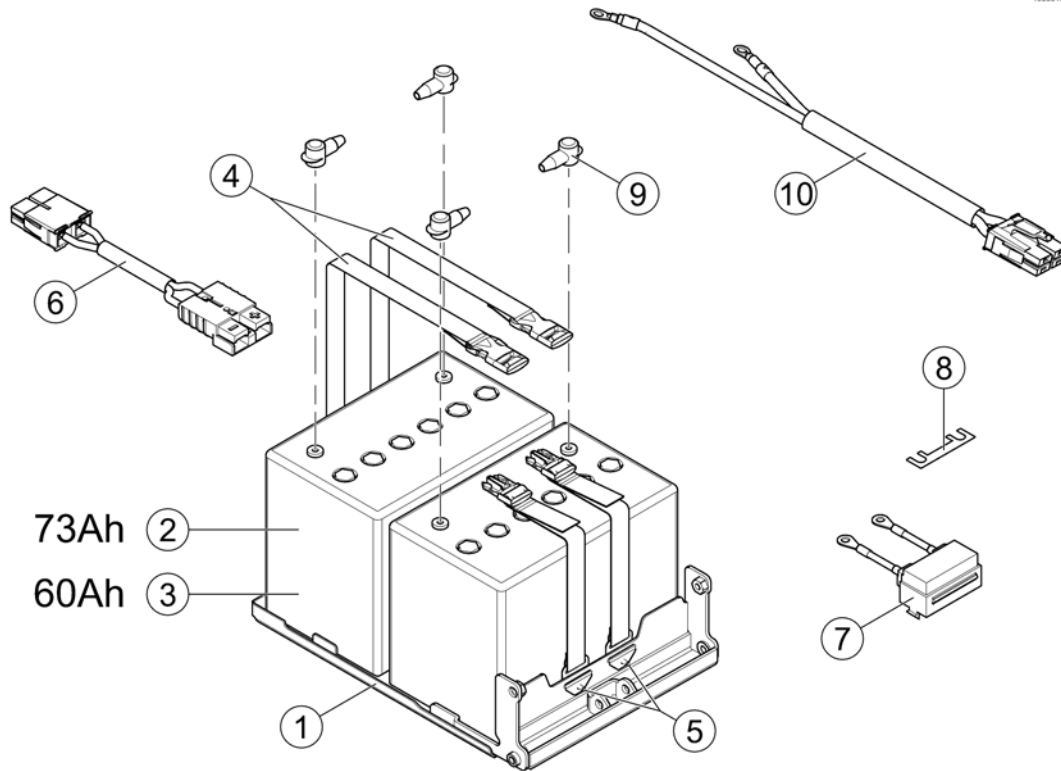
1563310.000



-	(Important advice! Replace batteries in pairs. Never mix: manufactures, date codes or technologies (Gel with AGM batteries))	1
-	SP1523150 (Strip, Self Adhesive Wear (2-1/2'Long x1/2"Wide))	1
1	SP1179723 (Kit, Battery Tray Assy - Standard Base)	1
2	P50016 BATTERI 70AH MK12-73 TROJAN	2
2	SP1587897 (Battery, 12V/73,6 Ah C20, Gel - preconditioned (Sweden only))	2
3	1432944 BATTERI 60AH MK12-60 TROJAN <i>Lowrider Option only!</i>	2
4	SP1179296 (Battery Strap Assy EU, GP24)	2
5	SP1523266 (Rivet (3/16 x7/16"))	4
6	SP1527780 BATTERIKABEL TDX	1
7	1432561 SIKRINGSHOLDER STORM3 KOMP.M 80A SIKRING	1
8	1433734 SIKRING 80A TIL BATTERI	1
9	P52996 POLSKO DEKSEL	4
10	- (Not available so far.)	0

(Battery, Tray & Wire Harness - Narrow Base)

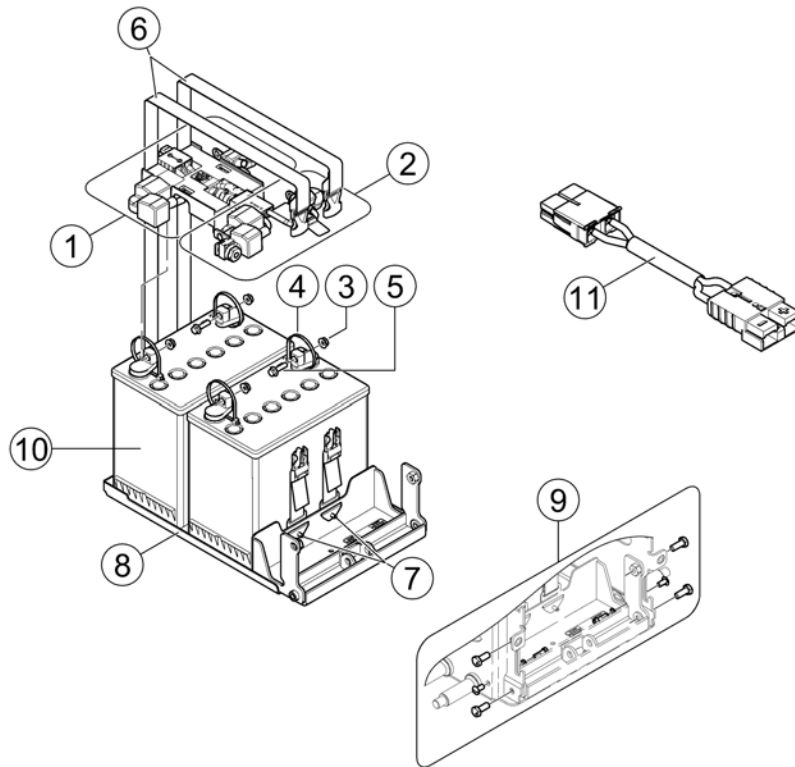
1583311.00



-	(Important advice! Replace batteries in pairs. Never mix: manufactures, date codes or technologies (Gel with AGM batteries))	1
-	SP1523150 (Strip, Self Adhesive Wear (2-1/2"Long x1/2"Wide))	1
1	SP1120363 (Kit, Battery Wiring w/Terminal Hardware, rear)	1
2	SP1120362 (Kit, Battery Wiring w/Terminal Hardware, front)	1
4	SP1522970 (Package, Cable Ties (11-1/2"Long))	1
5	SP1522972 (Screw, Hex Flange Whiz Lock (1/4-20x7/8"))	4
6	SP1144682 (Battery strap)	1
7	SP1523266 (Rivet (3/16 x7/16"))	4
8	SP1179722 (Kit, Battery Tray Assy - Narrow Base)	1
9	SP1179896 (Kit, Battery Tray Stop Hardware)	1
10	PG50989 BATTERI 50AH MK12-50 TROJAN	2
11	SP1527780 BATTERIKABEL TDX	1

(Chargers)


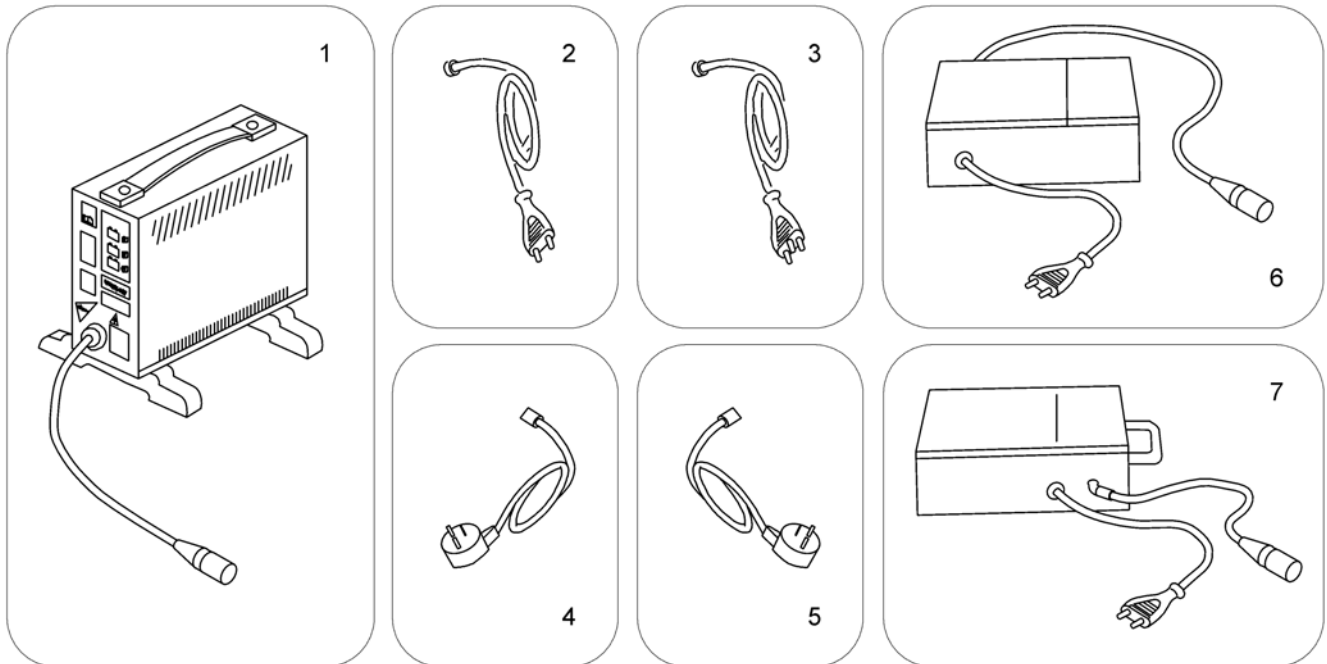
1589799 Iso




1	SPP51981	(Charger 8A without Power Cord) <i>[1343]</i>	1
2	P21079	(Power Cord Continental (Schuko GCB-0401))	1
3	P21080	(Charger Cable (UK))	1
4	1443841	(Power Cord Swiss (black 6A/250V))	1
5	SP1553849	(Power Cord Israel (IEC 320 C13 to SI-32))	1
6	SP1545511	LADER LADAC/DEFA 701239 MED PLUGG <i>[1389]</i>	1
7	5813987	(Charger 8A (Protectionclass IP54) incl. Power Cord Continental (Sweden only)) <i>[1388]</i>	1

(Lightsystem)

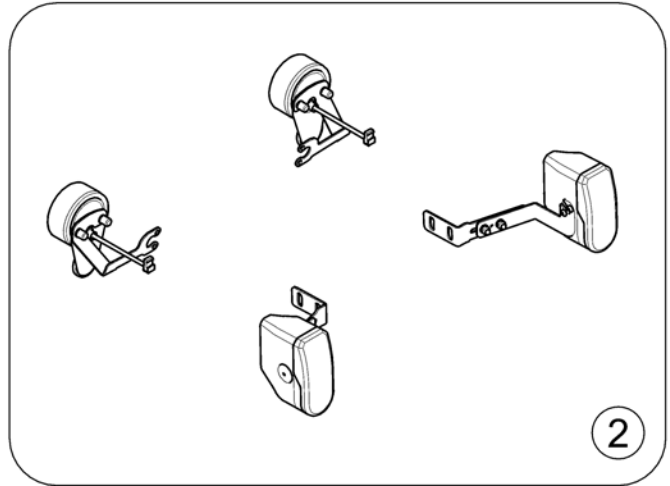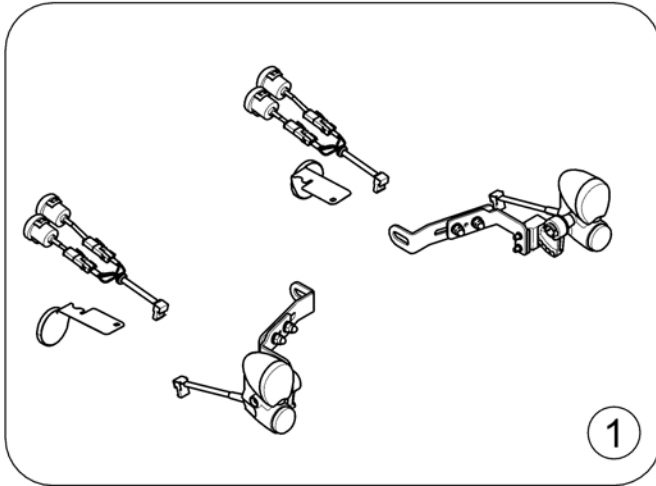1582964.iso




1	-	(Lightsystem LED) [1437]	1
2	-	(Lightsystem german Homologation) [1367]	1

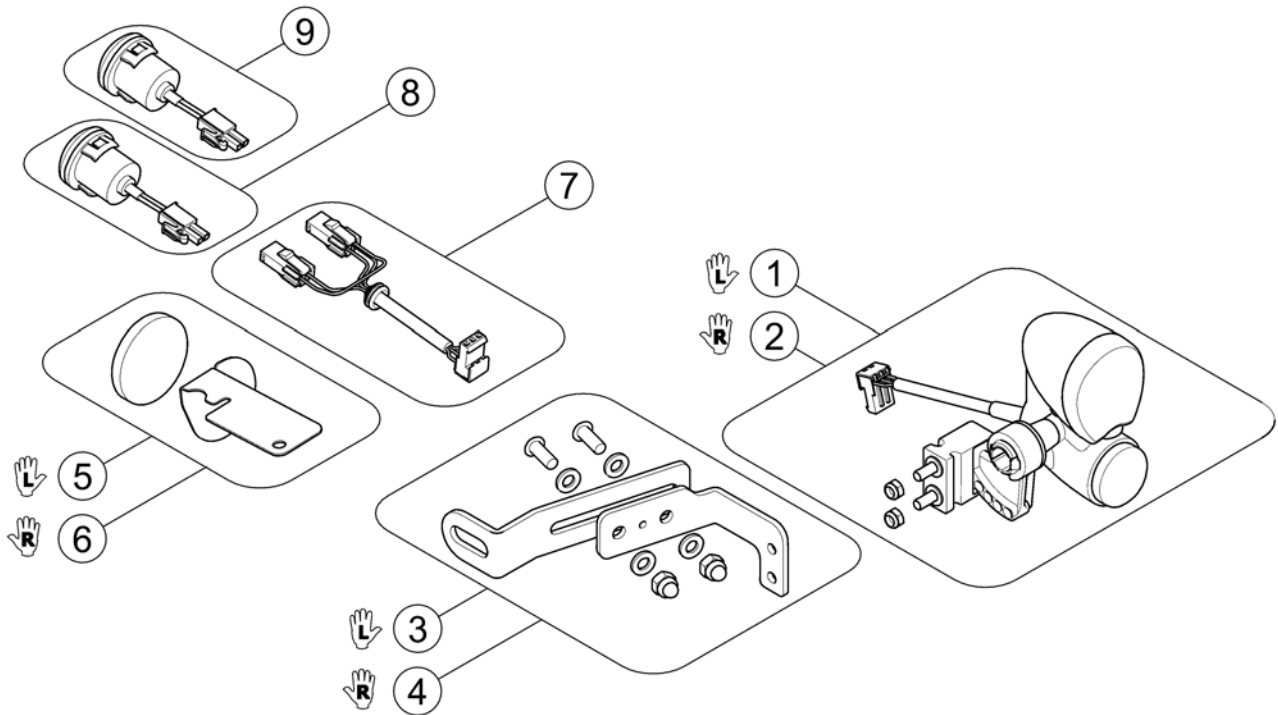


1	SP1536030	LYKT LED VE KPL U/FESTE TDX-SP2 LOW-RIDER	1
2	SP1576193	LED FRONTLYS TDX SP2 HØYRE SIDE	1
3	SP1583093	HOLDER LYKT LED VE TDX-SP2 LOWRIDER	1
4	SP1583092	HOLDER LYKT LED HØ TDX-SP2 LOWRIDER	1
5	SP1583095	(Reflector rear LH)	1
6	SP1583094	(Reflector rear RH)	1
7	SP1529712	BAKLYS CABLE STORM4	1
8	SP1530288	BAKLYS STORM4	1
9	SP1530287	BLINKLYS BAK STORM4	1

(Lightsystem german Homologation)

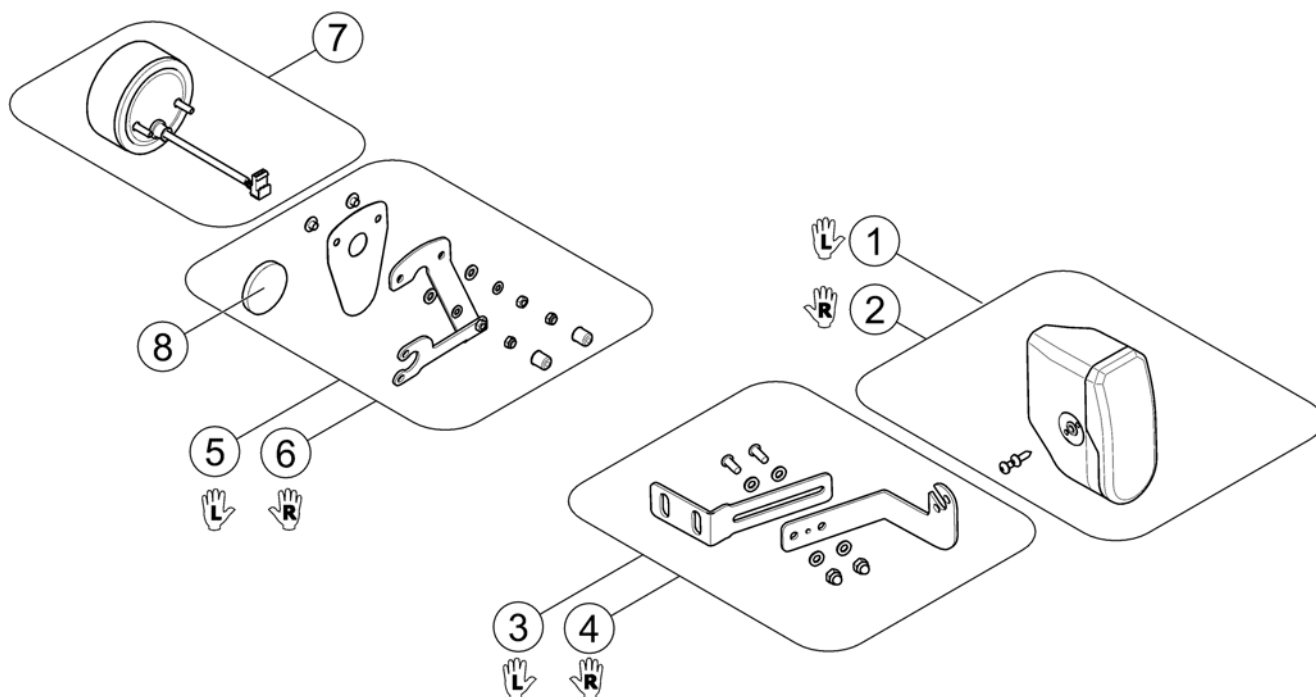1583077.00



-	P31108	LYSPÆRE 24V/10W	1
-	P31191	PÆRE 24V/5W BAKLYS STORM 3	1
-	SP1419271	(Light Bulb 24V/21W, yellow (Turn Signal front/rear))	1
!!!	-	(Rear Light Protector Kit)	0
1	1440426	FRONTLYS VE KPL THYPOON	1
2	1420412	FRONTLYS HØ KPL STORM3 TORNADO G50	1
3	SP1583088	LYSHOLDER FORAN VE TDX-SP2	1
4	SP1583087	LYSHOLDER FORAN HØ TDX-SP2	1
5	SP1583089	(Light Bracket rear LH)	1
6	SP1583090	(Light Bracket rear RH)	1
7	1442736	BAKLYS KPL STORM3	1
8	SP1541591	(Reflector, red)	1

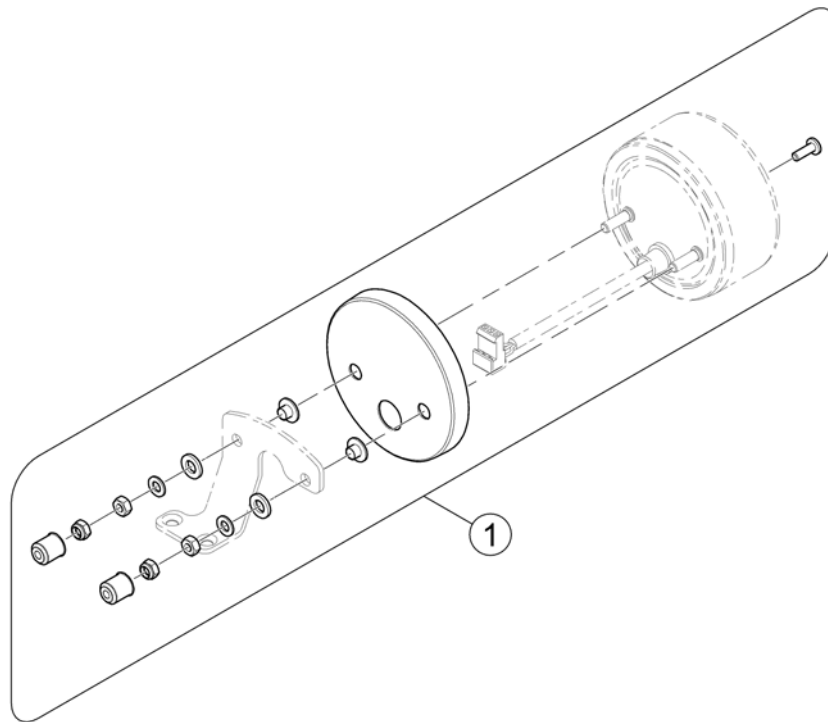


1 SP1581827

(Rear Light Protector Exchange Kit)

1

(Reflective Sign)


-	(Reflective Sign) <i>[1378]</i>	1
-	(Reflective Sign (Ultra Low Maxx Seat)) <i>[1378]</i>	1

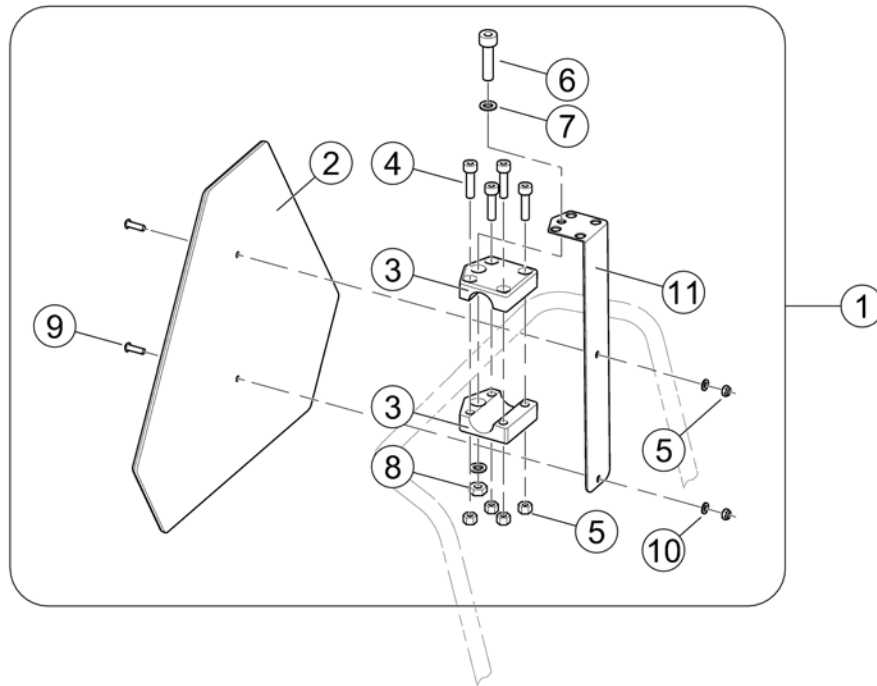

1	SP1430410	(Reflective Sign ECE-R 69, complete)	1
2	1430834	(Rear Marking Sign ECE-R 69)	1
3	Z12008A.P40	(Gripping clip with hexagon hole, RAL 9005 black)	2
4	D50004	(Cylinder Bolt M6x30-8.8 bk - DIN 912)	4
5	F20109	MUTTER 6 KT M6	6
6	D00102	SEKSKANT SKRUE M8X45	1
7	G00107	STOPPSKIVE 8,4X17X1,6	2
8	F00002	SEKSKANTMUTTER M8	1
9	D90187	(Fillister Head Screw M6x20-10.9 bk - DIN 7380)	2
10	G00115	SKIVE 6,4X12,5X1,6 MISTRAL	2
11	1430418-9005	(Mounting Bracket for Reflective Sign - RAL 9005 black)	1
12	SP1580671	(Bracket, Reflective Sign, cpl.)	1

(Reflective Sign (Ultra Low Maxx Seat))

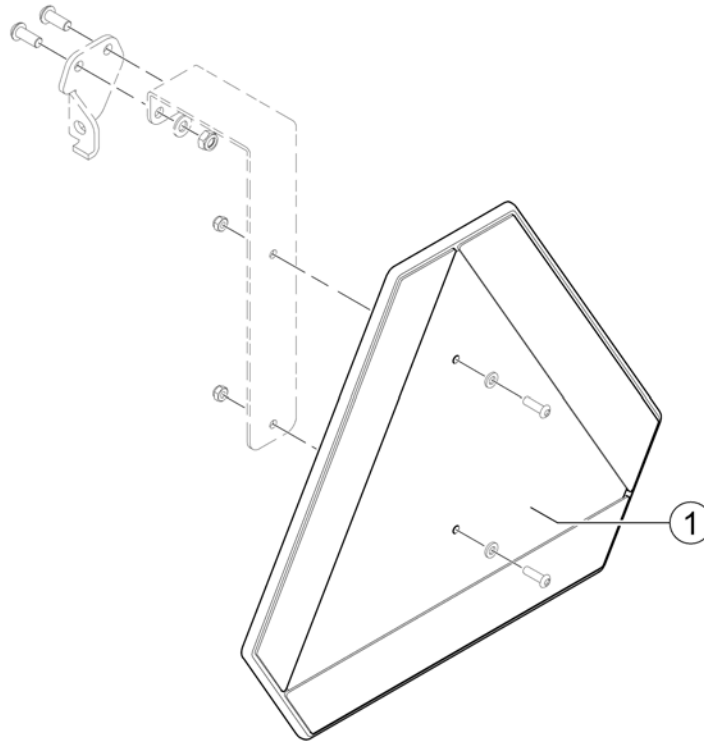
1508516.no



2	1430834	(Rear Marking Sign ECE-R 69)	1
1	1430834	(Rear Marking Sign ECE-R 69)	1

Invacare® TDX-SP2 Lowrider, TDX-SP2 NB | Sikkerhetsutstyr og tilbehør
(Adapter-Kit DAHL Docking Station - Standard Base)

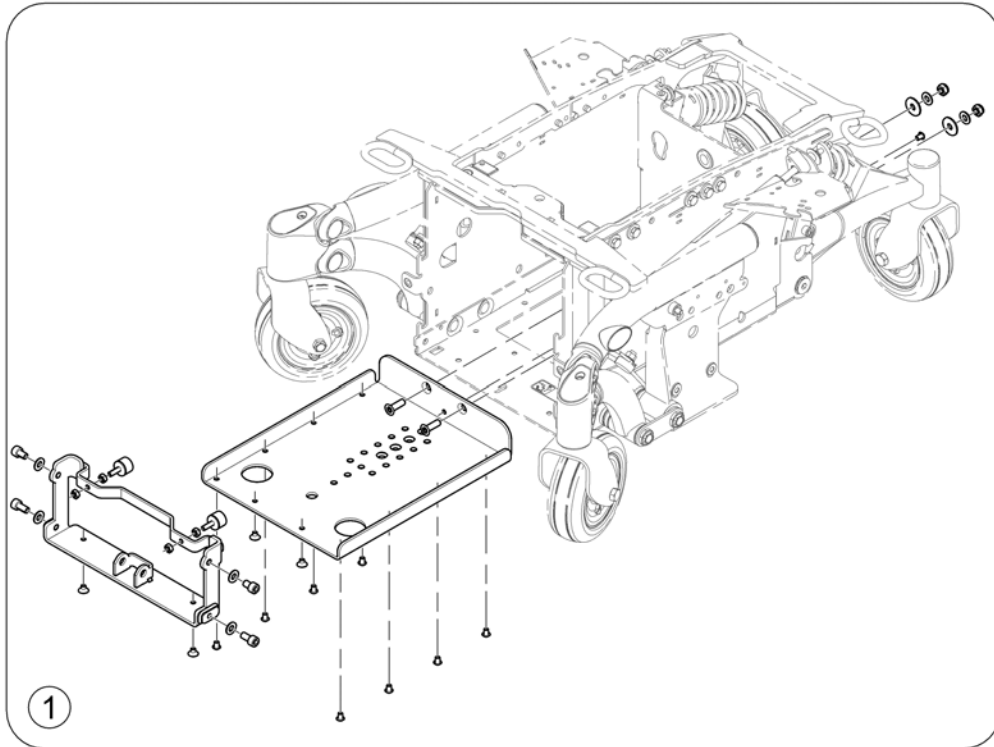1583706 iso




1 SP1564843

(Retrofit-Kit for DAHL Docking Station - Standard Base)

1

(Adapter-Kit DAHL Docking Station - Narrow Base)

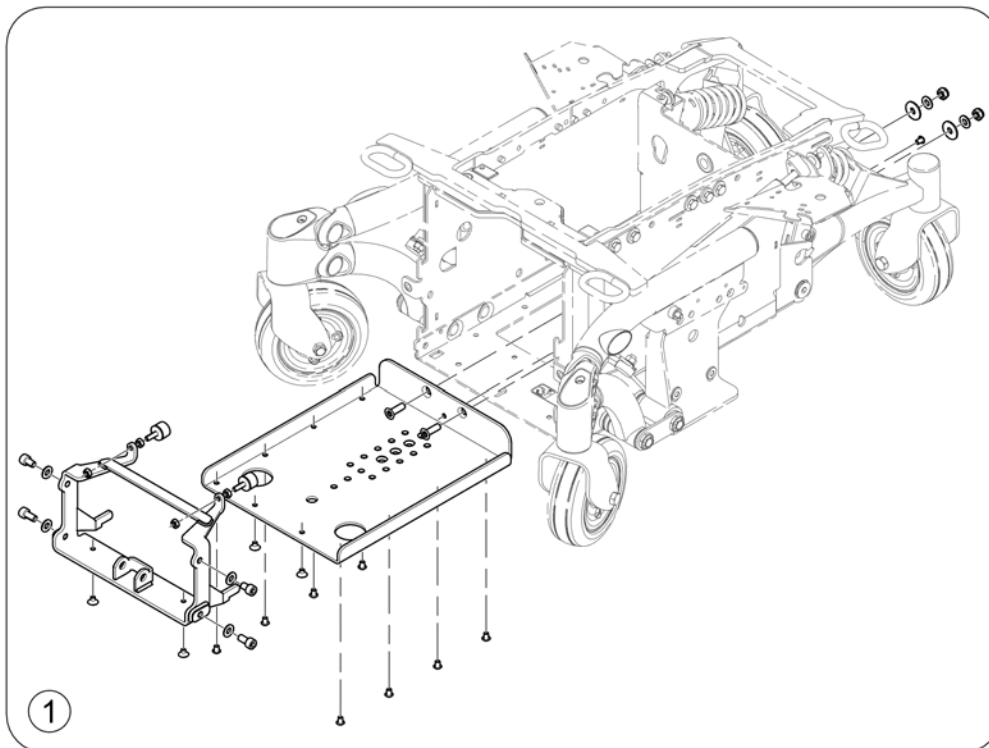1583707 no



1 SP1564807


(Retrofit-Kit for DAHL Docking Station - Narrow Base)

1


-	(Lap Belts (Modulite Seat)) <i>[1329/2113 for 2100]</i>	1
-	(Lap Belt (Ultra Low Maxx)) <i>[2113 for 2020]</i>	1
-	(Lap Belt (Max Seat)) <i>[1329 for 1145]</i>	1
-	(H-Belts (Modulite Seat)) <i>[1332 for 2100]</i>	1

(Lap Belts (Modulite Seat))

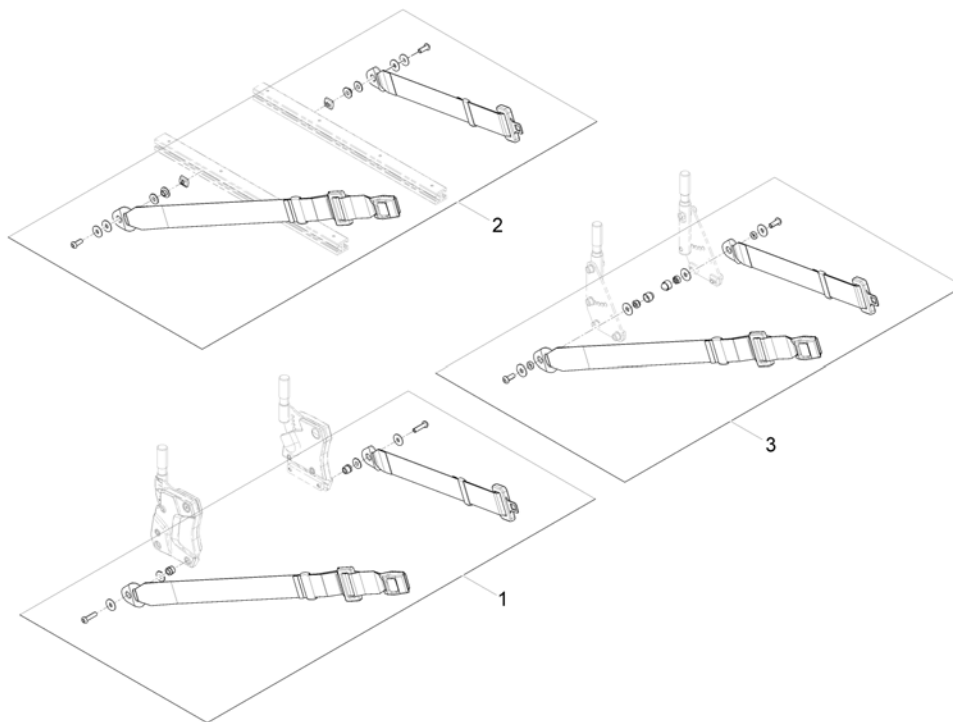1571836 iso



1	SP1570381	HOFTEBELTE MODULITE FESTET PÅ RYGG	1
2	SP1571183	HOFTEBELTE C TIL SKINNE MODULITE	1
3	1422845	HOFTEBELTE MED SNAP LÅS STORM 3/ DRAGON	1

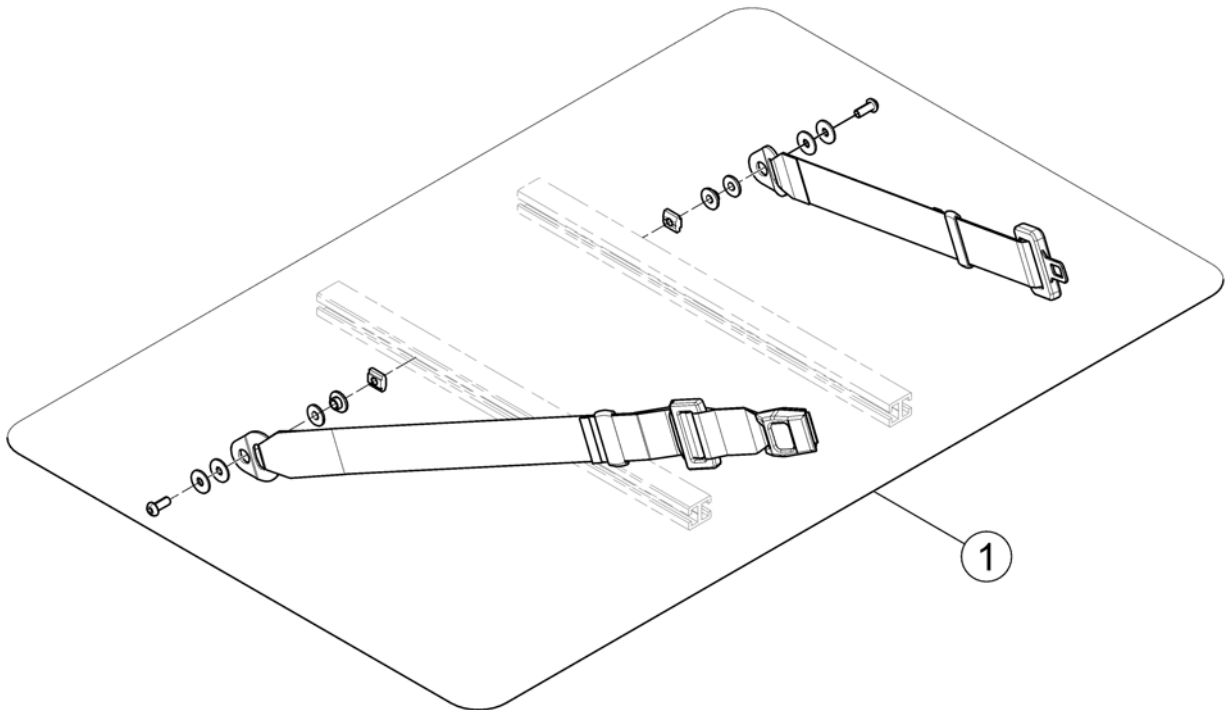


1 SP1571183

HOFTEBELTE C TIL SKINNE MODULITE

1

(Lap Belt (Max Seat))

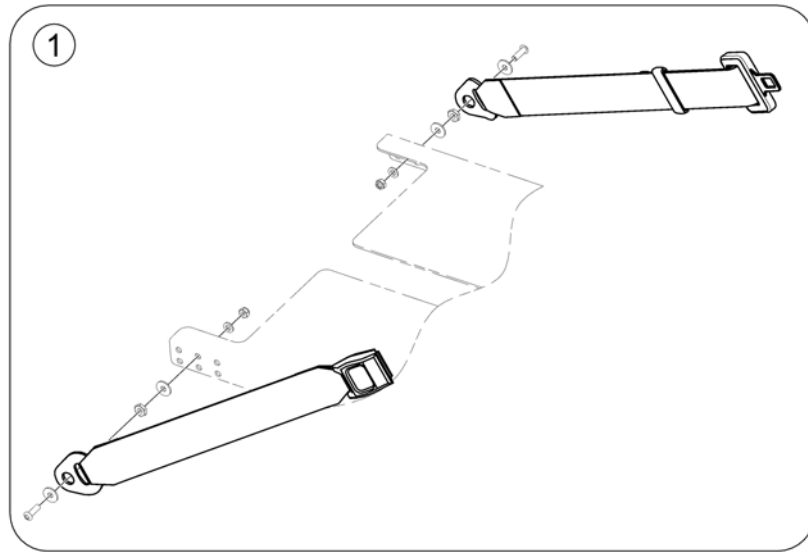
1582371.00



1 1551067


HOFTEBELTE, LANG MAX

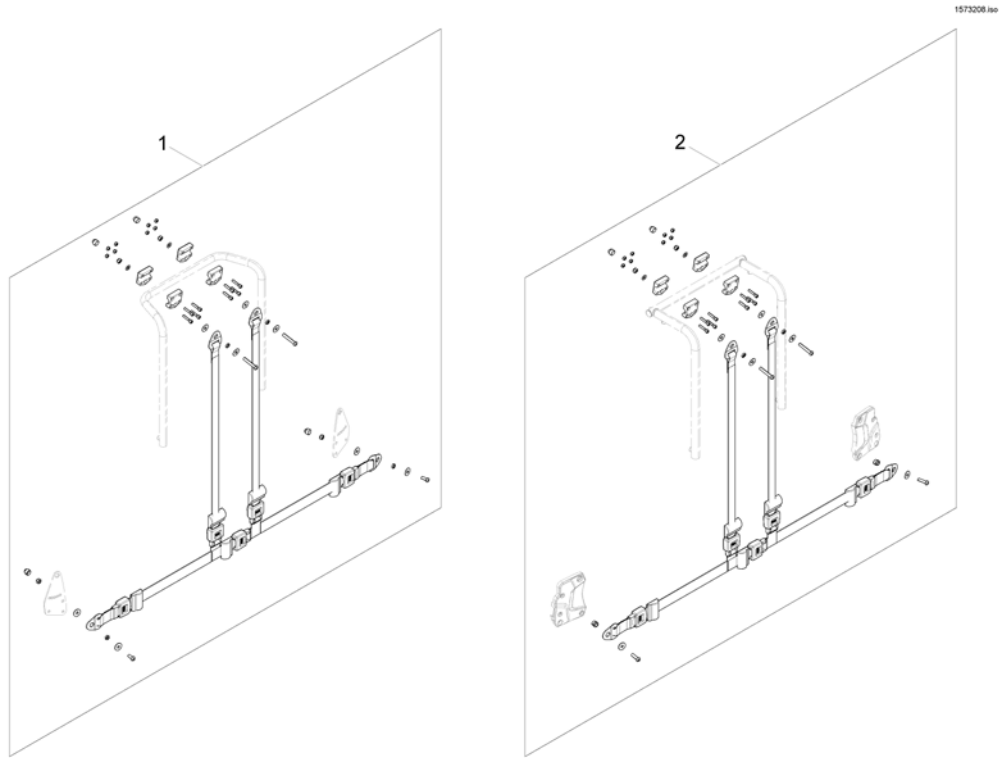1



1	SP1570385	(H-Belt for Backrest Frame)	1
2	SP1570380	(H-Belt for Back Posts)	1

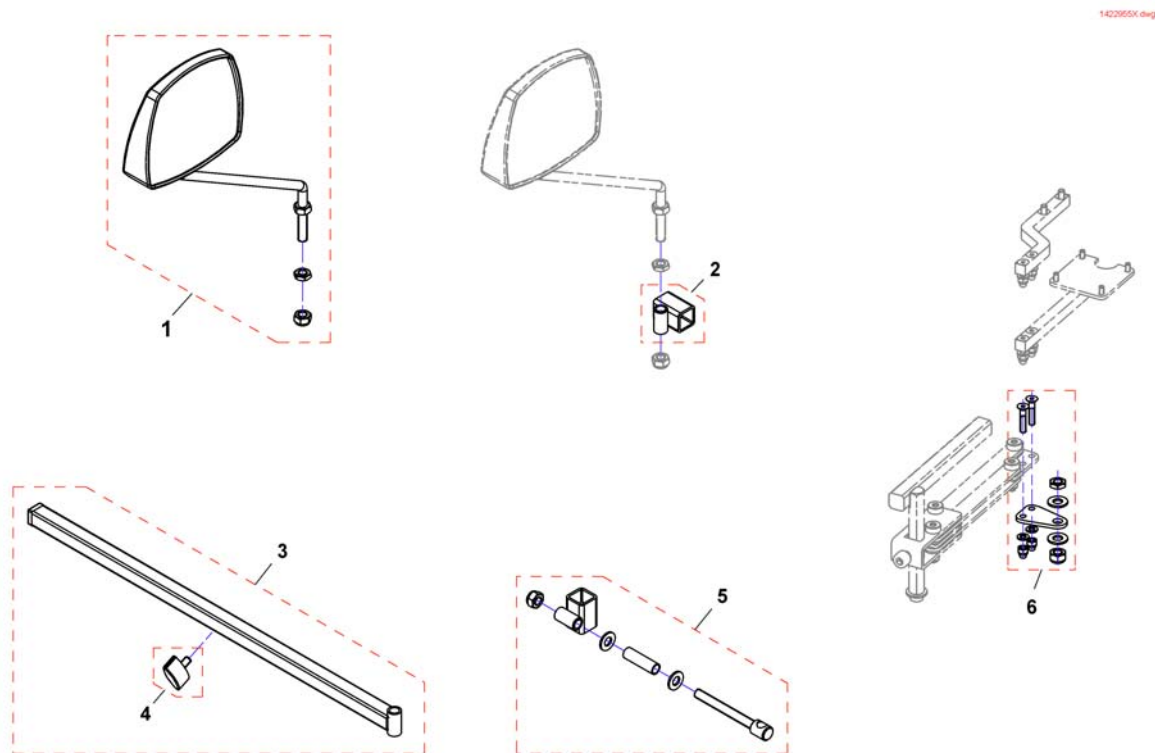(Rearview Mirror)




1	SP1422955	SPEIL KPL KUN SPEIL STORM 3	1
2	F276005.00	HOLDER SPEIL DRAGON STORM4 HUSK SP1527151	1
3	SP1423435	HOLDER SPEIL STORM3	1
4	L31001	VINGSKRUV M 6X10 GARANT G24 S	1
5	SP1422957	SPEIL FESTE BRAKETT STORM4 PÅ SWINGAWAY	1
6	SP1527151	FESTE FOR SPEILFESTE STORM4 TDX BORA	1


(Headrests)

-	(Head- / Neckrest)	1
-	(Headrest sidewise adjustable)	1
-	(Headrest with cheek support)	1
-	(Headrests - Ultra Low Maxx Seat)	1

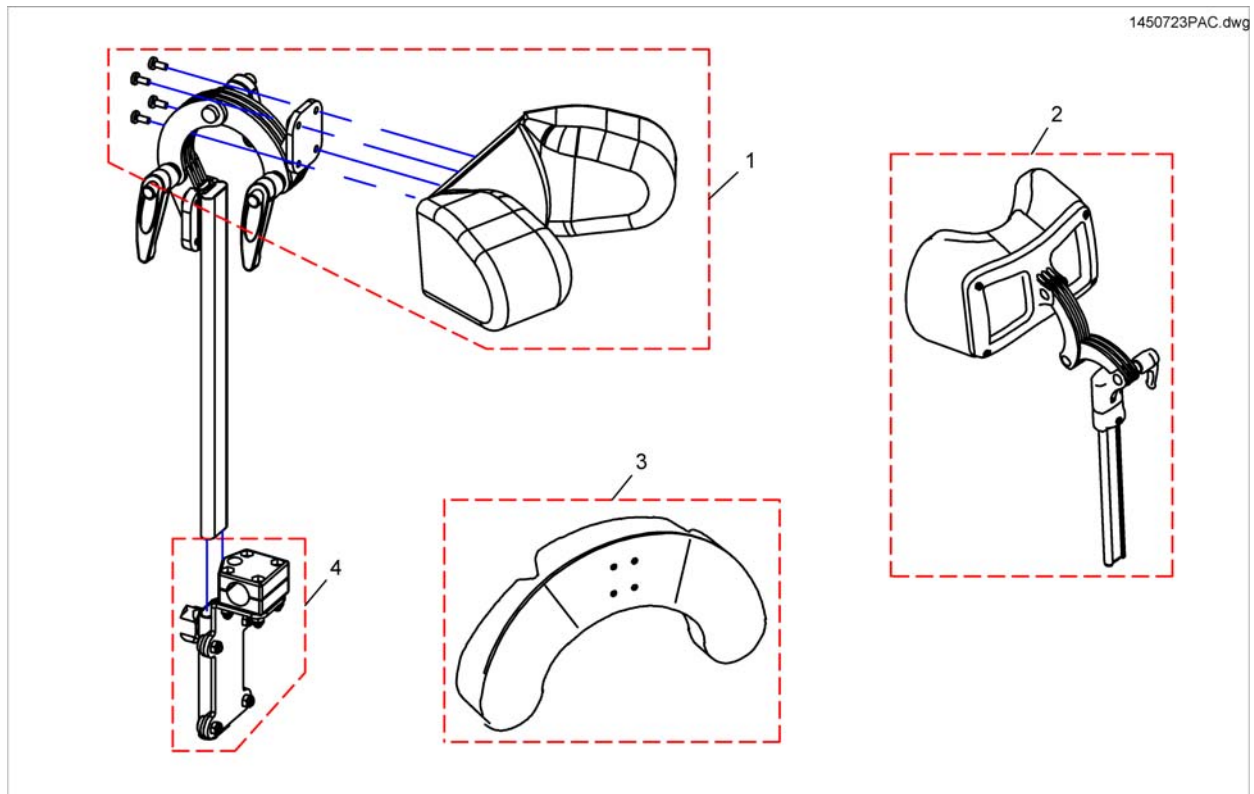(Head- / Neckrest)




1	1487404	HODESTØTTEPUTE KPL BORA	1
2	SP1527050	NAKKESTØTTE DYBDEJUSTERB STORM4 DARTEX	1
3	1444304	PUTE HAKESTØTTE BELLIS/A ZALEA	1
4	SP1422914	HOLDER HODESTØTTE RYGG STORM4 TDX	1

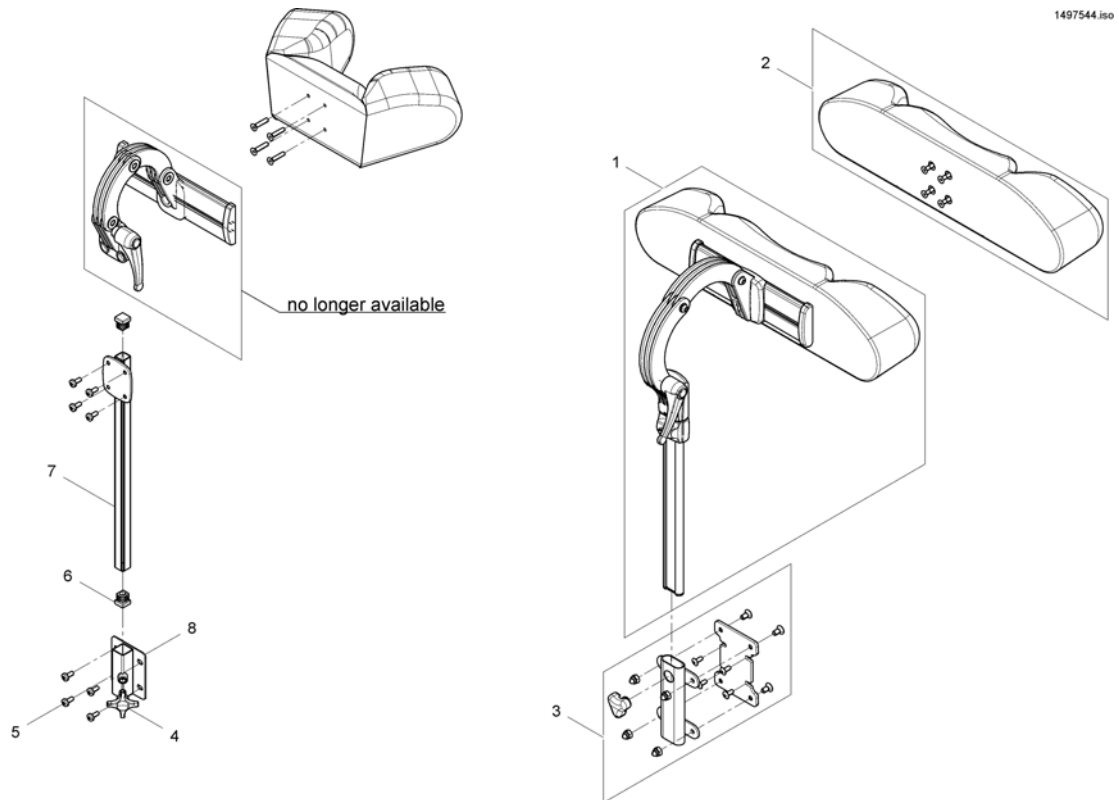


1	SP1530119	NAKKESTØTTE SIDEV JUST STORM4 EASYADAPT	1
2	SP1422865	HODESTØTTEPUTE U FESTE STORM4	1
2	1422866	(HEADREST PAD STD./CONTOUR GREY)	1
2	SP1522214	(HEADREST PAD STD./CONTOUR BLACK PLUSH)	1
3	SP1544680	MONTERINSKIT HODESTØTTE REA	1
4	256099	(Wing Screw M8 X 20 639)	1
5	D30197	(Lens head screw DIN7985-M05x12-4.8 zc)	1
6	L01055	ENDESTYKKE FESTEBR.	1
7	1451639-0205	(Headrest REA Angle Adjustment)	1
8	26002351	(Tube Headrest)	1

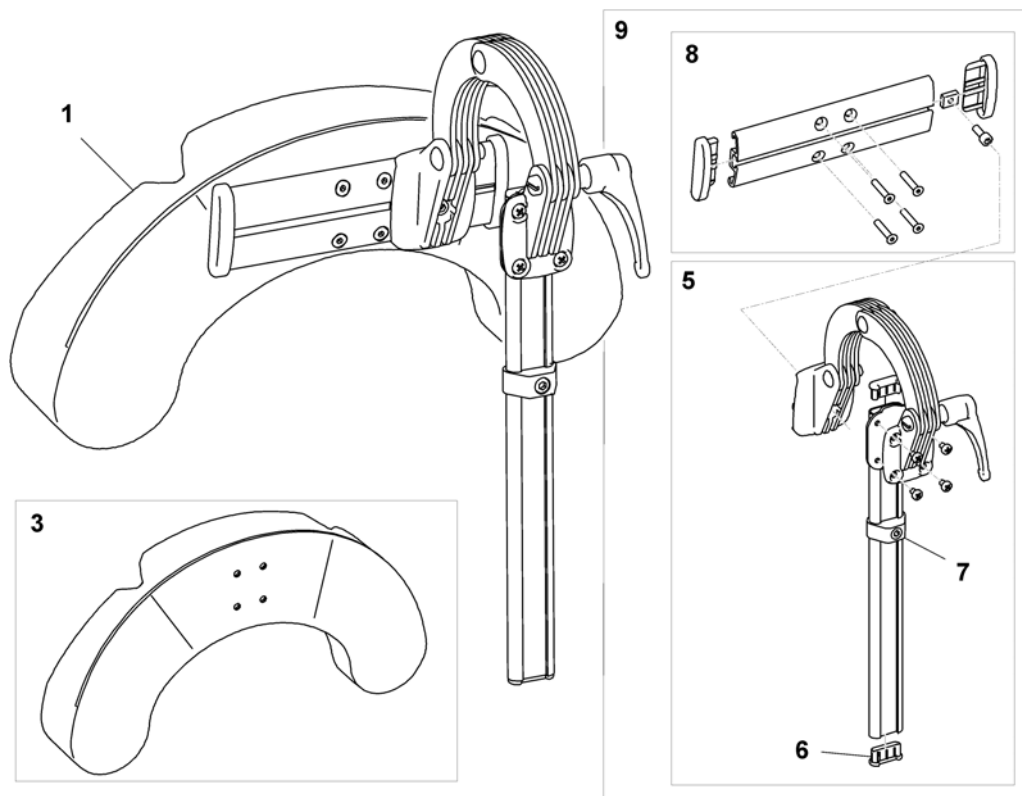(Headrest with cheek support)



1	1446560	HODE/SKINN STØTTE KPL TR26 MED SIDEFORFLYTNING	1
1	SP1541658	(Head/cheek support TR26 cpl. sidewise adjustable EasyAdapt)	1
3	1444304	PUTE HAKESTØTTE BELLIS/A ZALEA	1
5	SP1529572	HOLDER HODESTØTTE EASYADAPT	1
6	5716497	HYLSE MESSINGSKRUE FOR LOMME REA	1
7	SP1448661	(stop neck and headrest)	1
8	1445838	SKINNE HODESTØTTE KPL BE LLIS AZAL FOC SPI3	1
9	1485538	FESTE HODESTØTTE STORM4	1

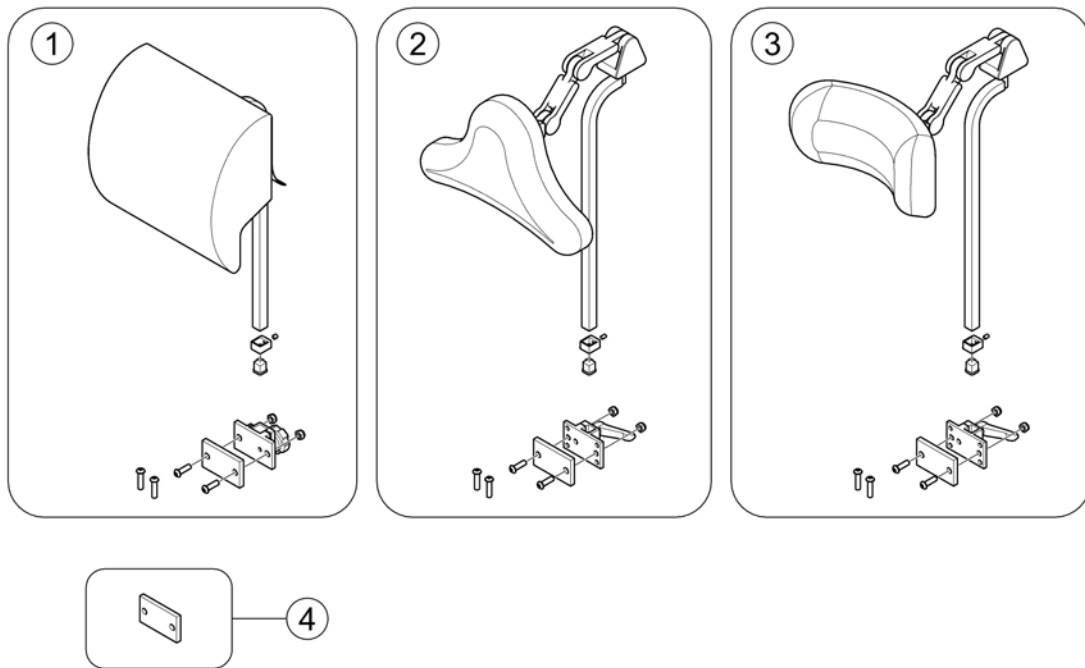
1588517.no



1	SPTRC1929	(Onyx Auto Style Headrest)	1
2	SPTRC1930	(All-Position Headrest)	1
3	SPTRC1931	(Comfort Plus Headrest)	1
4	SPTRX2685	(Headrest Shim Kit Posture Back)	1
4	SP1544188	(Headrest Shim Kit Deep Back SW 410-430mm)	1
4	SP1544187	(Headrest Shim Kit Deep Back SW 460-610mm)	1

(Lateral Supports)


-	(Lateral Supports - Modulite Seat)	1
-	(Lateral Support - Ultra Low Maxx Seat)	1

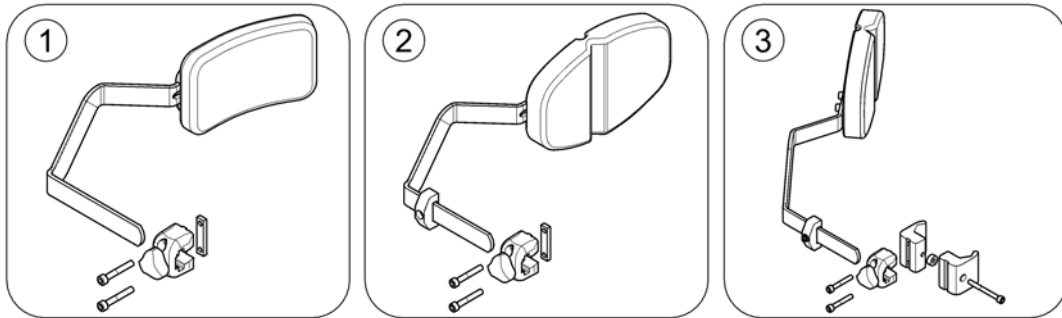

1	-	(Lateral supports, Flex3 backrest)	1
2	-	(Lateral supports, Flex3 backrest - multi adjustable)	1
3	-	(Lateral supports, Flex3 backrest (tension adjustable) - multi adjustable)	1

(Lateral supports, Flex3 backrest)

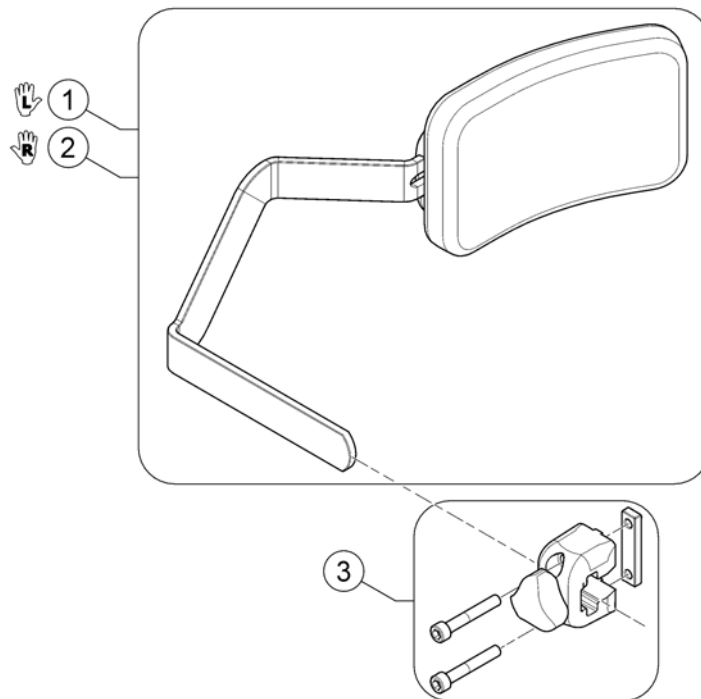
1007155.00



-	-	(Design change on 01. July 2016 - we suggest to change older supports on both sides.)	1
1	SP1552678	Sidestøtte forlenget ve	1
2	SP1552679	Sidestøtte forlenget hø	1
3	1543574	SIDESTØTTEFESTE KPL LAGUNA 2, CLEMATIS	1

(Lateral supports, Flex3 backrest - multi adjustable)

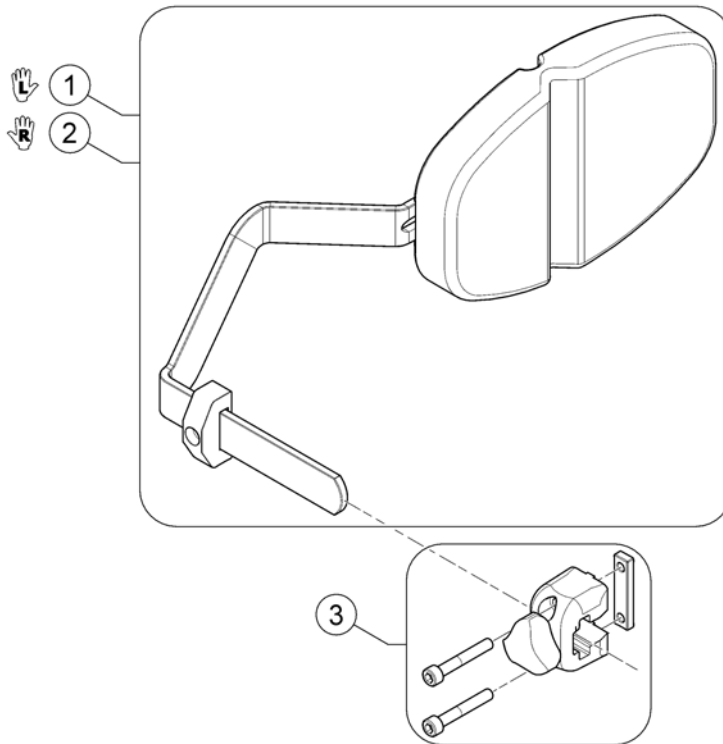
160208_00



-	-	(Design change on 01. July 2016 - we suggest to change older supports on both sides.)	1
1	SP1555692	Sidestøtte multi multi r eg forlenget ve	1
2	SP1555691	Sidestøtte multi multi r eg forlenget hø	1
3	1543574	SIDESTØTTEFESTE KPL LAGUNA 2, CLEMATIS	1

(Lateral supports, Flex3 backrest (tension adjustable) - multi adjustable)

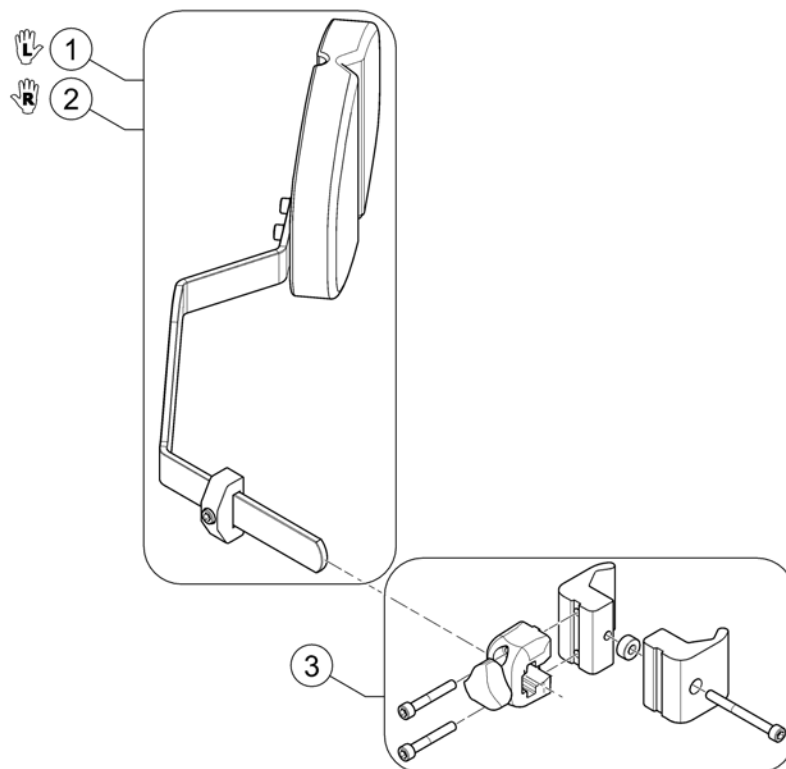
1607219.no



1	SP1562920	(Lateral Support, Flex3 back, tension adj., extended, multi-adj. - LH)	1
2	SP1562919	(Lateral Support, Flex3 back, tension adj., extended, multi-adj. - RH)	1
3	SP1554143	SIDESTØTTE FESTE FLEX3 TIL VELCRO RYGG	1

Invacare® TDX-SP2 Lowrider, TDX-SP2 NB | Setetilbehør
(Lateral Support - Ultra Low Maxx Seat)

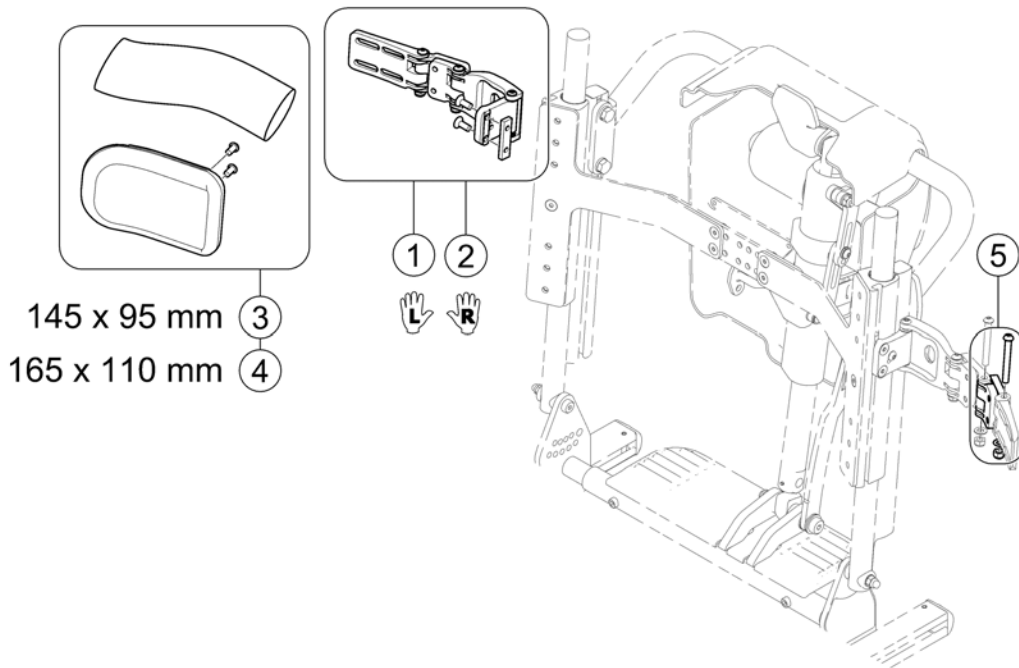1588596.iss




1	SPTRC1806	(Lateral Hardware Assy LH (w/o sleeve))	1
2	SPTRC1805	(Lateral Hardware Assy RH (w/o sleeve))	1
3	SPTRC2123-C03	(Lateral Pad Kit 145x95 mm & sleeve for hardware)	1
4	SPTRC2123-C04	(Lateral Pad Kit 165x110 mm & sleeve for hardware)	1
5	SPTRC2091	(Extension Kit, Lateral Support (Set of 2 pcs) incl. longer sleeves)	1

(Cane holder)

-	(Cane holder (Modulite: backrest frame))	1
-	(Cane holder (Modulite: back posts man recliner))	1
-	(Cane holder (Modulite: back posts pow recliner))	1
-	(Cane holder (MAX Seat))	1

(Cane holder (Modulite: backrest frame))

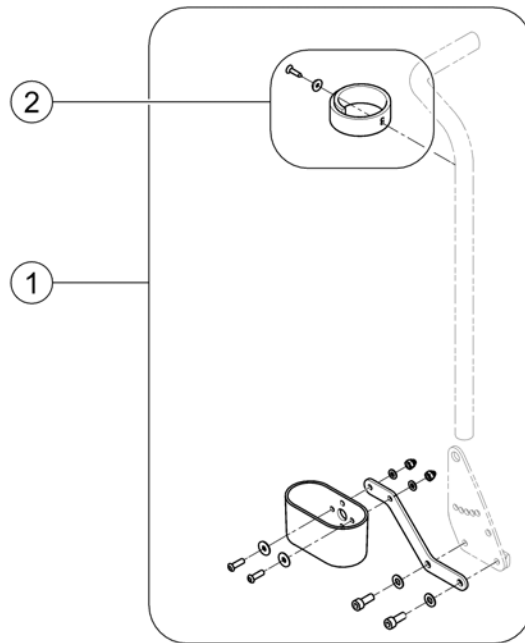
1571962.jpg



1	SP1570396	(Cane holder, standard backrest)	1
2	SP1420493	BÅND TIL KRYKKEHOLDER STORM4	1

(Cane holder (Modulite: back posts man recliner))

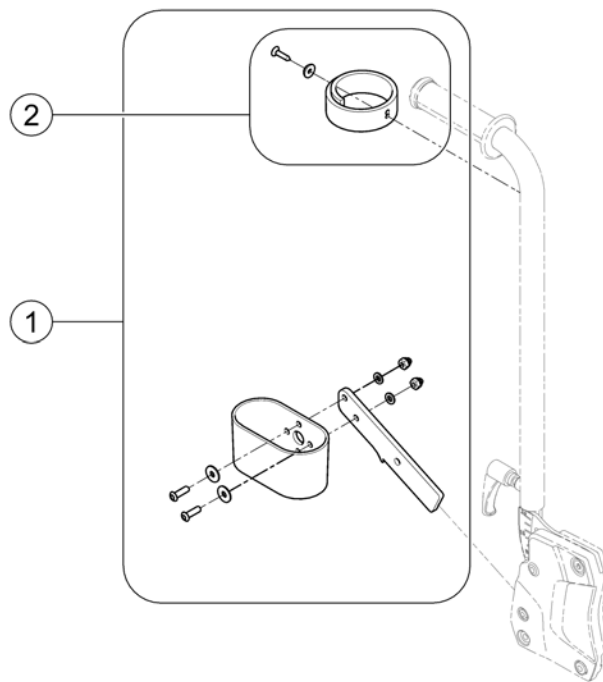
1571903.no



1	SP1570443	(Cane holder, manual backrest)	1
2	SP1420493	BÅND TIL KRYKKEHOLDER STORM4	1

(Cane holder (Modulite: back posts pow recliner))

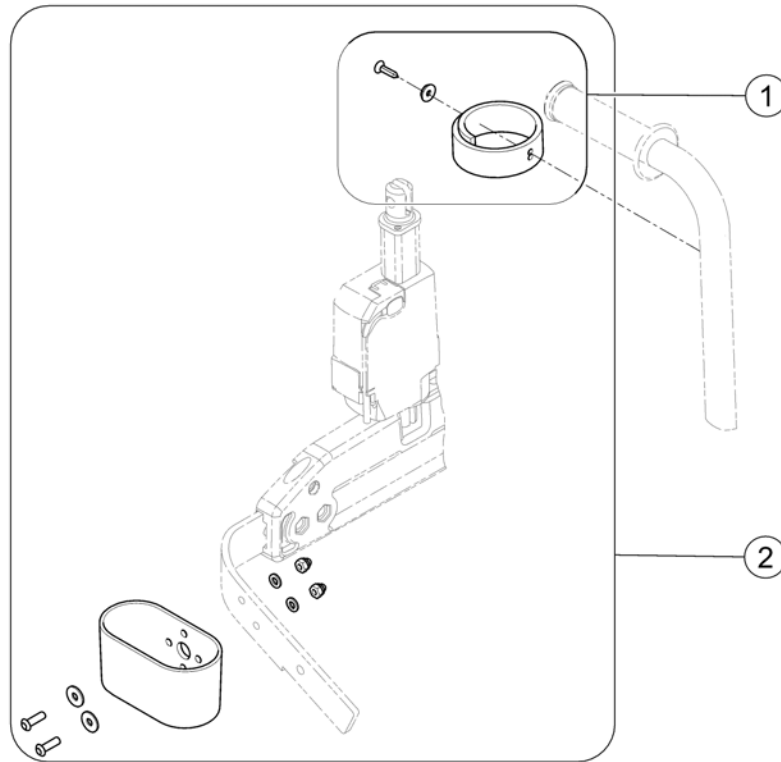
1571837.no





1	SP1420493	BÅND TIL KRYKKEHOLDER STORM4	1
2	SP1570382	KRYKKEHOLDER STORM4 STORM4 MODULITE	1

(Cane holder (MAX Seat))

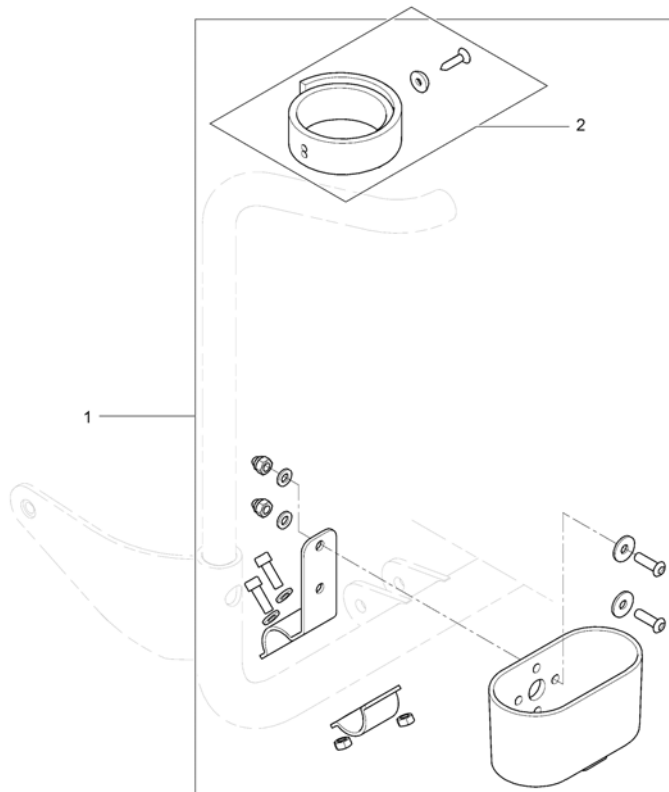1567255.iso




1	SP1543078	KRYKKEHOLDER STORM4 MAX	1
2	SP1420493	BÅND TIL KRYKKEHOLDER STORM4	1

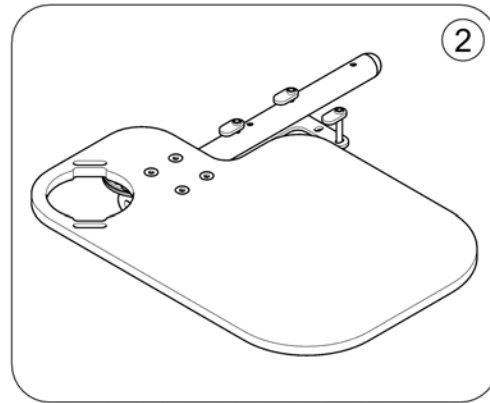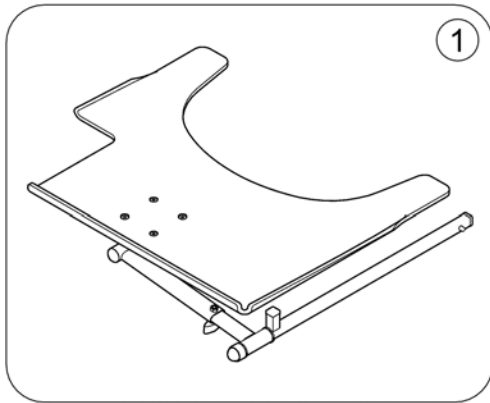


1	-	(Trays)	1
2	-	(Trays (Ultra Low Maxx Seat))	1

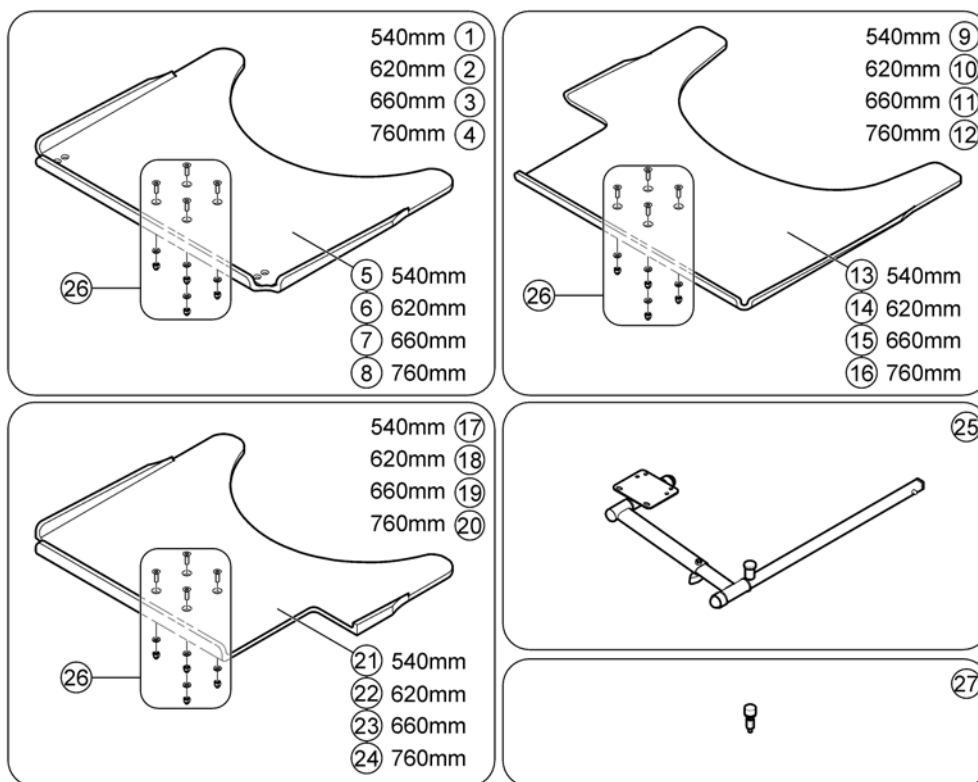(Trays)



1586046.iso


1	SP1548198	BORD UTEN UTSPARR SB38 STORM4	1
2	SP1548201	BORD UTEN UTSPARRING STORM4	1
3	SP1422954	BORD INKL FESTE KPL STORM UTEN UTSPARRING	1
4	SP1584644	(Swing-Out Tray, 760mm, Seatwidth 480mm to 580mm, no Cutout, cpl. with holder)	1
5	SP1535050	BORD 540MM SB33-43 UTEN UTSKJÆRING STORM4	1
6	SP1535010	BORD 620MM SB38-48 UTEN UTSKJÆRING STORM4	1
7	SP1535011	BORD 660MM SB44-53 UTEN UTSKJÆRING STORM4	1
8	SP1586049	(Tray table, 760mm, no Cutout)	1
9	SP1548138	(Swing-Out Tray, 540mm, Seatwidth 330mm to 430mm, Cutout RH, cpl. with holder)	1
10	SP1548199	BORD STYREBOKS HØ SIDE STORM4	1
11	SP1422953	BORD NEDFELLBART VE STYREBOKS HØ STORM3	1
12	SP1584663	(Swing-Out Tray, 760mm, Seatwidth 480mm to 580mm, Cutout RH, cpl. with holder)	1
13	SP1535051	BORD 540MM SB33-43 STORM4 UTSKJÆRING HØ	1
14	SP1535012	BORD 620 MM SB38-48 UTSK HØ STORM 4 UTEN FESTE	1
15	SP1535049	BORD 660MM SB 44-53 UTSK HØ SIDE UTEN FESTE	1
16	SP1586048	(Tray table, 760mm, Cutout RH)	1
17	SP1548197	(Swing-Out Tray, 540mm, Seatwidth 330mm to 430mm, Cutout LH, cpl. with holder)	1
18	SP1548200	BORD STYREBOKS VE SIDE STORM 4	1
19	SP1422952	BORD NEDFELLBART HØ STYREBOKS VE STORM3	1
20	SP1584662	(Swing-Out Tray, 760mm, Seatwidth 480mm to 580mm, Cutout LH, cpl. with holder)	1
21	SP1535009	BORD 540MM SB33-43 STORM4 UTSKJÆRING VE	1
22	SP1534974	BORD 620 MM SB38-48 UTSK VE STORM 4 UTEN FESTE	1
23	SP1534975	BORD 660MM SB38-48 UT VE STORM4 UTSKJÆRING VE	1
24	SP1586047	(Tray table, 760mm, Cutout LH)	1
25	SPZ177011	HOLDER FOR BORD SVINGBAR KPL STORM3/STORM4	1
26	Z177025	FESTE KIT BOR STORM4	1
27	SP1590090	(Tray Support bumper for Soft Armpad (increased height))	1

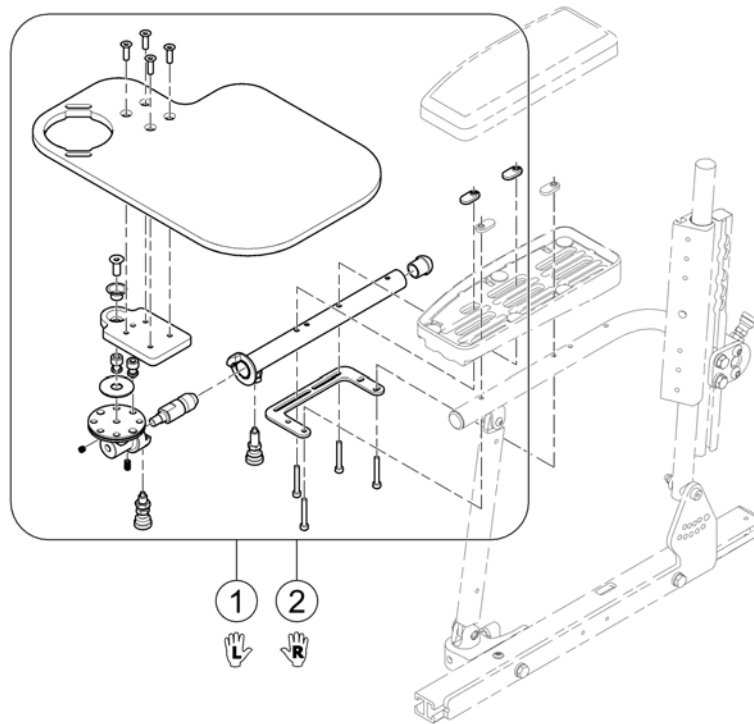

1	SPTRC2065	(Flip and fold away Tray LH)	1
2	SPTRC2066	(Flip and fold away Tray RH)	1

(Hour Counter)

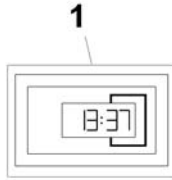1567259.iso




1 SP1541531 (Hour Counter)

1

(T-Nuts for Accessories Modulite & ULM Seat)

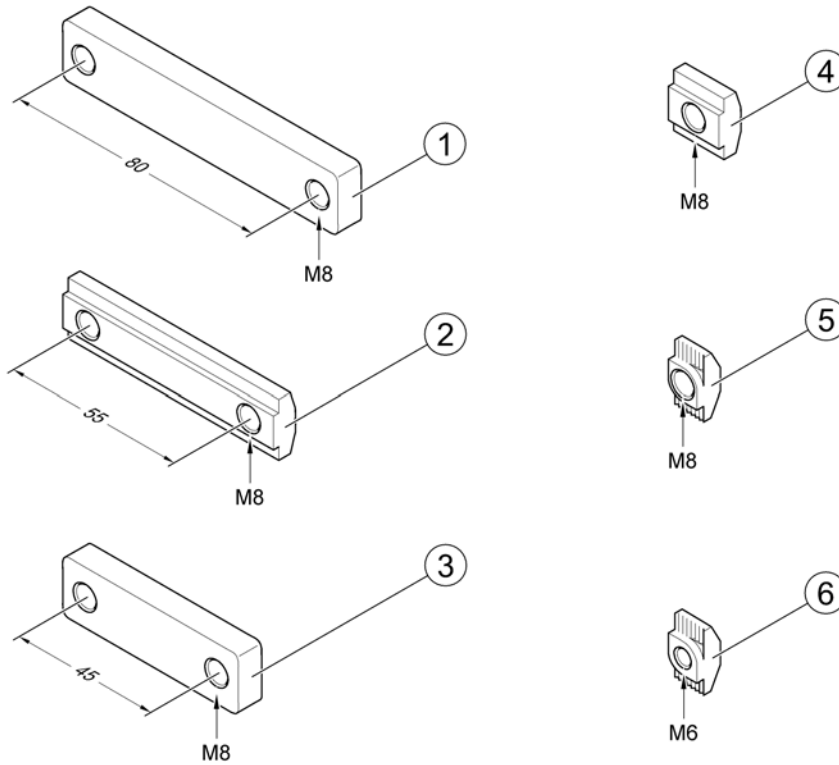157-4067 no



1	SP1556401	MUTTER SKINNE 80MM 2*M8 MODULITE	1
2	SP1572444	MUTTER SKINNE 55MM 2*M8 MODULITE	1
3	SP1555168	MUTTER SKINNE 45MM 2*M8 MODULITE	1
4	SP1557341	MUTTER SKINNE M8 BRED MODULITE	1
5	SP1573158	MUTTER SKINNE VRIDBAR M8 MODULITE	1
6	SP1573721	MUTTER SKINNE VRIDBAR M6 MODULITE	1



Yes, you can.®



Belgium & Luxemburg:

Invacare nv, Autobaan 22, B-8210 Loppem
Tel: (32) (0)50 83 10 10, Fax: (32) (0)50 83 10 11
belgium@invacare.com / www.invacare.be

Denmark:

Invacare A/S, Sdr. Ringvej 37, DK-2605 Brøndby
Tel: (45) (0)36 90 00 00, Fax: (45) (0)36 90 00 01
denmark@invacare.com / www.invacare.dk

Deutschland:

Invacare Aquatec GmbH, Alemannenstraße 10, D-88316 Isny
Tel: (49) (0)75 62 7 00 0, Fax: (49) (0)75 62 7 00 66
info@invacare-aquatec.com / www.invacare-aquatec.de

European Distributor Organisation:

Invacare, Kleiststraße 49, D-32457 Porta Westfalica
Tel: (49) (0)57 31 754 540, Fax: (49) (0)57 31 754 541
edo@invacare.com / www.invacare.eu.com

España:

Invacare SA, c/Areny s/n, Polígon Industrial de Celrà, E-17460 Celrà (Girona)
Tel: (34) (0)972 49 32 00, Fax: (34) (0)972 49 32 20
contactsp@invacare.com / www.invacare.es

France:

Invacare Poirier SAS, Route de St Roch, F-37230 Fondettes
Tel: (33) (0)2 47 62 64 66, Fax: (33) (0)2 47 42 12 24
contactfr@invacare.com / www.invacare.fr

Ireland:

Invacare Ireland Ltd, Unit 5 Seatown Business Campus, Seatown Road, Swords, County Dublin - Ireland
Tel: (353) 1 810 7084, Fax: (353) 1 810 7085
ireland@invacare.com / www.invacare.ie

Italia:

Invacare Mecc San s.r.l., Via dei Pini 62, I-36016 Thiene (VI)
Tel: (39) 0445 38 00 59, Fax: (39) 0445 38 00 34
italia@invacare.com / www.invacare.it

Nederland:

Invacare BV, Celsiusstraat 46, NL-6716 BZ Ede
Tel: (31) (0)318 695 757, Fax: (31) (0)318 695 758
nederland@invacare.com, csede@invacare.com / www.invacare.nl

Norge:

Invacare AS, Grensesvingen 9, Postboks 6230, Etterstad, N-0603 Oslo
Tel: (47) (0)22 57 95 00, Fax: (47) (0)22 57 95 01
norway@invacare.com, island@invacare.com / www.invacare.no

Österreich:

Invacare Austria GmbH, Herzog Odilostrasse 101, A-5310 Mondsee
Tel: (43) 6232 5535 0, Fax: (43) 6232 5535 4
info-austria@invacare.com / www.invacare.at

Portugal:

Invacare Lda, Rua Estrada Velha, 949, P-4465-784 Leça do Balio
Tel: (351) (0)225 1059 46/47, Fax: (351) (0)225 1057 39
portugal@invacare.com / www.invacare.pt

Sverige & Suomi:

Invacare AB, Fagerstagatan 9, S-163 91 Spånga
Tel: (46) (0)8 761 70 90, Fax: (46) (0)8 761 81 08
sweden@invacare.com, finland@invacare.com / www.invacare.se

Switzerland:

Invacare AG, Benkenstrasse 260, CH-4108 Witterswil
Tel: (41) (0)61 487 70 80, Fax: (41) (0)61 487 70 81
switzerland@invacare.com / www.invacare.ch

United Kingdom:

Invacare Limited, Pencoed Technology Park, Pencoed, Bridgend CF35 5AQ
Switchboard Tel: (44) (0)1656 776200, Fax: (44) (0)1656 776201
Customer services Tel: (44) (0)1656 776222, Fax: (44) (0)1656 776220
uk@invacare.com / www.invacare.co.uk