



Exceptional pressure re-distribution for heavier patients

Softform Bariatric provides a high level of patient comfort whilst offering effective pressure re-distribution. Suitable for use on all types of profiling beds, the mattress is essential for patients whose needs cannot be met with standard products. This high specification pressure redistributing mattress has been designed to meet the demands of Community and Acute healthcare environments. Practical and durable, this mattress is suitable for those heavier patients that are considered at 'Very High Risk'¹ of developing pressure ulcers.



The challenges of obese patients

Obese patients (a BMI >40, WHO 2000) have an increased risk of pressure ulcer development. Due to reduced vascularity and profusion of adipose tissue in obese patients wound healing is prolonged so preventative measures for Bariatric patients are essential. This mattress has been specifically designed to meet the needs of heavier patients and the challenges faced by carers and clinicians in managing their care.

Unique Glide System


The patented Glide design and split section at the backrest section of the mattress enables it to conform effectively with a profiling bed. This significantly reduces inappropriate and unwanted patient movement that can contribute to tissue damage caused by shear and friction forces, so reducing the need for unnecessary moving and handling, a difficult task for this group of patients.

Features and Options



Added Protection and Comfort

Multi stretch, breathable, waterproof cover with corner welds, reducing the opportunity for accumulation of dust and liquid residues.

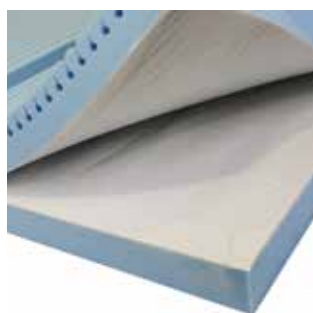

Added protection

High frequency welded waterproof seams protect the mattress core from contamination. The base of the **Softform Bariatric** is made from a toughened PU coated material which is designed to extend the life of the mattress. The vapor permeable fabric allows the skin to breathe and can reduce the effects of skin damage caused by profuse sweating, which can be a common problem for many obese patients.



The cut foam castillations distribute pressure evenly, reducing the risk of pressure ulcers for the bariatric patient who will likely have a limited range of movement. Many obese patients have difficulty lying flat due to breathing difficulties. The hinged backrest section of the **Softform Bariatric** mattress will be an advantage to many.

Technical data

Softform Bariatric	1150 - 1200 mm	Length: 2000 mm	Height: 152 mm	Weight: 15 kg
				Max. User Weight: 340 kg (53.5 Stone)

Softform Bariatric	Cleaning: 80°
	Bariatric use

Warranty 1 year

Invacare Limited
Pencoed Technology Park - Pencoed - Bridgend CF35 5AQ - United Kingdom
Tel: +44 1 656 776 222
Fax: +44 1656 776 220
E-mail: uk@invacare.com
Sales Order E-mail: ordersuk@invacare.com
www.invacare.co.uk

Invacare Ireland Ltd
Unit 5 - Seatown Business Campus Seatown Road - SWORDS - County Dublin - Ireland
Tel. +353 1 8107084
Fax +353 1 8107085
Email: Ireland@invacare.com
www.invacare.ie

© 2016 Invacare International Sàrl
All rights reserved.

All information quoted is believed to be correct at time of print. Invacare reserves the right to alter product specifications without prior consultation.

Softform Bariatric - UK - 02/2016 - 1582396

Colour



Blue

1/ Essential nursing care is pivotal in pressure ulcer prevention. This mattress will positively contribute to the outcome of a pressure ulcer prevention care plan. Education, clinical judgement and action based planning based on vulnerability are fundamental factors in the prevention of pressure ulcers. A range of assessment scales can be used as a formal method of assessing risk from pressure ulcer development, and should be used in conjunction with an informal assessment (informed nursing judgement). Informal assessment is considered to be of greater importance and clinical value.



Yes, you can.®