

Soft Tilt[®] 3

Effortless repositioning.
Maximum comfort.



Automated repositioning system
For medical profiling beds



Enhanced patient care



01



Efficient pressure care

Moves patient automatically from side to side to relieve pressure.

02



Single-handed care

Supports manual patient handling activities and helps reduce physical strain.

03




Sensory stimulation

Both wings elevate to envelope the patient, creating a calming and comfortable effect.

Soft Tilt 3 is designed to enhance patient care with automated repositioning functions that work in combination with medical profiling beds, for maximum patient and caregiver comfort.*

This intuitive bed option aims to support pressure management, sensory stimulation and single-handed care.

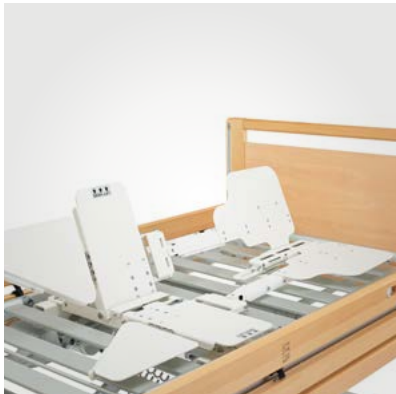

Wings elevate
40° max when automated.
73° max manually.



30° for envelopment.

08



Control box

Can be programmed to automatically turn patients.

07



04



Quality Rest and Sleep

Less caregiver interaction enables undisturbed sleep.

05



Easy to install

Secures to the mattress platform and aligns with bed profiling mechanisms.

06



Hand control

Intuitive hand control for manual use.

06

*Soft Tilt 3 fits Invacare beds, NordBed Optimo and Ultra (90 cm wide) and SB 755 (85 / 90 cm wide).

Technical specifications



Max angle in
manual mode
73°



Max angle in
automated mode
40°



Time values in
automated mode
0 – 180 minutes



Max user
weight
165 kg



Heaviest
part
12 kg



Max cleaning
temperature
80° (cover)

Colour



White

Please note: Soft Tilt 3 fits Invacare beds, NordBed Optimo and Ultra (90 cm wide) and SB 755 (85 / 90 cm wide).

Tel: +41 61 487 70 70
globalcomms@invacare.com
www.invacare.eu.com

Invacare International GmbH
Neuhofweg 51
CH-4147 Aesch BL
SWITZERLAND

© 2025 Invacare International GmbH. All rights reserved. All information quoted is believed to be correct at time of print. Invacare reserves the right to alter product specifications without prior consultation. EU_21102025

