


# Softform<sup>®</sup> Premier

**Pressure-redistributing mattress**

A high-quality mattress for comfort and pressure redistribution with a contoured insert pad.

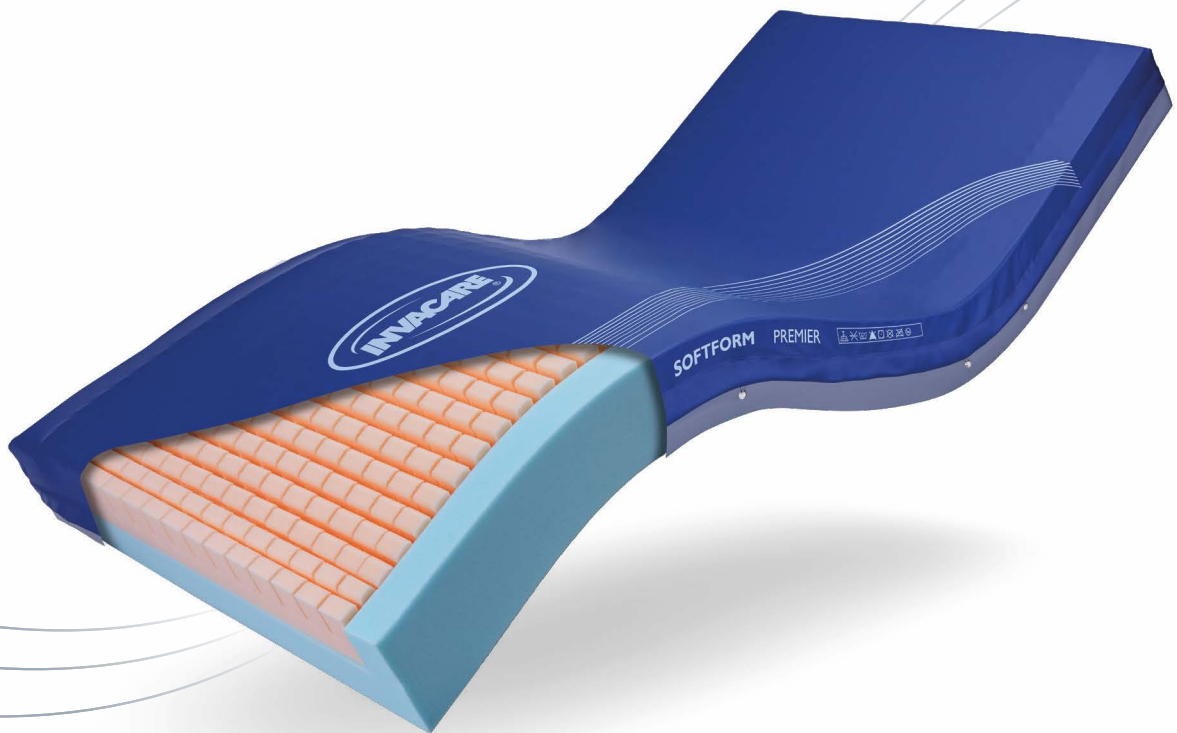

**Stability and pressure redistribution**

For personalised pressure injury management



**Yes, you can.<sup>®</sup>**

# Effective pressure injury prevention



## Pressure Redistribution

- › High quality foam offers a comfortable and effective pressure-redistributing support surface for patients at high to very high risk of developing pressure injuries.
- › Castellated foam helps redistribute weight over a greater area.
- › Non-glued construction allows good fit with profiled and electric beds.



## Shear & Friction Reduction

- › Castellated foam conforms to movement and body position, helping to reduce shear and friction.
- › The top surface of the mattress tapers to relieve pressure on the heels, which are particularly vulnerable to pressure injuries.
- › Two-way stretch, water-resistant, polyurethane cover is durable and helps to reduce shear and friction forces.



## Comfort

- › Castellations provide excellent comfort and air channels contributing to microclimate management.
- › Robust side walls and base facilitate stability during transfers and prevents the mattress from bottoming out.
- › Four-sided concealed welded zip reduces potential fluid ingress and allows easy detachment for inspection and laundering.

## Technical data ▼

For more information about this product, including the user manual, please visit our website: [www.invacare.co.uk](http://www.invacare.co.uk)

Sizes available ▼

|                      |         | 830 mm | 850 mm | 880 mm | 900 mm | 1050 mm | 1200 mm | 1370 mm* |
|----------------------|---------|--------|--------|--------|--------|---------|---------|----------|
| <b>140 mm height</b> | 2000 mm | -      | •      | •      | -      | -       | -       | -        |
|                      | 1900 mm | •      | •      | •      | •      | -       | •       | •        |
| <b>150 mm height</b> | 1970 mm | •      | •      | •      | •      | •       | •       | •        |
|                      | 2000 mm | •      | •      | •      | •      | -       | •       | •        |
| <b>200 mm height</b> | 2000 mm | -      | -      | •      | -      | -       | -       | -        |



Max. User Weight ▼

247.6 kg



Max. Cleaning Temperature ▼

95 °C

Other sizes may be available upon request. Please ask your local Invacare representative.

Manufactured to BS 7177 (Ignition Source 5); Foam Interior: BS 5852:1982 (Ignition Source 5); Cover BS 7175 (Ignition Source 5).

\*NB. The 1370 mm wide version mattress are not available with welded seams, only sewn seams

Essential nursing care is pivotal in pressure injury management. This mattress may positively contribute to the outcome of a pressure injury prevention plan that includes education, clinical judgment and action planning based on vulnerability, which are fundamental factors in the prevention of pressure injuries. A range of assessment scales can be used as a formal method of assessing risk from pressure injury development, and should be used in conjunction with an informal assessment (informed nursing judgement) as part of a comprehensive and individualised plan.

T: +44 1656 776222  
[uk@invacare.com](mailto:uk@invacare.com)  
[www.invacare.co.uk](http://www.invacare.co.uk)  
[www.invacare.ie](http://www.invacare.ie)

Invacare UK & Ireland  
 Unit 4 Pencoed Technology  
 Park  
 Bridgend CF35 5AQ



Yes, you can.®



© 2024 Invacare International GmbH. All rights reserved. All information quoted is believed to be correct at time of print. Invacare reserves the right to alter product specifications without prior consultation. Softform Premier - UK - 11/2024