


Softform[®]

Premier Original

Pressure redistributing mattress

High-quality and comfortable
pressure-redistributing support surface.

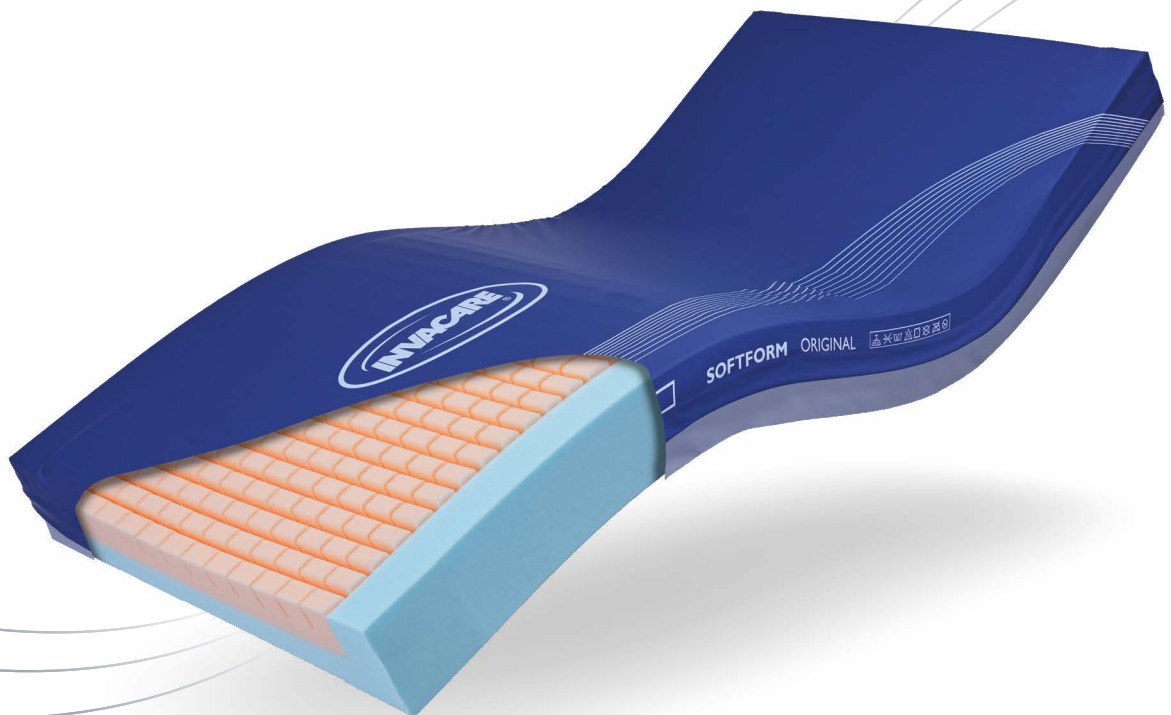

Stable and effective pressure redistribution

For personalised pressure injury management



Yes, you can.[®]

Effective pressure injury prevention



Pressure Redistribution

- › High-quality foam offers a comfortable and effective pressure-redistributing support surface for people at high to very high risk of developing pressure injuries.
- › Castellated foam helps redistribute weight over a greater area.
- › Non-glued construction allows good fit with profiled and electric beds.



Shear & Friction Reduction

- › Castellated foam conforms to movement and body position, helping to reduce shear and friction.
- › Two-way stretch, water-resistant, polyurethane cover is durable and helps to further reduce shear and friction forces.



Comfort

- › Castellations provide excellent comfort and air channels contributing to microclimate management.
- › Robust side walls and base facilitate stability during transfers and prevents the mattress from bottoming out.
- › Four-sided concealed welded zip reduces potential fluid ingress and allows easy detachment for inspection and laundering.

Technical data ▼

For more information about this product, including the user manual, please visit our website: www.invacare.co.uk



Sizes Available ▼



Max. User Weight ▼



Max. Cleaning Temperature ▼

150 mm height

2000 x 855 mm
2000 x 880 mm
2050 x 880 mm

247.6 kg

95 °C

Manufactured to: BS 7177 (Ignition Source 5).

Essential nursing care is pivotal in pressure injury management. This mattress may positively contribute to the outcome of a pressure injury prevention plan that includes education, clinical judgment and action planning based on vulnerability, which are fundamental factors in the prevention of pressure injuries. A range of assessment scales can be used as a formal method of assessing risk from pressure injury development, and should be used in conjunction with an informal assessment (informed nursing judgement) as part of a comprehensive and individualised plan.

T: +44 1656 776222
uk@invacare.com
www.invacare.co.uk
www.invacare.ie

Invacare UK & Ireland
Unit 4 Pencoed Technology
Park
Bridgend CF35 5AQ



Yes, you can.®



© 2024 Invacare International GmbH. All rights reserved. All information quoted is believed to be correct at time of print. Invacare reserves the right to alter product specifications without prior consultation. Softform Premier Original - UK - 11/2024