

P O W E R S E R I E S

TDX[®] SP2 Modulite



Highly responsive mid wheel drive power chair
with versatile Modulite seating



| TDX® SP2 Modulite

Optimal comfort and performance.

The **TDX®** SP2 Modulite is one of our highest performing mid wheel drive power chairs.


Quality and technology at its core, combined with our highly adaptable Modulite seating and advanced LiNX® controls to give optimal comfort and performance every time.



| TDX® SP2 Modulite

It's more than a base.

The **TDX® SP2** is built to perform and effortlessly take you where you want to go. Easy to control and manoeuvre to help give you more freedom and independence in any environment.

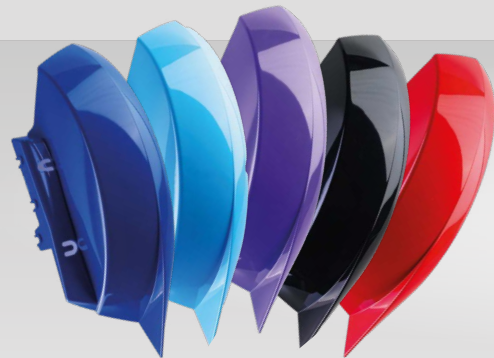

New shroud collection

Choose from a selection of exciting new colour shrouds for a stylish look.



Width **657 mm**



640 mm turning radius



Excellent manoeuvrability

The mid wheel drive makes navigation through tight spaces easier and gives you more controlled movement.



Patented SureStep® suspension

Provides a smooth and comfortable ride even over challenging surfaces.



Enhanced stability & safety

Designed to ensure all six wheels remain firmly on the ground with the integrated stability lock function.

| TDX® SP2 Modulite

New colours. Extra durability.

Add a fresh and stylish touch with a new shroud colour. The shroud's extra durability helps protect and withstand everyday challenges.



Obsidian Black

Ruby Red

Purple Passion

Peacock Blue

Arctic Blue

| TDX® SP2 Modulite

Extensive positioning options.

Our versatile Modulite seating system gives a personalised position for extra comfort and support.

Ranging from more basic seating requirements through to more complex postural needs, the modular design maximises positioning possibilities.



Complete solutions

Extensive range of seating options including backrests, seat cushions, armrests, headrests and many others.



Designed for safety & comfort

50° tilt, 300mm seat rise with 45° tilt bring comfort and security.



LNX legrest

Extends in both length and elevation, working with powered recline to help reduce shear.



Award-winning Matrix®

Personalised Matrix seating options to ensure an optimum fit every time.

Range of tilt options

Provide pressure management and maximised comfort.



300 mm of seat rise combined with up to **45°** tilt positioning

| TDX® SP2 Modulite

Modulite seating. Modular versatility.

Our most popular seating solution.
A modular approach to seating that maximises possibilities. Adjusts with ease for the perfect fit and optimised high levels of comfort.

Tech spec

- Seat width: 380 - 580 mm
- Seat depth: 410 - 560 mm
- Seat height: from 450 mm
- Armrest height: 245 - 390 mm
- Footrest length: 290 - 460 mm



Adjustable
seating
for a
perfect fit



| TDX® SP2 Modulite

Featuring advanced **LiNX**® technology.

LiNX® technology helps give a personalised and confident drive experience. The innovative technology learns habits and preferences over time making it a super intuitive system.



| TDX® SP2 Modulite

LiNX®

Discover more.



New Resolve joystick mount

- Huge range of achievable joystick positions for optimised comfort.
- Adjustable tension release.
- Optimised comfort and control.



3.5" Touchscreen

Impressive 3.5" colour touchscreen with joystick, speedometer and an odometer.



Specialist controls

Including compact remotes, head arrays and four-switch proximity arrays.



Built-in Bluetooth®

Enables connectivity to devices and environmental controls.



Control input options

A variety of control input options allow a personalised set up for maximum control and reassurance that the power chair will react consistently.




| TDX® SP2 Modulite

Technical data.

Further information about this product, including the user manual, please visit our website www.invacare.eu.com.

New colours

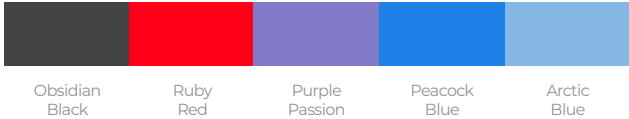

	Seat width	Seat depth	Seat height	Backrest height	Armrest height	Footrests	Seat angle
1-piece seat plate ▶	430 mm	410 - 510 mm	450 mm, 475 mm	480 mm	275 - 390 mm	290 - 460 mm	0° - 20°
Telescopic seat ▶	380 - 580 mm	410 - 560 mm	450 mm, 475 mm	480 - 610 mm	245 - 360 mm	290 - 460 mm	-5° - 50°
Max seat ▶	550 - 710 mm	560 - 630 mm	470 mm, 495 mm	570 - 680 mm	250 - 380 mm	360 - 480 mm	0° - 20°
	Backrest angle	Total height	Total length inc. footrests	Total length without footrests	Driving unit width	Max. user weight	Battery capacity
1-piece seat plate ▶	90° - 120°	1010 mm	1208 mm	980 mm	657 mm	180 kg	60 Ah or 73.5 Ah
Telescopic seat ▶	90° - 120°	1170 mm	1208 mm	980 mm	657 mm	180 kg	60 Ah or 73.5 Ah
Max seat ▶	90° - 115°	1125 mm	1240 mm	970 mm	657 mm	180 kg	60 Ah or 73.5 Ah
	Turning diameter	Max. climbable obstacle height	Max. safe slope	Speed	Driving range	Total weight	
1-piece seat plate ▶	1280 mm	70 mm	9° (15.8%)	3, 6, 10 or 12 km/h	up to 43.6 km	154 kg	
Telescopic seat ▶	1280 mm	70 mm	9° (15.8%)	3, 6, 10 or 12 km/h	up to 43.6 km	164 kg	
Max seat ▶	1280 mm	70 mm	9° (15.8%)	3, 6, 10 or 12 km/h	up to 43.6 km	170 kg	

| TDX® SP2 Modulite

Share your thoughts and help shape the future of power chairs.

We invite you to **share your feedback**.
Your insights will help us better serve
the power chair community.

Please scan and complete the form.



Thank you for letting us join your journey.

Invacare Power Series Team





Share your feedback.

SUBSCRIBE



My chair allows me to live
an independent life and
get outdoors everyday.

Ian Mackay
Invacare Ambassador

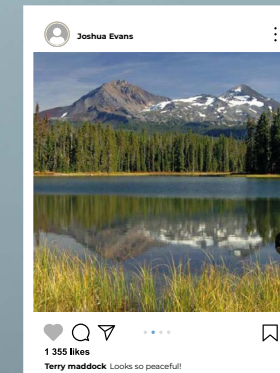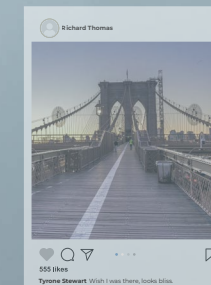

P O W E R S E R I E S

TDX SP2

AVIVA RX

AVIVA FX

Esprit Action

Kompas

Storm 4

Tel: +41 61 487 70 70
hqeuropa@invacare.com
www.invacare.eu.com

© 2026 Invacare International GmbH. All rights reserved.
All information quoted is believed to be correct at time of print.
Invacare reserves the right to alter product specifications without
prior consultation. EU - TDX SP2 Modulite - 22052026

