



ULTRALowmaxX
By Motion Concepts

Invacare® | **TDX2 Ultra**

Ultimate control. Ultimate comfort.

LINX®



The *Invacare TDX2 Ultra* is an elite, high performance, centre wheel drive power wheelchair equipped with the advanced Ultra Low Maxx seating system and revolutionary **LiNX** smart technology. Designed to be the new benchmark in the complex rehab arena, this powerchair is the perfect mix of functionality, customisation and style that adapts to the clients needs and desired lifestyle. Key features include:

- Enhanced stability
- Intuitive driving experience
- **LiNX** Smart Technology
- Range of alternative controls
- Ultimate posture support and shear management



Enhanced stability

The **TDX2 Ultra** is tried, tested and tough. Offering maximum stability and safety, the Stability Lock function ensures all six wheels stay firmly on the ground on uneven terrain. The patented **SureStep**® suspension technology enables a smooth ride while negotiating challenging surfaces. With larger batteries for a greater driving range, the **TDX2 Ultra** enables you to go further.

Intuitive driving experience

Invacare TDX2 Ultra customers can say “Yes” to mobility and “No” to compromise. With an incredibly tight turning radius, thanks to its Centre Wheel Drive, this remarkable wheelbase provides excellent manoeuvrability and an intuitive driving experience.

LiNX Smart Technology: Redefining Mobility

Invacare LiNX is our insight inspired control system with advanced technology that provides a superb driving experience and allows professionals to configure and tailor powerchairs quickly and more intuitively.

The smart control system enables programming and maintenance updates to be made wirelessly and in real time, plus it has straightforward access to statistical and diagnostic information, giving complete visibility to professionals.



Touchscreen technology

The REM400 remote has an impressive 3.5” colour touchscreen that requires minimal force to operate. Bringing modern technology to the powerchair market, this innovative interactive remote can be configured to a clients’ physical and cognitive ability.



Display only version

The REM500 is a display only version of the REM400, retaining the 3.5” interactive touch screen, and allows alternative controls to be used to operate the powerchair.



Built in Bluetooth

The REM400 remote allows individuals to be connected to 3 devices simultaneously and has mouse mover functionality.

Range of alternative controls

To maximise control and independence, the REM400 and REM500 allow a selection of alternative controls and joysticks to be used on the **TDX2 Ultra** powerchair:



Compact remote with LED display

Ideal for those who have a limited range of motion and reduced strength.



Extremity control joystick

Ideal for chin control with the outer shell being resistant to saliva.



MEC joystick

Ideal for those with muscle weakness or fine motor control.



Compact single switch joystick

Ideal for individuals with high muscle tone/uncontrolled movement.



Paediatric compact joystick

Ideal for midline mount.



Atom head array

Three proximity sensors mounted in the headrest with in-built Bluetooth.



Proton head array

Adjustable headrest with three proximity sensors and in-built Bluetooth.



Sip 'n' puff head array with lip switch

Combines simple sip and puff movements with head movements.



Four switch proximity array

Ideal for those with a defined range of motion and no force.



Attendant control unit

Gives carers access to certain chair functions.

IDC dual control

Allows access to the drive function, controlled by a speed handle with emergency stop feature.

Intelligent seating system

The Ultra Low Maxx seating system can combine up to 50° tilt and 170° of recline. Not only does this ensure effective redistribution of pressure for ischial relief and tissue reperfusion, it also warrants sufficient tilt for effective positioning. Not forgetting the importance of comfort, the ability to tilt to 50° allows the individuals to find a position of comfort optimal for them.



Ultimate recline technology

The **Extended Shear Reduction (ESR)** mechanism built into the Ultra Low Maxx system reduces the risk of shear which can cause tissue breakdown, prone to happen when the chair goes in and out of recline.

Never lose contact

Ultra Low Maxx has been specially designed to replicate the natural pivot points of the body, which ensures the client's head, shoulders, hips and feet don't slide or move when the chair is transitioning between positions.

Modular design for optimal support

The customisable seating system adapts to the client's body over time due to its varying width and depth options, its selection of backrests, headrests and armrests plus its built in rail, which allows for a range of positioning accessories to be easily mounted to support individual needs.





LiNX G-Trac® technology

G-Trac ensures the directions indicated by the control are followed by the chair.



SureStep® and Stability Lock

Allows easy obstacle climbing and guarantees a smooth downward transition.



Seat lifter

The 300mm seat lifter provides easy access to higher obstacles.



Hip and lateral supports

Available in multiple pad sizes and are easy to insert, adjust, move and operate.




Single sided fork

Minimises the front width of the base with no protruding screws or fork parts.



True centre wheel drive

Allows excellent manoeuvrability in narrow spaces and ensures an intuitive driving experience.

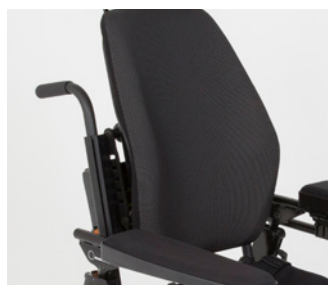

Multiple backrests

Offers multiple, distinct backrests providing varying levels of lateral trunk support.



USB charger

2 x USB 2.0 ports supply up to 1A each to keep devices charged. Rail mounting feature enables easy access.



DuraWatt motors

Quiet motors with sealed lifetime brushes reduce potential brush tolerance issues.



LED lighting system

Attractive, energy saving LED lights can be integrated in a well-protected position in the rear shroud.

Black tyres and castors

Black non-marking tyres and castors with the black rim gives the chair a stylish, modern look.



Anterior assist

5° or 10° options are available for individuals needing assistance with transfers or to enhance day to day reach and function.



Choice of armrests

Available with both cantilever style armrests, to ensure an open rail for positioning components, and with reclining armrests for a high level of support.



LNX powered centre mount legrests


Offers an active seating position and enables legs to be positioned at a 90° angle.

Headrest options

Different headrest styles available depending on the level of support needed.

Technical data

For more comprehensive information about this product, including the product's user manual, please visit www.invacare.co.uk

 <p>Seat width</p>	 <p>Seat depth</p>	 <p>Seat height</p>	 <p>Backrest height</p>	 <p>Armrest height</p>	 <p>Length of footrests</p>	 <p>Seat angle</p>
Invacare TDX2 Ultra 405 - 610 mm	405 - 580 mm	440, 460 and 480 mm ¹	355 - 505 mm ²	230 - 405 mm ³	270 - 474 mm ⁴	0° to 50° ⁵
 <p>Backrest angle</p>	 <p>Total height</p>	 <p>Total length inc. footrests</p>	 <p>Total length without footrests</p>	 <p>Driving unit width</p>	 <p>Max. user weight</p>	 <p>Battery capacity</p>
Invacare TDX2 Ultra 60° - 170°	1080 - 1535 mm	1120 mm ⁶	920 mm	650 mm	136kg	73.5 Ah
 <p>Turning diameter</p>	 <p>Max. climbable obstacle height</p>	 <p>Max. safe slope</p>	 <p>Motor capacity</p>	 <p>Speed</p>	 <p>Driving range</p>	 <p>Total weight</p>
Invacare TDX2 Ultra 1120 mm	75 mm	9° / 16 %	2 x 340 W	6, 8 and 10 km/h	up to 31 km ⁷	156kg ⁸ 169kg ⁹ 179kg



matrx
SEATING SERIES

Invacare **Module**



1. With a standard tilt
2. Depending on back support choice
3. Depending on armrest choice
4. Depending on legrest choice
5. 0°-45° for versions with lifter
-5° and -10° anterior mounting options
+5° pre tilt mounting option
6. Depending on legrest choice and seat position

7. Measured according to ISO 7176-4:2008 with 6 km/h and 73.5 Ah batteries; driving range will vary by configuration (battery size, speed).
8. TDX SP2 NB with tilt module
9. TDX SP2 with lifter/tilt module
10. Invacare Modulite™ seat components are cross compatible across the Invacare power wheelchair range and partly with the Invacare Rea® family.

Shroud and rim insert colours

Island Blue	Lights Out Matte - Black	RockStar Red	SpaceStation Silver
Invacare Blue	Wet Black	Grape JellyBean	

The LINX word, mark and logos are registered trademarks owned by Dynamic Controls and any use of such marks by Invacare is under license.

Don't forget, you can customise your **TDX SP2** by using the Invacare UNIQUE service.

TDX2 Ultra - UK - 07/2022
-1587754



Like us on Facebook:
Invacare UK



Follow us on Twitter:
@Think_Mobility



Follow us on YouTube:
Invacare Ltd - UK

Invacare Limited
Pencoed Technology Park - Pencoed
Bridgend CF35 5AQ - United Kingdom
Tel: +44 1656 776 222
Fax: +44 1656 776 220
E-mail: uk@invacare.com
Sales Order E-mail: ordersuk@invacare.com
www.invacare.co.uk

Invacare Ireland Ltd
Unit 5 - Seatown Business Campus
Seatown Road - SWORDS - County Dublin - Ireland
Tel. +353 1 8107084
Fax +353 1 8107085
Email: Ireland@invacare.com
www.invacare.ie

